

	Form, Engineering Build		Doc No.: FM-0066
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 1 of 5

I. Planning Information

Product Information				
Build No. / Lot No.	Project No.		Customer	
SLQ-80040411241E	00047		SILQ	
Part Number	Rev	DMR Status	Description	
SLQ-8004	A	☒ Draft ☐ Released	SILQ 2-way ClearTract 16 Fr x 10 ml Catheter, Non-Sterile, 10-Pack	
Comments				
☐ N/A The purpose of this build is to coat 16 Fr x 10 mL catheters on the 1-up machine with polymer solution provided by SILQ to generate samples so that SILQ can perform assessments on the polymer shelf life. 50 units are needed. Pathway Medtech does not need to fully build up the catheters to the drawing as the catheters will not be pouched or put in boxes. DMR is in draft state since the catheter is being coated with the 1-up coater.				

Build Requirements						
#	Purpose / Test Description	Test Protocol #	Notes	Sterile	Labeling	Qty Required
1	N/A	N/A	☒	☐	☐	50
Total (Sum)						50

Special Requirements / Instructions	
Notes / Comments	
☐ N/A Rather than use the method shown in the 1-up MP for rinsing the catheter, the following method will be used. Three glass tubes of water from SILQ's Culligan system will be filled and the catheter will be dunked once in each tube and then placed on the rack to dry. After visual inspection of the catheter once it has finished running through the 1-up machine, the catheters will be sent to SILQ to perform the CV and rinse test. Catheters will be sleeved but not pouched, labeled, or boxed. After sleeving the catheters, seal 25 in one poly bag and 25 in another.	
Sterilization	
☐ N/A Product will not be sterilized.	
Labeling	
☐ N/A Product will not be labeled.	

Recorded By		
Name	Signature	Date
Will Hail		04/11/2024

	Form, Engineering Build		Doc No.: FM-0066
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 2 of 5

II. Manufacturing Information

Material Traceability						
#	Component No.	Rev	Description	Lot No.	Exp. Date	Qty Needed
1	SLQ-6004	A	Catheter, 16 Fr x 10 mL, Silicone, NY Foley, 2-Way	20220902	N/A	50
2	HDXPLX-3000	N/A	Solution, SILQ Polymer	HYX-825021	N/A	A/R
3	SLQ-5000	A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	N/A	N/A	50
4	MP-SLQ-0001	B	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	N/A	N/A	N/A
5	MP-HDX-0000	A	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	N/A	N/A	N/A

Recorded By		
Name	Signature	Date
Will Hail		04/11/2024

Procedures & Personnel Traceability				
Document No.	Rev	Performed By / Signatures		
MP-SLQ-0001	B			N/A WH 04/16/24
MP-HDX-0000	A			N/A WH 04/16/24

☐ N/A	Equipment / Critical Parameters				
Equipment No./ Description	Cal Due Date	Setting	Value	Date / Time	Initials
SILQ's 1-up coater (provided by SILQ for this build)	N/A	Cure Time	120 seconds	04/16/2024 12:15 PM	W.H.

Manufacturing Activity Information		
Mfg Dates	Mfg Completion Date	Final Qty Mfg
04/16/2024 - 04/18/2024	04/18/2024	50

	Form, Engineering Build		Doc No.: FM-0066
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 3 of 5

Deviations	
☒ N/A	

Recorded By		
Name	Signature	Date
Will Hail		04/18/2024

III. Quality Information

Acceptance Activities							
#	Description	Qty Insp	Qty Acc	Qty Rej	Specification	Rev	Result (P/F)
1	Visually inspect the catheter for particulates and / or cured polymer solution	5	5	0	SLQ-8004	B	P
Comments							
☒ N/A							
Date Inspected		Disposition		NCR No.			
04/18/2024		ACCEPT		N/A			

Inspected By		
Name	Signature	Date
Will Hail		04/18/2024
Verified By		
Name	Signature	Date
Sean Tamayo		04/18/2024

Documentation Checklist / Review					
Place documents in order specified below:					
#	Document	Present	Completed	Meets Spec	N/A
1	Work Order (WO) Record		☒		
2	Line Clearance	☒	☒		☐
3	Product's Specification/Drawing & BOM	☒	☒	☒	
4	Labeling Generation/Reconciliation Record	☐	☐	☐	☒
5	Nonconforming Reports (NCR)	☐	☐		☒
6	Inspection Records	☐	☐	☐	☐

 pathway	Form, Engineering Build	Doc No.: FM-0066
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 4 of 5

7	Critical Process Parameter Record	☐	☐	☐	☒
---	-----------------------------------	--------------------------	--------------------------	--------------------------	-------------------------------------

Verified By		
Name	Signature	Date
WILL HALL		04/18/2024

	Form, Engineering Build		Doc No.: FM-0066
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 5 of 5

Release Authorization for Further Processing		
(e.g., terminally sterilization) – The individual signing below is certifying that they independently reviewed the Device History Record (DHR) and confirmed that all required documentation is included, records are completed and good documentation practices were followed, and specifications have been met (if not NCR is included)		
Authorized By		
Name	Signature	Date
WILL HAIL		04/18/2024

☒ N/A	Post Sterilization Acceptance Activities						
#	Description	Qty Insp	Qty Acc	Qty Rej	Specification	Rev	Result (P/F)
1	N/A						
2							
Comments							
☒ N/A							
Date Inspected		Disposition		NCR No.			
N/A		N/A		N/A			

Inspected By		
Name	Signature	Date
N/A	N/A	N/A
Verified By		
Name	Signature	Date
N/A	N/A	N/A

Final Documentation Checklist / Review					
Place documents in order specified below:					
#	Document	Present	Completed	Meets Spec	N/A
8	Sterilization Processing Record	☐	☐	☐	☒

Verified By		
Name	Signature	Date
N/A	N/A	N/A

	Line Clearance Form	Doc No.: FM-0057
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 1 of 1

Line Clearance - The individual performing the line clearance certifies that line clearance activities were conducted per WI-0014 prior to starting manufacturing activities. Quality management or designee verifies the line clearance was performed.

Actions		Performed by (Initials)	Date
1. Clear Area of previous WO Raw Materials.	☒ Completed	WH	04/16/2024
2. Verify kitted Raw Materials (RM) for: a. RM meet BOM revision level requirements b. RM have not Expired c. RM have no physical damage	a. ☒ Completed b. ☒ Completed c. ☒ Completed	WH	04/16/2024
3. Equipment is calibrated and properly functioning. <i>Note: Equipment identification and calibration due date will be documented in the work order.</i>	☐ Completed ☒ N/A	WH	04/16/2024
4. Proper documentation is present and available.	☒ Completed	WH	04/16/2024
5. Production Associates are trained and qualified to perform the process.	☒ Completed ☐ N/A	WH	04/16/2024
6. All work surfaces cleaned with IPA.	☒ Completed ☐ N/A	WH	04/16/2024
7. Verify Labels generated meet DMR requirements	☐ Completed ☒ N/A	WH	04/16/2024
8. Verify validated processes have been set to critical points within specification.	☐ Completed ☒ N/A	WH	04/16/2024
9. Verify Start-up activities within specification	☐ Completed ☒ N/A	WH	04/16/2024
Conclusion: Line clearance performed	☒ Completed	WH	04/16/2024

	Name	Signature	Date
Quality Management	B. Tanager	B.T.	04/16/2024



Date Inspected		Part No.		Rev.	Description				
4/18/2024		SLQ-8004		B	SILQ 2-Way Cleartract Foley Catheter, 16 Fr x 10 mL, Non-Sterile, 10-Pack				
Lot No.		Lot QTY		Approved Supplier(s)					
SLQ-80040411241E		50		Pathway					
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	B	B	B	B	B	B	B	B
	Lower Spec Limit (LSL)	N/A	N/A	0.300	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	N/A	0.020	N/A	N/A	N/A	N/A
	UOM	EA	EA	AU	AU	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Visual	Spectrophotometer	Spectrophotometer	Visual	Visual	Visual	Visual
	Asset No.	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
	Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size	5	5	N/A	N/A	N/A	N/A	N/A	N/A
Sample No.		Variable/ Attribute Data							
1	P	P	N/A	N/A	N/A	N/A	N/A	N/A	N/A
2	P	P	N/A	N/A	N/A	N/A	N/A	N/A	N/A
3	P	P	N/A	N/A	N/A	N/A	N/A	N/A	N/A
4	P	P	N/A	N/A	N/A	N/A	N/A	N/A	N/A
5	P	P	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Result Summary									
Pass QTY	5	5	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Fail QTY	0	0	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Result (P/F)	P	P	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Statistical Summary									



Min	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Comments									
Disposition:	ACCEPT						NCMR No:	N/A	
Inspected By:							Date:	4/18/2024	
Verified By:							Date:	4/18/2024	

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY PATHWAY PRIOR TO PRODUCTION.

2. DRAWING IS FOR REFERENCE ONLY. REFER TO THE RELEASED BILL OF MATERIALS FOR THE COMPLETE PARTS LIST.

3. UP UNTIL SEALING THE POUCH, PRODUCT SHALL BE PACKAGED IN A CLEANROOM ENVIRONMENT. POST SEALING THE POUCH, ASSEMBLY WILL CONTINUE OUTSIDE OF THE CLEANROOM ENVIRONMENT.
2. ASSEMBLY PROCEDURE
1. APPLY POLYMER SOLUTION TO THE CATHETER IN THE HATCHED AREA SHOWN. ALLOW THE CATHETER TO DRY OVERNIGHT.

2. INSERT THE CATHETER INTO A POLYETHYLENE SLEEVE SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CLOSED END.

3. PLACE THE POUCH LABEL ON THE CLEAR SIDE OF THE POUCH. WHEN VIEWING THE CLEAR SIDE OF THE POUCH FROM ABOVE, WITH THE CHEVRON TO THE RIGHT, THE LABEL TEXT SHALL READ LEFT TO RIGHT. WITH BEST EFFORT, APPLY THE LABEL SUCH THAT IT IS CENTERED.

4. INSERT THE SLEEVE / CATHETER INTO THE POUCH SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CHEVRON END OF THE POUCH.

5. SEAL THE POUCH.

6. PLACE 10 SLEEVED, POUCHED, SEALED, AND LABELED CATHETERS INTO ONE SHELF CARTON. APPLY TWO SHIPPER BOX LABELS TO THE LONG, WIDER SIDE OF THE CARTON. VERIFY THAT THE TEXT READS LEFT TO RIGHT WHEN THE OPENING OF THE CARTON IS TO THE LEFT. WITH BEST EFFORT, CENTER THE LABEL ON ITS RESPECTIVE FACE OF THE CARTON.

7. PLACE 10 SHELF CARTONS INTO ONE SHIPPER BOX. SEAL WITH PACKAGING TAPE AS NEEDED. PLACE A SHIPPER BOX LABEL IN THE UPPER LEFT CORNER ON TWO ADJACENT FACES, READING LEFT TO RIGHT.
3. DIMENSIONAL
1. () INDICATES LABEL REFERENCE DIMENSIONS.

2.  INDICATES CRITICAL DIMENSIONS.
4. CLEANING AND HANDLING
1. ALL PRODUCT, UP UNTIL INSERTING THE POUCHES INTO SHELF CARTONS, SHALL BE HANDLED BY GLOVED HANDS APPROPRIATE FOR A CLEANROOM ENVIRONMENT.
5. FUNCTIONAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
1. THE CATHETER'S ABSORBANCE VALUE AT 585 NANOMETERS SHALL BE GREATER THAN OR EQUAL TO 0.300 AU. REFER TO TM-HDX-0001.

2. THE CATHETER'S ABSORBANCE VALUE AT 440 NANOMETERS SHALL BE LESS THAN OR EQUAL TO 0.020 AU. REFER TO TM-HDX-0002.

3. THE PEEL STRENGTH FOR A 1 INCH WIDE SEALED SECTION OF THE STERILE BARRIER PACKAGING SHALL BE >= 1.00 LB WHEN TESTED PER TM-0003. IT IS NOT REQUIRED TO PLACE CATHETERS IN THESE TEST POUCHES, OR LABEL THEM. INSPECT N=3 SAMPLES, 1 AT THE BEGINNING, MIDDLE, AND END OF THE MANUFACTURING RUN.
6. VISUAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
1. THE CATHETER SHALL BE FREE OF CURED POLYMER AND DEBRIS.

2. THE CATHETER SHALL BE PROPERLY ORIENTED IN THE SLEEVE AND POUCH WITH THE LABEL APPLIED TO THE APPROPRIATE FACE AND IN THE CORRECT ORIENTATION.

3. INSIDE THE STERILE BARRIER, DEVICE IS TO BE FREE OF PARTICULATES WHEN VIEWED AT 18 INCHES AWAY UNDER AMBIENT ROOM LIGHTING AND WITHOUT THE USE OF MAGNIFICATION.


4. THE STERILE BARRIER PACKAGING SHALL NOT INCLUDE ANY OF THE FOLLOWING DEFECTS:

1. UNSEALED AREAS, NONHOMOGENEOUS OR UNSEALED AREAS, OVERSEALED AREAS, ANRRROW SEALS, CHANNELS, WRINKLES / FOLDOVERS / CRACKS, OR TEARS / PINHOLES.

5. THE SHELF CARTON SHALL HAVE THE CORRECT NUMBER OF CATHETERS INSIDE AND FACING THE CORRECT ORIENTATION. THE LABELS SHALL BE AFFIXED IN THE CORRECT LOCATION AND ORIENTATION.

6. THE SHIPPER BOX SHALL HAVE THE CORRECT NUMBER OF SHELF CARTONS. THE LABELS SHALL BE APPLIED IN THE CORRECT LOCATION AND ORIENTATION.

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	
	B	UPDATE NOTES BLOCK, SPLIT POUCH AND SHELF CARTON LABEL INTO 2 ITEMS, REPLACED SHELF CARTON	

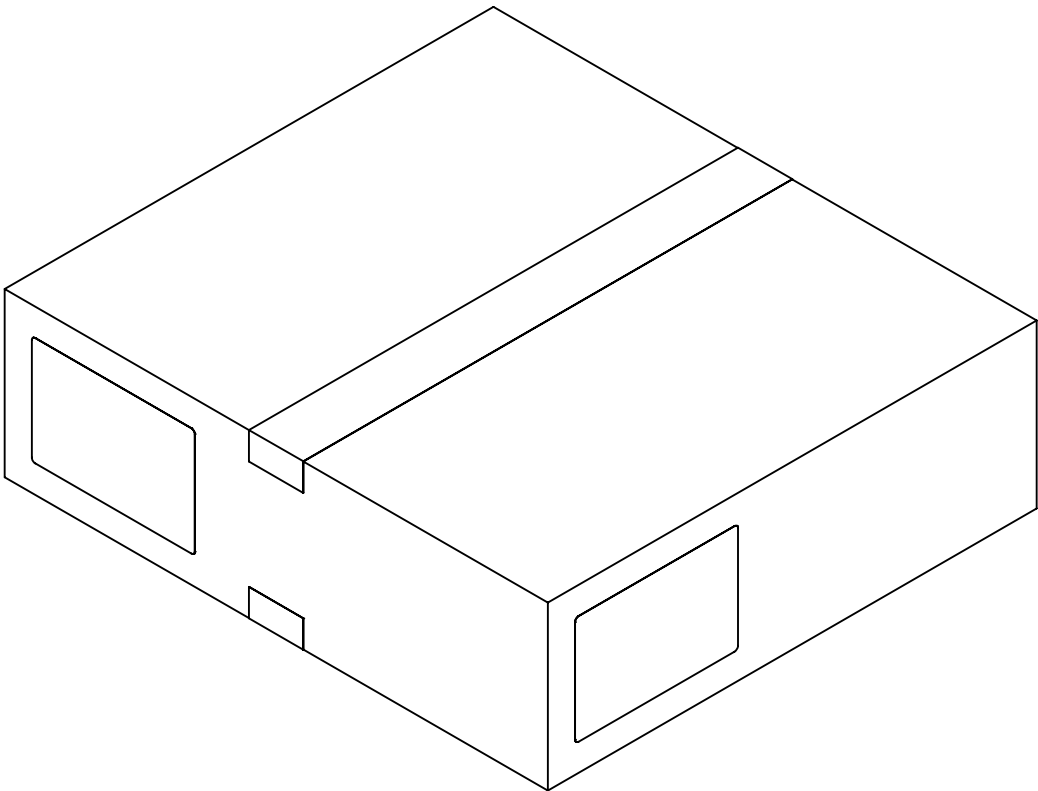

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.

DO NOT SCALE DRAWING




TITLE:
CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY, NON-STERILE, 10-PACK

SIZE	DWG. NO.	REV
B	SLQ-8004	B



SCALE: 1:5 SHEET 1 OF 3



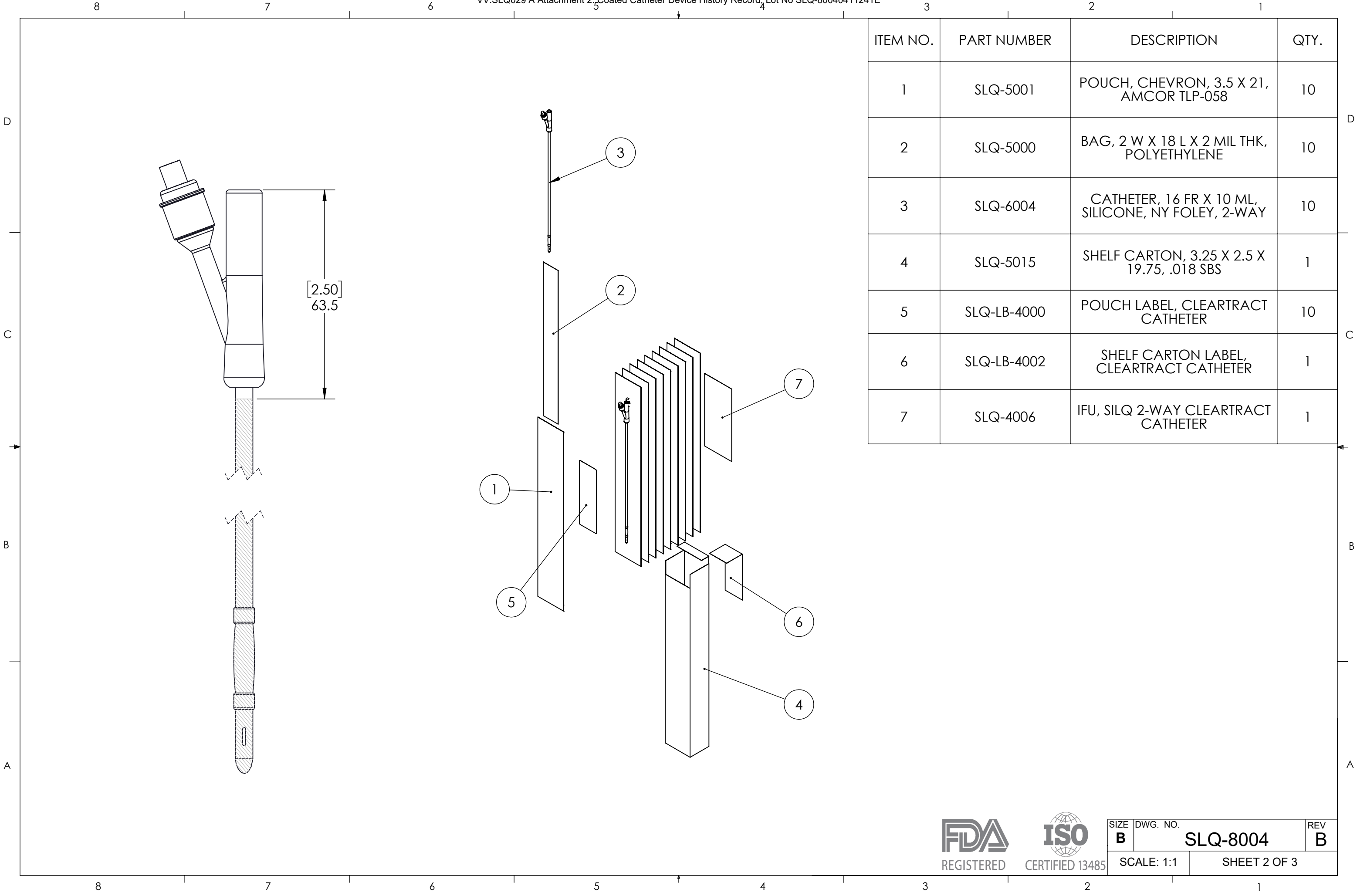

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
1	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	10
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	10
3	SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	10
4	SLQ-5015	SHELF CARTON, 3.25 X 2.5 X 19.75, .018 SBS	1
5	SLQ-LB-4000	POUCH LABEL, CLEARTRACT CATHETER	10
6	SLQ-LB-4002	SHELF CARTON LABEL, CLEARTRACT CATHETER	1
7	SLQ-4006	IFU, SILQ 2-WAY CLEARTRACT CATHETER	1

ITEM NO.	PARTNO	DESCRIPTION	QTY.
8	SLQ-5012	BOX, 20 IN. L X 18 IN. W X 6 IN. D, 200 LB TEST	1
9	PWY-5000	TAPE, PACKAGING, 2 MIL, 2 IN X 110 YDS, CLEAR	AR
10	SLQ-LB-4001	SHIPPER BOX LABEL, CLEARTRACT CATHETER	2

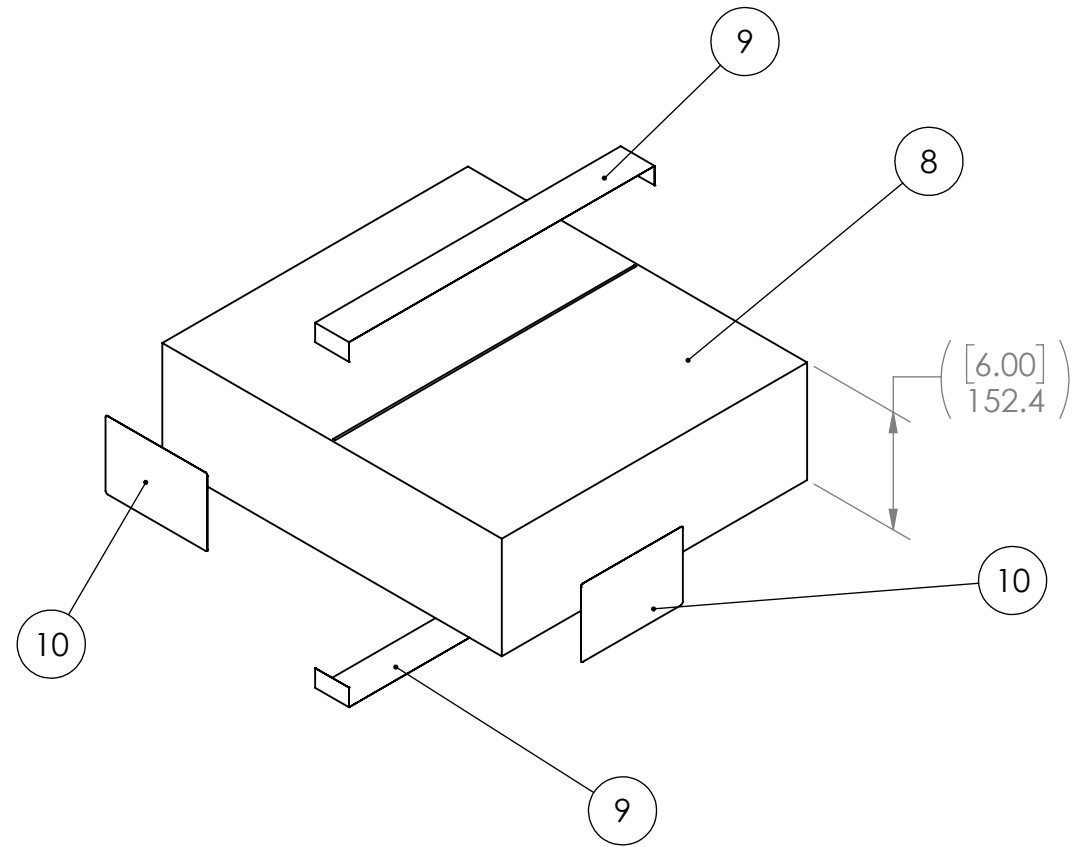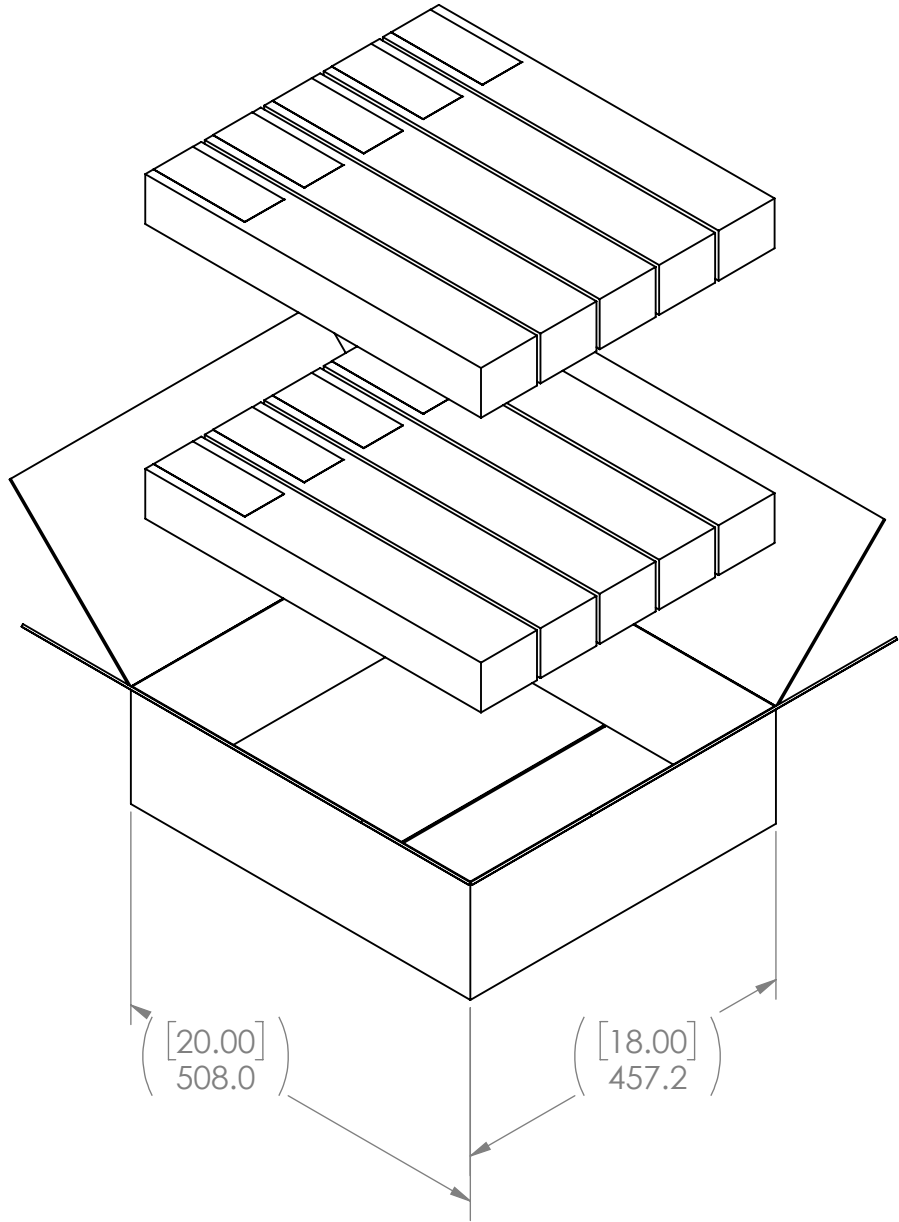

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 1 of 12

1.0 Purpose:

- 1.1 Document the process for the Hydrophilix Polymer Coating application to silicone products (i.e. catheters) utilizing the Hydrophilix UV Reactor.

2.0 Scope:

- 2.1 This procedure applies to the application of the Hydrophilix Polymer Coating to catheters utilizing the Hydrophilix UV Reactor.

3.0 References

- 3.1 WI-0012 Gowning and Prohibited Materials, Clean Room
- 3.2 WI-0014 Line Clearance
- 3.3 WI-0018 Gowning & Prohibited Mtl, For Assy/Packing of Silicone...

4.0 Responsibilities:

- 4.1 It is the responsibility of the individual(s) conducting these activities to comply with the requirements indicated in this procedure.
- 4.2 Clean Room must be thoroughly wiped down with diluted Clorox or 70% isopropyl alcohol and tacky roller should be used on all working surfaces to minimize the presence of lint/debris.

5.0 General Requirements:

- 5.1 Line clearance shall be conducted per WI-0014.
- 5.2 UV Safety Glasses are required at ALL TIMES when the Hydrophilix UV Reactor is in use.

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 2 of 12

6.0 Procedure – Hydrophilix UV Reactor Setup:

6.1 Setup the Hydrophilix UV Reactor as shown in Figure 1, Figure 2 and Figure 3.

1. Plug the Hydrophilix UV Reactor power cord into a 115V ± 10V outlet.
 2. Check to ensure that the 3 UV Light Switches are clockwise until knob stops and the Aluminum Cap is on top of the system.
 3. Toggle switch to the “ON” position (numbers on display panel will illuminate).
- Numbers on the “SV” display will start to countdown from 60 (sixty) seconds to 0 (zero) seconds.
- UV lights in the Reactor will automatically shut down when “SV” displays 0 (zero).

NOTE!

 - TOP ALUMINUM CAP SHOULD ONLY BE REMOVED WHEN THE “SV” DISPLAY INDICATES “0” (ZERO) TO MINIMIZE EXPOSURE TO UV LIGHT.

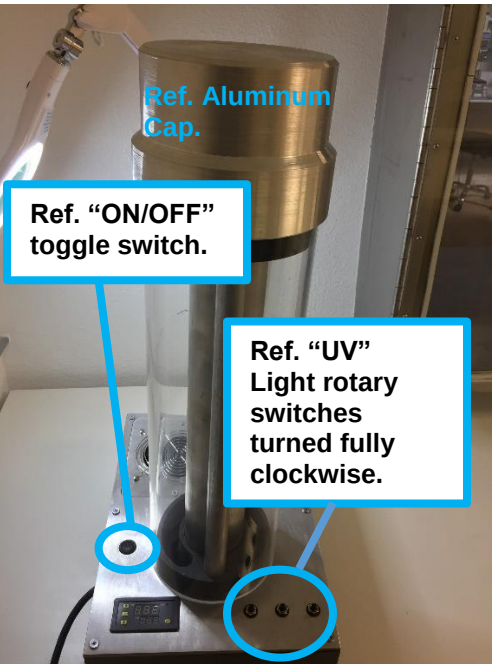


FIG. 1 – Hydrophilix UV Reactor

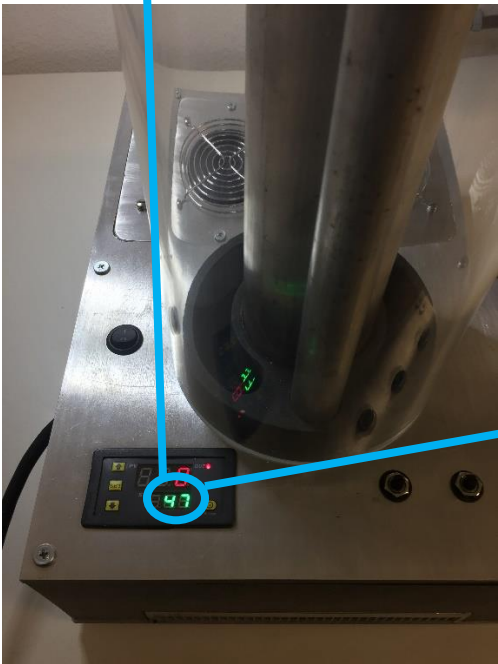


FIG. 2 – Count down display



FIG. 3 – UV Lights OFF when display “0” (Zero)

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 3 of 12

7.0 Procedure - Hydrophilix Polymer Coating Application:

7.1 Apply and UV cure the Hydrophilix Polymer Coating P/N HYDPLX-3000 to the 14Fr Silicone Catheter P/N HYDPLX-6000 or the 18Fr Silicone Catheter P/N HYDPLX-6001 as shown in Figure 4 through Figure 19.

1. Check to ensure that the cap on the Hydrophilix Polymer Solution P/N HYDPLX-3000 bottle is secured to the bottle.
2. Shake bottle to thoroughly mix the liquid contents in the bottle.

NOTE

- Check the bottle to ensure that thicker viscous material on the bottom of the bottle is not present.
- Check to ensure there is no separation between the liquid contents in the bottle. Mixture should have a milky white appearance.

1. Use a clean 10 mL syringe with attached 5-inch tube to draw approximately 8 mL of the Hydrophilix solution from the bottle as follows.
 - a. Make sure syringe plunger is pressed in until it bottoms out in syringe body.
 - b. Place the tubing that is attached to the syringe into the solution.
 - c. Slowly pull up on the plunger to draw the solution into the syringe (stop pulling when desired measurement).

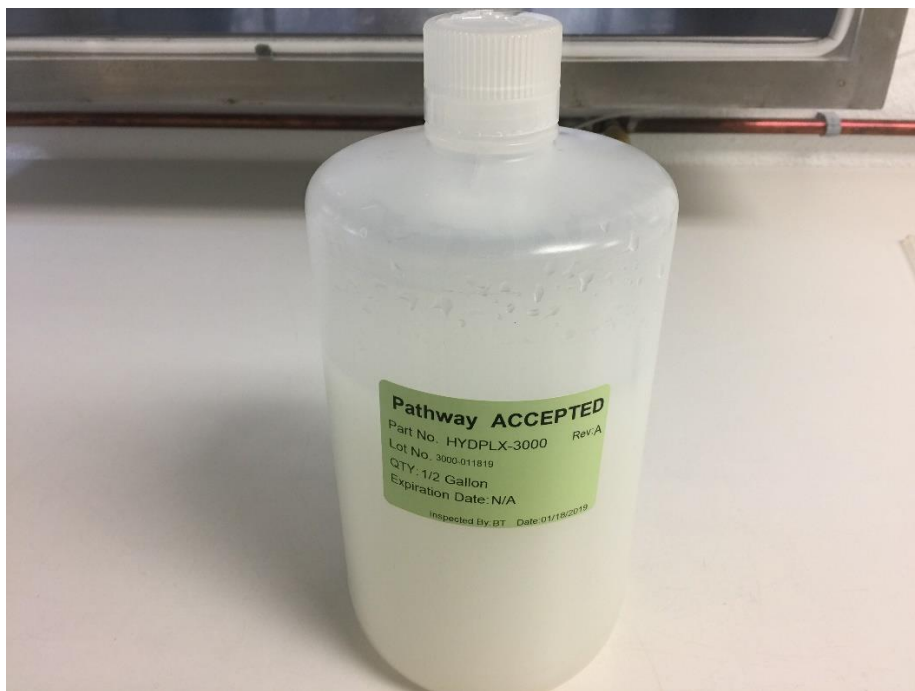


FIG. 4 - Mix Hydrophilix Polymer Solution

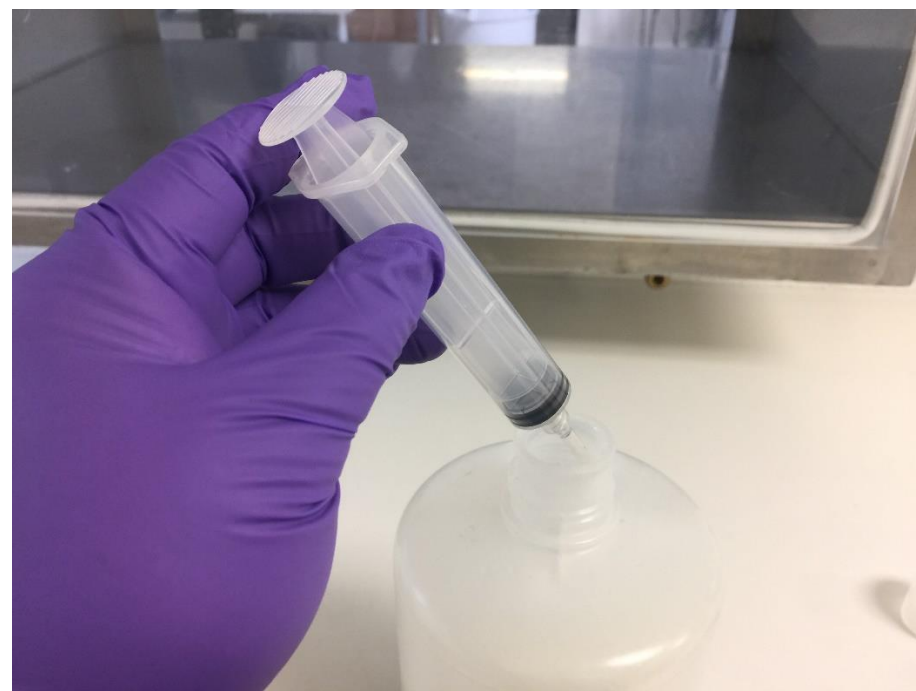


FIG. 5 – Drawing Solution into Syringe

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 4 of 12

1. Place/set the appropriate Quartz Tube (use larger diameter Quartz tubes for the 18 Fr Catheters) onto the holding rack.
2. Insert Syringe discharge tubing into the Quartz Tube.
3. Depress the syringe plunger to dispense the Hydrophilix Solution into the Quartz Tube.



FIG. 6 – Hydrophilix Solution dispensed into Quartz Tube

1. Fully insert the appropriate size Catheter P/N HYDPLX-6000 for 14 Fr or HYDPLX-6001 for 18 Fr into the Quartz tube.



FIG. 7 – Catheter inserted into Quartz Tube w/Solution

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 5 of 12

1. Check to ensure the “SV” display reading is at “0” (Zero).
2. Remove the Aluminum Cap from the top of the UV Reactor



FIG. 8 – Aluminum Cap removal

- View from the top of the UV Reactor
- Ref. Hole noted with **GREEN CIRCLE** is location in which the Quartz Tube is to be inserted.



FIG. 9 – View denoting Quartz tube hole

1. Insert the Quartz Tube w/Catheter and Hydrophilix Solution into the center hole on the top of the UV reactor (ref. Figure 9).
2. Check to ensure that the Quartz Tube is seated flush to the top of the UV Reactor (ref. Figure 11 and Figure 12).



FIG. 10 – Quartz Tube insertion

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 6 of 12

INCORRECT

1. Green Handle on the Quartz Tube is not seated flush (gap between bottom of the Green Handle and the black surface of the UV Reactor) to the surface of the UV Reactor.

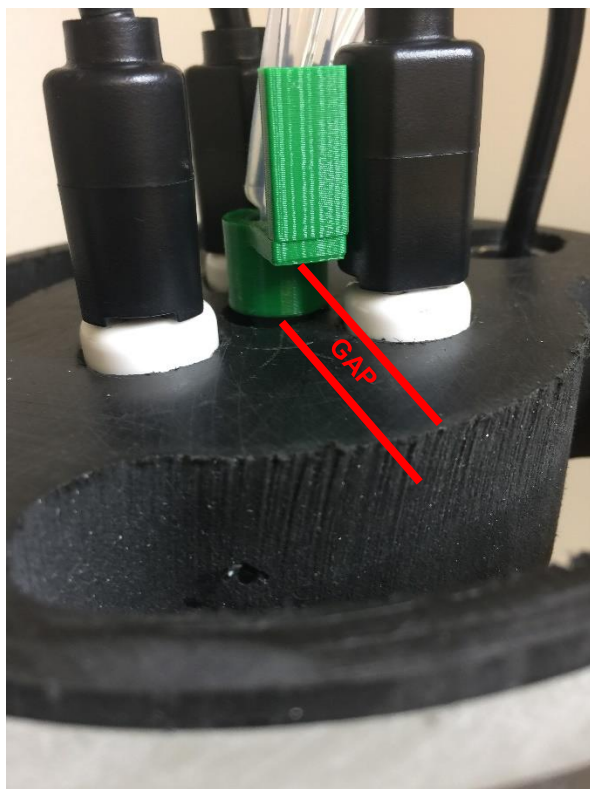


FIG. 11 – Quartz tube handle not seated

CORRECT

1. Green Handle on the Quartz Tube is seated flush (no gap between bottom of the Green Handle on the Quartz Tube and black surface of the UV Reactor) to the top surface of the UV Reactor



FIG. 12 – Quartz Tube handle seated flush

1. Place/set the Aluminum Cap on top of the UV Reactor.
2. Make sure that the Aluminum Cap is fully seated on top of the UV Reactor.



FIG. 13 – Aluminum Cap re-installed

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 7 of 12

1. Press the “Restart” switch on the display panel to begin the 1-minute UV Polymer coating process.



FIG. 14 – UV Reactor Start

- Display time will automatically countdown from “60, 59 58...9, 8, 7...” until 1-minute time frame has elapsed.



FIG. 15 – UV Reactor 1-minute countdown

1. Check to ensure that the time display “0”.

NOTE

- “0” (Zero) indicates that the UV Lights in the Reactor are “OFF” and it is safe to remove the Aluminum Cap.

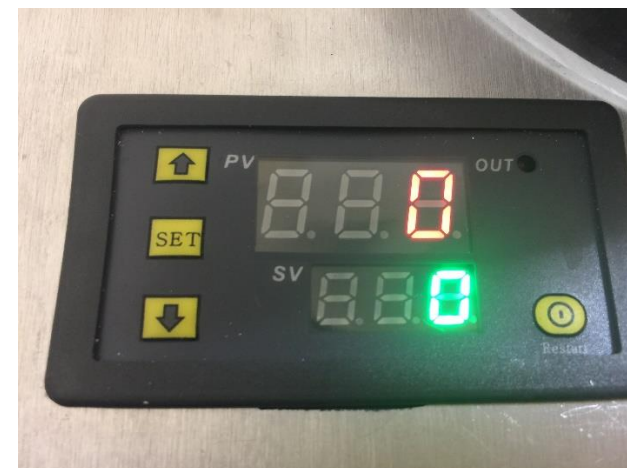


FIG. 16 – End of 1-minute countdown

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 8 of 12

1. Remove the Aluminum Cap from the top of the UV Reactor.



FIG. 17 – Quartz Tube removal

1. Check to ensure that the Hydrophilix Solution is light yellow/orange in color (this indicates that the UV lamps were “ON” and that the coating process was successful.

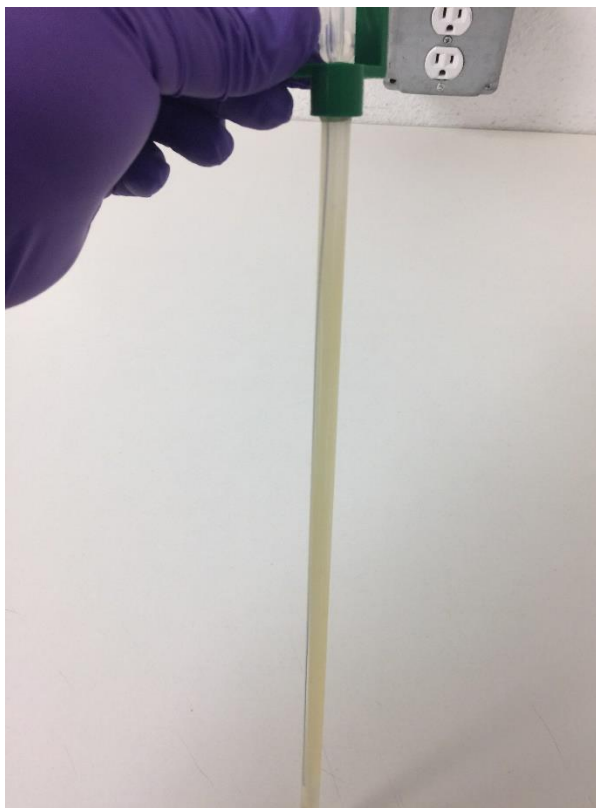


FIG. 18 – Solution color verification

1. Remove the Catheter from the Quartz Tube.

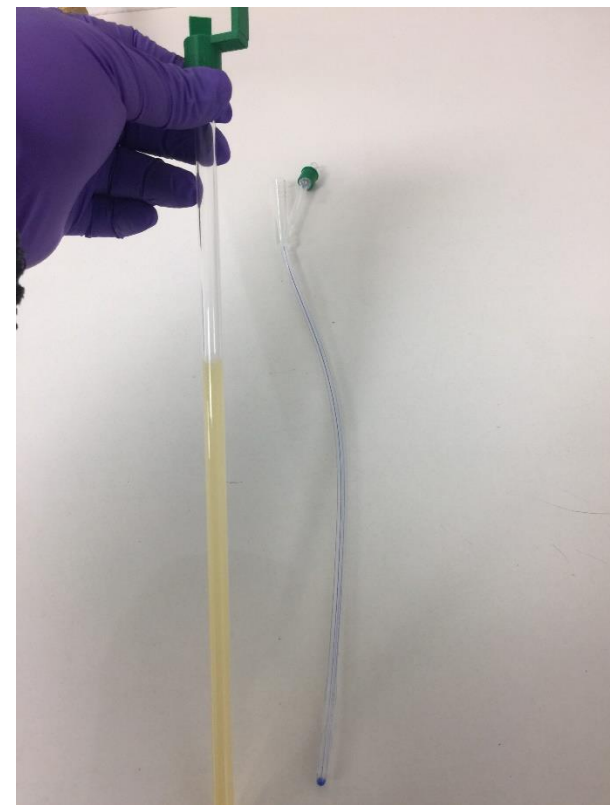


FIG. 19 – Catheter removal

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 9 of 12

8.0 Procedure – Catheter/Quartz Tube Cleaning:

8.1 Rinse the Polymer coated Catheter with De-ionized Water P/N HYDPLX-3001 as shown in Figure 20 and Figure 21.

1. Fill a clean 20 – 30 mL syringe with De-ionized water.

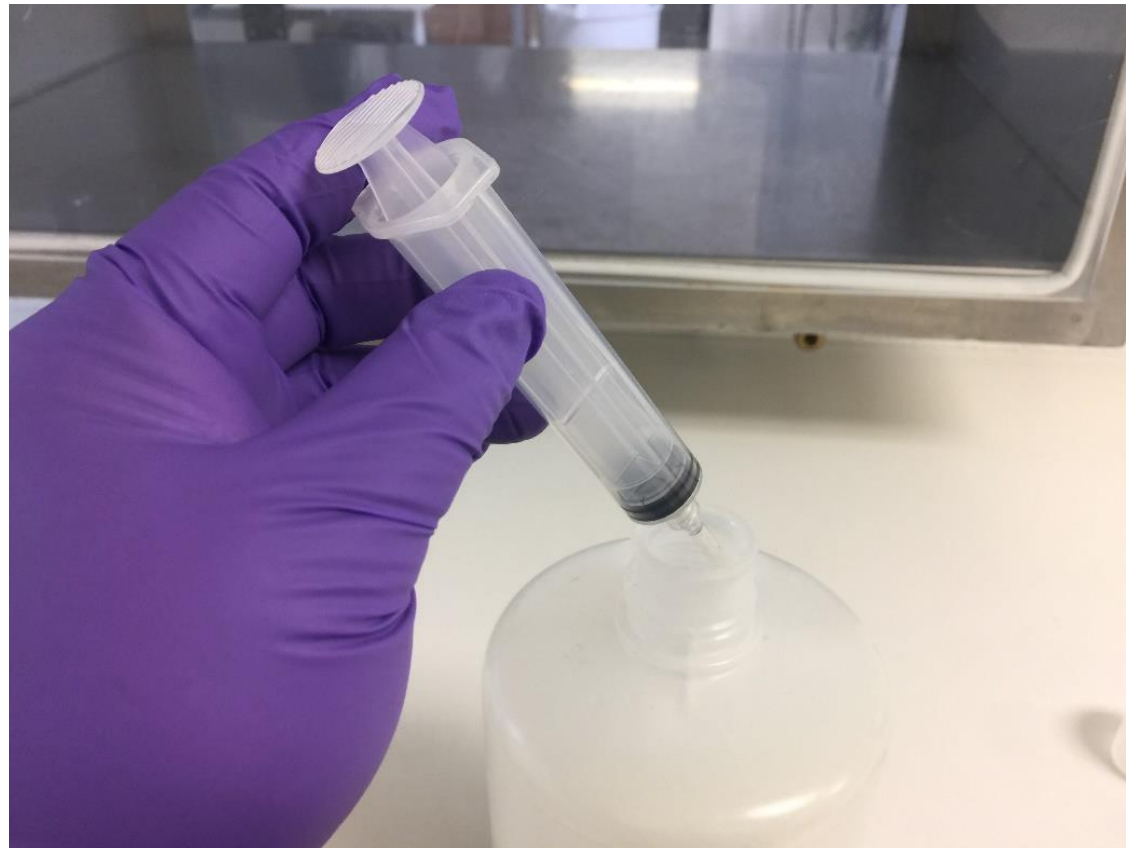


FIG. 20 – Fill Clean Syringe with De-ionized water

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 10 of 12

1. Use the de-ionized water to clean/rinse out the inside and outside of the Catheter (use designated liquid waste container to capture liquid).
2. Hang or place the clean Catheter onto a drying rack or equivalent and allow to air dry.

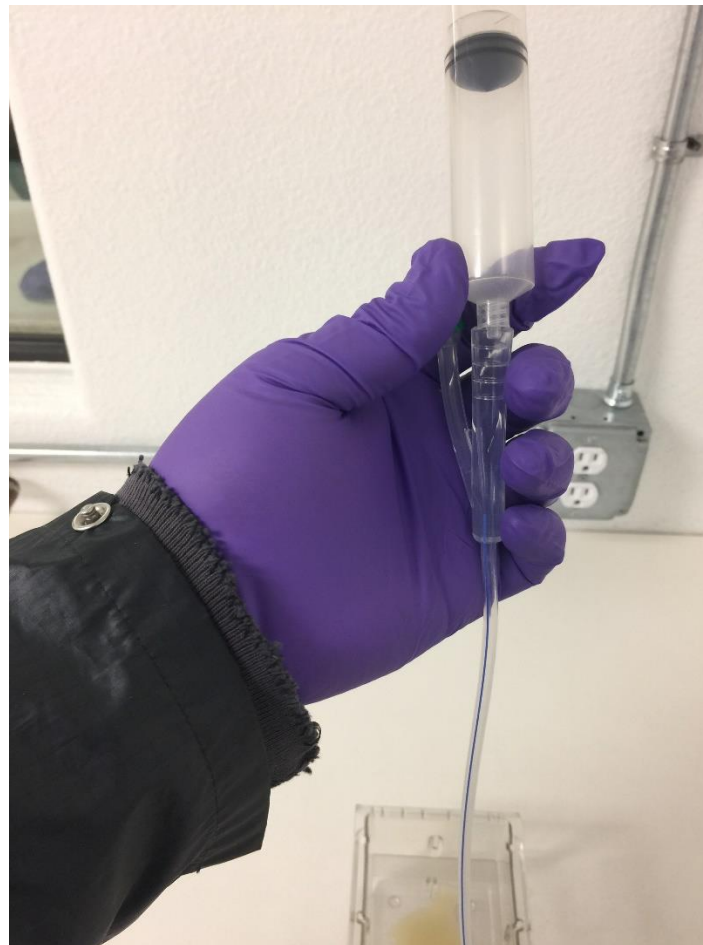


FIG. 21 – Catheter clean/rinse

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 11 of 12

9.0 Procedure – Quartz Tube Cleaning:

- 9.1 Clean/rinse the Quartz Tube using De-ionized Water P/N HYDPLX-3001 and a Pipe Cleaner as shown in Figure 22 and Figure 23.

1. Empty the used Hydrophilix Solution from the Quartz Tube into the designated liquid waste bin.
2. Discard used Hydrophilix Solution.

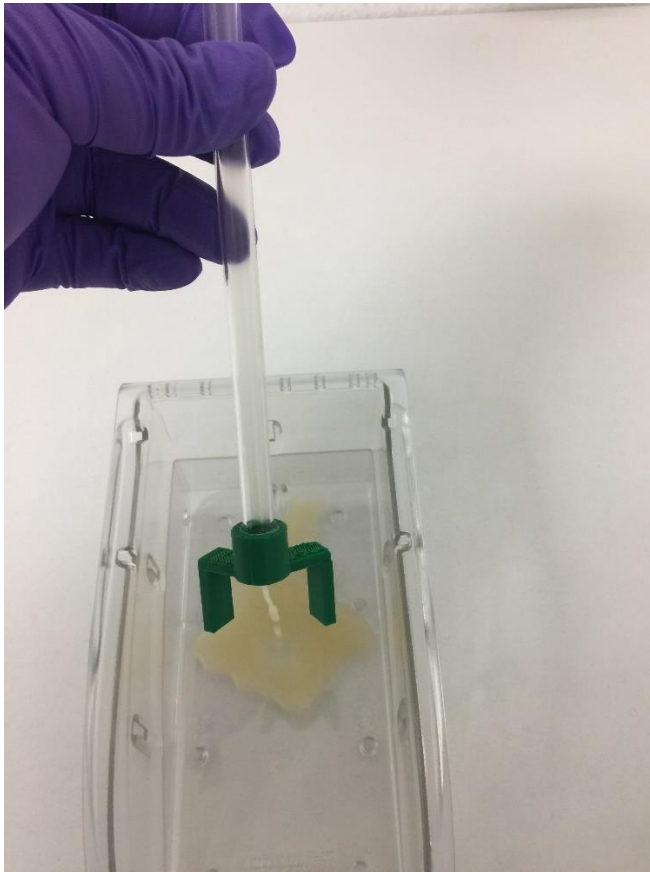


FIG. 22 – Empty Quartz Tube of used liquid

1. Use a pipe cleaner and De-ionized water to clean the inside of the Quartz Tube.
2. Discard used De-ionized water and pipe cleaner.

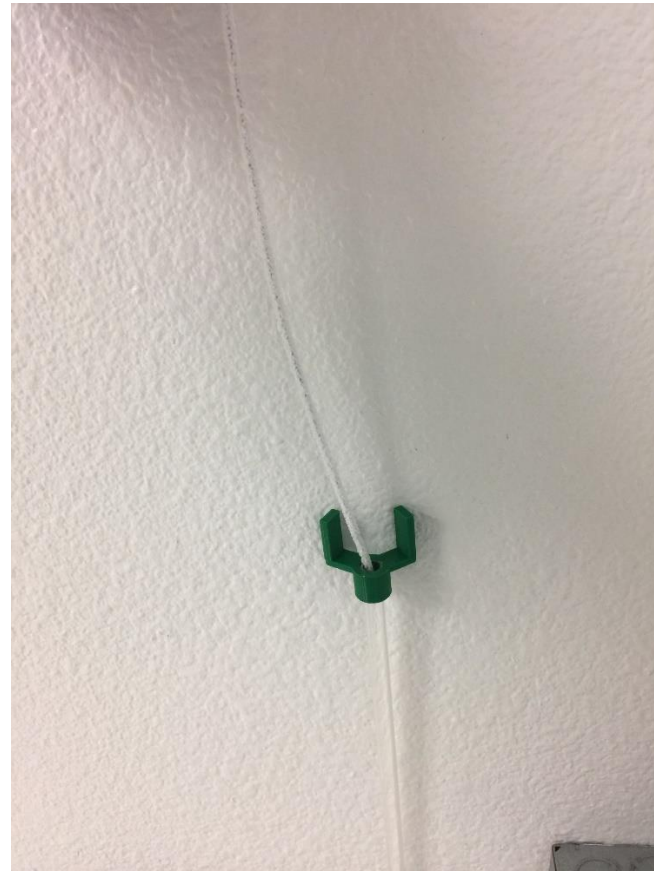


FIG. 23 – Use De-ionized Water and pipe cleaner

Pathway Innovation LLC.	Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing UV Reactor	Doc No.: MP- HDX-0000
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 12 of 12

10.0 Change History

Revision	Change Description
A	Initial Release

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 1 of 27

1.0 Purpose:

- 1.1 Document the process for the labeling and packaging of the Silq ClearTract Catheters.

2.0 Scope:

- 2.1 This procedure applies to the labeling and packaging of the Silq ClearTract Catheters.

3.0 References

- 3.1 WI-0012 Gowning & Prohibited Materials, Clean Room
- 3.2 WI-0018 Gowning & Prohibited Materials, for the Assy/Packing of Silicone or any Tacky/Sticky...
- 3.3 WI-0014 Line Clearance

4.0 Responsibilities:

- 4.1 It is the responsibility of the individual(s) conducting these activities to comply with the requirements indicated in this procedure.
- 4.2 Clean Room must be thoroughly wiped down with diluted Clorox or 70% isopropyl alcohol and tacky roller should be used on all working surfaces to minimize the presence of lint/debris.

5.0 General Requirements:

- 5.1 Line clearance shall be conducted per WI-0014.

6.0 Material and Equipment:

- 6.1 ACCU-SEAL Impulse Sealer or equivalent Validated Continuous Band Sealer

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 2 of 27

7.0 Procedure – Lot/Expiration Verification:

- 7.1 Verify that the information printed on the on the pouch, box, and shipper box label aligns with the Work Order as shown in Figure 1 Figure 2.

1. Check to ensure that the Catheter Fr size, Quantity, Lot Number, REF Number and Expiration Date that is printed on the Pouch and Box Label matches the recorded information on the Work Order. The pouch label will have a QTY of 1 and the box label will have a QTY of 10. Verify the lot number date reference is in the form MMDDYY and the expiration date is in the form YYYY-MM-DD (Ex: 2023-01-21). Verify item (01) is a fourteen-digit number representing the GTIN value. Verify item (10) matches the lot number. Verify item (17), which is tied to the expiration date is in the form YYMMDD. The below label is for a 16 Fr, 10 mL catheter expiring on January 21, 2023. After the first few labels, there is no need to keep verifying the information. View Appendix A for the pouch and box artwork that is applied to the SPT catheters



2-WAY 100% SILICONE CLEARTRACT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	1

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex



SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

R_x Only



SLQ-MMDDYYYY



XXXXXX



20YY-MM-DD



(01) #####
(10) SLQ-MMDDYYYY
(17) YYMMDD

SLQ-LB-4000 Rev. A

FIG. 1 – Pouch Label

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 3 of 27

QUANTITY	SIZE	 SILQ	R_x Only      
10	XX Fr, XX mL		 (01) ##### (10) SLQ-MMDDYYYY (17) YYMMDD

Symbols are defined in Instruction for Use (IFU)
 Store away from direct light exposure at room temperature
 Not made with natural rubber latex

 SILQ Technologies Corp.
 323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

LOT SLQ-MMDDYYYY

REF XXXXXX

 20YY-MM-DD

SLQ-LB-4002 Rev. A

FIG. 2 – Box Label

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 4 of 27

1. Check to ensure that the REF, LOT, Expiry Date printed on the Box Label matches the recorded information on the Work Order and match the pouch and 10-pack box labels. Verify the lot number date references and expiration dates are in the format below. Also refer to the image above for more detailed descriptions. View Appendix A for the shipper box artwork that is applied to SPT catheters



2-WAY 100% SILICONE CLEARTRACT CATHETER

REF

XXXXXX

LOT

SLQ-MMDDYYYY



20YY-MM-DD

QTY

100



SILQ Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4001 Rev. A

FIG. 3 – Shipping Box Label

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 5 of 27

8.0 Procedure – Inner Sleeve/Pouch Packaging and Sealing:

- 8.1 Insert the Polymer coated Silicone Catheter into the Catheter Sleeve and Pouch and seal as shown in Figure 4 through Figure 9.

1. Orient and apply the verified Pouch Label (ref. Figure 1) to the Polyethylene side of the Pouch P/N HDX-5001 such that the text reads away from the chevron end of the Pouch.

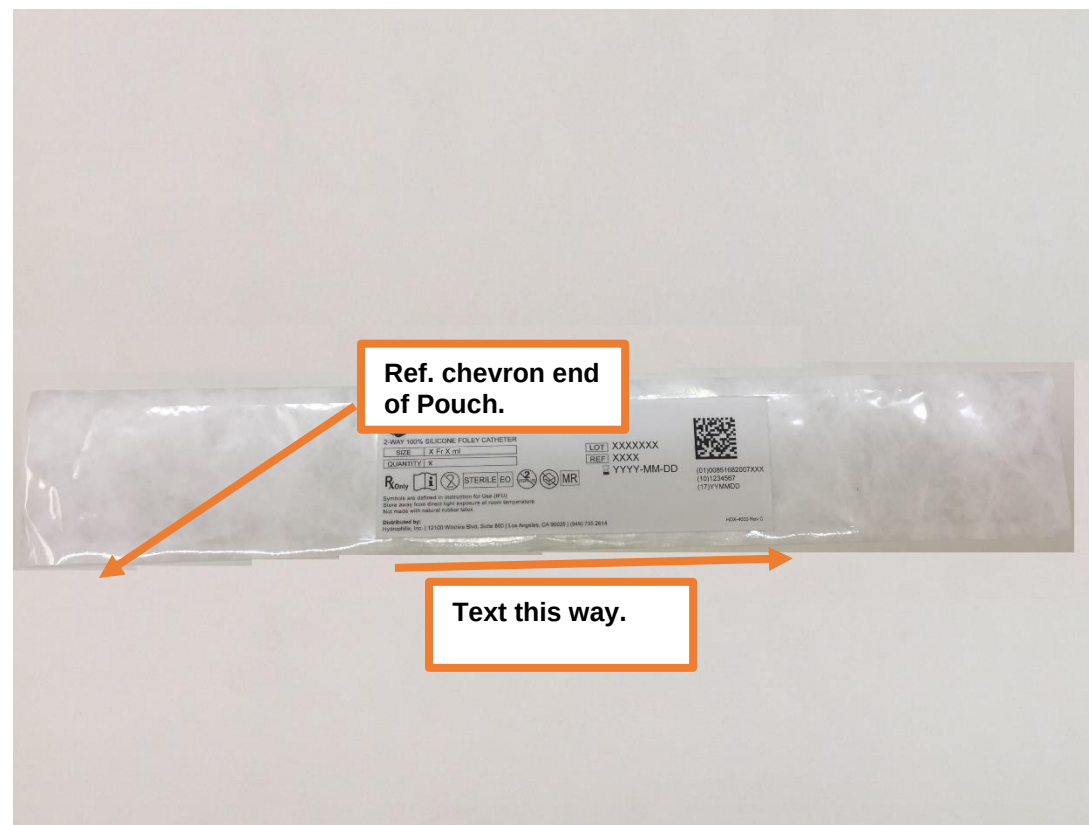


FIG. 4 – Label applied to Pouch

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
		Rev.: B
		Page 6 of 27

1. Visually inspect the Catheter (see Work Order to determine P/N) to ensure that the outer surface is free of lint/debris (silicone is tacky and has a tendency attract dust/lint).

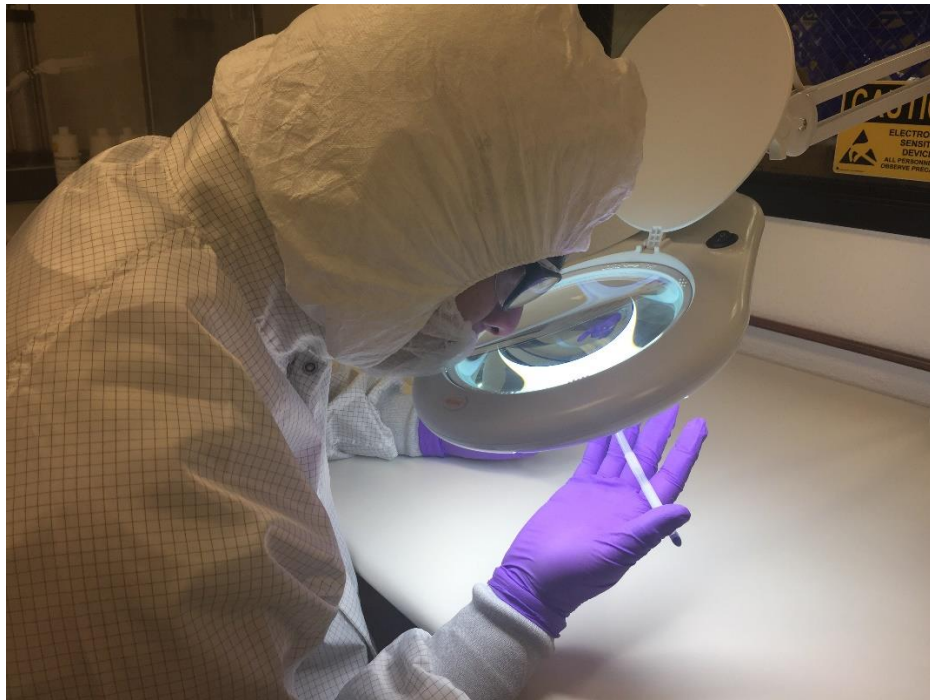


FIG. 5 – Catheter/Polyethylene Bag cleanliness verification

1. Check to ensure that the inside and outside of the 2" W x 18" L Polyethylene Bag is free of lint/debris (discard if necessary).
2. Insert the Catheter (tip end first) into the Polyethylene Bag until the tip of the catheter bottoms out to the closed end of the bag.



FIG. 6 – Catheter insertion to Polyethylene Bag

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
		Rev.: B
		Page 7 of 27

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

1. Visually inspect the sleeved Catheter to ensure that the outer and inner surface is free of lint/debris.
2. **Check inside of the Labeled Pouch to ensure it is free of any foreign matter/debris.**
3. **Insert sleeved catheter (closed end of sleeve) into the pouch until the sleeve bottoms out to the chevron end of the pouch.**



FIG. 7 – Sleeved catheter/pouch cleanliness verify & insert

1. **Visually inspect the pouched catheter to ensure that NO debris was introduced when the sleeved catheter was inserted.**



FIG. 8 – Pouched Catheter cleanliness inspection

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 8 of 27

1. Set/position the pouched catheter with the Tyvek (white side) facing up into the ACCU-SEAL Impulse Sealer.



FIG. 9 – View, sealing of pouched catheter utilizing impulse sealer

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: B
			Page 9 of 27

9.0 Procedure –10-Pack Inner Box and Outer Shipping Box Packaging and Labeling:

9.1 Assemble and insert the individually sealed Catheter Pouches into 10-Pack Inner Box into the Outer Shipping Box as shown in Figure 10 through Figure 15.

1. Place 10 catheters into the box with the text on the pouch label oriented as shown below.
2. Prior to boxing, check to ensure that there **ARE NO VISIBLE** signs of foreign matter (i.e. lint, hair....) in any of the sealed catheter pouches.

1. Insert the IFU into the loaded 10-pack box.

1. Close box flaps.
2. Apply the box label. It will fold over the blank edge of the box. Apply it such that the SILQ logo, text, size, and quantity are all vertical on the short edge of the box.
3. The rest of the label will be on the top face.
4. Verify the label is applied evenly and is free of large wrinkles.
5. If applied correctly, the label will not cover any pre-printed artwork.

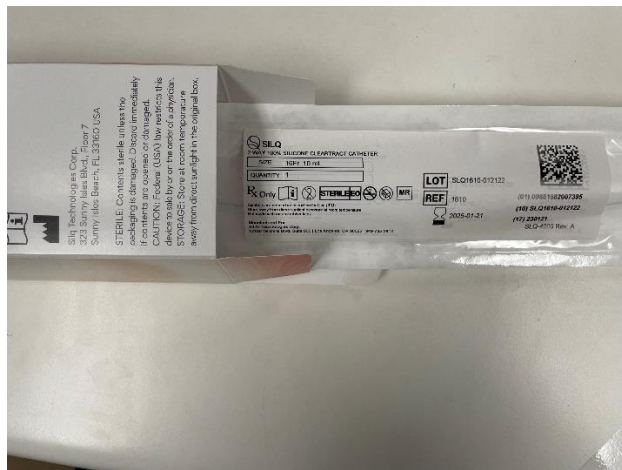


FIG. 10 – 10-Pack Inner Box Loading



FIG. 11 – IFU Loading



FIG. 12 – Box Label Application

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 10 of 27

1. Assemble the Shipping Box using Packaging Tape or equivalent.
2. Apply shipper box labels to adjacent sides of the box as shown.



FIG. 13 – Box and label application

1. Load a maximum of ten 10-Pack Catheter boxes into the shipping box.



FIG. 14 – 10-Pack Catheter Boxes load to shipping box

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 11 of 27

1. IF NECESSARY, FILL VOIDS INSIDE BOX with Bubble Wrap or equivalent.
2. Seal box using packing tape.
3. Route labeled box to Finished Goods Inventory (FGI) or shipping.



FIG. 15 – Seal shipping box and route to Finished Goods or Shipping

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 12 of 27

10.0 ACCU-SEAL Impulse Sealer (ES-0099-01) Equipment Setup:

- 10.1 Reference Figure A through Figure D for the ACCU-SEAL Model 5300-15-B Impulse Sealer Controls ([See Step 10.2 for validated operating ranges](#)).

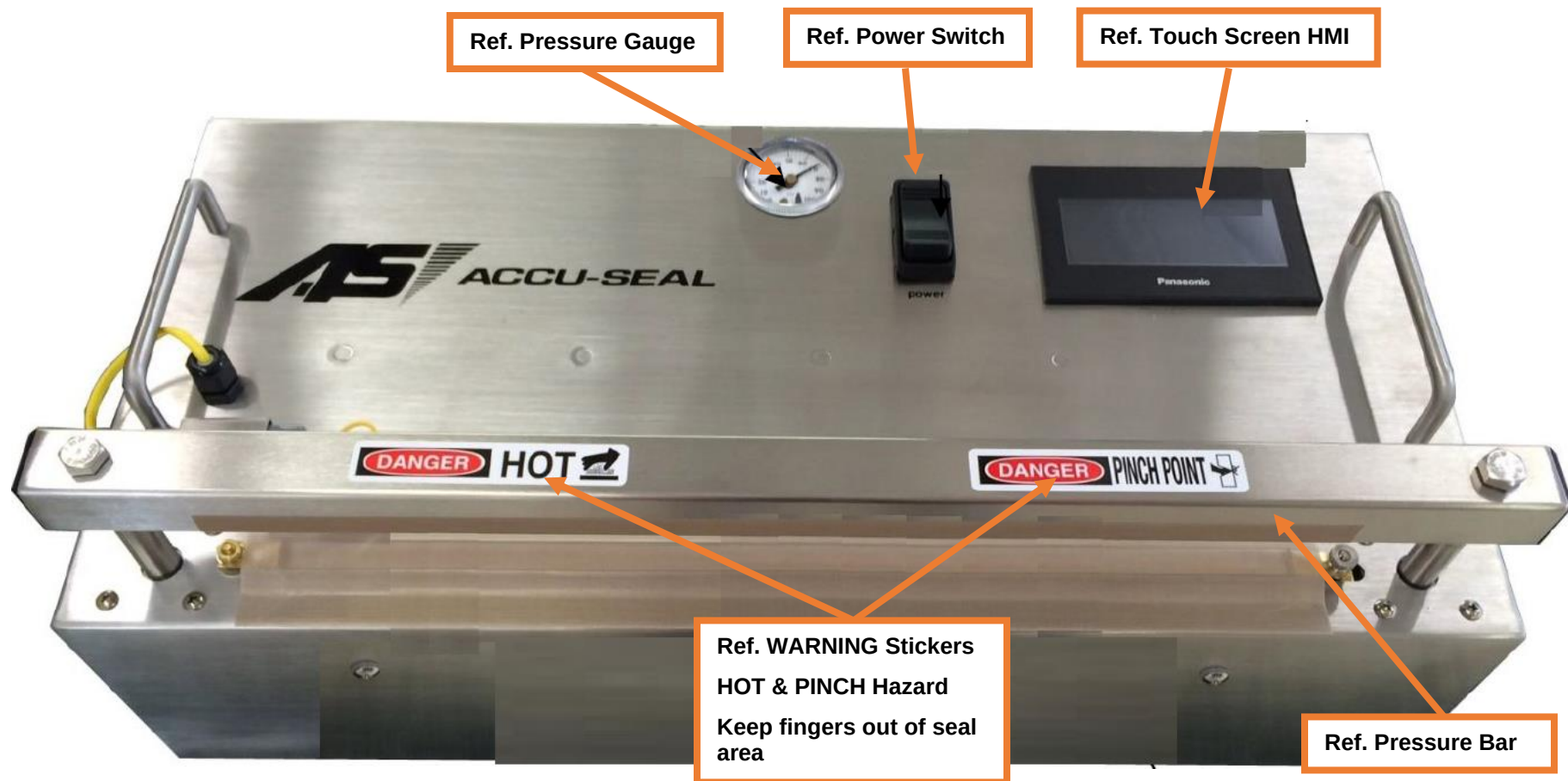


FIG. A – View from Front of Impulse Sealer

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 13 of 27

1. Connect Foot Pedal, Power cords to appropriate ports
2. Connect compressed air line to the "AIR IN" port.
3. Use Pressure Regulator to adjust pressure to the Validated Range of 55 psi.

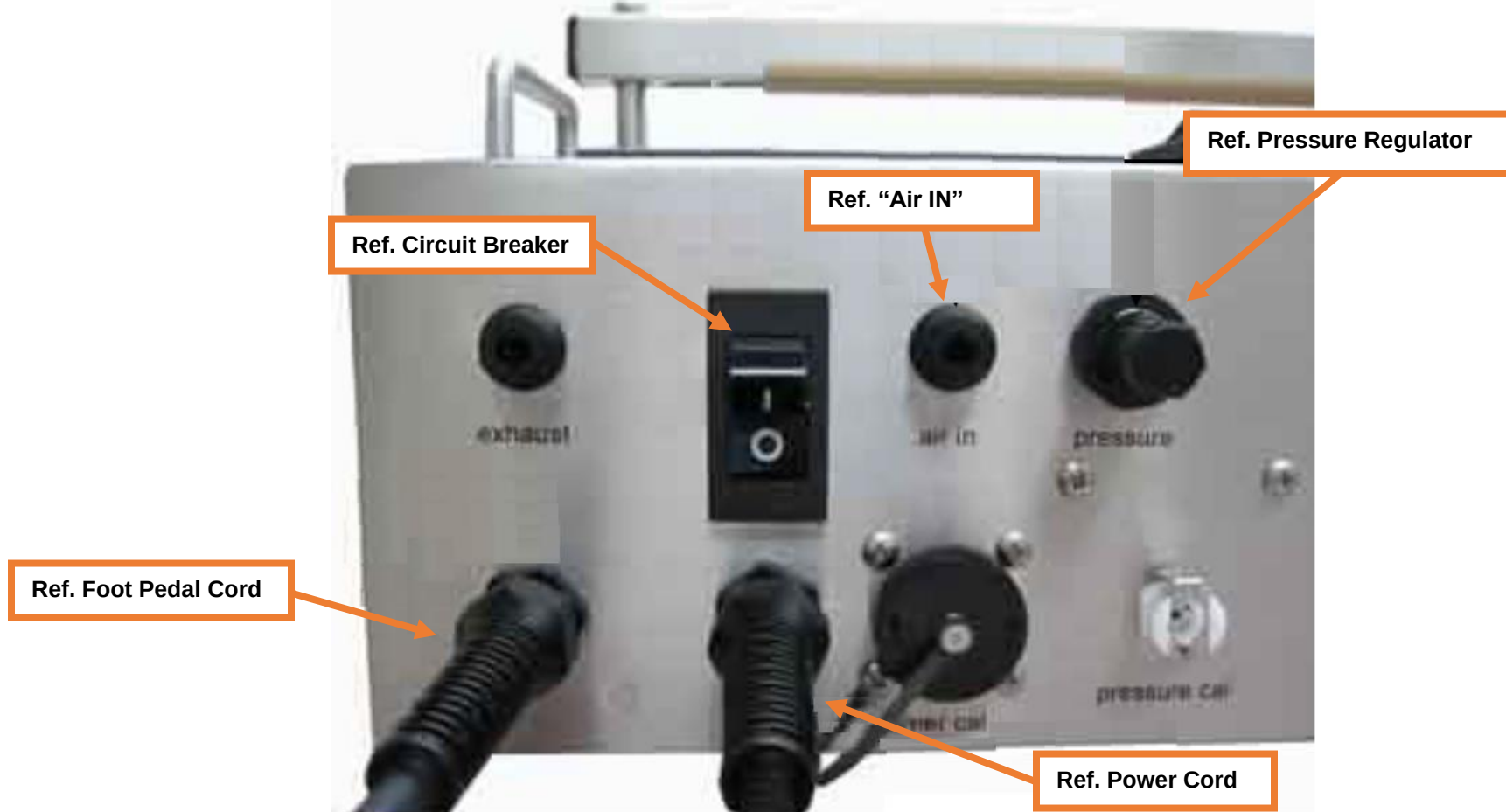


FIG. B – View, Back of Impulse Sealer

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 14 of 27

CONTROLS

Overview

The machine is programmed for two modes of control. Non-access code operation and Password Protected access. Non-access code operation only allows the operator to run the set seal parameters and to reset alarms.

The touch screen allows Authorized personnel, utilizing an Access code (password), access for editing all of the machine seal parameters, alarm settings, pressure settings, as well as the Calibration and Maintenance screens.

Controls and their Function:

- **Power Switch:** Turns power to the machine ON and OFF.
- **Foot Switch:** The Foot Switch is used to start the sealing cycle.
- **Circuit Breaker:** Breaks the Main Power to the machine when there is an overload condition.
- **Touch Screen HMI:** The Touch Screen HMI is used to Monitor, Set or Edit all of the machine parameters.

Controls Operation:

- Turn the machine ON by pressing the Main Power Switch located on top of the machine. When the Power is ON, the touch screen will be illuminated. To turn OFF Main Power, press the bottom side of the Main Power Switch.
- Current machine seal parameters are displayed on the main Run Screen. To edit or change machine, press the "Edit" button. When the Edit button is pressed, the "Enter Password" Screen will appear. The correct Access Code (password) must be entered to continue. When the correct code is entered you will have access to the functions screens and have the ability to change the machine parameters.

NOTE:

- Initial Set-Up Pass code has been programmed into the machine. See Initial Set up instruction in the maintenance Section. Use this code for first time access to set New password.

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 15 of 27

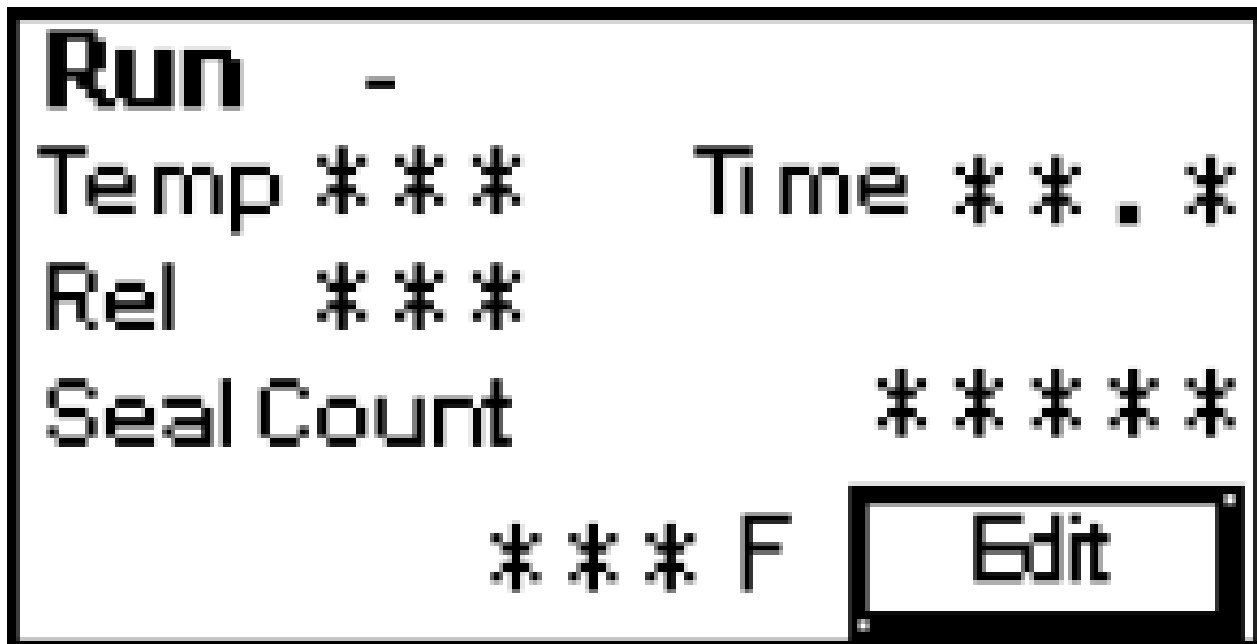


FIG. H – Touch Screen Functions



FIG. C – Edit Menu Screens

- **TEMP:** Displays set Seal Temperature
 - **REL:** Displays set Release Temperature
 - **TIME:** Displays set Seal Time
 - **SEAL COUNT:** Displays the number of heat seals completed
 - *****F:** Active temperature display. Display indicates the current temperature of the element.
 - [EDIT: Allows AUTHORIZED PERSONNEL, UTILIZING ACCESS CODE \(PASSWORD\), ACCESS TO ALL OF THE MACHINE SETTINGS:](#)
 - Seal Parameters
 - Change Password
 - Reset Seal Count
 - Alarms
 - Calibration Timer Test
- Edit Menu Screens**
- Display the functions that can be selected for Editing.

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 16 of 27

Seal Parameters

- This screen will allow you to change the seal parameters for:

- SEAL TEMP: Temperature held for the time specified by the Seat Time (Dwell).
- SEAL TIME: The amount of time the Seat Temperature is maintained.
- RELEASE: Temperature the seal bar releases and the cycle resets.

Other Function Screens

- Change Password
- Reset Seal Count
- Timer Test for Calibration
- Coyote Timer Setting

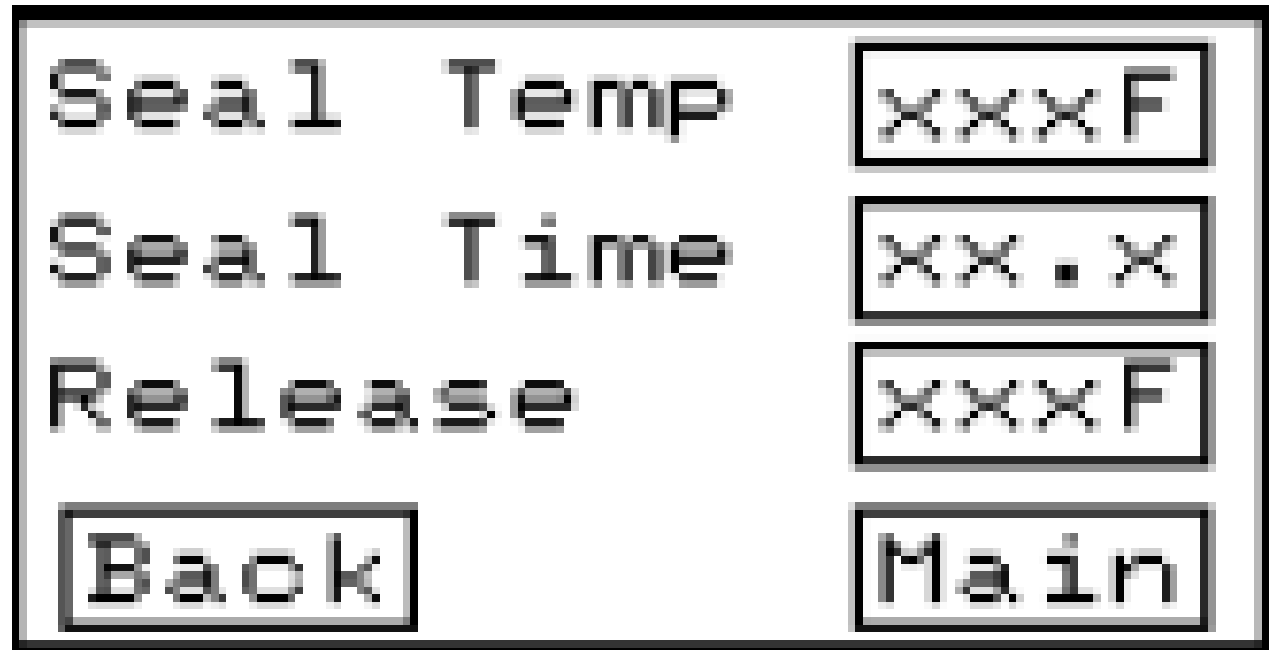


FIG. D – Seal Parameter Screens

Alarms

- JAW OBSTRUCITON: Indicates that for some reason the Seal Bar will not close completely.
- COLD SEAL: If the temperature does not reach the Seal Temperature within the set cold seal time, alarm will activate.
- OPEN SENSOR: Thermocouple is disconnected or has an open fault.
- OVERTEMPERATURE FAULT: Is a heater power control warning and indicates abnormal operation of the temperature control hardware.

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 17 of 27

ENTERING/EDITING MACHINE PARAMETERS: Any parameter that is enclosed in a box can be change. To change the parameter, press anywhere in the box of the parameter you want to change. This will then bring up a “keyboard” that will allow you to enter a new value. Once the desired value is entered, either the box of the parameter you are changing can be pressed or the “enter button on the Keyboard can be pressed to save the new value into memory. (Enter=the return arrow key) Parameter that are not enclosed in a box are for viewing only.

CHANGING ACCESS CODE: The Access Code can be changed, by pressing the “Change Password” button. To change the Access Code, press on the box with 00000. Enter a new Access Code. The number must be 5 digits and cannot start with ‘0’. The number you entered will show in the box. If the number is correct, Press the “ENTER” button to save the new Access Code to memory. The machine has been programmed with a Master Default Access Code in case your code is lost or forgotten. Please call if needed.

SEAL COUNT FUNCTION: The machine counts every seal that is started and successfully completed and displays this on the Main Screen. The machine will count up to 99,999 seals before it is automatically reset to “0”. The seal count can also be manually reset through the Seal Count Reset function.

PRESSURE: When the Regulator on the back of the machine is adjusted, the pressure gauge on the top of the machine will read the set pressure. Typical operating pressure should be above 50 psi.

JAW OBSTRUCTION ALARM: Indicates that for some reason the seal jaw will not close completely. If the alarm is activated, the machine will go into a lock-up mode. The machine will not operate until condition is corrected and the alarm is reset, by pressing the box that reads “ALARM RESET”, on the flashing red alarm screen. The Main Menu will now appear.

COLD SEAL ALARM: If the temperature does not reach the Seal Temperature within the set Cold Seal time, alarm will activate. If the alarm is activated, the machine will go into a lock-up mode. The machine will not operate until condition is corrected and the alarm is reset, by pressing the box that reads “ALARM RESET”, on the flashing red alarm screen. The Main Menu will now appear.

OPEN SENSOR ALARM: The machine continuously monitors the status of the machine temperature sensor. If the sensor is disconnected or has an open fault, the machine will alarm with the open sensor alarm screen. The sensor monitor is active whenever the machine is on so the alarm can occur even when the machine is in the resting state. The machine will not operate until condition is corrected and the alarm is reset, by pressing the box that reads “ALARM RESET”, on the flashing red alarm screen. The Main Menu will now appear.

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 18 of 27

COYOTE TIMER SETTING: The internal PLC monitors the health of the seal process to detect any abnormal operation of the temperature control hardware. This monitoring action is called a “coyote”. The coyote may be implemented in the hardware or software, or a combination of both. It is not a safety feature, but acts to detect abnormal temperature control operation and shut down the machine. If an abnormal temperature event occurs, the Over Temperature Fault Alarm will display on a flashing red screen. Contact Maintenance, cycle the power of the machine to reset the alarm. If the alarm occurs again, CALL ACCU-SEAL

NOTE: See Over Temperature Fault Alarm.

TIMER TEST: This screen is used during calibration for testing the timer function.

OPERATING PROCEDURES

Start UP

1. Verify that the Air supply is connected.
2. Turn the machine ON by pressing the Power Switch.
3. The Touch screen will light up ‘green’. When the boot-up sequence is complete, you will see the Main Run” Screen.
4. Verify that the nuts at each end of the element are tight.
5. View the Main Menu Screen and verify that the temperature and pressure readouts are reasonable.
6. Keeping your fingers away from the Seal Area, press and release the Foot Switch. Cycle the machine once or twice to warm it up.
7. Place the product to be sealed underneath the Seal Bar making sure the product inside the pouch is not protruding into the Seal Area.

CAUTION: Small products in the pouch that protrude into the Seal Area may be too small to cause a Jaw Obstruction Alarm, but could cause damage to element or a High temperature condition.

8. Keeping your fingers away from the Seal Area, press and release the Foot Switch. When the Foot Switch is pressed the machine will immediately begin the seal cycle.
9. When the seal cycle is complete the Seal Bar will open.

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 19 of 27

10. Remove the Product and inspect the seal.

11. At the end of each shift or production run, Press the Power switch to turn the machine OFF>

10.2 Validated Operating Range is as follows:

- **Amcor Pouch (SLQ-5001)**
- **Sealing Temperature** **270 – 280 °F**
- **Sealing Dwell time** **4.0 Fixed**
- **Sealing Pressure** **55 psi Fixed**

- **SteriPax Pouch (SLQ-5013)**
- **Sealing Temperature** **250 – 260 °F**
- **Sealing Dwell time** **3.0 Fixed**
- **Sealing Pressure** **50 psi Fixed**

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 20 of 27

11.0 Equipment Setup:

11.1 Reference Figure AA through EE for Emplex Model MPS6300 Continuous Band Sealer Controls.

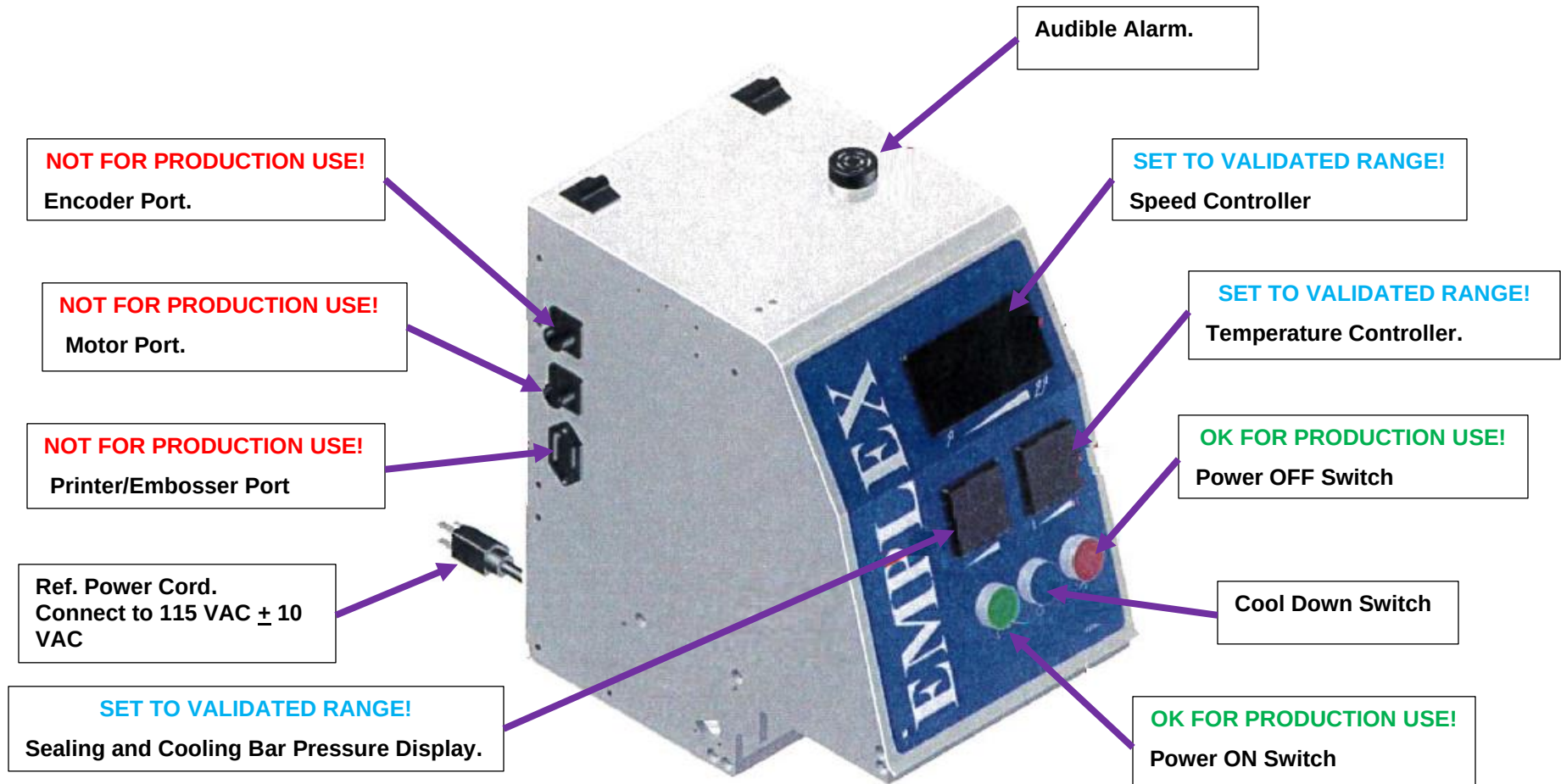


FIG. AA – Control Box Features

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 21 of 27

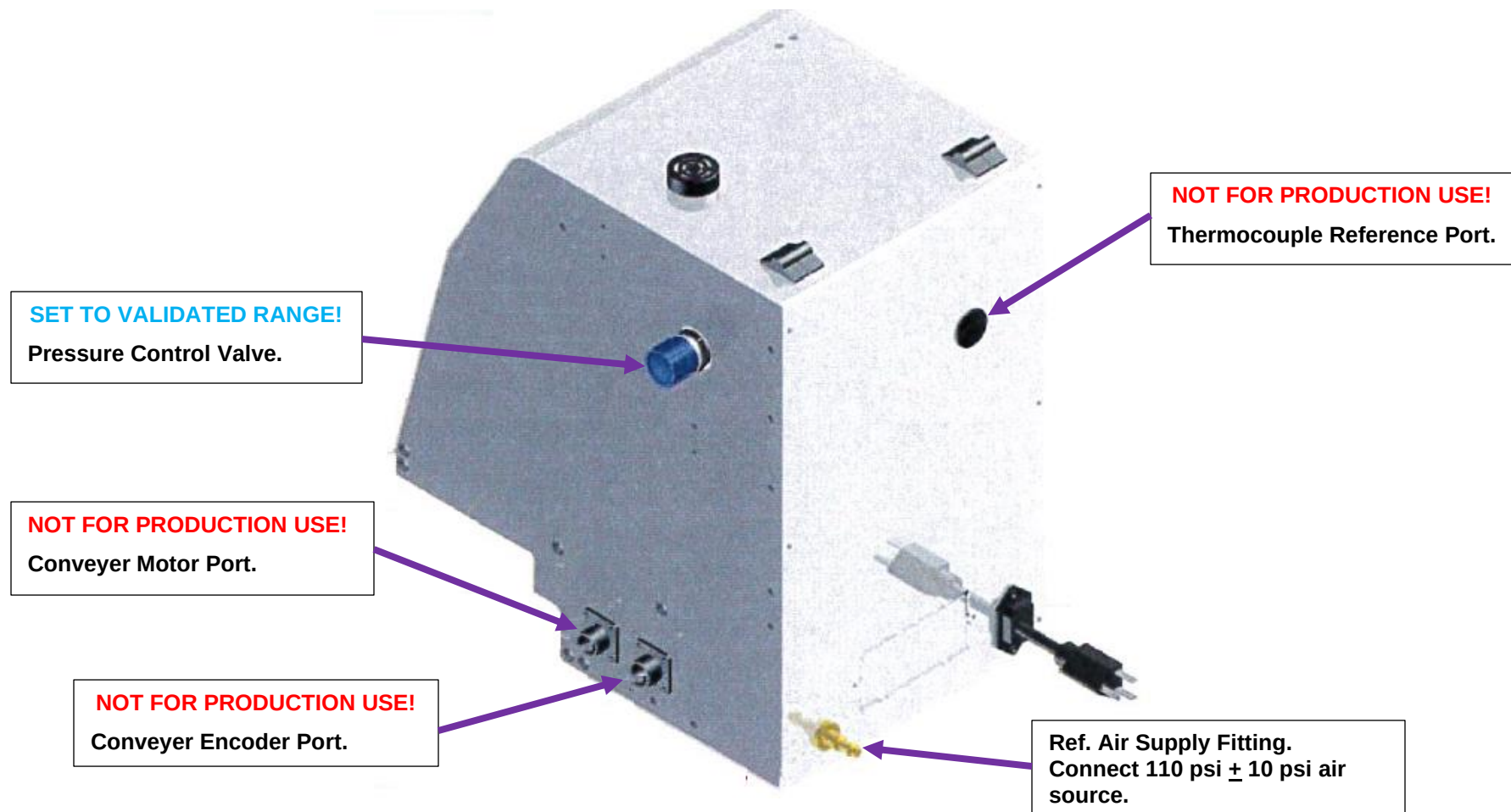


FIG. BB – Control Box Features continued

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 22 of 27

SETTING THE SEAL TEMPERATURE

1. Press "MODE" switch once.
2. Press UP ARROW (▲) or DOWN ARROW (▼) until Seal Temperature is within the Validated operating range.
3. Press "MODE" switch again to save.



FIG. CC – Seal Temperature Control

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 23 of 27

SETTING THE SEALSPEED

1. Press UP ARROW (▲) or DOWN ARROW (▼) until Seal Speed is within the Validated operating range.
2. Press "Enter" switch to save.

**FIG. DD – Set Seal Speed**

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 24 of 27

ADJUSTING SEALING AND COOLING BAR PRESSURE

1. Unlock the Pressure Control Valve by gently pulling the blue knob out from the control box.
2. Turn the Pressure Control Valve knob clockwise to increase pressure or counter-clockwise to decrease pressure (see Figure 12 for location of the Sealing and Cooling Bar pressure display) until Sealing and Cooling Bar Pressure is within the Validated operating range.
3. Once the desired setting is achieved, push the blue knob towards the control box to lock it in position.
4. The pressure display is identical to the temperature controller (ref. Figure 14). When the pressure regulator is adjusted the "Process Value" on the pressure display changes. If the "Process Value" deviates from the "Set Value" the Alarm will sound to alert the operator of a change in pressure.
5. Press "MODE" switch once.
6. Press UP ARROW (▲) or DOWN ARROW (▼) until "Set Value" matches the "Process Value" and is within the Validated operating range.
7. Press "MODE" switch again to save.

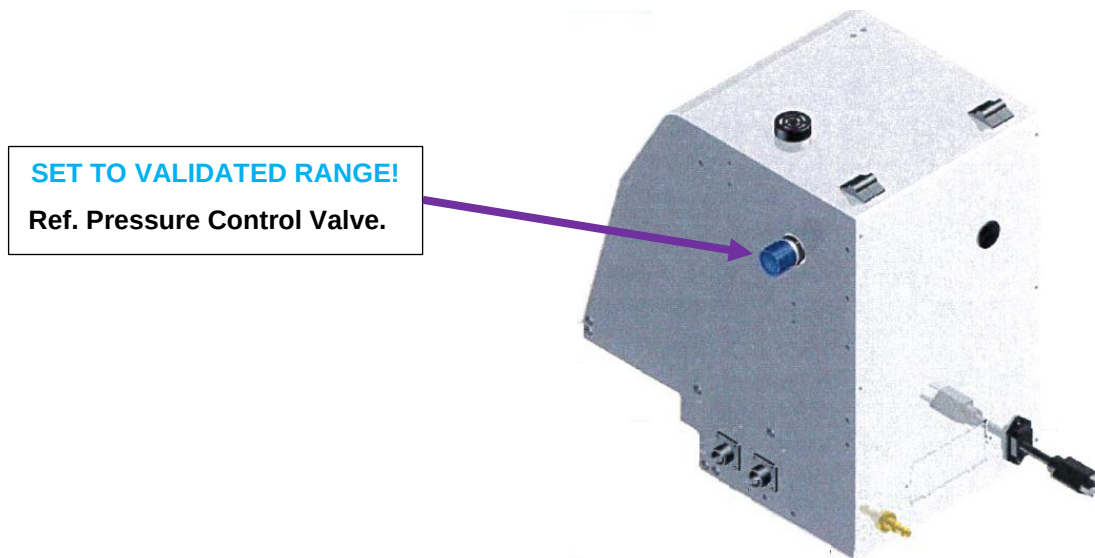


FIG. EE – Set Sealing and Cooling Bar Pressure

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: B
			Page 25 of 27

12.0 Appendix A: SPT Catheter Artwork



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	1

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex



SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160










 SLQ-MMDDYYYY
 XXXXXX
 20YY-MM-DD


 (01) #####
 (10) SLQ-MMDDYYYY
 (17) YYMMDD

SLQ-LB-4008 Rev. A

2-WAY 100% SILICONE CLEARTRACT SPT CATHETER	
SIZE	XX Fr, XX mL
QUANTITY	10













 (01) #####
 (10) SLQ-MMDDYYYY
 (17) YYMMDD

 SLQ-MMDDYYYY
 XXXXXX
 20YY-MM-DD

Symbols are defined in Instruction for Use (IFU)
 Store away from direct light exposure at room temperature
 Not made with natural rubber latex


 SILQ Technologies Corp.
 323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 26 of 27



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF XXXXXX

LOT SLQ-MMDDYYYY



20YY-MM-DD

QTY 100



SILQ Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters		Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 27 of 27

13.0 Change History

Revision	Change Description
A	Initial Release
B	Updated images to reflect new labels and shipper boxes. Added SteriPax pouch validated range to the sealer settings section.



INDIVIDUAL TRAINING RECORD

Doc No. FM-0004

Rev.	B
------	---

Page 1 of 1

[illegible]

This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled document. Users of this document must verify that they are using the most current revision of the controlled document before performing work.