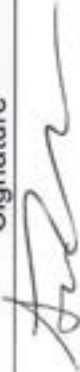


	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 1 of 7


Production Planning – Document the products identification information with the planned build quantity and any deviation(s) that may apply. Attach a copy of the device's specification/drawing/BOM and any applicable deviation(s) to the work order form prior to transfer to manufacturing.

Part No.	Rev	Lot No./ Serial No.	Expiration Date	Quantity	Deviation No.
SLQ-211810SPT	B	SLQ-11182024	2027-11-30	805/805	NA
Comments packs of 10					
SLQ 2-way Cleartract SPT 18 Fr x 10ml Catheter, Sterile, 10-Pack					
Name			Signature		Date
Sergio Tamayo					11/18/2024

Lot No/Serial Number	SLQ-11182024
----------------------	--------------

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in handwritten form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 2 of 7

Line Clearance - The individual performing the line clearance certifies that line clearance activities were conducted per WI-0014 prior to starting manufacturing activities. Quality management or designee verifies the line clearance was performed. *Note: The date line clearance occurred is considered the start date of manufacturing.*

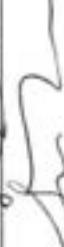Actions		Performed by	Date
1. Clear Area of previous WO Raw Materials.	☒ Completed.		11/18/2024
2. Verify kitted Raw Materials (RM) for: a. RM meet BOM revision level requirements. b. RM have not Expired. c. RM have no physical damage	☒ Completed. ☒ Completed. ☒ Completed.		11/18/2024
3. Equipment is calibrated and properly functioning. <i>Note: Equipment identification and calibration due date will be documented in the work order.</i>	☒ Completed. ☐ N/A		11/18/2024
4. Proper documentation is present and available.	☒ Completed.		11/18/2024
5. Production Associates are trained and qualified to perform the process.	☒ Completed. ☐ N/A		11/18/2024
6. All work surfaces cleaned with IPA.	☒ Completed. ☐ N/A		11/18/2024
7. Labels printed per Label Generation and Reconciliation Form.	☐ Completed. ① ☒ N/A		11/18/2024
8. Verify validated processes have been set to critical points within specification.	☒ Completed. ☐ N/A		11/18/2024
9. Verify Start-up activities are within specification	☒ Completed. ☐ N/A		11/18/2024
Conclusion: Line clearance performed	☒ Completed		11/18/2024
Quality Management			
Name		Signature	
Larry Dawin			
		Date	
		11/18/2024	

① Labels to be printed at a later date. See label reconciliation forms for print date. ST 11/18/2024

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission therefrom. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 3 of 7

Material/ Procedure Traceability (Check appropriate Box):

- ☐ – Identification information for the procedures and materials used correspond to the attached product's BOM and are captured on the attached Production Pick List (ref. P/N FM-0059).
- ☒ - Identification information for the procedures and materials used correspond to the attached product's BOM and are captured in the table below (N/A or delete table if not applicable).

BOM Item No.	Type	Rev	Material Designation	Lot No./Serial No.	Expiration Date	Rev	Material Designation	Lot No./Serial No.	Expiration Date
SLQ-8102	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-11182024	2027-11-30		☐ Primary ☐ Substitute		
1	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	20240401	N/A		☐ Primary ☐ Substitute		
2	☐ Document ☒ Material	B	☒ Primary ☐ Substitute	2302-01-2	02/23/2024, REF NCR 24-027		☐ Primary ☐ Substitute		
3	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
4	☐ Document ☒ Material	A	☐ Primary ☒ Substitute	A199054	N/A		☐ Primary ☐ Substitute		
5	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-11182024	2027-11-30		☐ Primary ☐ Substitute		
6	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	000072 011223	N/A		☐ Primary ☐ Substitute		
7	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-11182024	2027-11-30		☐ Primary ☐ Substitute		
8	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
9	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
10	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
11	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
12	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
13	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		

	Work Order Form		This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission therefrom. If any portion of this document is in handcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.	Doc No.: FM-0039
				Rev.: G
				Page 4 of 7

BOM Item No.	Type	Rev	Material Designation	Lot No./Serial No.	Expiration Date	Rev	Material Designation	Lot No./Serial No.	Expiration Date
14	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
15	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
16	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
SLQ-211810SPT	☐ Document ☒ Material	B	☒ Primary ☐ Substitute	SLQ-11182024	2027-11-30		☐ Primary ☐ Substitute		
1	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
2	☒ Document ☐ Material	E	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
3	☒ Document ☐ Material	G	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
4	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
5	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
6	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
7	☒ Document ☐ Material	F	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
8	☒ Document ☐ Material	C	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
9	☒ Document ☐ Material	F	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
10	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		

N/A ST 12/11/2024

Lot No./Serial Number	SLQ-11182024
-----------------------	--------------

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 5 of 7

Training Traceability – All individuals indicated in the following table have been trained to the referenced document(s) which support the manufacturing, inspection and test activities of the product identified in the work order. Leadership representative has certified that the individual(s) have received proper training prior to conducting manufacturing activities.

Document No.	Rev	Print Name	Signature
MP-SLQ-0001	B	James Romero	James Romero
		Francisco VARGAS	Francisco
		Mark Legaspi	Mark Legaspi
		Saverdy Pen	Saverdy Pen
MP-SLQ-0002	B	James Romero	James Romero
		Francisco VARGAS	Francisco
		Mark Legaspi	Mark Legaspi
		N/A	N/A 12/11/2024
Leadership Representative			
Name		Signature	Date
Sean Tamayo		Sean Tamayo	11/18/2024

~~Francisco V~~ ST 11/26/2024
 Note: Mark Legaspi added to the line on 11/26/2024 ST 11/26/2024
 Note: Saverdy Pen added to the line on 12/03/2024 ST 12/03/2024

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission therefrom. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 7 of 7

Documentation Checklist – If applicable, the following documentation shall be included upon completion of the manufacturing build. Document packet shall be organized in the order indicated below. Ensure that good documentation practices are applied per QP-003 "Quality Records". Check box(es) ☐ that apply to each requirement. The individual who signs has confirmed that all requirements have been met.

1. Work Order (WO) Record	☒ Present and Completed	
2. Product's Specification/Drawing & BOM (Arena PLM)	☒ Present	
3. Deviations	☐ Present	☒ N/A
4. Product Labeling & Generation/ Reconciliation Record	☐ Present, Completed and Results Meet Specification	☒ N/A
5. Non-Conforming Reports (NCR)	☒ Present and Completed	☐ N/A
6. Inspection Records	☒ Present, Completed and Results Meet Specification	
7. Critical Process Parameter Record	☒ Present, Completed and Results Meet Specification	☐ N/A
Name		Date
		01/07/2025

Release Authorization for Further Processing (e.g., terminal sterilization) – The individual signing below is certifying that they independently reviewed the Device History Record (DHR) and confirmed that all required documentation is included, records are complete, good documentation practices were followed, and specifications have been met (if not NCR is included)

Name	Signature	Date
Larry Daron		01/07/2025

Lot No/Serial Number	SLQ-11182024
----------------------	--------------

SLQ-211810SPT

In Production

SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack


Revision »B - ECO-001669 Status: Effective, Status Date: 05/14/2024 03:24:42 PM, Shared

Contains **10** first-level items, **50** line items, **41** unique items, **33** of which are shared.

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	△ t□ t□ t□ t□			
▶ 1 SLQ-8102 rev A	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack	Subassembly	In Prod	t□	3	1 EA	
2 WI-0006 rev E	Work Order Process	Work Instruction	In Prod	t□	2	1 each	
3 FM-0039 rev G	Form, Work Order	Form	In Prod	t□	2	1 each	
4 WI-0014 rev B	Line Clearance Process	Work Instruction	In Prod	t□	2	1 each	
5 MP-SLQ-0003 rev A	Manufacturing Process Flow, SILQ Foley Catheter	Manufacturing Procedure	In Prod	t□	2	1 EA	
6 SLQ-0001 rev A	Sterilization Load Shipping Pallet Configuration	Subassembly	In Prod	t□	2	1 EA	
▶ 7 QP-032 rev F	EO Sterilization Process Control	Quality System Procedure	In Prod	t□	7	1 each	
8 WI-0021 rev C	Work Instruction, Pyrogen (LAL) Evaluation Process	Work Instruction	In Prod	t□	2	1 each	
▶ 9 QP-010 rev F	Final Product Inspection	Quality System Procedure	In Prod	t□	6	1 each	
10 FM-0029 rev B	Form, Certificate of Conformance	Form	In Prod	t□	2	1 each	

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
 1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY SIQ PRIOR TO PRODUCTION.
 2. FINAL PRODUCT INSPECTION PER QP-010.
 3. PRODUCT EXPIRY DATE IS 36 MONTHS FROM THE DATE OF MANUFACTURE.
2. STERILIZATION
 1. PRODUCT IS STERILIZED VIA ETHYLENE OXIDE (ETO) STERILIZATION PER QP-002 CYCLE PWH-7374.
3. TESTING
 1. PYROGEN EVALUATION PER WH-0021 AT ≤ 0.5 EU/ML OR ≤ 20.0 EU/DEVICE

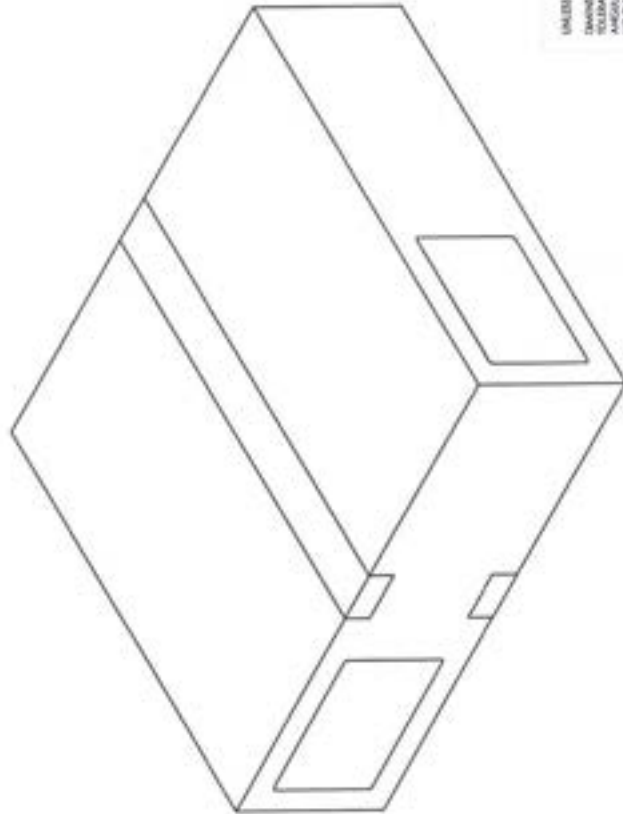

REVISIONS	
REV.	DESCRIPTION
A	INITIAL RELEASE
B	CHANGE PRODUCT EXPIRY DATE FROM 12 TO 36 MONTHS, REVISED FORMAT TO MATCH 1-018

pathway

TITLE: 2-WAY CLEARTRACT CATHETER, 18 FR X 10 ML, STERILE, 10-PACK
 SIZE: DWG. NO. B SLQ-211810SP
 REV B
 SHEET 1 OF 1

UNLESS OTHERWISE SPECIFIED:
 DIMENSIONS ARE IN MILLIMETERS
 DECIMALS ARE TO TWO PLACES
 ANGLES ARE TO NEAREST 5 DEGREES
 NO PLACE DECIMAL 1 0.5
 ONE PLACE DECIMAL 1 0.1
 TWO PLACE DECIMAL 1 0.05
 THREE PLACE DECIMAL 1 0.025

ASSEMBLY AND CONFORMANCE:
 THE INFORMATION CONTAINED HEREIN IS THE PROPERTY OF PATHWAY MEDICAL, INC. AND IS TO BE USED FOR THE PURPOSES OF THE INFORMATION CONTAINED HEREIN. NO PART OF THIS DOCUMENT IS TO BE REPRODUCED OR TRANSMITTED IN ANY FORM OR BY ANY MEANS, ELECTRONIC OR MECHANICAL, INCLUDING PHOTOCOPYING, RECORDING, OR BY ANY INFORMATION STORAGE AND RETRIEVAL SYSTEM, WITHOUT THE WRITTEN PERMISSION OF PATHWAY MEDICAL, INC.

Drawn By: B. SANDOVAL
 Drawn At: Scale



Printed by Sean Tamayo
Printed on 11/15/2024
Local time zone (GMT-08:00) Pacific Time (America/Los Angeles)
Page Arena > Workspace Items > Bill of Materials > Indented
All information contained in this document is proprietary and confidential.

SLQ-8102

In Production

SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack

Revision »A - ECO-001669 Status: Effective, Status Date: 05/14/2024 03:24:42 PM, Shared

Contains 16 first-level Items, 31 line Items, 24 unique Items, 23 of which are shared.

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE
1	SLQ-6002 rev A	Catheter, 18 Fr x 10 mL, Silicone, NY Foley, 2-Way	Fabricated Part	In Prod 3 10 EA
2	SLQ-3000 rev B	Solution, SILQ Polymer	Purchased Component	In Prod 4 80 ML
SUBSTITUTES				
	RANK#	ITEM NUMBER	ITEM NAME	NOTES
	1	HDX-3000 rev A	Solution, Hydrophilix Polymer	80 ML
3	SLQ-5000 rev A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	Packaging	In Prod 3 10 EA
4	SLQ-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	Packaging	In Prod 4 10 EA
SUBSTITUTES				
	RANK#	ITEM NUMBER	ITEM NAME	NOTES
	1	HDX-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	10 EA
5	SLQ-LB-4008 rev A	Label, Pouch, SILQ ClearTract SPT Catheter	Labeling	In Prod 2 10 EA

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE			
6	SLQ-5016 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS, Pre-Printed	Packaging	In Prod		3	1 EA
SUBSTITUTES							
	RANK\ITEM NUMBER	ITEM NAME		NOTES			
	1	SLQ-5015 rev A					
		Shelf Carton, 3.25 1 EA x 2.5 x 19.75, .018 SBS					
	2	SLQ-5014 rev A					
		Shelf Carton, 3.25 1 EA x 2.5 x 19.75, .024 SBS					
7	SLQ-LB-4009 rev A	Label, Box, SILQ ClearTract SPT Catheter	Labeling	In Prod		2	2 EA
8	SLQ-4007 rev A	Artwork, IFU, SILQ 2-Way ClearTract SPT Catheter	Labeling	In Prod		2	1 EA
9	SLQ-5012 rev A	Box, 20 in. L x 18 in. W x 6 in. D, 200 LB Test	Packaging	In Prod		2	0.1 EA
10	PWY-5000 rev A	Tape, Packaging, 2 Mil, 2 in x 110 yds, Clear	Packaging	In Prod		2	5.2 IN
11	SLQ-LB-4010 rev A	Label, Shipper Box, SILQ ClearTract SPT Catheter	Labeling	In Prod		2	0.2 EA
12	MP-SLQ-0002 rev B	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Manufacturing Procedure	In Prod		2	1 EA
13	FM-SLQ-0001 rev A	Form, Run Sheet for SILQ Automation Machine	Form	In Prod		2	1 EA
14	TM-HDX-0001 rev B	Test Method, UV Light Absorbance, Hydrophilix Catheters	Test Method	In Prod		2	1 each
15	TM-HDX-0002 rev A	Test Method, Rinse Test, Hydrophilix Catheter	Test Method	In Prod		2	1 each
16	MP-SLQ-0001 rev B	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	Manufacturing Procedure	In Prod		2	1 EA


NOTE: UNLESS OTHERWISE SPECIFIED

1. GENERAL
 1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY PATHWAY PRIOR TO PRODUCTION.
 2. DRAWING IS FOR REFERENCE ONLY. REFER TO THE RELEASED BILL OF MATERIALS FOR THE COMPLETE PARTS LIST.
 3. UP UNTIL SEALING THE POUCH, PRODUCT SHALL BE PACKAGED IN A CLEANROOM ENVIRONMENT. POST SEALING THE POUCH, ASSEMBLY WILL CONTINUE OUTSIDE OF THE CLEANROOM ENVIRONMENT.

2. ASSEMBLY PROCEDURE
 1. APPLY POLYMER SOLUTION TO THE CATHETER IN THE HATCHED AREA SHOWN. ALLOW THE CATHETER TO DRY OVERNIGHT.
 2. INSERT THE CATHETER INTO A POLYETHYLENE SLEEVE SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CLOSED END.
 3. PLACE THE POUCH LABEL ON THE CLEAR SIDE OF THE POUCH. WHEN VIEWING THE CLEAR SIDE OF THE POUCH FROM ABOVE, WITH THE CHEVRON TO THE RIGHT, THE LABEL TEXT SHALL READ LEFT TO RIGHT. WITH BEST EFFORT, APPLY THE LABEL SUCH THAT IT IS CENTERED.
 4. INSERT THE SLEEVE / CATHETER INTO THE POUCH SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CHEVRON END OF THE POUCH.
 5. SEAL THE POUCH.
 6. PLACE 10 SLEEVES, POUCHES, SEALED, AND LABELED CATHETERS INTO ONE SHELF CARTON. APPLY TWO SHIPPER BOX LABELS TO THE LONG, WIDER SIDE OF THE CARTON. VERIFY THAT THE TEXT READS LEFT TO RIGHT WHEN THE OPENING OF THE CARTON IS TO THE LEFT. WITH BEST EFFORT, CENTER THE LABEL ON ITS RESPECTIVE FACE OF THE CARTON.
 7. PLACE 10 SHELF CARTONS INTO ONE SHIPPER BOX. SEAL WITH PACKAGING TAPE AS NEEDED. PLACE A SHIPPER BOX LABEL IN THE UPPER LEFT CORNER ON TWO ADJACENT FACES, READING LEFT TO RIGHT.

3. DIMENSIONAL
 1. [] INDICATES LABEL REFERENCE DIMENSIONS.
 2. [] INDICATES CRITICAL DIMENSIONS.
4. CLEANING AND HANDLING
 1. ALL PRODUCT, UP UNTIL INSERTING THE POUCHES INTO SHELF CARTONS, SHALL BE HANDLED BY GLOVED HANDS APPROPRIATE FOR A CLEANROOM ENVIRONMENT.
5. FUNCTIONAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
 1. THE CATHETER'S ABSORBANCE VALUE AT 265 NANOMETERS SHALL BE GREATER THAN OR EQUAL TO 0.300 AU. REFER TO TM-HOK-0001.
 2. THE CATHETER'S ABSORBANCE VALUE AT 440 NANOMETERS SHALL BE LESS THAN OR EQUAL TO 0.020 AU. REFER TO TM-HOK-0002.
 3. THE PEEL STRENGTH FOR A 1 INCH WIDE SEALED SECTION OF THE STERILE BARRIER PACKAGING SHALL BE ≥ 1.00 LB WHEN TESTED PER TM-0003. IT IS NOT REQUIRED TO PLACE CATHETERS IN THESE TEST POUCHES, OR LABEL THEM. INSPECT N=3 SAMPLES, 1 AT THE BEGINNING, MIDDLE, AND END OF THE MANUFACTURING RUN.

6. VISUAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
 1. THE CATHETER SHALL BE FREE OF CURED POLYMER AND DEBRIS.
 2. THE CATHETER SHALL BE PROPERLY ORIENTED IN THE SLEEVE AND POUCH WITH THE LABEL APPLIED TO THE APPROPRIATE FACE AND IN THE CORRECT ORIENTATION.
 3. INSIDE THE STERILE BARRIER, DEVICE IS TO BE FREE OF PARTICULATES WHEN VIEWED AT 18 INCHES AWAY UNDER AMBIENT ROOM LIGHTING AND WITHOUT THE USE OF MAGNIFICATION.
 4. THE STERILE BARRIER PACKAGING SHALL NOT INCLUDE ANY OF THE FOLLOWING DEFECTS:
 1. UNSEALED AREAS, NONHOMOGENEOUS OR UNSEALED AREAS, OVERSEALED AREAS, ANHROW SEALS, CHANNELS, WRINKLES / FOLDOVERS / CRACKS, OR TEARS / PINHOLES.
 5. THE SHELF CARTON SHALL HAVE THE CORRECT NUMBER OF CATHETERS INSIDE AND FACING THE CORRECT ORIENTATION. THE LABELS SHALL BE AFFIXED IN THE CORRECT LOCATION AND ORIENTATION.
 6. THE SHIPPER BOX SHALL HAVE THE CORRECT NUMBER OF SHELF CARTONS. THE LABELS SHALL BE APPLIED IN THE CORRECT LOCATION AND ORIENTATION.

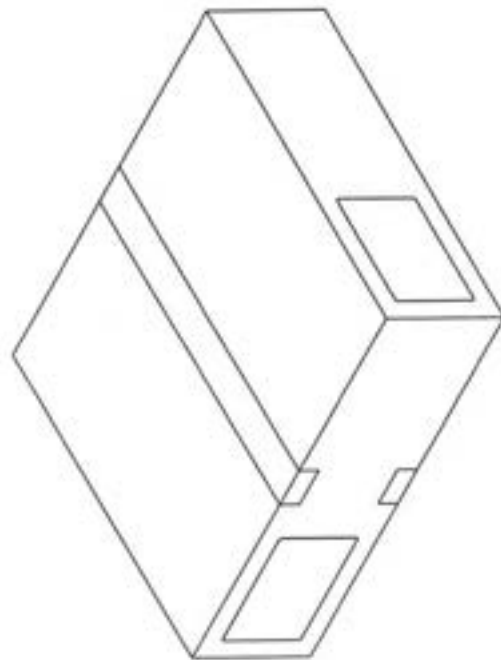

REVISIONS		
ZONE	REV.	DATE
	A	
INITIAL RELEASE		

pathway
New Product Introduction

TITLE: CATHETER, SILICONE, NY POLEY, NON-STERILE, 10-PACK, MASTER DRAWING
 SIZE: DWG. NO. B SLQ-81XX
 REV. A
 SCALE: 1:5
 SHEET 1 OF 3

UNLESS OTHERWISE SPECIFIED:
 DIMENSIONS ARE IN INCHES
 DIMENSIONS IN PARENTHESES ARE IN MILLIMETERS
 HOLE PLACEMENT SHALL BE TO DIMENSION
 ONE PLACE DECIMAL 15.3
 TWO PLACE DECIMAL 15.32
 THREE PLACE DECIMAL 15.325

PROPRIETARY AND CONFIDENTIAL
 THE INFORMATION CONTAINED HEREIN IS THE PROPERTY OF PATHWAY MEDICAL, INC. AND IS NOT TO BE REPRODUCED OR TRANSMITTED IN ANY FORM OR BY ANY MEANS, ELECTRONIC OR MECHANICAL, INCLUDING PHOTOCOPYING, RECORDING, OR BY ANY INFORMATION STORAGE AND RETRIEVAL SYSTEM, WITHOUT THE WRITTEN PERMISSION OF PATHWAY MEDICAL, INC.


ISO
9001
CERTIFIED

FDA
REGISTERED



DO NOT SCALE DRAWINGS

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
1	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	10
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	10
3	SLQ-60XX	SEE TABLE FOR PART SPECIFICATION	10
4	SLQ-5016	SHELF CARTON, 3.25 X 2.5 X 19.75, .018 SBS, PRE-PRINTED	1
5	SLQ-LB-4008	POUCH LABEL, CLEARTRACT SP CATHETER	10
6	SLQ-LB-4009	SHELF CARTON LABEL, CLEARTRACT SP CATHETER	1
7	SLQ-4007	IFU, SILQ 2-WAY CLEARTRACT SP CATHETER	1

PART NUMBER	DESCRIPTION	SLQ-81XX
SLQ-6000	CATHETER, 14 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8100
SLQ-6001	CATHETER, 16 FR X 5 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8101
SLQ-6002	CATHETER, 18 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8102
SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8104

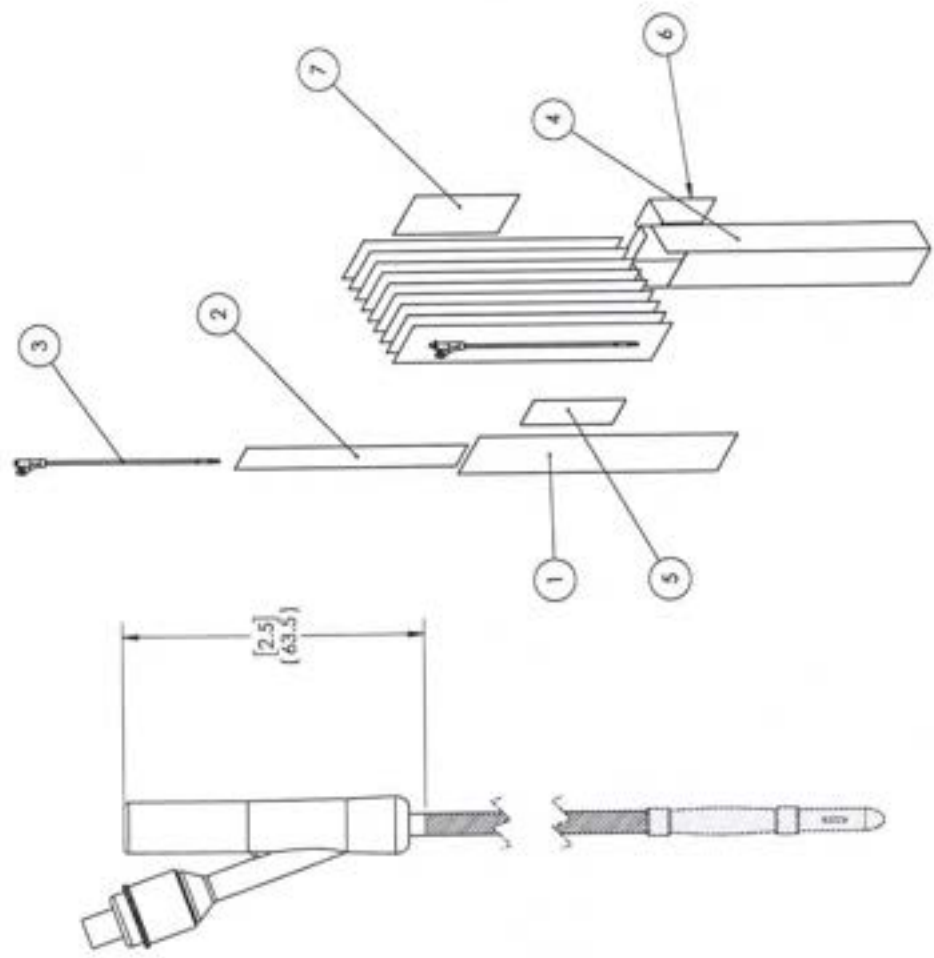

ITEM NO.	PARTNO	DESCRIPTION	QTY.
8	SLQ-5012	BOX, 20 IN. L X 18 IN. W X 6 IN. D, 200 LB TEST	1
9	PWY-5000	TAPE, PACKAGING, 2 MIL. 2 IN X 110 YDS, CLEAR	AR
10	SLQ-LB-4010	SHIPPER BOX LABEL, CLEARTRACT SP CATHETER	2

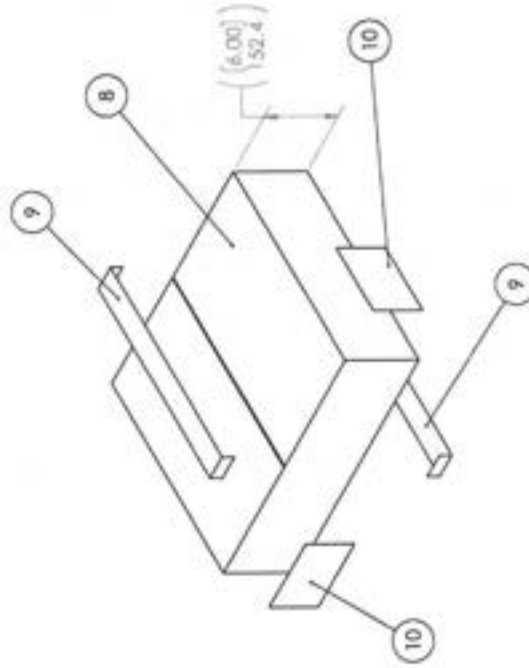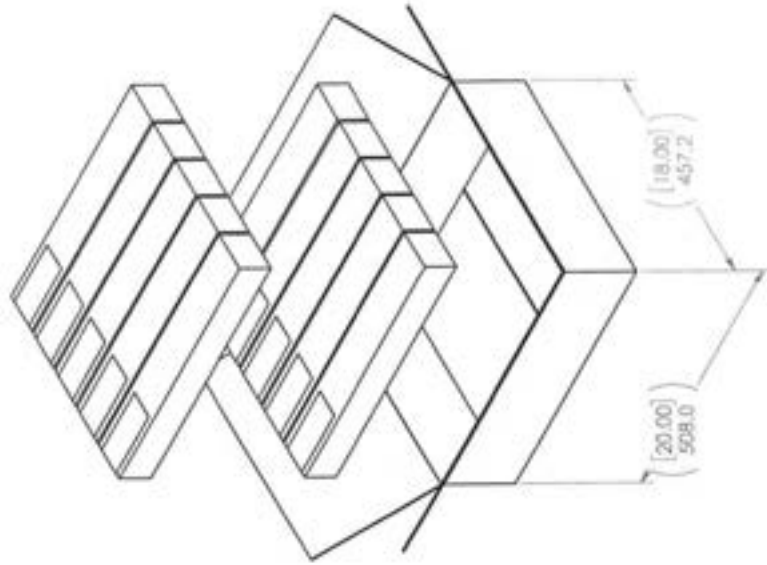

FDA
REGISTERED

ISO
CERTIFIED 13485

EDGE TEMPL NO. **8** **SLQ-81XX**

REV **A**

SCALE: 1:1 SHEET 3 OF 3

	Label Generation & Reconciliation				Doc No.: FM-0033	
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.					
	Rev.: D Page 1 of 2					

Label Information						
Part No.	SLQ-LB-4009	Rev.	A	Lot No.	SLQ-11182024 ☐ N/A	Label QTY
						10 ☐ N/A
Description	Label, Box, SILQ ClearTract SPT Catheter	UDI No.	10850041415168 ☐ N/A		Expiration Date:	2027-11-18 ☐ N/A 2027-11-30 ST 11/27/2024

Label Quantity Calculation	
Required Label Quantity (A)	81
2% Overage ($A \times 0.02 = B$)	3
QC Retain Labels	2
Quantity of Labels Required ($A + B + 2 = C$)	86
Quantity of Labels Printed ($C = D$)	86
Production's Initial	ST Date 11/27/2024

Label Verification	
What is the revision of the label per the BOM?	A
What is the revision on the printed label?	A
Is the printed information correct per specifications?	☒ Yes ☐ No
First and last label printed is affixed to back of form	☒ Yes ☐ No
Quantity of labels released to production (E)	84
QC's Initials	LS Date 12/3/2024



Label Reconciliation	
Quantity of labels used by production (F)	81
Quantity scrapped by production (G)	3
Quantity of unused labels returned to QC and scrapped (H)	0
($E - F - G - H = 0$)	0
Are all labels accounted for?	☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

1. Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
2. Review the appropriate workstations and surrounding areas for the missing labels.



Investigational Results
☒ N/A

	Label Generation & Reconciliation		Doc No.: FM-0033
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: D
			Page 2 of 2
QC's Initials			Date 12/11/2024

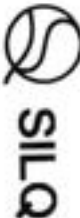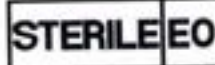QUANTITY 10	SIZE 18 Fr, 10 mL	2-WAY 100% SILICONE CLEARTRACT SPT CATHETER 	       
			 (01) 10850041415168 (10) SLQ-11182024 (17) 271130 Symbols are defined in Instruction for Use (IFU) Store away from direct light exposure at room temperature Not made with natural rubber latex  SILQ Technologies Corp. 323 Sunny Isles Blvd, 7th Floor Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

First

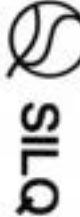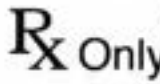QUANTITY 10	SIZE 18 Fr, 10 mL	2-WAY 100% SILICONE CLEARTRACT SPT CATHETER 	       
			 (01) 10850041415168 (10) SLQ-11182024 (17) 271130 Symbols are defined in Instruction for Use (IFU) Store away from direct light exposure at room temperature Not made with natural rubber latex  SILQ Technologies Corp. 323 Sunny Isles Blvd, 7th Floor Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

Last

	Label Generation & Reconciliation				Doc No.: FM-0033	
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.					
	Rev.: D					
						Page 1 of 2

Label Information						
Part No.	SLQ-LB-4008	Rev.	A	Lot No.	SLQ-11182024 ☐ N/A	Label QTY
						1 ☐ N/A
Description	Label, Pouch, SILQ ClearTract SPT Catheter	UDI No.	00850041415161 ☐ N/A		Expiration Date:	2022-11-16 ☐ N/A 2027-11-30 51 11/27/2024

Label Quantity Calculation	
Required Label Quantity (A)	805
2% Overage (A×0.02=B)	16
QC Retain Labels	2
Quantity of Labels Required (A+B+2=C)	823
Quantity of Labels Printed (C=D)	823
Production's Initial	ST
Date	11/27/2024

Label Verification	
What is the revision of the label per the BOM?	A
What is the revision on the printed label?	A
Is the printed information correct per specifications?	☒ Yes ☐ No
First and last label printed is affixed to back of form	☒ Yes ☐ No
Quantity of labels released to production (E)	821
QC's Initials	LO
Date	11/3/2024

Label Reconciliation	
Quantity of labels used by production (F)	805
Quantity scrapped by production (G)	0
Quantity of unused labels returned to QC and scrapped (H)	16
(E-F-G-H=0)	0
Are all labels accounted for?	☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

1. Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
2. Review the appropriate workstations and surrounding areas for the missing labels.

Investigational Results
☒ N/A

NOTES: UNLESS OTHERWISE SPECIFIED

- GENERAL
 - PRIOR TO PRINTING, VERIFY THE INFORMATION ON THE LABEL IS ACCURATE WHEN COMPARED AGAINST THE WORK ORDER.
 - PRIOR TO PRINTING, VERIFY THE LABEL HAS GONE THROUGH THE INSTALLATION QUALIFICATION (IQ) PROCESS FOR LABELS. REFER TO FM-0044.
 - LABEL CONTENT SHALL BE DEVELOPED PER W6-0009.

- MATERIAL
 - LABEL: STANDARD WHITE THERMAL TRANSFER
 - RIBBON: WAX / RESIN

- DIMENSIONS
 - (I) INDICATES LABEL REFERENCE DIMENSIONS.

- SAMPLING PLAN
 - IN ADDITION TO THE LABEL GENERATION AND RECONCILIATION FORM, FM-0033, INSPECT N+3 ADDITIONAL SAMPLES AT THE BEGINNING OF THE RUN.


- GS1 DATA MATRIX VALUES
 - (01): GLOBAL TRADE ITEM NUMBER (GTIN)
 - (10): BATCH OR LOT NUMBER
 - (17): EXPIRATION DATE

VISUAL INSPECTION - LABEL TEXT

- POST-PRINTING, ALL TEXT ON THE LABEL SHALL BE LEGIBLE.
- ALL VARIABLE LABEL CONTENT (LOT NUMBER, EXPIRATION DATE, ETC.) SHALL MATCH THE WORK ORDER.
- THE GS1 DATA MATRIX SHALL BE FREE OF MISSING PIXELS AND SHALL NOT HAVE SNAJGUES OR STRIAKS THROUGH IT. WHEN SCANNED WITH A BARCODE SCANNER, THE TEXT ADJACENT TO THE GS1 DATA MATRIX SHALL MATCH THE TEXT DISPLAYED BY THE SCANNER.
- (01): VERIFY THE GTIN IS DETAILED IN AN HUMAN READABLE FORMAT ADJACENT TO THE GS1 DATA MATRIX. THE GTIN SHALL BE A 14 DIGIT NUMBER.
- (10): VERIFY THE BATCH OR LOT NUMBER FORMAT IS IN COMPLIANCE WITH QP-019.
- (17): VERIFY THE EXPIRATION DATE IS IN THE FOLLOWING FORMAT: YYMMDD.
- VERIFY THE CATHETER STYLE AND CATHETER SIZE ARE REPRESENTATIVE OF A CONFIGURATION FOUND IN THE SKU TABLE.
- VERIFY THE GTIN IS FOUND IN THE SKU TABLE AND CORRESPONDS WITH THE CORRECT CATHETER SIZE AND STYLE.
- VERIFY THE REFERENCE NUMBER FOR THE CATHETER CORRESPONDS TO THAT FOUND IN THE SKU TABLE.
- VERIFY "QUANTITY" IS 10.

VISUAL INSPECTION - LABEL SYMBOLS

- POST-PRINTING, ALL SYMBOLS ON THE LABEL SHALL BE FREE OF SNAJGUES / MISSING ARTWORK.
- THE EXPIRATION DATE NEXT TO ITS CORRESPONDING SYMBOL SHALL BE IN THE FORMAT: YYYY-MM-DD.

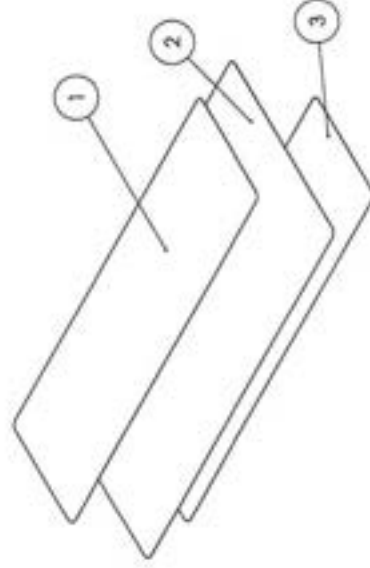

ITEM NO.	PARTNO	DESCRIPTION	QTY.
1	SLQ-AW-4009	LABEL BOX, SILQ CLEARTRACT SPT CATHETER	1
2	PWY-5007	RIBBON, THERMAL TRANSFER, WAX/RESIN, 4.3W	1
3	SLQ-5004	LABEL, 2.3125 IN. X 7.50 IN.	1

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

pathway
New Product Installation

TITLE:

LABEL, BOX, SILQ CLEARTRACT SPT CATHETER

SIZE DWG. NO.

B SLQ-LB-4009

SCALE: 1:1

SHEET 1 OF 2



REGISTERED CERTIFIED

DO NOT SCALE DRAWINGS

2

4

5

6

7

8

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
1. PRIOR TO PRINTING, VERIFY THE INFORMATION ON THE LABEL IS ACCURATE WHEN COMPARED AGAINST THE WORK ORDER.
2. PRIOR TO PRINTING, VERIFY THE LABEL HAS GONE THROUGH THE INSTALLATION QUALIFICATION (IQ) PROCESS FOR LABELS. REFER TO FM-0044.
3. LABEL CONTENT SHALL BE DEVELOPED PER WI-0009.

2. MATERIAL
 1. LABEL: STANDARD WHITE THERMAL TRANSFER
 2. RIBBON: WAX / RESIN

3. DIMENSIONS
1. ☐ INDICATES LABEL REFERENCE DIMENSIONS.

4. **SAMPLING PLAN**
1. IN ADDITION TO THE LABEL GENERATION AND RECONCILIATION FORM, FM-0033, INSPECT N-3 ADDITIONAL SAMPLES AT THE BEGINNING OF THE RUN.


5. GS1 DATA MATRIX VALUES
 1. (01): GLOBAL TRADE ITEM NUMBER (GTIN)
 2. (10): BATCH OR LOT NUMBER
 3. (17): EXPIRATION DATE

6. VISUAL INSPECTION - LABEL TEXT

1. POST-PRINTING, ALL TEXT ON THE LABEL SHALL BE LEGIBLE.
2. ALL VARIABLE LABEL CONTENT (LOT NUMBER, EXPIRATION DATE, ETC.) SHALL MATCH THE WORK ORDER.
3. THE GS1 DATA MATRIX SHALL BE FREE OF MISSING PIXELS AND SHALL NOT HAVE SMUDGES OR STRIAKS THROUGH IT. WHEN SCANNED WITH A BARCODE SCANNER, THE TEXT ADJACENT TO THE GS1 DATA MATRIX SHALL MATCH THE TEXT DISPLAYED BY THE SCANNER.
4. (01): VERIFY THE GTIN IS DETAILED IN AHUMAN READABLE FORMAT ADJACENT TO THE GS1 DATA MATRIX. THE GTIN SHALL BE A 14 DIGIT NUMBER.
5. (10): VERIFY THE BATCH OR LOT NUMBER FORMAT IS IN COMPLIANCE WITH QIP-019.
6. (17): VERIFY THE EXPIRATION DATE IS IN THE FOLLOWING FORMAT: YYMMDD. VERIFY THE CATHETER STYLE AND CATHETER SIZE ARE REPRESENTATIVE OF A CONFIGURATION FOUND IN THE SKU TABLE.
7. VERIFY THE GTIN IS FOUND IN THE SKU TABLE AND CORRESPONDS WITH THE CORRECT CATHETER SIZE AND STYLE.
8. VERIFY THE REFERENCE NUMBER FOR THE CATHETER CORRESPONDS TO THAT FOUND IN THE SKU TABLE.
9. VERIFY "QUANTITY" IS 1.

7. VISUAL INSPECTION - LABEL SYMBOLS
1. POST-PRINTING, ALL SYMBOLS ON THE LABEL SHALL BE FREE OF SMUDGES / MISSING ARTWORK.
2. THE EXPIRATION DATE NEXT TO ITS CORRESPONDING SYMBOL SHALL BE IN THE FORMAT: YYYY-MM-DD.

ITEM NO.	PARTNO	DESCRIPTION	QTY.
1	SLQ-AW-4008	LABEL, POUCH, SILQ CLEARTRACT SPT CATHETER	1
2	PWY-5007	RIBBON, THERMAL TRANSFER, WAX/RESIN, 4.3W	1
3	SLQ-5004	LABEL, 2.3125 IN. X 7.50 IN.	1

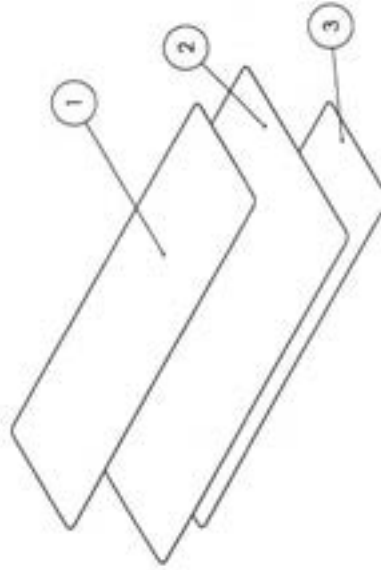

REVISIONS			DATE
DONE	REV.	DESCRIPTION	
	A	INITIAL RELEASE	

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGLES: $\pm 0.5^\circ$ MAX
HOLE PLACES: DECIMAL ± 0.13
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ± 0.13
HOLE PLACES: FRACTIONS ± 0.06

RESEARCH AND COMMERCE



TABLE 1. Continued

LABEL, POUCH, SILQ CLEARTRACT
SPT CATHETER

DATE	NAME: NO	REV
01-1-2000	01	01

B 510-18-4008

SCALE 1:1

THE UNIVERSITY OF CHICAGO PRESS


ISO 9001

REGISTERED

DEVICE WARNINGS / CONTRAINDICATIONS

CATHETER STYLE

CATHETER SIZE

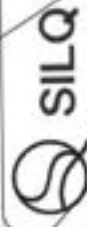

2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	1

Rx Only

Symbols are defined in Instruction for Use (IFU)
 Store away from direct light exposure at room temperature
 Not made with natural rubber latex

SILQ Technologies Corp.
 323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160



STERILE EO



LOT SLQ-MMDDYYYY

REF XXXXXX



(01) #####
 (10) SLQ-MMDDYYYY
 (17) YYMMDD



20"Y-MM-DD


SLQ-LB-4008 Rev. A

TRACEABILITY INFORMATION

COMPANY INFORMATION

CATHETER STYLE	CATHETER SIZE	GTIN	REF
2-WAY	16 Fr, 10 mL	00850041415154	211610SPT
2-WAY	14 Fr, 10 mL	00850041415130	211410SPT
2-WAY	16 Fr, 5 mL	00850041415147	211605SPT
2-WAY	18 Fr, 10 mL	00850041415161	211810SPT

	Label Generation & Reconciliation	Doc No.: FM-0033	
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: D
			Page 2 of 2
QC's Initials	<i>L</i>	Date	<i>12/11/2024</i>

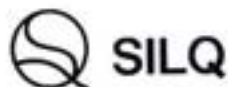

2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	18 Fr, 10 mL
------	--------------

QUANTITY	1
----------	---

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

 SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

R_x Only



STERILE EO



MR

LOT SLQ-11182024

REF 211810SPT

 2027-11-30



(01) 00850041415161
(10) SLQ-11182024
(17) 271130

SLQ-LB-4008 Rev. A

First



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	18 Fr, 10 mL
------	--------------

QUANTITY	1
----------	---

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

 SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

R_x Only



STERILE EO



MR

LOT SLQ-11182024

REF 211810SPT

 2027-11-30



(01) 00850041415161
(10) SLQ-11182024
(17) 271130

SLQ-LB-4008 Rev. A

Last

	Label Generation & Reconciliation				Doc No.: FM-0033	
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.					
	Rev.: D Page 1 of 2					

Label Information						
Part No.	SLQ-LB-4010	Rev.	A	Lot No.	SLQ-11182024 ☐ N/A	Label QTY
						100 ☐ N/A
Description	Label, Shipper Box, SLQ ClearTract SPT Catheter	UDI No.	10550041415166 ☒ N/A ST 12/10/2024		Expiration Date:	2027-11-18 ☐ N/A 2027-11-30 ST 12/10/2024

Label Quantity Calculation	
Required Label Quantity (A)	9
2% Overage (A×0.02=B)	1
QC Retain Labels	2
Quantity of Labels Required (A+B+2=C)	12
Quantity of Labels Printed (C=D)	12
Production's Initial	ST
Date	12/10/2024

Label Verification	
What is the revision of the label per the BOM?	A
What is the revision on the printed label?	A
Is the printed information correct per specifications?	☒ Yes ☐ No
First and last label printed is affixed to back of form	☒ Yes ☐ No
Quantity of labels released to production (E)	10
QC's Initials	Lp
Date	12/10/2024

Label Reconciliation	
Quantity of labels used by production (F)	9
Quantity scrapped by production (G)	0
Quantity of unused labels returned to QC and scrapped (H)	1
(E-F-G-H=0)	0
Are all labels accounted for?	☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

- Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
- Review the appropriate workstations and surrounding areas for the missing labels.

Investigational Results
☒ N/A

	Label Generation & Reconciliation	Doc No.: FM-0033
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>	Rev.: D
		Page 2 of 2
QC's Initials	<i>Lo</i>	Date
		<i>12/11/2024</i>



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF 211810SPT

LOT SLQ-11182024

 2027-11-30

QTY 100



SILO Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

First



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF 211810SPT

LOT SLQ-11182024

 2027-11-30

QTY 100



SILO Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

Last

Source:	☐ Receiving ☒ Start-Up Verification (QC) ☐ In-Process (QC) ☐ In-Process (Production)				Date:	05/14/2024
	☐ Final Product ☐ RMA (No.) ☐ Miscellaneous:					
Part No.:	SLQ-3000	Revision:	B	Lot No. / SN:	2302-01-2, 2302-01-1	
Description	Solution, SILQ Polymer	Supplier	SILQ Industries	Inspector	Larry Dacoron	

NONCONFORMITY(S) DETAIL(S)

#	Description	QTY Rejected	QTY Inspected	AQL	Specification	
					Procedure	Revision
1	Expired solution (02/22/2024 and 02/23/2024)	2	2	N/A	SLQ-3000	B

CONTAINMENT ACTIVITY

Contained QTY	Comments
2 barrels	N/A



INITIAL DISPOSITION

#	Disposition
1	☐ Increase Sampling ☐ Sort (See instructions below) ☐ Rework (See instructions below) ☒ Hold/Quarantine ☐ N/A

Note: Applies to Rework Dispositions Only	Risk Analysis Required?:	Risk Analysis:
	☐ Yes ☒ No	N/A

Additional Instructions:	N/A
--------------------------	-----

APPROVALS

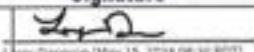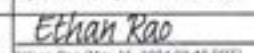Signature	Date	Signature	Date
☒ QA  Larry Dacoron (May 15, 2024 08:30 PDT)	05/15/24	☐ MFG	
☒ CUSTOMER  Ethan Rao (May 15, 2024 08:48 PDT)	05/15/24	☐	

SORT / REWORK COMPLETED BY

Signature	Date	Signature	Date
N/A	N/A	N/A	N/A

RE-INSPECTION

Description	QTY Rejected	QTY Inspected	AQL	Specification	
				Procedure	Revision
N/A	N/A	N/A	N/A	N/A	N/A

Inspected By:	N/A	Date	N/A
Signature			

FINAL DISPOSITION



#	Disposition	QTY
1	Meets Specification after Initial Disposition	0
	Use As Is (UAI)	2
	Customer Disposition Required?: ☒ Yes ☐ No	
	Return to Vendor (RTV)	0
	Scrap	0
Total:		2

Note: Applies to UAI Dispositions Only	Risk Analysis Required?:	Risk Analysis:
	☒ Yes ☐ No	Expired solutions may undergo chemical changes over time, leading to potential degradation of the product. Using expired solutions may yield inaccurate or unreliable data, undermining the validity of the study results.

Additional Instructions:	The suspension will be used in the manufacturing process of the product. Any product that will use the expired solution will be placed in quarantine after sterilization. After manufacturing of the product, the expired suspension will be subject to client's full chemical analysis QC process analogous to what was done upon initial validation of the suspension. An additional 50 16Fr catheters were manufactured using the same expired suspension, all of which will be tested after sterilization to further validate the suspension's acceptability. These 18Fr catheters will be subject to the standard sampling plan & QC process. If the results of the chemical analysis of the expired suspension, the testing of the 50 16Fr catheters, and the standard testing of this product are all acceptable, the product will be released from quarantine.
--------------------------	--

CAPA No. ☒ N/A	N/A
--	-----

APPROVALS

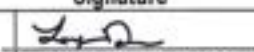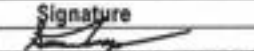Signature	Date	Signature	Date
☒ QA  Larry Dacoron (May 15, 2024 08:30 PDT)	05/15/24	☒ MFG  Sean Tamayo (May 15, 2024 08:30 PDT)	05/15/24

Pathway Medtech, LLC.		Nonconforming Record		NCR No.:	24-027	
☒ CUSTOMER	<u>Ethan Rao</u> Ethan Rao (May 15, 2024 08:48 PDT)	05/15/24	☐ N/A	N/A		N/A
Customer Comments / Special Disposition Instructions						
MATERIAL DESTRUCTION (APPLIES TO SCRAP ONLY)						☒ N/A
Performed by / Signature		Date		QTY Destroyed		
N/A		N/A		N/A		

NCR 24-027 - Expired solution

Final Audit Report

2024-05-15

Created:	2024-05-15
By:	Larry Dacoron (ldacoron@pathwaynpi.com)
Status:	Signed
Transaction ID:	CBJCHBCAABAAJBWGX-TJUPOJiKfnYn3DZGHlz0g4N4

"NCR 24-027 - Expired solution" History

-  Document created by Larry Dacoron (ldacoron@pathwaynpi.com)
2024-05-15 - 3:29:41 PM GMT
-  Document emailed to Ethan (ethanr@silq.tech) for signature
2024-05-15 - 3:29:45 PM GMT
-  Document emailed to Sean (stamayo@pathwaynpi.com) for signature
2024-05-15 - 3:29:45 PM GMT
-  Document emailed to Larry Dacoron (ldacoron@pathwaynpi.com) for signature
2024-05-15 - 3:29:45 PM GMT
-  Email viewed by Larry Dacoron (ldacoron@pathwaynpi.com)
2024-05-15 - 3:30:12 PM GMT
-  Document e-signed by Larry Dacoron (ldacoron@pathwaynpi.com)
Signature Date: 2024-05-15 - 3:30:28 PM GMT - Time Source: server
-  Email viewed by Sean (stamayo@pathwaynpi.com)
2024-05-15 - 3:35:08 PM GMT
-  Signer Sean (stamayo@pathwaynpi.com) entered name at signing as Sean Tamayo
2024-05-15 - 3:35:40 PM GMT
-  Document e-signed by Sean Tamayo (stamayo@pathwaynpi.com)
Signature Date: 2024-05-15 - 3:35:42 PM GMT - Time Source: server
-  Email viewed by Ethan (ethanr@silq.tech)
2024-05-15 - 3:47:49 PM GMT
-  Signer Ethan (ethanr@silq.tech) entered name at signing as Ethan Rao
2024-05-15 - 3:48:19 PM GMT



 Document e-signed by Ethan Rao (ethanr@silq.tech)
Signature Date: 2024-05-15 - 3:48:21 PM GMT - Time Source: server

 Agreement completed,
2024-05-15 - 3:48:21 PM GMT



pathway

Date Inspected		Part No.		Rev.		Description			
11/01/2024 - 11/11/2024		SLO-8102		A		SILQ 2-Way Cleartrac SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack			
Lot No.		Lot QTY		Rinsed catheter absorbance @440nm		Approved Supplier(s)			
SLO-11182024		60360.3 packs of 10		Ref TM-HDX-0001		Pathway			
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLO-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ±0.300 AU	Ref TM-HDX-0002	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text. 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text. 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box.
	Revision	A	A	A	A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	0.300	N/A	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	0.020	N/A	N/A	N/A	N/A	N/A
	UOM	EA	AU	AU	EA	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Spectrophotometer	Spectrophotometer	Visual	Visual	Visual	Visual	Visual
	Asset No	N/A	See Comments	See Comments	N/A	N/A	N/A	N/A	N/A
	Cal Due Date	N/A	See Comments	See Comments	N/A	N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size	11	11	11	11	11	11	11	7
Variable/ Attribute Data									
Sample No.		P	P	N/A	N/A	P	P	P	P
1		P	P	N/A	N/A	P	P	P	P
2		P	P	N/A	N/A	P	P	P	P
3		P	P	N/A	N/A	P	P	P	P
4		P	P	N/A	N/A	P	P	P	P
5		P	P	N/A	N/A	P	P	P	P
6		P	P	N/A	N/A	P	P	P	P
7		P	P	N/A	N/A	P	P	P	P
8		P	P	N/A	N/A	P	P	P	N/A
9		P	P	N/A	N/A	P	P	P	N/A
10		P	P	N/A	N/A	P	P	P	N/A
11		P	P	N/A	N/A	P	P	P	N/A



Result Summary												
Pass QTY	11	11	11	11	11	11	11	11	11	11	11	7
Fail QTY	0	0	0	0	0	0	0	0	0	0	0	0
Result (PIF)	P	P	P	P	P	P	P	P	P	P	P	P
Statistical Summary												
Min	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
StdDev	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Comments												
Rinse test and UV Absorbance Test were performed at SLQ. The rest of the visual inspections were performed at Pathway.												
Disposition:	ACCEPT					NCMR No:			N/A			
Inspected By:	<i>Lyndee S. Smith</i>					Date:			12/11/2024			
Verified By:	<i>Joe H.</i>					Date:			12/11/2024			



Date Inspected		Part No.		Rev.	Description				
12/10/2024		SLO-8102		A	SLO 2-Way Clearact SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack				
Lot No.		Lot QTY		Approved Supplier(s)					
SLO-11182024		603		Pathway					
In-Process Inspection	Catheter has been coated with SLO's polymer solution (SLO-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥ 0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤ 0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text. 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text. 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	A	B	A	A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	0.300	N/A	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	0.020	N/A	N/A	N/A	N/A	N/A
	UOM	EA	AU	AU	EA	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Spectrophotometer	Visual	Visual	Visual	Visual	Visual	Visual
	Asset No.	N/A	See Comments	N/A	See Comments	N/A	N/A	N/A	N/A
	Cal Due Date	N/A	See Comments	N/A	See Comments	N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size	N/A	11	N/A	N/A	N/A	N/A	N/A	N/A
Variable/ Attribute Data									
Sample No.									
1	N/A	N/A	0.412	0.000	N/A	N/A	N/A	N/A	N/A
2	N/A	N/A	0.436	0.003	N/A	N/A	N/A	N/A	N/A
3	N/A	N/A	0.461	0.006	N/A	N/A	N/A	N/A	N/A
4	N/A	N/A	0.471	0.002	N/A	N/A	N/A	N/A	N/A
5	N/A	N/A	0.464	0.002	N/A	N/A	N/A	N/A	N/A
6	N/A	N/A	0.481	0.004	N/A	N/A	N/A	N/A	N/A
7	N/A	N/A	0.432	0.001	N/A	N/A	N/A	N/A	N/A
8	N/A	N/A	0.455	0.000	N/A	N/A	N/A	N/A	N/A
9	N/A	N/A	0.440	0.003	N/A	N/A	N/A	N/A	N/A
10	N/A	N/A	0.465	0.004	N/A	N/A	N/A	N/A	N/A
11	N/A	N/A	0.429	0.002	N/A	N/A	N/A	N/A	N/A

Result Summary									
Pass QTY	N/A	N/A	11	11	N/A	N/A	N/A	N/A	N/A
Fail QTY	N/A	N/A	0	0	N/A	N/A	N/A	N/A	N/A
Result (P/F)	N/A	N/A	P	P	N/A	N/A	N/A	N/A	N/A
Statistical Summary									
Min	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
Comments									
Rinse test and UV Absorbance Test performed at SLQ. The rest of the visual inspections were performed at Pathway.									
Disposition:	ACCEPT				NCMR No:		N/A		
Inspected By:					Date:		1/9/2025		
Verified By:					Date:		01/09/2025		

Date Inspected	Part No.	Rev.	Description						
11/27/2024	SLQ-8102	A	SILQ 2-Way Cleartract SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack						
Lot No.	Lot QTY	ES-0059-01 Emplex Model MP56300-M Continuous Band Sealer - Validated Range: Temperature = 142 °C - 144 °C, Speed = 175 in/min., Pressure = 89 psi ES-0099-01 ACCU-SEAL Model 5300-15-B, Impulse Sealer - Validated Range: Temperature = 270 °F - 280 °F, Dwell = 4 sec., Pressure = 55 psi							
SLQ-11182024	603/60.3 packs of 10								
Specification	TM-0003								
Revision	C								
PRODUCT PREP FOR TEST									
* Seal one sample pouch P/N SLQ-5001 at the start, middle and end of each shift.									
* Cut two 1.00 inch wide x 2.00 inch long strips from the manufacturer sealed end of the pouch.									
* Execute peel strength test per TM-0003 and record results.									
Specification	≥1.0	≥1.0	142-144 °C	175 in/min	89 psi	270 - 280 °F	4.0	55.000	
Lower Spec Limit (LSL)	1.0	1.0	142.0	175.0	89	270.0	4.0	55	
Upper Spec Limit (USL)	N/A	N/A	144.0	175.0	89	280.0	4.0	55	
UOM	lb	lb	°C	in/min	psi	°F	seconds	psi	
Measuring Equipment Type	Force Gauge	Force Gauge	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	
Asset No.	ES-0028-03	ES-0028-03	N/A	N/A	N/A	ES-0099-01	ES-0099-01	ES-0099-01	
Cal Due Date	12/28/2024	12/28/2024	N/A	N/A	N/A	8/9/2024	8/9/2024	8/9/2024	
Sample ID	Top Strip	Bottom Strip	N/A	N/A	N/A	1	1	1	
Sample No.	Variable/ Attribute Data								
1: Start (Top Strip)	3.640	N/A	N/A	N/A	N/A	270.0	4.0	55.0	
2: Start (Bottom Strip)	N/A	5.000	N/A	N/A	N/A	N/A	N/A	N/A	
3: Middle (Top Strip)	5.740	N/A	N/A	N/A	N/A	N/A	N/A	N/A	
4: Middle (Bottom Strip)	N/A	4.400	N/A	N/A	N/A	N/A	N/A	N/A	
5: End (Top Strip)	7.180	N/A	N/A	N/A	N/A	N/A	N/A	N/A	
6: End (Bottom Strip)	N/A	4.780	N/A	N/A	N/A	N/A	N/A	N/A	
Result Summary									
Pass QTY	3	3	N/A	N/A	N/A	1	1	1	
Fail QTY	0	0	N/A	N/A	N/A	0	0	0	
Result (P/F)	P	P	N/A	N/A	N/A	P	P	P	
Statistical Summary									
Min	3.640	4.400	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	
Max	7.180	5.000	#REF!	#REF!	0.00	#REF!	#REF!	#REF!	
Mean	5.520	4.727	#REF!	#REF!	0.00	#REF!	#REF!	#REF!	
StDev	1.780	0.304	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	
Comments									
N/A									
Disposition:	ACCEPT					NCMR No: N/A			
Inspected By:						Date: 11/27/2024			

	Form, Run Sheet for SILQ Automation Machine		Doc No.: FM-SQL-0001
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission therefrom. If any portion of this document is in noncopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: A
			Page 1 of 1

Part Number	Revision	Description	Machine Asset No.
SLQ-8102	A	SILQ 2-Way ClearTract SPT 18Fr x 10ml Catheter, Non-Sterile, 10-Pack	ES-0091-01

Equipment Setup / Startup		Adjustments during Production		
Setting [Units]	Value	Value	Value	Value
Setup / Adjustment Date	11/18/2024	N/A ST 12/11/2024		
Setup / Adjustment Time	11:18:20 ①			
Fill Time [seconds]	9 seconds			
Confirm Cylinder Actuation All 4 pairs [P/F]	P			
UV Cure Time [seconds]	120 seconds			
Drain Time [seconds]	30 seconds			
Rinse Time [seconds]	6 seconds			
Final Drain [seconds]	30 seconds			
Number of Drain / Rinse Cycles [Count]	3			
Confirm Cylinder Actuation All 4 pairs [P/F]	P			
Size Setting	18FR			
Performed By Name	Sean Tawye			
Performed By Signature				
Verified By Name	Larry Daaron			
Verified By Signature				

Engineering / Lead Supervisor

QA

① 9:30 AM ST 11/18/2024

Date Inspected	Part No.	Rev.	Description							
01/21/2025	SLQ-211810SPT	B	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack							
Lot No.	Lot QTY									
SLQ-11182024	603/ 60.3 Packs of 10									
Identification of the product is traceable between the DHR, sterilization records, and laboratory report	DHR is completed and meets good documentation practices per GP-003	Sterilization processing record meets parameters of validated process	Sterility laboratory report meets acceptance criteria indicated in QP-010	Endotoxin laboratory report (LAL) meets acceptance criteria per V5-0021	No Damage to Outer Packing on Shipment back from Sterilizer	Verify Quantity Received back from Sterilizer matches Packing List				
Specification	QP-010 Revision	F	F	F	F	F	F			
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	20.0	N/A	N/A			
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A	N/A			
UOM	N/A	N/A	N/A	N/A	EU	N/A	N/A			
Measuring Equipment Type	Visual	Visual	Visual	Visual	Visual	Visual	Visual			
Asset No.	N/A	N/A	N/A	N/A	N/A	N/A	N/A			
Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A	N/A			
AQL	See Below	See Below	See Below	See Below	See Below	See Below	See Below			
Sample Size	DHR	DHR	Sterilization Processing Record	Sterility Laboratory Report	Endotoxin (LAL) Laboratory Report	Shipment	Shipment			
Sample No.	Variable/ Attribute Data									
1	P	P	P	P	F	P	P			
Result Summary										
Pass QTY	1	1	1	1	0	1	1			
Fail QTY	0	0	0	0	1	0	0			
Result (P/F)	P	P	P	P	F	P	P			
Statistical Summary										
Min	N/A	N/A	N/A	N/A	N/A	N/A	N/A			
Max	N/A	N/A	N/A	N/A	N/A	N/A	N/A			
Mean	N/A	N/A	N/A	N/A	N/A	N/A	N/A			
StDev	N/A	N/A	N/A	N/A	N/A	N/A	N/A			
Comments										
Disposition:	Reject						NCR No:	25-002		
Inspected By:	Larry DeCaron						Date:	01/21/2025		
Verified By:	Gustavo Siquiera						Date:	01/21/2025		



LGGS, Inc. Technical Report

SPONSOR: Pathway LLC
8779 Cottonwood Ave Suite 105
Santee, CA 92071

Lab No. 24-5200-00
P.O.No. P1975
Lot No. SLQ-11182024
Date Received: 12-24-24

Analysis: Limulus Amebocyte Lysate (LAL) Kinetic Turbidimetric Assay

Test Article: SLQ-211810SPT

Condition of Test Article at Receipt: Good

Procedure: Three (3) products were extracted following procedures outlined in test specification number, BEKT-PTW-001-00. One-tenth ml of the test and/or control solution and one-tenth of lysate solution, reconstituted per suppliers directions, were placed in sterile microplate wells and incubated at 37°C for a minimum of one hour in the Kinetic QCL Reader. The reaction time with the LAL/substrate was determined turbidimetrically and compared to a standard curve to determine endotoxin concentration. Positive Product Control (PPC) was simultaneously prepared to demonstrate that the test article did not significantly interfere with the LAL/substrate reaction.

Results:	Test Article Extract Dilution:	<u>2.0</u>
	PPC Testing:	<u>Acceptable</u>
	Test Article Extract:	<u>2.30 EU/ml; 92.0 EU/device</u>
	Date Tested:	<u>12-24-24</u>

Key to Results: Type of Product

Current FDA Requirements

1) Medical Device

1) Less than or equal to 0.5 EU/ml
Less than or equal to 20.0 EU/device

Comments: Test article as submitted does not conform to FDA and USP requirements for end-product release of medical devices.

Deviations: None

Completed: 12/26/2024 Tech RR / GS Approved by: 
Technical Manager (or deputy)

The conclusions and data herein apply only to the items tested and are not indicative of the quality or condition of apparently identical product. The report is submitted for the exclusive use of the sponsor. No use or mention of LGGS, Inc. in connection with any form of advertising or public announcement or reproduction of the report except in full may be made without written approval of LGGS, Inc.

STERIS: Ethylene Oxide Certificate Of Processing

Prepared PATHWAY INNOVATION LLC
For
Ethylene Oxide Run ID 1162-16205

<u>Product Code</u>	<u>Lot Number</u>	<u>Quantity (EA)</u>
PWI-7374 Full Cycle	SLQ-11182024	6
Customer ItemID: SLQ-211810SPT		
PWI-7374 Full Cycle	SLQ-11192024	2
Customer ItemID: SLQ-211410SPT		
PWI-7374 Full Cycle	SLQ-11202024	9
Customer ItemID: SLQ-211610SPT		

Preconditioning

<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>
Preconditioning Room	13-Dec-2024 8:07 AM	14-Dec-2024 4:27 AM	1220
			<u>Min</u> <u>Max</u>
Duration (min) Spec:			960 1920
Relative Humidity (%RH) Spec:			40.0 80.0
Actual:			47.0 65.0
Temperature (°F) Spec:			105.0 125.0
Actual:			107.0 117.0

Vessel Processing

<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>
Chamber 6	14-Dec-2024 4:37 AM	14-Dec-2024 1:32 PM	535
Cycle ID: PWI-7374	Check Value: N/A		
			<u>Min</u> <u>Max</u>
Duration (min) Spec:			1 1440

EO Drum Information

<u>Drum ID</u>	<u>Drum Serial Number</u>	<u>Drum Lot Number</u>	<u>EO Usage (LB)</u>
2188	E006254	UTLX931034K24	13.0
Total EO Usage:			13.0

EO Drum Information (Spec)

	<u>EO Usage (LB)</u>
Min	1.0
Max	100.0

Aeration

<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>	
Aeration Room	14-Dec-2024 1:42 PM	16-Dec-2024 1:57 PM	2895	
Total Aeration Time:			2895	
			<u>Min</u>	<u>Max</u>
Duration (min) Spec:			2880	4320
Aeration in Room Temperature (°F) Spec:			103.0	123.0
Actual:			110.0	114.0

EO Spec ID/Rev/Description: ODMS Spec #5 Rev 1 (PWI-7374 Rev. 1 Medical Product)

Product meets Customer specifications; zero nonconformities occurred during this run.

Customer Specification ID: 769

Signature Manifest

Reviewed and E-Signed By: **Shirley Carter (Senior Quality Technician)**

Date/Time E-signed: 18-Dec-2024 12:31 PM

Document Content Revision 1

Operating facilities are in compliance with applicable state and federal regulations (FDA, NRC, EPA, and OSHA) and provide services under a quality system which meets the requirements of FDA QSR, EN/ISO 13485, and in alignment with EN ANSI/AAMI/ISO 11135.

Processing Location

43425 Business Park Dr
Temecula CA 92590
Phone: 951-694-9340

*** End of Ethylene Oxide Certificate of Processing Document ***

STERIS: Receiving Report
PATHWAY INNOVATION LLC (8130)
Ethylene Oxide

Product Received From:
PATHWAY INNOVATION LLC
8779 Cottonwood Ave
Suite 105
Santee CA 92071
325

Carrier: NorthStar
Trailer: N/A
Bill of Lading: 4385521
Seal: None
PO(s): P1970
Product counts provided by Customer

Receipt Date/Time: 13-Dec-2024 6:38 AM
Pallet Count: 1
Unloaded By: Brian Neilson
Counted By: Diana Burns
Container Unload: ☐
Promise Date: 12-Jan-2025 6:38 AM

Receipt ID
1161-11952



Item [Pkg]	Lot Number	Quantity		
		Customer	STERIS	CD UOM
PWI-7374 Full Cycle [1] Customer Item ID: SLQ-211810SPT	SLQ-11182024	6	6	0 EA
PWI-7374 Full Cycle [1] Customer Item ID: SLQ-211410SPT	SLQ-11192024	2	2	0 EA
PWI-7374 Full Cycle [1] Customer Item ID: SLQ-211610SPT	SLQ-11202024	9	9	0 EA
Receipt Totals:		17	17	0

Inventory by Received Pallet

Pallet (SN)	Item ID [Pkg]	Quantity			Warehouse		Pallet Total
		Customer	STERIS	UOM	Location		
1 (92640) PWI-7374 Full Cycle [1] Customer Item ID: SLQ-211810SPT	SLQ-11182024	6	6	EA	N/A		1
1 (92640) PWI-7374 Full Cycle [1] Customer Item ID: SLQ-211410SPT	SLQ-11192024	2	2	EA	N/A		1
1 (92640) PWI-7374 Full Cycle [1] Customer Item ID: SLQ-211610SPT	SLQ-11202024	9	9	EA	N/A		1

STERIS: Receiving Report
PATHWAY INNOVATION LLC (8130)
Ethylene Oxide

Product Received From:

PATHWAY INNOVATION LLC
8779 Cottonwood Ave
Suite 105
Santee CA 92071
325

Carrier: NorthStar

Trailer: N/A

Bill of Lading: 4385521

Seal: None

PO(s): P1970

Product counts provided by Customer

Receipt ID
1161-11952



Receipt Date/Time: 13-Dec-2024 6:38 AM

Pallet Count: 1

Unloaded By: Brian Neilson

Counted By: Diana Burns

Container Unload: ☐

Promise Date: 12-Jan-2025 6:38 AM

Signature Manifest

Reviewed By: Diana Burns (Manager, Production)

Date/Time Signed: 13-Dec-2024 6:43 AM

Document Content Revision: 1

Processing Location

43425 Business Park Dr
Temecula CA 92590
Phone: 951-694-9340

GW1-08-004-01

Last Revised in Rel 2.0.0.0

Rel Date: 31-Mar-2017

Document ID: 56293

Page 2 of 2

PATHWAY INNOVATION LLC

8779 Cottonwood Ave
Suite 105
Santee, CA 92071
(619)415-0103

Packing List No: PL-SLQ-0052

Page No: 5 of 5

Ship To:

Company: STERIS AST
Attn: Diana Burns
Address: 43425 Business Park Drive
Temecula, CA 92590

Phone: 951-694-9340

Ship Date: 12/11/2024

POW P1970


Tracking No: TBD

Delivery Company: NorthStar Courier

Delivery Method: Courier

FOB: Santee, CA

[illegible]

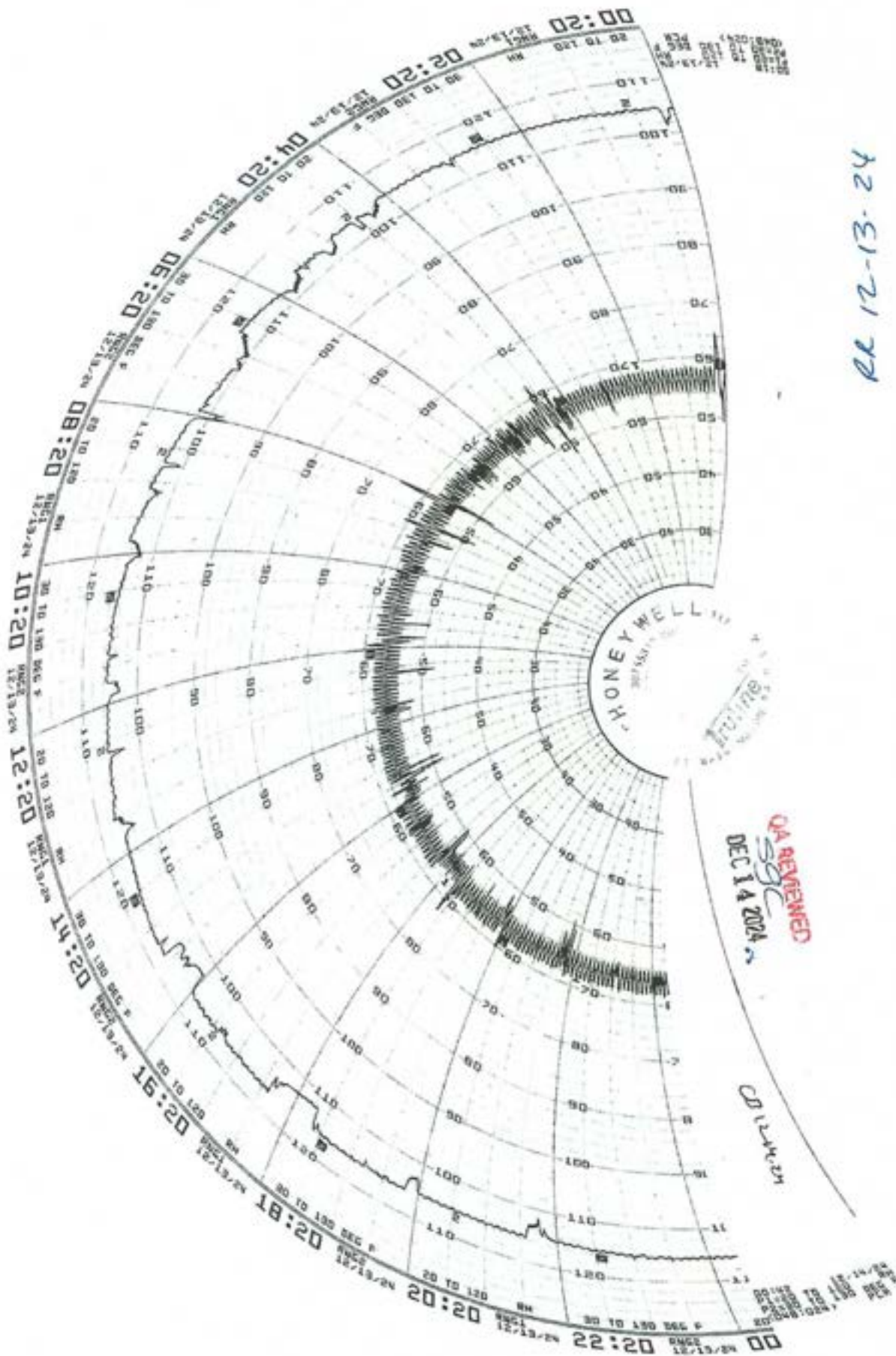

RR 12-13-24

STERIS: Ethylene Oxide Certificate Of Processing

Prepared PATHWAY INNOVATION LLC
For
Ethylene Oxide Run ID 1162-16205

Product Code	Lot Number	Quantity (EA)
PWI-7374 Full Cycle	SLQ-11182024	6
Customer ItemID:	SLQ-211810SPT	
PWI-7374 Full Cycle	SLQ-11192024	2
Customer ItemID:	SLQ-211410SPT	
PWI-7374 Full Cycle	SLQ-11202024	9
Customer ItemID:	SLQ-211610SPT	

Preconditioning

Equipment Name	Start Date	End Date	Duration (min)
Preconditioning Room	13-Dec-2024 8:07 AM	14-Dec-2024 4:27 AM	1220
Duration (min) Spec:		Min	Max
		960	1920
Relative Humidity (%RH) Spec:		40.0	80.0
		Actual:	47.0
		65.0	
Temperature (°F) Spec:		105.0	125.0
		Actual:	107.0
		117.0	

Vessel Processing

Equipment Name	Start Date	End Date	Duration (min)
Chamber 6	14-Dec-2024 4:37 AM	14-Dec-2024 1:32 PM	535
Cycle ID: PWI-7374	Check Value: N/A		
Duration (min) Spec:		Min	Max
		1	1440

EO Drum Information

Drum ID	Drum Serial Number	Drum Lot Number	EO Usage (LB)
2188	E006254	UTLX931034K24	13.0
Total EO Usage:			13.0

EO Drum Information (Spec)

EO Usage (LB)	
Min	1.0
Max	100.0

Aeration

Equipment Name	Start Date	End Date	Duration (min)	
Aeration Room	14-Dec-2024 1:42 PM	16-Dec-2024 1:57 PM	2895	
		Total Aeration Time:	2895	
			Min	Max
	Duration (min)	Spec:	2880	4320
	Aeration in Room Temperature (°F)	Spec:	103.0	123.0
		Actual:	110.0	114.0

EO Spec ID/Rev/Description ODMS Spec #5 Rev 1 (PWI-7374 Rev. 1 Medical Product)

Product meets Customer specifications; zero nonconformities occurred during this run.

Customer Specification ID: 769

Signature Manifest

Reviewed and E-Signed By: **Shirley Carter (Senior Quality Technician)**

Date/Time E-signed 18-Dec-2024 12:31 PM

Document Content Revision 1

Operating facilities are in compliance with applicable state and federal regulations (FDA, NRC, EPA, and OSHA) and provide services under a quality system which meets the requirements of FDA QSR, EN/ISO 13485, and in alignment with EN ANSI/AAMI/ISO 11135.

Processing Location

43425 Business Park Dr
Temecula CA 92590
Phone: 951-694-9340

*** End of Ethylene Oxide Certificate of Processing Document ***

STERIS: Receiving Report
PATHWAY INNOVATION LLC (8130)
Ethylene Oxide

Product Received From:

PATHWAY INNOVATION LLC
8779 Cottonwood Ave
Suite 105
San Jose CA 95071
325

Carrier: NorthStar

Trailer: N/A

Bill of Lading: 4385521

Seal: None

PO(s): P1970

Product counts provided by Customer

Receipt Date/Time: 13-Dec-2024 6:38 AM

Pallet Count: 1

Unloaded By: Brian Neilson

Counted By: Diana Burns

Container Unload: ☐

Promise Date: 12-Jan-2025 6:38 AM

Receipt ID
1161-11952



Item [Pkg]

PWI-7374 Full Cycle [1]

Customer Item ID: SLQ-211810SPT

PWI-7374 Full Cycle [1]

Customer Item ID: SLQ-211410SPT

PWI-7374 Full Cycle [1]

Customer Item ID: SLQ-211610SPT

Lot Number

SLQ-11182024

SLQ-11192024

SLQ-11202024

Quantity

Customer STERIS CD UOM

6 6 0 EA

2 2 0 EA

9 9 0 EA

Receipt Totals:

17 17 0

Inventory by Received Pallet

Pallet (SN) Item ID [Pkg]

1 (92640) PWI-7374 Full Cycle [1]

Customer Item ID: SLQ-211810SPT

1 (92640) PWI-7374 Full Cycle [1]

Customer Item ID: SLQ-211410SPT

1 (92640) PWI-7374 Full Cycle [1]

Customer Item ID: SLQ-211610SPT

Quantity

Customer STERIS UOM

6 6 EA

2 2 EA

9 9 EA

Warehouse

Location

N/A

N/A

N/A

Pallet

Total

1

1

1

STERIS: Receiving Report
PATHWAY INNOVATION LLC (8130)
Ethylene Oxide

Product Received From:

PATHWAY INNOVATION LLC
8779 Cottonwood Ave
Suite 105
Santee CA 92071
325

Carrier: NorthStar

Trailer: N/A

Bill of Lading: 4385521

Seal: None

PO(s): P1970

Product counts provided by Customer

Receipt Date/Time: 13-Dec-2024 6:38 AM

Pallet Count: 1

Unloaded By: Brian Neilson

Counted By: Diana Burns

Container Unload: ☐

Promise Date: 12-Jan-2025 6:38 AM

Receipt ID
1161-11952



Signature Manifest

Reviewed By: Diana Burns (Manager, Production)

Date/Time Signed: 13-Dec-2024 6:43 AM

Document Content Revision 1

Processing Location

43425 Business Park Dr
Temecula CA 92590
Phone: 951-694-9340

GW1-08-004-01

Last Revised in Rel 2.0.0.0

Rel Date: 31-Mar-2017

Document ID: 56293

Page 2 of 2

PATHWAY INNOVATION LLC

8779 Cottonwood Ave
Suite 105
Santee, CA 92071
(619)415-0103

Packing List No: PL-SLD-0052

Page No.: 1 of 1

Ship To:

Company: STERIS AST
Attn: Diana Burns
Address: 43425 Business Park Drive
Temecula, CA 92590

Phone: 951-694-9340

Ship Date: 12/11/2024

PO# P1970


Tracking No: TBD

Delivery Company: NorthStar Courier

Delivery Method: Courier

FOB: Santee, CA


[illegible]

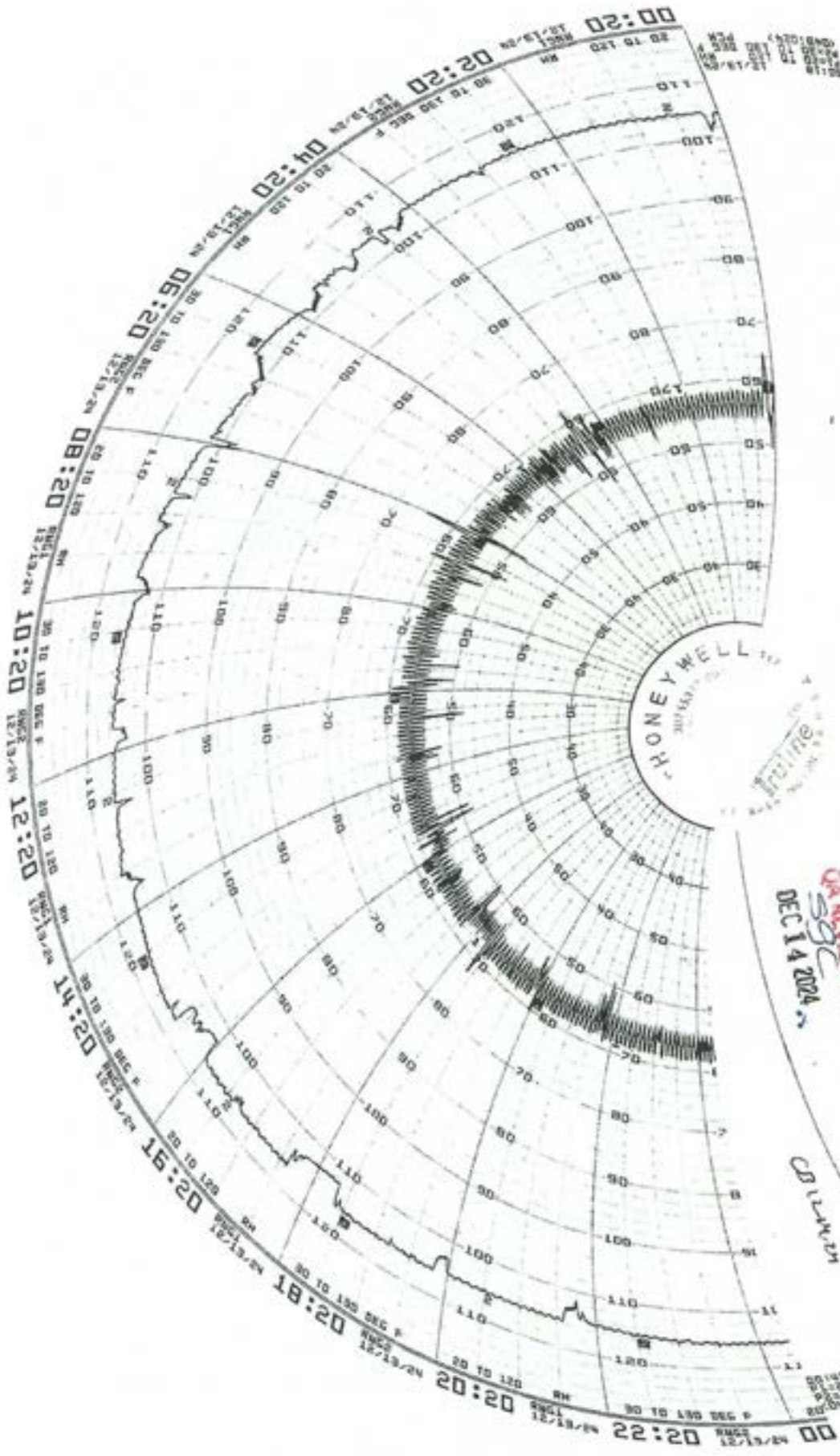

QA REVIEWED
DEC 14 2024

CD 11/13/24

RL 12-13-24

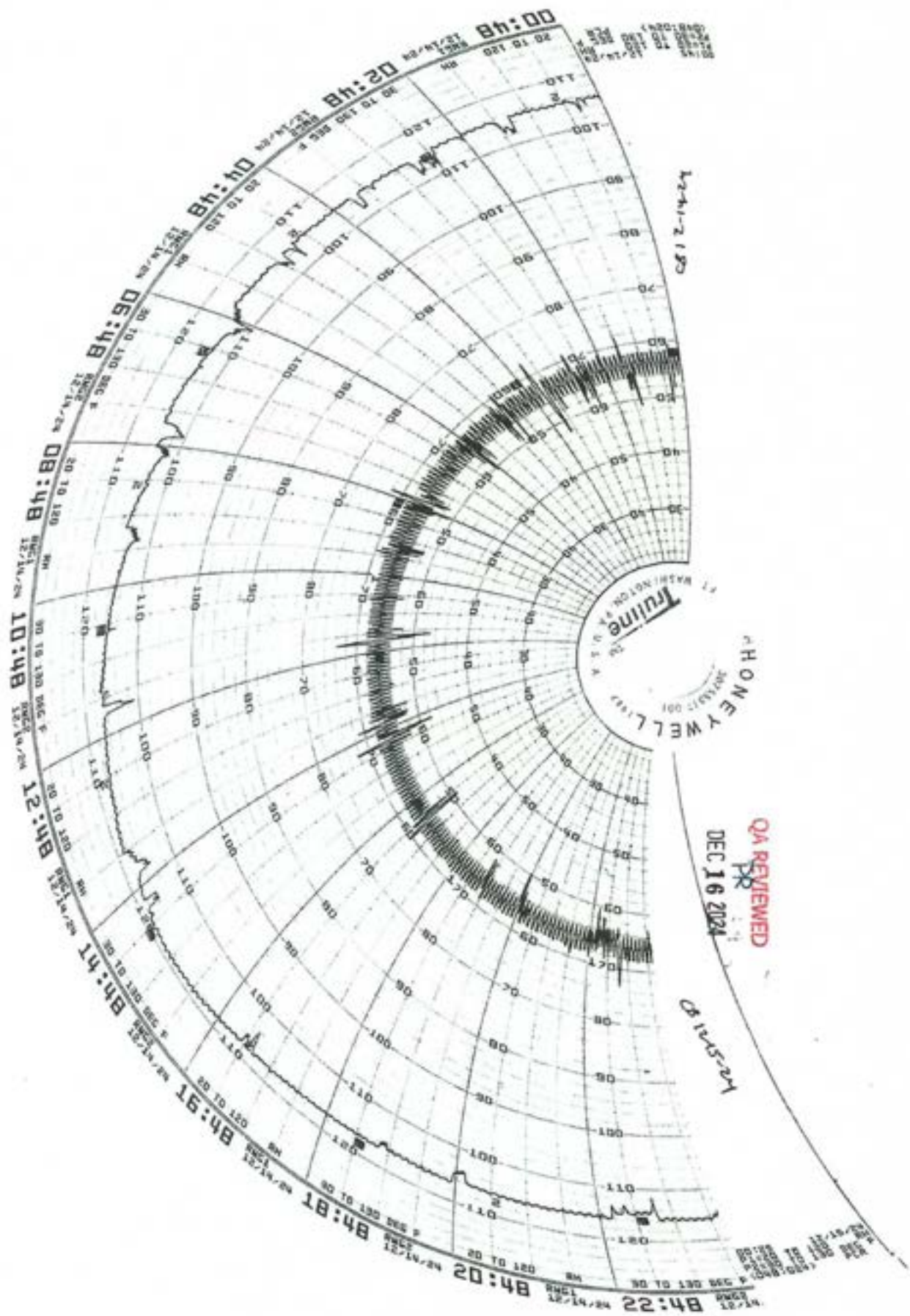

Customer Name Pathway FC
Load # PWI 1162-16205
Date Processed 12/14/2024 04:37:45
Pressure Units "HgA

Cycle Number 7374
Cycle Last Modified: 4/30/2018 1:56:49PM
Chamber # 6
Cycle Counter # 4,784
Temperature Units Fahrenheit

Phase Name

Initial Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 04:37:45	00:00:00	29.1	126
12/14/2024 04:42:45	00:05:00	21.7	124
12/14/2024 04:47:45	00:10:00	15.0	124
12/14/2024 04:52:45	00:15:00	8.2	123
12/14/2024 04:57:45	00:20:00	3.1	124
12/14/2024 04:59:54	00:22:08	2.0	124

Summary of Initial Evacuation

Pressure Start	29.1	Min. Pressure	2.0	Min Temp	123	Pressure Change	-27.1
Pressure End	2.0	Max Pressure	29.1	Max Temp	126	Phase Time	00:22:08

Stabilization

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 04:59:54	00:22:09	2.0	124
12/14/2024 05:01:53	00:24:09	2.1	125
12/14/2024 05:01:53	00:24:09	2.1	125

Summary of Stabilization

Pressure Start	2.0	Min. Pressure	2.0	Min Temp	124	Pressure Change	0.1
Pressure End	2.1	Max Pressure	2.1	Max Temp	125	Phase Time	00:02:00

Leak Check

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 05:01:53	00:24:09	2.1	125
12/14/2024 05:02:54	00:25:09	2.1	125
12/14/2024 05:03:53	00:26:09	2.1	125
12/14/2024 05:04:54	00:27:09	2.2	126
12/14/2024 05:05:54	00:28:09	2.2	126
12/14/2024 05:06:54	00:29:09	2.2	126

Summary of Leak Check

Pressure Start	2.1	Min. Pressure	2.1	Min Temp	125	Pressure Change	0.1
Pressure End	2.2	Max Pressure	2.2	Max Temp	126	Phase Time	00:05:00

Initial Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 05:06:54	00:29:09	2.2	126
12/14/2024 05:11:54	00:34:09	9.2	128
12/14/2024 05:16:54	00:39:09	16.8	128
12/14/2024 05:21:53	00:44:09	24.2	128
12/14/2024 05:24:27	00:46:42	28.0	128

Summary of Initial Nitrogen Wash

Pressure Start	2.2	Min Pressure	2.2	Min Temp	126	Pressure Change	25.8
Pressure End	28.0	Max Pressure	28.0	Max Temp	128	Phase Time	00:17:33

Initial Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 05:24:27	00:46:42	28.0	128
12/14/2024 05:29:27	00:51:42	21.3	124
12/14/2024 05:34:27	00:56:42	14.3	123
12/14/2024 05:39:27	01:01:42	7.3	123
12/14/2024 05:44:27	01:06:42	2.7	123
12/14/2024 05:45:53	01:08:08	2.0	123

Summary of Initial Evacuation

Pressure Start	28.0	Min Pressure	2.0	Min Temp	122	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	128	Phase Time	00:21:27

Humidification

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	Avg. % RH
12/14/2024 05:45:53	01:08:08	2.0	123	5
12/14/2024 05:46:53	01:09:08	2.5	124	21
12/14/2024 05:47:54	01:10:09	2.9	125	31
12/14/2024 05:48:05	01:10:20	3.0	125	31

Summary of Humidification

Pressure Start	2.0	Min Pressure	2.0	Min Temp	123	Pressure Change	1.0	Min RH	4
Pressure End	3.0	Max Pressure	3.0	Max Temp	125	Phase Time	00:02:11	Max RH	32

Humidity Dwell

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	Avg. % RH
12/14/2024 05:48:05	01:10:20	3.0	125	31
12/14/2024 05:53:05	01:15:20	3.0	125	35
12/14/2024 05:58:05	01:20:20	3.0	125	36
12/14/2024 06:03:04	01:25:20	3.0	125	35
12/14/2024 06:08:05	01:30:20	3.0	125	35
12/14/2024 06:13:05	01:35:20	3.0	126	34
12/14/2024 06:18:05	01:40:20	3.0	125	35
12/14/2024 06:23:05	01:45:20	3.1	125	35
12/14/2024 06:28:06	01:50:20	3.1	125	35
12/14/2024 06:33:06	01:55:20	3.1	125	35
12/14/2024 06:38:05	02:00:20	3.1	125	35
12/14/2024 06:43:05	02:05:21	3.1	125	35
12/14/2024 06:48:04	02:10:20	3.1	125	36

Summary of Humidity Dwell

Pressure Start	3.0	Min Pressure	3.0	Min Temp	125	Pressure Change	0.1	Min RH	31
Pressure End	3.1	Max Pressure	3.1	Max Temp	126	Phase Time	01:00:00	Max RH	36

Sterilant Injection

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	EO Drum 1 (lb)	EO Drum 2 (lb)	EO CONC
12/14/2024 06:48:04	02:10:20	3.1	125	631	0	0
12/14/2024 06:53:05	02:15:20	7.7	126	627	0	223
12/14/2024 06:58:05	02:20:20	12.8	125	620	0	522
12/14/2024 06:59:34	02:21:50	14.3	125	618	0	606

Customer Name Pathway FC
Load # PWI 1162-16205

STERIS
Applied Sterilization Technologies
Temecula, CA

Summary of Sterilant Injection

Pressure Start 3.1	Min Pressure 3.1	Min Temp 125	Drum 1 Start 631	Drum 2 Start 0
Pressure End 14.3	Max Pressure 14.3	Max Temp 126	Drum 1 End 618	Drum 2 End 0
			EIO Used 1 13	EIO Used 2 0
Pressure Change 11.2	MIN EO CONC 0.00			
Phase Time 00:11:30	MAX EO CONC 606.00			
			Total EO Weight 13	

Exposure

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	EO Drum 1 (lb)	EO Drum 2 (lb)	EO CONC
12/14/2024 06:59:34	02:21:50	14.3	125	618	0	606
12/14/2024 07:04:35	02:26:50	14.2	124	618	0	627
12/14/2024 07:09:35	02:31:50	14.1	125	618	0	622
12/14/2024 07:14:35	02:36:50	14.1	125	618	0	619
12/14/2024 07:19:35	02:41:50	14.3	125	618	0	628
12/14/2024 07:24:35	02:46:50	14.3	125	617	0	625
12/14/2024 07:29:35	02:51:50	14.3	125	617	0	622
12/14/2024 07:34:35	02:56:50	14.3	125	617	0	619
12/14/2024 07:39:35	03:01:50	14.3	125	617	0	619
12/14/2024 07:44:36	03:06:50	14.3	125	617	0	616
12/14/2024 07:49:36	03:11:51	14.2	125	617	0	616
12/14/2024 07:54:36	03:16:51	14.2	125	617	0	613
12/14/2024 07:59:35	03:21:51	14.2	125	617	0	612
12/14/2024 08:04:35	03:26:51	14.2	125	616	0	611
12/14/2024 08:09:36	03:31:51	14.2	125	616	0	609
12/14/2024 08:14:36	03:36:51	14.2	125	616	0	608
12/14/2024 08:19:36	03:41:51	14.2	125	616	0	608
12/14/2024 08:24:36	03:46:51	14.2	125	616	0	606
12/14/2024 08:29:37	03:51:51	14.2	125	616	0	606
12/14/2024 08:34:36	03:56:51	14.2	125	616	0	606
12/14/2024 08:39:36	04:01:51	14.2	125	616	0	605
12/14/2024 08:44:36	04:06:51	14.2	125	616	0	605
12/14/2024 08:49:37	04:11:51	14.2	125	616	0	603
12/14/2024 08:54:37	04:16:52	14.2	125	616	0	603
12/14/2024 08:59:37	04:21:52	14.2	125	616	0	603
12/14/2024 09:04:36	04:26:52	14.2	125	616	0	602
12/14/2024 09:09:36	04:31:52	14.2	125	615	0	603
12/14/2024 09:14:37	04:36:52	14.2	125	615	0	601
12/14/2024 09:19:37	04:41:52	14.2	125	615	0	601
12/14/2024 09:24:37	04:46:52	14.2	125	615	0	601
12/14/2024 09:29:37	04:51:52	14.2	125	615	0	599
12/14/2024 09:34:38	04:56:52	14.2	125	615	0	600
12/14/2024 09:39:37	05:01:52	14.2	125	615	0	600
12/14/2024 09:44:37	05:06:52	14.2	125	615	0	599
12/14/2024 09:49:37	05:11:52	14.2	125	615	0	600
12/14/2024 09:54:37	05:16:53	14.3	125	615	0	597
12/14/2024 09:59:38	05:21:53	14.3	125	615	0	597
12/14/2024 10:04:38	05:26:53	14.3	125	615	0	597
12/14/2024 10:09:38	05:31:53	14.3	125	615	0	597
12/14/2024 10:14:37	05:36:53	14.3	125	615	0	596
12/14/2024 10:19:38	05:41:53	14.3	125	615	0	595
12/14/2024 10:24:38	05:46:53	14.3	125	615	0	594
12/14/2024 10:29:38	05:51:53	14.3	125	614	0	594
12/14/2024 10:34:38	05:56:53	14.3	125	614	0	594
12/14/2024 10:39:39	06:01:53	14.3	125	614	0	593
12/14/2024 10:44:39	06:06:53	14.3	125	614	0	593
12/14/2024 10:49:38	06:11:53	14.3	125	614	0	593
12/14/2024 10:54:38	06:16:53	14.3	125	614	0	593
12/14/2024 10:59:35	06:21:50	14.3	125	618	0	592

Customer Name Pathway FC
Load # PWI 1162-16205

STERIS
Applied Sterilization Technologies
Temecula, CA

Summary of Exposure

Pressure Start 14.3	Min Pressure 14.1	Min Temp 124	Drum 1 Start 618	Drum 2 Start 0
Pressure End 14.3	Max Pressure 14.4	Max Temp 125	Drum 1 End 618	Drum 2 End 0
			EtO Used 1 0	EtO Used 2 0
Pressure Change 0.0	MIN EO CONC 592.00			
Phase Time 04:00:00	MAX EO CONC 630.00			
			Total EO Weight 0	

Sterilant Removal

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 10:59:35	06:21:50	14.3	125
12/14/2024 11:04:34	06:26:50	9.9	123
12/14/2024 11:09:35	06:31:50	5.3	122
12/14/2024 11:14:28	06:36:43	2.0	122

Summary of Sterilant Removal

Pressure Start 14.3	Min Pressure 2.0	Min Temp 122	Pressure Change -12.3
Pressure End 2.0	Max Pressure 14.3	Max Temp 125	Phase Time 00:14:53

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 11:14:28	06:36:43	2.0	122
12/14/2024 11:19:27	06:41:43	9.8	126
12/14/2024 11:24:28	06:46:43	17.0	126
12/14/2024 11:29:28	06:51:43	24.3	127
12/14/2024 11:31:59	06:54:15	28.0	127

Summary of Post Nitrogen Wash

Pressure Start 2.0	Min Pressure 2.0	Min Temp 122	Pressure Change 26.0
Pressure End 28.0	Max Pressure 28.0	Max Temp 127	Phase Time 00:17:32

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 11:31:59	06:54:15	28.0	127
12/14/2024 11:36:59	06:59:15	21.3	123
12/14/2024 11:42:00	07:04:15	14.3	122
12/14/2024 11:47:00	07:09:15	7.3	121
12/14/2024 11:52:00	07:14:15	2.7	121
12/14/2024 11:53:34	07:15:49	2.0	122

Summary of Post Evacuation

Pressure Start 28.0	Min Pressure 2.0	Min Temp 121	Pressure Change -26.0
Pressure End 2.0	Max Pressure 28.0	Max Temp 127	Phase Time 00:21:34

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 11:53:34	07:15:49	2.0	122
12/14/2024 11:58:34	07:20:49	9.1	126
12/14/2024 12:03:35	07:25:49	16.6	127
12/14/2024 12:08:34	07:30:49	24.1	127
12/14/2024 12:11:15	07:33:30	28.0	127

Summary of Post Nitrogen Wash

Pressure Start 2.0	Min Pressure 2.0	Min Temp 122	Pressure Change 26.0
Pressure End 28.0	Max Pressure 28.0	Max Temp 127	Phase Time 00:17:41

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 12:11:15	07:33:30	28.0	127
12/14/2024 12:16:15	07:38:30	21.3	123
12/14/2024 12:21:15	07:43:30	14.3	122
12/14/2024 12:26:15	07:48:30	7.3	121
12/14/2024 12:31:15	07:53:30	2.7	121
12/14/2024 12:32:43	07:54:58	2.0	121

Summary of Post Evacuation

Pressure Start	28.0	Min. Pressure	2.0	Min. Temp	121	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:21:28

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 12:32:43	07:54:58	2.0	121
12/14/2024 12:37:43	07:59:58	9.1	126
12/14/2024 12:42:42	08:04:58	16.6	127
12/14/2024 12:47:43	08:09:58	24.1	127
12/14/2024 12:50:21	08:12:36	28.0	127

Summary of Post Nitrogen Wash

Pressure Start	2.0	Min. Pressure	2.0	Min. Temp	121	Pressure Change	26.0
Pressure End	28.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:17:38

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 12:50:21	08:12:36	28.0	127
12/14/2024 12:55:21	08:17:36	21.3	123
12/14/2024 13:00:21	08:22:36	14.3	122
12/14/2024 13:05:21	08:27:36	7.3	121
12/14/2024 13:10:21	08:32:36	2.7	121
12/14/2024 13:11:44	08:33:59	2.0	121

Summary of Post Evacuation

Pressure Start	28.0	Min. Pressure	2.0	Min. Temp	120	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:21:23

Air Backfill

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
12/14/2024 13:11:44	08:33:59	2.0	121
12/14/2024 13:16:45	08:38:59	7.8	125
12/14/2024 13:21:44	08:44:00	14.1	125
12/14/2024 13:26:44	08:49:00	20.4	124
12/14/2024 13:31:44	08:54:00	26.5	124
12/14/2024 13:32:58	08:55:13	27.9	124

Summary of Air Backfill

Pressure Start	2.0	Min. Pressure	2.0	Min. Temp	121	Pressure Change	25.9
Pressure End	27.9	Max Pressure	27.9	Max Temp	125	Phase Time	00:21:13

Customer Name Pathway FC
Load # PWI 1162-16205

STERIS
Applied Sterilization Technologies
Temecula, CA

End of Cycle

Production Review



Signature


Date 12/14/2024

Quality Review



Signature


Date 12/18/24

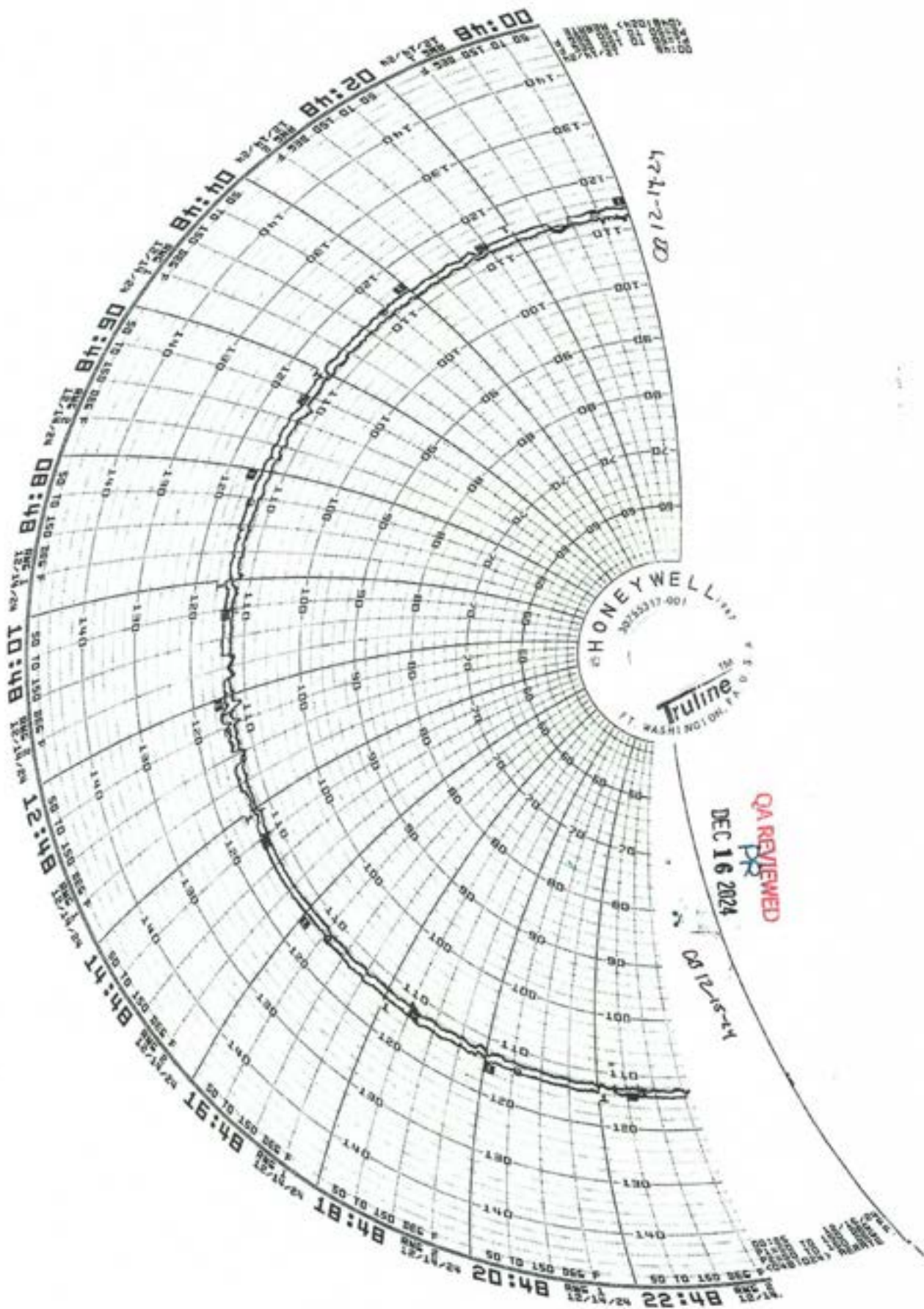


QA REVIEWED

DEC 16 2024

081215-24

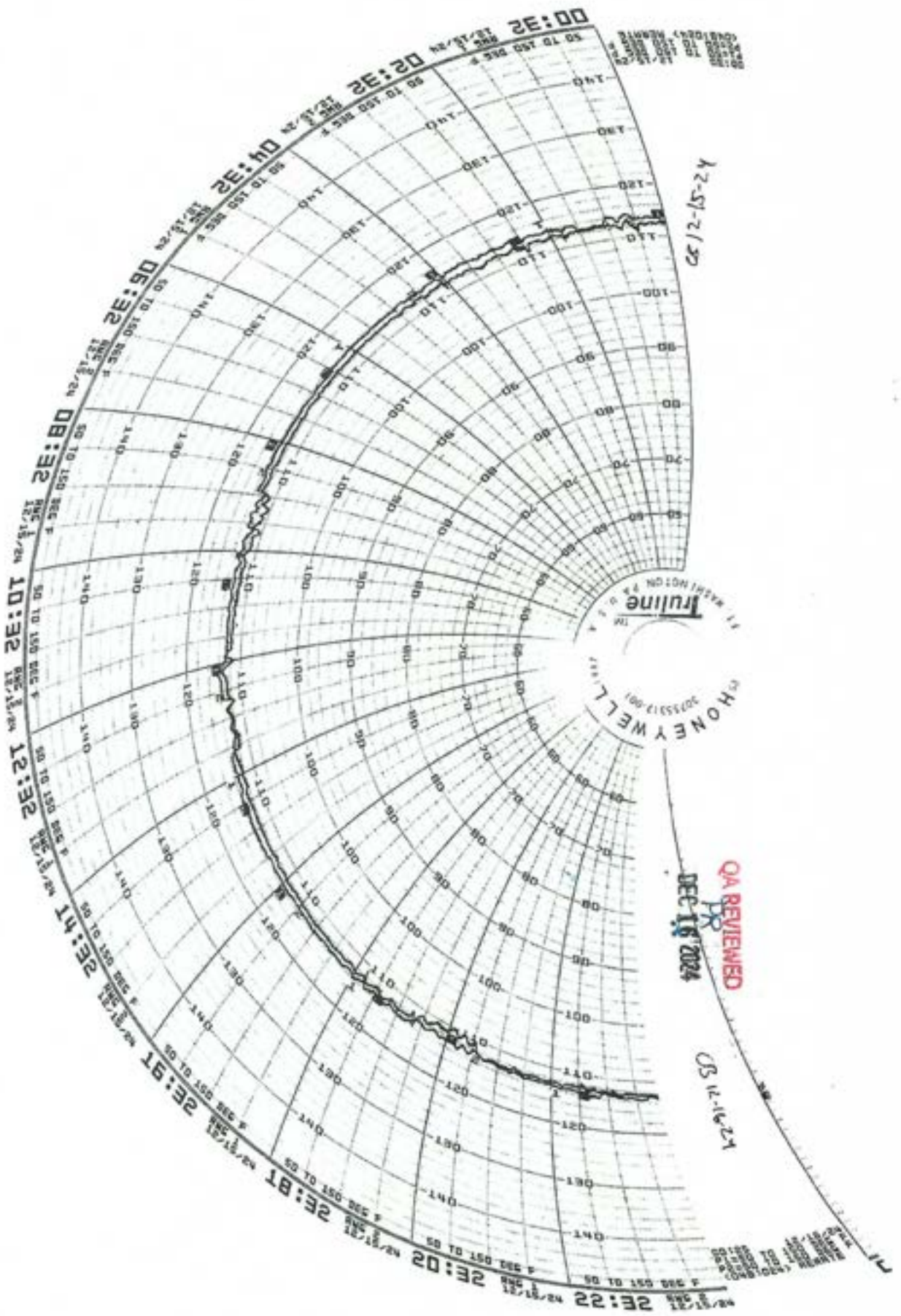

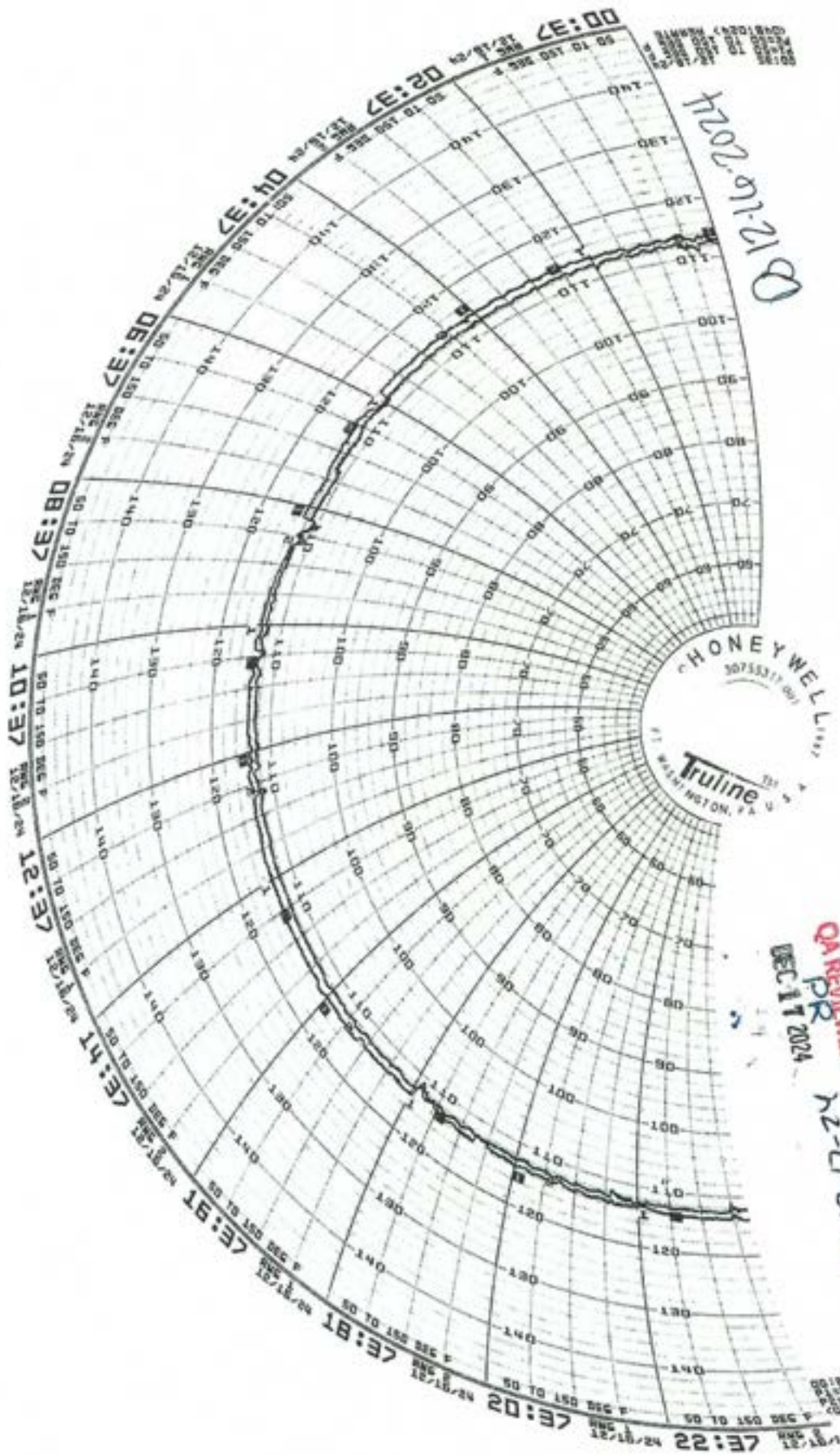


QA REVIEWED
PR
DEC 17 2024
H2-LI-21 20

PHASE 1 - PRE-INVESTIGATION			
NCR INITIATION			
NCR Owner (Full Name & Title):		Larry Dacoron, Sr. Quality Engineer	Date Issued: 01/21/2025
Supplier Name:		Ningbo Yinkmed	Qty Received: 603
Customer Name:		SILQ	Customer Approval Required? ☒ Yes ☐ No
Acceptance Criteria:		FDA Requirements: Less than or equal to 0.5 EU/ml Less than or equal to 20.0 EU/device	
Description of Nonconformance:		Test Article: SLQ-211810SPT Results: 92.0 EU/device	
NCR Source Type(s):		☐ Audits (internal & external) ☐ CAPAs ☐ Complaints ☐ Continuous Improvements ☐ Documentation Reviews ☐ Incoming Inspections ☒ Suppliers ☐ Management Reviews ☐ SCARs ☐ Other NCRs	
Source Reference #: (e.g., CAPA-24-001, Audit-24-001)		LGGS, Inc. Technical Report Lab No, 24-5200-00	
RISK ANALYSIS			
Product Risk Assessment (to be completed by Customer)			
Product Severity Rating:		☐ (1) No Patient Effect ☐ (2) Minor Injury ☐ (3) Moderate Injury ☒ (4) Serious Injury ☐ (5) Life Threatening Injury/Death/Disability	
Probability Rating:		☒ (1) Improbable ☐ (2) Remote ☐ (3) Occasional ☐ (4) Likely ☐ (5) Frequent	
Detection Rating:		☒ (1) Very High ☐ (2) High ☐ (3) Moderate ☐ (4) Low ☐ (5) Remote	
CORRECTION / CONTAINMENT			
Containment: Interim containment actions taken to limit the spread of the problem.		Material has been identified and quarantined.	
Interim Disposition:		☐ Inspect to New Revision ☐ Increase Sampling ☐ Sort-to-Accept-As-Is ☐ Sort-to-Scrap ☐ Sort-to-Return to Supplier ☐ Rework-to-Accept or Scrap ☒ Hold ☐ Use-as-Is	
PHASE 1 APPROVAL			
NCR Owner (Name/Title):	Larry Dacoron, Sr. Quality Engineer	Signature/Date:	 Larry Dacoron (11-Feb-2025 11:45 PST)
Customer (Name/Title):	Ethan Rao,	Signature/Date:	 Ethan Rao
Quality Assurance (Name/Title):	Aaron Rogers, Director of Quality and Regulatory	Signature/Date:	 Aaron Rogers (11-Feb-2025 11:48 PST)

PHASE 2 - INVESTIGATION	
Investigation Plan(s)	1. Equipment review: Identify specific equipment that may contribute to contamination 2. Catheters: Investigate the potential endotoxin contamination in catheters 3. Proprietary Solution: Focus on proprietary solution used with catheters
Investigation Result(s)	1. Tubing Issue: The tubing contains stagnant water that may remain for extended periods. Additionally, there is no planned maintenance (PM) for cleaning or replacing the tubing. 2. LAL Testing on Pre-Coated Catheters: LAL testing performed on catheters before they are coated with the proprietary solution showed that they passed the LAL test. 3. LAL Testing on Coated Catheters: LAL testing performed on catheters after coating with the proprietary solution also confirmed that they passed the LAL test.
Is a CAPA or SCAR required?	☐ Yes (refer to FM-0020, CAPA Form or FM-0087, SCAR Form & reference CAPA/SCAR #): ☐ No, provide justification: <i>The endotoxin issue is determine to be a one-time or isolated occurrence.</i>
Defect Codes: <i>check all that apply</i>	☐ Code 1 – Documentation (Missing documents, inadequate content, wrong Standard/Regulation) ☐ Code 2 – Device Defect (Any defect with actual product including damaged packaging) ☐ Code 3 - Data Entry Error (Manual or electronic information input) ☐ Code 4 – Expiration Date (Product expiration) ☐ Code 5 – Revision (Differences in identification of the “method of control”) ☐ Code 6 – Procedural (Procedures not being followed) ☐ Code 7 – Recurring NCR (Same NCR issue that was repeated) ☐ Code 8 – Tooling & Equipment Error (Expired calibration, unqualified tool, inadequate preventive maintenance) ☐ Code 9 – User Error (Incorrect product use/build by Operator) ☒ Code 10 – Requirement (Product requirement issues) ☐ Code 11 – Process Defect (NCR during manufacturing of product) ☐ Code 12 – GDP (Good Documentation Practice issues) ☐ Code 13 – Risk Analysis (Risk analysis evaluation issues) ☐ Code 14 – Facility/Environmental (Issues with Manufacturing Facility/Environment) ☐ Code 15 – Other:
Root Cause Codes: <i>check all that apply</i>	☐ Code 1 – Man (Inadequate Training) ☐ Code 2 – Man (Procedure Not Used/Used Improperly/Step Omitted) ☐ Code 3 – Man (Lack of Qualification, Demonstrated Proficiency) ☐ Code 4 – Method (Design/Inadequate Specifications) ☐ Code 5 – Method (Procedure Incomplete/Wrong/Ambiguous/Non-existent) ☐ Code 6 – Method (Manufacturing Process Controls) ☐ Code 7 – Method (Training Materials/Process) ☐ Code 8 – Method (Process Definition) ☐ Code 9 – Material (Defective Component/Material) ☐ Code 10 – Material (Incorrect Raw Material) ☐ Code 11 – Material (Expired/Damaged/Eroded) ☐ Code 12 – Material (Specification Incomplete/Wrong/Ambiguous) ☐ Code 13 – Machine (Defective Machine/Tool) ☐ Code 14 – Machine (Machine/Tool Capability) ☒ Code 15 – Machine (Inadequate Maintenance) ☐ Code 16 – Machine (Equipment Controls and Calibration) ☐ Code 17 – Environmental Controls (Temperature/Humidity) ☐ Code 18 – Environmental Controls (Facilities/Clean Room) ☐ Code 19 – Miscellaneous (Unable to determine Root Cause) ☐ Code 20 – Miscellaneous (Closed without Root Cause Determination) ☐ Code 21 - Miscellaneous (Recertification)

PHASE 2 - INVESTIGATION			
		☐ Code 22 - Miscellaneous (Remediation) ☐ Code 23 – Other (Design) ☐ Code 24 – Other (Labeling) ☐ Code 25 – Other (Manufacturing) ☐ Code 26 – Other (Not Confirmed) ☐ Code 27 – Other (User/Use Error) ☐ Code 28 – Other (Supplier Design) ☐ Code 29 – Other (Supplier Manufacture) ☐ Code 30 – Other (Supplier Labeling) ☐ Code 31 – Other (Wear and Tear) ☐ Code 32 – Other (Cause by Other) ☐ Code 33 – Other (Handling Damage) ☐ Code 34 – Other (Undeterminable) ☐ Code 35 – Other: <	
Corrective Action(s) identified to address the root cause(s): <i>Attach any supporting documents.</i>		Replaced the tubing that was holding stagnant water to prevent potential contamination	
Corrective Action(s) Reference: <i>ECO #, CAPA #, SCAR #, None</i>		None	
Is Rework required?		☐ Yes (refer to Rework Instructions below) ☒ No, provide justification and check the N/A Box below	
PHASE 2 APPROVAL			
NCR Owner (Name/Title):	Larry Dacoron, Sr. Quality Engineer	Signature/Date:	 Larry Dacoron (11-Feb-2025 11:45 PST)
MFG/Operations (Name/Title):	Sean Tamayo, Production Manager	Signature/Date:	 Sean Tamayo (11-Feb-2025 11:48 PST)
Engineering (Name/Title):	Will Hail, Sr. R&D Engineer	Signature/Date:	 Will Hail (11-Feb-2025 12:13 PST)
Quality Assurance (Name/Title):	Larry Dacoron, Sr. Quality Engineer	Signature/Date:	 Larry Dacoron (11-Feb-2025 11:45 PST)
Regulatory Affairs (Name/Title): ☐ N/A	Aaron Rogers, Director of Regulatory and Quality	Signature/Date:	 Aaron Rogers (11-Feb-2025 11:48 PST)
Customer (Name/Title): ☐ N/A	Ethan Rao	Signature/Date:	
REWORK/CORRECTION INSTRUCTIONS (☒ N/A if not applicable)			
Provide Rework Instructions (include pictures as required):			

REWORK OUTCOME (☐ N/A if not applicable)						
Acceptance Criteria (include document # and revision)	Description of NCR (material/product/process)	Part / Equipment / MP #	Part Rev	Part / Equipment / MP Description	Lot #	NCR Qty
FDA Requirements: Less than or equal to 0.5 EU/ml Less than or equal to 20.0 EU/device	Test Article: SLQ-211810SPT Results: 92.0 EU/device	SLQ-211810SPT	B	SILQ 2 way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack	SLQ-11182024	603
REWORK OUTCOME / FINAL APPROVAL						
Inspected By (Name/Title):	Larry Dacoron, Sr. Quality Engineer			Signature/Date :	 Larry Dacoron (11-Feb-2025 11:45 PST)	
MFG/Operations (Name/Title):	Sean Tamayo, Production Manager			Signature/Date :	 Sean Tamayo (11-Feb-2025 11:48 PST)	
Engineering (Name/Title):	Will Hail, Sr. R&D Engineer			Signature/Date :	 Will Hail (11-Feb-2025 12:13 PST)	
Quality Assurance (Name/Title):	Aaron Rogers, Director of Regulatory and Quality			Signature/Date :	 Aaron Rogers (11-Feb-2025 11:48 PST)	

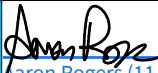PHASE 3 - CLOSURE			
Verify all information is accurate, required documentation is completed and attached?	☒ Yes ☐ No or N/A, provide Justification below		
Justification:	N/A		
Corrective action(s) is properly referenced and completed?	☒ Yes ☐ No or N/A, provide Justification below		
Justification:	N/A		
Affected material disposition is completed?	☒ Yes ☐ No or N/A, provide Justification below		
Justification:	N/A		
Notify Supplier/Customer of NCR:	☒ Yes, provide a copy of completed NCR to Supplier/Customer ☐ N/A		
NCR Log Verification:	☒ Verify NCR Log is updated		
PHASE 3 APPROVAL			
Quality Assurance (Name/Title):	Aaron Roger, Director of Regulatory and Quality	Signature/Date:	 Aaron Rogers (11-Feb-2025 15:48 PST)
Comments (if any):			
NCR CANCELLATION (☒ N/A if not applicable)			
Cancellation?	☐ Yes, provide justification below		
Justification:			
NCR CANCELLATION APPROVAL (☒ N/A if not applicable)			
Quality Assurance (Name/Title):		Signature/Date:	

Created:	2025-02-11 (Pacific Standard Time)
By:	Larry Dacoron (ldacoron@pathwaynpi.com)
Status:	Approved
Transaction ID:	CBJCHBCAABAAOTX9DgQXyCFTUSxyU00uf_jE7HPGsIX3
Number of Documents:	1
Document page count:	5
Number of supporting files:	0
Supporting files page count:	0

"NCR-25-002" History

-  Document created by Larry Dacoron (ldacoron@pathwaynpi.com)
2025-02-11 - 11:39:53 AM PST
-  Document emailed to Ethan Rao (ethanr@silq.tech) for approval
2025-02-11 - 11:43:26 AM PST
-  Document emailed to Larry Dacoron (ldacoron@pathwaynpi.com) for approval
2025-02-11 - 11:43:26 AM PST
-  Document emailed to Will Hail (whail@pathwaymedtech.com) for approval
2025-02-11 - 11:43:27 AM PST
-  Document emailed to stamayo@pathwaymedtech.com for approval
2025-02-11 - 11:43:27 AM PST
-  Document emailed to Aaron Rogers (arogers@pathwaymedtech.com) for approval
2025-02-11 - 11:43:27 AM PST
-  Agreement viewed by Larry Dacoron (ldacoron@pathwaynpi.com)
2025-02-11 - 11:43:27 AM PST
-  Larry Dacoron (ldacoron@pathwaynpi.com) has agreed to the terms of use and to do business electronically with Pathway Medtech, LLC
2025-02-11 - 11:45:10 AM PST
-  Document approved by Larry Dacoron (ldacoron@pathwaynpi.com)
Approval Date: 2025-02-11 - 11:45:10 AM PST - Time Source: server

-  Email viewed by Aaron Rogers (arogers@pathwaymedtech.com)
2025-02-11 - 11:45:32 AM PST
-  Aaron Rogers (arogers@pathwaymedtech.com) authenticated with Email OTP by verifying one-time code sent via email
2025-02-11 - 11:46:04 AM PST
-  Agreement viewed by Aaron Rogers (arogers@pathwaymedtech.com)
2025-02-11 - 11:46:05 AM PST
-  Email viewed by stamayo@pathwaymedtech.com
2025-02-11 - 11:46:59 AM PST
-  Aaron Rogers (arogers@pathwaymedtech.com) has agreed to the terms of use and to do business electronically with Pathway Medtech, LLC
2025-02-11 - 11:48:11 AM PST
-  Document approved by Aaron Rogers (arogers@pathwaymedtech.com)
Approval Date: 2025-02-11 - 11:48:11 AM PST - Time Source: server
-  stamayo@pathwaymedtech.com authenticated with Email OTP by verifying one-time code sent via email
2025-02-11 - 11:48:18 AM PST
-  Agreement viewed by stamayo@pathwaymedtech.com
2025-02-11 - 11:48:19 AM PST
-  Signer stamayo@pathwaymedtech.com entered name at signing as Sean Tamayo
2025-02-11 - 11:48:48 AM PST
-  Sean Tamayo (stamayo@pathwaymedtech.com) has agreed to the terms of use and to do business electronically with Pathway Medtech, LLC
2025-02-11 - 11:48:50 AM PST
-  Document approved by Sean Tamayo (stamayo@pathwaymedtech.com)
Approval Date: 2025-02-11 - 11:48:50 AM PST - Time Source: server
-  Email viewed by Ethan Rao (ethanr@silq.tech)
2025-02-11 - 12:00:55 PM PST
-  Ethan Rao (ethanr@silq.tech) authenticated with Email OTP by verifying one-time code sent via email
2025-02-11 - 12:01:23 PM PST
-  Agreement viewed by Ethan Rao (ethanr@silq.tech)
2025-02-11 - 12:01:24 PM PST

-  Ethan Rao (ethanr@silq.tech) has agreed to the terms of use and to do business electronically with Pathway Medtech, LLC
2025-02-11 - 12:01:42 PM PST
-  Document approved by Ethan Rao (ethanr@silq.tech)
Approval Date: 2025-02-11 - 12:01:42 PM PST - Time Source: server
-  Email viewed by Will Hail (whail@pathwaymedtech.com)
2025-02-11 - 12:12:24 PM PST
-  Will Hail (whail@pathwaymedtech.com) authenticated with Email OTP by verifying one-time code sent via email
2025-02-11 - 12:12:55 PM PST
-  Agreement viewed by Will Hail (whail@pathwaymedtech.com)
2025-02-11 - 12:12:56 PM PST
-  Will Hail (whail@pathwaymedtech.com) has agreed to the terms of use and to do business electronically with Pathway Medtech, LLC
2025-02-11 - 12:13:19 PM PST
-  Document approved by Will Hail (whail@pathwaymedtech.com)
Approval Date: 2025-02-11 - 12:13:19 PM PST - Time Source: server
-  Agreement completed.
2025-02-11 - 12:13:19 PM PST

NCR-25-002 -

Final Audit Report

2025-02-11

Created:	2025-02-11 (Pacific Standard Time)
By:	Larry Dacoron (ldacoron@pathwaympi.com)
Status:	Approved
Transaction ID:	CBJCHBCAABAAsPCSLWivxzDZ7GsDoxP-p7wgZzZXJTrq
Number of Documents:	1
Document page count:	8
Number of supporting files:	0
Supporting files page count:	0

"NCR-25-002 -" History

-  Document created by Larry Dacoron (ldacoron@pathwaympi.com)
2025-02-11 - 3:45:43 PM PST
-  Document emailed to Aaron Rogers (arogers@pathwaymedtech.com) for approval
2025-02-11 - 3:46:15 PM PST
-  Email viewed by Aaron Rogers (arogers@pathwaymedtech.com)
2025-02-11 - 3:48:01 PM PST
-  Aaron Rogers (arogers@pathwaymedtech.com) authenticated with Email OTP by verifying one-time code sent via email
2025-02-11 - 3:48:14 PM PST
-  Agreement viewed by Aaron Rogers (arogers@pathwaymedtech.com)
2025-02-11 - 3:48:15 PM PST
-  Aaron Rogers (arogers@pathwaymedtech.com) has agreed to the terms of use and to do business electronically with Pathway Medtech, LLC
2025-02-11 - 3:48:34 PM PST
-  Document approved by Aaron Rogers (arogers@pathwaymedtech.com)
Approval Date: 2025-02-11 - 3:48:34 PM PST - Time Source: server
-  Agreement completed.
2025-02-11 - 3:48:34 PM PST