

	Document Title: Suspension Processing Traveler	Revision: BA
	Document #: TR-C.SLQ001	Effective Date: 23 Sep 2022

General Production Information

Work Order #:	2502-01	NCMR #:	NCMR-001
Part Number:	C.SLQ001	Part Revision:	B
Lot #:	2502-01-1	Manufacture Date / Expiration Date:	Mfg date: 11 Feb 2026 Exp date: 11 Feb 2027
Build Purpose:	Commercial Product		
Additional Build Info:	N/A		

Material Traceability

Part Number	Description	Part Revision	Lot / Batch #	Expiration Date	Qty. Obtained	Obtained By: (Initial/Date)
M.SLQ001	Polymeric Dispersion Solution, Unpurified, CHEMPSBTC1010	B	MY 27B-FO	05 Feb 2026	13Kg	GA 11 Feb 2025

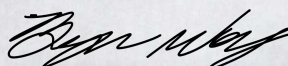
Equipment Traceability

Part Number	Description	Equipment ID	Calibration Due Date
E.SLQ006	Diafiltration System, KrosFlo KTF-2000	ST-005	01 Nov 2026
E.SLQ010	Floor Scale, Concentrate Reactor	ST-006	11 Dec 2025
E.SLQ012	Floor Scale, Suspension	ST-007	11 Dec 2025

Build Record

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
1	QM.SLQ049	Line clearance completed	☒ Yes ☐ No ☐ N/A	GA 11 Feb 2025
2	MP-C.SLQ001	Verify conductivity indicator on Portable Exchange Deionizer System (E.SLQ001) is green.	☒ Indicator is green	GA 11 Feb 2025

Reviewed By (Signature and Date):



29 Apr 2025

Batch/Lot No.: 2502-01-1

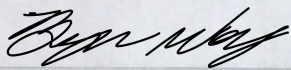
Silq Technologies Corporation

PROPRIETARY AND CONFIDENTIAL
This document is for reference only

Page 1 of 12

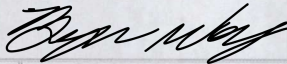
	Document Title: Suspension Processing Traveler	Revision: BA
	Document #: TR-C.SLQ001	Effective Date: 23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
32	MP-C.SLQ001	Drain liquid in diafiltration system/tubing into waste tank.	☒ Yes ☐ No	11 Feb 2025 CA
43	MP-C.SLQ001	Connect Retantate Retentate Lines and Feed Lines to the appropriate pumps. Route lines through the appropriate membrane modules, flow meters, and backpressure valves.	☒ Yes ☐ No	11 Feb 2025 MA
54	MP-C.SLQ001	Place Permeate Line into Peristaltic Pump; change Valve 1 to 'permeate line' position.	☒ Yes ☐ No	11 Feb 2025 CA
65	MP-C.SLQ001	Tare Concentrate Reactor Floor Scale. Fill Concentrate Reactor with ~30kg of high purity water (S.SLQ001).	☒ Floor scale tared Floor scale reading: 30.0 kg	11 Feb 2025 YR
76	MP-C.SLQ001	'Feed (remaining)' value on system software is within ± 0.2 kg of floor scale reading.	'Feed (Remaining)' value: 30.0 kg ☒ Within ± 0.2 kg of floor scale reading	11 Feb 2025 CA
87	MP-C.SLQ001	Manually run Pump 01 and Pump 02 until membrane modules and retentate lines are filled.	☒ Yes ☐ No	11 Feb 2025 CA
98	MP-C.SLQ001	Confirm no leaks are observed in all connections.	☒ No leaks observed	11 Feb 2025 CA
109	MP-C.SLQ001	Fill Concentrate Reactor until scale reads ~30kg.	Floor scale reading: 30.0 kg	11 Feb 2025
110	MP-C.SLQ001	Set parameters on main page of system software. • Starting Volume: 30L • DV Setpoint: 2.0 DV	☒ Parameters set	CA
121	MP-C.SLQ001	Click START PROCESS button to fill membrane modules with water. After membrane modules are filled with water, ensure Permeate Valve is set to vertical position and open the four clamps on permeate line.	☒ Membrane modules filled with water ☒ Permeate Valve set to vertical position ☒ Open four clamps on permeate line	
132	MP-C.SLQ001	Click Stop button when 'Perm. Total' value reads 60L on system software.	'Perm. Total' value: 60.4 kg	11 Feb 2025 YR

Reviewed By (Signature and Date):  29 Apr 2025	Batch/Lot No.: 2502-01-1
Silq Technologies Corporation	PROPRIETARY AND CONFIDENTIAL This document is for reference only
Page 2 of 12	

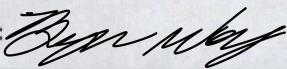
	Document Title: Suspension Processing Traveler	Revision: BA
	Document #: TR-C.SLQ001	Effective Date: 23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
143	MP-C.SLQ001	Mount Pneumatic Drum Stirrer into holder and turn on.	☒ Yes ☐ No	MR 11 Feb 2025
154	MP-C.SLQ001	Obtain container(s) of Unpurified Polymeric Dispersion Solution (M.SLQ001). Record traceability information in the Material Traceability table above.	☒ Container(s) obtained ☒ Material traceability recorded	
165	MP-C.SLQ001	Record volume of container(s).	Volume of container(s): 13 kg	
176	MP-C.SLQ001	Fill Concentrate Reactor with high purity water until the Designated Fill Value (DFV) equals the volume of container(s).	☒ Yes ☐ No Designated Fill Value (DFV): 33.1 L	MR 11 Feb 2025
187	MP-C.SLQ001	Calculate the Designated Operating Value (DOV). $DOV = 10 / 3 * DFV$.	Designated Operating Value (DOV): 110.4 L	MR 11 Feb 2025
198	MP-C.SLQ001	Slowly pour contents of Unpurified Polymeric Dispersion Solution container(s) into concentrate reactor. Ensure no gelation develops.	☒ Yes ☐ No	MR 11 Feb 2025
2019	MP-C.SLQ001	Fill concentrate reactor with high purity water until DOV is reached (within ± 2 kg). Record floor scale and system software 'Feed (Remaining)' value. Verify values are within ± 2 kg of DOV and each other.	Floor scale reading: 110.6 kg "Feed (Remaining)" value: 110.6 kg ☒ Values are within ± 2 kg of DOV and each other	
210	MP-C.SLQ001	Set parameters on main page of system software. - Starting Volume: Designated Operating Value - DV Setpoint: 10 DV	☒ Parameters set	
224	MP-C.SLQ001	Click START PROCESS button. Verify Pump 01 and Pump 02 are running	☒ Pump 01 and Pump 02 running	
232	MP-C.SLQ001	Run process until "Perm. Total" value reaches 10 times the value of the Designated Fill Value (see Step 17 above).	"Perm. Total" value: 1106.1 L	

Reviewed By (Signature and Date):  29 Apr 2025	Batch/Lot No.: 2502-01-1
Silq Technologies Corporation	PROPRIETARY AND CONFIDENTIAL This document is for reference only
Page 3 of 12	

	Document Title: Suspension Processing Traveler	Revision: BA
	Document #: TR-C.SLQ001	Effective Date: 23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
243	MP-C.SLQ001	Record weight of product tank (C.SLQ003) without the lid.	Product tank weight: 8.2 kg	11 Feb 2025 CR
254	MP-C.SLQ001	Transfer Suspension to plastic drum using the peristaltic pump via the product line (C.SLQ003).	☒ Yes ☐ No	11 Feb 2025 CR
265	MP-C.SLQ001	Take sample of suspension in 50mL tubes at the beginning, middle, and end of the transfer. Label each tube Aliquot A, Aliquot B, and Aliquot C, respectively.	☒ Yes ☐ No	11 Feb 2025 CR
276	MP-C.SLQ001	Place clean containers below Pump 01 and Pump 02. Disconnect retentate lines from pumps and allow suspension to drain into containers. Pour contents of the containers into the product tank.	☒ Yes ☐ No	11 Feb 2025 CR
287	MP-C.SLQ001	Record the initial weight of the suspension. Assumption: The initial weight (in kg) of the suspension is equal to the <u>weight of product tank subtracted from the weight of the suspension and product tank, initial volume (in L).</u>	Initial Suspension Weight: 131.1 kg	11 Feb 2025 CR
298	MP-C.SLQ001	Seal product tank.	☒ Yes ☐ No	11 Feb 2025 CR
3029	MP-C.SLQ001	Create product label for Suspension (L.SLQ005) indicating lot number and expiration date of suspension. Apply label to product tank.	☒ Yes ☐ No	
3130	MP-C.SLQ001	Place product tank in quarantine pending QC analysis of suspension.	☒ Yes ☐ No	
321	MP-C.SLQ001	Prepare KTF-2000 System for storage Fill concentrate reactor with ~30kg of high purity water and manually start running Pump 01 and Pump 02. Allow Membrane Modules 1 and 2 and Retentate Lines 1A, 1B, 2A, and 2B are filled. Allow pumps to run for an additional 10 minutes.	☒ Yes ☐ No	

Reviewed By (Signature and Date):  29 Apr 2025		Batch/Lot No.: 2502-01-1
Silq Technologies Corporation	PROPRIETARY AND CONFIDENTIAL This document is for reference only	Page 4 of 12



Document Title:
Suspension Processing Traveler

Revision:
BA

Document #:
TR-C.SLQ001

Effective Date:
23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
332	MP-C.SLQ001	Switch Permeate Line Valve to horizontal (closed) position and Waste Line Valve to vertical (open) position. Place Waste Line through peristaltic pump and into the concentrate reactor. Manually run the peristaltic pump to fully drain the contents of the concentrate reactor into the waste tank.	☒ Yes ☐ No	11 Feb 2025 LM
343	MP-C.SLQ001	Create Storage Solution. Combine 20L of high purity water + 1.5mL of bleach solution (S.SLQ007).	☒ Yes ☐ No	
354	MP-C.SLQ001	Pour Storage Solution into Concentrate Reactor. Close 4 Permeate Line Clamps and switch Permeate Line Valve to horizontal (closed) position and Waste Line Valve to vertical (open) position.	☒ Yes ☐ No	
365	MP-C.SLQ001	Manually run Pump 01 and Pump 02 until Membrane Modules 1 and 2 and Retentate Lines 1A, 1B, 2A, and 2B are filled.	☒ Yes ☐ No	
376	MP-C.SLQ001	Once filled, stop Pump 01 and Pump 02 and close the ball valves for Retentate Lines 1B and 2B at the Concentrate Reactor. Allow 5 minutes to elapse for settling and degassing. Ensure the Storage Solution has filled Membrane Modules 1 and 2 entirely. Storage Solution should at least fill Retentate Lines 1B and 2B up to the respective line's Flow Meter.	☒ Yes ☐ No	
387	MP-C.SLQ001	Wait 5 minutes to ensure the Membrane Modules are full and stable.	☒ Yes ☐ No	
398	MP-C.SLQ001	Manually run the Peristaltic Pump to fully drain the contents of the Concentrate Reactor into the Waste Tank.	☒ Yes ☐ No	
4039	MP-C.SLQ001	Power down Pump 03, Pump 04, concentrate reactor scale, suspension scale, and KTF system software.	☒ Yes ☐ No	
410	MP-C.SLQ001	Prepare the high purity water tank for storage Take a peristaltic pump with sterilized tubing to pump out any residual high purity water from the tank.	☒ Yes ☐ No	
421	MP-C.SLQ001	Allow high purity tank to dry for 24 hours before securing lid to the tank.	☒ Yes ☐ No	

Reviewed By (Signature and Date): *Byron M. King* 29 Apr 2025

Batch/Lot No.: 2502-01-1



Document Title:
Suspension Processing Traveler

Revision:
BA

Document #:
TR-C.SLQ001

Effective Date:
23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
432	MP-C.SLQ001	Clean containers using in the manufacturing process. Allow containers to dry and store after they have dried.	☒ Yes ☐ No	
443	QC-C.SLQ001	Execution of QC analysis of suspension samples. QC analysis report attached to DHR.	☒ Yes ☐ No	11 Feb 2025 UK
454	MP-C.SLQ001	Attach Suspension Quality Control Report and COA to DHR.	☒ Yes ☐ No	07 Mar 2025 UK
465	MP-C.SLQ001	Test results meet acceptance criteria.	☒ Yes ☐ No	07 Mar 2025 UK
476	MP-C.SLQ001	Dilution Procedure Calculate the Average Suspension Sample Concentration (g/L) by taking the average 'Residue Net Weight' of six (6) freeze dried suspension samples based off the 10 mL sample.	Average Suspension Sample Concentration: <u>8.55</u> g/L	07 Mar 2025 UK
48	MP-C.SLQ001	<u>Determine if Concentration Process is required. If Average Suspension Sample Concentration is <10g/L, Concentrate Process is required. If Average Suspension Sample Concentration is >10g/L, Dilution Procedure can proceed.</u>	<u>Is Concentration Process Required?</u> ☒ Yes, skip Steps 49-54 and proceed to Step A1. ☐ No, proceed to Step 49	07 Mar 2025 UK
497	MP-C.SLQ001	Calculate the volume (L) of high purity water to be added to the product tank. Volume = [Average Suspension Sample Concentration (g/L) * Initial Suspension Volume (L)]/[10 (g/L)] - Initial Suspension Volume (L)	Volume of high purity water to be added to product tank: _____ L	
5048	MP-C.SLQ001	Tare suspension floor scale and place product tank onto floor scale.	☐ Yes ☐ No	
5149	MP-C.SLQ001	Add volume of high purity water (from Step 35) to the product tank.	☐ Yes ☐ No	
5250	MP-C.SLQ001	Record the weight of product tank <u>with suspension</u> .	_____ kg	
531	MP-C.SLQ001	Calculate weight of suspension.	_____ kg	
542	MP-C.SLQ001	Update product label with weight of suspension.	☐ Yes ☐ No	

Reviewed By (Signature and Date):

[Signature]

29 Apr 2025

Batch/Lot No.: 2502-01-1

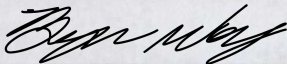
	Document Title: Suspension Processing Traveler	Revision: BA
	Document #: TR-C.SLQ001	Effective Date: 23Sep2022

Comments/Additional Notes:

Build Record Suspension Concentration Process

(Check one: ☐ Required, ☐ Not Required)

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
A1	QM.SLQ049	Line clearance completed	☒ Yes ☐ No ☐ N/A	JH MAR 2025 GF
A2	MP-C.SLQ001	Drain liquid in diafiltration system/tubing into waste tank.	☒ Yes ☐ No	JH MAR 2025 GF
A3	MP-C.SLQ001	Connect Retentate Lines and Feed Lines to the appropriate pumps. Route lines through the appropriate membrane modules, flow meters, and backpressure valves.	☒ Yes ☐ No	JH MAR 2025 GF
A4	MP-C.SLQ001	Place Permeate Line into Peristaltic Pump; change Valve 1 to 'permeate line' position.	☒ Yes ☐ No	JH MAR 2025 GF
A5	MP-C.SLQ001	Tare Concentrate Reactor Floor Scale. Fill Concentrate Reactor with ~30kg of high purity water (S.SLQ001).	☒ Floor scale tared Floor scale reading: <u>30.6 kg</u>	JH MAR 2025 GF
A6	MP-C.SLQ001	'Feed (remaining)' value on system software is within ±0.2kg of floor scale reading.	'Feed (Remaining)' value: <u>30.6 kg</u> ☒ Within ±0.2kg of floor scale reading	JH MAR 2025 GF
A7	MP-C.SLQ001	Manually run Pump 01 and Pump 02 until membrane modules and retentate lines are filled.	☒ Yes ☐ No	JH MAR 2025 GF
A8	MP-C.SLQ001	Confirm no leaks are observed in all connections.	☒ No leaks observed	JH MAR 2025 GF
A9	MP-C.SLQ001	Fill Concentrate Reactor until scale reads ~30kg.	Floor scale reading: <u>30.6 kg</u>	JH MAR 2025 GF
A10	MP-C.SLQ001	Set parameters on main page of system software. • Starting Volume: 30L • DV Setpoint: 2.0 DV	☒ Parameters set	
A11	MP-C.SLQ001	Click START PROCESS button to fill membrane modules with water. After membrane modules are filled with	☒ Membrane modules filled with water	

Reviewed By (Signature and Date):  29 Apr 2025

Batch/Lot No.: 2502-01-1



Document Title:
Suspension Processing Traveler

Revision:
BA

Document #:
TR-C.SLQ001

Effective Date:
23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
		water, ensure Permeate Valve is set to vertical position and open the four clamps on permeate line.	☒ Permeate Valve set to vertical position ☒ Open four clamps on permeate line	
A12	MP-C.SLQ001	Click Stop button when 'Perm. Total' value reads 60L on system software.	'Perm. Total' value: 60.0	
A13	MP-C.SLQ001	Mount Pneumatic Drum Stirrer into holder and turn on.	☒ Yes ☐ No	
A14	MP-C.SLQ001	Drain liquid in diafiltration system/tubing into waste tank.	☒ Yes ☐ No	
A15	MP-C.SLQ001	Calculate Concentration Permeate Required Value: Concentration Permeate Required = $\frac{[12 \text{ (g/L)} - \text{Average Suspension Sample Concentration (g/L)}] \times 12 \text{ (g/L)} \times \text{Initial Suspension Volume (L)}}{16 \text{ (L)}}$	Concentration Permeate Required (L): 53.7 L	24 MAR 2025 G
A16	MP-C.SLQ001	Tare Concentrate Reactor Floor Scale.	☒ Floor scale tared	
A17	MP-C.SLQ001	Transfer Suspension from plastic drum using the peristaltic pump via the product line (C.SLQ003) into Concentrate Reactor	☒ Yes ☐ No	24 MAR 2025 G
A18	MP-C.SLQ001	Manually run Pump 01 and Pump 02 until "Perm. Total" value reaches Concentration Permeate Required Value. Note: No High Purity Water will be added to Concentrate Reactor and concentrate Reactor Floor Scale reading will gradually decrease	"Perm. Total" value: 53.8 L	24 MAR 2025 G
A19	MP-C.SLQ001	Tare suspension floor scale and place product tank onto floor scale.	☒ Yes ☐ No	
A20	MP-C.SLQ001	Record weight of product tank (C.SLQ003) without the lid.	Product tank weight: 8.0 kg	
A21	MP-C.SLQ001	Transfer Suspension to plastic drum using the peristaltic pump via the product line (C.SLQ003).	☒ Yes ☐ No	24 MAR 2025 G
A22	MP-C.SLQ001	Take sample of suspension in 50mL tubes at the beginning, middle, and end of the transfer. Label each tube Aliquot A, Aliquot B, and Aliquot C, respectively.	☒ Yes ☐ No	

Reviewed By (Signature and Date): *[Signature]* 29 Apr 2025

Batch/Lot No.: 2502-01-1



Document Title:
Suspension Processing Traveler

Revision:
BA

Document #:
TR-C.SLQ001

Effective Date:
23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
A23	MP-C.SLQ001	Place clean containers below Pump 01 and Pump 02. Disconnect retentate lines from pumps and allow suspension to drain into containers. Pour contents of the containers into the product tank.	☒ Yes ☐ No	JH MWR 2025 G
A24	MP-C.SLQ001	Record the initial weight of the suspension. Assumption: The initial weight (in kg) of the suspension is equal to the weight of product tank subtracted from the weight of the suspension and product tank	Initial Suspension Weight: 63.2 kg	JH MWR 2025 G
A25	MP-C.SLQ001	Seal product tank.	☒ Yes ☐ No	JH MWR 2025 G
A26	MP-C.SLQ001	Create product label for Suspension (L.SLQ005) indicating lot number and expiration date of suspension. Apply label to product tank.	☒ Yes ☐ No	JH MWR 2025 G
A27	MP-C.SLQ001	Place product tank in quarantine pending QC analysis of suspension.	☒ Yes ☐ No	JH MWR 2025 G
A28	MP-C.SLQ001	<u>Prepare KTF-2000 System for storage</u> Fill concentrate reactor with ~30kg of high purity water and manually start running Pump 01 and Pump 02. Allow Membrane Modules 1 and 2 and Retentate Lines 1A, 1B, 2A, and 2B are filled. Allow pumps to run for an additional 10 minutes.	☒ Yes ☐ No	JH MWR 2025 G
A29	MP-C.SLQ001	Switch Permeate Line Valve to horizontal (closed) position and Waste Line Valve to vertical (open) position. Place Waste Line through peristaltic pump and into the concentrate reactor. Manually run the peristaltic pump to fully drain the contents of the concentrate reactor into the waste tank.	☒ Yes ☐ No	JH MWR 2025 G
A30	MP-C.SLQ001	Create Storage Solution. Combine 20L of high purity water + 1.5mL of bleach solution (S.SLQ007).	☒ Yes ☐ No	JH MWR 2025 G
A31	MP-C.SLQ001	Pour Storage Solution into Concentrate Reactor. Close 4 Permeate Line Clamps and switch Permeate Line Valve to horizontal (closed) position and Waste Line Valve to vertical (open) position.	☒ Yes ☐ No	JH MWR 2025 G

Reviewed By (Signature and Date):

[Signature]

29 Apr 2025

Batch/Lot No.: 2502-01-1



Document Title:
Suspension Processing Traveler

Revision:
BA

Document #:
TR-C.SLQ001

Effective Date:
23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
A32	MP-C.SLQ001	Manually run Pump 01 and Pump 02 until Membrane Modules 1 and 2 and Retentate Lines 1A, 1B, 2A, and 2B are filled.	☒ Yes ☐ No	JH AVR 2025 G
A33	MP-C.SLQ001	Once filled, stop Pump 01 and Pump 02 and close the ball valves for Retentate Lines 1B and 2B at the Concentrate Reactor. Allow 5 minutes to elapse for settling and degassing. Ensure the Storage Solution has filled Membrane Modules 1 and 2 entirely. Storage Solution should at least fill Retentate Lines 1B and 2B up to the respective line's Flow Meter.	☒ Yes ☐ No	JH AVR 2025 G
A34	MP-C.SLQ001	Wait 5 minutes to ensure the Membrane Modules are full and stable.	☒ Yes ☐ No	JH AVR 2025 G
A35	MP-C.SLQ001	Manually run the Peristaltic Pump to fully drain the contents of the Concentrate Reactor into the Waste Tank.	☒ Yes ☐ No	JH AVR 2025 G
A36	MP-C.SLQ001	Power down Pump 03, Pump 04, concentrate reactor scale, suspension scale, and KTF system software.	☒ Yes ☐ No	JH AVR 2025 G
A37	MP-C.SLQ001	<u>Prepare the high purity water tank for storage</u> Take a peristaltic pump with sterilized tubing to pump out any residual high purity water from the tank.	☒ Yes ☐ No	JH AVR 2025 G
A38	MP-C.SLQ001	Allow high purity tank to dry for 24 hours before securing lid to the tank.	☒ Yes ☐ No	JH AVR 2025 G
A39	MP-C.SLQ001	Clean containers using in the manufacturing process. Allow containers to dry and store after they have dried.	☒ Yes ☐ No	JH AVR 2025 G
A40	QC-C.SLQ001	Execution of QC analysis of suspension samples. QC analysis report attached to DHR.	☒ Yes ☐ No	JH AVR 2025 G
A41	MP-C.SLQ001	Attach Suspension Quality Control Report and COA to DHR.	☒ Yes ☐ No	JH AVR 2025 G
A42	MP-C.SLQ001	Test results meet acceptance criteria.	☒ Yes ☐ No	CPA 01 Apr 2025
A43	MP-C.SLQ001	<u>Dilution Procedure</u> Calculate the Average Suspension Sample Concentration (g/L) by taking the average 'Residue Net Weight' of six (6) freeze	Average Suspension Sample Concentration: <u>13.88</u> g/L	CR 01 Apr 2025

Reviewed By (Signature and Date): *[Signature]* 29 Apr 2025

Batch/Lot No.: 2502-01-1



Document Title:
Suspension Processing Traveler

Revision:
BA

Document #:
TR-C.SLQ001

Effective Date:
23Sep2022

Step	Procedure	Step Description	Status / Results	Performed By: (Initial/Date)
		dried suspension samples based off the 10 mL sample.		
A44	MP-C.SLQ001	Determine if Concentration Process is required. If Average Suspension Sample Concentration is <10g/L, Concentrate Process is required again. If Average Suspension Sample Concentration is >10g/L, Dilution Procedure can proceed.	Is Concentration Process Required? ☐ Yes, skip Steps A45-A50 and proceed to Step A1 and repeat. ☒ No, proceed to Step A45.	01 APR 2025 ZK
A45	MP-C.SLQ001	Calculate the volume (L) of high purity water to be added to the product tank. Volume = [Average Suspension Sample Concentration (g/L) * Initial Suspension Volume (L)]/[10 (g/L)] - Initial Suspension Volume (L)	Volume of high purity water to be added to product tank: <u>18.2</u> L	01 APR 2025 MK
A46	MP-C.SLQ001	Tare suspension floor scale and place product tank onto floor scale.	☒ Yes ☐ No	01 APR 2025 Q
A47	MP-C.SLQ001	Add volume of high purity water (from Step 35) to the product tank.	☒ Yes ☐ No	01 APR 2025 Q
A48	MP-C.SLQ001	Record the weight of product tank with suspension	<u>89.4</u> kg	01 APR 2025 Q
A49	MP-C.SLQ001	Calculate weight of suspension.	<u>81.2</u> kg	01 APR 2025 Q
A50	MP-C.SLQ001	Update product label with weight of suspension.	☒ Yes ☐ No	01 APR 2025 Q

End of procedure

Comments/Additional Notes:

Reviewed By (Signature and Date): Byron M. King 29 Apr 2025

Batch/Lot No.: 2502-01-1

	Document Title: Suspension Processing Traveler	Revision: BA
	Document #: TR-C.SLQ001	Effective Date: 23Sep2022

Part Description: Suspension

Part Number: C.SLQ001

Revision: B

Amount: 81.2kg

Lot Number: 2502-01-1

Expiration 13FEB2027



Silq Technologies Corp.
311 Otterson Dr., Suite 60
Chico, CA 95928



Reviewed By (Signature and Date):  29 Apr 2025		Batch/Lot No.: 2502-01-1
Silq Technologies Corporation	PROPRIETARY AND CONFIDENTIAL This document is for reference only	Page 12 of 12