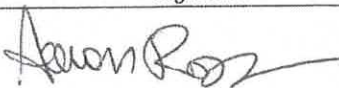
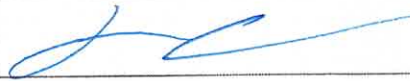


	Design Review		Doc No.: FM-0052
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.		Rev.: B
			Page 1 of 2

Project No.	Product Name	Part No./ Rev	Review Phase	Review Date
00017	Hydrophilix 100% Silicone Foley Catheter	HDX-8000 / B HDX-8001 / B HDX-8002 / B	Phase 3 Design Outputs	08/20/2020
Subject of Review				
#	Topic	Details	Reference Documents	
1	Action Items from Last Review (07/25/2019)	No action items	FILE-004280 rev 07/25/2019	
2	Component Specification	The documents were reviewed. The requirement is complete.	HDX-3000 rev A, HDX-4000 rev B, HDX-4001 rev B, HDX-5000 rev A, HDX-5001 rev A, HDX-5002 rev A, HDX-5003 rev A, HDX-5004 rev B, HDX-5005 rev A, HDX-5006 rev A, HDX-6000 rev A, HDX-6001 rev A, HDX-6002 rev A	
3	Preliminary BOM / Assembly Specification	The documents were reviewed. The requirement is complete.	HDX-8000 rev B, HDX-8001 rev B, HDX-8002 rev B	
4	Project Plan	No changes required to the project plan. The plan is current.	PL-0013 rev A	
5	Risk	No new risks discovered, nor changes required to the RMP or HA. The risk management plan is current. With no new risks identified the HA is complete.	RM-0018 rev A, RM-0019 rev A	
6	Review of Phase 3 Deliverables	Each document was reviewed and discussed. Project team in attendance has concluded that the deliverables are complete and satisfy the Design Control requirements in accordance with QP-007 rev E. It is agreed upon to freeze the design, close PHASE 3 and progress into PHASE 4.	Refer to #1-#2, QP-007 rev E.	
Action Items				
#	Description		Assigned To	Due Date
N/A	N/A		N/A	N/A

	Design Review	Doc No.: FM-0052
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>	Rev.: B
		Page 2 of 2

Project No.	Product Name	Part No./ Rev	Review Phase	Review Date
00017	Hydrophilix 100% Silicone Foley Catheter	HDX-8000 / B HDX-8001 / B HDX-8002 / B	Phase 3 Design Outputs	08/20/2020

Reviewers					
Functional Group	Required?	Name	Title	Signature	Date
Quality Assurance / Regulatory Affairs	☒	Aaron Rogers	Director of Regulatory and Quality		08/25/2020
R&D/ Engineering / Manufacturing	☒	Jonas Cochran	Managing Partner		08/25/2020
Marketing	☐				
Independent Reviewer	☒	Tom Cassou	Director of Business Development		8/25/2020
Other / Customer	☒	Brian McVerry	Chief Technology Officer		8.25.20

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	
Page 1 of 9		

Material Requirements				☐ N/A	
Product Number: Hydrophilix PolymerSoluti					
Dimensional Requirements				☒ N/A	
Visual Requirements				☒ N/A	
Label Content Specification				☒ N/A	
Certificate of Compliance (C of C)		☒ YES ☐ NO ☐ N/A	Certificate of Analysis (C of A)		☐ YES ☐ NO ☒ N/A
Material Compliance					
REACH Compliant		☐ YES ☐ NO ☒ N/A	ROHS Compliant		☐ YES ☐ NO ☒ N/A
Conflict Material		☐ YES ☐ NO ☒ N/A	FDA (Medical Grade) Compliant		☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)		☒ YES ☐ NO ☐ N/A	BPA Free		☐ YES ☐ NO ☒ N/A
Latex Free		☐ YES ☐ NO ☒ N/A	DEHP Free		☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free		☐ YES ☐ NO ☒ N/A	Documented evidence for compliance to incorporated material specification file		
Storage Requirements				☐ N/A	
Protect from heat and direct sunlight. Store in a cool location.					
Shelf-Life Requirements				☒ N/A	
Packaging Requirements				☒ N/A	
Supplier will notified Pathway Innovation in advance of any changes prior to the implementation of the change in the supplied material.					

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 2 of 9



Issue date 02/20/2019

Safety Data Sheet (SDS)

OSHA HazCom Standard 29 CFR 1910.1200(g) and GHS Rev 03.

Page 1/7

Reviewed on 02/20/2019

1 Identification

- Product Identifier
- Trade Name: Hydrophilix Polymer Solution
- Product Number: Hydrophilix PolymerSoluti
- Relevant identified uses of the substance or mixture and uses advised against:
- Product Description: Opaque Liquid
- Details of the Supplier of the Safety Data Sheet:
- Manufacturer/Supplier:
Hydrophilix, Inc.
12100 Wilshire Blvd., Suite 800
Los Angeles, CA 90025
ph: (310) 806-9202
email: Brian@Hydrophilix.com
website: Hydrophilix.com
- Emergency telephone number: (949) 735-2614

2 Hazard(s) Identification

- Classification of the substance or mixture:
The product does not need classification according to OSHA HazCom Standard 29 CFR paragraph (d) of §1910.1200(g) and GHS Rev 03.
- Label elements:
 - GHS label elements Non-Regulated Material
 - Hazard pictograms: Non-Regulated Material
 - Signal word: Non-Regulated Material
 - Hazard statements: Non-Regulated Material
 - Unknown acute toxicity:
This value refers to knowledge of known, established toxicological or ecotoxicological values.
0 % of the mixture consists of component(s) of unknown toxicity.
 - Classification system: NFPA/HMIS Definitions: 0-Least, 1-Slight, 2-Moderate, 3-High, 4-Extreme
 - NFPA ratings (scale 0 - 4)



Health = 0
Fire = 0
Reactivity = 0

HMIS-ratings (scale 0 - 4)



Health = 0
Fire = 0
Physical Hazard = 0

- Hazard(s) not otherwise classified (HNOC): None known

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 3 of 9



Issue date 02/20/2019

Safety Data Sheet (SDS)

OSHA HazCom Standard 29 CFR 1910.1200(g) and GHS Rev 03.

Page 2/7

Reviewed on 02/20/2019

Trade Name: Hydrophilix Polymer Solution

3 Composition/Information on Ingredients

- Chemical characterization: Mixtures
- Description: Non-Regulated Material
- Dangerous Components: Non-Regulated Material
- Additional information:

The exact percentages of the ingredients of this mixture are considered to be proprietary and are withheld in accordance with the provisions of paragraph (i) of §1910.1200 of 29 CFR 1910.1200 Trade Secrets.

4 First-Aid Measures

- Description of first aid measures
- General information: No special measures required.
- After inhalation: Supply fresh air; consult doctor in case of complaints.
- After skin contact: If skin irritation occurs, consult a doctor.
- After eye contact:

If eye irritation occurs, consult a doctor.

Rinse opened eye for several minutes under running water.

- After swallowing: If swallowed and symptoms occur, consult a doctor.
- Information for doctor
- Most important symptoms and effects, both acute and delayed: No further relevant information available.
- Indication of any immediate medical attention and special treatment needed:

No further relevant information available.

5 Fire-Fighting Measures

- Extinguishing media
- Suitable extinguishing agents: Use fire fighting measures that suit the environment.
- For safety reasons unsuitable extinguishing agents: No further relevant information.
- Special hazards arising from the substance or mixture: No further relevant information available.
- Advice for firefighters
- Special protective equipment for firefighters:

As in any fire, wear self-contained breathing apparatus pressure-demand (NIOSH approved or equivalent) and full protective gear to prevent contact with skin and eyes.

6 Accidental Release Measures

- Personal precautions, protective equipment and emergency procedures: Not required.
- Environmental precautions: Dilute with plenty of water.
- Methods and material for containment and cleaning up:

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 4 of 9



Issue date 02/20/2019

Safety Data Sheet (SDS)

OSHA HazCom Standard 29 CFR 1910.1200(g) and GHS Rev 03.

Page 3/7

Reviewed on 02/20/2019

Trade Name: Hydrophilix Polymer Solution

Absorb with liquid-binding material (i.e. sand, diatomite, acid binders, universal binders, sawdust).

Dispose of the collected material according to regulations.

· Reference to other sections:

See Section 7 for information on safe handling.

See Section 8 for information on personal protection equipment.

See Section 13 for disposal information.

· Protective Action Criteria for Chemicals

· **PAC-1:**

None of the ingredients are listed.

· **PAC-2:**

None of the ingredients are listed.

· **PAC-3:**

None of the ingredients are listed.

7 Handling and Storage

· Handling

· Precautions for safe handling: No special measures required.

· Information about protection against explosions and fires: No special measures required.

· Conditions for safe storage, including any incompatibilities

· Storage

· Requirements to be met by storerooms and receptacles:

Protect from heat and direct sunlight.

Store in a cool location.

· Information about storage in one common storage facility: Not required.

· Further information about storage conditions: None.

· Specific end use(s): No further relevant information available.

8 Exposure Controls/Personal Protection

· Additional information about design of technical systems: No further data; see section 7.

· Control parameters:

· Components with occupational exposure limits:

The product does not contain any relevant quantities of materials with critical values that have to be monitored at the workplace.

· Additional information: The lists that were valid during the creation of this SDS were used as basis.

· Exposure controls:

· Personal protective equipment

· General protective and hygienic measures:

The usual precautionary measures for handling chemicals should be followed.

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work	Page 5 of 9



Issue date 02/20/2019

Safety Data Sheet (SDS)

OSHA HazCom Standard 29 CFR 1910.1200(g) and GHS Rev 03.

Page 4/7

Reviewed on 02/20/2019

Trade Name: Hydrophilix Polymer Solution

- Breathing equipment: Not required.
- Protection of hands: Not required.
- Material of gloves: Not applicable.
- Penetration time of glove material: Not applicable.
- Eye protection:



Goggles recommended during refilling.

- Limitation and supervision of exposure into the environment: None

9 Physical and Chemical Properties

- Information on basic physical and chemical properties
- General Information
- Appearance:
 - Form: Liquid
 - Color: Opaque (white)
- Odour: Odorless
- Odor threshold: Not determined.
- pH-value @ 20 °C (68 °F): 4-6
- Change in condition
 - Melting point/Melting range: Not determined.
 - Boiling point/Boiling range: ≥ 100 °C (≥ 212 °F)
- Flash point: None
- Flammability (solid, gaseous): Not applicable.
- Ignition temperature: Not applicable
- Decomposition temperature: Not determined.
- Auto igniting: Product is not self-igniting.
- Danger of explosion: Product does not present an explosion hazard.
- Explosion limits:
 - Lower: Not determined.
 - Upper: Not determined.
- Vapor pressure @ 20 °C (68 °F): ≤ 23 hPa (≤ 17.3 mm Hg)
- Density @ 20 °C (68 °F): 1 g/cm³ (8.345 lbs/gal)
- Relative density: Not determined.
- Vapor density: Not determined.
- Evaporation rate: Not determined.

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 6 of 9



Issue date 02/20/2019

Safety Data Sheet (SDS)

OSHA HazCom Standard 29 CFR 1910.1200(g) and GHS Rev 03.

Page 5/7

Reviewed on 02/20/2019

Trade Name: Hydrophilix Polymer Solution

- Solubility in / Miscibility with:
Water: Fully miscible.
- Partition coefficient (n-octanol/water): Not determined.
- Viscosity:
Dynamic: Not determined.
Kinematic: Not determined.
- Solvent content:
Water: 99.0 %
VOC content: 0.00 %
0.0 g/l / 0.00 lb/gal
Solids content: 0.0 %
- Other information: No further relevant information available.

10 Stability and Reactivity

- Reactivity: The product is stable under normal conditions.
- Chemical stability: Stable under normal conditions.
- Thermal decomposition / conditions to be avoided: No decomposition if used according to specifications.
- Possibility of hazardous reactions: No dangerous reactions known.
- Conditions to avoid: No further relevant information available.
- Incompatible materials: No further relevant information available.
- Hazardous decomposition products: No dangerous decomposition products known.

11 Toxicological Information

- Information on toxicological effects:
- Acute toxicity:
- LD/LC50 values that are relevant for classification: No data available.
- Primary irritant effect:
- On the skin: No irritating effect.
- On the eye: No irritating effect.
- Additional toxicological information:

The product is not subject to classification according to internally approved calculation methods for preparations.

When used and handled according to specifications, the product does not have any harmful effects according to our experience and the information provided to us.

... Carcinogenic categories:

• IARC (International Agency for Research on Cancer):
None of the ingredients are listed.
• NTP (National Toxicology Program):
None of the ingredients are listed.
• OSHA-Ca (Occupational Safety & Health Administration):
None of the ingredients are listed.

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 7 of 9



Issue date 02/20/2019

Safety Data Sheet (SDS)

OSHA HazCom Standard 29 CFR 1910.1200(g) and GHS Rev 03.

Page 6/7

Reviewed on 02/20/2019

Trade Name: Hydrophilix Polymer Solution

12 Ecological Information

- Toxicity:
- Aquatic toxicity: No further relevant information available.
- Persistence and degradability: No further relevant information available.
- Behavior in environmental systems:
- Bioaccumulative potential: No further relevant information available.
- Mobility in soil: No further relevant information available.
- Additional ecological information:
- General notes: Generally not hazardous for water.
- Results of PBT and vPvB assessment:
- PBT: Not applicable.
- vPvB: Not applicable.
- Other adverse effects: No further relevant information available.

13 Disposal Considerations

- Waste treatment methods
- Recommendation:
Smaller quantities can be disposed of with household waste.
Observe all federal, state and local environmental regulations when disposing of this material.
- Uncleaned packaging
- Recommendation: Disposal must be made according to official regulations.
- Recommended cleansing agent: Water, if necessary with cleansing agents.

14 Transport Information

- UN-Number:
- DOT, ADR/ADN, ADN, IMDG, IATA Non-Regulated Material
- UN proper shipping name:
- DOT, ADR/ADN, ADN, IMDG, IATA Non-Regulated Material
- Transport hazard class(es):
- DOT, ADR/ADN, ADN, IMDG, IATA
- Class: Non-Regulated Material
- Packing group:
- DOT, ADR/ADN, IMDG, IATA Non-Regulated Material
- Environmental hazards: Not applicable.
- Special precautions for user: Not applicable.
- Transport in bulk according to Annex II of MARPOL73/78 and the IBC Code: Not applicable.
- UN "Model Regulation": Non-Regulated Material

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work	Page 8 of 9



Issue date 02/20/2019

Safety Data Sheet (SDS)

OSHA HazCom Standard 29 CFR 1910.1200(g) and GHS Rev 03.

Page 7/7

Reviewed on 02/20/2019

Trade Name: Hydrophilix Polymer Solution

15 Regulatory Information

- Section 355 (extremely hazardous substances):

None of the ingredients are listed.

- Section 313 (Specific toxic chemical listings):

None of the ingredients are listed.

- TSCA (Toxic Substances Control Act):

All ingredients are listed or exempt from listing.

- California Proposition 65:

- Chemicals known to cause cancer:

None of the ingredients are listed.

- Chemicals known to cause reproductive toxicity for females:

None of the ingredients are listed.

- Chemicals known to cause reproductive toxicity for males:

None of the ingredients are listed.

- Chemicals known to cause developmental toxicity:

None of the ingredients are listed.

- New Jersey Right-to-Know List:

None of the ingredients are listed.

- New Jersey Special Hazardous Substance List:

None of the ingredients are listed.

- Pennsylvania Right-to-Know List:

None of the ingredients are listed.

- Pennsylvania Special Hazardous Substance List:

None of the ingredients are listed.

- Carcinogenic categories:

- EPA (Environmental Protection Agency):

None of the ingredients are listed.

- TLV (Threshold Limit Value established by ACGIH):

None of the ingredients are listed.

- NIOSH-Ca (National Institute for Occupational Safety and Health):

None of the ingredients are listed.

- GHS label elements Non-Regulated Material

- Hazard pictograms: Non-Regulated Material

- Signal word: Non-Regulated Material

- Hazard statements: Non-Regulated Material

- National regulations:

None of the ingredients are listed.

- Chemical safety assessment: A Chemical Safety Assessment has not been carried out.

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-3000
	Description: Specification, Solution, Hydrophilix Polymer	Rev.: A
	This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work	Page 9 of 9



Issue date 02/20/2019

Safety Data Sheet (SDS)

OSHA HazCom Standard 29 CFR 1910.1200(g) and GHS Rev 03.

Page 7/7

Reviewed on 02/20/2019

Trade Name: Hydrophilix Polymer Solution

16 Other Information

The information and recommendations in this safety data sheet are, to the best of our knowledge, accurate as of the date of issue. Nothing herein shall be deemed to create warranty, expressed or implied, and shall not establish a legally valid contractual relationship. It is the responsibility of the user to determine applicability of this information and the suitability of the material or product for any particular purpose.

· Date of preparation / last revision: 02/20/2019 / -

· Abbreviations and acronyms:

ADR: The European Agreement concerning the International Carriage of Dangerous Goods by Road

ADN: The European Agreement concerning the International Carriage of Dangerous Goods by Inland Waterways

IMDG: International Maritime Code for Dangerous Goods

DOT: US Department of Transportation

IATA: International Air Transport Association

ACGIH: American Conference of Governmental Industrial Hygienists

EINECS: European Inventory of Existing Commercial Chemical Substances

ELINCS: European List of Notified Chemical Substances

CAS: Chemical Abstracts Service (division of the American Chemical Society)

NFPA: National Fire Protection Association (USA)

HMIS: Hazardous Materials Identification System (USA)

VOC: Volatile Organic Compounds (USA, EU)

LC50: Lethal concentration, 50 percent

LD50: Lethal dose, 50 percent

PBT: Persistent, Bioaccumulative and Toxic

vPvB: very Persistent and very Bioaccumulative

NIOSH: National Institute for Occupational Safety and Health

OSHA: Occupational Safety & Health Administration

TLV: Threshold Limit Value

PEL: Permissible Exposure Limit

REL: Recommended Exposure Limit

· * Data compared to the previous version altered.

SDS created by MSDS Authoring Services www.msdsauthoring.com +1-877-204-9106

Change History

Revision	Change Description
A	Initial release

Date Inspected	Part No.			Rev.	Description											
	HDX-3000			A	Solution, Polymer, Hydrophilix											
Lot No.		Lot QTY		Approved Supplier(s)												
				Hydrophilix												
Specification	No Damage to Outer Packaging.	Container is sealed and is not leaking	C of C is present with supplier's inspection data that complies to the specification.													
Revision			N/A													
Lower Spec Limit (LSL)	N/A	N/A	N/A													
Upper Spec Limit (USL)	N/A	N/A	N/A													
UOM	N/A	N/A	N/A													
Measuring Equipment Type	N/A	N/A	N/A													
Asset No.	N/A	N/A	N/A													
Cal Due Date	N/A	N/A	N/A													
AQL	N/A	N/A	N/A													
Sample Size	Shipment	Shipment	Shipment													
Sample No.	Variable/ Attribute Data															
1																
2																
3																
4																
5																
Result Summary																
Pass QTY																
Fail QTY																
Result (P/F)																
Statistical Summary																
Min	N/A	N/A	0.000	0.000	0.000	0.000	0.000	0.000	0.000	0.000						
Max	N/A	N/A	0.000	0.000	0.000	0.000	0.000	0.000	0.000	0.000						
Mean	N/A	N/A	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!						
StDev	N/A	N/A	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!						
Comments																
Disposition:							NCMR No:									
Inspected By:							Date:									

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-5000
	Description: Bag, 2" W x 18" L x 2 Mil Thk, Polyethylene	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 1 of 2

Material Requirements ☐ N/A			
- Primary: Supplied by Customer - Secondary: ULINE P/N S-1426			
Dimensional Requirements ☐ N/A			
2" W x 18" L x 2 Mil Thk			
Visual Requirements ☒ N/A			
Label Content Specification ☒ N/A			
Certificate of Compliance (C of C)	☐ YES ☐ NO ☒ N/A	Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A
Material Compliance			
REACH Compliant	☐ YES ☐ NO ☒ N/A	ROHS Compliant	☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES ☐ NO ☒ N/A	FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A	BPA Free	☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES ☐ NO ☒ N/A	DEHP Free	☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A	Documented evidence for compliance to incorporated material specification file	
Storage Requirements ☒ N/A			
Shelf-Life Requirements ☒ N/A			
Packaging Requirements ☒ N/A			
Supplier will notified Pathway Innovation in advance of any changes prior to the implementation of the change in the supplied material.			

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-5000
	Description: Bag, 2" W x 18" L x 2 Mil Thk, Polyethylene	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 2 of 2



2 x 18" 2 Mil Industrial Poly Bags



[Enlarge](#)

OVER 290 SIZES IN STOCK

Most popular thickness for shipping parts, clothing or protecting food.

Meets FDA and USDA specifications for food contact.
100% virgin, high clarity polyethylene film.

SOLD IN CARTON QUANTITIES

MODEL NO.	OUTSIDE W x L	LBS./ CTN.	PRICE PER CARTON				QTY./ CTN.	ADD TO CART	
			1	5	10	20+			
S-1426	2 x 18"	4	\$35	\$33	\$31	\$29	1,000	1	ADD

Change History

Revision	Change Description
A	Initial release

Date Inspected	Part No.			Rev.	Description															
	HDX-5000			A	Bag, 2" W x 18" L x 2 Mil Thk, Polyethylene															
Lot No.		Lot QTY		Approved Supplier(s)																
				ULINE																
Specification	No Damage to Outer Packaging	Verify parts are free of foreign matter or contamination																		
Revision	N/A	N/A																		
Lower Spec Limit (LSL)	N/A	N/A																		
Upper Spec Limit (USL)	N/A	N/A																		
UOM	N/A	N/A																		
Measuring Equipment Type	Visual	Visual																		
Asset No.	N/A	N/A																		
Cal Due Date	N/A	N/A																		
AQL	N/A	4.0																		
Sample Size	Shipment																			
Sample No.	Variable/ Attribute Data																			
1																				
2																				
3																				
4																				
5																				
6																				
7																				
8																				
9																				
10																				
11																				
12																				
13																				
14																				
15																				
16																				
17																				
18																				
19																				
20																				
21																				
22																				
23																				
24																				
25																				
26																				
27																				
28																				
29																				
Result Summary																				
Pass QTY																				
Fail QTY																				
Result (P/F)																				
Statistical Summary																				
Min	N/A	N/A	N/A	N/A	N/A	N/A														
Max	N/A	N/A	N/A	N/A	N/A	N/A														
Mean	N/A	N/A	N/A	N/A	N/A	N/A														
StDev	N/A	N/A	N/A	N/A	N/A	N/A														
Comments																				
Disposition:									NCR No:											
Inspected By:									Date:											

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-5001
	Description: Specification, Pouch, Tyvek Breathable, 3.5" x 21", Amcor TLP-058	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	
		Page 1 of 3

Material Requirements				☐ N/A			
Amcor 3.5" x 21" Pouch, P/N TLP-058 3.5X21							
Dimensional Requirements				☐ N/A			
3.5" x 21"							
Visual Requirements				☐ N/A			
No damage to outer packaging box							
Label Content Specification				☐ N/A			
Certificate of Compliance (C of C)	☒ YES	☐ NO	☐ N/A	Certificate of Analysis (C of A)	☐ YES	☐ NO	☒ N/A
Material Compliance							
REACH Compliant	☐ YES	☐ NO	☒ N/A	ROHS Compliant	☐ YES	☐ NO	☒ N/A
Conflict Material	☐ YES	☐ NO	☒ N/A	FDA (Medical Grade) Compliant	☐ YES	☐ NO	☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES	☐ NO	☒ N/A	BPA Free	☐ YES	☐ NO	☒ N/A
Latex Free	☐ YES	☐ NO	☒ N/A	DEHP Free	☐ YES	☐ NO	☒ N/A
Animal Derive Material (ADM) Free	☐ YES	☐ NO	☒ N/A	Documented evidence for compliance to incorporated material specification file			
Storage Requirements				☒ N/A			
Shelf-Life Requirements				☒ N/A			
Packaging Requirements				☒ N/A			
Supplier will notified Pathway Innovation in advance of any changes prior to the implementation of the change in the supplied material.							

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-5001
	Description: Specification, Pouch, Tyvek Breathable, 3.5" x 21", Amcor TLP-058	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 2 of 3

Poly Structure

Advanced Breathable Pouches for Medical Applications

Breathability and performance suitable for a wide variety of medical applications. Selection of laminates to meet your needs for tensile strength, tear strength, or puncture resistance. Available with coated and uncoated Tyvek® and different gauges of material to meet sterility applications. Extensive line of laminates allow for the proper choice of materials. Incorporates an adhesive lamination of polyester and polyethylene. Suitable for radiation and EO sterilization for wide range of medical devices. Uncoated Tyvek provides breathability and microbial barrier. Film portion of pouch has excellent clarity and strength properties. Consistent peelable seal is feature of pouch.

Enquire



Applications

Healthcare

Medical

- Minimally Invasive Devices
- Catheters, Guidewires & Stents

Technical Specifications

Delivery

- Product delivered in boxes. Delivery process depending on order quantity.

Overview

Features

- Suitable for EtO and Irradiation Sterilisation
- Available as Uncoated Tyvek
- Clean fiber tear free opening

Customer benefits

- Breathability and performance suitable for a wide variety of medical applications
- Selection of laminates to meet your needs for tensile strength, tear strength, or puncture resistance;

End consumer benefits

- Easy Opening - Easy and clean opening
- Available in various sizes
- product clarity

End consumer types

- Hospitals
- Medical-Device-Manufacturers

Materials

- Tyvek

Continents available

- North America
- South America

Marketing benefits

- Consistent peelable seal reduces mishandling errors

Minimum order Units per year

- 1 box

Product safety

- Comprehensive Factbook available with aging data

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-5001
	Description: Specification, Pouch, Tyvek Breathable, 3.5" x 21", Amcor TLP-058	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 3 of 3

Change History

Revision	Change Description
A	Initial release

Date Inspected	Part No.			Rev.	Description					
	HDX-5001			A	Pouch, Chevron, 3.5" W x 21.0" L, Amcor TLP-058					
Lot No.		Lot QTY		Approved Supplier(s)						
				Amcor Flexibles						
Specification	No Damage to Outer Packaging	Verify there are no voids/channels in the 3 sealed sides	Verify there is no pinholes, tears, or voids in the tyvek and poly material	Verify there are no creases on the unsealed side	Verify parts are free of foreign matter or contamination	C of C is present with supplier's inspection data that complies to the specification.				
Revision	N/A	N/A	N/A	N/A	N/A	N/A				
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A				
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A				
UOM	N/A	N/A	N/A	N/A	N/A	N/A				
Measuring Equipment Type	Visual	Visual	Visual	Visual	Visual	N/A				
Asset No.	N/A	N/A	N/A	N/A	N/A	N/A				
Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A				
AQL	N/A	2.5	2.5	2.5	2.5	N/A				
Sample Size	Shipment		0	0	0	Shipment				
Sample No.	Variable/ Attribute Data									
1										
2										
3										
4										
5										
6										
7										
8										
9										
10										
11										
12										
13										
14										
15										
16										
17										
18										
19										
20										
21										
22										
23										
24										
25										
26										
27										
28										
29										
30										
31										
32										
33										
34										
35										
Result Summary										
Pass QTY										
Fail QTY										
Result (P/F)										
Statistical Summary										
Min	N/A	N/A	N/A	N/A	N/A	N/A				
Max	N/A	N/A	N/A	N/A	N/A	N/A				
Mean	N/A	N/A	N/A	N/A	N/A	N/A				
StDev	N/A	N/A	N/A	N/A	N/A	N/A				
Comments										
Disposition:							NCR No:			
Inspected By:							Date:			

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-5002
	Description: Specification, Box, 10-Pack Inner, 20.75 L x 4.12 W x 2.37 H	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 1 of 1

Material Requirements ☐ N/A			
Box 10-Pack Inner, 20.75 L x 4.12 W x 2.37 H (Primary: Supplied by PSM Secondary: Supplied by SquarePeg)			
Dimensional Requirements ☐ N/A			
20.75 L x 4.12 W x 2.37 (REF)			
Visual Requirements ☐ N/A			
No damage to outer packaging box			
Label Content Specification ☐ N/A			
Certificate of Compliance (C of C)	☐ YES ☐ NO ☒ N/A	Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A
Material Compliance			
REACH Compliant	☐ YES ☐ NO ☒ N/A	ROHS Compliant	☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES ☐ NO ☒ N/A	FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A	BPA Free	☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES ☐ NO ☒ N/A	DEHP Free	☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A	Documented evidence for compliance to incorporated material specification file	
Storage Requirements ☒ N/A			
Shelf-Life Requirements ☒ N/A			
Packaging Requirements ☒ N/A			
Supplier will notified Pathway Innovation in advance of any changes prior to the implementation of the change in the supplied material.			

Change History

Revision	Change Description
A	Initial release

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-5003
	Description: Specification, Box, 22 in. L X 14 in. W X 12 in. H, 275 LB. Test Double Wall	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 1 of 2

Material Requirements ☐ N/A			
• Uline Model No. S-18166 or equivalent			
Dimensional Requirements ☐ N/A			
• 22" L X 14" W X 12" H (REF)			
Visual Requirements ☒ N/A			
•			
Label Content Specification ☒ N/A			
•			
Certificate of Compliance (C of C)	☐ YES ☐ NO ☒ N/A	Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A
Material Compliance			
REACH Compliant	☐ YES ☐ NO ☒ N/A	ROHS Compliant	☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES ☐ NO ☒ N/A	FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A	BPA Free	☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES ☐ NO ☒ N/A	DEHP Free	☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A	Documented evidence for compliance to incorporated material specification file	
Storage Requirements ☒ N/A			
•			
Shelf-Life Requirements ☒ N/A			
•			
Packaging Requirements ☒ N/A			
•			
Supplier will notify Pathway Innovation in advance of any changes prior to the implementation of the change in the supplied material.			



3200 Premium Wax/Resin Ribbons

The Zebra 3200 high-performance wax/resin thermal transfer ribbon is a special ink formulation for use with Zebra thermal transfer printers. This ribbon provides outstanding print quality combined with chemical and solvent resistance. By providing smear and scratch resistance on Zebra coated paper and matte synthetic label stocks, this ribbon results in excellent staying power and durability in harsh environments.

UL/ cUL / CSA Certifications with the following materials:

- Z-Xtreme 4000T White and Silver (UL, cUL, CSA)
- Z-Xtreme 4000T High-Tack White & Silver (UL, cUL, CSA)
- Z-Xtreme 2000T White, Silver & Clear (UL/cUL)
- Z-Supreme 3000T White (UL)



TECHNICAL SPECIFICATIONS

Attribute	Specification
Color	Black
Ink Formulation	Wax/Resin
Base Film	Polyester
Base Film Caliper	4.5 microns
Total Caliper	8.6 microns +/- 10%
Scanning Capabilities	IR and Visible Light

Recommended Storage Conditions: 72° F (22° C) at 50% RH

Shelf Life 1 year

For testing of this product, please order sample ribbon 03200BK08005.

SUGGESTED APPLICATIONS

- Labels exposed to solvents or chemicals



- Labels subject to abrasion or friction



- Outdoor labeling applications in extreme environments



- Chemical drum labeling applications



- Labels exposed to water, steam, alkali or acid solutions



For more information, visit www.zebra.com/supplies

Product Performance and Suitability: The information contained in this document is to be used for guidance only and is not intended for use in setting specifications. All purchasers of Zebra products shall be solely responsible for independently determining if the product conforms to all requirements of their unique application.

NA and Corporate Headquarters | +1 800 423 0442 | inquiry4@zebra.com

©2016 ZIH Corp and/or its affiliates. All rights reserved. ZEBRA and the stylized Zebra head are trademarks of ZIH Corp, registered in many jurisdictions worldwide. All other trademarks are the property of their respective owners. 10/2016

Pathway Innovation LLC.	Raw Material Specification	Part No.: HDX-5005
	Description: Tape, Packaging	Rev.: A
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Page 1 of 2

Material Requirements ☐ N/A			
• Clear Polypropylene w/acrylic adhesive			
Dimensional Requirements ☐ N/A			
• REF 2" x 55 yds x 2 mil or 2" x 110 yds x 2 mil			
Visual Requirements ☒ N/A			
•			
Label Content Specification ☒ N/A			
•			
Certificate of Compliance (C of C)	☐ YES ☐ NO ☒ N/A	Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A
Material Compliance			
REACH Compliant	☐ YES ☐ NO ☒ N/A	ROHS Compliant	☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES ☐ NO ☒ N/A	FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A	BPA Free	☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES ☐ NO ☒ N/A	DEHP Free	☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A	Documented evidence for compliance to incorporated material specification file	
Storage Requirements ☒ N/A			
•			
Shelf-Life Requirements ☒ N/A			
•			
Packaging Requirements ☒ N/A			
•			
Supplier will notified Pathway Innovation in advance of any changes prior to the implementation of the change in the supplied material.			

Change History

Revision	Change Description
A	Initial release

Pathway Innovation LLC.	Raw Material Specification				Part No.: HDX-5005
	Description: Tape, Packaging				Rev.: A
	This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work				Page 2 of 2

Uline Industrial Tape - 2 Mil, 2" x 110 yds, Clear

2.0 MIL



Our best selling tape. Strong solvent acrylic adhesive sticks to virtually any surface.

DURABLE: Industrial Tape stands up to heat, humidity and cold. **Application range 0°F to 140°F.**

RELIABLE: Long shelf life.

APPEARANCE: Crystal clear - Won't yellow. Makes packages look their best.

± Free Offer

Enlarge

20 LBS. PER INCH TENSILE STRENGTH

MODEL NO.	COLOR	MIL	SIZE	ROLLS/ CASE	LBS./ CASE	PRICE PER ROLL					ADD TO CART		
						18	36	72	144	288+			
S-423	Clear	2.0	2" x 110 yds.	36	22	\$3.25	\$1.69 per roll					18	ADD

[Additional Info](#)

[Email Page](#)

[Shopping Lists](#)

[Request a Catalog](#)

DIMENSIONS:

- Core: 3"
- Diameter: 4 1/2"

FEATURES:

- Rolls have a paper starter tab.

MATERIAL:

- Backing: Polypropylene film
- Latex free

SHELF LIFE:

- 1 year

TEMPERATURE:

- Performance: -60 to 140 F

USAGE:

- Takes approximately 40 minutes to reach maximum bond.

Availability: [In Stock](#)

Unit Weight: 0.56 lbs.

[SDS Sheet](#)

[Catalog Page 769](#)

Date Inspected	Part No.			Rev.	Description		
	HDX-6000			A	Catheter, 14 Fr x 10ml Silicone		
Lot No.		Lot QTY		Approved Supplier(s)			
				Ningbo Yingmed			
Specification	No Damage to Outer Packaging and Catheters are double bagged.	Minimum identification is present (P/N, Rev., Lot No.).	C of C is present with supplier's inspection data that complies to the specification.	Parts are clean and free of contamination	Confirm print on catheter matches configuration "14 Fr 4.7 mm 10ml"	TM-0013 Foley Catheter Leak Test	
Revision	N/A	N/A	N/A	N/A			
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A	
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A	
UOM	N/A	N/A	N/A	N/A	N/A	N/A	
Measuring Equipment Type	N/A	N/A	N/A	N/A	N/A	N/A	
Asset No.	N/A	N/A	N/A	N/A	N/A	N/A	
Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A	
AQL	N/A	N/A	N/A	4.0	4.0	1.0	
Sample Size	Shipment	Shipment	Shipment				
Sample No.	Variable/ Attribute Data						
1							
2							
3							
4							
5							
6							
7							
8							
9							
10							
11							
12							
13							
14							
15							
16							
17							
18							
19							
20							
21							
22							
23							
24							
25							
26							
27							
28							
29							
Result Summary							
Pass QTY							
Fail QTY							
Result (P/F)							
Statistical Summary							
Min	N/A	N/A	N/A	N/A	N/A	N/A	
Max	N/A	N/A	N/A	N/A	N/A	N/A	
Mean	N/A	N/A	N/A	N/A	N/A	N/A	
StDev	N/A	N/A	N/A	N/A	N/A	N/A	
Comments							
Disposition:						NCMR No:	
Inspected By:						Date:	

ITEM NO.	DESCRIPTION	MATERIAL	HARDNESS	COLOR	QTY.	NOTES
1	TUBING	SILICONE	80 SHORE A	NATURAL/XRAY	1	X-RAY COLOR IS BLUE
2	BALLOON CUFF	SILICONE	25 SHORE A	NATURAL	1	ORDINARY BALLOON TYPE
3	TIP	SILICONE	50 SHORE A	BLUE	1	TIP HAS 33% CONCENTRATION OF BARIUM SULFATE
4	3-WAY FUNNEL	SILICONE	57 SHORE A	NATURAL	1	N/A
5	CHECK VALVE	PP/SS/PS/SILICONE	70 SHORE A	WHITE	1	BESPAK MED
6	SLEEVE	ABS	80 SHORE A	ORANGE	1	N/A

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
N/A	A	INITIAL RELEASE	2019-06-10

D

D

C


C

B

B

A

A

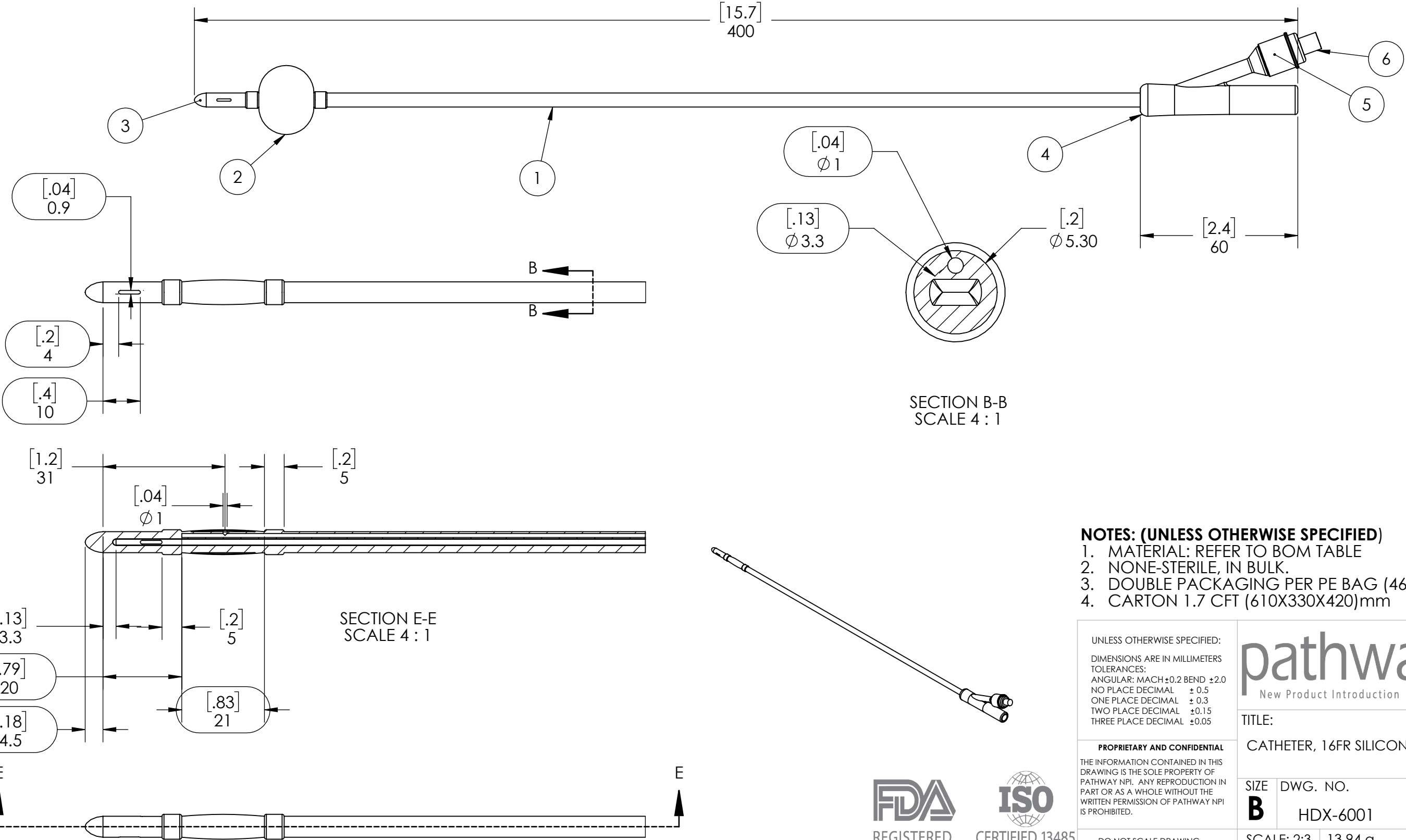

- NOTES: (UNLESS OTHERWISE SPECIFIED)**
- 1. MATERIAL: REFER TO BOM TABLE
 - 2. NONE-STERILE, IN BULK.
 - 3. DOUBLE PACKAGING PER PE BAG (460X690)mm
 - 4. CARTON 1.7 CFT (610X330X420)mm

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.



TITLE:
CATHETER, 16FR SILICONE, 5 mL

SIZE	DWG. NO.	REV
B	HDX-6001	A

SCALE: 2:3 13.94 g SHEET 1 OF 1

DO NOT SCALE DRAWING



Date Inspected	Part No.			Rev.	Description		
	HDX-6001			A	Catheter, 16 Fr x 5ml Silicone		
Lot No.		Lot QTY		Approved Supplier(s)			
				Ningbo Yingmed			
Specification	No Damage to Outer Packaging and Catheters are double bagged.	Minimum identification is present (P/N, Rev., Lot No.).	C of C is present with supplier's inspection data that complies to the specification.	Parts are clean and free of contamination	Confirm print on catheter matches configuration "16 Fr 5.3 mm 5ml"	TM-0013 Foley Catheter Leak Test	
Revision	N/A	N/A	N/A	N/A			
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A	
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A	
UOM	N/A	N/A	N/A	N/A	N/A	N/A	
Measuring Equipment Type	N/A	N/A	N/A	N/A	N/A	N/A	
Asset No.	N/A	N/A	N/A	N/A	N/A	N/A	
Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A	
AQL	N/A	N/A	N/A	4.0	4.0	1.0	
Sample Size	Shipment	Shipment	Shipment				
Sample No.	Variable/ Attribute Data						
1							
2							
3							
4							
5							
6							
7							
8							
9							
10							
11							
12							
13							
14							
15							
16							
17							
18							
19							
20							
21							
22							
23							
24							
25							
26							
27							
28							
29							
Result Summary							
Pass QTY							
Fail QTY							
Result (P/F)							
Statistical Summary							
Min	N/A	N/A	N/A	N/A	N/A	N/A	
Max	N/A	N/A	N/A	N/A	N/A	N/A	
Mean	N/A	N/A	N/A	N/A	N/A	N/A	
StDev	N/A	N/A	N/A	N/A	N/A	N/A	
Comments							
Disposition:						NCMR No:	
Inspected By:						Date:	

ITEM NO.	DESCRIPTION	MATERIAL	HARDNESS	COLOR	QTY.	NOTES
1	TUBING	SILICONE	80 SHORE A	NATURAL/XRAY	1	X-RAY COLOR IS BLUE
2	BALLOON CUFF	SILICONE	25 SHORE A	NATURAL	1	ORDINARY BALLOON TYPE
3	TIP	SILICONE	50 SHORE A	BLUE	1	TIP HAS 33% CONCENTRATION OF BARIUM SULFATE
4	3-WAY FUNNEL	SILICONE	57 SHORE A	NATURAL	1	N/A
5	CHECK VALVE	PP/SS/PS/SILICONE	70 SHORE A	WHITE	1	BESPAK MED
6	SLEEVE	ABS	80 SHORE A	RED	1	N/A

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
N/A	A	INITIAL RELEASE	2019-06-10

D

D

C


C

B

B

A

A

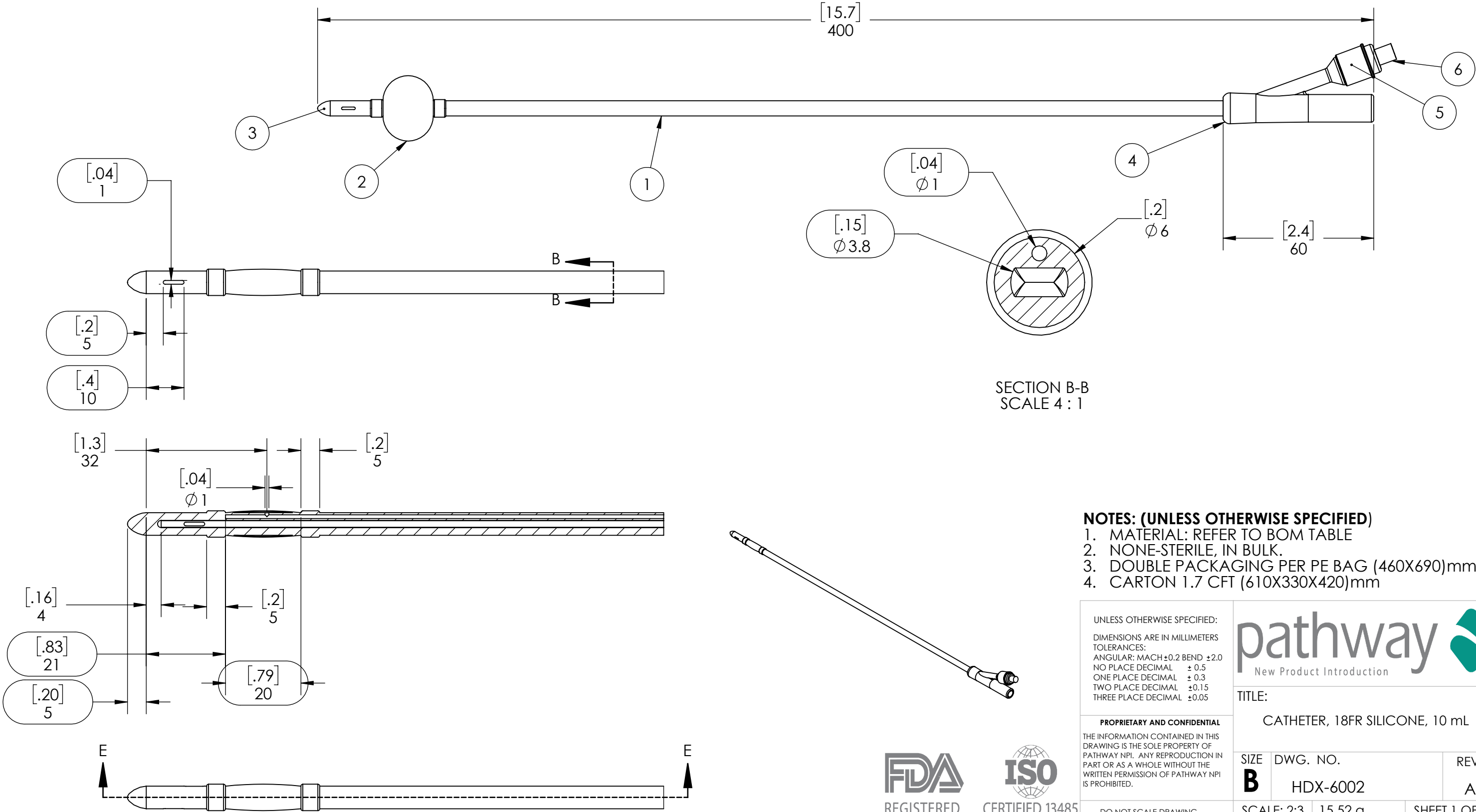

NOTES: (UNLESS OTHERWISE SPECIFIED)
1. MATERIAL: REFER TO BOM TABLE
2. NONE-STERILE, IN BULK.
3. DOUBLE PACKAGING PER PE BAG (460X690)mm
4. CARTON 1.7 CFT (610X330X420)mm

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.

DO NOT SCALE DRAWING



TITLE:
CATHETER, 18FR SILICONE, 10 mL

SIZE B	DWG. NO. HDX-6002	REV A
------------------	----------------------	----------

SCALE: 2:3	15.52 g	SHEET 1 OF 1
------------	---------	--------------


Last Saved: Friday, July 31, 2020 2:24:06 PM
SOFTWARE: SW 2018 SPS




Date Inspected	Part No.			Rev.	Description		
	HDX-6002			A	Catheter, 18 Fr x 10ml Silicone		
Lot No.		Lot QTY		Approved Supplier(s)			
				Ningbo Yingmed			
Specification	No Damage to Outer Packaging and Catheters are double bagged.	Minimum identification is present (P/N, Rev., Lot No.).	C of C is present with supplier's inspection data that complies to the specification.	Parts are clean and free of contamination	Confirm print on catheter matches configuration "18 Fr 6.0 mm 10ml"	TM-0013 Foley Catheter Leak Test	
Revision	N/A	N/A	N/A	N/A			
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A	
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A	
UOM	N/A	N/A	N/A	N/A	N/A	N/A	
Measuring Equipment Type	N/A	N/A	N/A	N/A	N/A	N/A	
Asset No.	N/A	N/A	N/A	N/A	N/A	N/A	
Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A	
AQL	N/A	N/A	N/A	4.0	4.0	1.0	
Sample Size	Shipment	Shipment	Shipment				
Sample No.	Variable/ Attribute Data						
1							
2							
3							
4							
5							
6							
7							
8							
9							
10							
11							
12							
13							
14							
15							
16							
17							
18							
19							
20							
21							
22							
23							
24							
25							
26							
27							
28							
29							
Result Summary							
Pass QTY							
Fail QTY							
Result (P/F)							
Statistical Summary							
Min	N/A	N/A	N/A	N/A	N/A	N/A	
Max	N/A	N/A	N/A	N/A	N/A	N/A	
Mean	N/A	N/A	N/A	N/A	N/A	N/A	
StDev	N/A	N/A	N/A	N/A	N/A	N/A	
Comments							
Disposition:						NCMR No:	
Inspected By:						Date:	

NOTES:(UNLESS OTHERWISE SPECIFIED)

- 1 LABEL ARTWORK PRINTED TOWARDS OPEN END OF POUCH. ITEM 6.
- 2 IFU ITEM 5 TO BE INSERTED INTO 10-PACK BOX ITEM 7.
- 3 SEAL SHIPPING BOX WITH PACKING TAPE OR EQUIVALENT.
- 4 REFERENCE ITEM 4 FOR POUCH, 10-PACK BOX AND SHIPPING BOX LABEL ARTWORK REQUIREMENTS.
- 5 APPLY POLYMER SOLUTION PART NO. HDX-3000 TO ITEM 9 ON HATCHED AREA AS SHOWN.
- 6 BALLOON VOLUME 10mL.

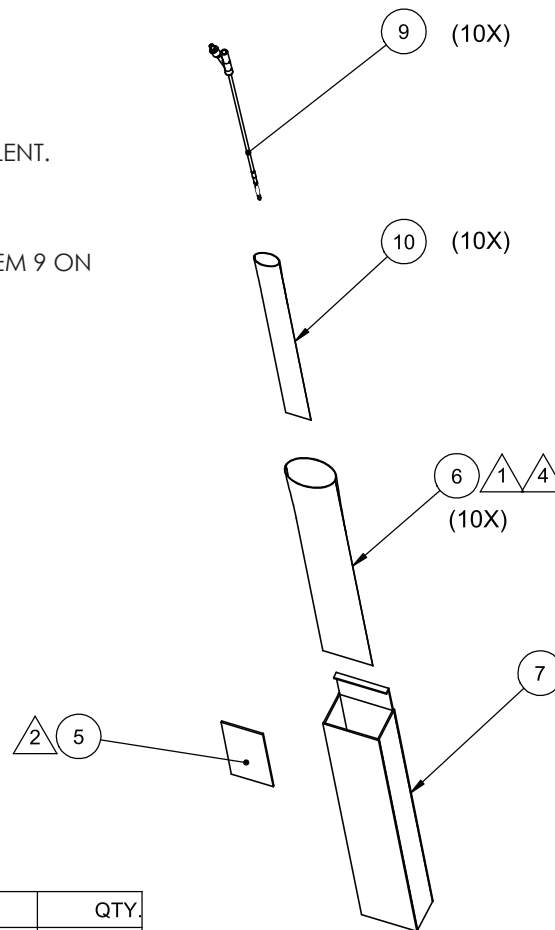

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
N/A	A	INITIAL RELEASE	06/10/2019
N/A	B	ADDED PRODUCTION SPECIFICATIONS ON PG 3.	03/18/2020

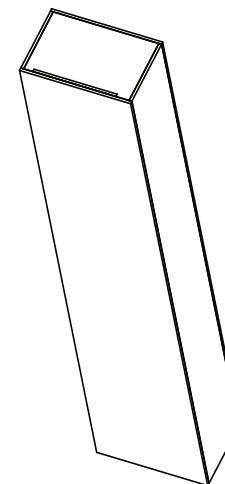

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
4	HDX-4000	ARTWORK, LABEL, HDX FOLEY CATHETER (POUCH AND INNER CARTON)	3
5	HDX-4001	ARTWORK,IFU,HDX FOLEY CATHETER	1
6	HDX-5000	POUCH, CHEVRON, 3.5X21, AMCRO TLP-058	1
7	HDX-5001	BOX, 10-PACK, CATHETER	2
8	HDX-5002	BOX	1
9	HDX-6000	CATHETER, 14FR SILICONE	1
10	HDX-5003	SLEEVE	1
11	HDX-5004	LABEL	3

UNLESS OTHERWISE SPECIFIED:
 DIMENSIONS ARE IN MILLIMETERS
 TOLERANCES:
 ANGULAR: MACH ± 0.2 BEND ± 2.0
 NO PLACE DECIMAL ± 0.5
 ONE PLACE DECIMAL ± 0.3
 TWO PLACE DECIMAL ± 0.15
 THREE PLACE DECIMAL ± 0.05

PROPRIETARY AND CONFIDENTIAL
 THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.




TITLE:
BNS, 14 FR CATHETER

SIZE DWG. NO.
B HDX-8000


SCALE: 1:25 g

SHEET 1 OF 3

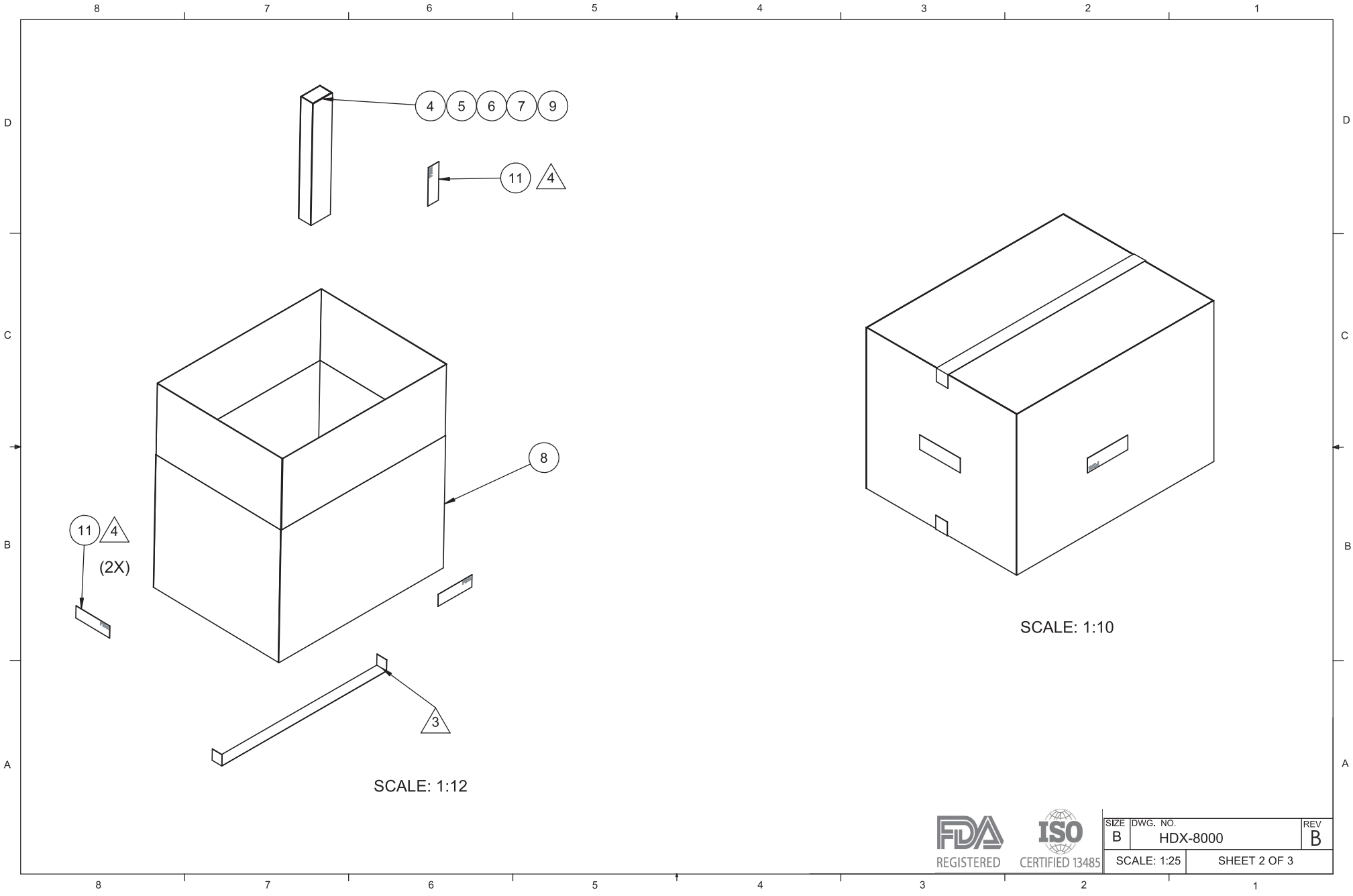


SCALE: 1:10

ITEM NO.	DESCRIPTION	FUNCTION
1	TUBING	CATHETER BODY
2	BALLOON CUFF	RETENTION BALLOON
3	TIP	SMOOTH INSERTION
4	2-WAY FUNNEL	FLUID DRAINAGE
5	CHECK VALVE	BALLOON INFLATION/DEFLATION
6	SLEEVE	SYRINGE INTERFACE

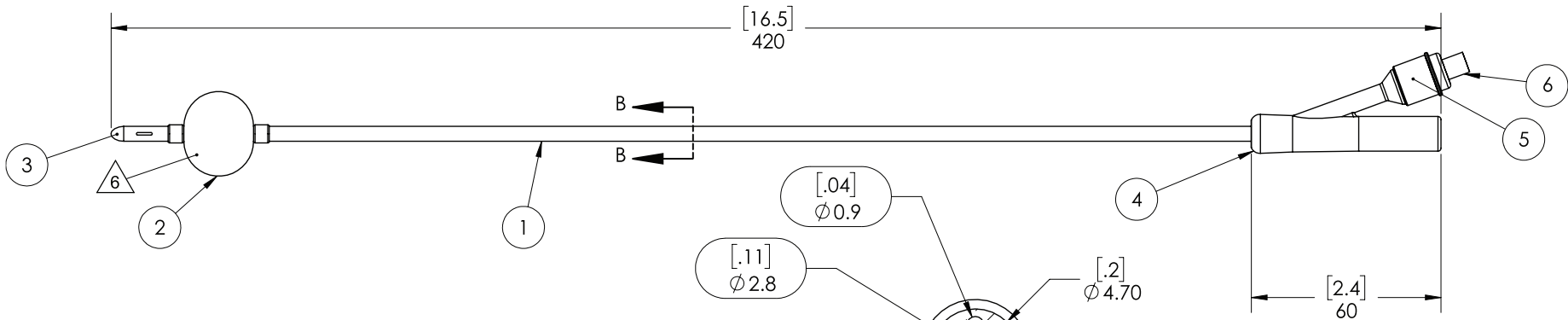


SECTION B-B
SCALE 4 : 1

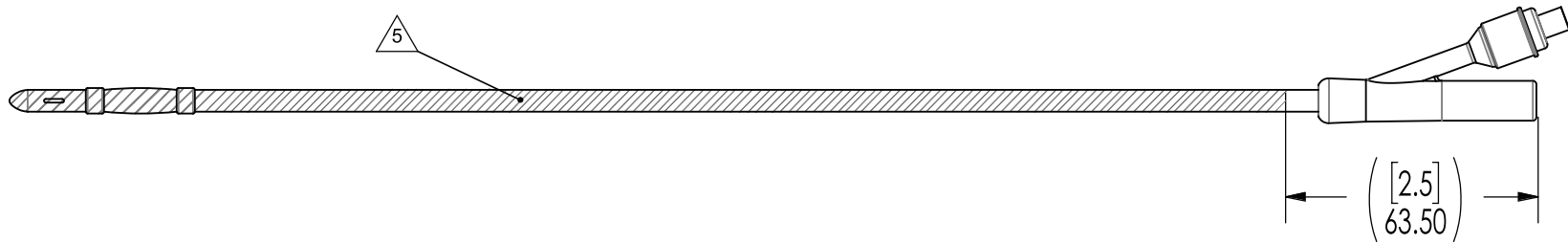


Date Inspected	Part No.				Rev.	Description			
	HDX-8000				B	Hydrophilix 14 Fr x 10 ml Catheter, Non-Sterile, 10-Pack			
Lot No.		Lot QTY			Approved Supplier(s)				
Specification	Catheter has been coated with Polymer Solution P/N HDX-3000. MP-HDX-0000	Catheter has been rinsed and free of debris. MP-HDX-0000	Coated and rinsed Catheter meets the requirements of the Spectrophotometer Test. TM-HDX-0000	P/N, Lot No. and Exp date on labels match Work Order requirements. MP-HDX-0000	Catheter is properly oriented in the sleeve and pouch. Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. MP-HDX-0001	Correct quantity in 10-Pack Box. MP-HDX-0001			
Revision		0		0		0			
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A			
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A			
UOM	N/A	N/A	N/A	N/A	N/A	N/A			
Measuring Equipment Type	Visual	Visual	Spectrophotometer	Visual	Visual	Visual			
Asset No.	N/A	N/A		N/A	N/A	N/A			
Cal Due Date	N/A	N/A		N/A	N/A	N/A			
AQL	6.5	6.5	6.5	6.5	6.5	6.5			
Sample Size									
Sample No.	Variable/ Attribute Data								
1									
2									
3									
4									
5									
6									
7									
8									
9									
10									
11									
12									
13									
14									
15									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	N/A	N/A	N/A	N/A	N/A	N/A	0.000	0.000	0.000
Max	N/A	N/A	N/A	N/A	N/A	N/A	0.000	0.000	0.000
Mean	N/A	N/A	N/A	N/A	N/A	N/A	#DIV/0!	#DIV/0!	#DIV/0!
StDev	N/A	N/A	N/A	N/A	N/A	N/A	#DIV/0!	#DIV/0!	#DIV/0!
Comments									
Disposition:						NCMR No.:			
Inspected By:						Date:			

NOTES:(UNLESS OTHERWISE SPECIFIED)

- 1 LABEL ARTWORK PRINTED TOWARDS OPEN END OF POUCH. ITEM 6.
- 2 IFU ITEM 5 TO BE INSERTED INTO 10-PACK BOX ITEM 7.
- 3 SEAL SHIPPING BOX WITH PACKING TAPE OR EQUIVALENT.
- 4 REFERENCE ITEM 4 FOR POUCH, 10-PACK BOX AND SHIPPING BOX LABEL ARTWORK REQUIREMENTS.
- 5 APPLY POLYMER SOLUTION PART NO. HDX-3000 TO ITEM 9 ON HATCHED AREA AS SHOWN.
- 6 BALLOON VOLUME 5mL.

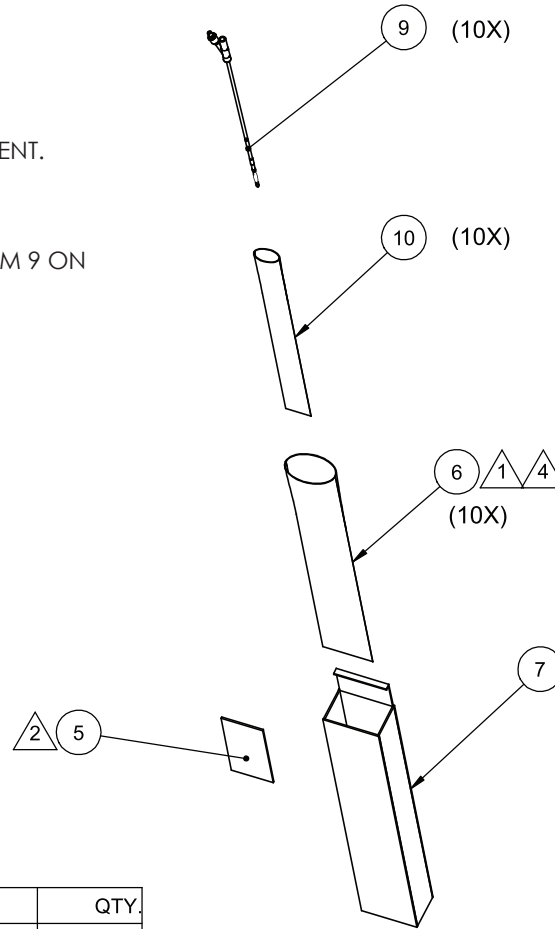

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
N/A	A	INITIAL RELEASE	06/10/2019
N/A	B	ADDED PRODUCTION SPECIFICATIONS ON PG 3.	03/18/2020

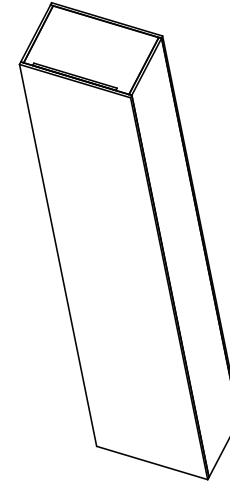

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
4	HDX-4000	ARTWORK, LABEL, HDX FOLEY CATHETER (POUCH AND INNER CARTON)	3
5	HDX-4001	ARTWORK,IFU,HDX FOLEY CATHETER	1
6	HDX-5000	POUCH, CHEVRON, 3.5X21, AMCRO TLP-058	1
7	HDX-5001	BOX, 10-PACK, CATHETER	2
8	HDX-5002	BOX	1
9	HDX-6001	CATHETER, 16FR SILICONE	1
10	HDX-5003	SLEEVE	1
11	HDX-5004	LABEL	3

UNLESS OTHERWISE SPECIFIED:
 DIMENSIONS ARE IN MILLIMETERS
 TOLERANCES:
 ANGULAR: MACH ± 0.2 BEND ± 2.0
 NO PLACE DECIMAL ± 0.5
 ONE PLACE DECIMAL ± 0.3
 TWO PLACE DECIMAL ± 0.15
 THREE PLACE DECIMAL ± 0.05

PROPRIETARY AND CONFIDENTIAL
 THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.





TITLE:
BNS, 16 FR CATHETER

SIZE DWG. NO.
B HDX-8001

REV
B

DO NOT SCALE DRAWING SCALE: 1:25 g SHEET 1 OF 3

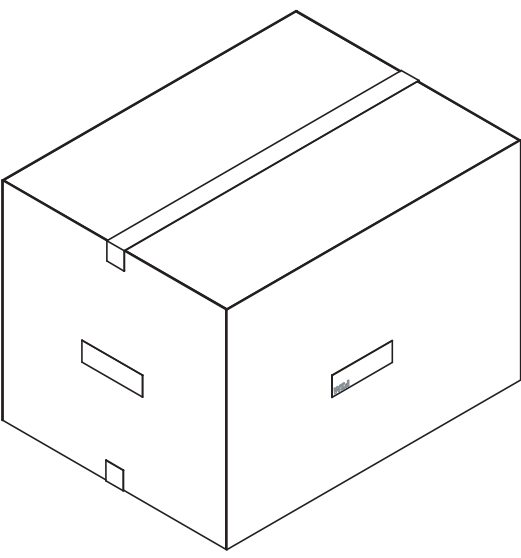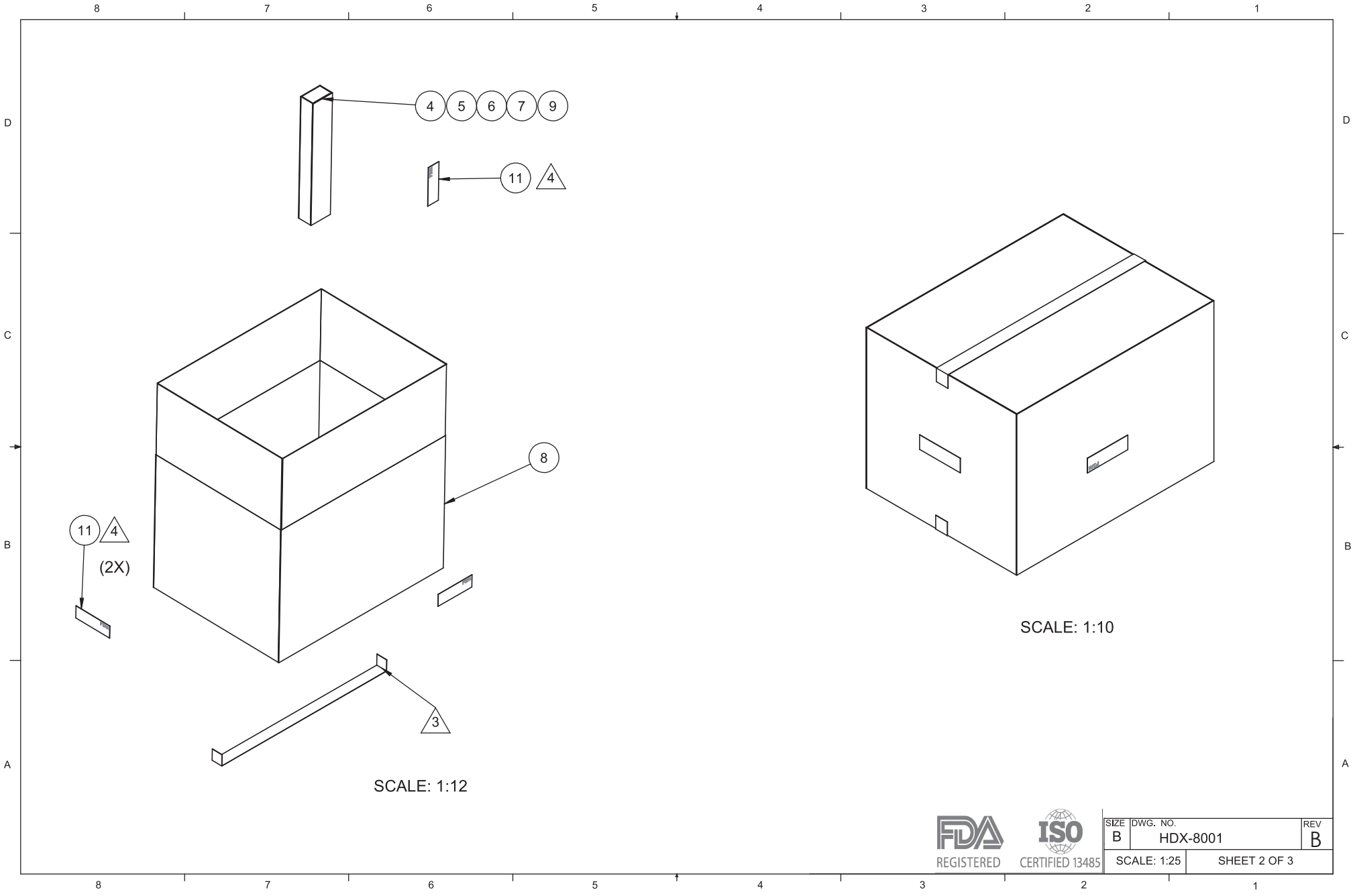

SCALE: 1:10

ITEM NO.	DESCRIPTION	FUNCTION
1	TUBING	CATHETER BODY
2	BALLOON CUFF	RETENTION BALLOON
3	TIP	SMOOTH INSERTION
4	2-WAY FUNNEL	FLUID DRAINAGE
5	CHECK VALVE	BALLOON INFLATION/DEFLATION
6	SLEEVE	SYRINGE INTERFACE

D

C

B


A

D

C

B


A

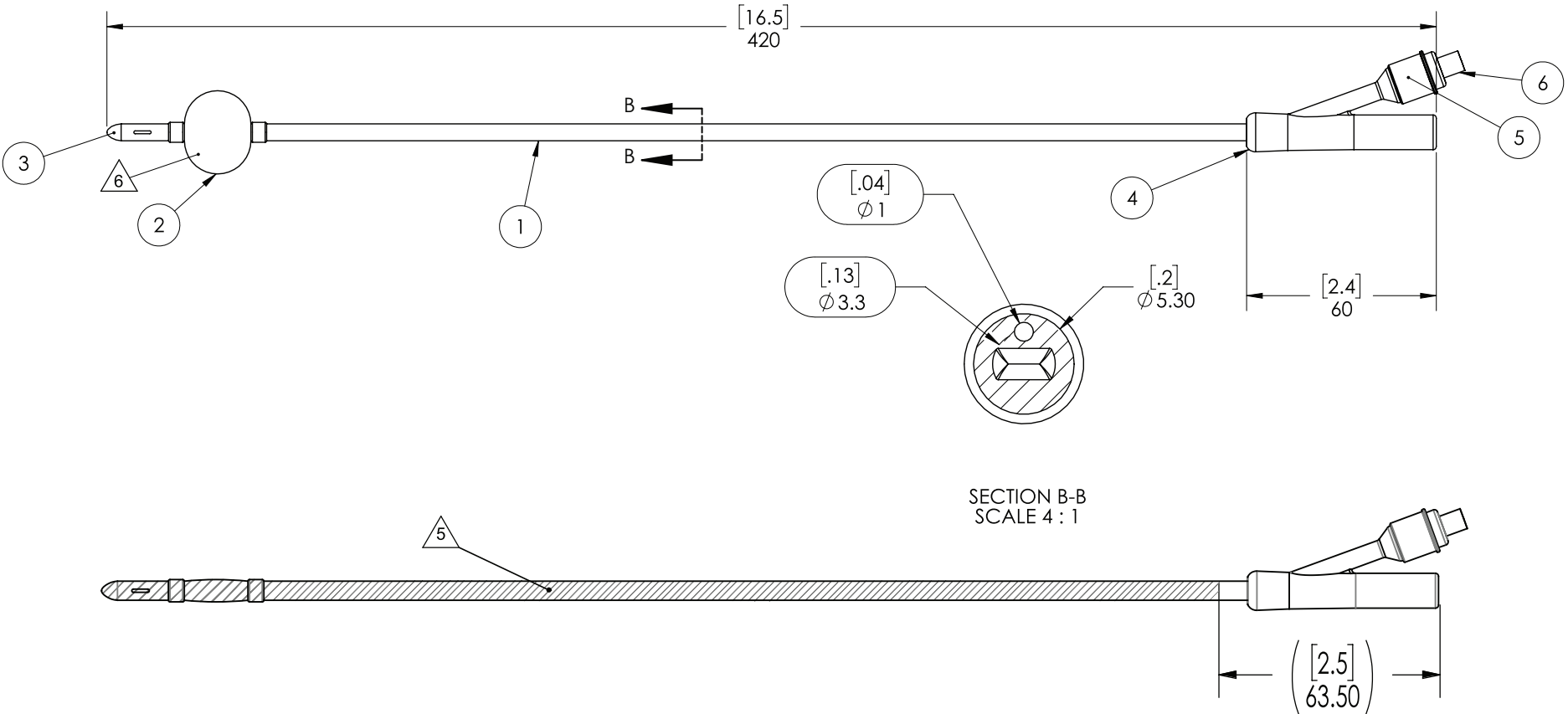


Date Inspected	Part No.				Rev.	Description			
	HDX-8001				B	Hydrophilix 16 Fr x 5ml Catheter, Non-Sterile, 10-Pack			
Lot No.		Lot QTY			Approved Supplier(s)				
Specification	Catheter has been coated with Polymer Solution P/N HDX-3000. MP-HDX-0000	Catheter has been rinsed and free of debris. MP-HDX-0000	Coated and rinsed Catheter meets the requirements of the Spectrophotometer Test. TM-HDX-0000	P/N, Lot No. and Exp date on labels match Work Order requirements. MP-HDX-0000	Catheter is properly oriented in the sleeve and pouch. Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. MP-HDX-0001	Correct quantity in 10-Pack Box. MP-HDX-0001			
Revision		0		0		0			
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A			
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A			
UOM	N/A	N/A	N/A	N/A	N/A	N/A			
Measuring Equipment Type	Visual	Visual	Spectrophotometer	Visual	Visual	Visual			
Asset No.	N/A	N/A		N/A	N/A	N/A			
Cal Due Date	N/A	N/A		N/A	N/A	N/A			
AQL	6.5	6.5	6.5	6.5	6.5	6.5			
Sample Size									
Sample No.	Variable/ Attribute Data								
1									
2									
3									
4									
5									
6									
7									
8									
9									
10									
11									
12									
13									
14									
15									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	N/A	N/A	N/A	N/A	N/A	N/A	0.000	0.000	0.000
Max	N/A	N/A	N/A	N/A	N/A	N/A	0.000	0.000	0.000
Mean	N/A	N/A	N/A	N/A	N/A	N/A	#DIV/0!	#DIV/0!	#DIV/0!
StDev	N/A	N/A	N/A	N/A	N/A	N/A	#DIV/0!	#DIV/0!	#DIV/0!
Comments									
Disposition:						NCMR No.:			
Inspected By:						Date:			

NOTES:(UNLESS OTHERWISE SPECIFIED)

- 1 LABEL ARTWORK PRINTED TOWARDS OPEN END OF POUCH. ITEM 6.
- 2 IFU ITEM 5 TO BE INSERTED INTO 10-PACK BOX ITEM 7.
- 3 SEAL SHIPPING BOX WITH PACKING TAPE OR EQUIVALENT.
- 4 REFERENCE ITEM 4 FOR POUCH, 10-PACK BOX AND SHIPPING BOX LABEL ARTWORK REQUIREMENTS.
- 5 APPLY POLYMER SOLUTION PART NO. HDX-3000 TO ITEM 9 ON HATCHED AREA AS SHOWN.
- 6 BALLOON VOLUME 10mL.

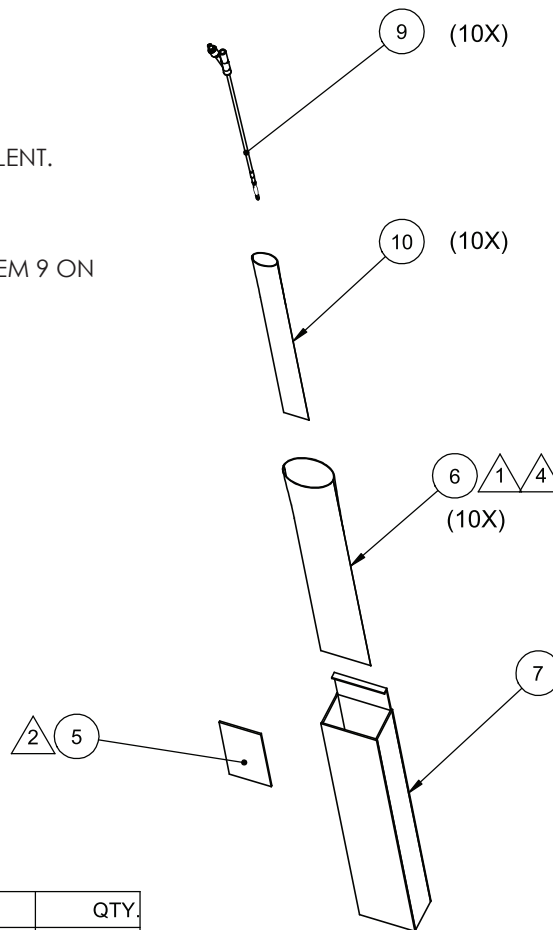

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
N/A	A	INITIAL RELEASE	06/10/2019
N/A	B	ADDED PRODUCTION SPECIFICATIONS ON PG 3.	03/18/2020

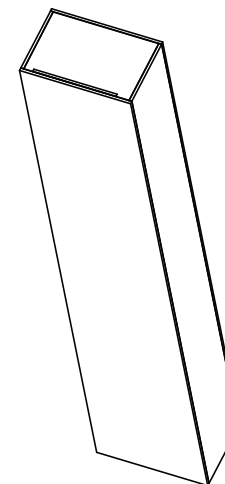

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
4	HDX-4000	ARTWORK, LABEL, HDX FOLEY CATHETER (POUCH AND INNER CARTON)	3
5	HDX-4001	ARTWORK,IFU,HDX FOLEY CATHETER	1
6	HDX-5000	POUCH, CHEVRON, 3.5X21, AMCRO TLP-058	1
7	HDX-5001	BOX, 10-PACK, CATHETER	2
8	HDX-5002	BOX	1
9	HDX-6002	CATHETER, 18FR SILICONE	1
10	HDX-5003	SLEEVE	1
11	HDX-5004	LABEL	3

UNLESS OTHERWISE SPECIFIED:
 DIMENSIONS ARE IN MILLIMETERS
 TOLERANCES:
 ANGULAR: MACH ± 0.2 BEND ± 2.0
 NO PLACE DECIMAL ± 0.5
 ONE PLACE DECIMAL ± 0.3
 TWO PLACE DECIMAL ± 0.15
 THREE PLACE DECIMAL ± 0.05

PROPRIETARY AND CONFIDENTIAL
 THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.




TITLE:
BNS, 18 FR CATHETER

SIZE DWG. NO.
B HDX-8002

REV
B

DO NOT SCALE DRAWING SCALE: 1:25 g SHEET 1 OF 3

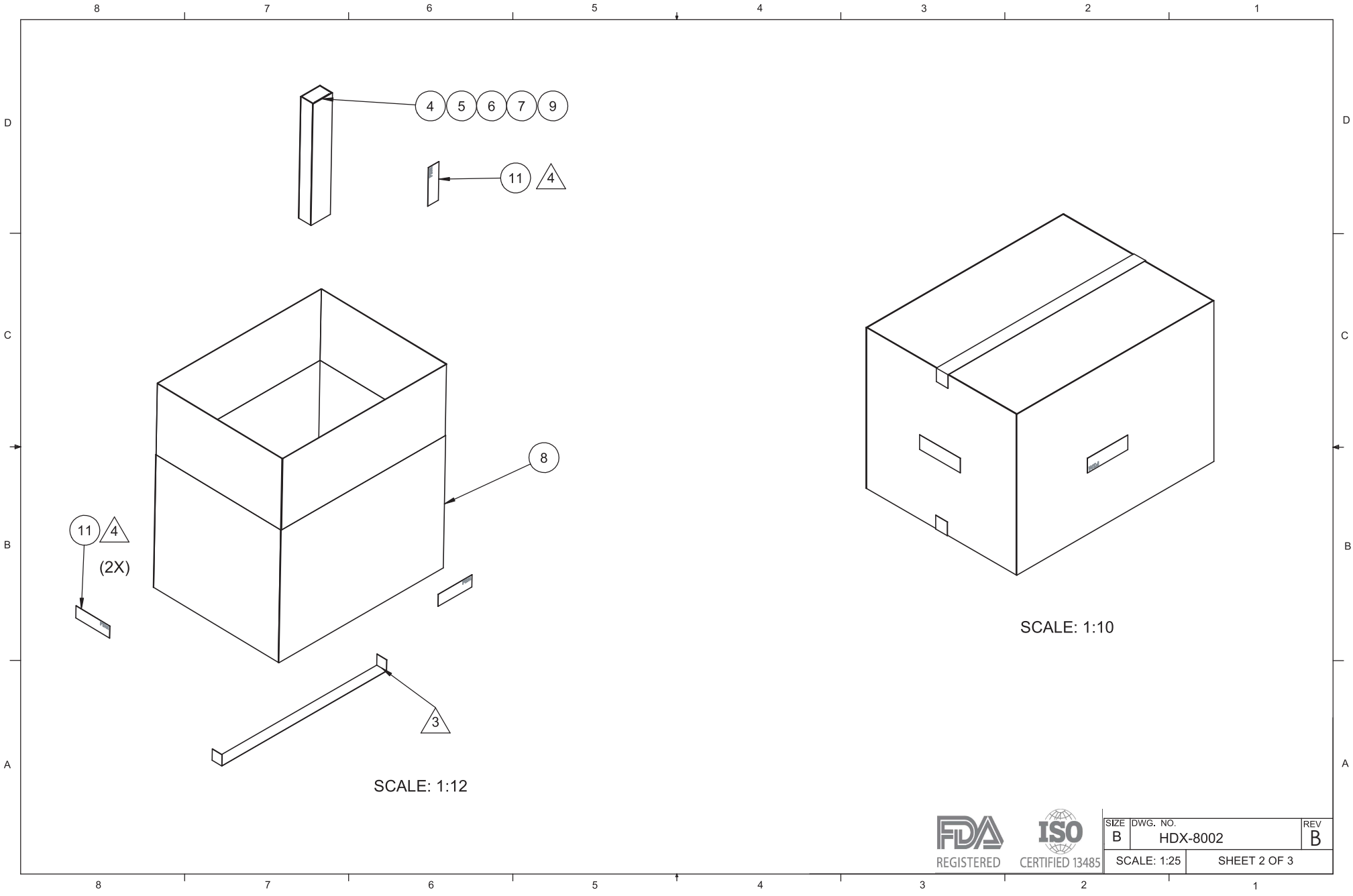

SIZE	DWG. NO.	REV
B	HDX-8002	B
SCALE: 1:25		SHEET 2 OF 3

ITEM NO.	DESCRIPTION	FUNCTION
1	TUBING	CATHETER BODY
2	BALLOON CUFF	RETENTION BALLOON
3	TIP	SMOOTH INSERTION
4	2-WAY FUNNEL	FLUID DRAINAGE
5	CHECK VALVE	BALLOON INFLATION/DEFLATION
6	SLEEVE	SYRINGE INTERFACE

D

C

B


A

D

C

B

A

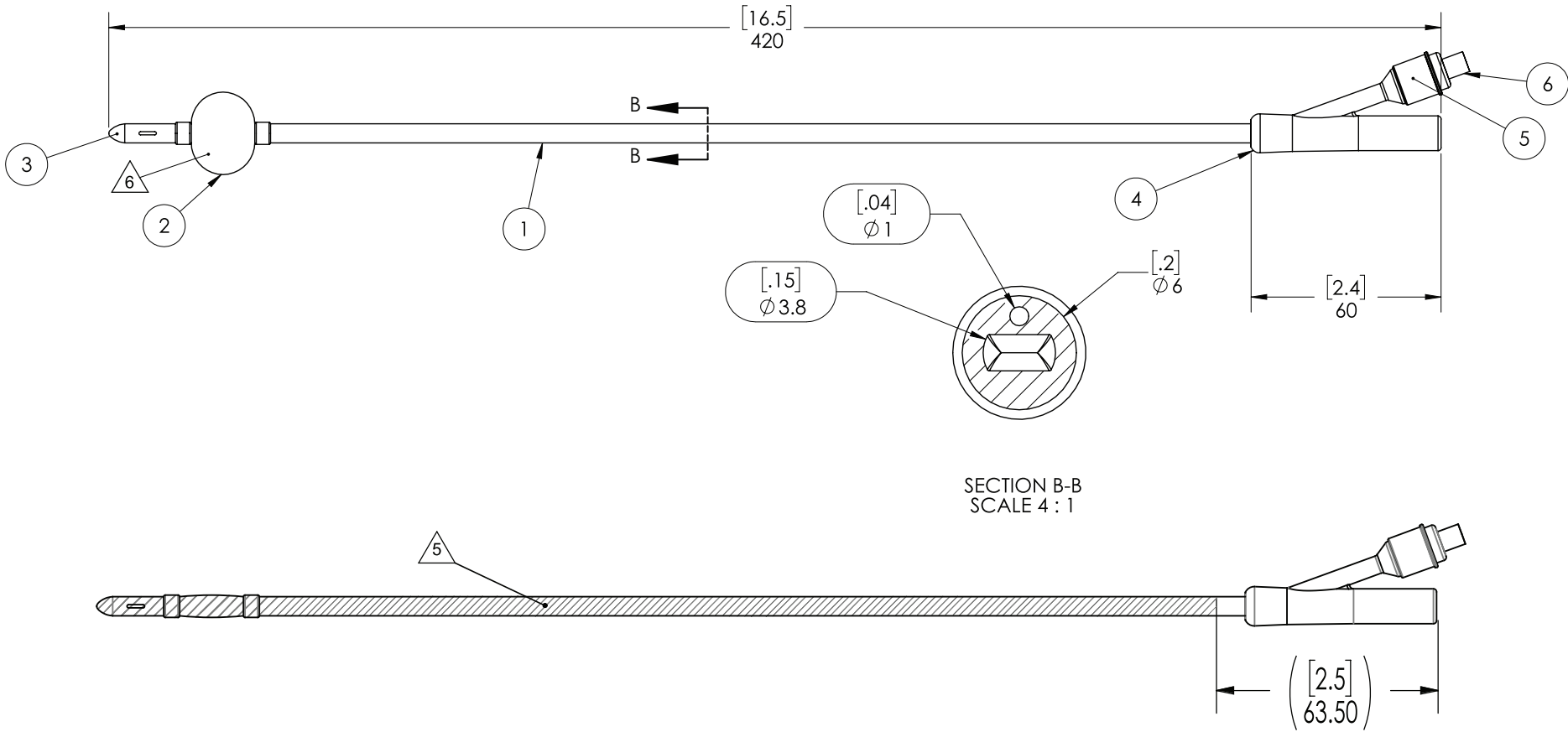

Date Inspected	Part No.				Rev.	Description			
	HDX-8002				B	Hydrophilix 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack			
Lot No.		Lot QTY			Approved Supplier(s)				
Specification	Catheter has been coated with Polymer Solution P/N HDX-3000. MP-HDX-0000	Catheter has been rinsed and free of debris. MP-HDX-0000	Coated and rinsed Catheter meets the requirements of the Spectrophotometer Test. TM-HDX-0000	P/N, Lot No. and Exp date on labels match Work Order requirements. MP-HDX-0000	Catheter is properly oriented in the sleeve and pouch. Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. MP-HDX-0001	Correct quantity in 10-Pack Box. MP-HDX-0001			
Revision		0		0		0			
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A			
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A			
UOM	N/A	N/A	N/A	N/A	N/A	N/A			
Measuring Equipment Type	Visual	Visual	Spectrophotometer	Visual	Visual	Visual			
Asset No.	N/A	N/A		N/A	N/A	N/A			
Cal Due Date	N/A	N/A		N/A	N/A	N/A			
AQL	6.5	6.5	6.5	6.5	6.5	6.5			
Sample Size									
Sample No.	Variable/ Attribute Data								
1									
2									
3									
4									
5									
6									
7									
8									
9									
10									
11									
12									
13									
14									
15									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	N/A	N/A	N/A	N/A	N/A	N/A	0.000	0.000	0.000
Max	N/A	N/A	N/A	N/A	N/A	N/A	0.000	0.000	0.000
Mean	N/A	N/A	N/A	N/A	N/A	N/A	#DIV/0!	#DIV/0!	#####
StDev	N/A	N/A	N/A	N/A	N/A	N/A	#DIV/0!	#DIV/0!	#####
Comments									
Disposition:						NCMR No.:			
Inspected By:						Date:			