

Pathway Innovation LLC.	Template, Design Change Assessment		Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: A
			Page 1 of 2

Project No. / Name	00047 / Hydrophilix 16Fr X 5 – 10 ml with 3 Year Shelf Life Study				
Product Design Specification / Rev	PS-0006 / B			Date	11-30-2021
Change Type	Product: ☒ Line Extension ☐ Change	Manufacturing: ☐ Change	Rationale	1) This is a product line extension by adding a version of the catheter with a 5 – 10 ml balloon. This line extension is not a material change but a change to the volume labeled on the balloon from 5ml to a range 5 – 10 ml. The following catalog number will be added: HDX-160510 2) The shelf life will be extended from 1 year to 3 years.	

Rationale = Describe the justification for the change type and explain the drivers behind the change.

The following table provides an analysis of the potential impact proposed design changes may have with the design requirements described in affected product design specification.

TABLE 1: Change and Requirement Identification						
ID	Type	Component	Change Description	Qualified Part / Rev	Proposed Part / Rev	Requirements Affected
1	A/B/C	Catheter	Balloon volume changed from 5 ml to range of 5 ml – 10 ml. No material changes. The shelf life will be extended to 3 years.	HDX-6001 / A	HDX-6006 / A	#2,3,4 - Labeling #11 - Balloon Integrity #12 - Balloon Volume Maintenance #13 - Non interconnect #14 - Deflation reliability #15 - Flow rate #16 - Balloon Traction #18 - Balloon & Shaft Size (DIMS) #21 - Connection type #22 - No of Lumens #23 - French size offerings #24 – Shelf Life
2	A/C	Box / Label	Change from 16Fr 5ml to 16F 5ml – 10 ml, change from (22x14x12) in. box to (20x18x6) in. box.	HDX-4000 / C HDX-5003 / A	HDX-4002 / A HDX-5012 / A	#2,3,4 – Labeling #26 – Transportation/Packaging
3	A/C	IFU	Change from 16Fr 5ml to 16F 5ml – 10 ml Name change from Hydrophilix to SILQ	HDX-4001 / C	HDX-4006 / A	#2,3,4 – Labeling

Type: A = Form, B = Fit, C = Function, Performance

The following table provides an impact assessment for the proposed design changes described in Table 1.

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Internal / External Resources	None
Products Affected	HDX-160510
Customers Affected	US Market
Distribution Channels	US Market

Pathway Innovation LLC.	Template, Design Change Assessment	Doc No.: T-007
		Rev.: A
		Page 2 of 2

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Physical/ Functional (e.g., Bench Testing) Assessment	Design verification will be performed on modified catheter. The product specification (PS-0006) will be revised to include this report in the trace matrix section.
Shelf-Life Assessment	All materials that construct the 2-way Catheter have been qualified with a twelve (12) month shelf life. Design verification will be performed on all performance aspects to extend the shelf life to 3 years.
Dimensional (e.g., 1 st Article) Assessment	Device proximal tip, balloon and shaft will be confirmed. The product specification (PS-0006) will be revised to include this report in the trace matrix section.
Sterilization Assessment	There is no impact to the sterilization validation performed because there were no sterile barrier or pallet configuration changes that affect sterility. A documented assessment to adopt the product into the validated cycle will be performed.
Biocompatibility Assessment	There is no impact to the biocompatibility evaluations performed because there were no material changes. No further action is required.
Packaging Assessment	There is no impact to the validated packaging configuration. The outer shipping box change only removes the dead space and optimizes the size. This has no impact to the design. Simulated distribution and packaging testing will be performed to validate the extended shelf life. The product specification (PS-0006) will be revised to include this report in the trace matrix section.
Regulatory Assessment	The guidance document “Deciding when to submit a 510(k) for a change to an existing device (1500054)” was used in this assessment, which the attached decision trees were completed. A 510(k) submission is not required. No further action is required.
Quality System Assessment	None

Contains Nonbinding Recommendations

Note that the first question is whether the change is being made with the intent to significantly improve the safety or effectiveness of the device, for example, to significantly improve clinical outcomes, to mitigate a known risk, in response to adverse events, etc. (Figure 1 – Main Flowchart). If so, the change likely could significantly affect safety or effectiveness and submission of a new 510(k) is likely required. If not, you should continue to follow the logic scheme shown in Figure 1, below.

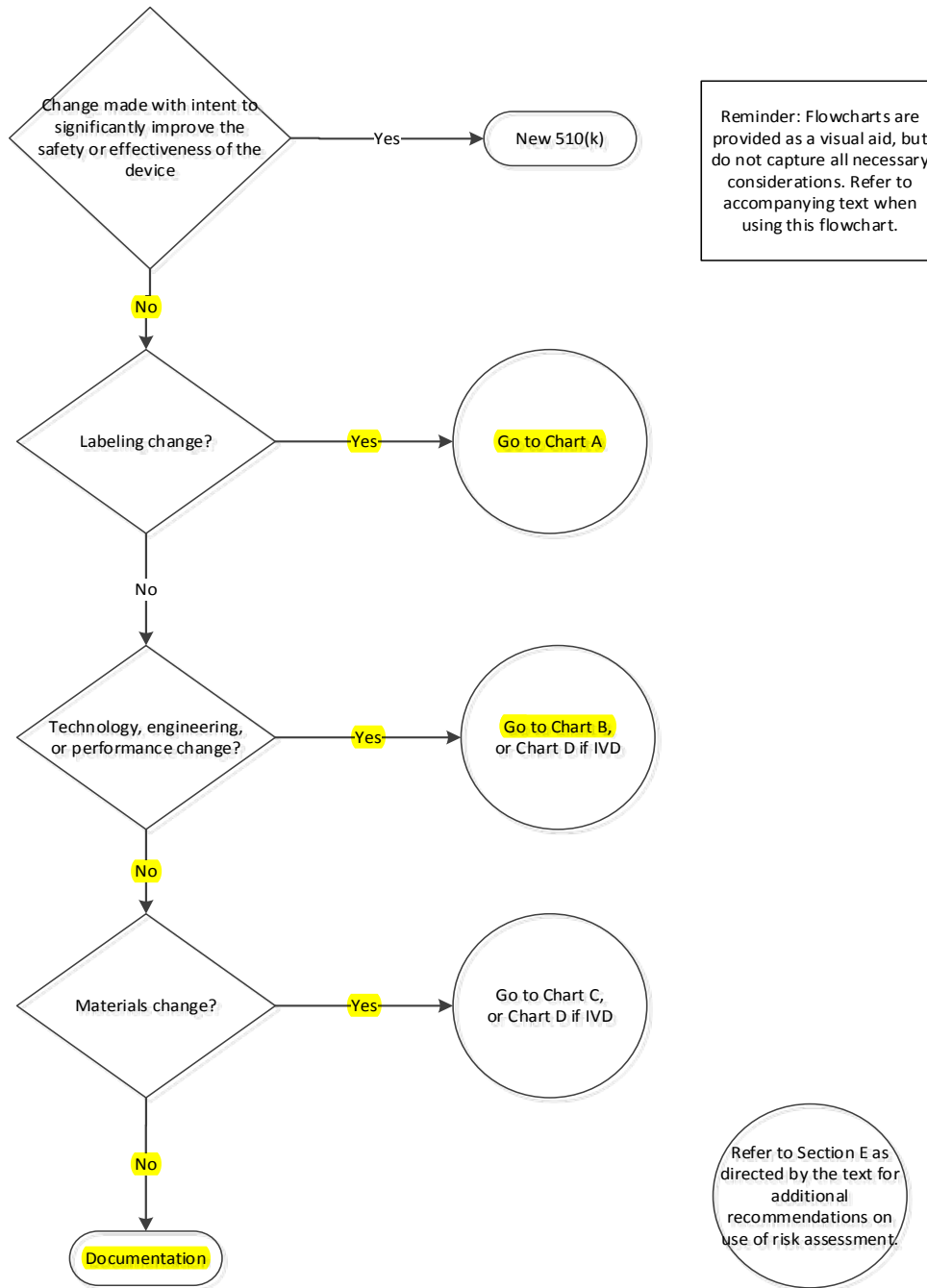


Figure 1 - Main Flowchart

Contains Nonbinding Recommendations

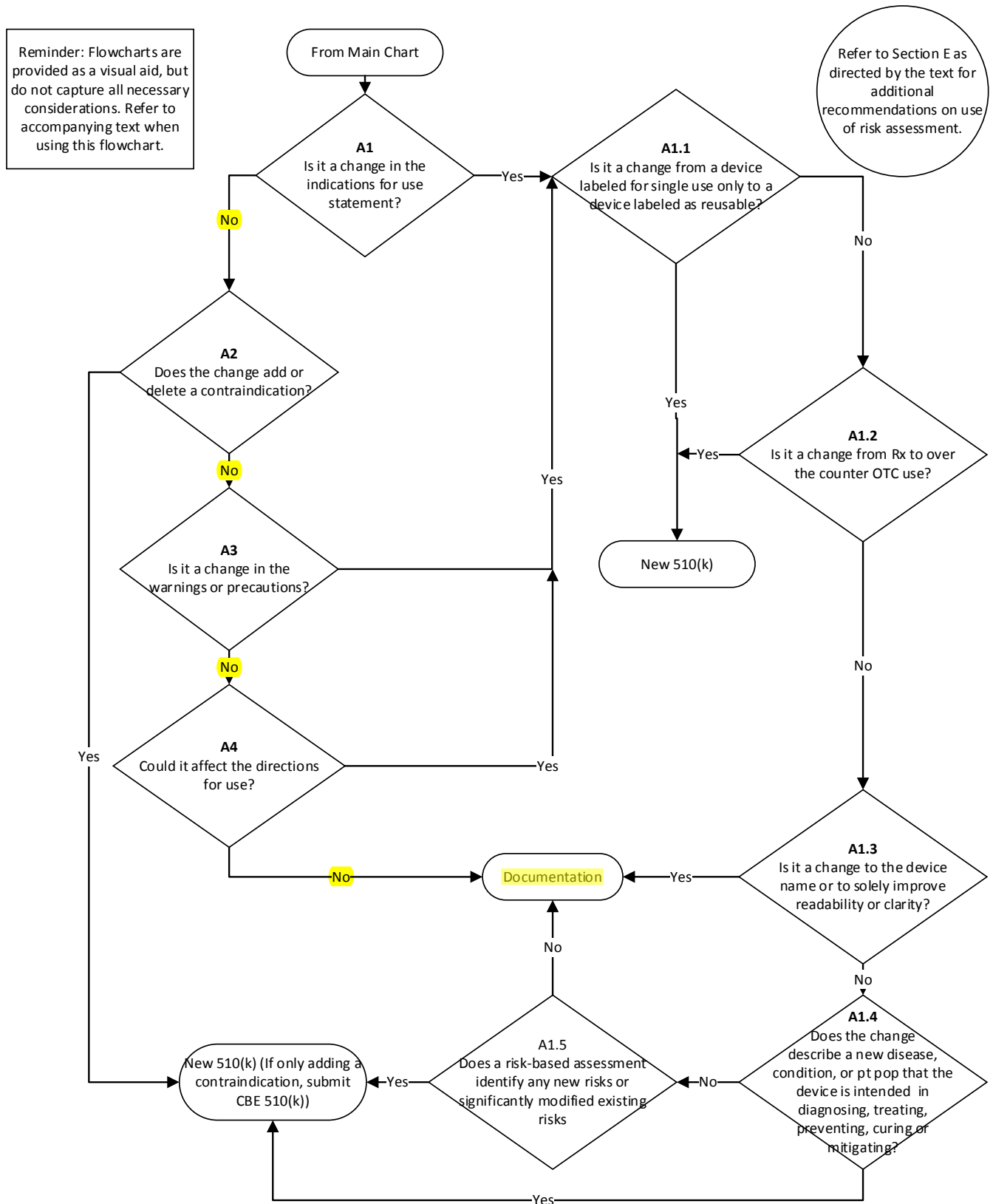


Figure 2 - Flowchart A: Labeling Changes

Contains Nonbinding Recommendations

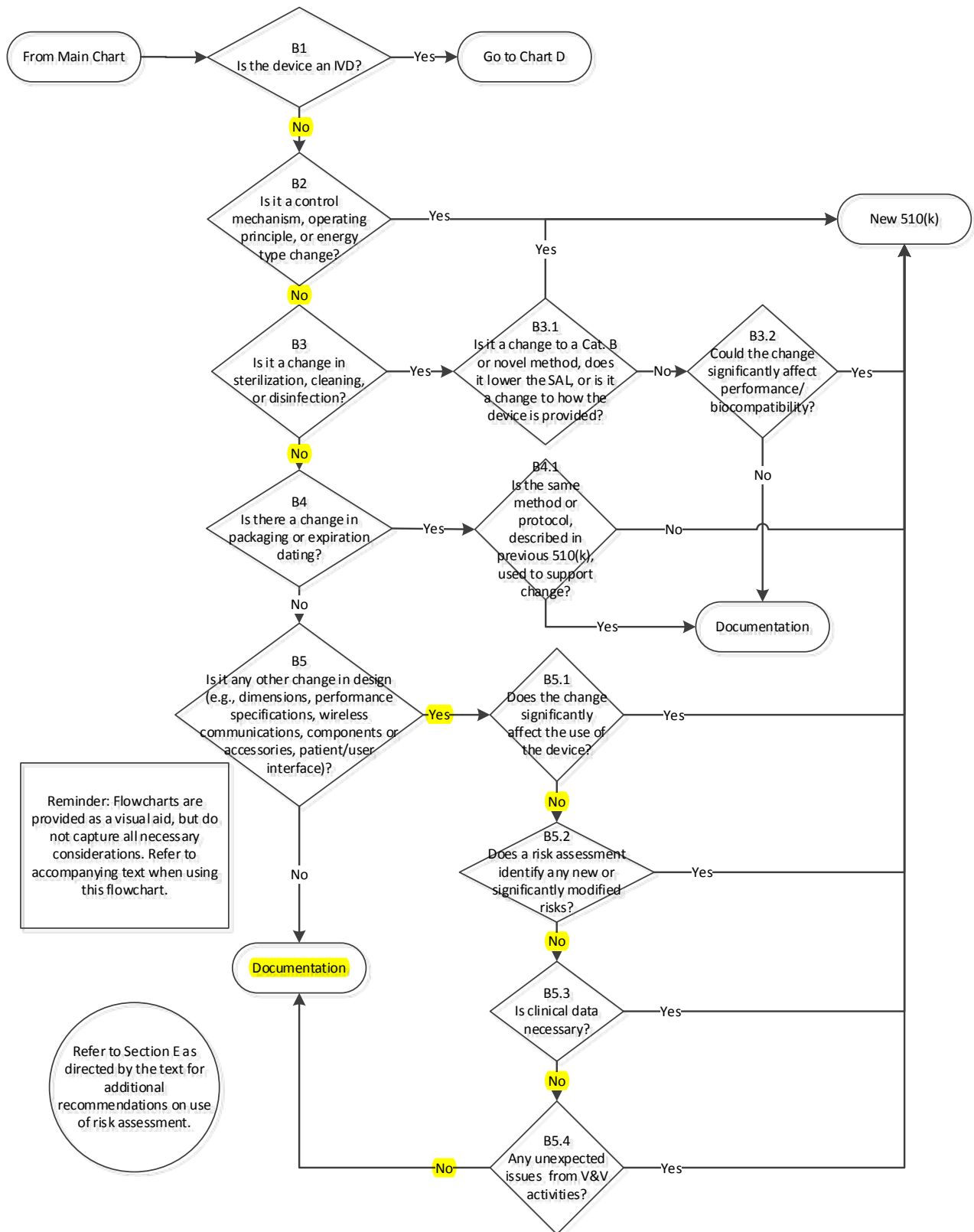


Figure 3 - Flowchart B: Technology, Engineering, and Performance Changes

Contains Nonbinding Recommendations

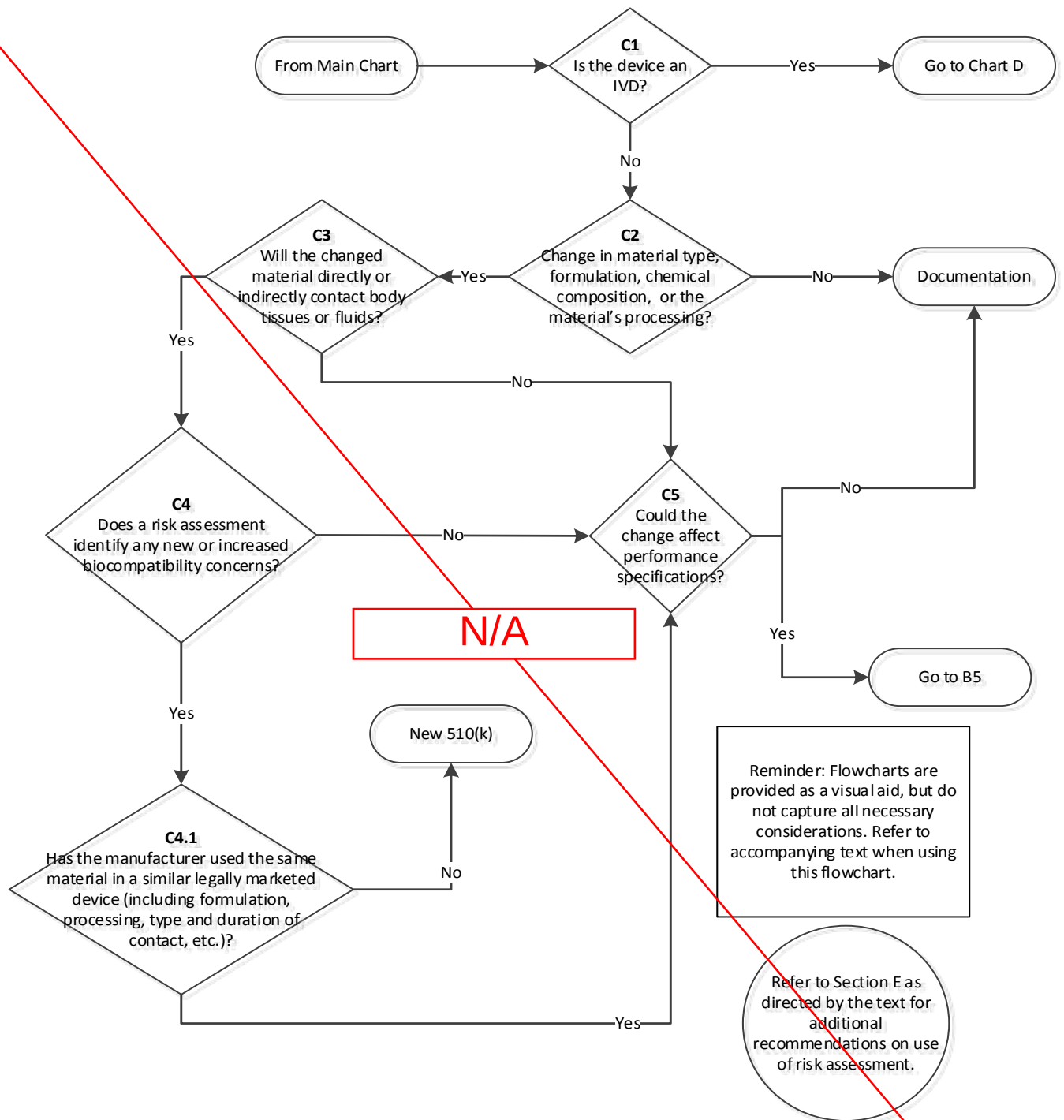


Figure 4 - Flowchart C: Materials Changes