

	Design Review		Doc No.: FM-0052
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: B
			Page 1 of 2

Project No.	Product Name	Part No./ Rev	Review Phase	Review Date
00049	SILQ ClearTract 16 Fr 10ml	SLQ-1610 / A	Line Extension Commercial Release	03/11/2022
Subject of Review				
#	Topic	Details	Reference Documents	
1	Change Assessment / Plan	The document was reviewed and complete.	FILE-008918 03/01/2022	
2	Device Master Record (DMR)	The documents were reviewed and complete. - Assembly Specifications / BOM: SLQ-1610 rev A, SLQ-8004 rev A - Components: SLQ-3000 rev A, SLQ-5000 rev A, SLQ-5001 rev A, SLQ-5002 rev A, SLQ-5004 rev A, SLQ-5005 rev A, SLQ-5006 rev A, SLQ-50011 rev A, SLQ-5012 rev A, SLQ-6004 rev A - Inspection Forms: SLQ-8004, SLQ-6004, SLQ-5001, SLQ-5000, SLQ-3000 - Labeling: SLQ-4000 rev A, SLQ-4001 rev A, SLQ-4006 rev A - Manufacturing Procedures: MP-SLQ-0001 rev A, MP-SLQ-0002 rev A	Documents are noted in "Details"	
3	Product Performance	No product performance testing was required. The 1610 is equivalent to 160510.	Refer to Change Assessment #1	
4	Sterilization Adoption	There is no impact to the sterilization validation performed because there were no sterile barrier or pallet configuration changes that affect sterility. No further action is required. The 1610 is equivalent to 160510.	Refer to Change Assessment #1	
5	Labeling Content Review	The previous labeling content review applies. The only changes to the labeling were renumbering their document numbers from HDX-4XXX to SLQ-4XXX. There were no content changes.	FILE-008308 12/01/2021, ECO-000899	
6	Risk Management Plan	The document was reviewed and up to date.	RM-0018 rev C	
7	Product Specification/Trace Matrix	No new requirements. Document is up to date. Trace matrix is up to date and complete.	PS-0006 rev D	
8	Risk – Hazard Analysis	The document was reviewed and up to date. No new risks discovered. Complete.	RM-0019 rev C	
9	Risk – DFMEA	The document was reviewed and up to date. No new risks discovered.	PS-0020 rev C	
10	Process Validation	The documents were reviewed. - MVP: TP-HDX-0007 rev B, TR-HDX-0007 rev A. - Mfg Line Evaluation – TP-HDX-0008 rev B, TR-HDX-0008 rev B Lot 2 and 3 were completed and report updated.	Documents are noted in "Details"	
11	Risk – PFMEA	The document was reviewed and up to date. No new risks discovered.	RM-0021 rev C	
12	Market Clearance	Letter to File Memo was reviewed and complete.	FILE-008919 03/02/2022	

	Design Review		Doc No.: FM-0052
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: B
			Page 2 of 2

Project No.	Product Name		Part No./ Rev	Review Phase	Review Date
00049	SILQ ClearTract 16 Fr 10ml		SLQ-1610 / A	Line Extension Commercial Release	03/11/2022
13	Review of Line Extension Deliverables	Each document was reviewed and discussed. Project team in attendance has concluded that the deliverables are complete and satisfy the Design Control requirements in accordance with QP-007 rev F for this design change (line extension). It is agreed upon to release the 1610 for commercial distribution and close the project design history file..		Refer to #1-#13, QP-007 rev F.	
Action Items					
#	Description			Assigned To	Due Date
N/A	N/A			N/A	N/A
Reviewers					
Functional Group	Required?	Name	Title	Signature	Date
Customer - SILQ Technologies Corp	☒	Brian McVerry	CTO		
Project Manager	☒	Michael Klumph	Project Manager		
Research & Development / Manufacturing	☒	Arthur Deptala	CEO, Managing Partner		
Quality Assurance / Regulatory Affairs	☒	Aaron Rogers	Director of Regulatory and Quality		
Quality Assurance	☒	Larry Dacoron	Sr. Quality Engineer		
Independent Reviewer	☒	Karim Natafji	Project Manager (not directly responsible for this project)		

Pathway Innovation LLC.	Template, Design Change Assessment		Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: A
			Page 1 of 2

Project No. / Name	00049 / SILQ 16Fr X 10 ml Product Introduction (US)				
Product Design Specification / Rev	PS-0006 / C			Date	03-01-2022
Change Type	Product: ☒ Line Extension ☐ Change	Manufacturing: ☐ Change	Rationale	1) This is a product line extension by adding a version of the catheter with a 10 ml balloon. This line extension is not a material change but a change to the volume labeled on the balloon of 10 ml. The following catalog number will be added: SLQ-1610	

Rationale = Describe the justification for the change type and explain the drivers behind the change.

The following table provides an analysis of the potential impact proposed design changes may have with the design requirements described in affected product design specification.

TABLE 1: Change and Requirement Identification						
ID	Type	Component	Change Description	Qualified Part / Rev	Proposed Part / Rev	Requirements Affected
1	A/B/C	Catheter	Balloon volume changed from a 5 ml – 10 ml range to 10 ml. No material changes.	HDX-6006 / A	SLQ-6004 / A	#2,3,4 - Labeling #11 - Balloon Integrity #12 - Balloon Volume Maintenance #13 - Non interconnect #14 - Deflation reliability #15 - Flow rate #16 - Balloon Traction #18 - Balloon & Shaft Size (DIMS) #21 - Connection type #22 - No of Lumens #23 - French size offerings #24 – Shelf Life
2	A/C	Box / Label	Change from 16Fr 5ml – 10 ml to 16F 10 ml	HDX-4002 / A	SLQ-4000 / A	#2,3,4 – Labeling
3	A/C	IFU	Change from 16Fr 5ml – 10 ml to 16F 10 ml	HDX-4006 / A	SLQ-4006 / A	#2,3,4 – Labeling

Type: A = Form, B = Fit, C = Function, Performance

The following table provides an impact assessment for the proposed design changes described in Table 1.

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Internal / External Resources	None
Products Affected	SLQ-1610
Customers Affected	US Market
Distribution Channels	US Market
Physical/ Functional (e.g., Bench Testing) Assessment	Design verification was performed on 16Fr 5 ml -10 ml catheter which is equivalent to the catheter labeled as 16Fr 10 ml. There is no impact to the design verification performed because the proposed model performs within the tested specification. The product specification (PS-0006) will be revised to include this model in the trace matrix section. No further action is required.

Pathway Innovation LLC.	Template, Design Change Assessment		Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: A
			Page 2 of 2

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Shelf-Life Assessment	All materials that construct the 2-way Catheter have been qualified with a twelve (12) month shelf life. No further action is required.
Dimensional (e.g., 1 st Article) Assessment	Device proximal tip, balloon and shaft was confirmed on 16Fr 5 ml and 16Fr 5 ml -10 ml catheter which is equivalent to the catheter labeled as 16Fr 10 ml. There is no impact to the Device proximal tip, balloon and shaft because the proposed model performs within the tested specification. The product specification (PS-0006) will be revised to include this model in the trace matrix section. No further action is required.
Sterilization Assessment	There is no impact to the sterilization validation performed because there were no sterile barrier or pallet configuration changes that affect sterility. No further action is required.
Biocompatibility Assessment	There is no impact to the biocompatibility evaluations performed because there were no material changes. No further action is required.
Packaging Assessment	There is no impact to the validated packaging configuration. No further action is required.
Regulatory Assessment	The guidance document “Deciding when to submit a 510(k) for a change to an existing device (1500054)” was used in this assessment, which the attached decision trees were completed. A 510(k) submission is not required. No further action is required.
Quality System Assessment	None

Contains Nonbinding Recommendations

Note that the first question is whether the change is being made with the intent to significantly improve the safety or effectiveness of the device, for example, to significantly improve clinical outcomes, to mitigate a known risk, in response to adverse events, etc. (Figure 1 – Main Flowchart). If so, the change likely could significantly affect safety or effectiveness and submission of a new 510(k) is likely required. If not, you should continue to follow the logic scheme shown in Figure 1, below.

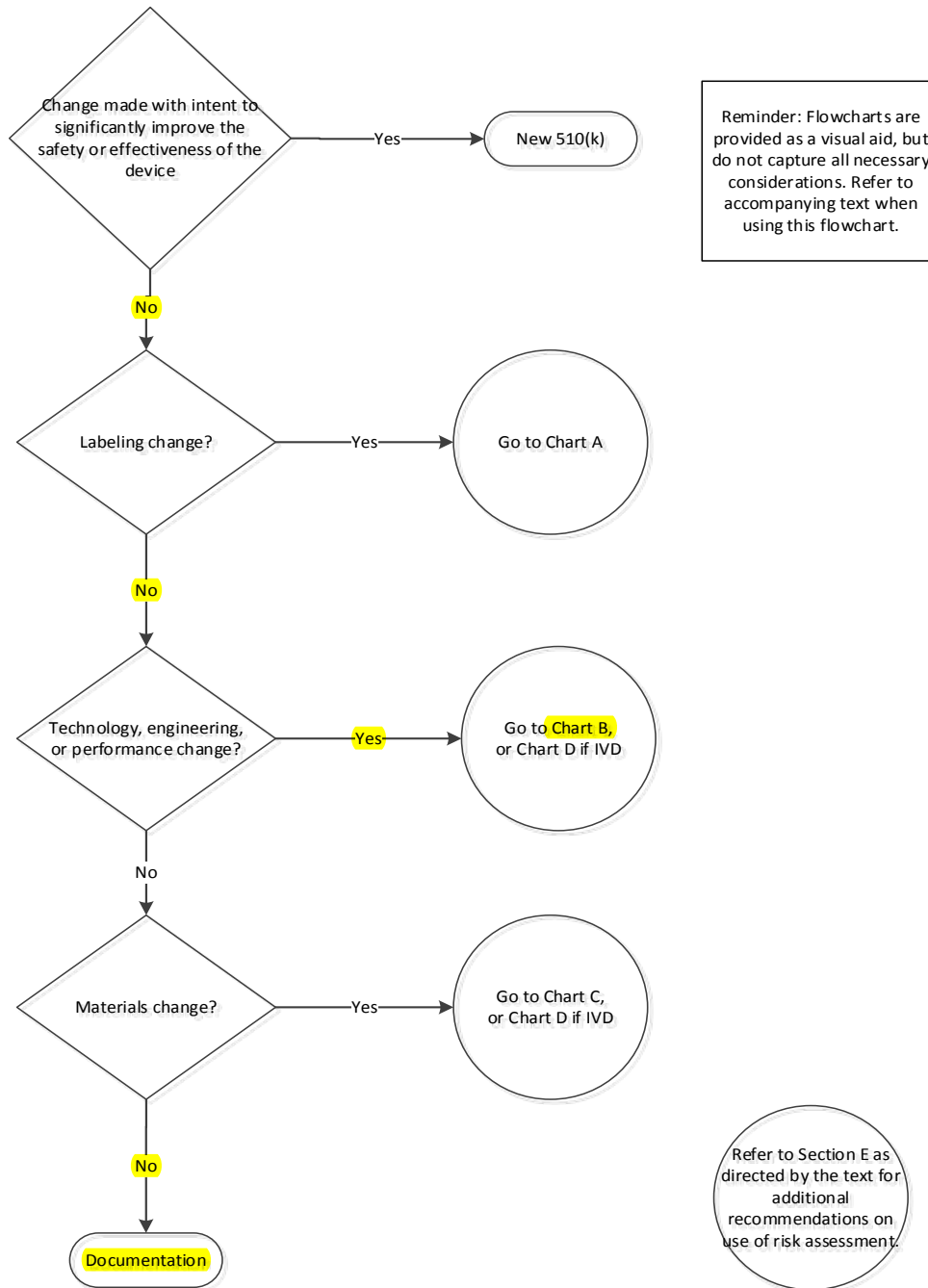


Figure 1 - Main Flowchart

Contains Nonbinding Recommendations

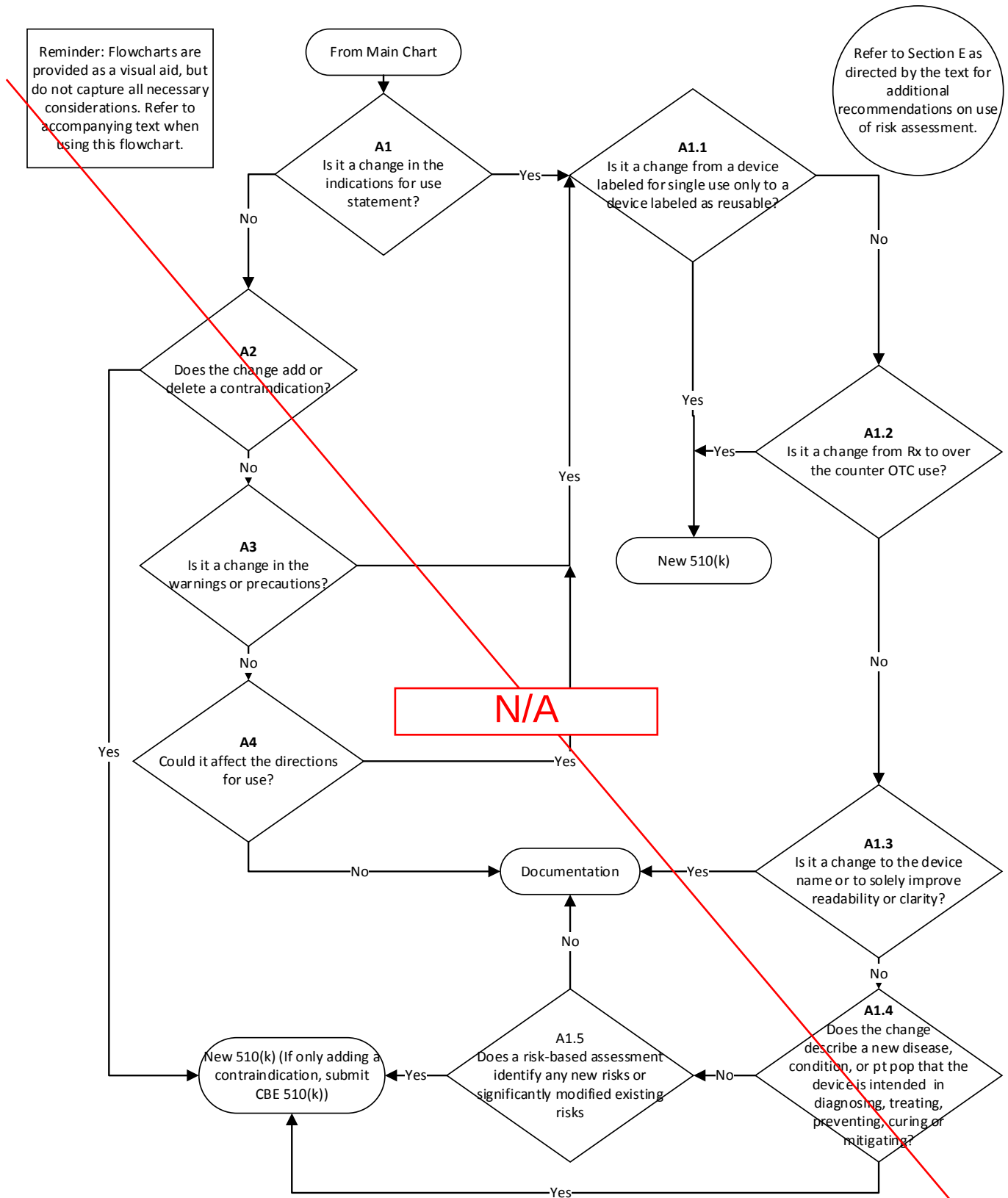


Figure 2 - Flowchart A: Labeling Changes

Contains Nonbinding Recommendations

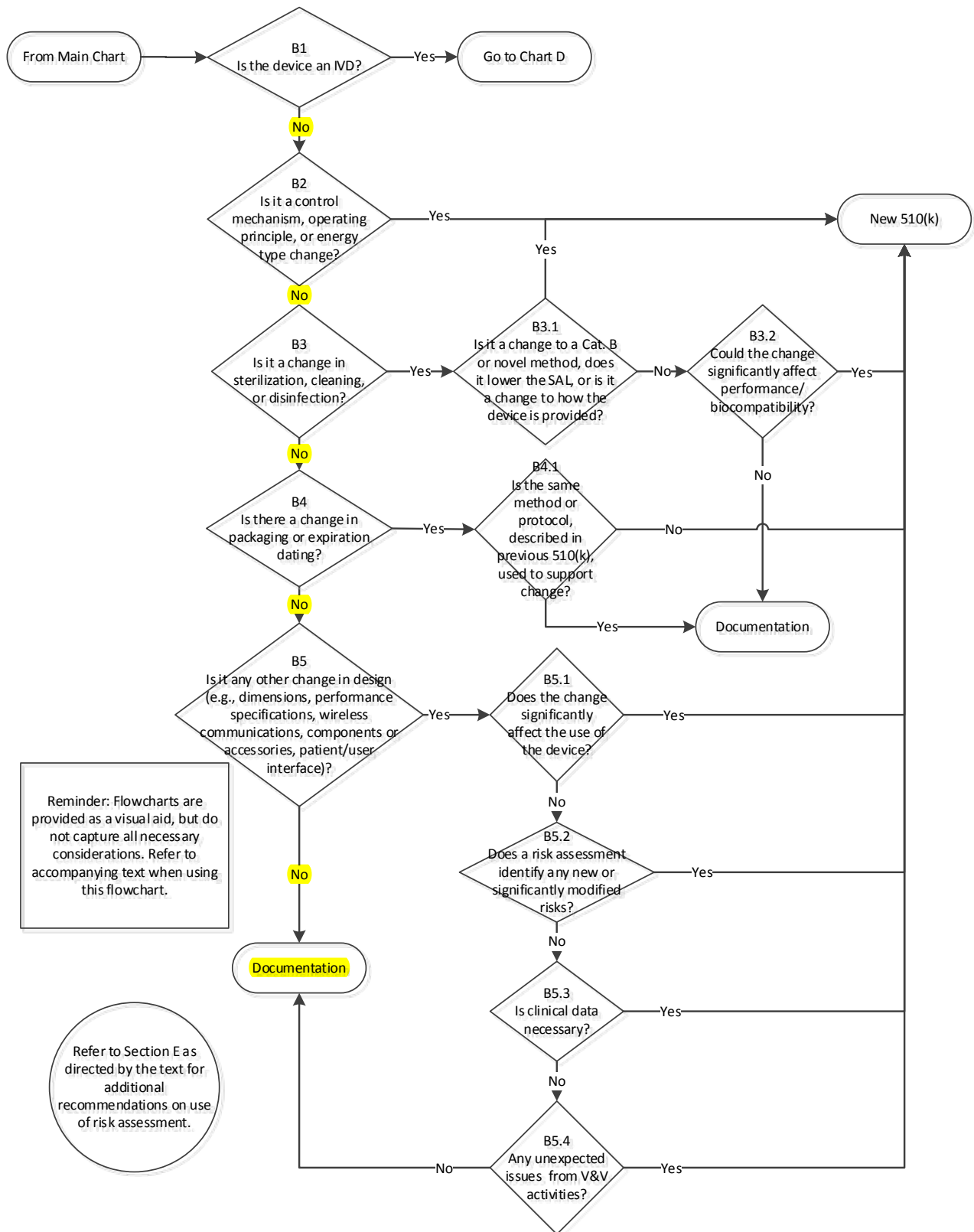


Figure 3 - Flowchart B: Technology, Engineering, and Performance Changes

Contains Nonbinding Recommendations


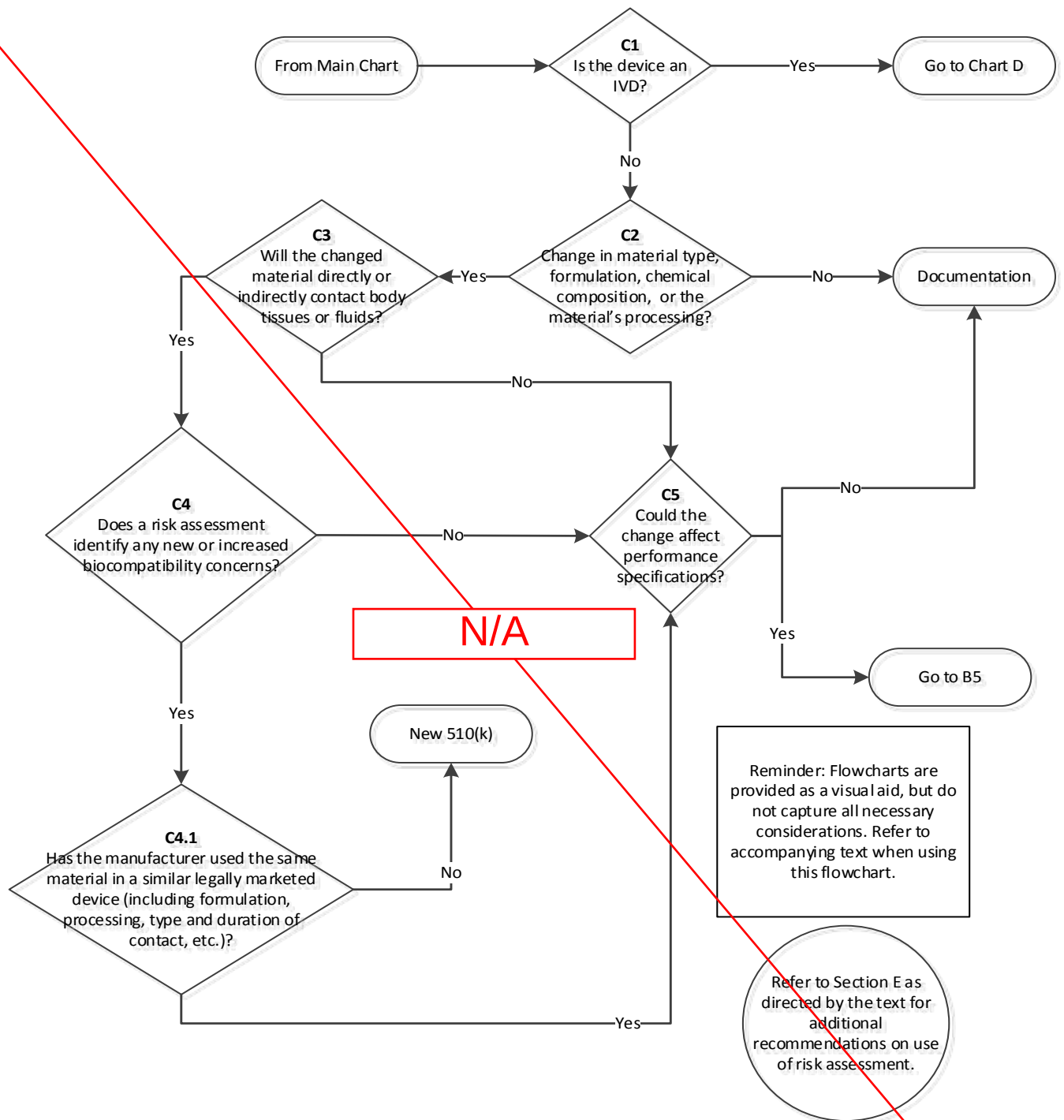


Figure 4 - Flowchart C: Materials Changes

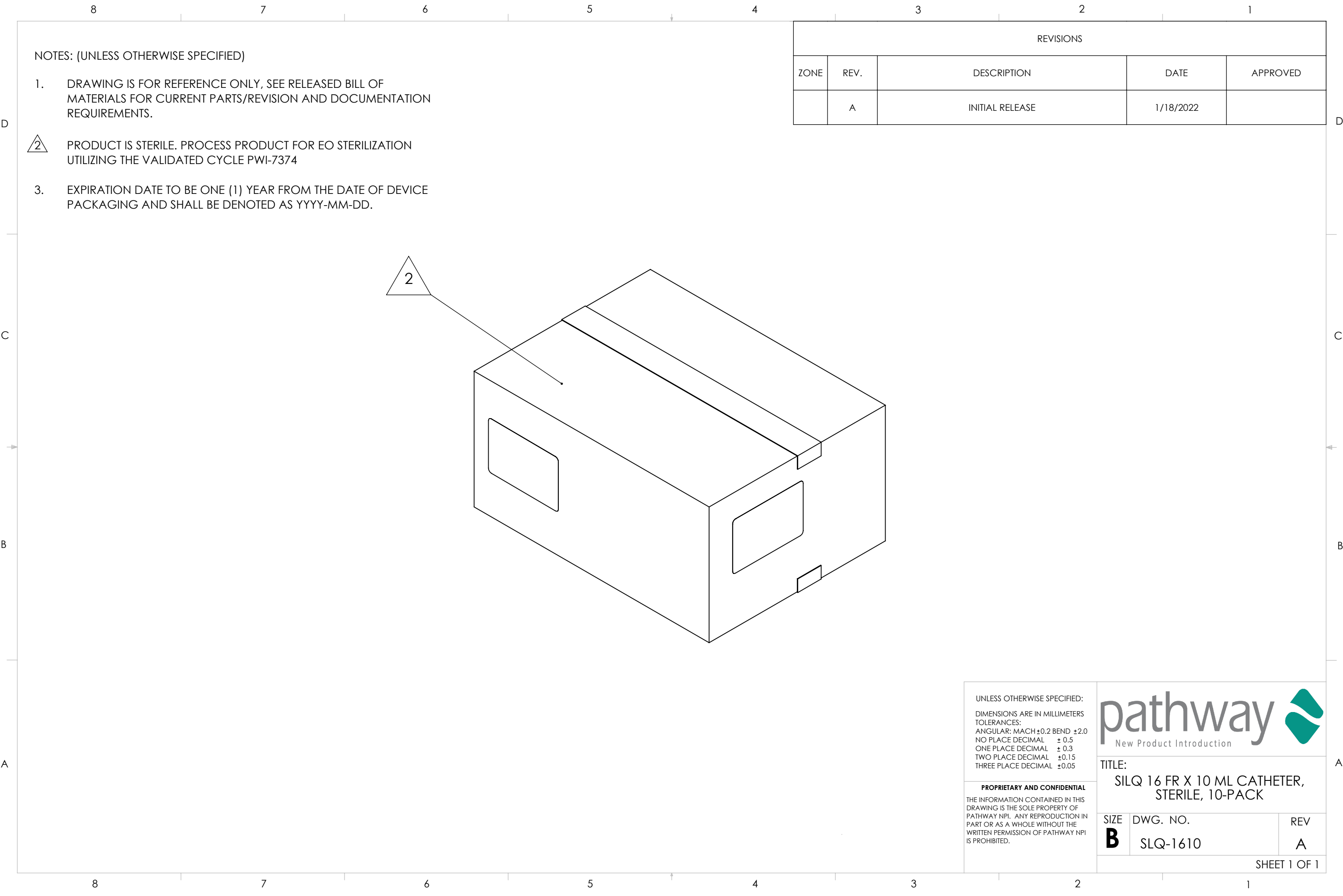

- NOTES: (UNLESS OTHERWISE SPECIFIED)
- 1. DRAWING IS FOR REFERENCE ONLY, SEE RELEASED BILL OF MATERIALS FOR CURRENT PARTS/REVISION AND DOCUMENTATION REQUIREMENTS.
 - 2. PRODUCT IS STERILE. PROCESS PRODUCT FOR EO STERILIZATION UTILIZING THE VALIDATED CYCLE PWI-7374
 - 3. EXPIRATION DATE TO BE ONE (1) YEAR FROM THE DATE OF DEVICE PACKAGING AND SHALL BE DENOTED AS YYYY-MM-DD.

REVISIONS				
ZONE	REV.	DESCRIPTION	DATE	APPROVED
	A	INITIAL RELEASE	1/18/2022	

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.



TITLE:
SILQ 16 FR X 10 ML CATHETER,
STERILE, 10-PACK

SIZE	DWG. NO.	REV
B	SLQ-1610	A

Pathway LLC

Printed by Aaron Rogers
Printed on March 01, 2022
Local time zone (GMT-08:00) Pacific Time (America/Los_Angeles)
Page Arena > Workspace Items > Bill of Materials > Indented
All information contained in this document is proprietary and confidential.

SLQ-1610

In Production

SILQ 2-way Cleartract 16 Fr x 10 ml Catheter, Sterile, 10-Pack

Revision »A - ECO-000907 Status: Effective, Status Date: 02/02/2022 04:50:37 PM, Unshared

Contains 4 first-level Items, 34 line Items, 33 unique Items, 2 of which are shared.

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	      	
1 FM-0029 rev A	Form, Certificate of Conformance	Form	In Prod	2	1 each
2 QP-010 rev F	Final Product Inspection	Quality System Procedure	In Prod	6	1 each
3 QP-032 rev C	EO Sterilization Process Control	Quality System Procedure	In Prod	4	1 each
4 SLQ-8004 rev A	SILQ 2-way Cleartract 16 Fr x 10 ml Catheter, Non-Sterile, 10-Pack	Final Assembly	In Prod	4	1 EA
1 SLQ-6004 rev A	Catheter, 16 Fr x 10 mL, Silicone, NY Foley, 2-Way Fabricated Part		In Prod	3	10 EA
2 SLQ-3000 rev A	Solution, SILQ Polymer	Purchased Component	In Prod	4	80 ML
SUBSTITUTES					
	RANK  ITEM NUMBER	ITEM NAME		NOTES	
	1 HDX-3000 rev A	Solution, Hydrophilix Polymer	80 ML		
3 SLQ-5005 rev A	Ribbon, Thermal Transfer, Resin/Wax, 4.3" W	Packaging	In Prod	2	91.2 IN
4 SLQ-5004 rev A	Label, 2.3125 in. X 7.50 in.	Packaging	In Prod	2	12 EA
5 SLQ-4000 rev A	Artwork, Label, Pouch/Box, SILQ 2-Way ClearTract Catheter	Labeling	In Prod	1	12 EA
6 SLQ-5011 rev A	Label, 4 in. x 6 in., Standard White Thermal Transfer	Packaging	In Prod	2	0.2 EA
7 SLQ-4001 rev A	Artwork, Label, Box, SILQ 2-Way ClearTract Catheter	Labeling	In Prod	1	0.2 EA

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE			
8 SLQ-4006 rev A	Artwork, IFU, SILQ 2-Way ClearTract Catheter	Labeling	In Prod	7	1 EA	
9 SLQ-5000 rev A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	Packaging	In Prod	3	10 EA	
10 SLQ-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	Packaging	In Prod	4	10 EA	
SUBSTITUTES						
	RANK  ITEM NUMBER	ITEM NAME 	NOTES			
	1 HDX-5001 rev A	Pouch, Chevron, 3.5 10 EA X 21, Amcor TLP-058				
	2 PSM-5000 rev A	Pouch, Chevron, 3.5 10 EA X 21, Amcor TLP-058				
11 SLQ-5002 rev A	Box, 10-Pack, Inner, 20.75 L x 4.12 W x 2.37 H, PSM	Packaging	In Prod	2	1 EA	
SUBSTITUTES						
	RANK  ITEM NUMBER	ITEM NAME 	NOTES			
	1 HDX-5002 rev A	Box, 10-Pack, Inner, 1 EA 20.75 L x 4.12 W x 2.37 H, PSM				
12 SLQ-5012 rev A	Box, 20 in. L x 18 in. W x 6 in. D, 200 LB Test	Packaging	In Prod	2	0.1 EA	
13 SLQ-5006 rev A	Tape, Packaging, 2 Mil, 2" x 110 yds, Clear	Packaging	In Prod	2	5.2 IN	
SUBSTITUTES						
	RANK  ITEM NUMBER	ITEM NAME 	NOTES			
	1 PWY-5000 rev A	Tape, Packaging, 2 5.2 IN Mil, 2 in x 110 yds, Clear				
14 MP-SLQ-0002 rev A	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Manufacturing Procedure	In Prod	2	1 EA	
15 FM-SLQ-0001 rev A	Form, Run Sheet for SILQ Automation Machine	Form	In Prod	2	1 EA	
16 TM-HDX-0001 rev A	Test Method, UV Light Absorbance, Hydrophilix Catheters	Test Method	In Prod	2	1 each	
17 TM-HDX-0002 rev A	Test Method, Rinse Test, Hydrophilix Catheter	Test Method	In Prod	2	1 each	
18 MP-0016 rev A	Manufacturing Procedure, Label Generation	Manufacturing Procedure	In Prod	2	1 each	Max 12 month product shelf life
19 FM-0033 rev C	Form, Label Generation and Reconciliation	Form	In Prod	2	4 each	

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE		
20	IP-0002 rev B	Inspection Procedure, Product Labeling and Packaging	Inspection Procedure	In Prod	1	1 each
21	MP-SLQ-0001 rev A	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	Manufacturing Procedure	In Prod	2	1 EA



Date Inspected	Part No.			Rev.	Description		
	SLQ-3000			A	Solution, Polymer, SILQ		
Lot No.		Lot QTY		Approved Supplier(s)			
				SILQ Technologies Corp. / Hydrophilix			
Specification	No Damage to Outer Packaging.	Verify container is sealed and is not leaking.	Verify container is labeled with Part No, Lot No, Quantity	C of A is present with supplier's inspection data that complies to the specification.			
Revision			N/A	N/A			
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A			
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A			
UOM	N/A	N/A	N/A	N/A			
Measuring Equipment Type	N/A	N/A	N/A	N/A			
Asset No.	N/A	N/A	N/A	N/A			
Cal Due Date	N/A	N/A	N/A	N/A			
AQL	N/A	N/A	N/A	N/A			
Sample Size	Shipment	Shipment	Shipment	Shipment			
Sample No.	Variable/ Attribute Data						
1							
2							
3							
4							
5							
Result Summary							
Pass QTY							
Fail QTY							
Result (P/F)							
Statistical Summary							
Min	N/A	N/A	0.000	0.000	0.000	0.000	0.000
Max	N/A	N/A	0.000	0.000	0.000	0.000	0.000
Mean	N/A	N/A	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!
StDev	N/A	N/A	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!
Comments							
Disposition:						NCMR No:	
Inspected By:						Date:	
Verified By:						Date:	

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-3000
	Description: Solution, SILQ Polymer		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 1

Material Requirements				☐ N/A
SILQ Polymer Solution (Proprietary)				
Traceability	☒ LOT ☐ SERIAL ☐ NONE			
Dimensional Requirements				☒ N/A
Visual Requirements				☒ N/A
Label Content Specification				☐ N/A
Product to be labeled with part number, lot number and quantity				
Certificate of Compliance (C of C)	☐ YES ☐ NO ☐ N/A	Certificate of Analysis (C of A)	☒ YES ☐ NO ☐ N/A	
Material Compliance				
REACH Compliant	☐ YES ☐ NO ☒ N/A	ROHS Compliant	☐ YES ☐ NO ☒ N/A	
Conflict Material	☐ YES ☐ NO ☒ N/A	FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A	
Material Safety Data Sheet (MSDS)	☒ YES ☐ NO ☐ N/A	BPA Free	☐ YES ☐ NO ☒ N/A	
Latex Free	☐ YES ☐ NO ☒ N/A	DEHP Free	☐ YES ☐ NO ☒ N/A	
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A	Documented evidence for compliance to incorporated material specification file		
Storage Requirements				☐ N/A
Protect from heat and direct sunlight. Store in a cool location.				
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☒ N/A
Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.				

Change History

Revision	Change Description
A	Initial Release



2-WAY 100% SILICONE CLEARTRACT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	1



Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

Manufactured for:
SILQ Technologies Corp.
12100 Wilshire Blvd, Suite 800 | Los Angeles, CA 90025 | 949-735-2614

LOT	SLQ1234-MMDDYY
REF	XXXXXX
	20YY-MM-DD



(01) 11223344556675
(10) XXXXXX
(17) YYMMDD



2-WAY 100% SILICONE CLEARTRACT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	10



Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

Manufactured for:
SILQ Technologies Corp.
12100 Wilshire Blvd, Suite 800 | Los Angeles, CA 90025 | 949-735-2614

LOT	SLQ1234-MMDDYY
REF	XXXXXX
	20YY-MM-DD



(01) 11223344556675
(10) XXXXXX
(17) YYMMDD



2-WAY 100% SILICONE CLEARTRACT CATHETER

REF

XXXXXX

LOT

SLQ1234-MMDDYY



20YY-MM-DD

QTY

100



SILQ Technologies
Corp.

12100 Wilshire Blvd, Suite 800 | Los Angeles, CA 90025 | 424-309-8523

SLQ-4001 Rev. A



R_x Only

2-Way 100% Silicone ClearTract Catheter

Balloon Inflation Volume:

XX Fr X-X ml

Please read all instructions before using this device

Not made with natural rubber latex

Foley Catheter is MR Safe.

Contraindications

This device is contraindicated in the presence of conditions which create unacceptable risk during catheterization.

Intended Use

Urological catheter is intended for drainage of the urinary tract. Catheterization is accomplished by inserting the catheter through the urethra and into the bladder. Intended population is adults and pediatrics.

Directions

Installation & Inflation

1. Prepare application site in a sterile fashion.
2. Visually inspect the catheter for any signs of damage or deterioration prior to use.
3. If a trial inflation is conducted, follow steps 7 & 8 on the sterile field then evacuate the sterile water with the syringe. Ensure the balloon is empty prior to insertion.
4. Apply medical grade lubricant to distal end of the catheter.
5. Gently advance the catheter through the urethra until distal tip is beyond the bladder neck.
6. Confirm the position of the distal tip. The return of urine typically indicates the distal tip is in the bladder.
7. Fill Luer-tipped syringe with sterile water and connect to the inflation valve.
8. Fill the catheter balloon to the labeled volume. Do NOT exceed labeled volume.
9. Slowly withdraw the catheter until resistance is felt. This is a sign the balloon is resting at the bladder neck.

Deflation & Removal

1. Gently connect Luer-tipped syringe to the inflation valve. Gently aspirate while allowing the pressure in the balloon to assist filling the syringe with water occupying the balloon. Remove and reconnect the syringe in the inflation valve if there is no or slow deflation. Rapid deflation could collapse the lumen and prevent adequate deflation of the balloon.
2. If during the removal process the balloon ruptures, ensure all fragments are carefully removed from patient.
3. Discard after use.



Precautions

1. Sterile technique shall be applied during installation and removal of the catheter.
2. It is recommended to be replaced after 7 days of use, however, shall not exceed 30 days.

Cautions

1. Federal law restricts the device to sale or by the order of a physician.
2. Never clamp the catheter. Use catheter plug or valve if required.

Warnings

1. Single Use ONLY. Do NOT reuse. Reuse may cause infection, illness, injury or death.
2. Read the instructions carefully before using the device.
3. Do not use this device beyond the labeled expiration date.
4. Do not use the device if the packaging is punctured or damaged.
5. Do not force catheter as injury to the urethra can occur.
6. Do not use needle.
7. Do not exceed recommended capacities.

Potential Adverse Events

With all Foley Catheters, the following potential adverse events could occur:

1. Urinary tract infection
2. Perforation or laceration of the urethra
3. Perforation or laceration of the bladder
4. Abrasions

Symbol Description

	Catalog Number
	Lot Number
	Use by
	Consult instructions for use
	Single Use
	Sterilized by Ethylene Oxide
	Do not resterilize
	Do not use if package is damaged
	Prescription Use Only (US)
	Foley Catheter is MR Safe

Manufactured for
SILQ Technologies Corp.
12100 Wilshire Blvd, # 800
Los Angeles, CA 90025
(949) 735-2614



Date Inspected	Part No.			Rev.	Description		
	SLQ-5000			A	Bag, 2" W x 18" L x 2 Mil Thk, Polyethylene		
Lot No.		Lot QTY		Approved Supplier(s)			
				Ningbo Yingmed			
Specification	No Damage to Outer Packaging.	Verify parts are free of foreign matter or contamination.					
Revision							
Lower Spec Limit (LSL)	N/A	N/A					
Upper Spec Limit (USL)	N/A	N/A					
UOM	N/A	N/A					
Measuring Equipment Type	N/A	N/A					
Asset No.	N/A	N/A					
Cal Due Date	N/A	N/A					
AQL	N/A	4.0					
Sample Size	Shipment						
Sample No.	Variable/ Attribute Data						
1							
2							
3							
4							
5							
6							
7							
8							
9							
10							
11							
12							
13							
14							
15							
16							
17							
18							
19							
20							
21							
22							
23							
24							
25							
26							
27							
28							
29							
Result Summary							
Pass QTY							
Fail QTY							
Result (P/F)							
Statistical Summary							
Min	N/A	N/A					
Max	N/A	N/A					
Mean	N/A	N/A					
StDev	N/A	N/A					
Comments							
Disposition:						NCMR No:	
Inspected By:						Date:	
Verified By:						Date:	

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-5000
	Description: Bag, Polyethylene		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 1

Material Requirements				☐ N/A
SLQ-5000: Supplied by Ningbo Yingmed with catheter shipments				
Traceability	☐ LOT ☐ SERIAL ☒ NONE			
Dimensional Requirements				☐ N/A
None.				
Visual Requirements				☐ N/A
No foreign matter or contamination				
Label Content Specification				☐ N/A
None.				
Certificate of Compliance (C of C)	☐ YES	☐ NO	☒ N/A	Certificate of Analysis (C of A) ☐ YES ☐ NO ☒ N/A
Material Compliance				
REACH Compliant	☐ YES	☐ NO	☒ N/A	ROHS Compliant ☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES	☐ NO	☒ N/A	FDA (Medical Grade) Compliant ☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES	☐ NO	☒ N/A	BPA Free ☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES	☐ NO	☒ N/A	DEHP Free ☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES	☐ NO	☒ N/A	Documented evidence for compliance to incorporated material specification file
Storage Requirements				☒ N/A
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☐ N/A
Product to be single bagged inside shipping container at a minimum. Double bagging is preferred.				
Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.				

Change History

Revision	Change Description
A	Initial Release



Date Inspected	Part No.			Rev.	Description		
	SLQ-5001			A	Pouch, Chevron, 3.5" W x 21.0" L		
Lot No.		Lot QTY		Approved Supplier(s)			
AMCOR							
	No Damage to Outer Packaging.	Verify parts are free of foreign matter or contamination.	Verify there are no voids/ channels in the 3 seals	Verify there are no pinholes, tears, or voids in the tyvek and poly materials	Verify there are no creases in the unsealed (open) side	Verify product is labeled with Part No, Lot No, Quantity	C of C is present with supplier's inspection data that complies to the specification.
Specification							
Revision							
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A	N/A
UOM	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Measuring Equipment Type	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Asset No.	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A	N/A
AQL	N/A	2.5	2.5	2.5	2.5	N/A	N/A
Sample Size	Shipment					Shipment	Shipment
Sample No.	Variable/ Attribute Data						
1							
2							
3							
4							
5							
6							
7							
8							
9							
10							
11							
12							
13							
14							
15							
16							
17							
18							
19							
20							
21							
22							
23							
24							
25							
26							
27							
28							
29							
Result Summary							
Pass QTY							
Fail QTY							
Result (P/F)							
Statistical Summary							
Min	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Comments							
Disposition:						NCMR No:	
Inspected By:						Date:	
Verified By:						Date:	

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-5001
	Description: Pouch, Tyvek, 3.5" x 21", Chevron		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 2

Material Requirements				☐ N/A
AMCOR 125666A, HS Pouch, UNP TYV 1059B//PET/PE RLA-203 (TLP-058)				
Traceability	☒ LOT ☐ SERIAL ☐ NONE			
Dimensional Requirements				☐ N/A
3.5 x 21"				
Visual Requirements				☐ N/A
- No foreign matter or contamination - No voids or channels in seals - No pinholes, tears or voids in Tyvek / poly materials - No creases on the unsealed (open) side				
Label Content Specification				☐ N/A
Product to be labeled with part number, lot number and quantity				
Certificate of Compliance (C of C)	☒ YES ☐ NO ☐ N/A		Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A
Material Compliance				
REACH Compliant	☐ YES ☐ NO ☒ N/A		ROHS Compliant	☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES ☐ NO ☒ N/A		FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A		BPA Free	☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES ☐ NO ☒ N/A		DEHP Free	☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A		Documented evidence for compliance to incorporated material specification file	
Storage Requirements				☒ N/A
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☐ N/A
Product to be double bagged inside shipping container.				
Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.				

Change History

Revision	Change Description
A	Initial Release

	RAW MATERIAL SPECIFICATION	
	Description: Pouch, Tyvek, 3.5" x 21", Chevron	Part No.: SLQ-5001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 2 of 2

Poly Structure

Advanced Breathable Pouches for Medical Applications

Breathability and performance suitable for a wide variety of medical applications. Selection of laminates to meet your needs for tensile strength, tear strength, or puncture resistance. Available with coated and uncoated Tyvek® and different gauges of material to meet sterility applications. Extensive line of laminates allow for the proper choice of materials. Incorporates an adhesive lamination of polyester and polyethylene. Suitable for radiation and EO sterilization for wide range of medical devices. Uncoated Tyvek provides breathability and microbial barrier. Film portion of pouch has excellent clarity and strength properties. Consistent peelable seal is feature of pouch.

[Enquire](#)



Applications

Healthcare

Medical

- Minimally Invasive Devices
- Catheters, Guidewires & Stents

Technical Specifications

Delivery

- Product delivered in boxes. Delivery process depending on order quantity.

Overview

Features

- Suitable for EtO and Irradiation Sterilisation
- Available as Uncoated Tyvek
- Clean fiber tear free opening

Customer benefits

- Breathability and performance suitable for a wide variety of medical applications
- Selection of laminates to meet your needs for tensile strength, tear strength, or puncture resistance.

End consumer benefits

- Easy Opening - Easy and clean opening
- Available in various sizes
- product clarity

End consumer types

- Hospitals
- Medical-Device-Manufacturers

Materials

- Tyvek

Continents available

- North America
- South America

Marketing benefits

- Consistent peelable seal reduces mishandling errors

Minimum order Units per year

- 1 box

Product safety

- Comprehensive Factbook available with aging data

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-5002
	Description: Box, 10-Pack Inner, 20.75 L x 4.12 W x 2.37 H		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 1

Material Requirements				☐ N/A
SquarePeg: Box, White, 20.75 L x 4.12 W x 2.37 H				
Traceability	☐ LOT ☐ SERIAL ☒ NONE			
Dimensional Requirements				☐ N/A
20.75 L x 4.12 W x 2.37 H (REF)				
Visual Requirements				☒ N/A
Label Content Specification				☒ N/A
Certificate of Compliance (C of C)	☐ YES ☐ NO ☒ N/A	Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A	
Material Compliance				
REACH Compliant	☐ YES ☐ NO ☒ N/A	ROHS Compliant	☐ YES ☐ NO ☒ N/A	
Conflict Material	☐ YES ☐ NO ☒ N/A	FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A	
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A	BPA Free	☐ YES ☐ NO ☒ N/A	
Latex Free	☐ YES ☐ NO ☒ N/A	DEHP Free	☐ YES ☐ NO ☒ N/A	
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A	Documented evidence for compliance to incorporated material specification file		
Storage Requirements				☒ N/A
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☒ N/A
Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.				

Change History

Revision	Change Description
A	Initial Release

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-5004
	Description: Label, 2.3125 in x 7.5 in		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 2

Material Requirements				☐ N/A
AMAZON: House Label P/N HL-99019-006, White Thermal Direct Label				
Traceability	☐ LOT ☐ SERIAL ☒ NONE			
Dimensional Requirements				☐ N/A
2-5/16" x 7-1/2" (REF)				
Visual Requirements				☒ N/A
Label Content Specification				☒ N/A
Certificate of Compliance (C of C)	☐ YES	☐ NO	☒ N/A	Certificate of Analysis (C of A) ☐ YES ☐ NO ☒ N/A
Material Compliance				
REACH Compliant	☐ YES	☐ NO	☒ N/A	ROHS Compliant ☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES	☐ NO	☒ N/A	FDA (Medical Grade) Compliant ☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES	☐ NO	☒ N/A	BPA Free ☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES	☐ NO	☒ N/A	DEHP Free ☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES	☐ NO	☒ N/A	<i>Documented evidence for compliance to incorporated material specification file</i>
Storage Requirements				☒ N/A
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☒ N/A
<i>Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.</i>				

Change History

Revision	Change Description
A	Initial Release

	RAW MATERIAL SPECIFICATION	
	Description: Label, 2.3125 in x 7.5 in	Part No.: SLQ-5004
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	
		Rev.: A
		Page 2 of 2



Compatible with DYMO 99019 1-Part Internet Postage Labels (2-5/16" x 7-1/2") - BPA Free! (6 Rolls; 150 Labels per Roll)


Brand: HouseLabels

★★★★☆ 32 ratings | 3 answered questions

Amazon's Choice for "99019 labels"

Price: \$25.99 & FREE Shipping. Details & FREE Returns

- These labels are the easiest internet postage labels to apply to any shipment as only ONE label, not multiple labels, are required for each package. Print from your favorite shipping programs, such as Endicia.com, Stamps.com, PayPal, or eBay. See details "For Endicia/DAZzle Users" below in the Products Description section.
- Compatible with DYMO Printers: 320, 400, 400 Duo, 400 Twin Turbo, 450, 450 Duo, 450 Turbo, 450 Twin Turbo, 4XL, SE300, SE450 CoStar LabelWriters: XL, XL+, BC, BC+, TURBO, SE200, SE250, ACII 200 / ASCII 250
- All labels sold by HouseLabels are tested to ensure that they fall within printer specifications Labels are smooth and bright white for professional results Labels are smudge-proof as well as water and oil resistant All temperature, strong permanent adhesive is excellent for all general purpose labeling applications Each roll is individually packed into a protective black polybag and then placed in a mini-carton to protect it from dirt, sunlight and scuffing.

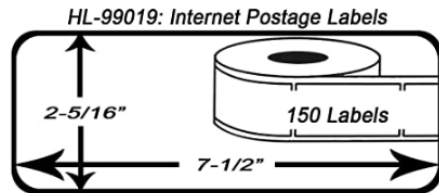

Roll over image to zoom in

GOLDISTOCK

Get your high quality



Heavier Tyvek Wristbands 7.5 Mil
- Goldistock™ Top 20* 200 Ct.
Variety Pack- ¼" Arm Bands...

Product information

Technical Details

Manufacturer	HouseLabels
Brand Name	HouseLabels
Item Weight	3.9 pounds
Product Dimensions	6 x 3.1 x 8 inches
Item model number	HL-99019-006
Is Discontinued By Manufacturer	No
Color	White
Size	2-5/16" x 7-1/2"
Manufacturer Part Number	DY-99019-6

Additional Information

ASIN	B017C7MDGM
Customer Reviews	★★★★☆ 32 ratings 4.5 out of 5 stars
Best Sellers Rank	#21,159 in Office Products (See Top 100 in Office Products) #3 in Internet Postage Labels
Date First Available	August 4, 2012

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-5005
	Description: Ribbon, ZEBRA 3200 Wax-Resin		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 2

Material Requirements				☐ N/A
- ZEBRA P/N: 03200GS11007 - Amazon P/N: 03200GS11007				
Traceability	☐ LOT ☐ SERIAL ☒ NONE			
Dimensional Requirements				☐ N/A
4.3" x 243' (REF)				
Visual Requirements				☒ N/A
Label Content Specification				☒ N/A
Certificate of Compliance (C of C)	☐ YES ☐ NO ☒ N/A		Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A
Material Compliance				
REACH Compliant	☐ YES ☐ NO ☒ N/A		ROHS Compliant	☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES ☐ NO ☒ N/A		FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A		BPA Free	☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES ☐ NO ☒ N/A		DEHP Free	☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A		Documented evidence for compliance to incorporated material specification file	
Storage Requirements				☒ N/A
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☒ N/A
Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.				

Change History

Revision	Change Description
A	Initial Release

 pathway	RAW MATERIAL SPECIFICATION								Part No.: SLQ-5005			
	Description: Ribbon, ZEBRA 3200 Wax-Resin								Rev.: A			
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>								Page 2 of 2			

For use with all 4" desktop printers

3200 Wax/Resin Ribbons

A high-performance wax/resin ribbon for use with Zebra thermal transfer printers. Provides outstanding barcode print quality. Chemical resistant when printed on Zebra paper and synthetic label stocks. Smear- and scratch-resistant, this ribbon offers excellent staying power and durability in harsh environments. UL Recognized when used with Zebra approved labels.

Core: 0.5" ID

Width		Length		Part #	Quantity Per		Outer Box	Stocking Locations				
inches	mm	feet	m		Carton	Outer Box	Weight	Price/roll	WI	GA	TX	CA
2.52"	64mm	243'	74m	03200GS06407	12	12	1 lbs.	\$4.12	X	X	X	X
3.31"	84mm	243'	74m	03200GS08407	12	12	2 lbs.	\$5.39	X	X	X	X
4.30"	110mm	243'	74m	03200GS11007	12	12	2 lbs.	\$7.06	X	X	X	X

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-5006
	Description: Tape, Packaging		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 2

Material Requirements				☐ N/A
ULINE P/N S-423, Tape, Clear				
Traceability	☐ LOT ☐ SERIAL ☒ NONE			
Dimensional Requirements				☐ N/A
2 mil, 2 in x 110 yds				
Visual Requirements				☒ N/A
Label Content Specification				☒ N/A
Certificate of Compliance (C of C)	☐ YES ☐ NO ☒ N/A	Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A	
Material Compliance				
REACH Compliant	☐ YES ☐ NO ☒ N/A	ROHS Compliant	☐ YES ☐ NO ☒ N/A	
Conflict Material	☐ YES ☐ NO ☒ N/A	FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A	
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A	BPA Free	☐ YES ☐ NO ☒ N/A	
Latex Free	☐ YES ☐ NO ☒ N/A	DEHP Free	☐ YES ☐ NO ☒ N/A	
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A	Documented evidence for compliance to incorporated material specification file		
Storage Requirements				☒ N/A
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☒ N/A
Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.				

Change History

Revision	Change Description
A	Initial Release

	RAW MATERIAL SPECIFICATION				Part No.: SLQ-5006
	Description: Tape, Packaging				Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>				Page 2 of 2

Uline Industrial Tape - 2 Mil, 2" x 110 yds, Clear

2.0 MIL



Our best selling tape. Strong solvent acrylic adhesive sticks to virtually any surface.

DURABLE: Industrial Tape stands up to heat, humidity and cold. **Application range 0°F to 140°F.**

RELIABLE: Long shelf life.

APPEARANCE: Crystal clear - Won't yellow. Makes packages look their best.

Free Offer

 [Enlarge](#)

20 LBS. PER INCH TENSILE STRENGTH

MODEL NO.	COLOR	MIL	SIZE	ROLLS/CASE	LBS./CASE	PRICE PER ROLL					ADD TO CART
						18	36	72	144	288+	
S-423	Clear	2.0	2" x 110 yds.	36	22	\$3.25	\$1.69 per roll				18 ADD

 [Additional Info](#)

 [Email Page](#)

 [Shopping Lists](#)

[Request a Catalog](#)

DIMENSIONS:

- Core: 3"
- Diameter: 4 1/2"

FEATURES:

- Rolls have a paper starter tab.

MATERIAL:

- Backing: Polypropylene film
- Latex free

SHELF LIFE:

- 1 year

TEMPERATURE:

- Performance: -60 to 140 F

USAGE:

- Takes approximately 40 minutes to reach maximum bond.

Availability: [In Stock](#)

Unit Weight: 0.56 lbs.

[SDS Sheet](#)

[Catalog Page 769](#)

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-5011
	Description: Label, 4 in x 6 in, White		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 2

Material Requirements				☐ N/A
ONLINE LABELS Product #RL786TT or equivalent, White Thermal Transfer Label				
Traceability	☐ LOT ☐ SERIAL ☒ NONE			
Dimensional Requirements				☐ N/A
4 in x 6 in (REF) 3 in Core, 8 in OD (REF)				
Visual Requirements				☒ N/A
Label Content Specification				☒ N/A
Certificate of Compliance (C of C)	☐ YES ☐ NO ☒ N/A		Certificate of Analysis (C of A)	☐ YES ☐ NO ☒ N/A
Material Compliance				
REACH Compliant	☐ YES ☐ NO ☒ N/A		ROHS Compliant	☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES ☐ NO ☒ N/A		FDA (Medical Grade) Compliant	☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES ☐ NO ☒ N/A		BPA Free	☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES ☐ NO ☒ N/A		DEHP Free	☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES ☐ NO ☒ N/A		Documented evidence for compliance to incorporated material specification file	
Storage Requirements				☒ N/A
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☒ N/A
<i>Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.</i>				

Change History

Revision	Change Description
A	Initial Release

	RAW MATERIAL SPECIFICATION	
	Description: Label, 4 in x 6 in, White	Part No.: SLQ-5011
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 2 of 2




 My Account
  Cart


Sheet Labels
 Roll Labels
 Printed Labels
 Design & Templates
 Help

 Live Chat
 |
  1-888-575-2235

Home > Labels on Rolls > 3" Core Roll Labels > RL786TT

Product RL786TT - 4" x 6" Thermal Roll Labels - 3" Core / 8" Outer Diameter

★★★★★ (4.9 based on 9 reviews for Standard White Thermal Transfer) | 5 Questions & Answers

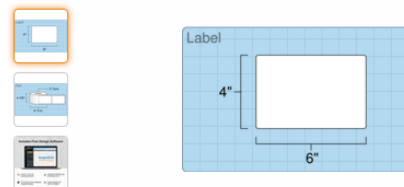

Product Specs for RL786 (detailed specs)

Label Size: 4" x 6"
 Core Diameter: 3"
 Outer Diameter: 8"
 Labels per Roll: 1,000

Measurement display: Decimals (Inches)

[View Label Templates for RL786](#)

[Printing Instructions for RL786TT](#)

1 Select a material:

[Help me find a material](#)

- ☒ Standard White Thermal Transfer (ribbon required)
- ☐ Standard White Direct Thermal
- ☐ Removable Thermal Transfer (ribbon required)
- ☐ Removable Direct Thermal

2 Select a quantity:

- ☐ Enter a quantity not listed below.
- ☐ 1 roll / 1,000 Labels (\$13.95 / roll) **\$13.95**
- ☐ 4 rolls / 4,000 Labels (\$12.49 / roll) **\$49.95**
- ☐ 12 rolls / 12,000 Labels (\$11.91 / roll) **\$142.95**

	RAW MATERIAL SPECIFICATION		Part No.: SLQ-5012
	Description: Box, 20 in. L x 18 in W x 6 in H, 200 LB Test		Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Page 1 of 2

Material Requirements				☐ N/A
ULINE P/N S-4903, Corrugate Box				
Traceability	☐ LOT ☐ SERIAL ☒ NONE			
Dimensional Requirements				☐ N/A
20 x 18 x 6 in.				
Visual Requirements				☒ N/A
Label Content Specification				☒ N/A
Certificate of Compliance (C of C)	☐ YES	☐ NO	☒ N/A	Certificate of Analysis (C of A) ☐ YES ☐ NO ☒ N/A
Material Compliance				
REACH Compliant	☐ YES	☐ NO	☒ N/A	ROHS Compliant ☐ YES ☐ NO ☒ N/A
Conflict Material	☐ YES	☐ NO	☒ N/A	FDA (Medical Grade) Compliant ☐ YES ☐ NO ☒ N/A
Material Safety Data Sheet (MSDS)	☐ YES	☐ NO	☒ N/A	BPA Free ☐ YES ☐ NO ☒ N/A
Latex Free	☐ YES	☐ NO	☒ N/A	DEHP Free ☐ YES ☐ NO ☒ N/A
Animal Derive Material (ADM) Free	☐ YES	☐ NO	☒ N/A	<i>Documented evidence for compliance to incorporated material specification file</i>
Storage Requirements				☒ N/A
Shelf-Life Requirements				☒ N/A
Packaging Requirements				☒ N/A
<i>Supplier will notify Pathway LLC in advance of any changes prior to the implementation of the change in the supplied material.</i>				

Change History

Revision	Change Description
A	Initial Release

 pathway	RAW MATERIAL SPECIFICATION					Part No.: SLQ-5012
	Description: Box, 20 in. L x 18 in W x 6 in H, 200 LB Test					Rev.: A
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>					Page 2 of 2


20 x 18 x 6" Corrugated Boxes



 Enlarge

Perfect for apparel, picture frames, books, scales, spools of wire and electronic equipment.

- Boxes ship flat for freight and storage savings.
- Larger quantities, call for quote.

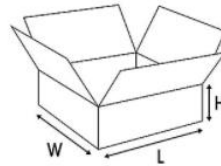

200 LB. TEST

SOLD IN BUNDLE QUANTITIES

MODEL NO.	INSIDE DIMENSIONS L x W x H	PRICE PER BOX					BUNDLE/ BALE QTY.	LBS./ BNDL.	IN STOCK SHIPS TODAY	
		25	100	250	500	1,000+				
S-4903	20 x 18 x 6"	\$2.32	\$2.18	\$2.05	\$1.90	\$1.83	25 / 250	41 lbs.	25	ADD

NOTES: (UNLESS OTHERWISE SPECIFIED)

1.

GENERAL

1.

ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY SILQ PRIOR TO PRODUCTION.

2.

PARTS SHALL ARRIVE NON-STERILE.

3.

UNLESS OTHERWISE SPECIFIED, INSPECT TO AQL 4.0.
2.

MATERIAL

1.

MATERIAL: SEE PAGE 2.

2.

FINISH: NONE
3.

DIMENSIONAL

1.

INDICATES CRITICAL DIMENSION.

2.

() INDICATES REFERENCE DIMENSION.
4.

DIMENSIONAL INSPECTION

1.

DIMENSIONS APPLY AFTER FINISHING AND CLEANING PROCESS.

2.

ALL MEASUREMENTS SHALL BE TAKEN WITH DEVICES HOLDING A CURRENT CALIBRATION CERTIFICATE.
5.

VISUAL INSPECTION

1.

WHEN VIEWED FROM 18 INCHES AWAY, WITHOUT THE USE OF MAGNIFICATION AND UNDER AMBIENT ROOM LIGHTING, THE PART SHALL BE FREE FROM PARTICULATES.

2.

WHEN VIEWED FROM 18 INCHES AWAY, WITHOUT THE USE OF MAGNIFICATION AND UNDER AMBIENT ROOM LIGHTING ALONG THE SHAFT AND AT THE TIP, THE PART SHALL BE FREE OF MOLDING DEFECTS SUCH AS FLASH, EXCESS MATERIAL, SINKS, ETC.

3.

AT THE CONNECTOR, THE SIZE OF THE CATHETER SHALL BE INDICATED AS FOLLOWS: "16 FR, 5.3MM. THE BALLON VOLUME SHALL BE INDICATED AS FOLLOWS: "10ML/CC." THE PRINTING SHALL BE LEGIBLE AND FREE FROM MISSING INFORMATION, SMUDGES, AND DUPLICATES.
6.

PACKAGING

1.

CATHETERS SHALL COME IN A POLYETHYLENE SLEEVE WITH A MAXIMUM OF 1 CATHETER PER SLEEVE.

2.

CATHETERS SHALL BE SEALED IN AN OUTER BAG WITH A MAXIMUM OF 500 UNITS PER BAG.
7.

LABELING


1.

OUTER BAGS SHALL BE LABELLED WITH PART NUMBER, REVISION, AND LOT NUMBER.
8.

CERTIFICATIONS

1.

CERTIFICATE OF CONFORMANCE (COFC) REQUIRED FOR EACH UNIQUE LOT NUMBER.

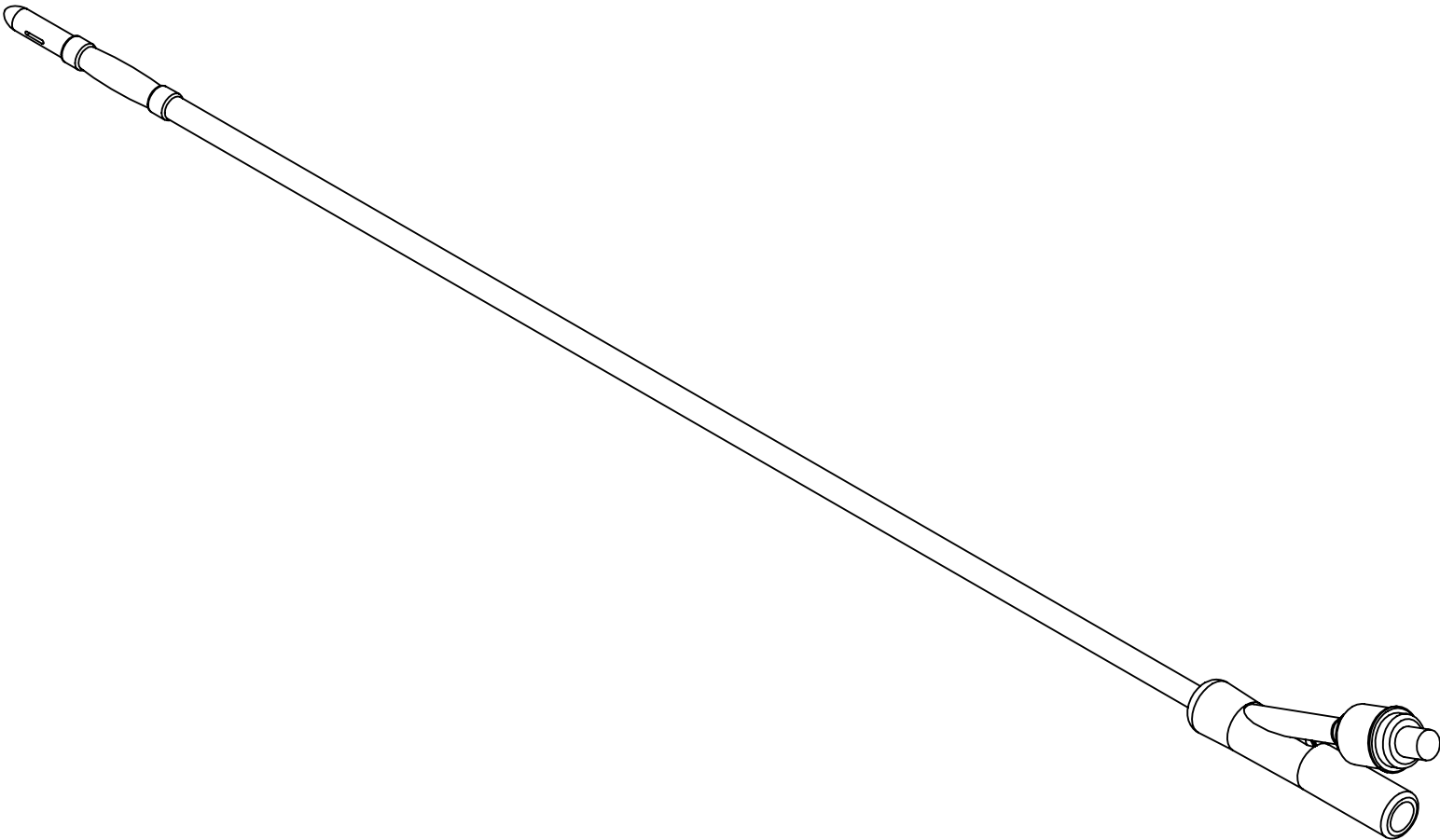

REVISIONS				
ZONE	REV.	DESCRIPTION	DATE	APPROVED
	A	INITIAL RELEASE		

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.








TITLE:
CATHETER, 16 FR X 10 ML,
SILICONE, NY FOLEY, 2-WAY

SIZE	DWG. NO.	REV
B	SLQ-6004	A
SCALE: 2:3		13.82 g
SHEET 1 OF 2		

DO NOT SCALE DRAWING

ITEM NO.	DESCRIPTION	MATERIAL	HARDNESS	COLOR	QTY.	NOTES
1	TUBING	SILICONE	80 SHORE A	NATURAL/XRAY	1	X-RAY COLOR IS BLUE
2	BALLOON CUFF	SILICONE	25 SHORE A	NATURAL	1	ORDINARY BALLOON TYPE
3	TIP	SILICONE	50 SHORE A	BLUE	1	TIP HAS 33% CONCENTRATION OF BARIUM SULFATE
4	2-WAY FUNNEL	SILICONE	57 SHORE A	NATURAL	1	N/A
5	CHECK VALVE	PP/SS/PS/SILICONE	70 SHORE A	WHITE	1	BESPAK MED
6	SLEEVE	ABS	80 SHORE A	ORANGE	1	N/A

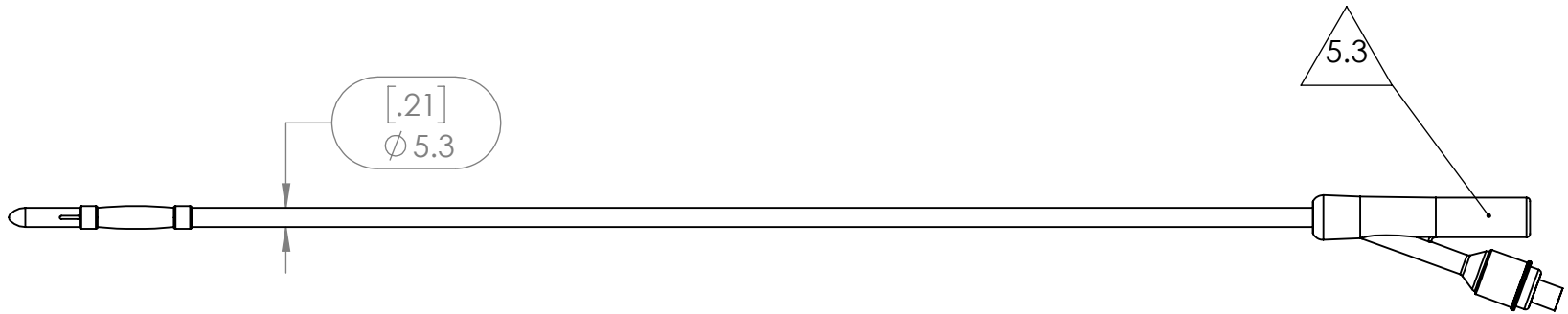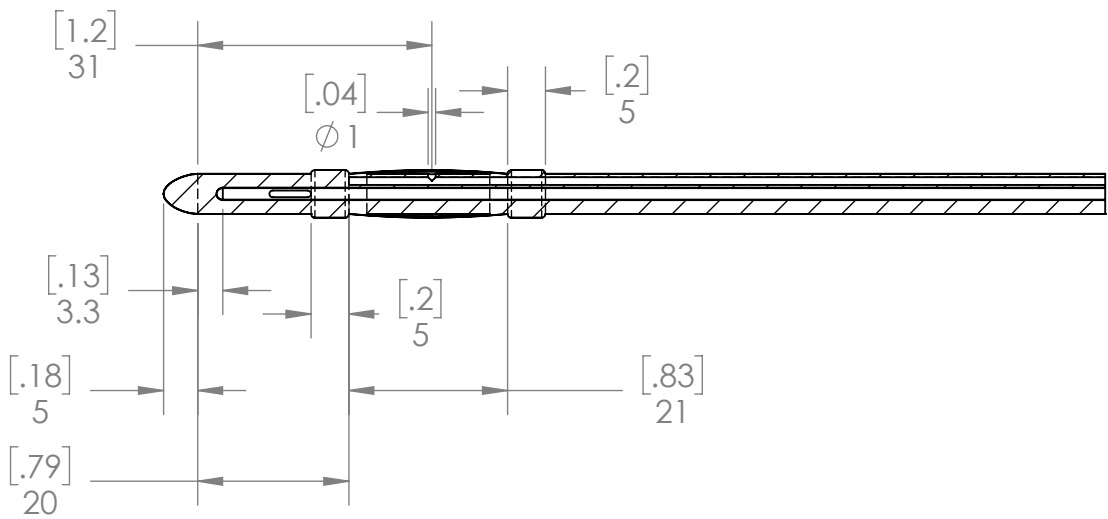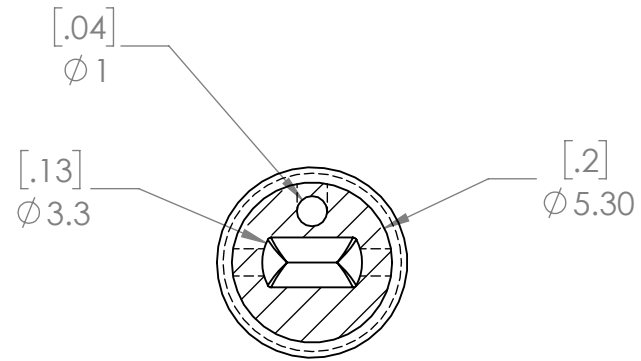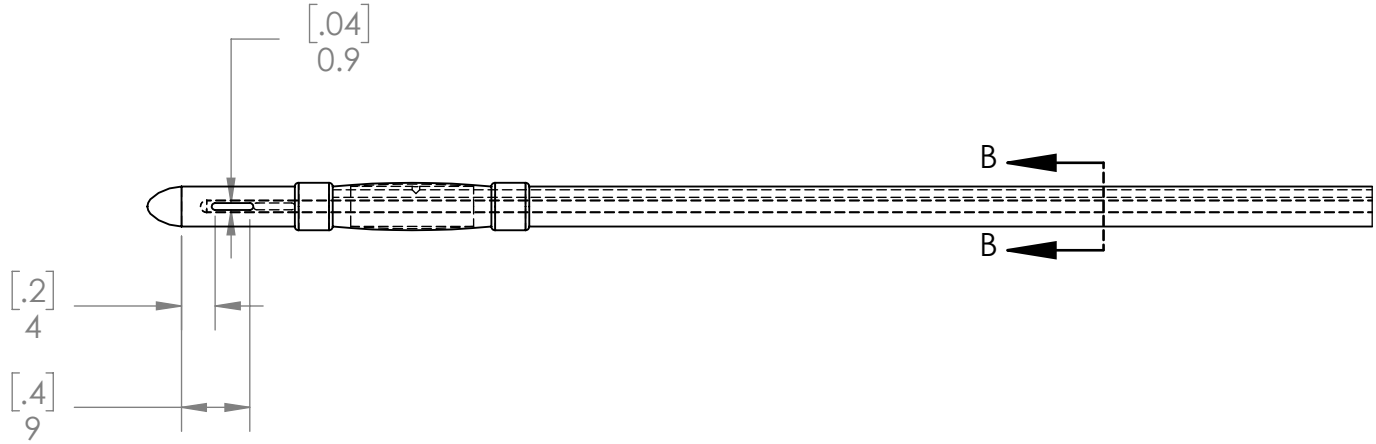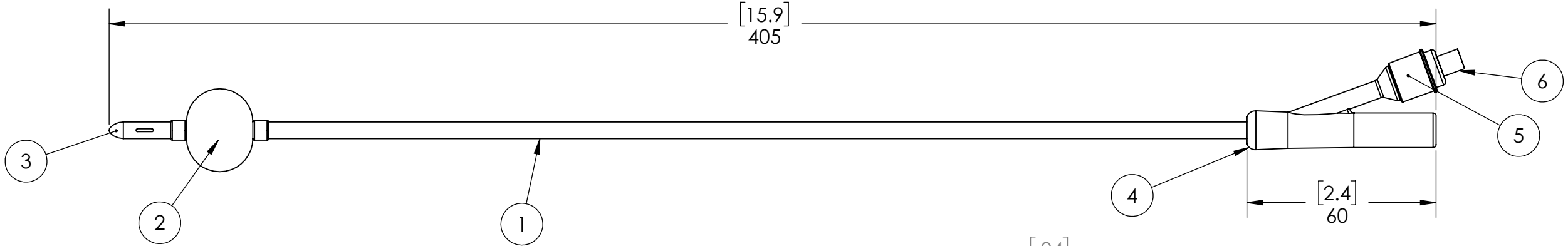

SIZE B	DWG. NO. SLQ-6004	REV A
SCALE: 2:3	SHEET 2 OF 2	




Date Inspected	Part No.			Rev.	Description							
	SLQ-6004			A	2-Way Foley Catheter, 16 Fr x 10 mL							
Lot No.		Lot QTY		Approved Supplier(s)								
Ningbo Yingmed												
Specification	No Damage to Outer Packaging and Catheters are double bagged.	Catheters shall come in a polyethylene sleeve with a maximum of 1 catheter per sleeve and the catheters shall be sealed in an outer bag with a maximum of 500 units per bag.	C of C is present with supplier's inspection data that complies to the specification.	When viewed from 18 inches away, without the use of magnification and under ambient room lighting, the part shall be free from particulates.	When viewed from 18 inches away, without the use of magnification and under ambient room lighting along the shaft and at the tip, the part shall be free from molding defects such as flash, excess material, sinks, etc.	At the connector, the size of the catheter shall be indicated as follows: "16 Fr, 5.3mm" The balloon volume shall be indicated as follows: "10 ml/cc." The printing shall be legible and free from missing information, smudges, and duplicates.	5.3 ± 0.3	TM-0013 Foley Catheter Leak Test				
Revision	N/A	N/A	N/A	N/A	A	A	A	A				
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	N/A	N/A	5.0	N/A				
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A	5.6	N/A				
UOM	N/A	N/A	N/A	N/A	N/A	N/A	mm	N/A				
Measuring Equipment Type	N/A	N/A	N/A	N/A	N/A	N/A	Calipers	N/A				
Asset No.	N/A	N/A	N/A	N/A	N/A	N/A		N/A				
Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A		N/A				
AQL	N/A	N/A	N/A	4.0	4.0	4.0	4.0	4.0				
Sample Size	Shipment	Shipment	Shipment									
Sample No.	Variable/ Attribute Data											
1												
2												
3												
4												
5												
6												
7												
8												
9												
10												
11												
12												
13												
14												
15												
16												
17												
18												
19												
20												
21												
22												
23												
24												
25												



26								
27								
28								
29								
Result Summary								
Pass QTY								
Fail QTY								
Result (P/F)								
Statistical Summary								
Min	N/A	N/A	N/A	N/A	N/A	N/A	0.000	N/A
Max	N/A	N/A	N/A	N/A	N/A	N/A	0.000	N/A
Mean	N/A	N/A	N/A	N/A	N/A	N/A	#DIV/0!	N/A
StDev	N/A	N/A	N/A	N/A	N/A	N/A	#DIV/0!	N/A
Comments								
Disposition:						NCMR No:		
Inspected By:						Date:		
Verified By:						Date:		

NOTES: (UNLESS OTHERWISE SPECIFIED)

1. DRAWING IS FOR REFERENCE ONLY. SEE RELEASED BILL OF MATERIALS FOR CURRENT PARTS/REVISION AND DOCUMENTATION REQUIREMENTS.
2. EXPIRATION DATE TO BE ONE (1) YEAR FROM THE DATE OF DEVICE PACKAGING AND SHALL BE DENOTED AS YYYY-MM-DD.
3. ARTWORK PRINTED ON LABEL WITH RIBBON.
4. VERIFY LABEL INFORMATION IS COMPLIANT WITH THE WORK ORDER.
5. POUCH LABEL APPLICATION SHALL BE CENTERED IN BOTH DIRECTIONS ON CLEAR SIDE OF POUCH WITH PEEL (CHEVRON) END ON THE RIGHT SIDE AND THAT LABEL TEXT IS ORIENTED IN THE SAME DIRECTION.
6. INNER BOX LABEL APPLICATION (2X) SHALL BE CENTERED IN BOTH DIRECTIONS AND LOCATED ON THE FRONT AND BACK FACES AND THAT LABEL TEXT IS ORIENTED IN THE SAME DIRECTION.
7. OUTER BOX LABEL APPLICATION (2X) SHALL BE LOCATED IN THE UPPER LEFT CORNER ON TWO ADJACENT SIDES OF THE BOX.
8. APPLY POLYMER COATING TO CATHETER ONLY IN THE HATCHED AREA.
9. PLACE CATHETER INTO SLEEVE AND INTO POUCH.
10. INSIDE THE STERILE BARRIER, DEVICE IS TO BE FREE OF PARTICULATES WHEN OBSERVED AT APPROXIMATELY 18 INCHES AWAY, UNDER AMBIENT ROOM LIGHTING, AND WITHOUT THE USE OF MAGNIFICATION.
11. STERILE BARRIER PACKAGING (INCLUDE SEALED AREA) SHALL NOT EXHIBIT THE FOLLOWING VISUAL DEFECTS:
 - UNSEALED AREAS
 - NONHOMOGENEOUS OR UNDERSEALED AREAS
 - OVERSEALED AREAS
 - NARROW SEALS
 - CHANNELS
 - WRINKLES/FOLDOVERS/CRACKS
 - TEARS/PINHOLES
 - THE PEEL STRENGTH FOR A 1 INCH WIDE SEALED SECTION OF THE STERILE BARRIERPACKAGING SHALL BE >= 1 LB WHEN TESTED PER TM-0003.
 - PACKAGED DEVICES (QTY 10) SHALL BE NEATLY STACKED INTO THE INNER BOX.
12. NEATLY FOLD IFU IN THIRDS AND PLACE INSIDE INNER BOX.
13. SECURE OUTER BOX CLOSED WITH PACKAGING TAPE AS NECESSARY

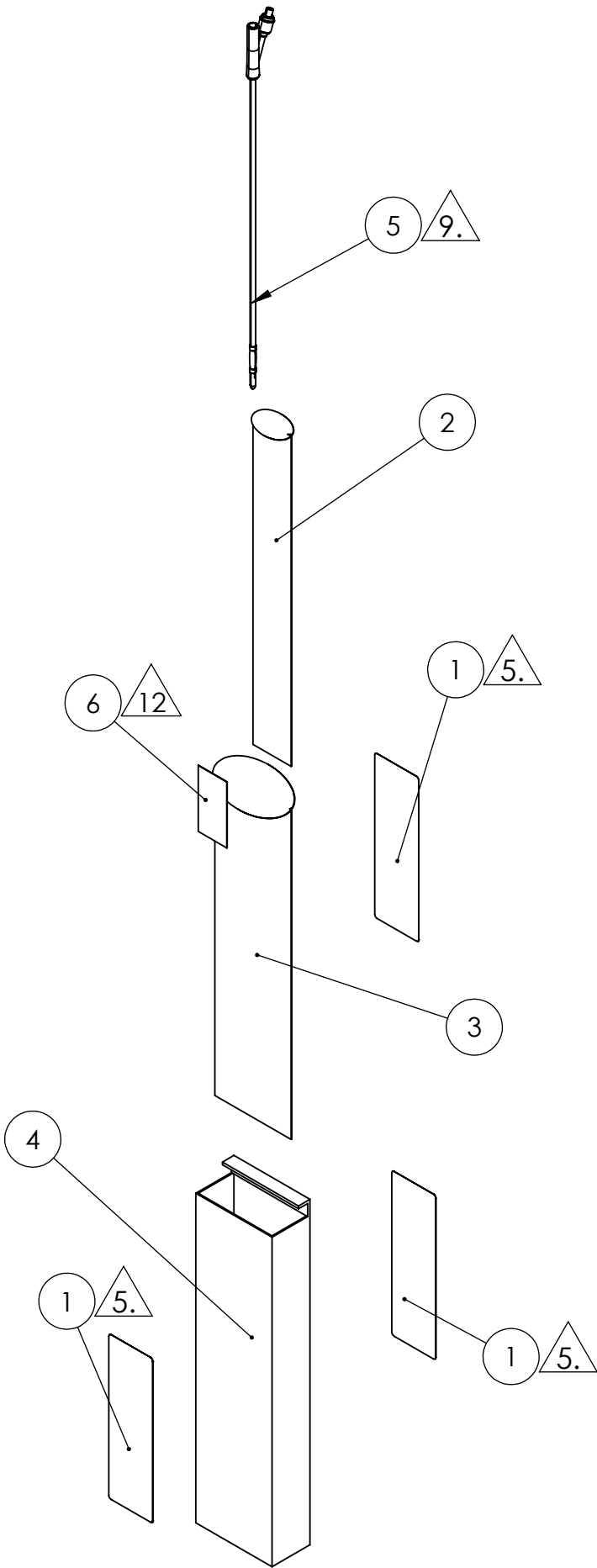

REVISIONS			
ZONE	REV	DESCRIPTION	DATE
N/A	A	INITIAL RELEASE	

ITEM NO.	PARTNO	DESCRIPTION	QTY.
1	SLQ-4000	ARTWORK, LABEL, POUCH/BOX, SILQ CLEARTRACT CATHETER	3
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	1
3	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	1
4	SLQ-5002	BOX, 10-PACK, INNER, 20.75 L X 4.12 W X 2.37 H, PSM	1
5	SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	1
6	SLQ-4005	ARTWORK, IFU, SILQ 2-WAY CLEARTRACT CATHETER	1

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.



TITLE:
SILQ 16 FR X 10 ML CATHETER, NON-STERILE, 10-PACK

SIZE	DWG. NO.	REV
B	SLQ-8004	A




DO NOT SCALE DRAWING

SCALE: 1:1

g


SHEET 1 OF 3

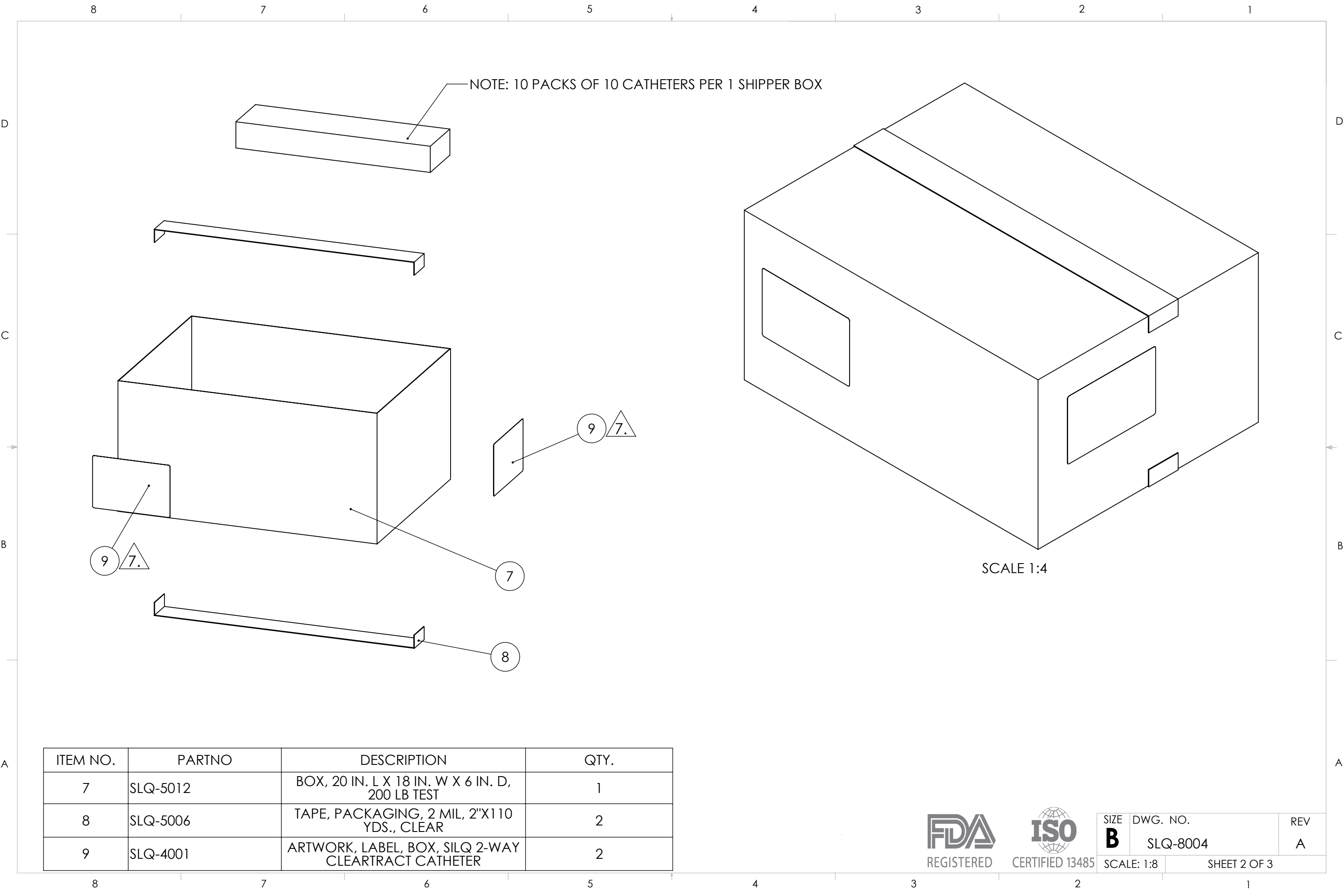

ITEM NO.	PARTNO	DESCRIPTION	QTY.
7	SLQ-5012	BOX, 20 IN. L X 18 IN. W X 6 IN. D, 200 LB TEST	1
8	SLQ-5006	TAPE, PACKAGING, 2 MIL, 2"X110 YDS., CLEAR	2
9	SLQ-4001	ARTWORK, LABEL, BOX, SILQ 2-WAY CLEARTRACT CATHETER	2



SIZE	DWG. NO.	REV
B	SLQ-8004	A
SCALE: 1:8		SHEET 2 OF 3

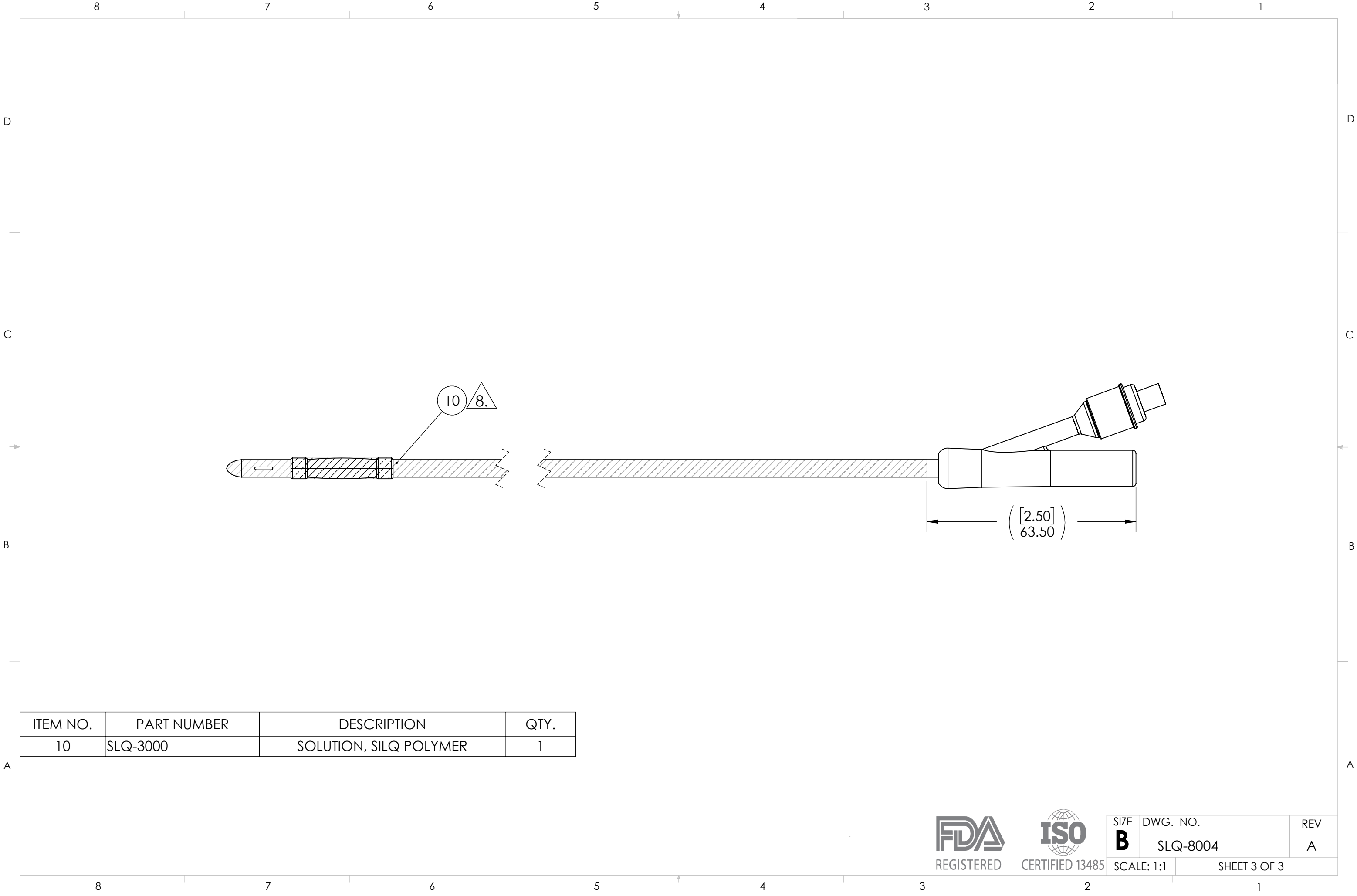



Date Inspected		Part No.		Rev.	Description				
		SLQ-8004		A	SILQ 2-Way Cleartract Foley Catheter, 16 Fr x 10 mL, Non-Sterile, 10-Pack				
Lot No.		Lot QTY		Approved Supplier(s)					
				Pathway					
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	A	A		A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	N/A	0.300	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	N/A	0.020	N/A	N/A	N/A	N/A
	UOM	EA	EA	AU	AU	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Visual	Spectrophotometer	Spectrophotometer	Visual	Visual	Visual	Visual
	Asset No.	N/A	N/A			N/A	N/A	N/A	N/A
	Cal Due Date	N/A	N/A			N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size								
Sample No.		Variable/ Attribute Data							
1									
2									
3									
4									
5									
6									
7									
8									
9									
10									
11									
12									
13									
14									



15									
16									
17									
18									
19									
20									
21									
22									
23									
24									
25									
26									
27									
28									
29									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
Comments									
Disposition:							NCMR No:		
Inspected By:							Date:		
Verified By:							Date:		

	Form, Run Sheet for SILQ Automation Machine		Doc No.: FM-SQL-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 1 of 1

Part Number	Revision	Description	Machine Asset No.

Equipment Setup / Startup			Adjustments during Production			
Setting [Units]		Value	Value	Value	Value	Value
Engineering / Lead Supervisor	Setup / Adjustment Date					
	Setup / Adjustment Time					
	Fill Time [seconds]					
	Confirm Cylinder Actuation All 4 pairs [P/F]					
	UV Cure Time [seconds]					
	Drain Time [seconds]					
	Rinse Time [seconds]					
	Final Drain [seconds]					
	Number of Drain / Rinse Cycles [Count]					
	Confirm Cylinder Actuation All 4 pairs [P/F]					
	Size Setting					
	Performed By Name					
	Performed By Signature					
QA	Verified By Name					
	Verified By Signature					

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 1 of 24

1.0 Purpose:

- 1.1 Document the process for the labeling and packaging of the Silq ClearTract Catheters.

2.0 Scope:

- 2.1 This procedure applies to the labeling and packaging of the Silq ClearTract Catheters.

3.0 References

- 3.1 WI-0012 Gowning & Prohibited Materials, Clean Room
- 3.2 WI-0018 Gowning & Prohibited Materials, for the Assy/Packing of Silicone or any Tacky/Sticky...
- 3.3 WI-0014 Line Clearance

4.0 Responsibilities:

- 4.1 It is the responsibility of the individual(s) conducting these activities to comply with the requirements indicated in this procedure.
- 4.2 Clean Room must be thoroughly wiped down with diluted Clorox or 70% isopropyl alcohol and tacky roller should be used on all working surfaces to minimize the presence of lint/debris.

5.0 General Requirements:

- 5.1 Line clearance shall be conducted per WI-0014.

6.0 Material and Equipment:

- 6.1 ACCU-SEAL Impulse Sealer or equivalent Validated Continuous Band Sealer

7.0 Procedure – Lot/Expiration Verification:

7.1 Verify that the information printed on the on the pouch/10 Pack Box Label as well as the Shipping Box Label matches those recorded on the Work Order as shown in Figure 1 Figure 2.

1. Check to ensure that the Catheter Fr size, Quantity, Lot Number, REF Number and Expiration Date that is printed on the Pouch and Box Label matches the recorded information on the Work Order. The pouch label will have a QTY of 1 and the box label will have a QTY of 10. Verify the lot number date reference is in the form MMDDYY and the expiration date is in the form YYYY-MM-DD (Ex: 2023-01-21). Verify item (01) is a fourteen-digit number representing the GTIN value. Verify item (10) matches the lot number. Verify item (17), which is tied to the expiration date is in the form YYMMDD. The below label is for a 16 Fr, 10 mL catheter expiring on January 21, 2023. After the first few labels, there is no need to keep verifying the information.

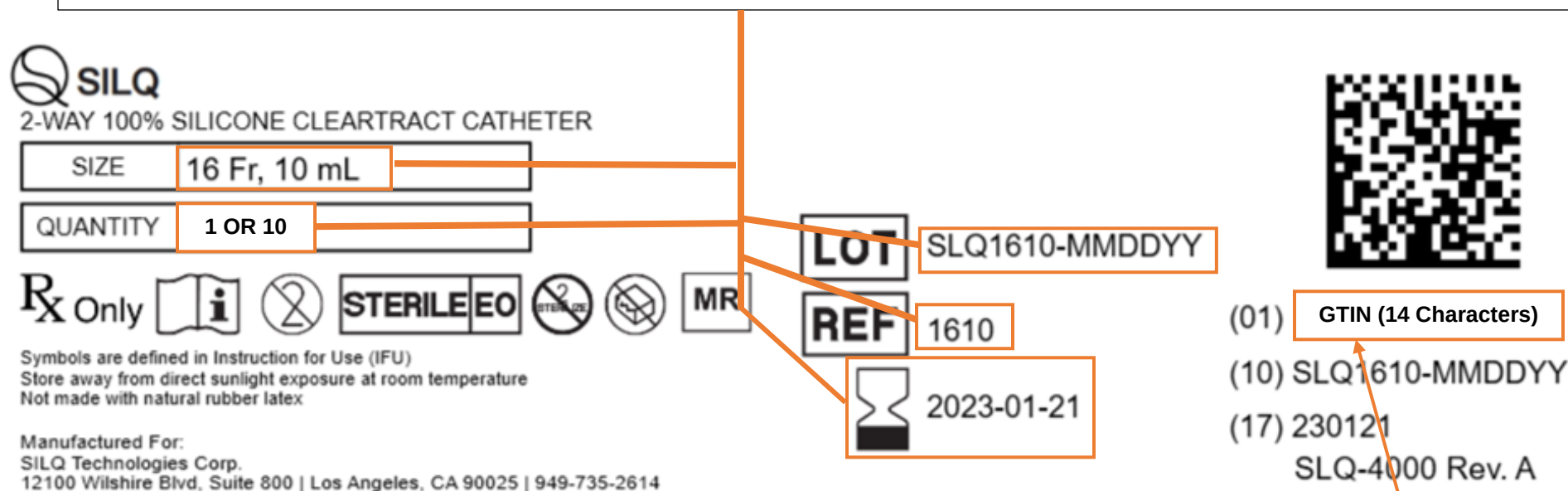


FIG. 1 – Catheter Fr size, Qty, REF and Lot Number/Expiration Date Verification

Single Pouch: starts with a "0"
10-Pack Box: starts with a "1"

1. Check to ensure that the REF, LOT, Expiry Date printed on the Box Label matches the recorded information on the Work Order and match the pouch and 10-pack box labels. Verify the lot number date references and expiration dates are in the format below. Also refer to the image above for more detailed descriptions.



2-WAY 100% SILICONE CLEARTRACT CATHETER

REF

XXXXXX

LOT

SLQ1234-MMDDYY



20YY-MM-DD

QTY

100


 SILQ Technologies Corp.

12100 Wilshire Blvd, Suite 800 | Los Angeles, CA 90025 | 424-309-8523 SLQ-4001 Rev. A

FIG. 2 – Shipping Box Label

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters <i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Doc No.: MP-SLQ-0001 Rev.: A Page 4 of 24
--	--	---

8.0 Procedure – Inner Sleeve/Pouch Packaging and Sealing:

- 8.1 Insert the Polymer coated Silicone Catheter into the Catheter Sleeve and Pouch and seal as shown in Figure 3 through Figure 8.

1. Orient and apply the verified Pouch Label (ref. Figure 1) to the Polyethylene side of the Pouch P/N HDX-5001 such that the text reads away from the chevron end of the Pouch.

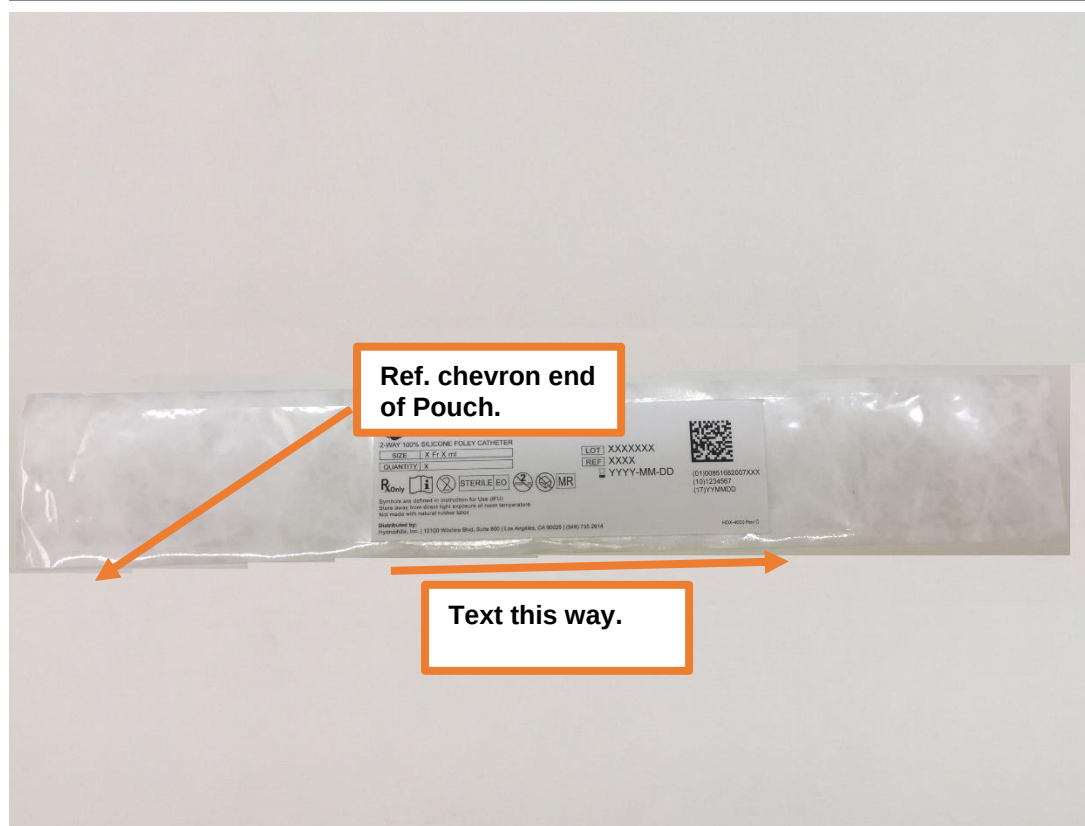


FIG. 3 – Label applied to Pouch



Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters

Doc No.: MP-SLQ-0001

Rev.: A

Page 5 of 24

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

1. Visually inspect the Catheter (see Work Order to determine P/N) to ensure that the outer surface is free of lint/debris (silicone is tacky and has a tendency attract dust/lint).

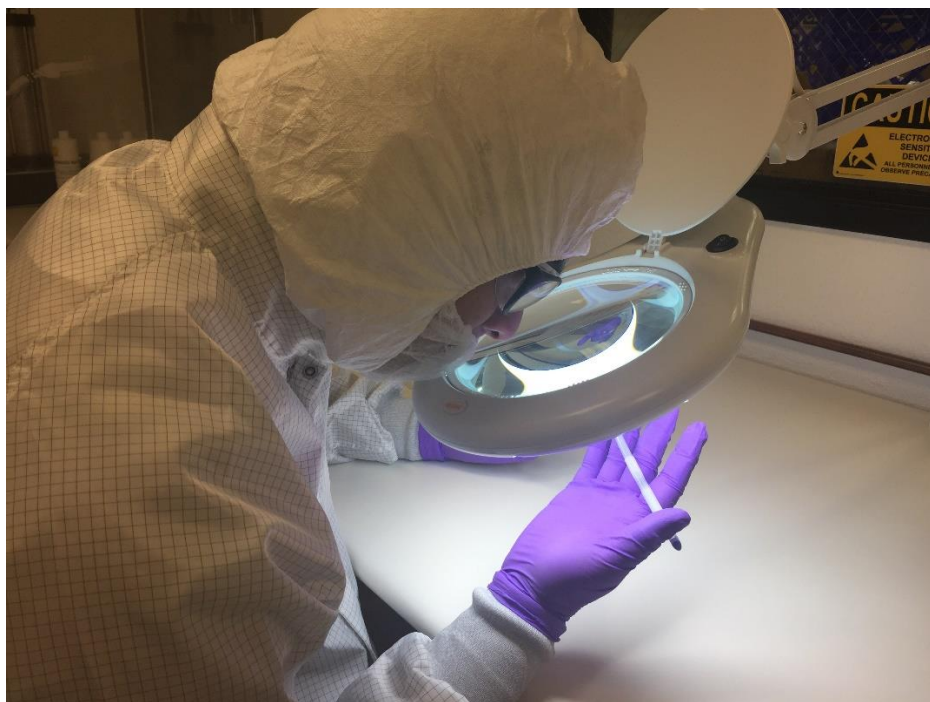


FIG. 4 – Catheter/Polyethylene Bag cleanliness verification

1. Check to ensure that the inside and outside of the 2" W x 18" L Polyethylene Bag is free of lint/debris (discard if necessary).
2. Insert the Catheter (tip end first) into the Polyethylene Bag until the tip of the catheter bottoms out to the closed end of the bag.



FIG. 5 – Catheter insertion to Polyethylene Bag



Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters

Doc No.: MP-SLQ-0001

Rev.: A

Page 6 of 24

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

1. Visually inspect the sleeved Catheter to ensure that the outer and inner surface is free of lint/debris.
2. **Check inside of the Labeled Pouch to ensure it is free of any foreign matter/debris.**
3. **Insert sleeved catheter (closed end of sleeve) into the pouch until the sleeve bottoms out to the chevron end of the pouch.**



FIG. 6 – Sleeved catheter/pouch cleanliness verify & insert

1. **Visually inspect the pouched catheter to ensure that NO debris was introduced when the sleeved catheter was inserted.**



FIG. 7 – Pouched Catheter cleanliness inspection

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters <i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Doc No.: MP-SLQ-0001 Rev.: A Page 7 of 24
--	--	---

1. Set/position the pouched catheter with the Tyvek (white side) facing up into the ACCU-SEAL Impulse Sealer.



FIG. 8 – View, sealing of pouched catheter utilizing impulse sealer

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 8 of 24

9.0 Procedure –10-Pack Inner Box and Outer Shipping Box Packaging and Labeling:

9.1 Assemble and insert the individually sealed Catheter Pouches into 10-Pack Inner Box into the Outer Shipping Box as shown in Figure 9 through Figure 14.

1. Fold in the bottom flaps of the 10-Pack Box.
2. Orient and apply the verified 10-Pack Box Label (ref. Figure 1) such that it is centered and reads towards the open end of the box.
3. Flip the box over and apply another 10-pack box label. Verify this label also has text reading towards the open end of the box.

1. Prior to boxing, check to ensure that there ARE NO VISIBLE signs of foreign matter (i.e. lint, hair....) in any of the sealed catheter pouches.
2. Orient and neatly place 10 (ten) of the sealed Catheter Pouches (Chevron end first) into the labeled 10-Pack box.

1. Insert the IFU P/N HDX-4004 into the loaded 10-Pack Box.
2. Close box flaps.



FIG. 9 – 10-Pack Inner Box Assy/Label



FIG. 10 – 10-Pack Inner Box Loading



FIG. 11 – IFU insertion into 10-Pack Box



Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters

Doc No.: MP-SLQ-0001

Rev.: A

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Page 9 of 24

1. Assemble the Shipping Box using Packaging Tape or equivalent.
2. Apply shipper box labels to adjacent sides of the box as shown.



FIG. 12 – Box and label application

1. Load a maximum of ten 10-Pack Catheter boxes into the shipping box.



FIG. 13 – 10-Pack Catheter Boxes load to shipping box

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters <i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Doc No.: MP-SLQ-0001 Rev.: A Page 10 of 24
--	--	--

1. IF NECESSARY, FILL VOIDS INSIDE BOX with Bubble Wrap or equivalent.
2. Seal box using packing tape.
3. Route labeled box to Finished Goods Inventory (FGI) or shipping.



FIG. 14 – Seal shipping box and route to Finished Goods or Shipping

10.0 ACCU-SEAL Impulse Sealer (ES-0099-01) Equipment Setup:

- 10.1 Reference Figure A through Figure D for the ACCU-SEAL Model 5300-15-B Impulse Sealer Controls ([See Step 10.2 for validated operating ranges](#)).

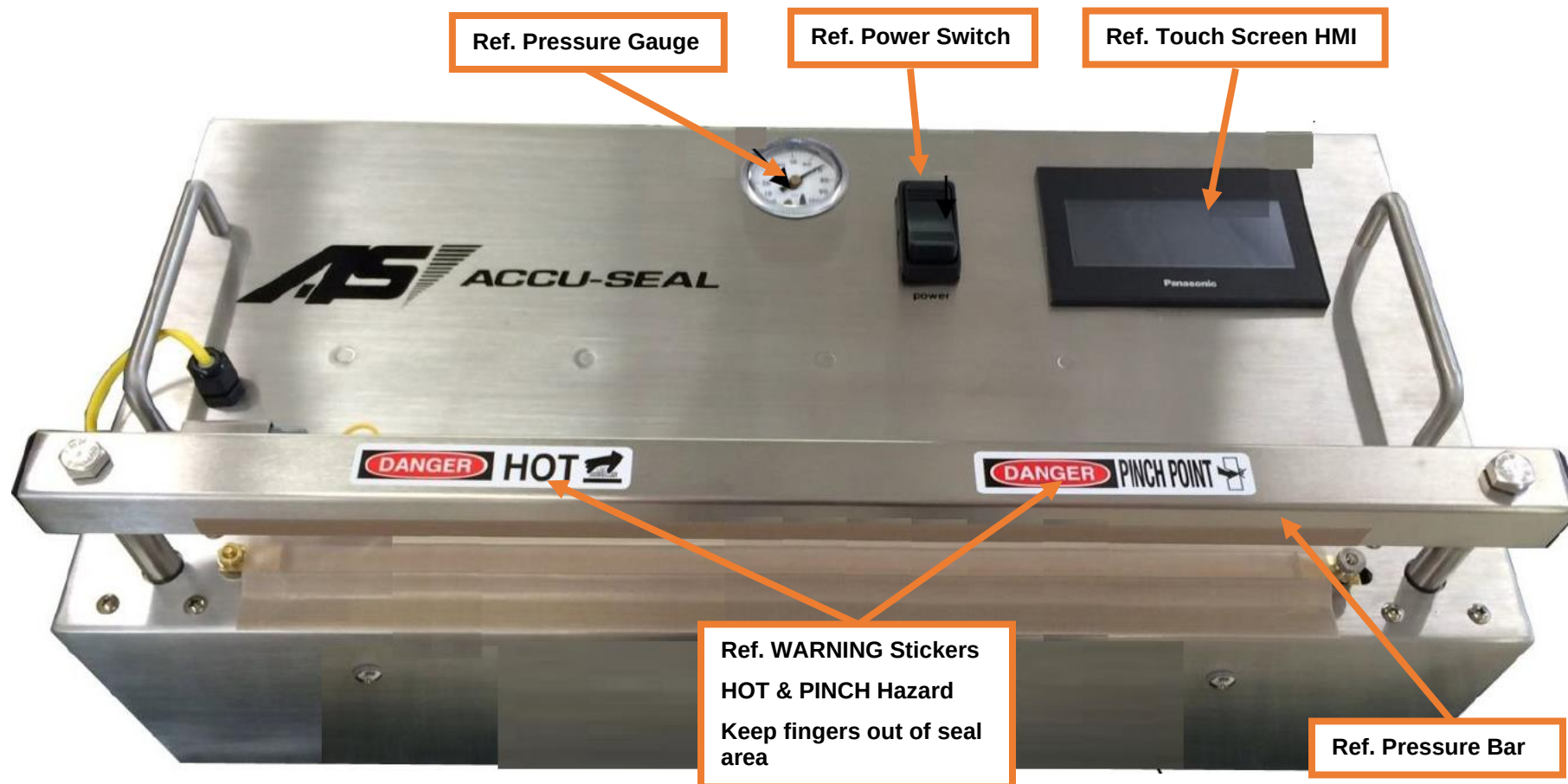


FIG. A – View from Front of Impulse Sealer

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters <i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Doc No.: MP-SLQ-0001 Rev.: A Page 12 of 24
--	--	--

1. Connect Foot Pedal, Power cords to appropriate ports
2. Connect compressed air line to the "AIR IN" port.
3. Use Pressure Regulator to adjust pressure to the Validated Range of 55 psi.

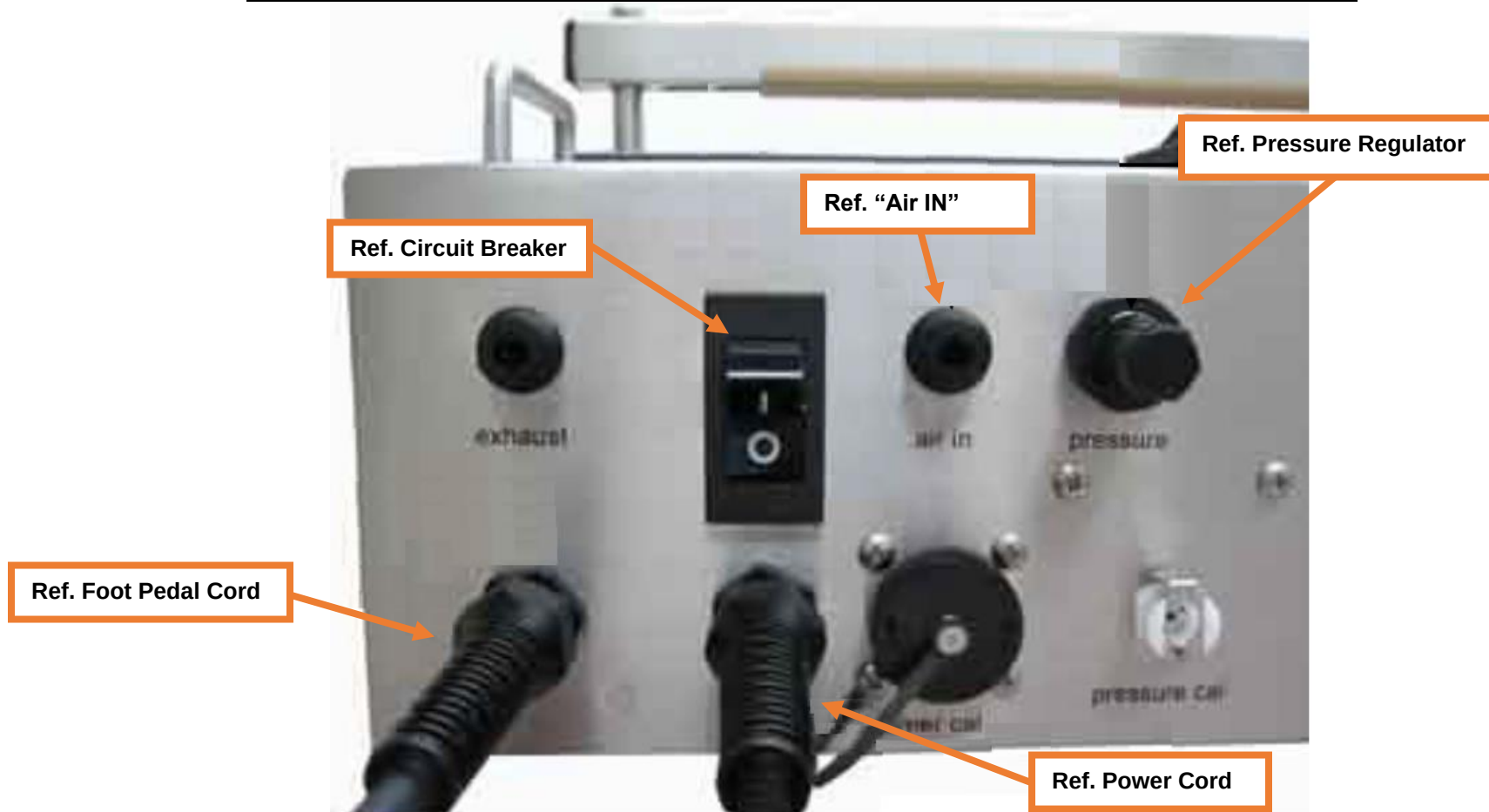


FIG. B – View, Back of Impulse Sealer

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 13 of 24

CONTROLS

Overview

The machine is programmed for two modes of control. Non-access code operation and Password Protected access. Non-access code operation only allows the operator to run the set seal parameters and to reset alarms.

The touch screen allows Authorized personnel, utilizing an Access code (password), access for editing all of the machine seal parameters, alarm settings, pressure settings, as well as the Calibration and Maintenance screens.

Controls and their Function:

- **Power Switch:** Turns power to the machine ON and OFF.
- **Foot Switch:** The Foot Switch is used to start the sealing cycle.
- **Circuit Breaker:** Breaks the Main Power to the machine when there is an overload condition.
- **Touch Screen HMI:** The Touch Screen HMI is used to Monitor, Set or Edit all of the machine parameters.

Controls Operation:

- Turn the machine ON by pressing the Main Power Switch located on top of the machine. When the Power is ON, the touch screen will be illuminated. To turn OFF Main Power, press the bottom side of the Main Power Switch.
- Current machine seal parameters are displayed on the main Run Screen. To edit or change machine, press the "Edit" button. When the Edit button is pressed, the "Enter Password" Screen will appear. The correct Access Code (password) must be entered to continue. When the correct code is entered you will have access to the functions screens and have the ability to change the machine parameters.

NOTE:

- Initial Set-Up Pass code has been programmed into the machine. See Initial Set up instruction in the maintenance Section. Use this code for first time access to set New password.

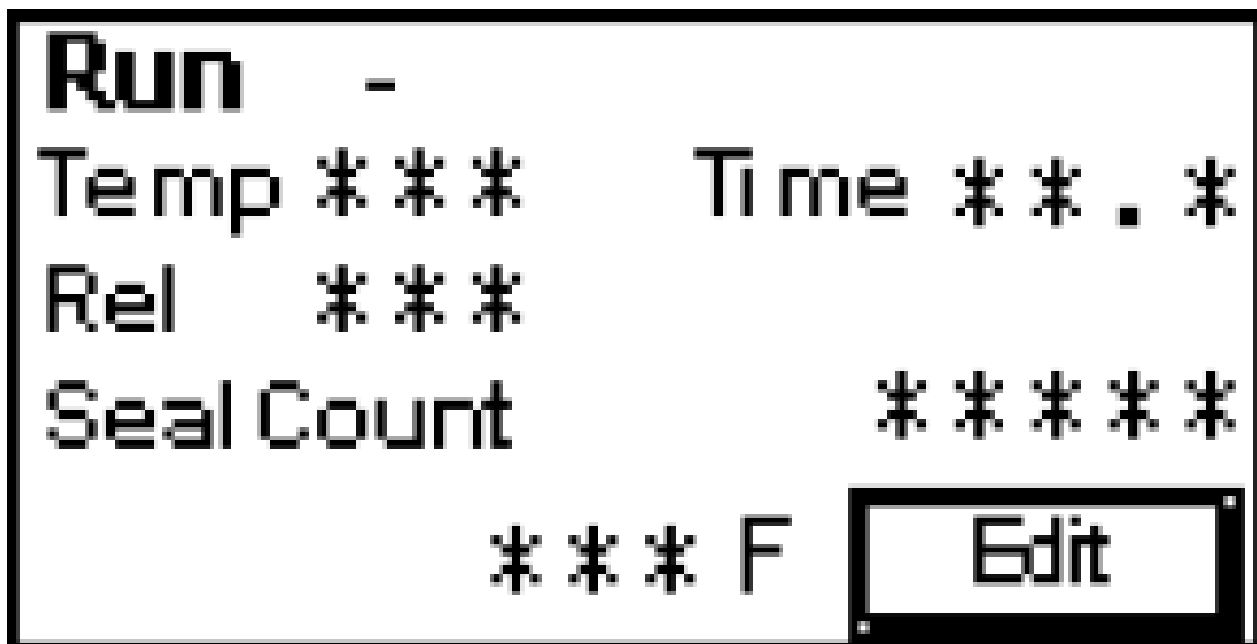


FIG. H – Touch Screen Functions

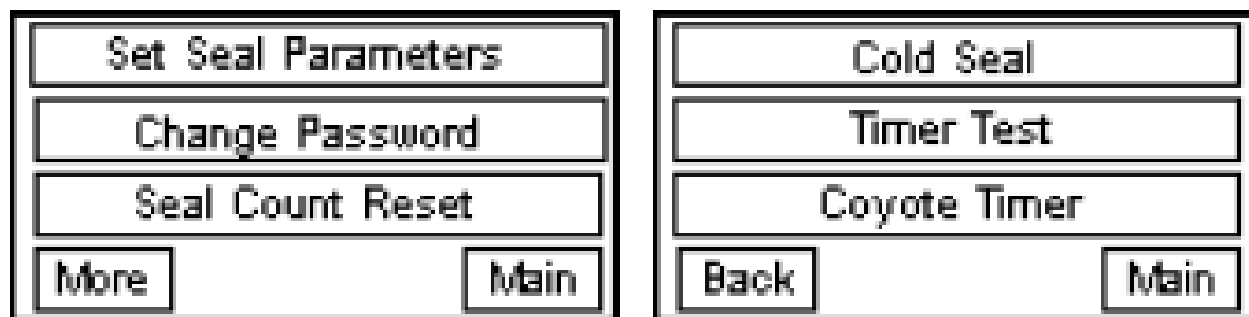


FIG. C – Edit Menu Screens

- **TEMP:** Displays set Seal Temperature
 - **REL:** Displays set Release Temperature
 - **TIME:** Displays set Seal Time
 - **SEAL COUNT:** Displays the number of heat seals completed
 - *****F:** Active temperature display. Display indicates the current temperature of the element.
 - [EDIT: Allows AUTHORIZED PERSONNEL, UTILIZING ACCESS CODE \(PASSWORD\), ACCESS TO ALL OF THE MACHINE SETTINGS:](#)
 - Seal Parameters
 - Change Password
 - Reset Seal Count
 - Alarms
 - Calibration Timer Test
- Edit Menu Screens**
- Display the functions that can be selected for Editing.



Seal Parameters

- This screen will allow you to change the seal parameters for:

- SEAL TEMP: Temperature held for the time specified by the Seal Time (Dwell).
- SEAL TIME: The amount of time the Seal Temperature is maintained.
- RELEASE: Temperature the seal bar releases and the cycle resets.

Other Function Screens

- Change Password
- Reset Seal Count
- Timer Test for Calibration
- Coyote Timer Setting

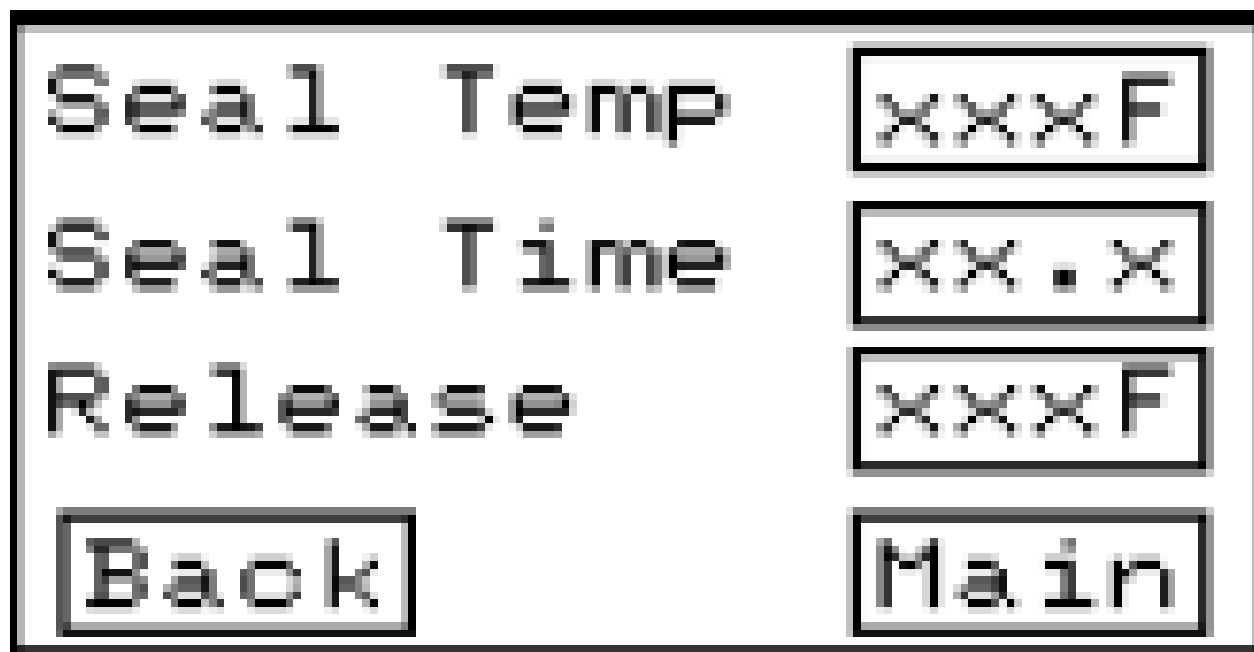


FIG. D – Seal Parameter Screens

Alarms

- JAW OBSTRUCITON: Indicates that for some reason the Seal Bar will not close completely.
- COLD SEAL: If the temperature does not reach the Seal Temperature within the set cold seal time, alarm will activate.
- OPEN SENSOR: Thermocouple is disconnected or has an open fault.
- OVERTEMPERATURE FAULT: Is a heater power control warning and indicates abnormal operation of the temperature control hardware.

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 16 of 24

ENTERING/EDITING MACHINE PARAMETERS: Any parameter that is enclosed in a box can be change. To change the parameter, press anywhere in the box of the parameter you want to change. This will then bring up a “keyboard” that will allow you to enter a new value. Once the desired value is entered, either the box of the parameter you are changing can be pressed or the “enter button on the Keyboard can be pressed to save the new value into memory. (Enter=the return arrow key) Parameter that are not enclosed in a box are for viewing only.

CHANGING ACCESS CODE: The Access Code can be changed, by pressing the “Change Password” button. To change the Access Code, press on the box with 00000. Enter a new Access Code. The number must be 5 digits and cannot start with ‘0’. The number you entered will show in the box. If the number is correct, Press the “ENTER” button to save the new Access Code to memory. The machine has been programmed with a Master Default Access Code in case your code is lost or forgotten. Please call if needed.

SEAL COUNT FUNCTION: The machine counts every seal that is started and successfully completed an displays this on the Main Screen. The machine will count up to 99,999 seals before it is automatically reset to “0”. The seal count can also be manually reset through the Seal Count Reset function.

PRESSURE: When the Regulator on the back of the machine is adjusted, the pressure gauge on the top of the machine will read the set pressure. Typical operating pressure should be above 50 psi.

JAW OBSTRUCTION ALARM: Indicates that for some reason the seal jaw will not close completely. If the alarm is activated, the machine will go into a lock-up mode. The machine will not operate until condition is corrected and the alarm is reset, by pressing the box that reads “ALARM RESET”, on the flashing red alarm screen. The Main Menu will now appear.

COLD SEAL ALARM: If the temperature does not reach the Seal Temperature within the set Cold Seal time, alarm will activate. If the alarm is activated, the machine will go into a lock-up mode. The machine will not operate until condition is corrected and the alarm is reset, by pressing the box that reads “ALARM RESET”, on the flashing red alarm screen. The Main Menu will now appear.

OPEN SENSOR ALARM: The machine continuously monitors the status of the machine temperature sensor. If the sensor is disconnected or has an open fault, the machine will alarm with the open sensor alarm screen. The sensor monitor is active whenever the machine is on so the alarm can occur even when the machine is in the resting state. The machine will not operate until condition is corrected and the alarm is reset, by pressing the box that reads “ALARM RESET”, on the flashing red alarm screen. The Main Menu will now appear.



COYOTE TIMER SETTING: The internal PLC monitors the health of the seal process to detect any abnormal operation of the temperature control hardware. This monitoring action is called a "coyote". The coyote may be implemented in the hardware or software, or a combination of both. It is not a safety feature, but acts to detect abnormal temperature control operation and shut down the machine. If an abnormal temperature event occurs, the Over Temperature Fault Alarm will display on a flashing red screen. Contact Maintenance, cycle the power of the machine to reset the alarm. If the alarm occurs again, CALL ACCU-SEAL

NOTE: See Over Temperature Fault Alarm.

TIMER TEST: This screen is used during calibration for testing the timer function.

OPERATING PROCEDURES

Start UP

1. Verify that the Air supply is connected.
2. Turn the machine ON by pressing the Power Switch.
3. The Touch screen will light up 'green'. When the boot-up sequence is complete, you will see the Main Run" Screen.
4. Verify that the nuts at each end of the element are tight.
5. View the Main Menu Screen and verify that the temperature and pressure readouts are reasonable.
6. Keeping your fingers away from the Seal Area, press and release the Foot Switch. Cycle the machine once or twice to warm it up.
7. Place the product to be sealed underneath the Seal Bar making sure the product inside the pouch is not protruding into the Seal Area.

CAUTION: Small products in the pouch that protrude into the Seal Area may be too small to cause a Jaw Obstruction Alarm, but could cause damage to element or a High temperature condition.

8. Keeping your fingers away from the Seal Area, press and release the Foot Switch. When the Foot Switch is pressed the machine will immediately begin the seal cycle.
9. When the seal cycle is complete the Seal Bar will open.

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 18 of 24

10. Remove the Product and inspect the seal.

11. At the end of each shift or production run, Press the Power switch to turn the machine OFF>

10.2 Validated Operating Range is as follows:

- **Sealing Temperature** **270 – 280 °F**
- **Sealing Dwell time** **4.0 Fixed**
- **Sealing Pressure** **55 psi Fixed**

11.0 Equipment Setup:

11.1 Reference Figure AA through EE for Emplex Model MPS6300 Continuous Band Sealer Controls.

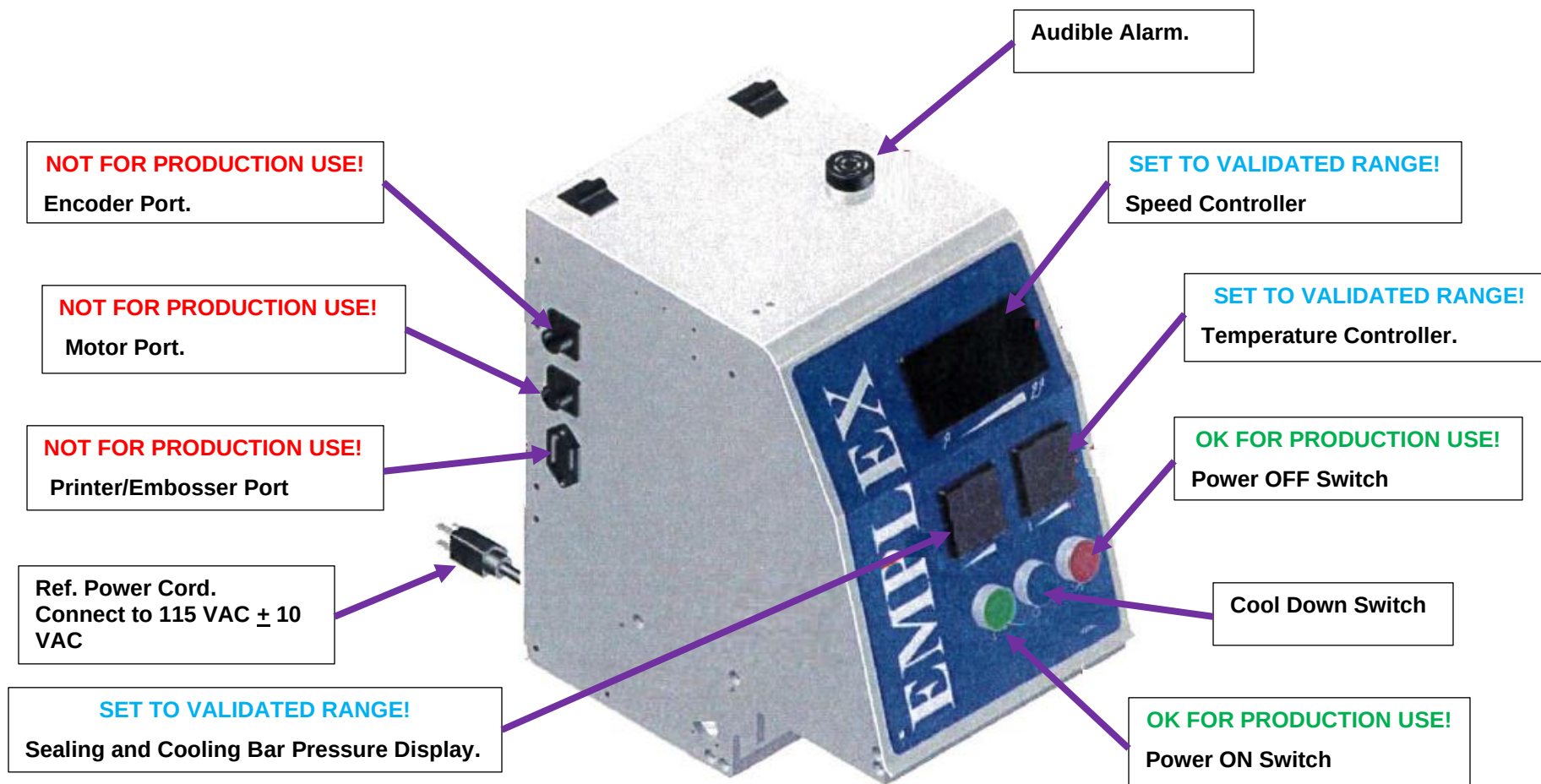


FIG. AA – Control Box Features

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 20 of 24

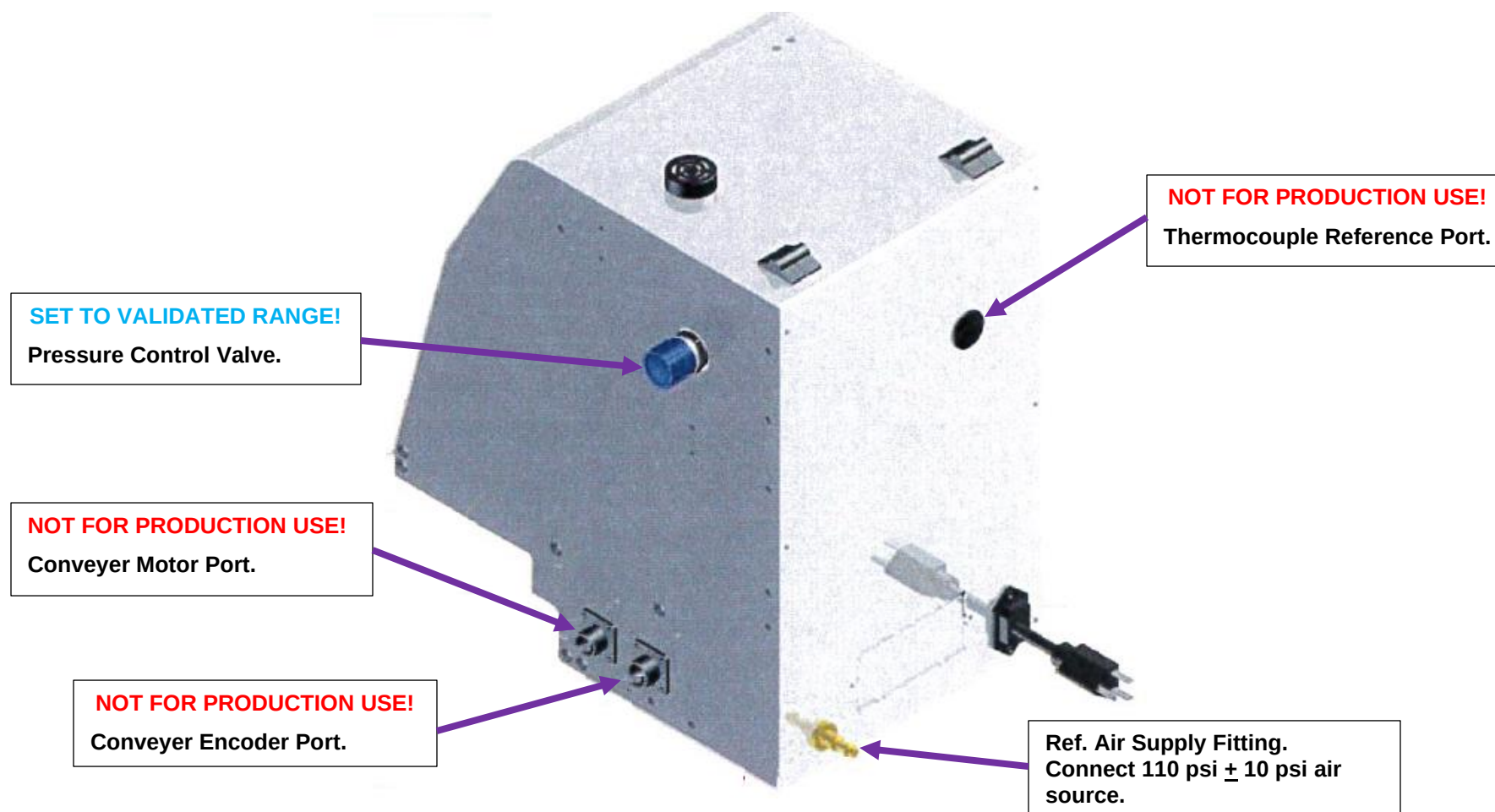


FIG. BB – Control Box Features continued

SETTING THE SEAL TEMPERATURE

1. Press "MODE" switch once.
2. Press UP ARROW (▲) or DOWN ARROW (▼) until Seal Temperature is within the Validated operating range.
3. Press "MODE" switch again to save.



FIG. CC – Seal Temperature Control

SETTING THE SEALSPEED

1. Press UP ARROW (▲) or DOWN ARROW (▼) until Seal Speed is within the Validated operating range.
2. Press "Enter" switch to save.



FIG. DD – Set Seal Speed

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 23 of 24

ADJUSTING SEALING AND COOLING BAR PRESSURE

1. Unlock the Pressure Control Valve by gently pulling the blue knob out from the control box.
2. Turn the Pressure Control Valve knob clockwise to increase pressure or counter-clockwise to decrease pressure (see Figure 12 for location of the Sealing and Cooling Bar pressure display) until Sealing and Cooling Bar Pressure is within the Validated operating range.
3. Once the desired setting is achieved, push the blue knob towards the control box to lock it in position.
4. The pressure display is identical to the temperature controller (ref. Figure 14). When the pressure regulator is adjusted the "Process Value" on the pressure display changes. If the "Process Value" deviates from the "Set Value" the Alarm will sound to alert the operator of a change in pressure.
5. Press "MODE" switch once.
6. Press UP ARROW (▲) or DOWN ARROW (▼) until "Set Value" matches the "Process Value" and is within the Validated operating range.
7. Press "MODE" switch again to save.

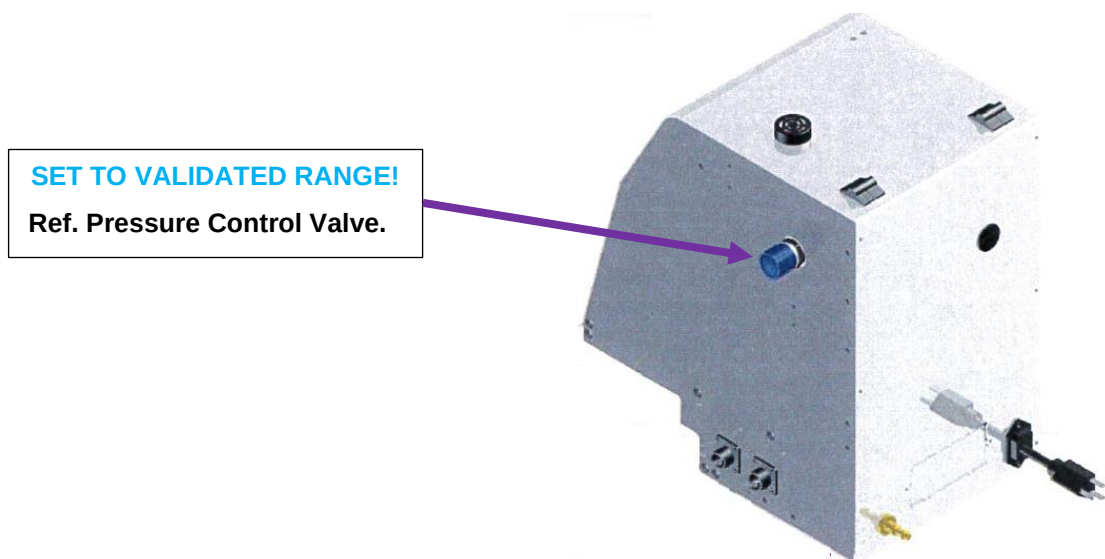


FIG. EE – Set Sealing and Cooling Bar Pressure

	Manufacturing Procedure, Labeling and Packaging, Silq ClearTract Catheters	Doc No.: MP-SLQ-0001
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 24 of 24

12.0 Change History

Revision	Change Description
A	Initial Release

Manufacturing Procedure – SILQ ClearTract Catheter Surface Treatment Using the Automation Machine




Note: This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, “For Reference Only”, unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work

	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Doc No.: MP-SLQ-0002
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 2 of 27

1. Purpose:
 1. Document the process for surface treating a SILQ ClearTract Catheter using the automation machine.
2. Scope:
 1. This procedure applies to surface treating of the SILQ ClearTract Catheter using the automation machine
3. References
 1. WI-0012 – Gowning and Prohibited Materials, Cleanroom
 2. WI-0014 - Line Clearance
 3. WI-0018 – Gowning & Prohibited Materials, for the Assy/Packing of Silicone...
4. Responsibilities:
 1. It is the responsibility of the individual(s) conducting these activities to comply with the requirements indicated in this procedure.
5. General Requirements:
 1. Line clearance shall be conducted per WI-0014 and PRIOR to any manufacturing activities.
6. Material and Equipment
 1. Automation Machine – ES-0091-01
 2. SILQ Polymer Solution Tank
 3. SILQ Clean Water & Waste Tanks
 4. Non-Latex Gloves
 5. Drying Rack
 6. Stopwatch, Calibrated


	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Doc No.: MP-SLQ-0002
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 3 of 27

Manufacturing Flow Chart

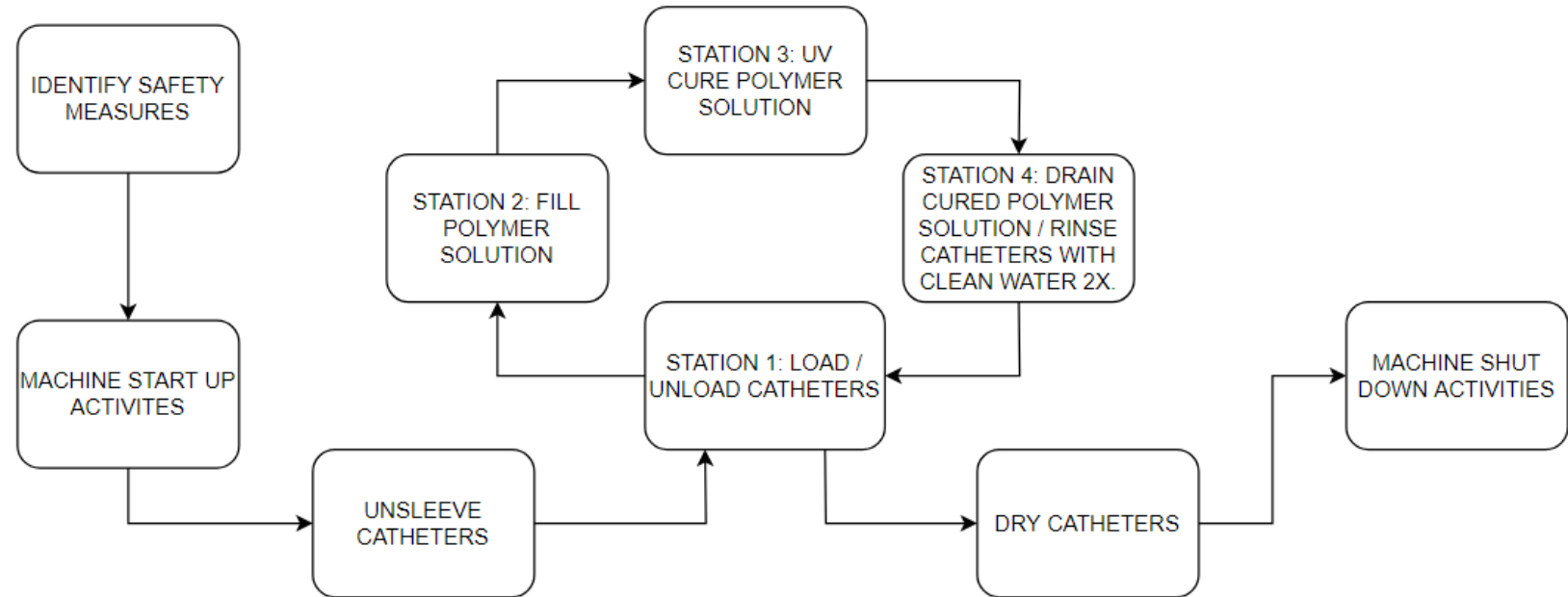

	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Doc No.: MP-SLQ-0002
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 4 of 27

Manufacturing Line Layout

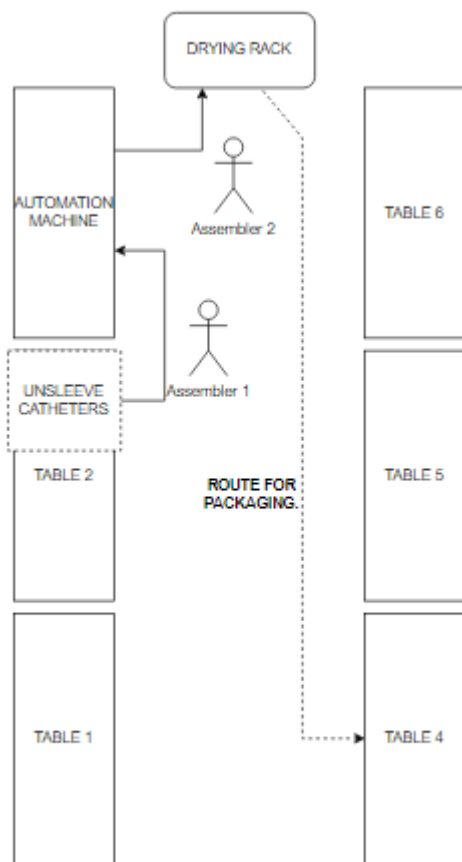



Procedural Notes

- The figure to the left depicts the manufacturing area set up for surface treating the catheters.
- It is recommended to use two assemblers.
- This procedure can be done under a different manufacturing line layout by the direction of the lead engineer or production supervisor.
- Manufacturing shall be done in a cleanroom (Room 104) environment at Pathway, located in Santee, CA, 92071.


	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Doc No.: MP-SLQ-0002
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 5 of 27

IDENTIFY SAFETY MEASURES

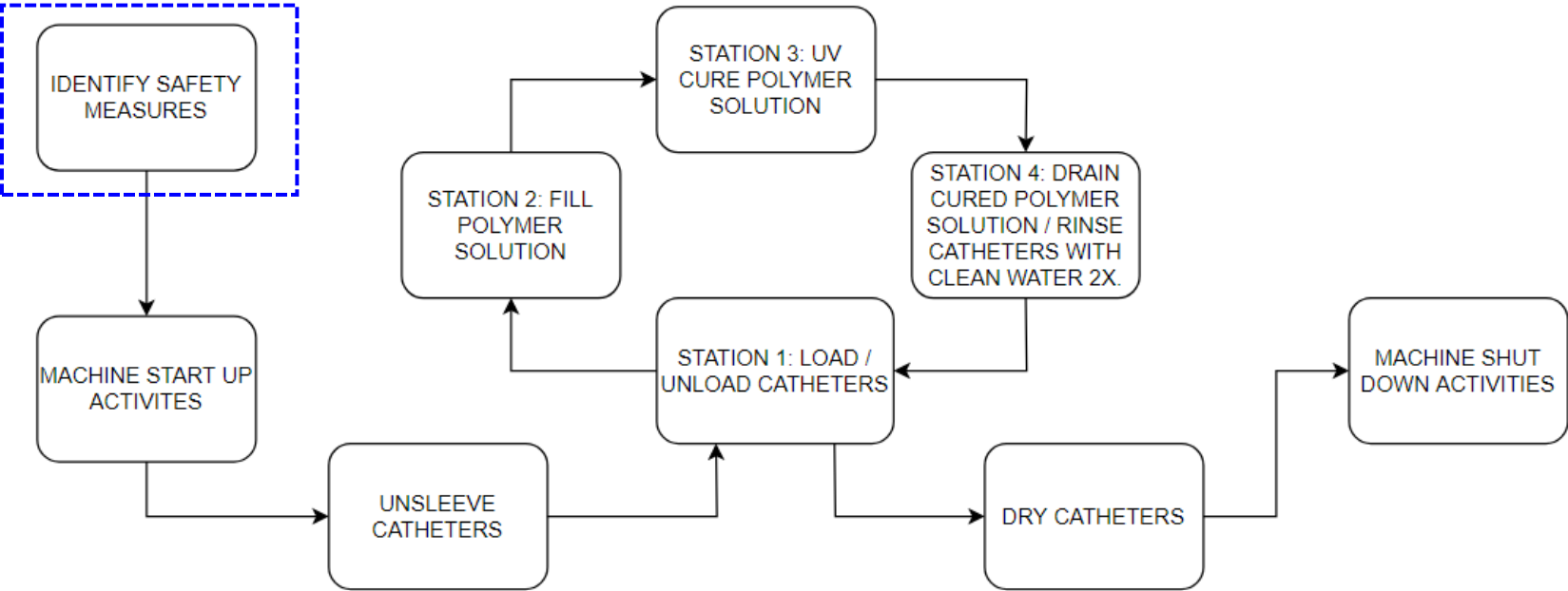

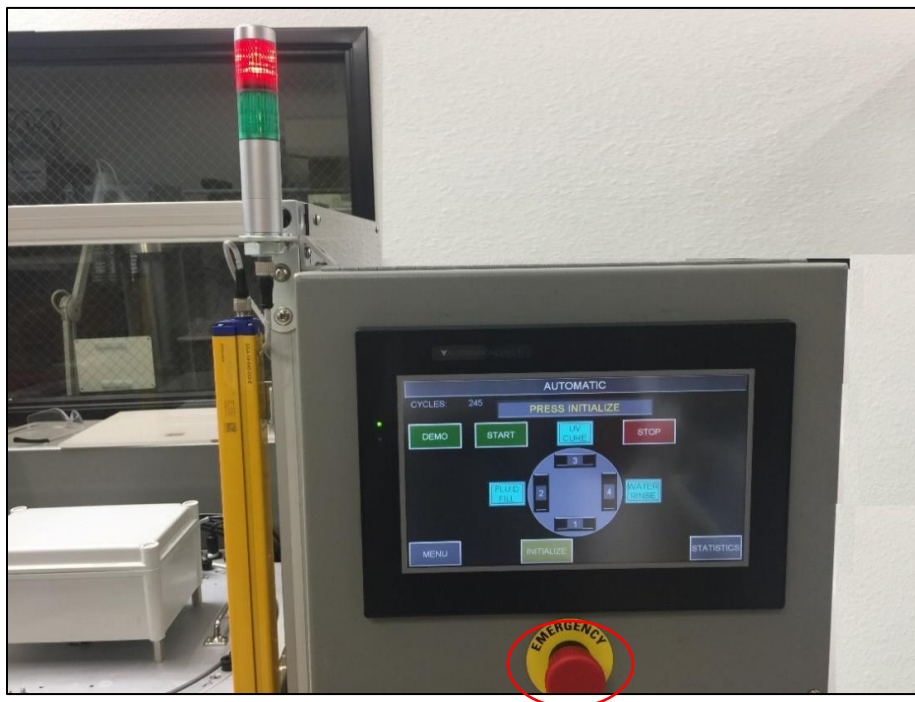


- Identify the emergency stop button. This button halts all machine activities.
- If pressed, consult the line supervisor or lead engineer before resuming manufacturing activities.

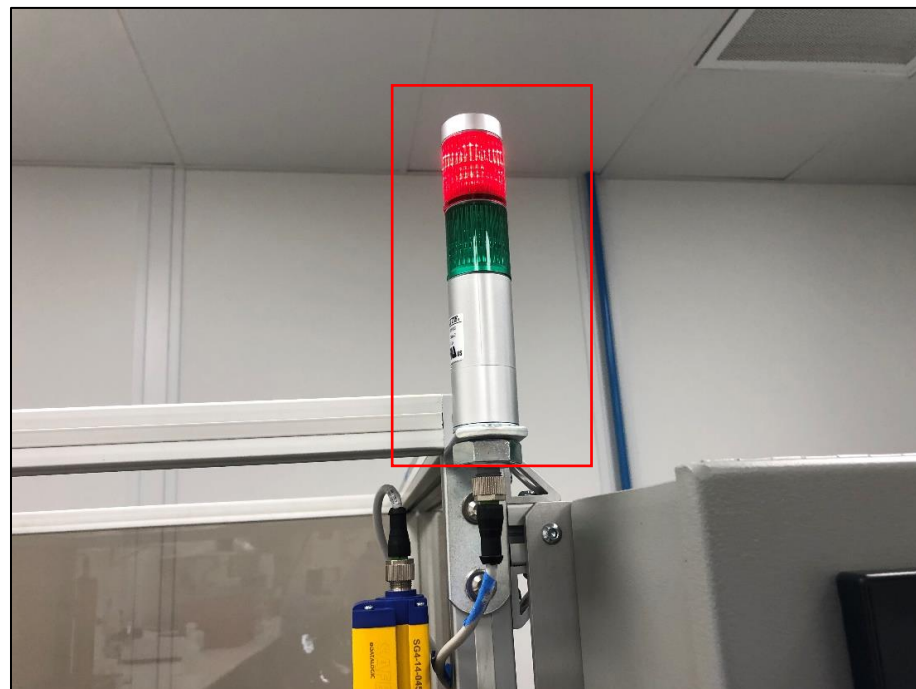


- Identify the light tower. The red light means that the machine is not currently going through any of its activities.

IDENTIFY SAFETY MEASURES

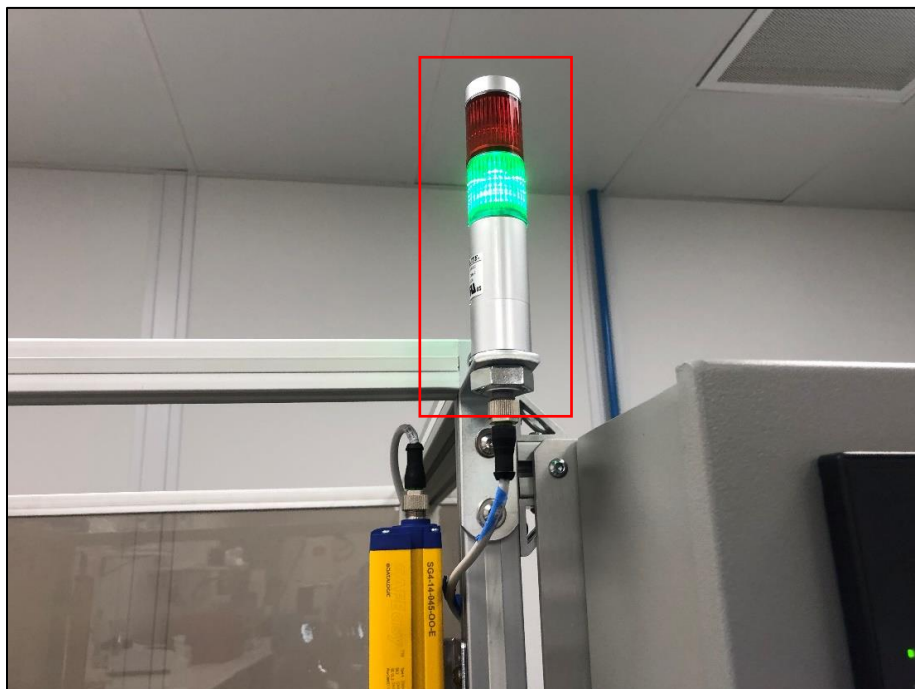

- Identify the green light on the light tower.
- The green light means the machine is performing one or some of the following activities.
 - Indexing
 - Filling
 - UV Curing
 - Draining / Rinsing

IDENTIFY SAFETY MEASURES

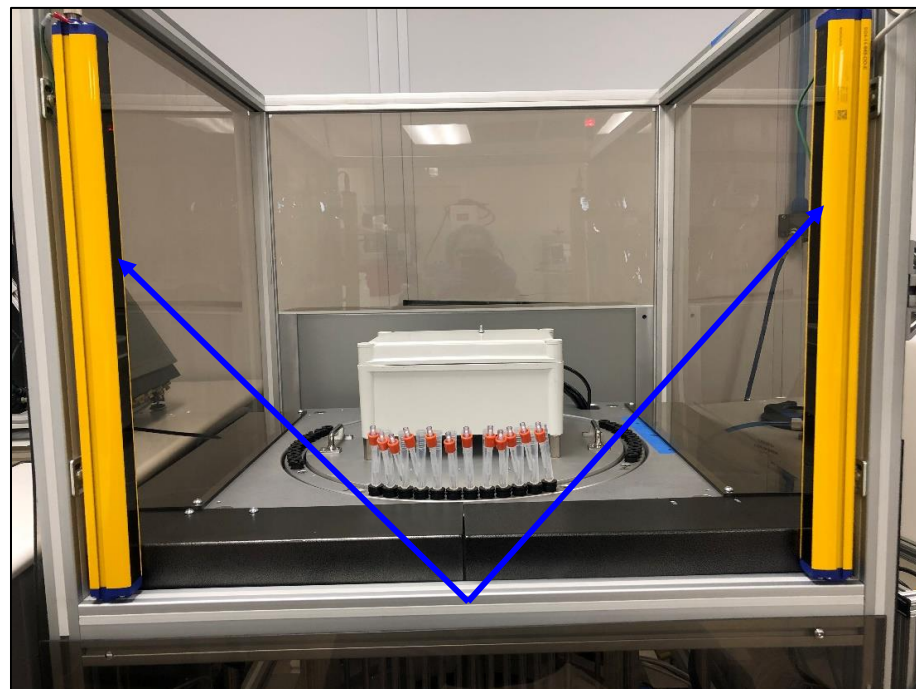

- Identify the motion sensors. These prevent the user from inserting hands into the machine while the machine is indexing.
- Note: when the machine is done indexing, the user can insert hands into the machine to retrieve catheters and insert new ones.
- Wait at least 3 full seconds after the machine has stopped indexing before inserting hands into the machine.



Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine


Doc No.: MP-SLQ-0002

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work

Rev.: A


Page 8 of 27

MACHINE START UP ACTIVITIES

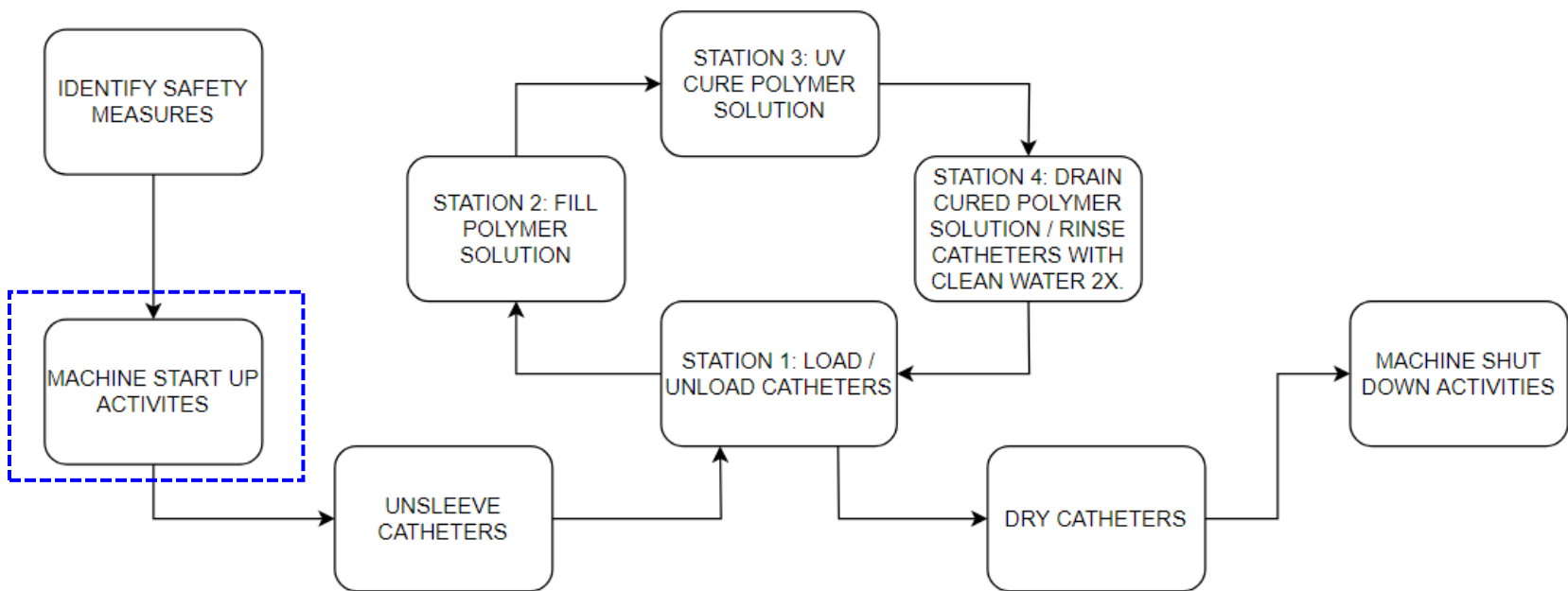




- On the rights side of the machine, locate the power switch.
- Turn it from off to on.

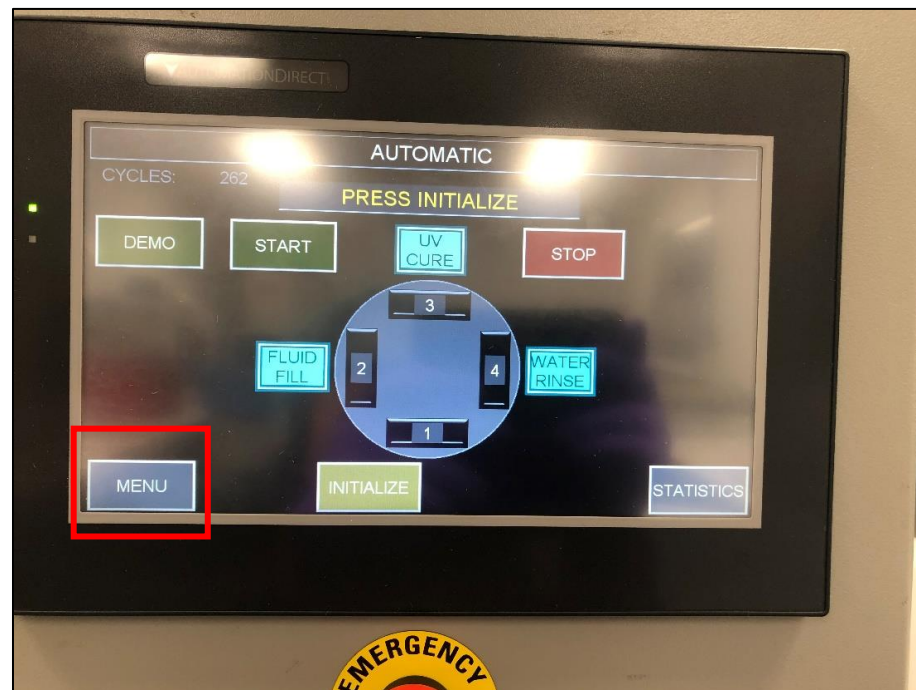


- After the machine starts up, it will be in the "automatic" mode.
- Press "MENU" to access the machine settings

MACHINE START UP ACTIVITIES

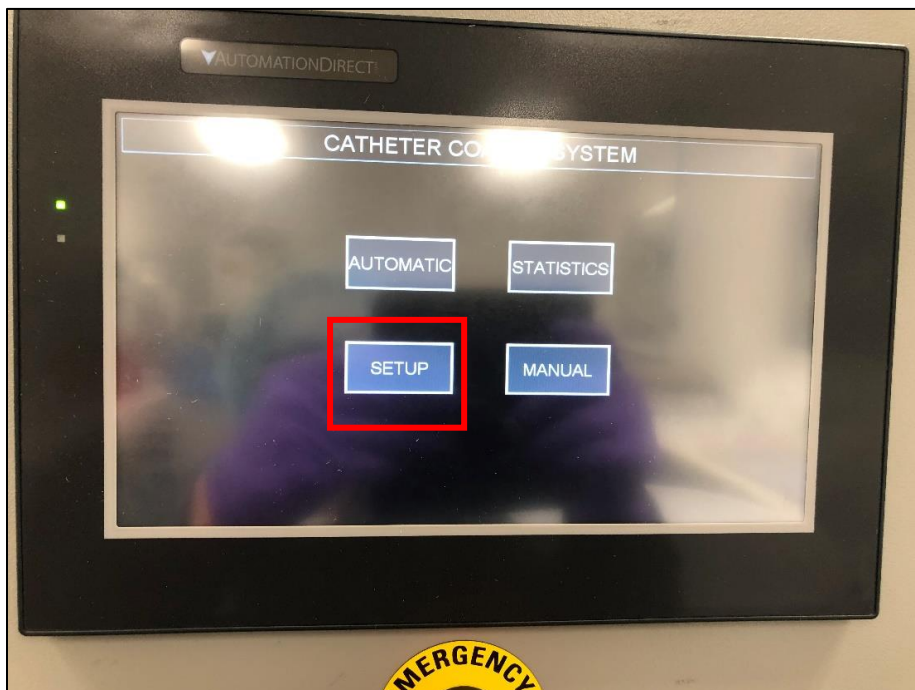


- Press "SETUP".


MACHINE START UP ACTIVITIES

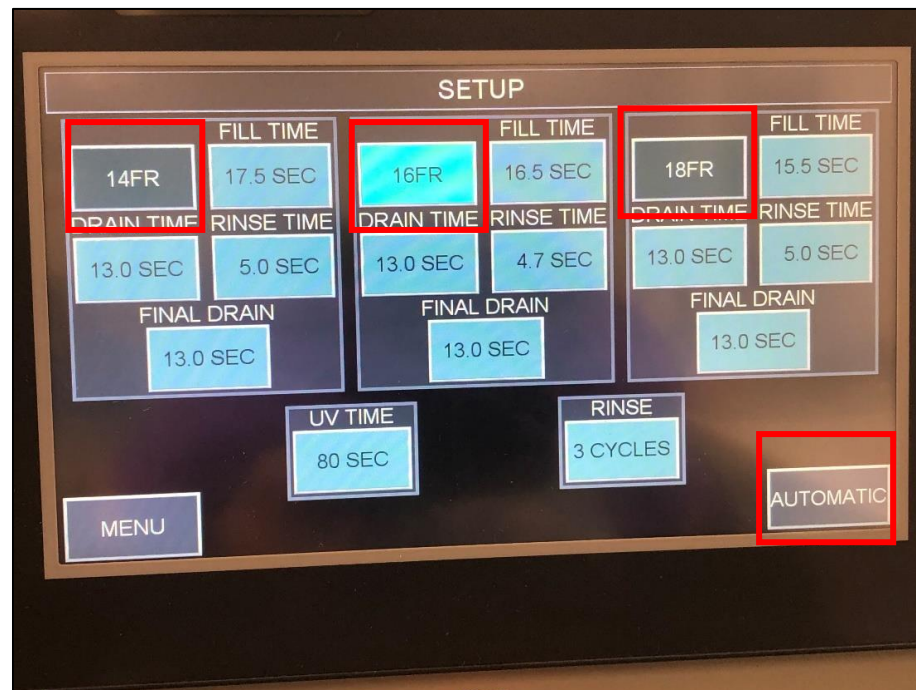

- Select the catheter size that will be used in the machine.
- Do not change any of the parameters.
- Record the parameters on the attachment, catheter run sheet **FM-SLQ-0001**.
- Press "Automatic" once finished
- After machine start up, and if during manufacturing activities, the settings need to change, consult the line supervisor or lead engineer.

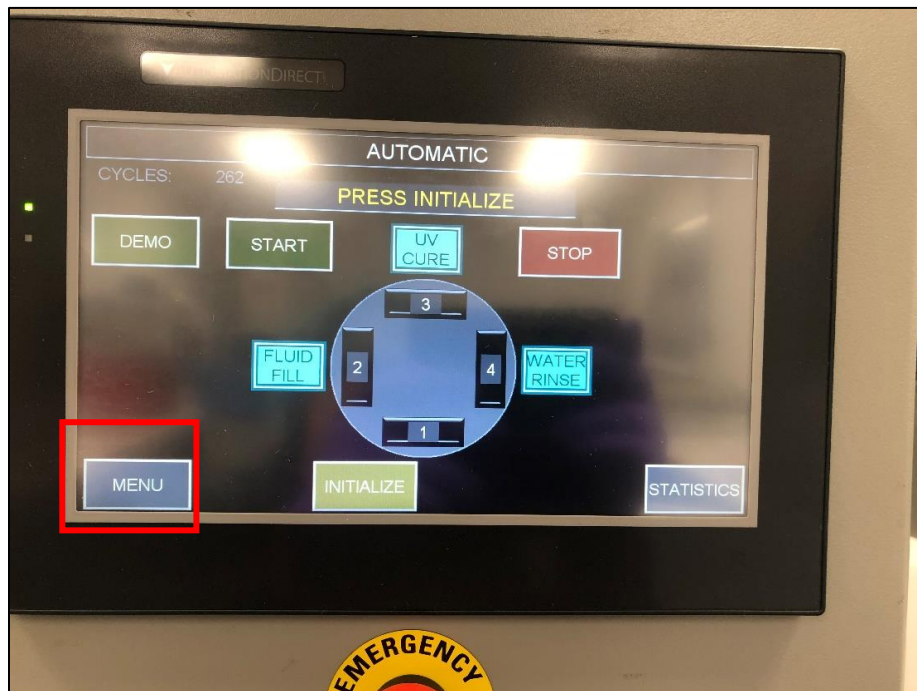


- Select "MENU."

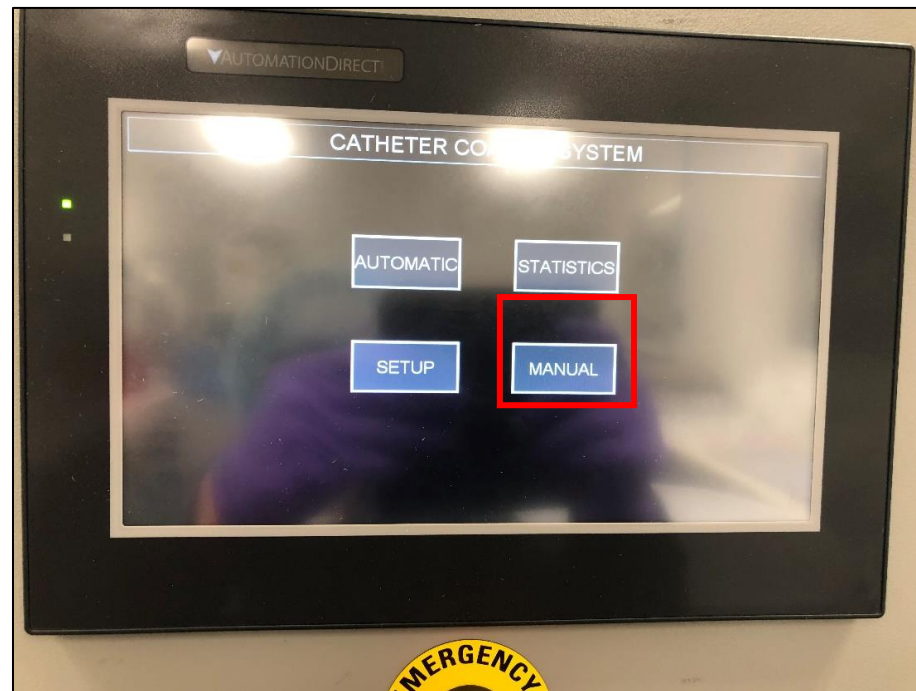


- Select "MANUAL." The automation machine is now out of automatic mode.

MACHINE START UP ACTIVITIES

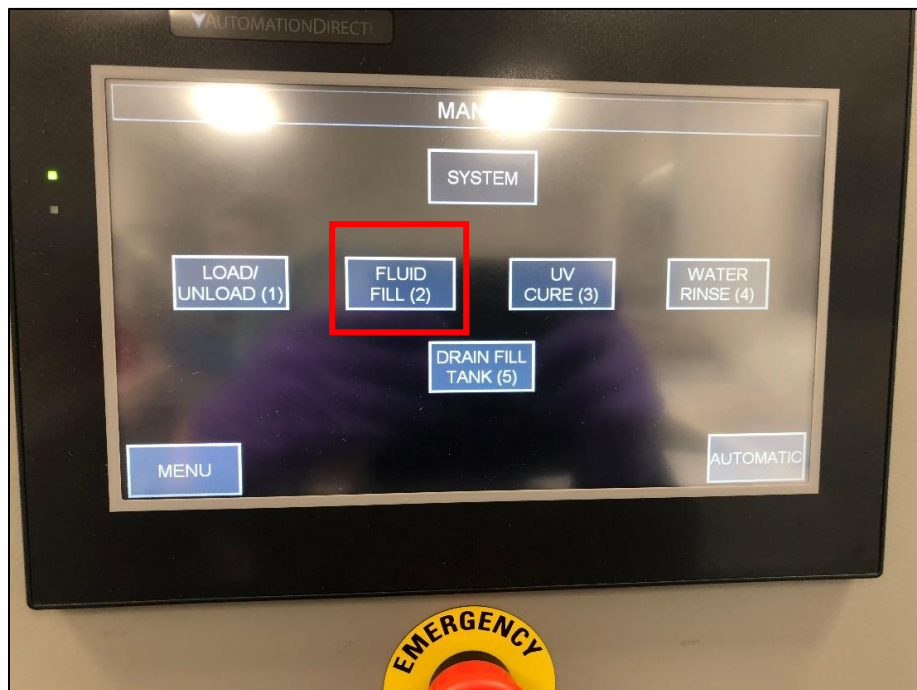

- Select "FLUID FILL (2)."

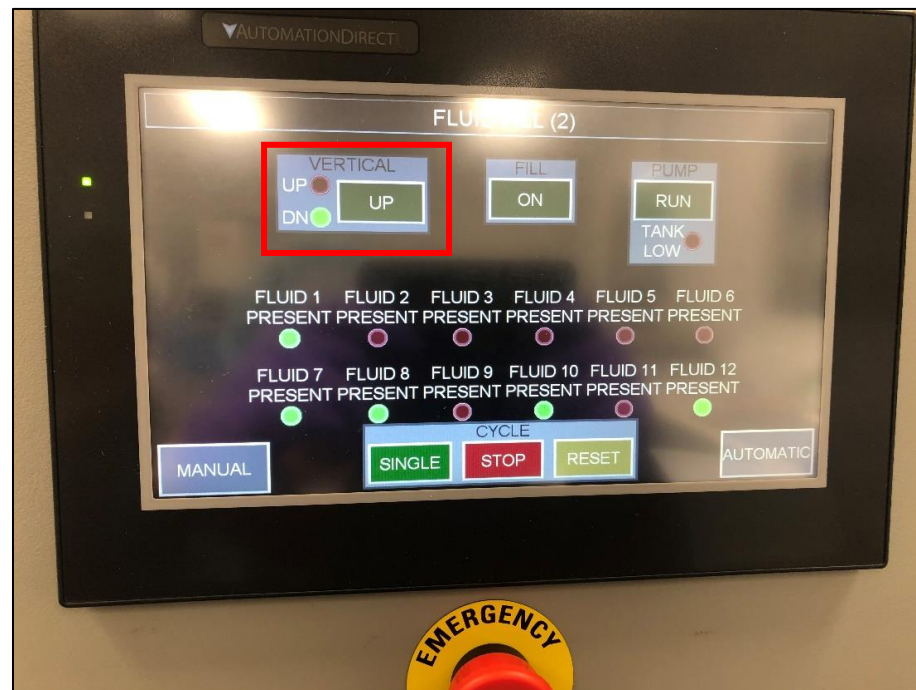


- Pressing the "UP" button will actuate the cylinder needed to fill the machine with the polymer solution.
- There will be an audible *click* coming from the cylinder moving into position.

MACHINE START UP ACTIVITIES




- If the cylinder actuates correctly, the "UP" light will activate to green.
- Press "DOWN" to actuate the cylinder back to its lower position.
- Repeat actuation 4 more times and verify proper function for each run.
- Press "MANUAL"


MACHINE START UP ACTIVITIES

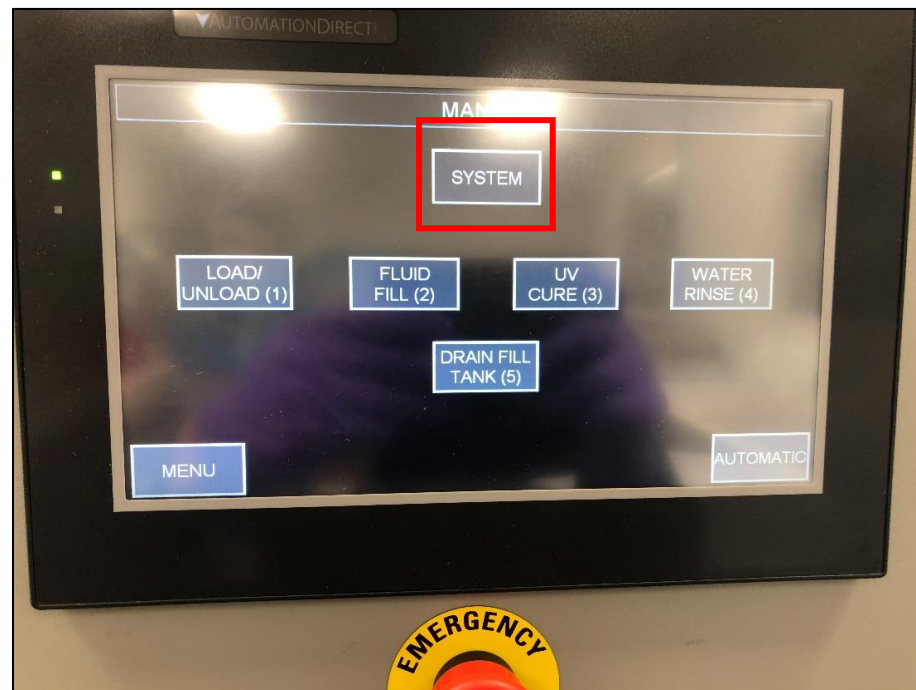

- Press "SYSTEM"

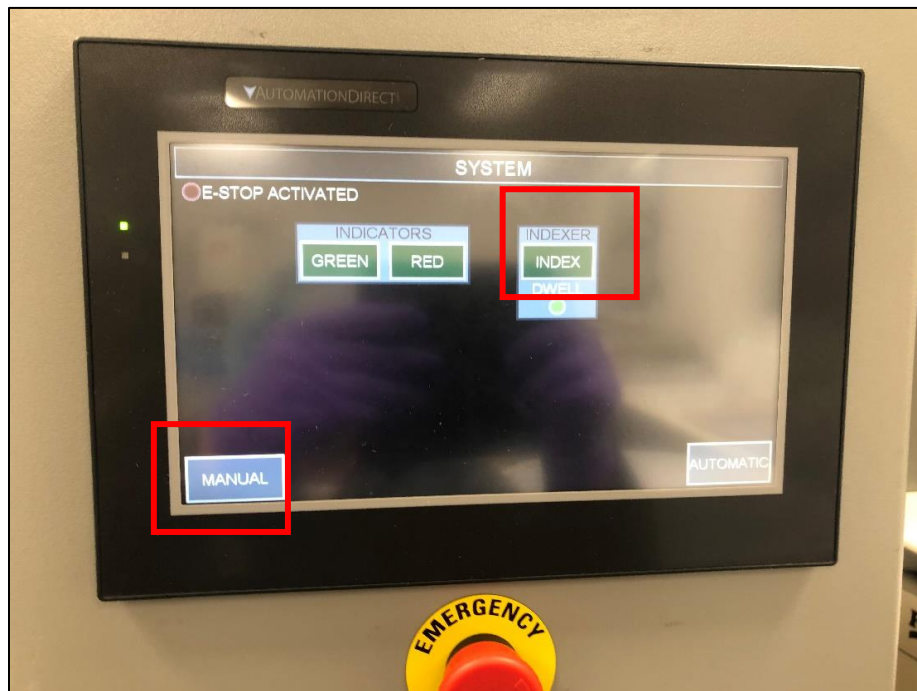


- Press "INDEX" to rotate the section to station 3. Wait 3 seconds.
- Press "INDEX" again to rotate the section to station 4.
- Press "MANUAL."

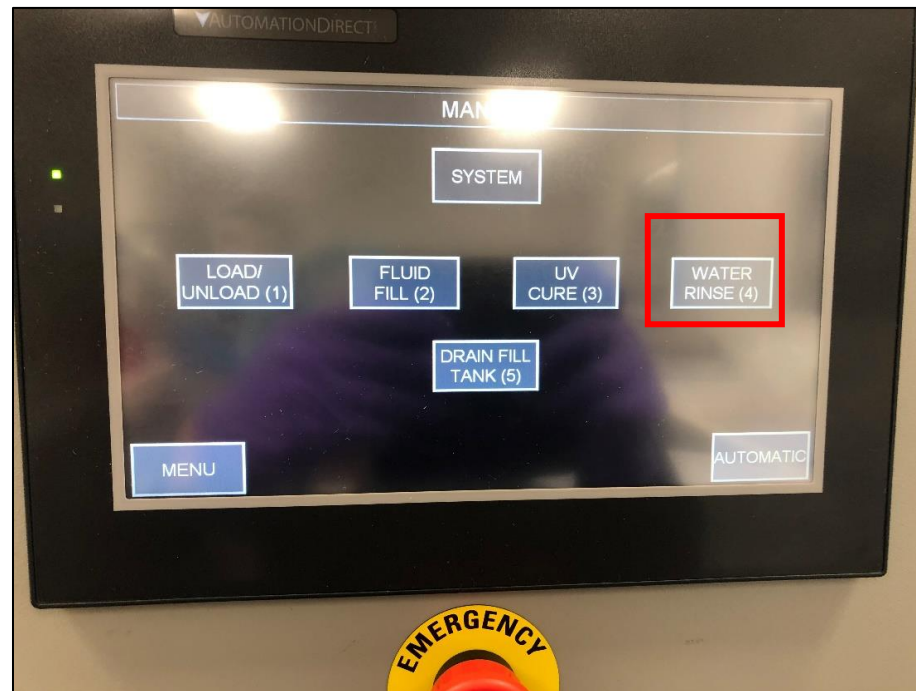

- Press "WATER RINSE (4)."


MACHINE START UP ACTIVITIES

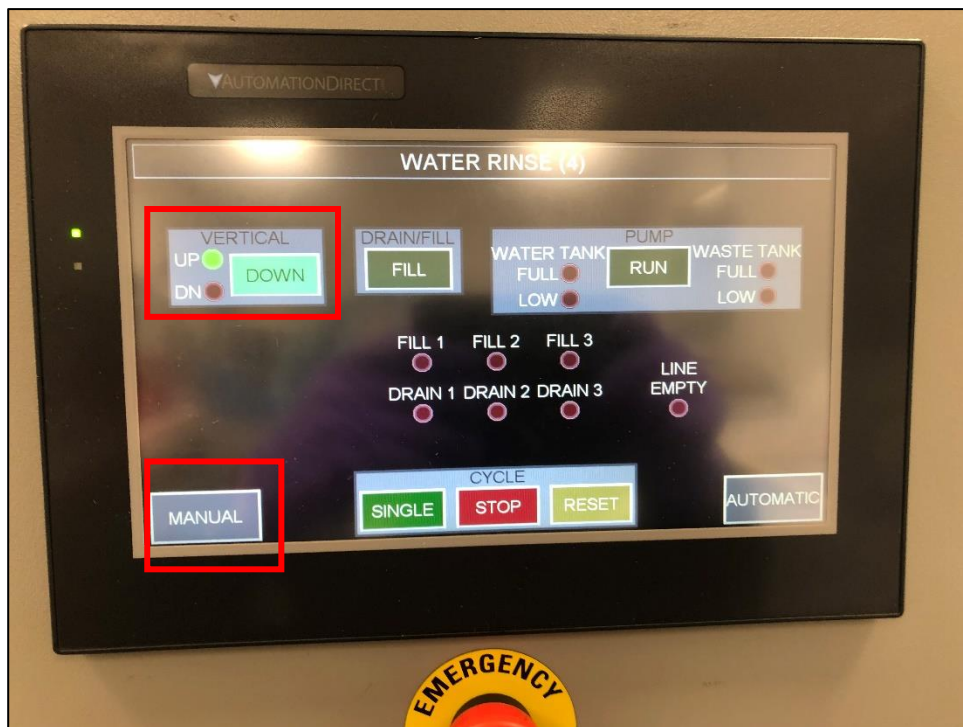

- Repeat the cylinder actuation steps from the fill side.
- Verify the cylinder correctly actuates 5 times.
- Press "MANUAL"

- There are 4 pairs of fittings on the automation machine that need to interface correctly with the cylinders at both the fill side and water rinse side.
- Using the methods already outlined. Verify the all 4 pairs of fittings interface correctly with both cylinders. Follow the flow detailed below.
 - Verify fitting 1 interfaces 5 times with the fill side cylinder. Index fitting 1 to the rinse side. Verify it interfaces 5 times with the rinse side cylinder. Record pass or fail on the run sheet.
 - Verify fitting 2 interfaces 5 times with the fill side cylinder. Index fitting 2 to the rinse side. Verify it interfaces 5 times with the rinse side cylinder. Record pass or fail on the run sheet.
 - Verify fitting 2 interfaces 5 times with the fill side cylinder. Index fitting 1 to the rinse side. Verify it interfaces 5 times with the rinse side cylinder. Record pass or fail on the run sheet.
 - Verify fitting 2 interfaces 5 times with the fill side cylinder. Index fitting 1 to the rinse side. Verify it interfaces 5 times with the rinse side cylinder. Record pass or fail on the run sheet.


MACHINE START UP ACTIVITIES

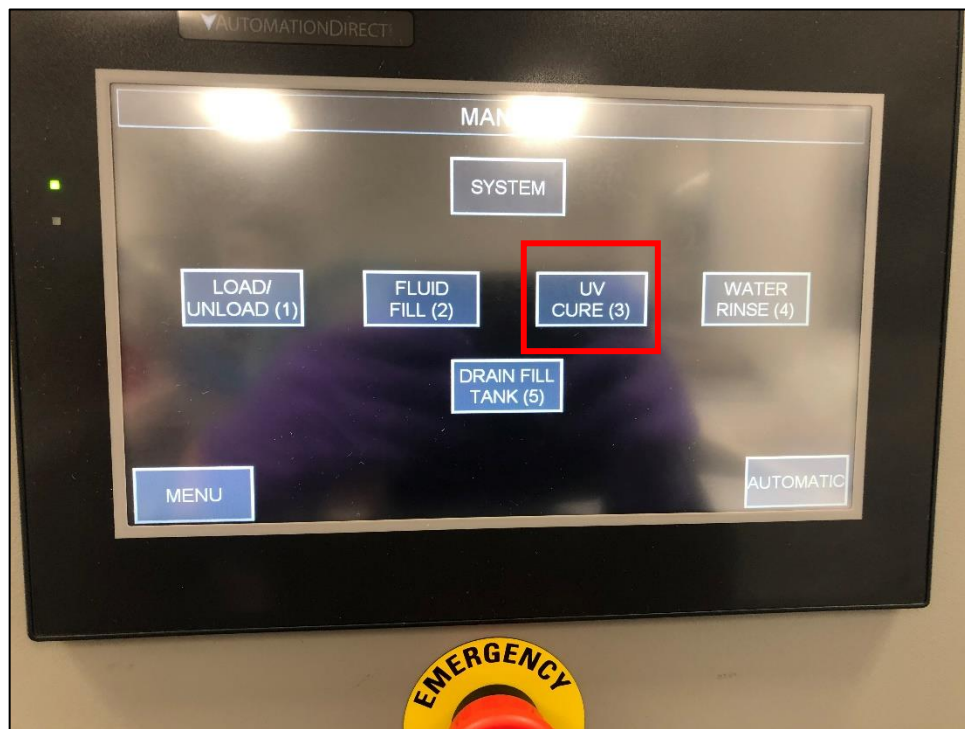


- Press "UV Cure (3)."
- Activate the UV lamps by pressing "ON."
- Let the lamps heat up for approximately 10 minutes before cycling any catheters through the machine.

MACHINE START UP ACTIVITIES

SURFACE TREATMENT ACTIVITIES

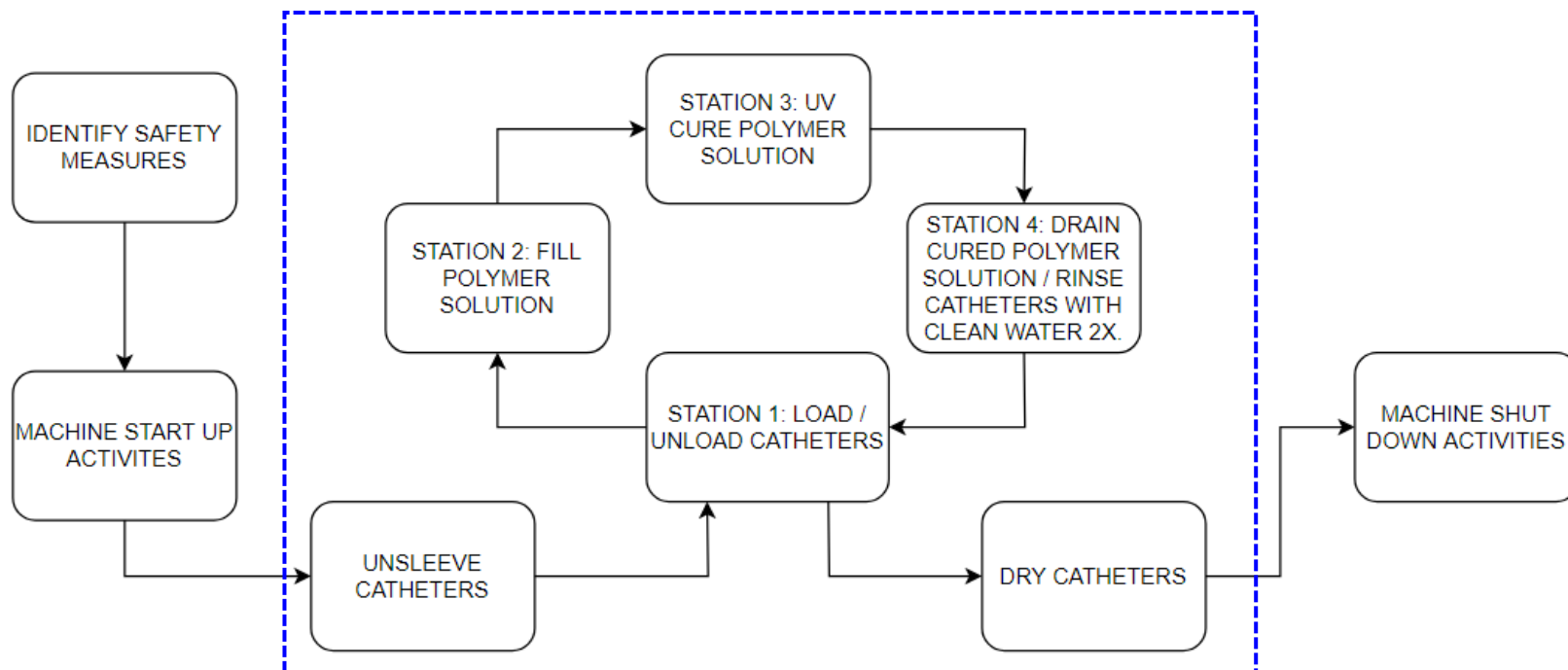


This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work

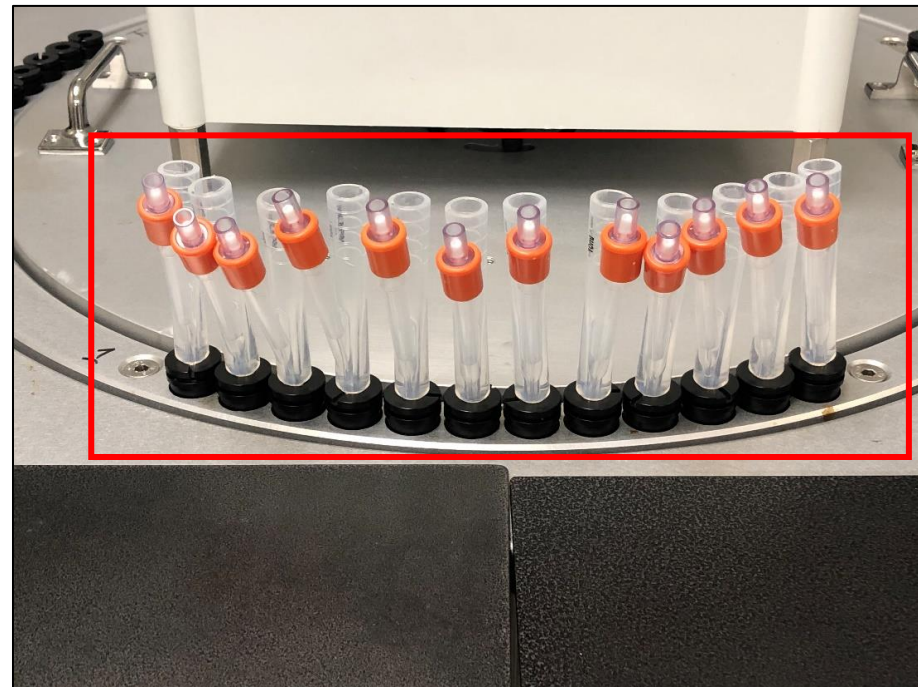


- Remove a catheter from its polyethylene sleeve.
- Discard the sleeve.
- Insert the catheter into the far-left spot.
- Repeat for 11 more catheters, moving one spot to the right, until the entire set is full.

- Rotate the catheters in the orientation shown above.

SURFACE TREATMENT ACTIVITIES

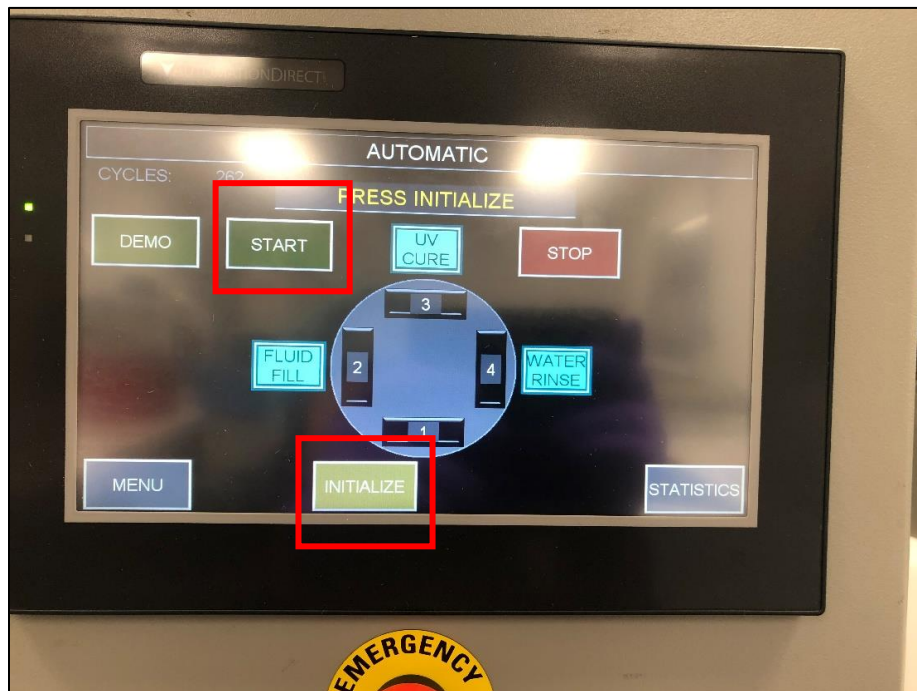


- Press "INITIALIZE" then press "START."

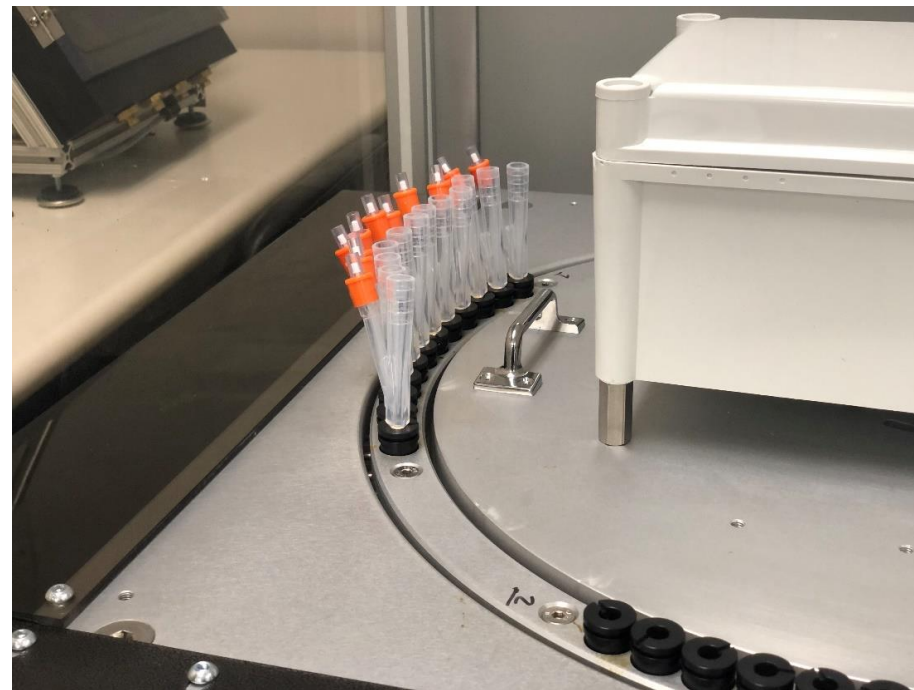


- The first set of catheters will rotate to the fill station.

SURFACE TREATMENT ACTIVITIES

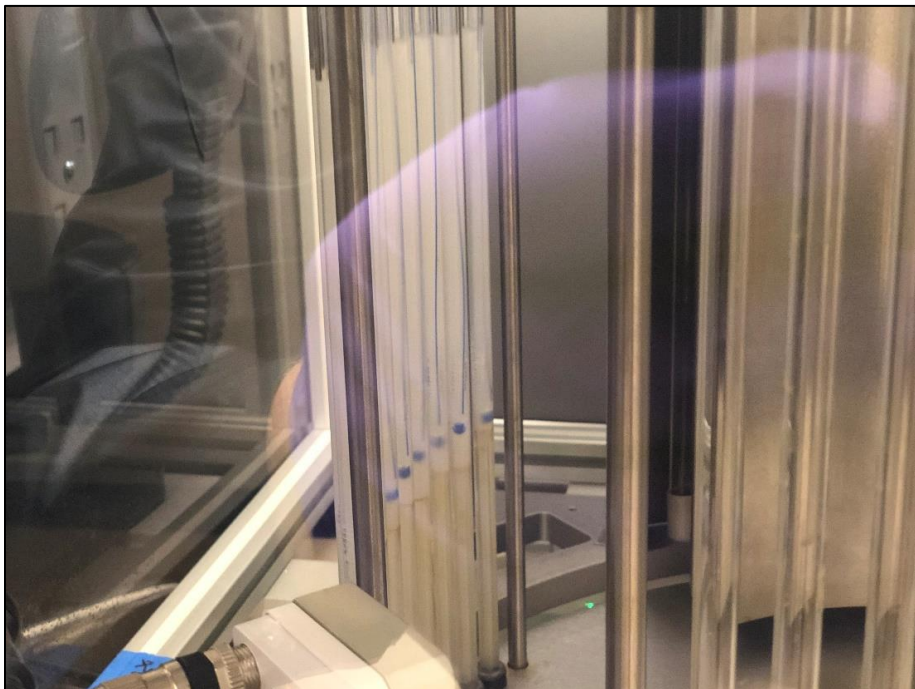

- Polymer solution will start filling the tubes.

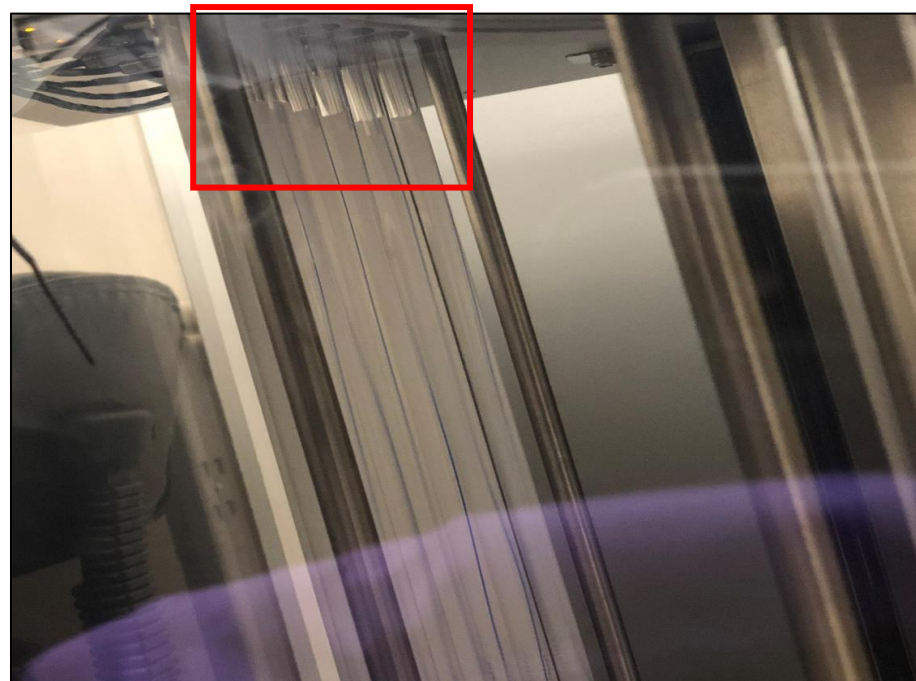


- The solution will stop filling at approximately the height shown above.
- While the solution is filling, remove 12 more catheters from their sleeves and insert them into the next set of 12 fittings.

SURFACE TREATMENT ACTIVITIES



- Once the next set of 12 is in the machine, press "START."
- The first set of catheters will rotate to the UV cure station while the second set of catheters will rotate to the fill station.

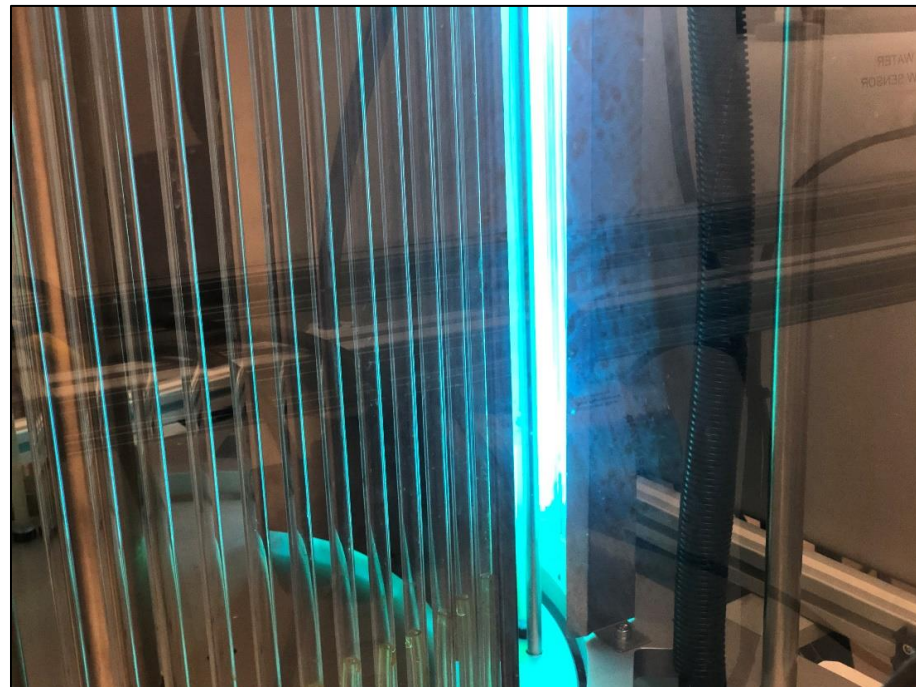


- While the first set of catheters is curing, insert a third set of catheters into the machine in the same way as the before sets.

SURFACE TREATMENT ACTIVITIES



- Once the UV curing step is complete.
- Press "START" on the machine to rotate the first batch of catheters to the drain/rinse station.
- The pumps will drain the cured solution into the waste tanks.

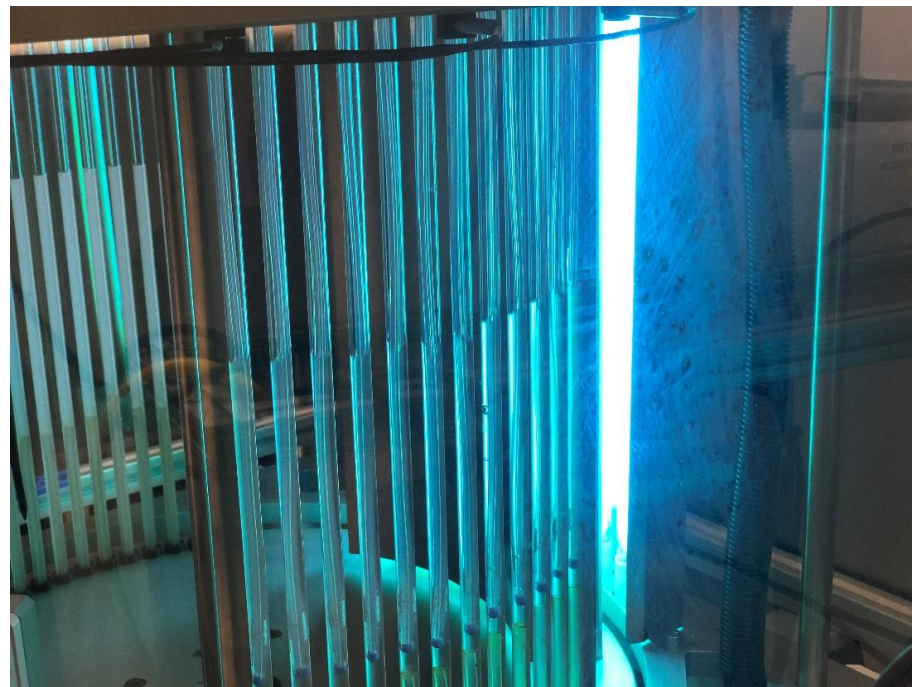

- The pumps will then push clean water into the tubes to rinse residual polymer solution. The pumps will then drain this water into the waste tanks.
- This process will repeat 2 or 3 times to ensure catheters fully rinsed.
- During this time, insert the fourth set of catheters into the machine.


SURFACE TREATMENT ACTIVITIES

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work



- After all machine activities are complete, press "START", and the fourth set of catheters will rotate to the fill station and first set will rotate to the load / unload station.
- Remove all 12 catheters and load them onto the drying racks as shown to the right. Insert a new set of 12 catheters after removing the finished set. From now, each machine rotation will involve unloading catheters and replacing with a new set.
- Do not place more than 500 catheters on the drying rack.
- Each section on the rack can hold 13 rows of 4 catheters.

SURFACE TREATMENT ACTIVITIES

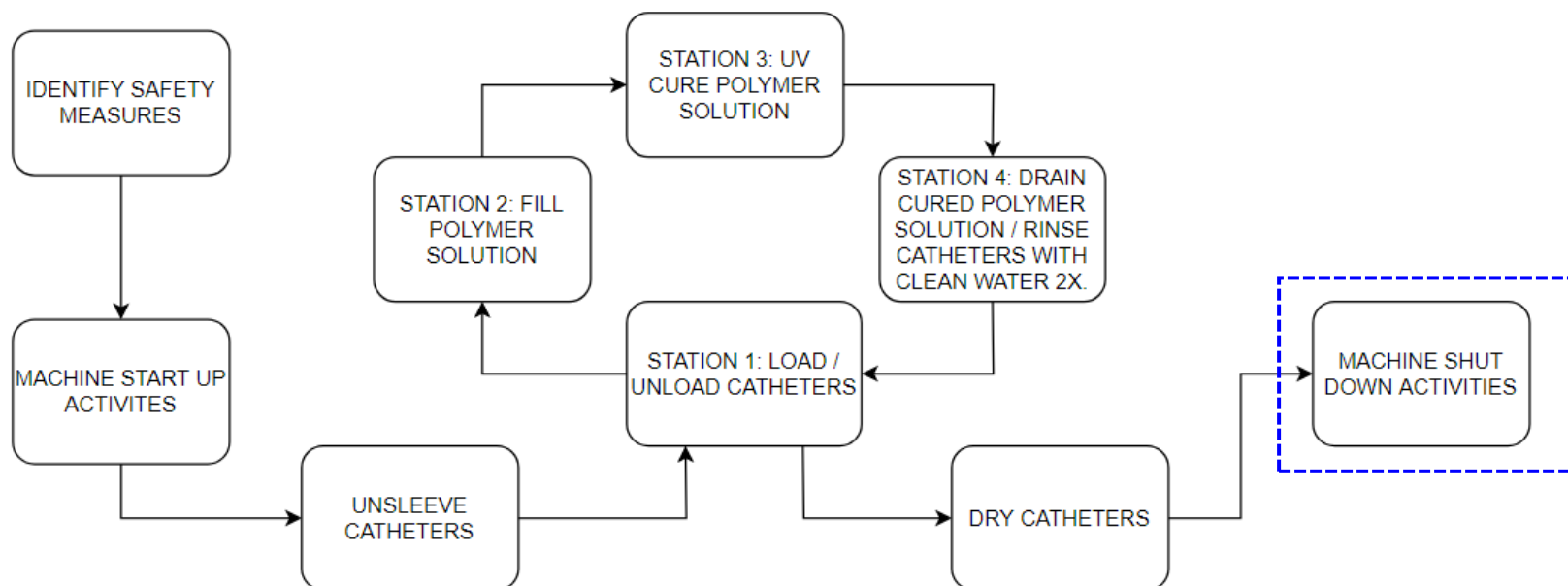

	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Doc No.: MP-SLQ-0002
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 25 of 27

- Perform the following as machine shut down activities
 - Fill each set of tubes to the top with clean water and drain it completely out three times.
 - Keep the power switch at "ON" and the mode in "AUTOMATIC."
 - Clean all surfaces with 70% IPA and a lint free wipe.
 - Remove the front panel from the automation machine and clean all interior surfaces with 70% IPA and a lint free wipe. Reinstall front panel.
 - If the waste tank is near the upper sensor, turn the power off, disconnect all fluid lines going to the waste and clean water tanks. Disconnect all power cables as well.
 - Remove the cart from the clean room to drain all waste fluids.
 - Also, if the clean water tank is nearing the lower sensor, replenish its supply with RO or DI water to just below the upper sensor.
 - Bring the cart back into the gowning room. Clean all surfaces with 70% IPA and lint free wipes.
 - Place the cart back into its position in the clean room. Connect the fluid lines and the power cables. Turn the power switch to "ON" and the mode to "AUTOMATIC."
- Perform these additional shut down activities
 - Clear all discarded polyethylene sleeves.
 - Leave catheters in the drying rack to dry overnight

MACHINE SHUT DOWN ACTIVITIES

	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Doc No.: MP-SLQ-0002
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 26 of 27

- After overnight drying, route the catheters for packaging and labeling.
 - Ensure the catheter run sheet is present.
- Return all equipment back their designated location.
- Return all unused materials back to their designated location or to the line supervisor.
- Clean all work surfaces with 70% Isopropyl Alcohol (IPA), clean cloths and/or lint free wipes.

- END OF MANUFACTURING PROCEDURE-

Revision	Change Description
A	Initial Release

	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Doc No.: MP-SLQ-0002
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled, "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work</i>	Rev.: A
		Page 27 of 27

APPENDIX A: Table of Approved Automation Machine Settings

Setting	Unit	Lower Limit	Upper Limit
Fill Time	Seconds	10	30
UV Cure Time	Seconds	110	130
Drain Time	Seconds	10	30
Rinse Time	Seconds	5	20
Final Drain	Seconds	10	30
Number of Drain / Rinse Cycles	N/A	2	3

Pathway LLC

Printed by Aaron Rogers

Printed on March 01, 2022

Local time zone (GMT-08:00) Pacific Time (America/Los_Angeles)

Page Arena > Changes > Summary

All information contained in this document is proprietary and confidential.

ECO-000899

Initial Release of SILQ components

Effective as of 01/28/2022 01:39:37 PM

Basic Information

Category	Engineering Change Order
Change Number	ECO-000899
Title	Initial Release of SILQ components
Description	These are SILQ (SLQ-) components for 16Fr 10ml catheter.
Associated Files	1 Reference File(s) + 0 Implementation File(s)
Routings	Drawing/Specification, Labeling / Promotional Material
Approval Deadline	None specified
Effectivity	01/28/2022
Expiration Date	N/A (This is a permanent Change.)
Lifecycle Status	Effective
Lifecycle Status Date	01/28/2022 01:39:37 PM
Creator	Aaron Rogers
Created on	01/13/2022 05:23:44 PM
Submitter	Aaron Rogers
Submitted on	01/28/2022 01:37:06 PM

Regulatory Assessment

QMS Impacted?	No
Medical Device Impacted?	Yes
Regulatory Impact	Project#00049

Details

Additional Attributes

Material affected? No

Risk Assessment N/A



0005 **R_x Only**

2-Way 100% Silicone ClearTract Catheter

Balloon Inflation Volume: 16 Fr 5 ml – 10 ml

Please read all instructions before using this device

Not made with natural rubber latex

Foley Catheter is MR Safe.

Contraindications

This device is contraindicated in the presence of conditions which create unacceptable risk during catheterization.

Intended Use

Urological catheter is intended for drainage of the urinary tract. Catheterization is accomplished by inserting the catheter through the urethra and into the bladder. Intended population is adults and pediatrics.

Directions

Installation & Inflation

1. Prepare application site in a sterile fashion.
2. Visually inspect the catheter for any signs of damage or deterioration prior to use.
3. If a trial inflation is conducted, follow steps 7 & 8 on the sterile field then evacuate the sterile water with the syringe. Ensure the balloon is empty prior to insertion.
4. Apply medical grade lubricant to distal end of the catheter.
5. Gently advance the catheter through the urethra until distal tip is beyond the bladder neck.
6. Confirm the position of the distal tip. The return of urine typically indicates the distal tip is in the bladder.
7. Fill Luer-tipped syringe with sterile water and connect to the inflation valve.
8. Fill the catheter balloon to the labeled volume. Do NOT exceed labeled volume.
9. Slowly withdraw the catheter until resistance is felt. This is a sign the balloon is resting at the bladder neck.

Deflation & Removal

1. Gently connect Luer-tipped syringe to the inflation valve. Gently aspirate while allowing the pressure in the balloon to assist filling the syringe with water occupying the balloon. Remove and reconnect the syringe in the inflation valve if there is no or slow deflation. Rapid deflation could collapse the lumen and prevent adequate deflation of the balloon.
2. If during the removal process the balloon ruptures, ensure all fragments are carefully removed from patient.
3. Discard after use.



Precautions

0016
0017

1. Sterile technique shall be applied during installation and removal of the catheter.
2. It is recommended to be replaced after 7 days of use, however, shall not exceed 30 days.

Cautions

0005

1. Federal law restricts the device to sale or by the order of a physician.
2. Never clamp the catheter. Use catheter plug or valve if required.

Warnings

1. Single Use ONLY. Do NOT reuse. Reuse may cause infection, illness, injury or death.
2. Read the instructions carefully before using the device.
3. Do not use this device beyond the labeled expiration date.
4. Do not use the device if the packaging is punctured or damaged.
5. Do not force catheter as injury to the urethra can occur.
6. Do not use needle.
7. Do not exceed recommended capacities.

0021
0025

Potential Adverse Events

With all Foley Catheters, the following potential adverse events could occur:

1. Urinary tract infection
2. Perforation or laceration of the urethra
3. Perforation or laceration of the bladder
4. Abrasions

0004
0007

Symbol Description

	Catalog Number
	Lot Number
	Use by
	Consult instructions for use
	Single Use
	Sterilized by Ethylene Oxide
	Do not resterilize
	Do not use if package is damaged
	Prescription Use Only (US)
	Foley Catheter is MR Safe

Manufactured for
SILQ Technologies Corp.
12100 Wilshire Blvd, # 800
Los Angeles, CA 90025
(949) 735-2614

0601
0002



2-WAY 100% SILICONE CLEARTRACT CATHETER

SIZE	16 Fr, 5 - 10 mL
QUANTITY	10

0025

Rx Only



STERILIZED



Symbols are defined in Instruction for Use (IFU)

Store away from direct light exposure at room temperature

Not made with natural rubber latex

Manufactured for:

SILQ Technologies Corp.
12100 Wilshire Blvd., Suite 800 | Los Angeles, CA 90025 | 949-735-2814

0010



(01) 11223344556675

(10) 165010

(17) YYMMDD

0011
0012

HDX-4002
Rev.B

0001
0002



2-WAY 100% SILICONE CLEARTRACT CATHETER

SIZE	16 Fr, 5 - 10 mL
QUANTITY	1

Rx Only

STERILE



Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

LOT	HDX1234-567890
REF	160510
	20YY-MM-DD



(01) 11223344556675
(10) 165010
(17) YYMMDD

Manufactured for:
SILQ Technologies Corp.
12100 Wilshire Blvd, Suite 800 | Los Angeles, CA 90025 | 949-735-2614

HDX-4002
Rev.B



2-WAY 100% SILICONE CLEARTRACT CATHETER

REF 160510

LOT HDX-MMDDYYYY

 20YY-MM-DD

QTY 100

 **SILQ Technologies Corp.** 12100 Wilshire Blvd, Suite 800 | Los Angeles, CA 90025 | 949-735-2614 HDX-4003 Rev. A

	Label Content Design Requirements	Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: C
		Page 2 of 7

- 5.1 It is the responsibility of the individual that designs product labels and labeling to comply with this procedure.

6.0 General Requirements

- 6.1 Paper stock and printing materials shall be selected based on its ability to maintain legibility throughout the labeling life (e.g., processing, storage, handling, and distribution).
- 6.2 Adhesive, if applicable, shall be selected based on its ability to maintain adherence throughout the labeling life (e.g., processing, storage, handling, and distribution).
- 6.3 Additional labeling requirements maybe required based on product specific regulations or guidance documents.

7.0 Labeling Requirements

#	Reference	Requirement Description	Device / Label / IFU	Markets
0001	21CFR801.1 21CFR812.5	Device label shall include the name of manufacturer, packer, or distributor. <i>Abbreviations for "Incorporated", etc. may be used. Qualifying phrase to distinguish relationship "Manufactured for ____", etc. may be used.</i>	Label	US
0002	21CFR801.1 21CFR812.5	Device label shall include street address, city, state, and zip code of manufacturer, packer, or distributor. <i>Street address may be omitted if listed in current directory.</i>	Label	US
0003	21CFR801.15	The predominant language shall be English.	Label / IFU	US
0004	21CFR801.15	Use of symbols shall include either adjacent text in all language(s) or a glossary of symbols.	Label / IFU	US
0005	21CFR801.15 21CFR801.109	The symbol statement "Rx only" or "[rx] only" or the statement "Caution: Federal law restricts this device to sale by or on the order of a ____" (i.e. "physician", etc.) <i>Except for surgical instruments.</i> "Only" needs to immediately follow "Rx." However, the symbol statement "Rx only" does not necessarily need to be bracketed in quotation marks, and the word "only" may appear in upper or lower case letters, for example, Rx only, Rx Only, or Rx ONLY.	Label / IFU	US
0006	21CFR801.109	The label of the device includes the method of its application or use. <i>Except for surgical instruments.</i>	Label / IFU	US
0007	21CFR801.15	If a glossary of symbols is used, it shall contain: - Each symbol used in the labeling - Title and designation number of standard(s) - Title of symbol and reference number (if any) - The symbols meaning	IFU	US
0008	21CFR801.15	If any representation appears in a foreign language, then all labeling requirements shall also appear in that foreign language.	Label / IFU	US
0009	21CFR801.18	Dates shall be expressed in the format YYYY-MM-DD (i.e. 4-digit year, 2-digit month, 2-digit day, each separated by a hyphen).	Label / IFU	US
0010	21CFR801.20 MDR Annex I 23.2	Label shall include a unique device identifier (UDI).	Label	US



Label Content Design Requirements

Doc No.: WI-0009

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 3 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0011	21CFR801.40	UDI must be presented in two forms: - Plain-text - Automatic Identification and Data Capture (AIDC)	Label	US
0012	21CFR801.40	UDI includes the following information: - Device identifier, which is a fixed portion that identifies the specific version or model of a device and the labeler of the device - Production identifier, which is a variable portion that identifies one or more of the following: lot number, serial number, expiration date, manufacturing date	Label	US
0013	21CFR801.45	Device meant to be reprocessed must be directly marked with UDI in either or both: - Plain-text - Automatic Identification and Data Capture (AIDC) or alternative technology <i>Except when it interferes with safety and efficacy, or not technically feasible.</i>	Device	US
0014	21CFR801.5	Instruction for Use includes statements of all conditions, purposes or uses for which such device is intended, prescribed or recommended.	IFU	US
0015	21CFR801.5	Instruction for Use includes Quantity of dose for: - each intended use - persons of different ages and physical conditions	IFU	US
0016	21CFR801.5 21CFR801.109	Instruction for Use includes frequency of administration or application.	IFU	US
0017	21CFR801.5 21CFR801.109	Instruction for Use includes duration of administration or application.	IFU	US
0018	21CFR801.5 21CFR801.109	Instruction for Use includes time of administration or application in relation to meals, onset of symptoms or other time factors.	IFU	US
0019	21CFR801.5 21CFR801.109	Instruction for Use includes route or method of administration or application.	IFU	US
0020	21CFR801.5	Instruction for Use includes preparation for use.	IFU	US
0021	21CFR801.109	Instruction for Use includes indications for Use, Effects, Hazards, Contraindications, side effects, precautions.	IFU	US
0022	21CFR801.109	Instruction for Use includes date of issuance or latest revision.	IFU	US
0023	21CFR812.5	An Investigational Device Exemption (IDE) label includes the quantity of contents.	Label	US
0024	21CFR812.5	An IDE label includes "CAUTION--Investigational device. Limited by Federal (or United States) law to investigational use."	Label	US
0025	21CFR812.5	The label or other labeling shall describe all relevant contraindications, hazards, adverse effects, interfering substances or devices, warnings, and precautions.	Label / IFU	US
0026	21CFR812.5	The labeling of an investigational device shall not bear any statement that is false or misleading in any particular and shall not represent that the device is safe or effective for the purposes for which it is being investigated.	Label / IFU	US



Label Content Design Requirements

Doc No.: WI-0009

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 4 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0027	21CFR812.5	An investigational device shipped solely for research on or with laboratory animals shall bear on its label the following statement: "CAUTION--Device for investigational use in laboratory animals or other tests that do not involve human subjects."	Label	US
0028	MDR Annex I 23.1	Information provided (medium, format, content, legibility, location) is appropriate for the device and healthcare practitioner.	Label / IFU	EU
0029	MDR Annex I 23.1	Where practical, information shall appear on the device itself.	Device	EU
0030	MDR Annex I 23.1	Information on device labeling shall take the form of symbols, where appropriate. In areas for which no standards exist, the symbols and colors shall be described in the documentation supplied with the device.	Label / IFU	EU
0031	MDR Annex I 23.1	Labels shall be provided in human-readable format and supplemented with machine-readable information.	Label	EU
0032	MDR Annex I 23.1	Each device shall be accompanied by Instructions for Use (IFU). <i>Except for class I and IIa devices that don't require information for safe use.</i>	IFU	EU
0033	MDR Annex I 23.1, 23.2, 23.3	Device name shall appear on the device label, sterile barrier packaging and IFU.	Label / IFU	EU
0034	MDR Annex I 23.1, 23.2, 23.3	Details needed to identify the device / device description shall appear on the device label, sterile barrier packaging.	Label	EU
0035	MDR Annex I 23.1, 23.2	Packaging contents shall appear on the device label.	Label	EU
0036	MDR Annex I 23.1, 23.2, 23.4	Intended for use statement shall appear on the device label.	Label	EU
0037	MDR Annex I 23.1, 23.2, 2.3, 23.4	Name and address of the manufacturer shall appear on the device label, sterile barrier packaging, IFU.	Label / IFU	EU
0038	MDR Annex I 23.2	Name and Address of authorized representative shall appear on the device label. <i>Except where the manufacturer resides within the EU.</i>	Label	EU
0039	MDR Annex I 23.2, 23.4	Where applicable, device label and IFU shall include an indication device contains: - Medicinal substance, human blood, plasma derivative - Tissues, cells or derivatives of human origin - Tissues, cells or derivatives of animal origin	Label / IFU	EU
0040	MDR Annex I 23.2, 23.4	Where applicable, device label and IFU shall include an indication device contains: - Substances which are carcinogenic, mutagenic or toxic to reproduction - Substances having endocrine-disrupting properties	Label / IFU	EU
0041	MDR Annex I 23.2	The label shall include the lot or serial number of the device along with the appropriate symbol.	Label	EU
0042	MDR Annex I 23.2, 23.3	The sterile barrier packaging shall include the date of manufacture (YYYY-MM).	Label	EU
0043	MDR Annex I 23.2, 23.3	The label and sterile barrier packaging shall include the expiration date (YYYY-MM), if applicable. Otherwise the label shall include the date of manufacture.	Label	EU
0044	MDR Annex I 23.2	The label shall include the UDI number.	Label	EU

N/A

N/A



Label Content Design Requirements

Doc No.: WI-0009

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 5 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0045	MDR Annex I 23.2, 23.4	The label and IFU shall indicate any storage and / or handling conditions that apply.	Label / IFU	EU
0046	MDR Annex I 23.3	The sterile barrier packaging shall be indicated it is the sterile barrier.	Label	EU
0047	MDR Annex I 23.2, 23.3, 23.4	The label, sterile barrier packaging, and IFU shall indicate its method and state of sterility that applies, with the appropriate symbol.	Label / IFU	EU
0048	MDR Annex I 23.2	The label shall include any immediate warnings or precautions that apply.	Label	EU
0049	MDR Annex I 23.2, 23.4	The label and IFU shall indicate the device is for single use with the appropriate symbol, if that applies.	Label / IFU	EU
0050	MDR Annex I 23.2	The label shall indicate the device is a reprocessed single-use device, and the number of reprocessing cycles performed any limitations that apply.	Label	EU
0051	MDR Annex I 23.2	The label and sterile barrier packaging shall indicate the device is "custom-made" if applicable.	Label	EU
0052	MDR Annex I 23.2	The label shall indicate the device is a medical device.	Label	EU
0053	MDR Annex I 23.2, 23.3	The label and sterile barrier packaging shall indicate the device is "exclusively for clinical investigation" if applicable.	Label	EU
0054	MDR Annex I 23.2, 23.4	The label and IFU shall indicate the device is composed of / combination of substances intended to be introduced / applied.	Label / IFU	EU
0055	MDR Annex I 23.2	The label shall indicate the serial number for implantable devices, if applicable.	Label	EU
0056	MDR Annex I 23.3	The sterile barrier packaging shall include an instruction to check the IFU if the sterile packaging is damaged or unintentionally opened before use.	Label	EU
0057	MDR Annex I 23.4	Instructions for Use includes device's intended purpose with a clear specification of indications, contra-indications, the patient target group or groups, and of the intended users, as appropriate.	IFU	EU
0058	MDR Annex I 23.4	Instructions for Use includes where applicable, a specification of the clinical benefits to be expected.	IFU	EU
0059	MDR Annex I 23.4	Instructions for Use includes where applicable, links to the summary of safety and clinical performance.	IFU	EU
0060	MDR Annex I 23.1, 23.4	Instructions for Use includes the performance characteristics of the device.	IFU	EU
0061	MDR Annex I 23.4	Instructions for Use includes where applicable, information allowing the healthcare professional to verify if the device is suitable and select the corresponding software and accessories.	IFU	EU
0062	MDR Annex I 23.4	Instructions for Use includes where applicable, information allowing the healthcare professional to verify if the device is suitable and select the corresponding software and accessories.	IFU	EU
0063	MDR Annex I 23.1, 23.4	Instructions for Use includes any residual risks, limitations, contra-indications, precautions, warnings and any undesirable side-effects, including information to be conveyed to the patient in this regard.	IFU	EU

N/A

N/A



Label Content Design Requirements

Doc No.: WI-0009

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 6 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0064	MDR Annex I 23.4	Instructions for Use includes specifications the user requires to use the device appropriately (e.g. degree of accuracy).	IFU	EU
0065	MDR Annex I 23.4	Instructions for Use includes details of any preparatory treatment or handling of the device before it is ready for use or during its use, such as sterilisation, final assembly, calibration, etc., including the levels of disinfection required to ensure patient safety and all available methods for achieving those levels of disinfection.	IFU	EU
0066	MDR Annex I 23.4	Instructions for Use includes any requirements for special facilities, or special training, or particular qualifications of the device user and/or other persons.	IFU	EU
0067	MDR Annex I 23.4	Instructions for Use includes the information needed to verify whether the device is properly installed and is ready to perform safely and as intended.	IFU	EU
0068	MDR Annex I 23.4	Instructions for Use includes instructions in the event of the sterile packaging being damaged or unintentionally opened before use, if supplied sterile.	IFU	EU
0069	MDR Annex I 23.4	Instructions for Use includes appropriate instructions for sterilization, if applicable.	IFU	EU
0070	MDR Annex I 23.4	Instructions for Use includes if the device is reusable, information on the appropriate processes for allowing reuse, including cleaning, disinfection, packaging and, where appropriate, the validated method of re-sterilization.	IFU	EU
0071	MDR Annex I 23.4	Instructions for Use includes an indication, if appropriate, that a device can be reused only if it is reconditioned under the responsibility of the manufacturer to comply with the general safety and performance requirements.	IFU	EU
0072	MDR Annex I 23.4	Instructions for Use includes information on known characteristics and technical factors known to the manufacturer that could pose a risk if the device were to be re-used. Applied to single use devices.	IFU	EU
0073	MDR Annex I 23.4	Instructions for Use includes information to identify other devices and equipment in order to obtain a safe combination, and / or any known restrictions. Applies to devices intended for use together with other devices and/or general-purpose equipment.	IFU	EU
0074	MDR Annex I 23.4	Instructions for Use includes detailed information as to the nature, type and where appropriate, the intensity and distribution of the emitted radiation and the means of protecting the patient, user, or other person from unintended radiation during use of the device. Applies to devices that emit radiation for medical purposes.	IFU	EU
0075	MDR Annex I 23.4	Instructions for Use includes information that allows the user and/or patient to be informed of any warnings, precautions, contra- indications, measures to be taken and limitations of use regarding the device.	IFU	EU
0076	MDR Annex I 23.4	Instructions for Use includes the overall qualitative and quantitative information on the materials and substances to which patients can be exposed with implantable devices.	IFU	EU
0077	MDR Annex I 23.4	Instructions for Use includes warnings or precautions to be taken in order to facilitate the safe disposal of the device, its accessories and the consumables used with it, if any.	IFU	EU

N/A

N/A

	Label Content Design Requirements		Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 7 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0078	MDR Annex I 23.4	Instructions for Use includes the circumstances in which the user should consult a healthcare professional, if intended for use by lay persons.	IFU	EU
0079	MDR Annex I 23.4	Instructions for Use includes date of issue of the instructions for use or, if they have been revised, date of issue and identifier of the latest revision of the instructions for use.	IFU	EU
0080	MDR Annex I 23.4	Instructions for Use includes notice to the user and/or patient that any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.	IFU	EU
0081	MDR Annex I 23.4	Instructions for Use includes information to be supplied to the patient with an implanted device.	IFU	EU
0082	MDR Annex I 23.4	Instructions for Use includes, for devices that incorporate electronic programmable systems, including software, or software that are devices in themselves, minimum requirements concerning hardware, IT networks characteristics and IT security measures, including protection against unauthorized access, necessary to run the software as intended.	IFU	EU
0083	MDR Article 20#5, Annex V	CE conformity marking shall consist of the initials "CE" followed by the identification number of the notified body responsible for the device's conformity assessment. The various components of the CE marking must have substantially the same vertical dimension, which may not be less than 5 mm.	Label / IFU	EU
0084	FDA (1768)	If a manufacturer elects to include a statement in medical product labeling indicating that natural rubber latex or synthetic derivatives of natural rubber latex were not used as materials in the manufacture of their medical product and container, FDA recommends the use of the statement "Not made with natural rubber latex." If this statement is made without any qualification, it would apply to the medical product, its container, and any packaging.	Label / IFU	US
0085	ASTM F2503	For devices which are brought into the magnetic resonance (MR) environment, permanently mark the item as MR Safe, MR Conditional, or MR Unsafe using the icons with color or black and white options. Refer to ASTM F2503 TABLE 1 and TABLE 2 for details.	Device / Label / IFU	US

8.0 Records

8.1 Label Content documents will be processed per QP-002. Refer to QP-003 for storage location and retention requirements.

9.0 Flowchart – N/A

10.0 Change History

Revision	Change Description
A	Initial release. This procedure will replace procedure QP-7.09
B	Update header, references, and format of requirements (table).
C	Remove reference to FDA 1150 Guidance

	Label Content Design Requirements		Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 1 of 7

1.0 Purpose

- 1.1 The purpose of this document is to indicate the requirements for labels and labeling associated with medical devices.

2.0 Scope

- 2.1 This applies to labeling that is considered the device labels (primary and secondary packages) and the Instructions for Use (IFU)
- 2.2 The requirements indicated in this document are associated with regulatory requirements associated with FDA and EU.
- 2.3 The procedure does not apply to promotions materials.

3.0 References

- 3.1 21 CFR 801 – Labeling
- 3.2 21 CFR 812 – Investigation Device Exemption
- 3.3 ASTM 2503-20 – Standard Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment
- 3.4 ISO 15223-1 – Medical devices – Symbols to be used with medical device labels, labeling, and information to be supplied – Part 1: General requirements
- 3.5 FDA (1768) Guidance – Recommendations for Labeling Medical Products to Inform Users that the Product or Product Container is not Made with Natural Rubber Latex, December 2, 2014
- 3.6 QP-002 – Document Control
- 3.7 QP-003 – Quality Records

4.0 Definitions

- 4.1 Label – Display of written, printed, or graphic matter upon the immediate container.
- 4.2 Labeling – All labels and other written, printed, or graphic matter upon any article or any of its containers or wrapper, or accompanying such article at any time while a device is held for sale after shipment or delivery for shipment in interstate commerce.
- 4.3 Investigational device – A device that is the object of a clinical investigation or research involving one or more human subjects to determine its safety and effectiveness.
- 4.4 Unique Device Identifier (UDI) – An identifier that adequately identifies a device through its distribution.

5.0 Responsibilities

	Risk Management Plan – SILQ Foley Catheter		Doc No.: RM-0018
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 1 of 2

1.0 Purpose

- 1.1 The purpose of this document is to indicate the requirements of the risk management plan that is associated with the development of the SILQ ClearTract Foley Catheter.

Note: The product name SILQ ClearTract Catheter and Hydrophilix Foley Catheter is used interchangeably. The company name Hydrophilix and SILQ are used interchangeably as the name recently was changed from Hydrophilix to SILQ.

2.0 Scope

- 2.1 The document applies for the design projects associated with the product introduction of the SILQ ClearTract Foley Catheter family of products.
- 2.2 This document will identify the activities planned for risk management.

3.0 References

- 3.1 ISO 14971:2019 – Medical devices – Application of risk management to medical devices
- 3.2 Project#00017 – Hydrophilix Foley Catheter product introduction (US)
- 3.3 Project#00047 – SILQ ClearTract Foley Catheter Product Line & Shelf Life Extension (US)
- 3.4 Project#00049 – SILQ ClearTract Foley Catheter Product Line Extension (US)
- 3.5 QP-002 – Document Control
- 3.6 QP-004 – Risk Management
- 3.7 QP-007 – Design Control
- 3.8 QP-011 – Post Market Surveillance

4.0 Responsibilities & Authorities

- 4.1 The risk management team will include members of Pathway. Team members will include individuals representing Research & Development (R&D), Manufacturing, Marketing, Regulatory Affairs (RA), and Quality Assurance (QA).
- 4.2 Individuals that participate in risk management activities will have knowledge of and experience with the device (or similar device) and how it is manufactured.
- 4.3 The risk management team will be required to complete all activities indicated in this document. Team members from Pathway will be required to review all content of risk management deliverables and agree on its conclusions with documented approval.

	Risk Management Plan – SILQ Foley Catheter		Doc No.: RM-0018
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 2 of 2

5.0 Product Description

- 5.1 The SILQ ClearTract Foley Catheter is intended for use to drain the urinary tract.

6.0 Product Life-Cycle Phases

6.1 Design

- 6.1.1 In the appropriate phase of the design control process, per QP-007, a Hazard Analysis (HA) and Design Failure Mode Effects Analysis (DFMEA) will be developed per, QP-004, as input to the design of the device.
- 6.1.2 Criteria of risk acceptability is identified in QP-004.
- 6.1.3 If additional risk(s) are identified during the progression through the design control phases, especially in verification and validations, risk management deliverables will be reviewed and updates appropriately.

6.2 Process Validation

- 6.2.1 In the "Transfer to Manufacturing" phase of the design control process per QP-007, a Process Failure Mode Effects Analysis (PFMEA) will be developed as a deliverable process validation. Criteria of risk acceptability is identified in QP-004.

6.3 Production & Post-Production

- 6.3.1 Information collected from production and post-production activities will be evaluated for risk as part of the post-market surveillance process per QP-011. Risk management report will be review and newly identified risk or changes to existing risk will be incorporated.
- 6.3.2 Revised information to the risk management reports will be reviewed with R&D, RA, and QA to determine if action shall be taken to the product design to mitigate identified risk. Action shall be incorporated in a change control project per QP-007.

7.0 Verification Activities

- 7.1 Activities associated for risk controls (initial & mitigation) will be documented in the risk management reports.

8.0 Risk Management File

- 8.1 The risk management file consists of this plan and the risk management report (i.e., Hazard Analysis, DFMEA, & PFMEA). The risk management plan will be maintained as a controlled document per QP-002 and associated in the product's Design History File (DHF).

9.0 Change History

Revision	Change Description
A	Initial release
B	Update name change "Hydrophilix" to "SILQ". Add ref to Project#00047
C	Add references to Project#00049

	Product Specification / Traceability Matrix – SILQ Foley Products		Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: D
			Page 1 of 7

1.0 Purpose

- 1.1 To document the design requirements for the SILQ ClearTract Foley Catheter products and provide traceability between the risk management reports (i.e, Hazard Analysis, Design Failure Mode Effects Analysis) and Verification & Validation deliverables.

2.0 Scope

- 2.1 SILQ Foley Catheter family of products

3.0 General Information

Table 1: General Information	
Attribute	Description
Device Name	SILQ 100% Silicone 2-Way Foley Catheter, 14Fr/10ml SILQ 100% Silicone 2-Way Foley Catheter, 16Fr/5ml SILQ 2-Way 100% Silicone ClearTract Catheter, 16Fr/5-10ml SILQ 2-Way 100% Silicone ClearTract Catheter, 16Fr/10ml SILQ 100% Silicone 2-Way Foley Catheter, 18Fr/10ml
Device's Intended Use / Indications for Use Statement	Urological catheter is intended for drainage of the urinary tract. Catheterization is accomplished by inserting the catheter through the urethra and into the bladder. Intended population is adults and pediatrics.
User Types	Healthcare practitioners / Clinical user (e.g. nurses, doctors)
User Environment	Hospital
Patient Population	Any patient population
Region or Country	US
Agency	FDA
Device Class	II
Regulation - Description	21 CFR 876.5130 - Urological catheter and accessories
Product Code	EZL - Catheter, Retention Type, Balloon

4.0 Requirements and Trace Matrix

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: D
								Page 2 of 7

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Biological	N/A	The device shall meet biocompatibility requirement for a surface device with mucosal membrane contact for a prolong exposure >24 to 30 days, which the following evaluations are required: 1. Cytotoxicity 2. Sensitization 3. Irritation or Intracutaneous Reactivity 4. Material Mediated Pyrogen	Device	1	1, 7, 9	HA#8.1, 8.2, 8.3, 9.1	TR-0071	N/A
Labeling	N/A	Information shall be presented in symbol if the symbol is recognized by the regulatory agency(s)	Labeling	2	4, 14	HA#12.1, 12.2, 13.1 DFMEA#6, 8, 10	N/A	SLQ-4000 SLQ-4001 SLQ-4006
Labeling	N/A	Product labeling (including direction for use) shall comply to FDA regulations on labeling. Specific requirement for devices with antimicrobial agents is indicated in guidance document (14).	Labeling	3	8, 13	N/A	N/A	SLQ-4000 SLQ-4001 SLQ-4006
Labeling	N/A	Product labeling shall include the following information associated with the device description: 1. Outer Diameter (French Size) 2. Balloon Size (cc) 3. Effective Shaft Size 4. Type of Coating (if applicable) 5. Amount/ concentration of antimicrobial agent (if applicable)	Labeling	4	7, 13	HA#11.3 DFMEA#8	N/A	SLQ-4000 SLQ-4001 SLQ-4006
Material	N/A	Connector shall be made of rigid or semi-rigid materials.	Connector	5	12	HA#10.1, 10.2 DFMEA#4	N/A	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6004
Material	N/A	The device shall be constructed on medical grade silicone material.	Device	6	N/A	HA#8.1, 8.2, 8.3, 8.4 DFMEA#1, 2, 3	N/A	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6004

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: D
								Page 3 of 7

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Packaging	N/A	The device will be packaged in a tyvek/poly pouch., which 10 pouched devices will be packaged in a box.	Packaged Device	7	N/A	HA#7.1, 14.1 DFMEA#5	N/A	HDX-8000 HDX-8001 HDX-8002 HDX-8006 SLQ-8004
Packaging	N/A	The peel strength for 1" MFG sealed section of the sterile barrier packaging shall be ≥ 1 lb.	Sterile Barrier	8	3	HA#7.1, 14.1 DFMEA#5	TR-0074	N/A
Packaging	N/A	The sterile barrier packaging (include sealed area) shall not exhibit the following visual defects: 1. Unsealed Areas 2. Nonhomogeneous or Undersealed Areas 3. Over-sealed Areas 4. Narrow Seals 5. Channels 6. Wrinkles/Fold-overs/Cracks 7. Tears/Pinholes	Sterile Barrier	9	3	HA#7.1, 14.1 DFMEA#5	TR-0074	N/A
Packaging	N/A	The sterile barrier packaging shall not exhibit a constant stream of bubbles when submerge in water and internal pressurized to a determine test pressure. Balloon Integrity (Resistance to Rupture) - The balloon shall be inflated easily with distilled or deionized water to labeled volume without showing any evidence of breakage throughout a 24 hour period.	Sterile Barrier Device	10 11	3 7	HA#7.1, 14.1 DFMEA#5	TR-0074 TR-0072	N/A N/A
Performance	N/A							
Performance	N/A	Balloon Volume Maintenance - The balloon shall maintain its volume throughout a 24 hour period.	Device	12	7	HA#10.1 DFMEA#3	TR-0072 TR-HDX-0001	N/A

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: D
								Page 4 of 7

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Performance	N/A	Connector shall be non-interconnectables, which does not appear to provide a secure connection when forcefully assembled and shall easily disengage.	Connector	13	12	HA#10.1 DFMEA#4	TR-0072 TR-HDX-0001	N/A
Performance	N/A	Deflation Reliability (Failure to Deflate) - The balloon shall deflate to within 4 French sizes of the label French size within 15 min or be otherwise manipulated to effect drainage within this time period.	Device	14	7	HA#10.1, 10.4 DFMEA#4	TR-0072 TR-HDX-0001	N/A
Performance	N/A	Flow Rate though Drainage Lumen - Shall have a minimum average flow rate of 100 cm3/min.	Device	15	7	HA#10.2 DFMEA#1	TR-0072 TR-HDX-0001	N/A
Performance	N/A	Inflated Balloon Response to Traction - The entire, inflated balloon of catheters shall not pass into or through the funnel barrel neck (28 Fr x 6 cm) when a 1 lb weight is affixed from the distal catheter end and subjected to static (2 minutes) and dynamic load (2 ft drop).	Device	16	7	HA#10.1 DFMEA#3	TR-0072 TR-HDX-0001	N/A
Performance	N/A	The device shall be easily inserted through the urethra until the deflated balloon passes through the urethral sphincter and in a manner that prevents urethral trauma.	Device	17	N/A	HA#10.3 DFMEA#1	TR-0072 TR-HDX-0001	N/A
Physical	N/A	Balloon & Shaft Size - The proximal tip, the balloon, and the shaft (for 10 in. distal to the balloon) shall meet the requirements on size and tolerances as is relates to the labeled French size per the tolerances indicated below: Tip = ± 1 Fr Shaft = ± 1 Fr Balloon = + 4 Fr (Uninflated and Deflated after Immersion)	Device	18	7	HA#10.3 DFMEA#1	TR-0072 TR-HDX-0001	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6004

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: D
								Page 5 of 7

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Physical	N/A	The balloon size shall be 5ml, 10ml or 5-10ml.	Device	19	N/A	HA#10.3 DFMEA#1	TR-0072 TR-HDX-0001	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6004
Physical	N/A	The device shall be able to connect into an industry standard or equivalent urine bag to facilitate sterile and easy urine collection.	Device	21	N/A	HA#10.1 DFMEA#4	TR-0072 TR-HDX-0001	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6004
Physical	N/A	The device shall be equipped with 2 lumens, which have the following purposes: 1. Drainage 2. Balloon Inflation	Device	22	N/A	HA#10.1, 10.2 DFMEA#1, 2	N/A	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6004
Physical	N/A	The product offering will be French Sizes (Outer Diameter) of 14 Fr (0.184"), 16 Fr (0.210"), and 18 Fr (0.236").	Device	23	N/A	HA#10.1 DFMEA#1	See Req#18	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6004
Stability	N/A	The device shall have a final shelf-life period of 3 years, but initially be qualified for shorter shelf-life periods (e.g., 1 year).	Sterile, Packaged Device	24	N/A	HA#7.1, 8.4 DFMEA#1, 2, 3, 4, 5	TR-0073	N/A

	Product Specification / Traceability Matrix – SILQ Foley Products						Doc No.: PS-0006
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.						Rev.: D
							Page 6 of 7

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Sterilization	N/A	Device shall be sterile to a sterility assurance level (SAL) of 10 ⁻⁶ when terminally sterilized via Ethylene Oxide (EO).	Device	25	2	HA#7.1	TR-0075 TR-HDX-0009	N/A
Performance	N/A	Packaged devices shall maintain its sterility and functionality through transportation	Sterile, Packaged Device	26	5, 6	HA#14.1 DFMEA#5	TR-0074	N/A
User	N/A	A healthcare provider must be able to easily and quickly implement device.	Device	27	N/A	HA#13.1	N/A	SLQ-4000 SLQ-4006

5.0 Reference Documents

Table 3: Reference Documents				
Source No	Organization	Document No	Edition	Title
1	ANSI/ AAMI/ ISO	10993-1	2009	Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management process
2	ANSI/ AAMI/ ISO	11135	2014	Sterilization of health care products – Ethylene Oxide – Requirements for development, validation and routine control of a sterilization process for medical devices
3	ISO	11607-1	2019	Packaging for terminally sterilized medical devices – Part 1: Requirements for materials, sterile barrier systems and packaging systems
4	ANSI/ AAMI/ ISO	15223-1	2016	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied
5	ASTM	D4169	2016	Standard Practice For Performance Testing Of Shipping Containers And Systems

	Product Specification / Traceability Matrix – SILQ Foley Products		Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: D
			Page 7 of 7

Table 3: Reference Documents				
Source No	Organization	Document No	Edition	Title
6	ASTM	D4332	2014	Standard Practice for Conditioning Containers, Packages, or Packaging Components for Testing.
7	ASTM	F623	2019	Standard Performance Specification for Foley Catheter
8	FDA	21 CFR 801	4/15/2018	Title 21 - Food and Drugs Chapter I – Food and Drug Administration Department of Health Services; Subchapter H – Medical Devices; Part 801 Labeling
9	FDA	FDA-2013-D-0350	6/16/2016	FDA Guidance: Use of International Standard ISO 10993, "Biological Evaluation of Medical Devices Part 1: Evaluation and Testing within a risk management process"
11	ISO	14971	2019	Medical Devices - Application of Risk Management to Medical Devices
12	ISO	80369-1	2018	Small-Bore Connectors For Liquids And Gases In Healthcare Applications - Part 1: General Requirements.
13	FDA	N/A	N/A	Guidance for the Content of Premarket Notifications for Conventional and Antimicrobial Foley Catheters
14	FDA	81 FR 38911	6/15/2016	Use of Symbols in Labeling
15	N/A	RM-0019	C	Risk Management Report, SILQ ClearTract Foley Catheter Hazard Analysis (HA)
16	N/A	RM-0020	C	Risk Management Report, SILQ ClearTract Foley Catheter Design Failure Mode and Effects Analysis (DFMEA)

6.0 Change History

Revision	Change Description
A	Initial Release
B	Update format, indications for use
C	Update format; name change Hydrophilix to SILQ; add 16Fr 5-10ml; general updates (spelling); correction to #22 (2-Way); remove #20 (Device is 2-way, not a 3-way), remove ref doc#10
D	Update references to HA and DFMEA, component specs SLQ-6004, labeling SLQ-4000, SLQ-4001, SLQ-4006; add 16Fr/10ml

1.0 Purpose

- 1.1 The purpose of this document is to identify the associated risk with the intended use of the SILQ ClearTract Foley Catheter.

2.0 Scope

- 2.1 The product numbers associated with the SILQ ClearTract Foley Catheter is indicated in the following table:

Table 1: SILQ ClearTract Catheter Products	
Part No.	Description
HDX-1410	2-way 100% Silicone ClearTract Catheter, 14 Fr/10ml
HDX-1605	2-way 100% Silicone ClearTract Catheter, 16 Fr/5ml
HDX-1810	2-way 100% Silicone ClearTract Catheter, 18 Fr/10ml
HDX-160510	2-way 100% Silicone ClearTract Catheter, 16 Fr/5-10ml
SLQ-1610	2-Way 100% Silicone ClearTract Catheter, 16Fr/10ml

- 2.2 The intended use of this device is the following:

- 2.2.1 Urological catheter is intended for drainage of the urinary tract. Catheterization is accomplished by inserting the catheter through the urethra and into the bladder. Intended population is adults and pediatrics.

3.0 References

- 3.1 QP-004 – Risk Management

4.0 Risk Tables Definitions

- 4.1 # – Enumerate each Hazard Category with a sequential number
- 4.2 Hazard Category – Describes the general grouping of a hazard(s) / hazardous situation(s) that is being evaluated for risk
- 4.3 PH# – Enumerate each Potential Hazard with a sequential number, using the format {Hazard Category #}.{PH#}, (e.g. 1.1, 1.2, 1.3)
- 4.4 Potential Hazard (PH) – Source of harm
- 4.5 Potential Harm – Injury or damage to the health of people, or damage to property or the environment
- 4.6 Severity (S) – Measure of the possible consequences of a hazard.
- 4.7 Potential Failure Mode (FM) – The potential consequences of the failure
- 4.8 Risk Controls - Process through which decisions are reached and protective measures are implemented for reducing risk to, or maintaining risk within specified levels.

	Hazard Analysis (HA) – Risk Management Report, SILQ ClearTract Foley Catheter	Doc No.: RM-0019
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: C
		Page 2 of 9

4.8.1 Risk Controls¹ – Initial risk controls that are incorporated.

4.8.2 Risk Controls² – Secondary risk controls, also known as “Mitigation”.

4.9 Verification – Evidence source of linkage that initial (i.e., Risk Controls¹) risk controls have implemented, verified, or validated.

4.10 Probability – Likelihood of the hazard occurring

4.10.1 Probability¹ (P¹) – The probability assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls¹)

4.10.2 Probability² (P²) – The probability assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls²)

4.11 Risk Rating – Residual risk determination based on the concept of ‘As Far As Possible (AFAP)’ approach. The AFAP approach provides three (3) possible outcomes, which are 1) the residual risk exceeds the risk acceptability (RED), 2) the residual risk is acceptable to the risk acceptability criteria (GREEN), and 3) the residual risk is between the two (2) states indicated and additional controls are suggested to lower the risk to the lowest practical level, without economically considerations, based on management decision (YELLOW).

4.11.1 Risk Rating¹ (R¹) – This calculation is post the implementation of initial risk controls (i.e., Risk Controls¹)

4.11.2 Risk Rating² (R²) – This calculation is post the implementation of the secondary risk controls (i.e., Risk Controls²)

5.0 Identify Characteristics that could Impact on Safety

Table 2: : Questions from ISO-14971:2012 Annex C w/ Answers		
#	ISO-14971:2012 Annex C Question	Answer
1	What is the intended use and how is the medical device to be used?	<i>Refer to 2.2</i>
2	Is the medical device intended to be implanted?	<i>No</i>
3	Is the medical device intended to be in contact with the patient or other persons?	<i>Yes</i>
4	What materials or components are utilized in the medical device or are used with, or are in contact with, the medical device?	<i>Silicone</i>
5	Is energy delivered to or extracted from the patient?	<i>No</i>
6	Are substances delivered to or extracted from the patient?	<i>Yes, urine</i>
7	Are biological materials processed by the medical device for subsequent re-use, transfusion or transplantation?	<i>No</i>
8	Is the medical device supplied sterile or intended to be sterilized by the user, or are other microbiological controls applicable?	<i>Yes, supplied sterile by EO.</i>
9	Is the medical device intended to be routinely cleaned and disinfected by the user?	<i>No</i>
10	Is the medical device intended to modify the patient environment?	<i>No</i>
11	Are measurements taken?	<i>Yes, measurements of urine.</i>

Table 2: : Questions from ISO-14971:2012 Annex C w/ Answers

#	ISO-14971:2012 Annex C Question	Answer
12	Is the medical device interpretative?	No
13	Is the medical device intended for use in conjunction with other medical devices, medicines or other medical technologies?	The device has been designed to be used with a urine collection bag.
14	Are there unwanted outputs of energy or substances?	No
15	Is the medical device susceptible to environmental influences?	No
16	Does the medical device influence the environment?	No
17	Are there essential consumables or accessories associated with the medical device?	No
18	Is maintenance or calibration necessary?	No
19	Does the medical device contain software?	No
20	Does the medical device have a restricted shelf-life?	Yes, refer to #23
21	Are there any delayed or long-term use effects?	No
22	To what mechanical forces will the medical device be subjected?	The device is subjected to external tensile forces imposed by the urine content in the collection bag, Internal pressure exerted by urine within the lumen, external compressive forces of the bladder and urethra, internal forces from the filled balloon.
23	What determines the lifetime of the medical device?	Product and sterile barrier performance
24	Is the medical device intended for single use?	Yes
25	Is safe decommissioning or disposal of the medical device necessary?	Yes
26	Does installation or use of the medical device require special training or special skills?	No
27	How will information for safe use be provided?	Instructions for use
28	Will new manufacturing processes need to be established or introduced?	No
29	Is successful application of the medical device critically dependent on human factors such as the user interface?	No
29.1	Can the user interface design features contribute to use error?	No
29.2	Is the medical device used in an environment where distractions can cause use error?	No
29.3	Does the medical device have connecting parts or accessories?	Yes, urine collection bag, Luer-lock syringe
29.4	Does the medical device have a control interface?	No
29.5	Does the medical device display information?	No
29.6	Is the medical device controlled by a menu?	No
29.7	Will the medical device be used by persons with special needs?	No
29.8	Can the user interface be used to initiate user actions?	No
30	Does the medical device use an alarm system?	No

Table 2: : Questions from ISO-14971:2012 Annex C w/ Answers

#	ISO-14971:2012 Annex C Question	Answer
31	In what way(s) might the medical device be deliberately misused?	None
32	Does the medical device hold data critical to patient care?	No
33	Is the medical device intended to be mobile or portable?	Yes
34	Does the use of the medical device depend on essential performance?	No

6.0 Risk Acceptability Determination

- 6.1 The criteria of risk acceptability will be determined with a Risk Priority Number (RPN) and on the concept of 'As Far As Possible' (AFAP) approach. Refer to the Risk Rating definition in §4.0.
- 6.2 There are two (2) components used in determining risk acceptability, which are Severity and Probability. Refer to the following tables:

Table 3: Severity Ratings	
Index	Description
1	No patient effect
2	Minor injury that resolves without treatment
3	Moderate injury that may require treatment
4	Serious injury requiring medical intervention
5	Life threatening injury, disability, death

Table 4: Probability Ratings	
Index	Description
1	Improbable - <0.001%
2	Remote – 0.001% to .01%
3	Occasional – 0.01% to 0.1%
4	Likely – 0.1% to 1%
5	Frequent - >1%

6.3 Risk acceptability will be determined per the following risk index:

Table 5: Risk Rating Table						
Severity	5	5	10	15	20	25
	4	4	8	12	16	20
	3	3	6	9	12	15
	2	2	4	6	8	10
	1	1	2	3	4	5
		1	2	3	4	5
		Probability				

7.0 Risk Assessment Tables

Table 6: SILQ ClearTract Foley Catheter													
#	Hazard Category	PH#	Potential Hazard	Potential Harm	S	Potential Failure Mode	Risk Controls ¹	Verification	P ¹	R ¹	Risk Controls ²	P ²	R ²
1	Acoustic Energy		The medical device does not have any acoustic energy hazards.										
2	Electromagnetic Energy		The medical device does not have any electromagnetic energy hazards.										
3	Mechanical Energy		The medical device does not have any mechanical energy hazards.										
4	Potential (stored) Energy		The medical device does not have any potential energy hazards.										
5	Radiation Energy		The medical device does not have any radiation energy hazards.										
6	Thermal Energy		The medical device does not have any thermal energy hazards.										
7	Biological	7.1	Biological Contamination	Fever/ infection to patient	3	Device is not provided sterile/ insufficient sterilization process	Device will be terminally sterilized with validated method (Design)	Sterilization validation TR-0075, TR-HDX-0009	1	3	N/A	N/A	N/A
						Device is reused	Warnings statement for single-use (Information)	IFU SLQ-4006	2	3	N/A	N/A	N/A
						Inadequate sterile barrier system	Sterile barrier system will be designed for material that are commonly used for single-use medical device (Design)	Packaging Validation TR-0074	1	3	N/A	N/A	N/A
							Warning will be incorporated for inspecting packaging prior to use as well not using the device if packaging is damaged (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A



Hazard Analysis (HA) – Risk Management Report, SILQ ClearTract Foley Catheter

Doc No.: RM-0019

This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 6 of 9

Table 6: SILQ ClearTract Foley Catheter

#	Hazard Category	PH#	Potential Hazard	Potential Harm	S	Potential Failure Mode	Risk Controls ¹	Verification	P ¹	R ¹	Risk Controls ²	P ²	R ²
						Use beyond expiration date	Warning statement (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
		7.2	Biological Contamination	Fever/ infection to clinician	3	Proper handling and disposal are not conducted during clinical application	Warning statement (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
8	Chemical	8.1	Exposure to known allergens	Allergic reaction	3	Materials (e.g., latex, rubber) that is known to cause allergic response	Material will be selected that do not cause allergic response (Design) / Biocompatibility Testing	Design verification – biocompatibility TR-0071	1	3	N/A	N/A	N/A
		8.2	Exposure to toxic materials / False positives in the blood work evaluation / red blood cell damage	Toxicity	3	Materials leach during clinical application	Materials will be selected that are non-toxic (Design)	Biocompatibility Test Report TR-0071	1	3	N/A	N/A	N/A
						Excessive ethylene oxide residuals from sterilization process	Device will be released for distribution with acceptable residual levels (Design)	Sterilization Validation TR-0075, TR-HDX-0009					
		8.3	Exposure to pyrogens	Fever / pyrogenic response	3	Contacting material causes pyrogenic response	Material selection (Design)	Design verification – biocompatibility (Material mediated pyrogenicity) TR-0071	1	3	N/A	N/A	N/A
		8.4	Procedure complication – device inoperative	Delay in treatment	3	Material degrades over usable life	Material selection (Design)	Stability evaluation TR-0073	1	3	N/A	N/A	N/A
							Sterilization method selection (Design)						
9	Immunological	9.1	Exposure to non-biocompatible materials	Cell Damage	3	Patient contact materials are not biocompatibility	Material selection (Design)	Biocompatibility Test Report TR-0071	1	3	N/A	N/A	N/A
10	Function (Operational)	10.1	Procedure complication – device inoperative	Delay in treatment	2	Balloon Burst – Material or dimensional issues (wall thickness)	Material selection (Design) Dimensional Assessment (Design)	Design verification – Balloon Integrity (Resistance to Rupture) TR-0072, TR-HDX-0001	1	2	N/A	N/A	N/A
					2	Balloon Burst – Material or dimensional issues (wall thickness)	Material selection (Design) Dimensional Assessment (Design)	Design verification – Balloon Integrity (Resistance to Rupture) TR-0072, TR-HDX-0001	1	2	N/A	N/A	N/A
					1	Balloon deflates or does not inflates – Valve malfunctions or leaks (wall thickness issue)	Valve selection (Design) Dimensional Assessment (Design)	Design verification – Balloon Volume Maintenance TR-0072, TR-HDX-0001	1	1	N/A	N/A	N/A
					1	Balloon does not inflate – Valve malfunctions	Valve selection (Design)	Design verification TR-0072, TR-HDX-0001	1	1	N/A	N/A	N/A
					2	Catheter breaks under normal tensile stresses – Material or dimensional issues	Valve selection (Design) Dimensional Assessment (Design)	Design verification – Inflated Balloon Response to Traction TR-0072, TR-HDX-0001	1	2	N/A	N/A	N/A
					1	Connection issues – Syringe (inflation,	Connection will be standardized to what is	Design verification TR-0072, TR-HDX-0001	1	1	N/A	N/A	N/A



Hazard Analysis (HA) – Risk Management Report, SILQ ClearTract Foley Catheter

Doc No.: RM-0019

This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 7 of 9

Table 6: SILQ ClearTract Foley Catheter

#	Hazard Category	PH#	Potential Hazard	Potential Harm	S	Potential Failure Mode	Risk Controls ¹	Verification	P ¹	R ¹	Risk Controls ²	P ²	R ²
						irrigation); Collection bag (drainage)	incorporated with marketed catheters (Design)						
		10.2	Procedure complication – downgraded performance	Delay in treatment	2	Poor drainage – Insufficient holes and lumen diameter, Lumen collapses (Materials or dimensional issues)	Material selection (Design) Dimensional Assessment (Design)	Design verification – Flow Rate though Drainage Lumen TR-0072, TR-HDX-0001	1	2	N/A	N/A	N/A
					2	Poor irrigation - Insufficient holes and lumen diameter, Lumen collapses (Materials or dimensional issues)	Material selection (Design) Dimensional Assessment (Design)	Design verification TR-0072, TR-HDX-0001	1	2	N/A	N/A	N/A
		10.3	Blunt / rough surface	Urinary tract damage	3	Tip geometry – Damage during insertion	Tip design will be standardized to what is incorporated with marketed catheters (Design)	First Article TR-0072, TR-HDX-0001	1	3	N/A	N/A	N/A
					3	Irrigation pressure too high – Insufficient hole geometry	Valve selection (Design) Dimensional Assessment (Design)	Design verification TR-0072, TR-HDX-0001	1	3	N/A	N/A	N/A
		10.4	Procedure complication – Difficult to remove	Extended procedure / Urinary tract damage	3	Balloon will not deflate – Valve malfunctions	Valve selection (Design)	Design Verification - Deflation Reliability (Failure to Deflate) TR-0072, TR-HDX-0001	1	3	N/A	N/A	N/A
11	User Error (Operational)	11.1	Procedure complication – device inoperative	Delay in treatment	1	Catheter dislodges for patient – insufficient balloon inflation	Warning and limitations statements (Information) Define balloon fill volume	IFU SLQ-4006	2	2	N/A	N/A	N/A
		11.2	Partial / leaking connection to drain bag	Rash/Infection – Skin	2	Inadequate connection with collection bag	Product instruction for connection	IFU SLQ-4006	2	4	N/A	N/A	N/A
		11.3	Unskilled/ untrained clinician	Patient discomfort	2	Clinicians – does not read some or all information in IFU	Product labeling will clearly state to read the IFU prior to performing clinical application (Information)	Product Labeling SLQ-4000, SLQ-4001, SLQ-4006	2	4	N/A	N/A	N/A
							Warning and limitations statements (Information)	IFU SLQ-4006					
					3	Balloon not deflated or partially deflated during removal	Warning and limitations statements (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
					3	Excessive tensile – Catheter not secure, weight of collection bag pulls on catheter	Warning and limitations statements (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
		2	Balloon burst - excessive balloon inflation	Warning and limitations statements (Information) Define balloon fill volume	IFU SLQ-4006	2	4	N/A	N/A	N/A			



Hazard Analysis (HA) – Risk Management Report, SILQ ClearTract Foley Catheter

Doc No.: RM-0019

This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 8 of 9

Table 6: SILQ ClearTract Foley Catheter

#	Hazard Category	PH#	Potential Hazard	Potential Harm	S	Potential Failure Mode	Risk Controls ¹	Verification	P ¹	R ¹	Risk Controls ²	P ²	R ²
					2	Excessive irrigation pressure	Warning and limitations statements (Information)	IFU SLQ-4006	2	4	N/A	N/A	N/A
		11.4	Reuse of single use device	Fever/infection to patient	3	Device was reused	Warning and limitations statements (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
		11.5	Use of a compromised device	Fever/infection to patient	3	Sterile barrier system was not inspected for damage	Warning and limitations statements (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
		11.6	Use of expired device	Fever/infection to patient	3	Device used after expiration date or outside biocompatibility duration limits	Warning and limitations statements (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
		11.7	Mishandled during disposal	Fever/infection to clinician	3	Occurs during disposal process	Warning and limitations statements (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
12	Labeling (Product Labeling)	12.1	Loss of traceability	Cannot identify the device or trace back to MFG	1	Incomplete or ambiguous instructions for use	Label and labeling will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review FILE-004318, FILE-008308	1	1	N/A	N/A	N/A
		12.2	Use of expired device / reuse of single use device	Fever/infection to patient	3	Incomplete or ambiguous instructions for use	Label and labeling will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review FILE-004318, FILE-008308	1	3	N/A	N/A	N/A
13	Operating Instructions (Product Labeling)	13.1	Procedural Complication / misuse	Patient discomfort	3	Incomplete or ambiguous instructions for use	IFU will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review FILE-004318, FILE-008308	1	3	N/A	N/A	N/A
14	Environmental	14.1	Deteriorated device / Sterile barrier failure	Patient health is compromised	3	Stored outside the prescribe environmental parameters	Warning statement (Information)	IFU SLQ-4006	2	6	N/A	N/A	N/A
						Packaging is not adequate to protect the device through transportation	Packaging will be designed to protect the device and its sterile barrier through transportation stresses (Design)	Design Verification (Distribution Evaluation) TR-0074	1	3	N/A	N/A	N/A

8.0 Risk/ Benefit Assessment

8.1 The following table indicates the risk/benefit analysis for the residual risks that were greater than the risk acceptability limit after all applicable controls/mitigations were applied:

Table 7: Risk/Benefit Assessment for Residual Risk above Acceptability Limit for Process Hazard Assessment								
Table No.	Hazard Category	PH#	Potential Hazard / Harm	Potential Failure Mode	S	P ²	Risk/Benefit Analysis	Does the benefit outweigh the residual risk?
N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A

8.2 There was no residual risk identified that was above the risk acceptability limit per §6.3.

9.0 Conclusion

9.1 All Identified risks regardless of acceptability were reduced as far as possible (AFAP). All risks, including residual risks and regardless of their dimension, were reduced “as far as possible” (AFAP). Risks have been reduced as low as possible by the following methods:

- Design
- Appropriate testing to ensure design criteria were met
- Protective measures (critical-to-quality)
- Manufacturing and quality control measures

9.2 Based on a review of the risk-management process for the SILQ Foley Catheter described herein, the risk management plan based on ISO 14971:2019 has been implemented and the individual, cumulative and overall residual risk have been reduced as low as possible regardless of the magnitude and the benefit of using this device will outweigh risks posed by the device during usage.

10.0 Change History

Revision	Change Description
A	Initial release
B	Add HDX-160510, general updates, add verification references
C	Add SLQ-1610, add column Potential Harm to Table 6, general updates



Design Failure Mode Effects Analysis (FMEA) – Risk Management Report SILQ ClearTract Foley Catheter

Doc No.: RM-0020

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 1 of 8

1.0 Purpose

- 1.1 The purpose of this document is to identify the associated risk with the design of the components and their interactions with each other.

2.0 Scope

- 2.1 The product numbers associated with the SILQ ClearTract Foley Catheter is indicated in the following table:

Table 1: SILQ ClearTract Catheter Products	
Part No.	Description
HDX-1410	2-way 100% Silicone ClearTract Catheter, 14 Fr/10ml
HDX-1605	2-way 100% Silicone ClearTract Catheter, 16 Fr/5ml
HDX-1810	2-way 100% Silicone ClearTract Catheter, 18 Fr/10ml
HDX-160510	2-way 100% Silicone ClearTract Catheter, 16 Fr/5-10ml
SLQ-1610	2-Way 100% Silicone ClearTract Catheter, 16Fr/10ml

3.0 References

- 3.1 QP-004 – Risk Management

4.0 Risk Tables Definitions (ordered left to right)

- 4.1 # – Enumerate each Description with a sequential number
- 4.2 Description – Describes the operation that is being evaluated for risk
- 4.3 Potential Failure Mode (FM) – The potential consequences of the failure
- 4.4 Severity (S) – Measure of the possible consequences of a hazard.
- 4.5 Potential Effect of Failure – All conceivable effects of the potential failure mode
- 4.6 Potential Cause of Failure – All conceivable causes of the potential failure mode
- 4.7 Risk Controls - Process through which decisions are reached and protective measures are implemented for reducing risk to, or maintaining risk within specified levels.
- 4.7.1 Risk Controls¹ – Initial risk controls that are incorporated.
- 4.7.2 Risk Controls² – Secondary risk controls, also known as "Mitigation".



Design Failure Mode Effects Analysis (FMEA) – Risk Management Report SILQ ClearTract Foley Catheter

Doc No.: RM-0020

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 2 of 8

- 4.8 Verification – Evidence source of linkage that initial (i.e., Risk Controls¹) risk controls have implemented, verified, or validated.
- 4.9 Probability – Likelihood of the hazard occurring
- 4.9.1 Probability¹ (P¹) – The probability assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls¹)
- 4.9.2 Probability² (P²) – The probability assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls²)
- 4.10 Detection - The ability to screen out a hazard through verification activities (e.g., inspection).
- 4.10.1 Detection¹ (D¹) – The detection assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls¹)
- 4.10.2 Detection² (D²) – The detection assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls²)
- 4.11 Frequency – A 'Frequency Score' will be calculated using the ratings from probability and detection components. The frequency equation is the following: ***Frequency*** = $\sqrt{Px D}$ **Note:** The value calculated should be rounded up.
- 4.11.1 Frequency¹ (F¹) – The frequency score calculated using Risk Controls¹ and Detection¹
- 4.11.2 Frequency² (F²) – The frequency score calculated using Risk Controls² and Detection²
- 4.12 Risk Rating – Residual risk determination based on the concept of 'As Far As Possible' (AFAP) approach. The AFAP approach provides three (3) possible outcomes, which are 1) the residual risk exceeds the risk acceptability (RED), 2) the residual risk is acceptable to the risk acceptability criteria (GREEN), and 3) the residual risk is between the two (2) states indicated and it can be an option to lower the risk to the lowest practical level without economic considerations (YELLOW).
- 4.12.1 Risk Rating¹ (R¹) – This calculation is post the implementation of initial risk controls (i.e., Risk Controls¹)
- 4.12.2 Risk Rating² (R²) – This calculation is post the implementation of the secondary risk controls (i.e., Risk Controls²)

5.0 Risk Acceptability Determination

- 5.1 The criteria of risk acceptability will be determined with a Risk Priority Number (RPN) and on the concept of 'As Far As Possible (AFAP)' approach. Refer to the Risk Rating definition in §4.0.
- 5.2 There are three (3) components used in determining risk acceptability, which are Severity, Probability, and Detection. Refer to the following tables:

Table 2: Severity Ratings	
Index	Description
1	No patient effect
2	Minor injury that resolves without treatment



Design Failure Mode Effects Analysis (FMEA) – Risk Management Report SILQ ClearTract Foley Catheter

Doc No.: RM-0020

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 3 of 8

Table 2: Severity Ratings

Index	Description
3	Moderate injury that may require treatment
4	Serious injury requiring medical intervention
5	Life threatening injury, disability, death

Table 3: Probability Ratings

Index	Description
1	Improbable - <0.001%
2	Remote – 0.001% to .01%
3	Occasional – 0.01% to 0.1%
4	Likely – 0.1% to 1%
5	Frequent - >1%

Table 4: Detection Ratings

Index	Description
1	Very High – Current controls are almost certain to detect the failure mode.
2	High – High likelihood that current controls will detect the failure mode.
3	Moderate – Moderate likelihood that current controls will detect the failure mode.
4	Low – Low likelihood that current controls will detect the failure mode.
5	Remote – Remote likelihood that current controls will detect the failure mode.

5.3 Risk acceptability will be determined per the following risk index:

Table 5: Risk Rating Table

Severity	5	4	3	2	1
	5	4	3	2	1
	10	8	6	4	2
	15	12	9	6	3
	20	16	12	8	4
	25	20	15	10	5
=Square Root (Probability X Detection)					

Note: The purpose of using the square root of the probability and detection scores is to deemphasize these two (2) components, which places more emphasize on the severity component.

#	Component Description - Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
1	Catheter tube transfers urine from patient into drainage bag	Poor drainage / occluded lumen	2	Patient discomfort – device performance	Insufficient holes, lumen diameter, wall thickness	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (flow test)	1	2	2	4	N/A	N/A	N/A	N/A	N/A
		Catheter breaks	2	Procedural Complication – device failure	Wall thickness	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (inflated balloon response to traction)	1	2	2	4	N/A	N/A	N/A	N/A	N/A
					Performance degrades over time	Material Selection (Design)	Stability Evaluation TR-0073									
		Oversized tube	3	Procedural Complication – Urethral damage	Wall thickness	Design follows standard ISO 20969	TR-0072, TR-HDX-0001 (Balloon size and shaft size)	1	2	2	6	N/A	N/A	N/A	N/A	N/A
2	Catheter tube allows for irrigation	Poor irrigation/ occluded lumen	1	Procedural Complication – device failure	Insufficient holes, lumen diameter, wall thickness	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (flow test)	1	2	2	4	N/A	N/A	N/A	N/A	N/A
3	Balloon inflates to secure catheter in place	Balloon leaks / deflates	2	Procedural Complication – Catheter dislodges	Wall thickness	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (Balloon Volume Maintenance)	1	2	2	4	N/A	N/A	N/A	N/A	N/A
					Performance degrades over time	Material Selection (Design)	Stability Evaluation TR-0073									
		Balloon won't inflate	1	Procedural Complication – No patient harm	Wall thickness	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (Balloon Volume Maintenance)	1	2	2	2	N/A	N/A	N/A	N/A	N/A
		Balloon bursts	2	Procedural Complication	Wall thickness	Design follows standard ASTM F623, Fill volume	TR-0072, TR-HDX-0001 (Balloon integrity – resistance to rupture)	1	2	2	4	N/A	N/A	N/A	N/A	N/A
Performance degrades over time	Material Selection (Design)				Stability Evaluation TR-0073											

Table 6: Product Risk Assessment

#	Component Description - Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²		
4	Valve allows inflation and deflation of balloon and connects with a syringe	Leaking Valve – fail to seal	2	Procedural Complication - Catheter dislodges	Valve malfunction	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (Balloon Volume Maintenance)	1	2	2	4	N/A	N/A	N/A	N/A	N/A		
					Performance degrades over time	Material Selection (Design)	Stability Evaluation TR-0073											
		Balloon won't deflate – valve fails to open	3	Procedural Complication – difficult to remove	Valve malfunction	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (Deflation Reliability – failure to deflate)	1	2	2	4	N/A	N/A	N/A	N/A	N/A		
					Performance degrades over time	Material Selection (Design)	Stability Evaluation TR-0073											
		Syringe won't connect to valve	1	Complications in manufacturing assembly process	Not designed to Luer standards	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (Verification Non-Interconnectable Connector)	1	1	1	1	N/A	N/A	N/A	N/A	N/A		
		Valve disconnects from catheter	1	Procedural Complication – device failure	Inadequate design	Design follows standard ASTM F623	TR-0072, TR-HDX-0001 (Verification Non-Interconnectable Connector)	1	2	2	2	N/A	N/A	N/A	N/A	N/A		
		5	Pouch – Provide sterile barrier	Breach in sterile barrier	3	Non-sterile device	Inadequate seal	Material selection (Device)	Design verification (seal integrity) TR-0074	1	1	1	3	N/A	N/A	N/A	N/A	N/A
							Inadequate material (e.g., sterilization process, product weight, sharp edges)		Design verification (distribution evaluation) TR-0074									
Performance degrades over time	Stability Evaluation TR-0073																	



Design Failure Mode Effects Analysis (FMEA) – Risk Management Report

SILQ ClearTract Foley Catheter

Doc No.: RM-0020

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 6 of 8

Table 6: Product Risk Assessment

#	Component Description - Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
6	Information for Use (DFU) – Provides device information	Complication with device use/ intended purpose	3	Patient harm	Incomplete or ambiguous instructions for use	IFU will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review FILE-004318, FILE-008308	1	1	1	3	N/A	N/A	N/A	N/A	N/A
7	Inner Carton – Protects and orients the pouched devices	Packaging is damaged	3	Device/ sterile barrier damage	Inadequate design	Material/ card weight selection (Design) Geometry is design to handle anticipate load (Design)	Design verification (distribution evaluation) TR-0074	1	1	1	3	N/A	N/A	N/A	N/A	N/A
8	Label, Inner Carton – Provides device information	Identification issues	1	Nuisance in preparation for the clinical application	Incomplete or ambiguous product information	Label and labeling will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review FILE-004318, FILE-008308	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Device use outside specification	3	Patient harm	Inadequate disclosure of specifications/ limitation (e.g., single-use, expiry)	Label and labeling will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review FILE-004318, FILE-008308	1	1	1	3	N/A	N/A	N/A	N/A	N/A
9	Shipper Case – Protects packaged devices	Packaging is damaged	3	Device/ sterile barrier damage	Inadequate design	Material/ card weight selection (Design) Geometry is design to handle anticipate load (Design)	Design verification (distribution evaluation) TR-0074	1	1	1	3	N/A	N/A	N/A	N/A	N/A



Design Failure Mode Effects Analysis (FMEA) – Risk Management Report SILQ ClearTract Foley Catheter

Doc No.: RM-0020

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 7 of 8

Table 6: Product Risk Assessment

#	Component Description - Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
10	Label, Shipper Case – Provides device information	Identification issues	1	Nuisance in preparation for the clinical application	Incomplete or ambiguous product information	Label and labeling will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review FILE-004318, FILE-008308	1	1	1	1	N/A	N/A	N/A	N/A	N/A

7.0 Risk/ Benefit Assessment

7.1 The following table indicates the risk/benefit analysis for the residual risks that were greater than the risk acceptability limit after all applicable controls/mitigations were applied:

Table 7.1 – Risk/Benefit Assessment for Residual Risk above Acceptability Limit for Process Hazard Assessment

Table No.	Process Description	Section No.	Potential Failure Mode	Potential Effect of Failure	Potential Cause of Failure	S	F ²	Risk/Benefit Analysis	Does the benefit outweigh the residual risk?
N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A

7.2 There was no residual risk identified that was above the risk acceptability limit per §5.3.

8.0 Conclusion

8.1 Risks have been reduced as low as possible by the following methods:

- Design the product with inherent safety
- Protective measures/ mechanisms to the product of the device.
- Provide information for safety

Based on a review of the risk-management process for the SILQ ClearTract Foley Catheter described herein, the risk management plan based on ISO 14971:2019 has been implemented and the individual, cumulative and overall residual risk have been reduced as low as possible regardless of the magnitude and the benefit of using this device will outweigh risks posed by the device during usage.

	Design Failure Mode Effects Analysis (FMEA) – Risk Management Report SILQ ClearTract Foley Catheter		Doc No.: RM-0020
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 8 of 8

9.0 Change History

Revision	Change Description
A	Initial release
B	Add HDX-160510, general updates, add verification references
C	Add SLQ-1610

	Master Validation Plan, SILQ ClearTract Catheter	Doc No.: TP-HDX-0007
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 1 of 3

1.0 Purpose

- 1.1 The Master Validation Plan (MVP) will document the validation requirements for processes that will be used to manufacture and verify the finished product of the SILQ ClearTract Catheter.
 - 1.1.1 Note: The product name SILQ ClearTract Catheter and Hydrophilix Foley Catheter is used interchangeably. The company name Hydrophilix and SILQ are used interchangeably as the name recently was changed from Hydrophilix to SILQ.

2.0 Scope

- 2.1 This document applies to the manufacturing and verification processes for the SILQ ClearTract Foley Catheter.

3.0 References

- 3.1 Project#00047 – Hydrophilix Foley Catheter Product Line & Shelf Life Extension (US)
- 3.2 Project#00049 – SILQ Foley Catheter Product Line (US)
- 3.3 DMR HDX-800X – Device Master Record, Hydrophilix Foley Catheter
- 3.4 DMR SLQ-800X – Device Master Record, SILQ Foley Catheter
- 3.5 FM-0055 – Form, Installation Qualification Checklist
- 3.6 IP-0002 – Inspection Procedure, Product Labeling and Packaging
- 3.7 MP-0016 – Manufacturing Procedure, Label Generation
- 3.8 MP-HDX-0001 – Manufacturing Procedure, Hydrophilix Foley Catheter, Labeling and Packaging
- 3.9 MP-HDX-0004 – Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing Automation Machine
- 3.10 MP-SLQ-0001 – Manufacturing Procedure, SILQ Foley Catheter, Labeling and Packaging
- 3.11 MP-SLQ-0002 – Manufacturing Procedure, SILQ Polymer Coating Application Utilizing Automation Machine
- 3.12 QP-002 – Document Control
- 3.13 QP-003 – Quality Records
- 3.14 QP-012 – Process Validation
- 3.15 QP-030 – Statistical Techniques
- 3.16 TM-HDX-0001 – Test Method, UV Light Absorbance, Hydrophilix Catheters
- 3.17 TM-HDX-0002 – Test Method, Rinse Test, Hydrophilix Catheter
- 3.18 TP-0104 – Test Protocol, Measurement System Analysis for UV Light Absorbance Test, Hydrophilix Products (TM-HDX-0001)
- 3.19 TP-0106 – Test Protocol, Operation Qualification (OQ), Sealing Process, Impulse Sealer

 pathway	Master Validation Plan, SILQ ClearTract Catheter		Doc No.: TP-HDX-0007
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 2 of 3

- 3.20 TP-HDX-0004 – Performance Qualification (PQ), Sealing Process, Impulse Sealer, Hydrophilix
- 3.21 TP-HDX-0005 – Test Protocol, Operation Qualification (OQ), Hydrophilix Surface Treatment Process
- 3.22 TP-HDX-0006 – Test Protocol, Performance Qualification (PQ), Hydrophilix Surface Treatment Process
- 3.23 TP-HDX-0008 – Test Protocol, Manufacturing Line Evaluation (Process Validation), Hydrophilix Foley Catheter
- 3.24 TP-SAL-0003 – Test Protocol, Operation Qualification (OQ), SteriPax Pouch Sealing Process Continuous Band Sealer
- 3.25 TP-SAL-0004 – Test Protocol, Performance Qualification (PQ), SteriPax Pouch Sealing Process Emplex Model MPS6300-M
- 3.26 WI-0013 – Work Instruction, Gage Repeatability & Reproducibility Evaluation

4.0 Responsibilities

- 4.1 Production/Manufacturing will be responsible for execution of this plan.

5.0 Sampling Requirements

- 5.1 Based on the activities being conducted, sampling will be selected in conjunction of QP-012, QP-030, and WI-0013. Sampling requirements will be documented with the protocol.

6.0 General Requirements

- 6.1 Protocols generated to support the MVP will be reviewed and approved per QP-002 prior to starting associated activities. Report will be accompanied with the appropriate attachments gathered during the activity then will be reviewed and approved per QP-002.
- 6.2 Good documentation practice will be followed when generating records per QP-003.
- 6.3 All requirements indicated in the approved protocol must be addressed.
 - 6.3.1 If a requirement is not addressed, rationale for not completing the requirement will be documented in the deviation section of the report.
- 6.4 Acceptance criteria will be indicated within the protocol.
- 6.5 Individuals executing any portion of the test protocol must read and understand the portions of the protocol that applies, and a training record will be completed prior to performing the activity.
- 6.6 Results from verifications activities will be collected and recorded on the provided data collection sheet that is released with the protocol. Results shall confirm indicated acceptance criteria or exceptional conditions must be identified, documented, investigated, and addressed.

7.0 Validation Approach

- 7.1 Validation activities associated to Installation Qualification (IQ), Operation Qualification (OQ), and Process Qualification (PQ) will be conducted per QP-012.

	Master Validation Plan, SILQ ClearTract Catheter		Doc No.: TP-HDX-0007
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 3 of 3

7.2 Measurement System Analysis (MSA) for verification activities will be conducted per WI-0013.

7.3 The following table provides the details of the manufacturing processes:

Table 1: Process Validation Requirements							
#	Process	Doc#	Equipment	Validation Required?	Validation Activity	Guidance Procedure	Validation Doc #
1	Labeling (Printing)	MP-0016 IP-0002	Label Printer	No 100% Verifiable	IQ	QP-012	FM-0055
2	Surface Treatment	MP-HDX-0004 MP-SLQ-0002	Automated Coater (ES-0091-01)	Yes	IQ / OQ / PQ	QP-012	FM-0055 TP-HDX-0005 TP-HDX-0006
3	UV Light Absorbance Test	TM-HDX-0001	Scenstar VIS Spectrometer (ES-0092-01)	Yes	IQ / MSA	QP-012 WI-0013	FM-0055 TP-0104
4	Rinse Test	TM-HDX-0002					
5	Sleeve Loading	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A
6	Pouch Label Application	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A
7	Pouch Loading	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A
8	Pouch Sealing (Sterile Barrier System)	MP-HDX-0001 MP-SLQ-0001	Pouch Sealer (ES-0059-01)	Yes	IQ /OQ /PQ	QP-012	FM-0055 TP-SAL-0003 TP-SAL-0004
			Pouch Sealer (ES-0099-01)				FM-0055 TP-0106 TP-HDX-0004
9	Box Assembly & Label Application	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A
10	Packaging	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A
11	Entire Process	MP-HDX-0001, MP-HDX-0004 MP-SLQ-0001, MP-SLQ-0002	N/A	Yes	PQ	QP-012	TP-HDX-0008

7.4 Revalidation requirements are indicated in QP-012.

8.0 Attachments

8.1 None

9.0 Change History

Revision	Change Description
A	Initial release
B	Add Project#49, Add references to SLQ numbered procedures

	Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter	Doc No.: TP-HDX-0008
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 1 of 4

1.0 Purpose

- 1.1 This protocol provides the requirements for conducting a manufacturing line evaluation for the SILQ ClearTract Foley Catheter treatment, labeling and packaging process to address the process validation requirement.

Note: The product name SILQ ClearTract Catheter and Hydrophilix Foley Catheter is used interchangeably. The company name Hydrophilix and SILQ are used interchangeably as the name recently was changed from Hydrophilix to SILQ.

2.0 Scope

- 2.1 Applies to the SILQ ClearTract Foley Catheter.

3.0 References

- 3.1 Project#00047 – Hydrophilix Foley Catheter Product Line & Shelf Life Extension (US)
- 3.2 Project#00049 – SILQ Foley Catheter Product Line (US)
- 3.3 DMR HDX-800X – Device Master Record, Hydrophilix Foley Catheter
- 3.4 DMR SLQ-800X – Device Master Record, SILQ Foley Catheter
- 3.5 QP-002 – Document Control
- 3.6 QP-003 – Quality Records
- 3.1 QP-009 – Receiving Inspection
- 3.2 QP-010 – Final Product Inspection
- 3.3 QP-012 – Process Validation
- 3.4 MP-HDX-0001 – Manufacturing Procedure, Hydrophilix Foley Catheter, Labeling and Packaging
- 3.5 MP-HDX-0004 – Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing Automation Machine
- 3.6 MP-SLQ-0001 – Manufacturing Procedure, SILQ Foley Catheter, Labeling and Packaging
- 3.7 MP-SLQ-0002 – Manufacturing Procedure, SILQ Polymer Coating Application Utilizing Automation Machine
- 3.8 TP-HDX-0007 – Test Protocol, Master Validation Plan, SILQ ClearTract Foley Catheter
- 3.9 WI-0004 – Work Instruction, Arena PLM - Deviations

4.0 Responsibilities

- 4.1 It is the responsibility of the individual conducting the evaluation to comply with requirements indicated in this procedure.

5.0 Sampling Requirements

	Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter	Doc No.: TP-HDX-0008
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 2 of 4

- 5.1 A minimum of three (3) runs is required to fulfill the manufacturing line evaluation.

6.0 Material and Equipment

- 6.1 The equipment utilized are described within manufacturing, test, and inspection procedures referenced in §3.0.
- 6.2 The material utilized are indicated in the DMRs referenced in §3.0.

7.0 Acceptance Criteria

- 7.1 Device History Record (DHR) Review
- 7.1.1 DHR review requirements indicated in QP-010 shall be met.

8.0 Procedure

8.1 General Requirements

- 8.1.1 The product manufactured in the manufacturing line evaluation is intended for Human Use distribution. The product will be quarantined immediately after processing. The product will only be released from quarantine, if all acceptance criteria indicated in the protocol are met and documented in an interim report for an individual lot, or a final report to disposition all three lots.
- 8.1.2 The PQ protocol and Device Master Record (DMR), which will include the BOM, Specification (Component, Sub-Assembly, and Finish Goods), Manufacturing Assembly, Test, and Inspection Procedures, will be approved per QP-002 prior to starting the PQ.
- 8.1.3 Validation activities will be performed by personnel trained to the protocol and/or all referenced documents. Training records will be attached to the final report and incorporate in the employees training file.
- 8.1.4 Good documentation practice will be followed when generating records per QP-003.
- 8.1.5 Materials used in the PQ have met receiving inspection requirement per QP-009.
- 8.1.6 Work orders will be generated with the PQ lots, which a unique lot number will be assigned.
- 8.1.7 Proper identification, segregation, line clearance will be conducted to eliminate potential mix-up
- 8.1.8 All equipment used will be calibrated (if applicable). Equipment asset numbers and calibration due dates will be documented on the appropriate data collection sheets.
- 8.1.9 Any deviation to the approved protocol will be documented in the final report.
- 8.1.10 Changes to Manufacturing Documentation
- 8.1.10.1 Deviations in the manufacturing documentation, as it relates to the process, during these builds, shall be made per WI-0004.

8.2 Manufacture of Process Validation Lots

- 8.2.1 Manufacture the product per released manufacturing assembly procedure referenced in §3.0

	Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter	Doc No.: TP-HDX-0008
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 3 of 4

8.2.2 Conduct verification activities per released test and inspection procedure referenced in §3.0.

8.3 The Device History Record (DHR)

8.3.1 The DHR will be reviewed per QP-010. Result of this review will be documented in the final report.

	Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter	Doc No.: TP-HDX-0008
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 4 of 4

9.0 Records

- 9.1 The records generated from this evaluation will be incorporated into a report, which is be retained per QP-002.

10.0 Attachments

- 10.1 None

11.0 Change History

Revision	Change Description
A	Initial release
B	Add Project#49, add SLQ numbered procedures

	Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter	Doc No.: TR-HDX-0008
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 1 of 4

1.0 Purpose

- 1.1 This report provides the results of conducting a manufacturing line evaluation for the SILQ ClearTract Foley Catheter treatment, labeling and packaging process to address the process validation requirement.

Note: The product name SILQ ClearTract Catheter and Hydrophilix Foley Catheter is used interchangeably. The company name Hydrophilix and SILQ are used interchangeably as the name recently was changed from Hydrophilix to SILQ.

2.0 Scope

- 2.1 Applies to the SILQ ClearTract Foley Catheter.

3.0 References

- 3.1 DMR HDX-160510 – Device Master Record, Hydrophilix Foley Catheter 16 FR x 5 ml – 10 ml, Sterile, 10-Pack
- 3.2 DMR SLQ-1610 – Device Master Record, SILQ Foley Catheter 16 FR x 10 ml, Sterile, 10-Pack
- 3.3 Project#00047 – Hydrophilix Foley Catheter Product Line & Shelf Life Extension (US)
- 3.4 Project#00049 – SILQ Foley Catheter Product Line (US)
- 3.5 MP-HDX-0001 rev C – Manufacturing Procedure, Hydrophilix Foley Catheter, Labeling and Packaging
- 3.6 MP-HDX-0004 rev B – Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing Automation Machine
- 3.7 MP-SLQ-0001 rev A – Manufacturing Procedure, SILQ Foley Catheter, Labeling and Packaging
- 3.8 MP-SLQ-0002 rev A – Manufacturing Procedure, SILQ Polymer Coating Application Utilizing Automation Machine
- 3.9 QP-002 rev G – Document Control
- 3.10 QP-003 rev F – Quality Records
- 3.11 QP-009 rev D – Receiving Inspection
- 3.12 QP-010 rev F – Final Product Inspection
- 3.13 QP-012 rev C – Process Validation
- 3.14 TP-HDX-0008 rev B – Test Protocol, Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter

4.0 Acceptance Criteria

- 4.1 Device History Record (DHR) Review
 - 4.1.1 DHR review requirements indicated in QP-010 shall be met.

	Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter	Doc No.: TR-HDX-0008
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 2 of 4

5.0 Deviations

5.1 None

6.0 Results & Discussion

6.1 Prior to execution, the following items were completed:

6.1.1 The Device Master Record was approved and released per:.

Table 1: DMR Summary	
DMR No.	ECO#
HDX-8006	ECO-000822
SLQ-8004	ECO-000907

6.1.2 The PQ protocol was approved and released per ECO-000859. It was updated and released per ECO-000921.

6.1.3 All production and quality control associates that took part in these manufacturing builds were trained and training records were generated. Refer to the attached training records (8.1).

6.1.4 Good documentation practice was followed when generating records per QP-003.

6.1.5 All material used for the manufacturing lots were inspected and released per QP-009.

6.1.6 Unique lot numbers were assigned and documented per work order.

6.1.7 Proper identification, segregation, line clearance was conducted and documented per work order.

6.1.8 All equipment used was calibrated (if applicable). Equipment asset numbers and calibration due dates were documented on the appropriate data collection sheets.

6.1.9 There were no deviations to the protocol.

6.2 Manufactured Lot Summary

6.2.1 The following table indicates the general information of the manufacturing lots produced during the manufacturing line evaluation:

Table 2: Manufactured Lot Summary				
#	Part No.	Rev	Lot No.	Build Qty
1	HDX-160510	A	HDX-8006110821	350
2	SLQ-1610	A	SLQ1610-012022	609
3	SLQ-1610	A	SLQ1610-012122	483

6.2.2 Three lots shall be released. This completes the manufacturing line evaluation.

6.3 The Device History Record (DHR) Review

	Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter	Doc No.: TR-HDX-0008
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 3 of 4

6.3.1 The DHR was reviewed per QP-010 and documented in the work order. The DHR met the acceptance criteria.

	Manufacturing Line Evaluation, SILQ ClearTract Foley Catheter	Doc No.: TR-HDX-0008
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: B
		Page 4 of 4

7.0 Conclusion

- 7.1 Based on the results, the lots successfully met all the requirements per TP-HDX-0008 rev B. One HDX-160510 lot and two SLQ-1610 lots will be released from quarantine to commercial distribution.

8.0 Attachments

- 8.1 Training Records

9.0 Change History

Revision	Change Description
A	Initial release
B	Add references to Project#00049, DMR SLQ-8XXX, MP-SLQ-0001, MP-SLQ-0002, Add Table 1, Add lots 2 and 3 to Table 2.

	Master Validation Execution, SILQ ClearTract Catheter	Doc No.: TR-HDX-0007
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 1 of 5

1.0 Purpose

- 1.1 This report provides the results of executing the master validation plan for the SILQ ClearTract Foley Catheter.

Note: The product name SILQ ClearTract Catheter and Hydrophilix Foley Catheter is used interchangeably. The company name Hydrophilix and SILQ are used interchangeably as the name recently was changed from Hydrophilix to SILQ.

2.0 Scope

- 2.1 Applies to the SILQ ClearTract Foley Catheter.

3.0 References

- 3.1 TP-HDX-0007 rev B – Test Protocol, Master Validation Plan, SILQ ClearTract Catheter
- 3.2 FM-0055 – Form, Installation Qualification Checklist
- 3.3 IP-0002 rev B – Inspection Procedure, Product Labeling and Packaging
- 3.4 MP-0016 rev A – Manufacturing Procedure, Label Generation
- 3.5 MP-HDX-0001 rev C – Manufacturing Procedure, Hydrophilix Foley Catheter, Labeling and Packaging
- 3.6 MP-HDX-0004 rev B – Manufacturing Procedure, Hydrophilix Polymer Coating Application Utilizing Automation Machine
- 3.7 MP-SLQ-0001 rev A – Manufacturing Procedure, SILQ Foley Catheter, Labeling and Packaging
- 3.8 MP-SLQ-0002 rev A – Manufacturing Procedure, SILQ Polymer Coating Application Utilizing Automation Machine
- 3.9 QP-002 rev G – Document Control
- 3.10 QP-003 rev F – Quality Records
- 3.11 QP-012 rev C – Process Validation
- 3.12 TM-HDX-0001 rev A – Test Method, UV Light Absorbance, Hydrophilix Catheters
- 3.13 TM-HDX-0002 rev A – Test Method, Rinse Test, Hydrophilix Catheter
- 3.14 TP-0104 rev A – Test Protocol, Measurement System Analysis for UV Light Absorbance Test, Hydrophilix Products (TM-HDX-0001)
- 3.15 TP-0106 rev A – Test Protocol, Operation Qualification (OQ), Sealing Process, Impulse Sealer
- 3.16 TP-HDX-0004 rev A – Performance Qualification (PQ), Sealing Process, Impulse Sealer, Hydrophilix
- 3.17 TP-HDX-0005 rev A – Test Protocol, Operation Qualification (OQ), Hydrophilix Surface Treatment Process

	Master Validation Execution, SILQ ClearTract Catheter	Doc No.: TR-HDX-0007
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 2 of 5

- 3.18 TP-HDX-0006 rev A – Test Protocol, Performance Qualification (PQ), Hydrophilix Surface Treatment Process
- 3.19 TP-HDX-0008 rev B – Test Protocol, Manufacturing Line Evaluation (Process Validation), Hydrophilix Foley Catheter
- 3.20 TP-SAL-0003 rev A – Test Protocol, Operation Qualification (OQ), SteriPax Pouch Sealing Process Continuous Band Sealer
- 3.21 TP-SAL-0004 rev A – Test Protocol, Performance Qualification (PQ), SteriPax Pouch Sealing Process Emplex Model MPS6300-M
- 3.22 TR-0090 rev A – Test Report, Installation Qualification (IQ) Checklist, Zebra Printer GX430t
- 3.23 TR-0122 rev A – Test Report, Measurement System Analysis for UV Light Absorbance Test, Hydrophilix Products (TM-HDX-0001)
- 3.24 TR-0126 rev A – Test Report, Operation Qualification (OQ), Sealing Process, Impulse Sealer
- 3.25 TR-0129 rev A – Test Report, ES-0099-01 Installation Qualification (IQ) ACCU-SEAL SenocorpWhite Impulse Sealer
- 3.26 TR-0137 rev A – Test Report, ES-0092-01 Installation Qualification (IQ) Scenstar VIS Spectrometer
- 3.27 TR-0139 rev A – Test Report, Installation Qualification, IQ, ZEBRA GX430T Label Printer
- 3.28 TR-0140 rev A – Test Report, Installation Qualification, IQ, SQUIX 4M Label Printer
- 3.29 TR-HDX-0004 rev A – Test Report, Performance Qualification (PQ), Sealing Process, Impulse Sealer, SILQ
- 3.30 TR-HDX-0005 rev A – Test Report, Operation Qualification (OQ), Hydrophilix Surface Treatment Process
- 3.31 TR-HDX-0006 rev A – Test Report, Performance Qualification (PQ), Hydrophilix Surface Treatment Process
- 3.32 TR-HDX-0008 rev B – Test Report, Manufacturing Line Evaluation (Process Validation), SILQ ClearTract Foley Catheter
- 3.33 TR-SAL-0000 rev A – Test Report, Installation Qualification (IQ), Plexpack Continuous Band Sealer
- 3.34 TR-SAL-0003 rev A – Test Report, Operation Qualification (OQ), SteriPax Pouch Continuous Band Sealer
- 3.35 TR-SAL-0004 rev A – Test Report, Performance Qualification (PQ), SteriPax Pouch Sealing Process Emplex Model MPS6300-M
- 3.36 WI-0013 rev A – Work Instruction, Gage Repeatability & Reproducibility Evaluation

	Master Validation Execution, SILQ ClearTract Catheter	Doc No.: TR-HDX-0007
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 3 of 5

4.0 Acceptance Criteria

4.1 Confirm MVP requirements have been addressed.

- 4.1.1 If a requirement is not addressed, rationale for not completing the requirement will be documented in the deviation section of the report.

5.0 Deviations

5.1 None

6.0 Results & Discussion

6.1 The following general requirements were met during the execution of the master validation plan:

- 6.1.1 Protocols generated to support the MVP were reviewed and approved per QP-002 prior to starting associated activities.
- 6.1.2 Acceptance criteria were indicated within the protocols.
- 6.1.3 All requirements indicated in their approved protocols were addressed.
- 6.1.4 Good documentation practices were followed when generating records per QP-003.
- 6.1.5 Individuals executing any portion of the test protocol read and understood the portions of the protocol that applied, and training records were completed prior to performing the activity.
- 6.1.6 Results from verifications activities were collected and recorded on data collection sheets. Results confirmed indicated acceptance criteria was met.
- 6.1.7 Reports were accompanied with their attachments, and were reviewed and approved per QP-002.

 pathway	Master Validation Execution, SILQ ClearTract Catheter		Doc No.: TR-HDX-0007
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 4 of 5

6.2 The following table provides a summary of the master validation execution:

Table 1: Master Validation Execution Summary								
#	Process	Doc#	Equipment	Validation Required?	Validation Activity	Guidance Procedure	Validation Doc #	Validation Report #
1	Labeling (Printing)	MP-0016 IP-0002	Label Printer	No 100% Verifiable	IQ	QP-012	FM-0055	TR-0090 TR-0139 TR-0140
2	Surface Treatment	MP-HDX-0004 MP-SLQ-0002	Automated Coater (ES-0091-01)	Yes	IQ / OQ / PQ	QP-012	FM-0055 TP-HDX-0005 TP-HDX-0006	TR-0117 TR-HDX-0005 TR-HDX-0006
3	UV Light Absorbance Test	TM-HDX-0001	Scenstar VIS Spectrometer (ES-0092-01)	Yes	IQ / MSA	QP-012 WI-0013	FM-0055 TP-0104	TR-0137 TR-0122
4	Rinse Test	TM-HDX-0002						
5	Sleeve Loading	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A	N/A
6	Pouch Label Application	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A	N/A
7	Pouch Loading	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A	N/A
8	Pouch Sealing (Sterile Barrier System)	MP-HDX-0001 MP-SLQ-0001	Pouch Sealer (ES-0059-01)	Yes	IQ /OQ /PQ	QP-012	FM-0055 TP-SAL-0003 TP-SAL-0004	TR-SAL-0000 TR-SAL-0003 TR-SAL-0004
			Pouch Sealer (ES-0099-01)				FM-0055 TP-0106 TP-HDX-0004	TR-0129 TR-0126 TR-HDX-0004
9	Box Assembly & Label Application	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A	N/A
10	Packaging	MP-HDX-0001 MP-SLQ-0001	N/A	No 100% Verifiable	N/A	N/A	N/A	N/A
11	Entire Process	MP-HDX-0001, MP-HDX-0004 MP-SLQ-0001, MP-SLQ-0002	N/A	Yes	PQ	QP-012	TP-HDX-0008	TR-HDX-0008

	Master Validation Execution, SILQ ClearTract Catheter		Doc No.: TR-HDX-0007
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 5 of 5

7.0 Conclusion

- 7.1 Based on the results, the master validation met all the requirements per TP-HDX-0007 rev B.

8.0 Attachments

- 8.1 None

9.0 Change History

Revision	Change Description
A	Initial release

1.0 Purpose

- 1.1 The purpose of this document is to identify the associated risk with the manufacturing processes of the SILQ ClearTract Foley Catheter.

2.0 Scope

- 2.1 The product numbers associated with the SILQ ClearTract Foley Catheter is indicated in the following table:

Table 1: SILQ ClearTract Catheter Products	
Part No.	Description
HDX-1410	2-way 100% Silicone ClearTract Catheter, 14 Fr/10ml
HDX-1605	2-way 100% Silicone ClearTract Catheter, 16 Fr/5ml
HDX-1810	2-way 100% Silicone ClearTract Catheter, 18 Fr/10ml
HDX-160510	2-way 100% Silicone ClearTract Catheter, 16 Fr/5-10ml
SLQ-1610	2-Way 100% Silicone ClearTract Catheter, 16Fr/10ml

3.0 References

- 3.1 QP-004 – Risk Management

4.0 Risk Tables Definitions

- 4.1 # – Enumerate each Description with a sequential number
- 4.2 Description – Describes the operation that is being evaluated for risk
- 4.3 Potential Failure Mode (FM) – The potential consequences of the failure
- 4.4 Severity (S) – Measure of the possible consequences of a hazard.
- 4.5 Potential Effect of Failure – All conceivable effects of the potential failure mode
- 4.6 Potential Cause of Failure – All conceivable causes of the potential failure mode
- 4.7 Risk Controls - Process through which decisions are reached and protective measures are implemented for reducing risk to, or maintaining risk within specified levels.
 - 4.7.1 Risk Controls¹ – Initial risk controls that are incorporated.
 - 4.7.2 Risk Controls² – Secondary risk controls, also known as "Mitigation".

- 4.8 Verification – Evidence source of linkage that initial (i.e., Risk Controls¹) risk controls have implemented, verified, or validated.
- 4.9 Probability – Likelihood of the hazard occurring
- 4.9.1 Probability¹ (P¹) – The probability assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls¹)
- 4.9.2 Probability² (P²) – The probability assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls²)
- 4.10 Detection - The ability to screen out a hazard through verification activities (e.g., inspection).
- 4.10.1 Detection¹ (D¹) – The detection assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls¹)
- 4.10.2 Detection² (D²) – The detection assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls²)
- 4.11 Frequency – A 'Frequency Score' will be calculated using the ratings from probability and detection components. The frequency equation is the following: ***Frequency*** = $\sqrt{Px D}$ **Note:** The value calculated should be rounded up.
- 4.11.1 Frequency¹ (F¹) – The frequency score calculated using Risk Controls¹ and Detection¹
- 4.11.2 Frequency² (F²) – The frequency score calculated using Risk Controls² and Detection²
- 4.12 Risk Rating – Residual risk determination based on the concept of 'As Far As Possible' (AFAP) approach. The AFAP approach provides three (3) possible outcomes, which are 1) the residual risk exceeds the risk acceptability (RED), 2) the residual risk is acceptable to the risk acceptability criteria (GREEN), and 3) the residual risk is between the two (2) states indicated and it can be an option to lower the risk to the lowest practical level without economic considerations (YELLOW).
- 4.12.1 Risk Rating¹ (R¹) – This calculation is post the implementation of initial risk controls (i.e., Risk Controls¹)
- 4.12.2 Risk Rating² (R²) – This calculation is post the implementation of the secondary risk controls (i.e., Risk Controls²)

5.0 Risk Acceptability Determination

- 5.1 The criteria of risk acceptability will be determined with a Risk Priority Number (RPN) and on the concept of 'As Far As Possible (AFAP)' approach. Refer to the Risk Rating definition in §4.0.
- 5.2 There are three (3) components used in determining risk acceptability, which are Severity, Probability, and Detection. Refer to the following tables:

Table 2: Severity Ratings	
Index	Description
1	No patient effect
2	Minor injury that resolves without treatment

Table 2: Severity Ratings

Index	Description
3	Moderate injury that may require treatment
4	Serious injury requiring medical intervention
5	Life threatening injury, disability, death

Table 3: Probability Ratings

Index	Description
1	Improbable - <0.001%
2	Remote – 0.001% to .01%
3	Occasional – 0.01% to 0.1%
4	Likely – 0.1% to 1%
5	Frequent - >1%

Table 4: Detection Ratings

Index	Description
1	Very High – Current controls are almost certain to detect the failure mode.
2	High – High likelihood that current controls will detect the failure mode.
3	Moderate – Moderate likelihood that current controls will detect the failure mode.
4	Low – Low likelihood that current controls will detect the failure mode.
5	Remote – Remote likelihood that current controls will detect the failure mode.

5.3 Risk acceptability will be determined per the following risk index:

Table 5: Risk Rating Table

Severity	5	5	10	15	20	25
	4	4	8	12	16	20
	3	3	6	9	12	15
	2	2	4	6	8	10
	1	1	2	3	4	5
		1	2	3	4	5
=Square Root (Probability X Detection)						

Note: The purpose of using the square root of the probability and detection scores is to deemphasize these two (2) components, which places more emphasize on the severity component.

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
1	Receiving Inspection / Quality control	Accept Non-Conforming Material	2	Procedural Complication	Operator Error	Process Instructions	QP-009	1	1	1	2	N/A	N/A	N/A	N/A	N/A
			1		Equipment Error	Calibration Program	QP-020	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Preventative Maintenance Program	QP-015									
		Reject Conforming Material	1	None	Operator Error	Process Instructions	QP-009	1	1	1	1	N/A	N/A	N/A	N/A	N/A
			1		Equipment Error	Calibration Program	QP-020	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Preventative Maintenance Program	QP-015									
2	Line Clearance / Area Preparation	Incorrect Material	1	Non-Functional Device	Operator Error	Process Instructions Process Settings/Line Clearance Form DMR	WI-0014 FM-0039	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Material Expired	1	Non-Functional Device	Operator Error	DMR	FM-0039	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Improper Material Handling	1	Damaged Materials	Operator Error	Process Instructions	WI-0014	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Out of Calibration Equipment	3	Non-Functional Device	Operator Error	Process Instructions	WI-0014	1	1	1	3	N/A	N/A	N/A	N/A	N/A
						Equipment Qualification	IQ Report (TR)									
		Untrained Operators	1	Compromised Device	Inadequate Training	Process Instructions	WI-0014	1	1	1	1	N/A	N/A	N/A	N/A	N/A
Training Program	QP-0017															



Process Failure Mode Effects Analysis (FMEA) – Risk Management Report

SILQ ClearTractc Foley Catheter

Doc No.: RM-0021

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 5 of 9

Table 6.0 – SILQ ClearTractc Foley Catheter

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
3	Catheter tube forming / curing	Incomplete / Malformed Tube	2	Non-Functional Device / Poor flow	Operator Error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	2	2	4	N/A	N/A	N/A	N/A	N/A
4	Balloon forming / dipping	Incomplete / Malformed Balloon	3	Non-Functional Device / Balloon Burst	Operator Error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	2	2	6	N/A	N/A	N/A	N/A	N/A
5	Valve installation	Loose valve / Valve disconnects	1	Non-Functional Device	Operator error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	1	1	1	N/A	N/A	N/A	N/A	N/A
6	Catheter Size Printing	Incorrect print information	2	Urethral damage / over inflation	Operator error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	1	1	2	N/A	N/A	N/A	N/A	N/A
		Illegible	2		Equipment	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	2	1	2	6	N/A	N/A	N/A	N/A	N/A
7	Labeling	Incorrect print information	1	Potential traceability issue	Operator error	Process instruction	MP-0016	2	1	2	2	Label inspection per IP-0002	1	1	1	1
		Illegible	1	Potential traceability issue	Equipment error	Process instruction	MP-0016	2	1	2	2	Label inspection per IP-0002	1	1	1	1



Process Failure Mode Effects Analysis (FMEA) – Risk Management Report

SILQ ClearTractc Foley Catheter

Doc No.: RM-0021

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 6 of 9

Table 6.0 – SILQ ClearTractc Foley Catheter

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
8	Surface Treatment	Incomplete / No Treatment	1	No effect	Equipment error	Process Validation	TR-HDX-0007 (MVP)	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Start-up Verification Process Instructions	MP-SLQ-0002 TM-HDX-0001 TM-HDX-0002									
						In-Process Inspection	Inspection Plans (Form) HDX-8XXX SLQ-8XXX									
					Operator error	Process instruction	MP-HDX-0004									
						In-Process Inspection	Inspection Plans (Form) HDX-8XXX SLQ-8XXX									
9	Apply Label to Pouch	Misaligned or wrinkled	1	Cosmetic issue	Operator Error	Process Instructions	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	N/A
10	Pouch Loading	Incorrect product configuration	1	Non-Functional Device	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	2	2	2	4	N/A	N/A	N/A	N/A	N/A
11	Tray Sealing (Sterile Barrier)	Weak Seal	4	Breach in sterile barrier	Equipment error	Process validation Inspection process	TR-HDX-0007 (MVP) Inspection Plans (Form) HDX-8XXX SLQ-8XXX	1	1	1	4	N/A	N/A	N/A	N/A	N/A
					Operator error – incorrect sealing parameters	Process instruction Inspection process	TR-0074 Inspection Plans (Form) HDX-8XXX SLQ-8XXX	1	1	1	4	N/A	N/A	N/A	N/A	N/A
					Contamination	Environmental controls	QP-018	1	1	1	4	N/A	N/A	N/A	N/A	N/A

Table 6.0 – SILO ClearTractc Foley Catheter

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²					
12	Assemble Box	Incorrect / Weak box	1	Damage product	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	N/A					
	Apply Label to Box	Misaligned or wrinkled	1	Cosmetic issue	Operator Error	Process Instructions	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	N/A					
13	Packaging	Incorrect packaging configuration	2	Damage product	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	2	2	2	4	Packaging inspection per IP-0002	1	1	1	2					
14	Sterilization	Product not processed / mix-up	4	Potential sterility issue	Operator error	Final Product Inspection	QP-010	1	1	1	4	N/A	N/A	N/A	N/A	N/A					
		Cycle Process Excursion	4	Potential sterility issue	Operator error	Supplier;s SOP	WI-00863 GWI-08-004-09	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1					
						Supplier Qualification Program	QP-005														
					Equipment error	Process Validation	TR 0075, TR-HDX-0009	1	1	1	4						Final Product Inspection per QP-010	1	1	1	1
						EO Program	QP-032														
		No BI's placed	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1					
						Supplier's SOP	WI-00939														
		Incorrect BI's used	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1					
						Supplier's SOP	WI-00939														
		Improper Material Handling with Lab	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1					
Supplier's SOP	WI-00939																				
Supplier Qualification Program	QP-005																				



Process Failure Mode Effects Analysis (FMEA) – Risk Management Report SILQ ClearTract Foley Catheter

Doc No.: RM-0021

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: C

Page 8 of 9

Table 6.0 – SILQ ClearTractc Foley Catheter

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
15	Final Product Inspection	Incomplete DHR	1	Potential traceability issue	Operator error	Process instruction	WI-0006 QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Incorrect DHR	1	Potential traceability issue	Operator error	Process instruction	WI-0006 QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		No CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Incomplete CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Incorrect CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A

	Process Failure Mode Effects Analysis (FMEA) – Risk Management Report SILQ ClearTract Foley Catheter							Doc No.: RM-0021
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: C
								Page 9 of 9

7.0 Risk/ Benefit Assessment

7.1 The following table indicates the risk/benefit analysis for the residual risks that were greater than the risk acceptability limit after all applicable controls/mitigations were applied:

Table 7.1 – Risk/Benefit Assessment for Residual Risk above Acceptability Limit for Process Hazard Assessment									
Table No.	Process Description	Section No.	Potential Failure Mode	Potential Effect of Failure	Potential Cause of Failure	S	F ²	Risk/Benefit Analysis	Does the benefit outweigh the residual risk?
N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A

7.2 There was no residual risk identified that was above the risk acceptability limit per §5.3.

8.0 Conclusion

8.1 Risks have been reduced as low as possible by the following methods:

- Design the product with inherent safety
- Appropriate testing to ensure design criteria were met
- Protective measures/ mechanisms to the product of the device. (critical-to-quality)
- Manufacturing and quality control measures
- Provide information for safety

Based on a review of the risk-management process for the SILQ ClearTract Foley Catheter described herein, the risk management plan based on ISO 14971:2019 has been implemented and the individual, cumulative and overall residual risk have been reduced as low as possible regardless of the magnitude and the benefit of using this device will outweigh risks posed by the device during usage.

9.0 Change History

Revision	Change Description
A	Initial release
B	Add HDX-160510, general updates, add verification references, MP-HDX-0000 to MP-HDX-0004
C	Add SLQ-1610, format updates to tables in §5.0, add sterilization process #14, update verification references

Memo

To: Letter to File K192034

From: Aaron Rogers

Date: March 2, 2022

Re: Project#00049, Modification to SILQ ClearTract Foley Catheter, 16Fr 10ml

PURPOSE

The purpose of this letter to file is to document the regulatory requirements associated with modifications to the SILQ ClearTract Foley Catheter 16 Fr 10ml (aka HDX 100% Silicone 2-Way Foley Catheter) (i.e., whether these modifications trigger the need for a 510(k), or can be implemented with supporting documentation) on behalf of SILQ Technologies Corp. Pathway LLC. is providing regulatory consult.

SCOPE

The scope of this letter is the design history file for Project#00049, SILQ ClearTract Foley Catheter Product Line & Shelf Life Extension (US), which introduces the 16Fr 10ml product model.

BACKGROUND

Following the FDA guidance for Industry and Food and Drug Administration Staff, "Deciding When to Submit a 510(k) for a Change to an Existing Device", issued on October 25, 2017, only a letter to file was deemed necessary for this design change. The following rationale is provided:

1. The change was not made with intent to significantly improve the safety or effectiveness of the device.
2. The change does not affect the existing labeling:
 - a. No changes to the Indication for Use.
 - b. No changes to the contraindication.
 - c. No changes to the warning and/or precautions.
 - d. No changes to the directions for use.
3. The change does not affect the existing technology, engineering, and/or performance in the following aspect:



Pathway, LLC.
8779 Cottonwood Ave, Suite 105
Santee, CA 92079
T +1 619 415 0103
www.pathwaynpi.com

- a. No changes to the control mechanism.
 - b. No changes to the device operating principle.
 - c. No changes to the sterilization process.
 - d. No changes to expiration dating.
 - e. No changes to the components.
 - f. No changes to the performance specifications.
4. No change was made to the material formulation.

Refer the Attachment 1 – MAIN FLOWCHART and Attachment 2 – FLOWCHART B: Technology, Engineering, and Performance Changes per FDA guidance, “Deciding When to Submit a 510(k) for a Change to an Existing Device” (October 25, 2017).

DESIGN CHANGE

A change was made to the technology, engineering, and/or performance with respect to the catheter component:

- The balloon volume changed from a range of 5ml – 10ml to 10ml. This new version was released as a new product model on ECO-000889.

The design changes were captured through the design control process per QP-007 of the Quality System. Activities and requirements for implementation of this design change was documented on the design change assessment per Project#00049.

CONCLUSION

As described above, and in accordance to the FDA guidance document, the device modification does not alter the safety or effectiveness of the product within scope. Therefore, this review has determined that the design change may be implemented with documentation to the Design History File.

ATTACHMENT

- Attachment 1 – MAIN FLOWCHART
- Attachment 3 – FLOWCHART B: Technology, Engineering, and Performance Changes

A handwritten signature in purple ink, reading "Aaron Rogers".

Aaron Rogers
Director of Regulatory and Quality



Pathway, LLC.
8779 Cottonwood Ave, Suite 105
Santee, CA 92079
T +1 619 415 0103
[www. pathwaynpi.com](http://www.pathwaynpi.com)

ATTACHMENT 1

Contains Nonbinding Recommendations

Note that the first question is whether the change is being made with the intent to significantly improve the safety or effectiveness of the device, for example, to significantly improve clinical outcomes, to mitigate a known risk, in response to adverse events, etc. (Figure 1 – Main Flowchart). If so, the change likely could significantly affect safety or effectiveness and submission of a new 510(k) is likely required. If not, you should continue to follow the logic scheme shown in Figure 1, below.

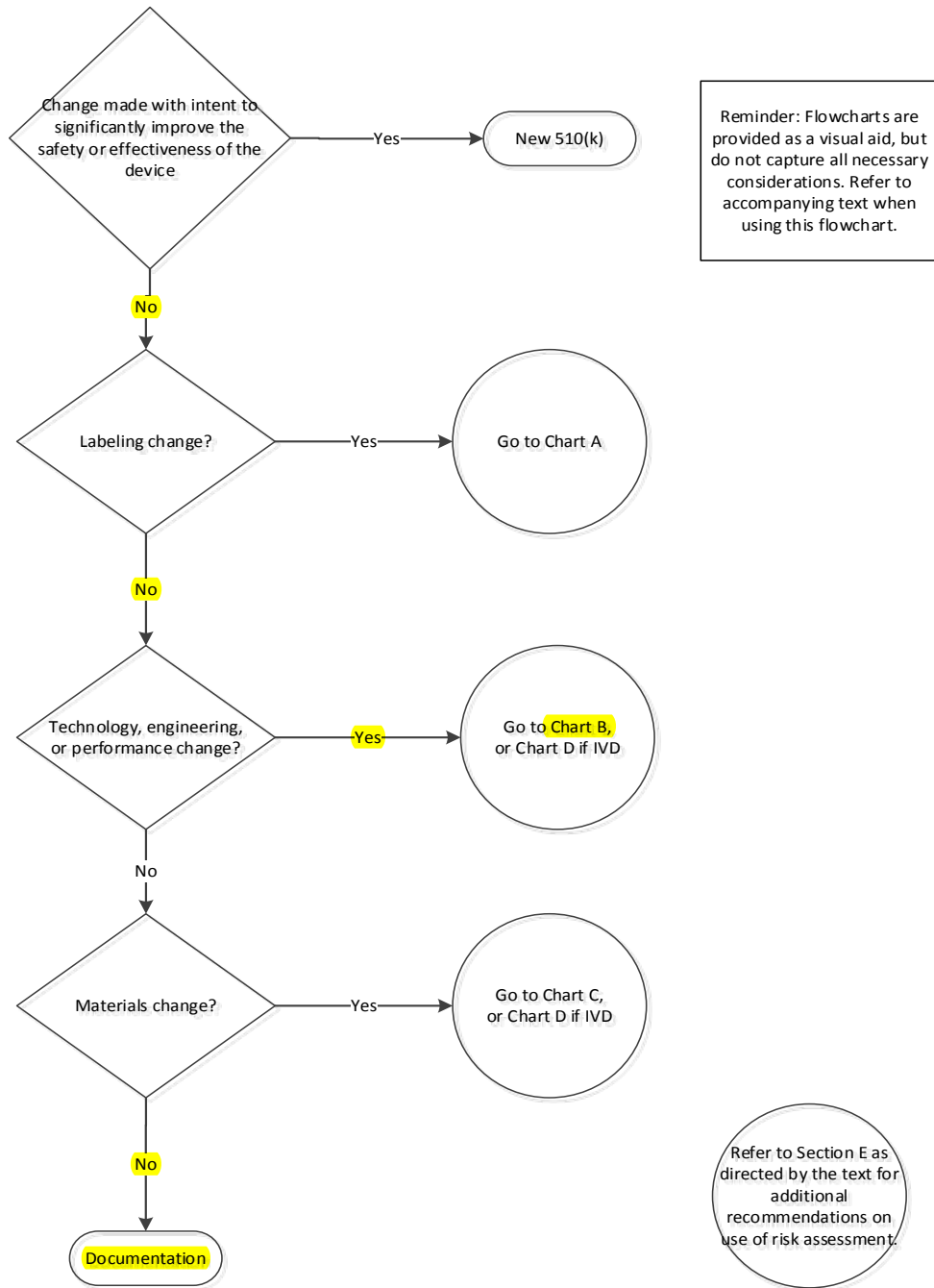


Figure 1 - Main Flowchart



Pathway, LLC.
8779 Cottonwood Ave, Suite 105
Santee, CA 92079
T +1 619 415 0103
[www. pathwaynpi.com](http://www.pathwaynpi.com)

ATTACHMENT 2

Contains Nonbinding Recommendations

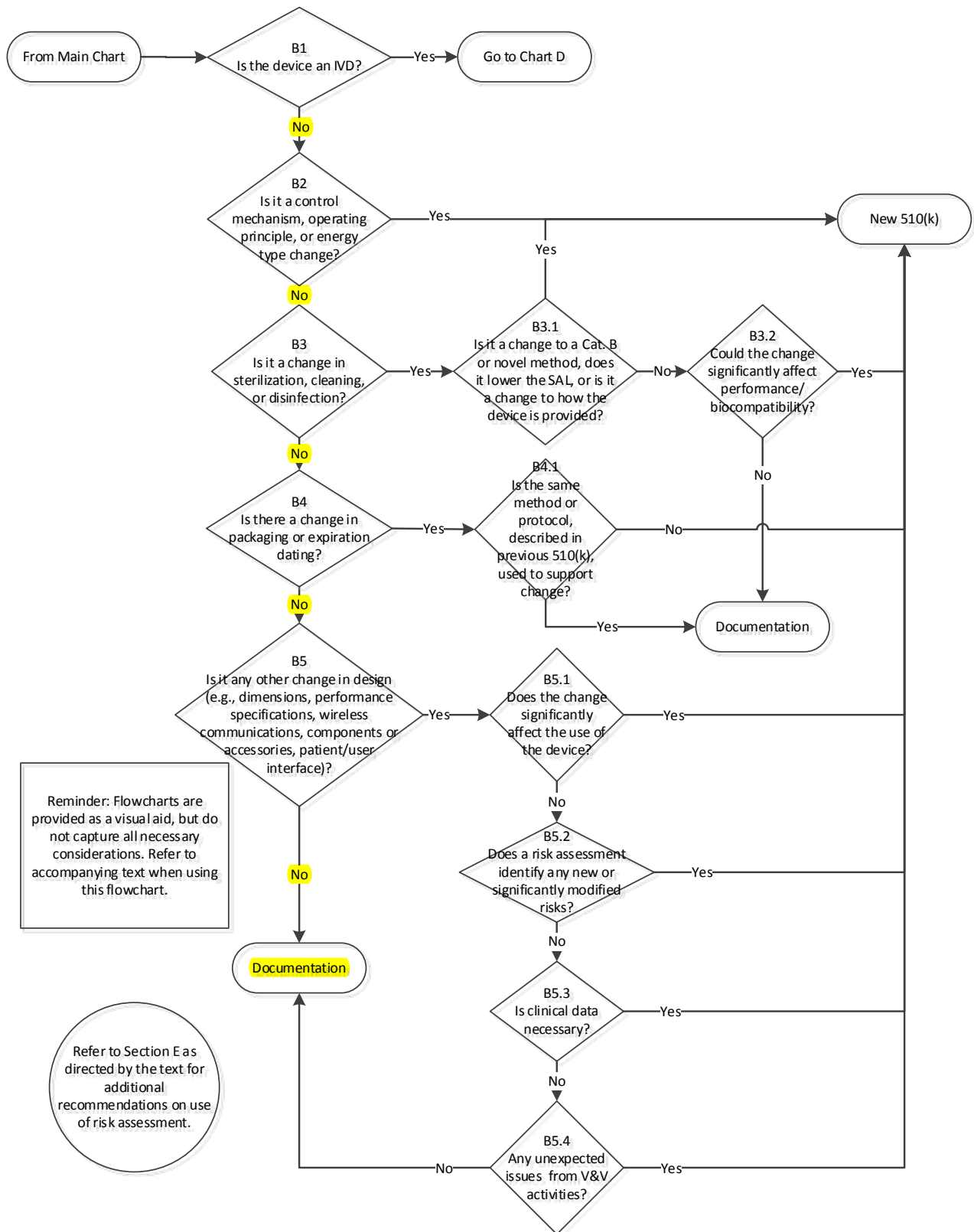


Figure 3 - Flowchart B: Technology, Engineering, and Performance Changes