

Pathway Innovation LLC.	Template, Design Change Assessment		Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: A
			Page 1 of 2

Project No. / Name	00049 / SILQ 16Fr X 10 ml Product Introduction (US)				
Product Design Specification / Rev	PS-0006 / C			Date	03-01-2022
Change Type	Product: ☒ Line Extension ☐ Change	Manufacturing: ☐ Change	Rationale	1) This is a product line extension by adding a version of the catheter with a 10 ml balloon. This line extension is not a material change but a change to the volume labeled on the balloon of 10 ml. The following catalog number will be added: SLQ-1610	

Rationale = Describe the justification for the change type and explain the drivers behind the change.

The following table provides an analysis of the potential impact proposed design changes may have with the design requirements described in affected product design specification.

TABLE 1: Change and Requirement Identification						
ID	Type	Component	Change Description	Qualified Part / Rev	Proposed Part / Rev	Requirements Affected
1	A/B/C	Catheter	Balloon volume changed from a 5 ml – 10 ml range to 10 ml. No material changes.	HDX-6006 / A	SLQ-6004 / A	#2,3,4 - Labeling #11 - Balloon Integrity #12 - Balloon Volume Maintenance #13 - Non interconnect #14 - Deflation reliability #15 - Flow rate #16 - Balloon Traction #18 - Balloon & Shaft Size (DIMS) #21 - Connection type #22 - No of Lumens #23 - French size offerings #24 – Shelf Life
2	A/C	Box / Label	Change from 16Fr 5ml – 10 ml to 16F 10 ml	HDX-4002 / A	SLQ-4000 / A	#2,3,4 – Labeling
3	A/C	IFU	Change from 16Fr 5ml – 10 ml to 16F 10 ml	HDX-4006 / A	SLQ-4006 / A	#2,3,4 – Labeling

Type: A = Form, B = Fit, C = Function, Performance

The following table provides an impact assessment for the proposed design changes described in Table 1.

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Internal / External Resources	None
Products Affected	SLQ-1610
Customers Affected	US Market
Distribution Channels	US Market
Physical/ Functional (e.g., Bench Testing) Assessment	Design verification was performed on 16Fr 5 ml -10 ml catheter which is equivalent to the catheter labeled as 16Fr 10 ml. There is no impact to the design verification performed because the proposed model performs within the tested specification. The product specification (PS-0006) will be revised to include this model in the trace matrix section. No further action is required.

Pathway Innovation LLC.	Template, Design Change Assessment		Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: A
			Page 2 of 2

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Shelf-Life Assessment	All materials that construct the 2-way Catheter have been qualified with a twelve (12) month shelf life. No further action is required.
Dimensional (e.g., 1 st Article) Assessment	Device proximal tip, balloon and shaft was confirmed on 16Fr 5 ml and 16Fr 5 ml -10 ml catheter which is equivalent to the catheter labeled as 16Fr 10 ml. There is no impact to the Device proximal tip, balloon and shaft because the proposed model performs within the tested specification. The product specification (PS-0006) will be revised to include this model in the trace matrix section. No further action is required.
Sterilization Assessment	There is no impact to the sterilization validation performed because there were no sterile barrier or pallet configuration changes that affect sterility. No further action is required.
Biocompatibility Assessment	There is no impact to the biocompatibility evaluations performed because there were no material changes. No further action is required.
Packaging Assessment	There is no impact to the validated packaging configuration. No further action is required.
Regulatory Assessment	The guidance document “Deciding when to submit a 510(k) for a change to an existing device (1500054)” was used in this assessment, which the attached decision trees were completed. A 510(k) submission is not required. No further action is required.
Quality System Assessment	None

Contains Nonbinding Recommendations

Note that the first question is whether the change is being made with the intent to significantly improve the safety or effectiveness of the device, for example, to significantly improve clinical outcomes, to mitigate a known risk, in response to adverse events, etc. (Figure 1 – Main Flowchart). If so, the change likely could significantly affect safety or effectiveness and submission of a new 510(k) is likely required. If not, you should continue to follow the logic scheme shown in Figure 1, below.

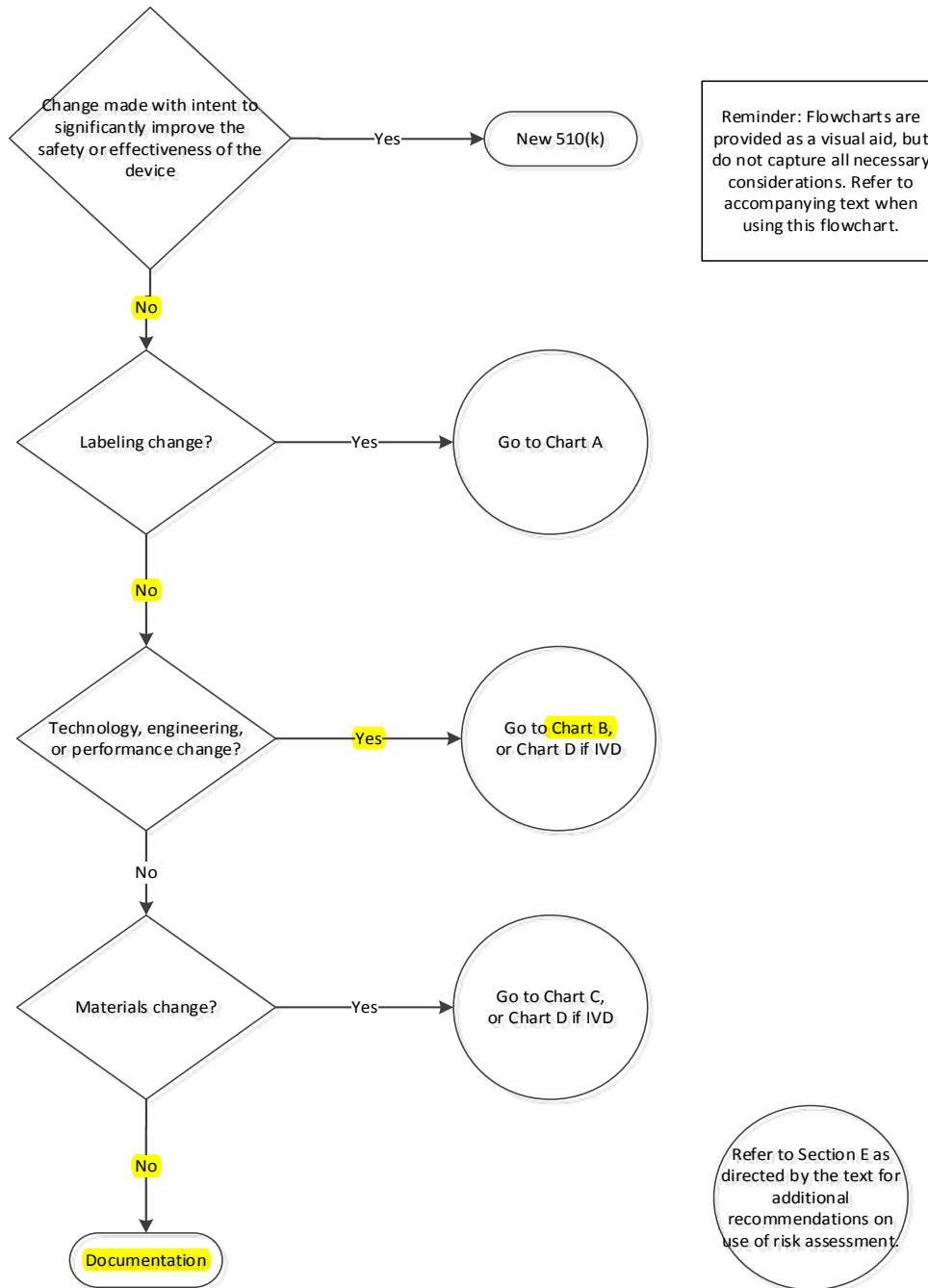


Figure 1 - Main Flowchart

Contains Nonbinding Recommendations

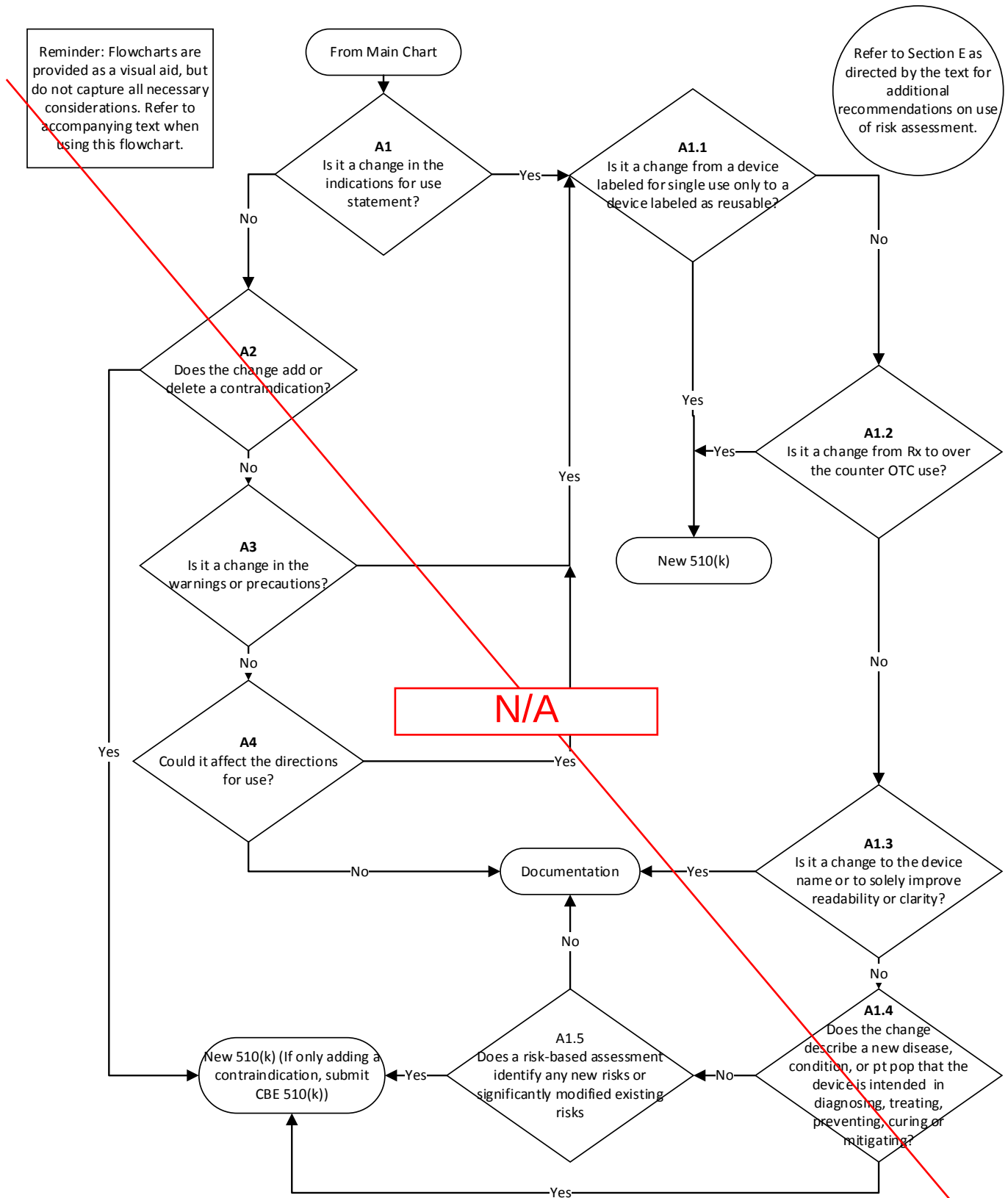


Figure 2 - Flowchart A: Labeling Changes

Contains Nonbinding Recommendations

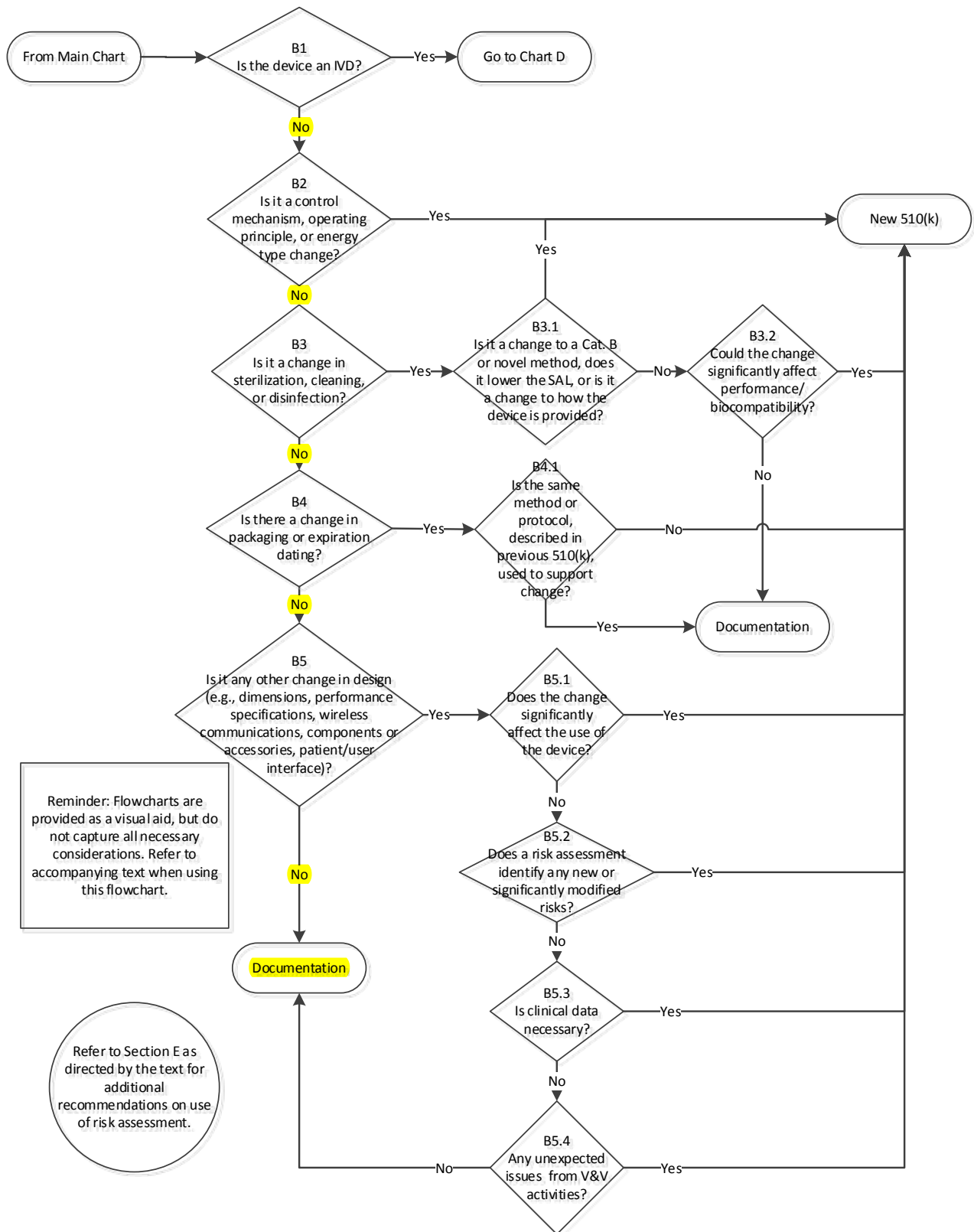


Figure 3 - Flowchart B: Technology, Engineering, and Performance Changes

Contains Nonbinding Recommendations

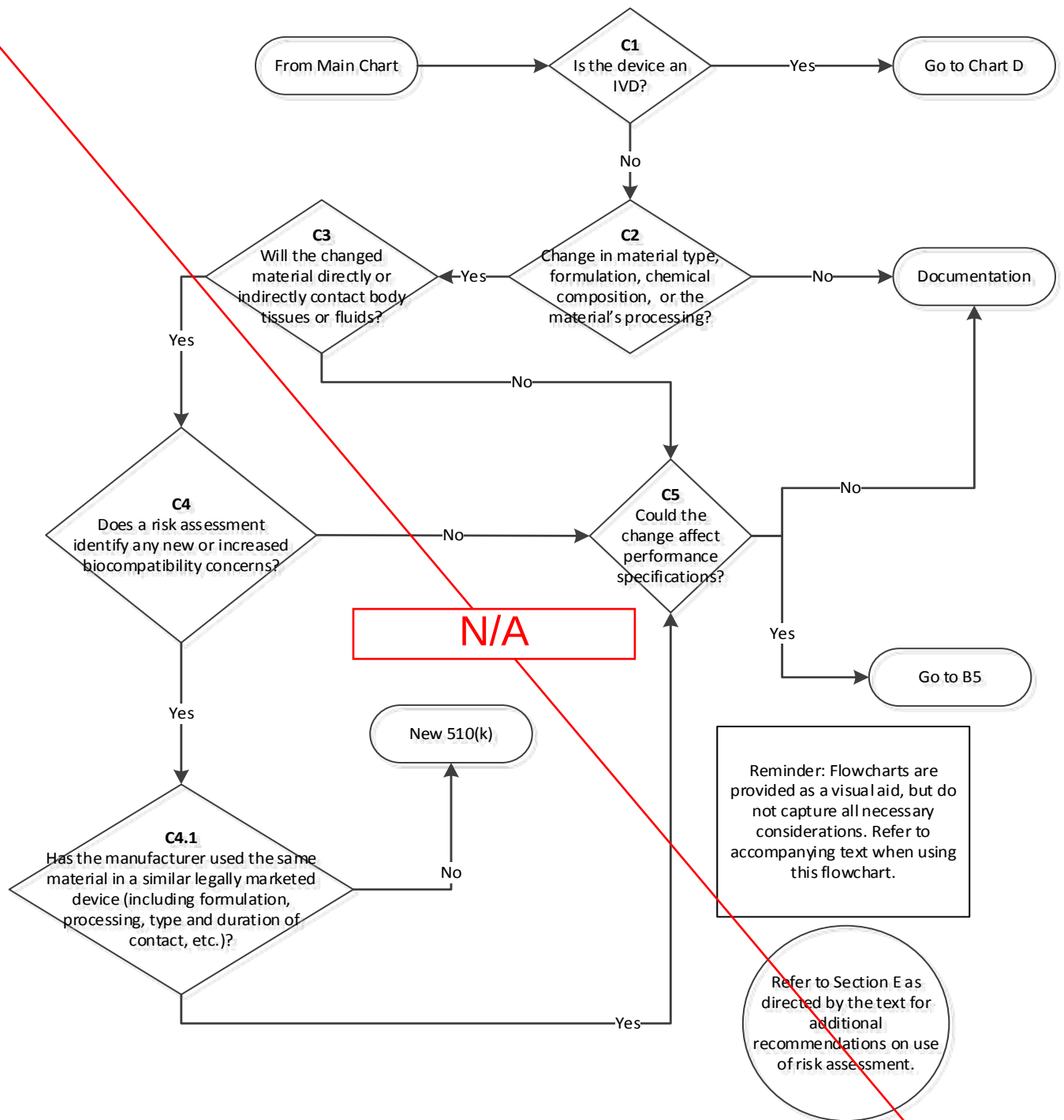


Figure 4 - Flowchart C: Materials Changes