

#1 CHANGE ASSESSMENT

Pathway Innovation LLC.	Template, Design Change Assessment		Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: A
			Page 1 of 2

Project No. / Name	00053 / SILQ ClearTract Foley Catheter Suprapubic Indications for Use				
Product Design Specification / Rev	PS-0006 / D			Date	05-11-2023
Change Type	Product: ☒ Line Extension ☐ Change	Manufacturing: ☐ Change	Rationale	This is a line extension to create new product models with expanded indications for use that includes suprapubic placement of the catheter. New SPT models will include 14Fr10ml, 16Fr5ml, 16Fr10ml, 18Fr10ml.	

Rationale = Describe the justification for the change type and explain the drivers behind the change.

The following table provides an analysis of the potential impact proposed design changes may have with the design requirements described in affected product design specification.

TABLE 1: Change and Requirement Identification						
ID	Type	Component	Change Description	Qualified Part / Rev	Proposed Part / Rev	Requirements Affected
1	A/B/C	IFU	Intended Use changes from: Urological catheter is intended for drainage of the urinary tract. Catheterization is accomplished by inserting the catheter through the urethra and into the bladder. Intended population is adults and pediatrics Intended Use changes to: Urological catheter is intended for drainage of the urinary tract. Catheterization is accomplished by inserting the catheter through the urethra and into the bladder, but also suprapubic placement or by nephrostomy. Intended population is adults and pediatrics	SLQ-4006 / B	SLQ-4007 / A	#1 – Biocompatibility: New endpoint identified (Genotoxicity) #2,3,4 - Labeling
2	A/B/C	Labels / IFU	Update Manufacturer address from: Manufactured for: SILQ Technologies Corp. 12100 Wilshire Blvd, Suite 800 Los Angeles, CA 90025 949-735-2614 Update Manufacturer address to: Manufactured by: SILQ Technologies Corp. 323 Sunny Isles Blvd., 7th Floor Sunny Isles Beach, FL 33160 (424) 309-8523	SLQ-4000 / A SLQ-4001 / A SLQ-4006 / A	SLQ-LBL-4008 / A SLQ-LBL-4009 / A SLQ-4007 / A	#2,3,4 – Labeling

Type: A = Form, B = Fit, C = Function, Performance

The following table provides an impact assessment for the proposed design changes described in Table 1.

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Internal / External Resources	NAMSA – biocompatibility testing
Products Affected	SLQ-211410SPT, SLQ-211605SPT, SLQ-211610SPT, SLQ-211810SPT
Customers Affected	US Market
Distribution Channels	US Market

Pathway Innovation LLC.	Template, Design Change Assessment	Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>	Rev.: A
		Page 2 of 2

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Physical/ Functional (e.g., Bench Testing) Assessment	There are no physical or functional changes to the product. PS-0006 will be updated with reference to design output documents for new model configurations.
Shelf-Life Assessment	All materials that construct the 2-way Catheter have been qualified with a twelve (12) month shelf life. No further action is required.
Dimensional (e.g., 1 st Article) Assessment	There are no dimensional changes to the product. No dimensional verification is required. No further action is required.
Sterilization Assessment	There is no impact to the sterilization validation performed because there were no sterile barrier or pallet configuration changes that affect sterility. No further action is required.
Biocompatibility Assessment	Expanding the indications for use will require genotoxicity testing. The biocompatibility summary report will be amended with additional testing results.
Packaging Assessment	There is no impact to the validated packaging configuration. No further action is required.
Regulatory Assessment	The guidance document “Deciding when to submit a 510(k) for a change to an existing device (1500054)” was used in this assessment, which the attached decision trees were completed. A 510(k) submission is required.
Quality System Assessment	No changes required to the QMS. No further action is required.
Risk Assessment	New risks for modified indications for use and route of administration to be documented and addressed by SILQ in their RMF. SILQ RMF will be updated with reference to new verification activities. RMF is the responsibility of SILQ. No further action is required by Pathway.
Manufacturing Assessment	No changes are required to the manufacturing process. No further action is required.

Contains Nonbinding Recommendations

Note that the first question is whether the change is being made with the intent to significantly improve the safety or effectiveness of the device, for example, to significantly improve clinical outcomes, to mitigate a known risk, in response to adverse events, etc. (Figure 1 – Main Flowchart). If so, the change likely could significantly affect safety or effectiveness and submission of a new 510(k) is likely required. If not, you should continue to follow the logic scheme shown in Figure 1, below.

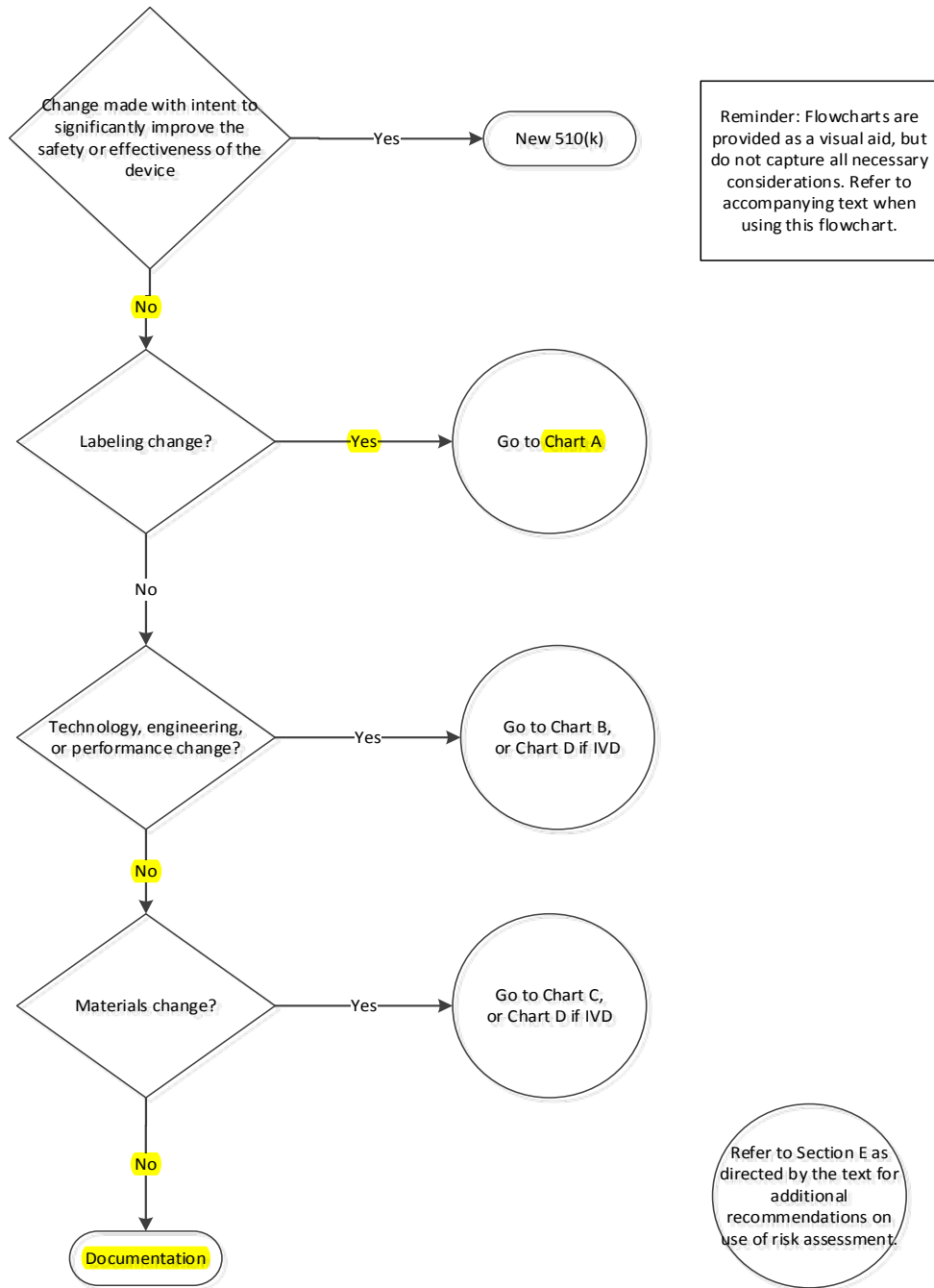


Figure 1 - Main Flowchart

Contains Nonbinding Recommendations

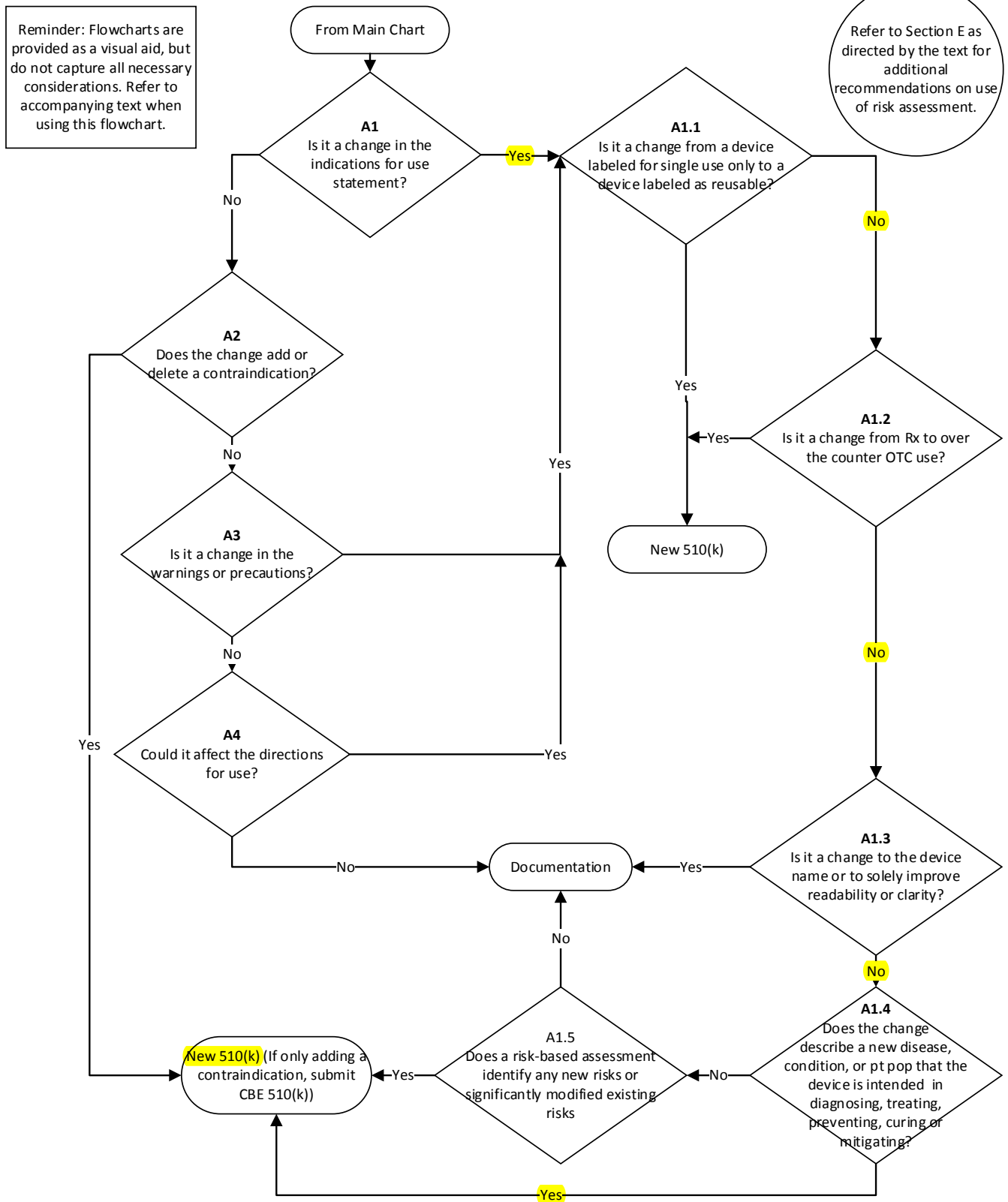


Figure 2 - Flowchart A: Labeling Changes

Contains Nonbinding Recommendations

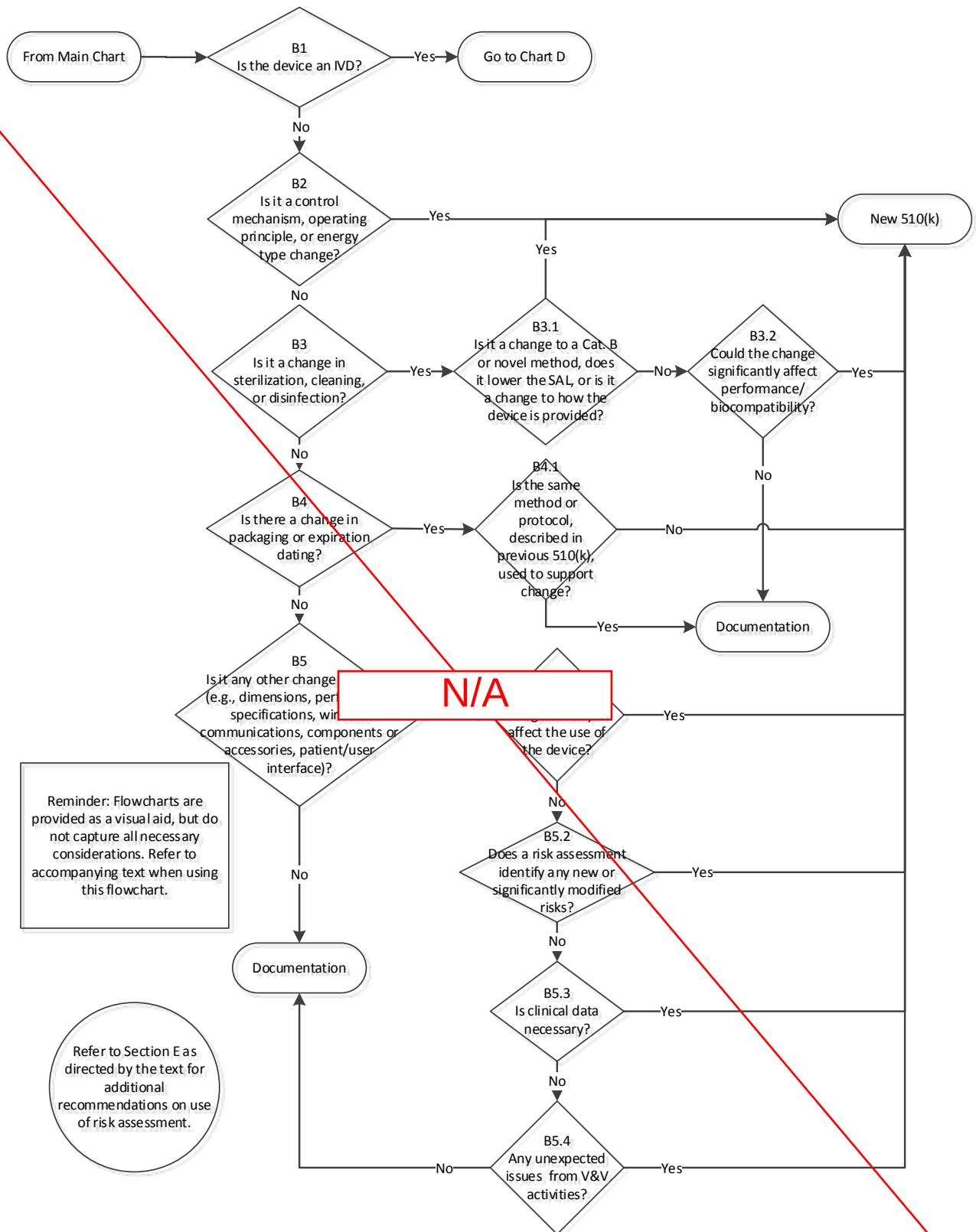


Figure 3 - Flowchart B: Technology, Engineering, and Performance Changes

Contains Nonbinding Recommendations

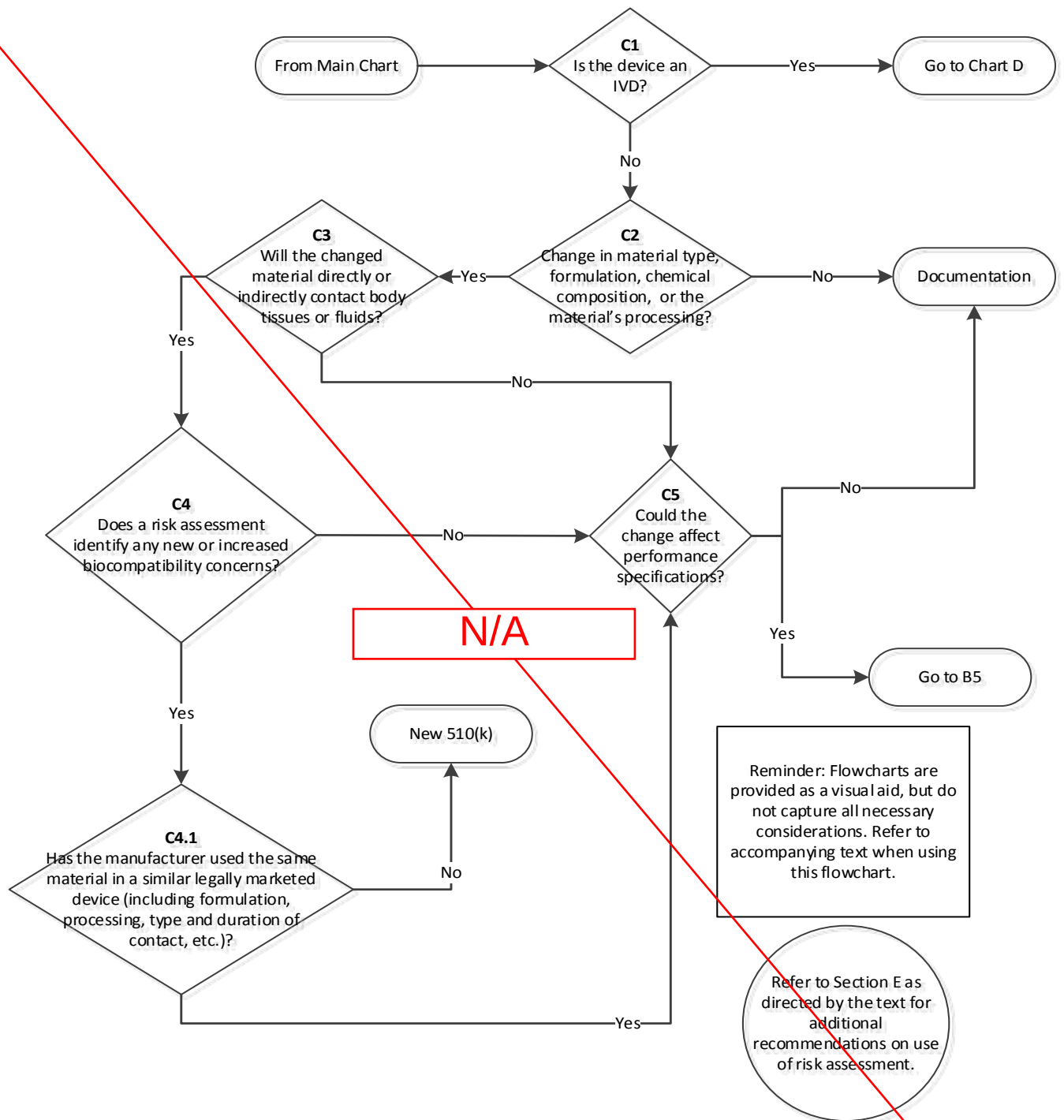


Figure 4 - Flowchart C: Materials Changes

#2 BIOCOMPATIBILITY REPORT

	Test Report, Biocompatibility, SILQ ClearTract Catheter		Doc No.: TR-SLQ-0012
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 1 of 7

1.0 Purpose

- 1.1 This report summarizes the additional results of the biocompatibility testing performed on the SILQ ClearTract Catheter which is to support an expanded indications for use and changes to the route of administration for suprapubic placement.
- 1.2 Controls associated with risks per §1.1 are summarized below

Table 1: Risk Assessment for expanded indications for use			
Device Change	Risk	Verification / Validation Method(s)	
Indications for Use change to add suprapubic placement	Immunological / Adverse tissue response	Biocompatibility evaluation in agreement with recommendations in CDRH's Biocompatibility Guidance. Based on our risk management procedures, biocompatibility testing was performed for new endpoint. Nature of body contact changes from Surface Device / Mucosal Membrane to External Communicating / Tissue.	Genotoxicity (ISO 10993-3:2014) using: Mouse Lymphoma Assay - 4 and 24 Hr Treatments, 2 extracts Ames Bacterial Reverse Mutation Assay - 2 extracts, 5 Strains
Route of administration change to include suprapubic placement	Immunological / Pyrogenic Response / Fever	Sterility evaluation in agreement with recommendations in CDRH's Sterility Information Guidance. Based on our risk management procedures, bacterial endotoxin testing was performed.	USP <161> Medical Devices – Bacterial Endotoxin and Pyrogen Tests and USP <85> Bacterial Endotoxins Test: by Kinetic Turbidimetric Assay

- 1.3 The following table summarizes the design requirement fulfilled by the biocompatibility testing provided in this report.

Table 2: Design Requirements to Design Verification Activities		
Requirement # / Type	Product Requirement Description	Verification Activity
1.1 / Biological	Patient contacting surfaces and materials shall meet biocompatibility requirements for an external communicating device with tissue/bone/dentin for prolonged contact duration (>24hr and <30days) per ISO 10993-1.	Biocompatibility Testing per ISO 10993-1

2.0 Scope

- 2.1 Applies to the SILQ ClearTract Catheters as defined in Project#00053 and PS-0006.

	Test Report, Biocompatibility, SILQ ClearTract Catheter		Doc No.: TR-SLQ-0012
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 2 of 7

- 2.2 Biocompatibility testing was performed by NAMSA located at 6750 Wales Road, Northwood, Ohio 43619.
- 2.3 Biocompatibility testing was performed by LGGS, Inc. located at 3335 E. Miraloma Avenue, Suite 140, Anaheim, CA 92806.
- 2.4 Based on the original clinical application, the device was categorized as a surface-contacting/ mucosal membrane with a prolong contact duration (> 24 hours to 30 days) per ISO 10993-1. The biocompatibility test program completed and reported in TR-0071 for device category was:
 - Cytotoxicity
 - Sensitization
 - Irritation or Intracutaneous Reactivity
 - Systemic Injection
 - Material mediated pyrogen
 - Intramuscular Implantation
 - Subacute Intraperitoneal systemic toxicity
 - Sub chronic Intravenous systemic toxicity
- 2.5 Based on the present clinical application, the device is categorized as an external communicating device with tissue/bone/dentin for prolonged contact duration (>24hr to 30 days) per ISO 10993-1. In agreement with recommendations in CDRH's Biocompatibility Guidance and CDRH's current guidance for *Submission and Review of Sterility Information in Premarket Notification (510(k)) Submissions for Devices Labeled as Sterile, January 21, 2016*, the following additional testing was performed and results provided in this report.
 - Genotoxicity Testing
 - Bacterial Endotoxin Testing

3.0 References

- 3.1 Project#00053 – SILQ ClearTract Foley Catheter Suprapubic Indications for Use
- 3.2 PS-0006 rev E – Product Design Specification and Trace Matrix, SILQ ClearTract Catheter
- 3.3 QP-002 – Document Control
- 3.4 Test Methods / Standards
 - 3.4.1 ISO 10993-3:2014 – Biological evaluation of medical devices - Part 3: Tests for genotoxicity, carcinogenicity and reproductive toxicity (FDA 2-228)
 - 3.4.2 ANSI AAMI ST72:2019 – Bacterial endotoxins - Test methods, routine monitoring, and alternatives to batch testing (FDA 14-541)
 - 3.4.3 USP <161> Medical Devices-Bacterial Endotoxin and Pyrogen Tests (FDA 14-564)
 - 3.4.4 USP <85> Bacterial Endotoxin Tests (FDA 14-570)
- 3.5 TR-0071 – Test Report, Biocompatibility, HDX Foley Catheter

Note: The product name SILQ ClearTract Catheter and Hydrophilix (HDX) Foley Catheter is used interchangeably. The company name Hydrophilix and SILQ are used interchangeably as the name was changed from Hydrophilix to SILQ.

	Test Report, Biocompatibility, SILQ ClearTract Catheter		Doc No.: TR-SLQ-0012
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 3 of 7

4.0 Testing Summary

Biocompatibility evaluations were conducted by NAMSA (Northwood, OH) in accordance with the provisions of 21 CFR 58, Good Laboratory Practice for Nonclinical Laboratory Studies. Test samples consisted of appropriate surface area/ mass of materials, as determined by the contract test laboratory. Test samples were fully manufactured and sterilized prior to being forwarded to the contract test laboratory. Samples were prepared using fully intact devices to ensure all tissue contacting surfaces were evaluated per the relevant ISO standard.

The outcomes of the biocompatibility evaluations, test methods, and results are listed in **Table 3**. All test results satisfied the biocompatibility requirements of the relevant ISO-10993-3 standard for genotoxicity.

Bacterial endotoxin testing (BET) was conducted by LGGS, Inc. (Anaheim, CA) in accordance with 21 CFR Part 820 Good Manufacturing Practice, as applicable to a testing laboratory. Test samples were fully manufactured devices and terminally sterilized via Ethylene Oxide (EO) prior to being forwarded to the contract test laboratory.

Table 3: Biocompatibility Testing Summary			
Test / Endpoints	Description (Method)	Result	Lab Report
Genotoxicity	Genotoxicity (ISO 10993-3:2014) using: Mouse Lymphoma Assay - 4 and 24 Hr Treatments, 2 extracts Ames Bacterial Reverse Mutation Assay - 2 extracts, 5 Strains	Mouse Lymphoma Assay: Under the conditions of this study, the test article was considered negative (non-mutagenic, nonclastogenic) because the test article induced mutant frequency (IMF) did not exceed the Global Evaluation Factor of 90×10^{-6} at Relative Total Growth (RTG) values $\geq 10\%$. Study Outcome: Pass (Non-mutagenic) Ames Bacterial Reverse Mutation Assay: A healthy lawn with no apparent thinning or increase in the size of micro-colonies was observed for all test article plates; therefore these background lawns were considered normal and were given a score of '1'. No precipitate was noted on the test article plates for each strain. Under the conditions of this study and for all strains tested, the mean number of revertant colonies for the test article was less than 2-fold over the mean number of revertant colonies for the concurrent negative control, the number of revertant colonies was within spontaneous negative control reversion ranges or Test Facility historical ranges, and the background lawn was considered normal; therefore, the test article was considered Non-Mutagenic. Study Outcome: Pass (Non-mutagenic)	Attachment 11.1
Bacterial Endotoxins	USP <161> Medical Devices – Bacterial Endotoxin and Pyrogen	N=10 / LOT Lot#1: 12.4 EU/device; 0.31 EU/ml	Attachment 11.2

	Test Report, Biocompatibility, SILQ ClearTract Catheter		Doc No.: TR-SLQ-0012
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 4 of 7

Table 3: Biocompatibility Testing Summary			
Test / Endpoints	Description (Method)	Result	Lab Report
	Tests and USP <85> Bacterial Endotoxins Test: by Kinetic Turbidimetric Assay	Lot#2: 12.4 EU/device; 0.31 EU/ml Lot#3: 13.2 EU/device; 0.33 EU/ml	

5.0 Test Plan, Biocompatibility (ISO 10993-1:2018)

- 5.1 Biocompatibility testing was performed by NAMSA located at 6750 Wales Road, Northwood, Ohio 43619.
- 5.2 LAL/BET testing was performed by LGGS, Inc. located at 3335 E. Miraloma Avenue, Suite 140, Anaheim, CA 92806.
- 5.3 Tissue contacting materials according to ISO 10993-1 are summarized in the following table:

Table 4: Tissue Contacting Materials				
Component	Material	Category	Contact	Contact Duration
Catheter: Funnel / Shaft / Balloon, Tip	Medical Grade Silicone YM827 Barium Sulfate	Externally communicating	Tissue / Bone / Dentin	Prolonged (> 24 hours to 30 days)
SILQ Surface Treatment	Zwitterionic polymer			

- 5.4 The biological endpoints and test methods are summarized in the following table:

Table 5: Biocompatibility Testing Summary	
Test / Endpoints	Description (Method)
Genotoxicity	Genotoxicity (ISO 10993-3:2014) using: Mouse Lymphoma Assay - 4 and 24 Hr Treatments, 2 extracts Ames Bacterial Reverse Mutation Assay - 2 extracts, 5 Strains
Bacterial Endotoxins	USP <161> Medical Devices – Bacterial Endotoxin and Pyrogen Tests and USP <85> Bacterial Endotoxins Test: by Kinetic Turbidimetric Assay

5.5 Test Device

The test device consisted of the SILQ Catheter.

6.0 Sampling Requirements

- 6.1 Sampling requirements for each biocompatibility evaluation are summarized in the following table:

Table 6: Biocompatibility Testing Sample Requirements		
Test / Endpoints	Description (Method)	Sample Size
Genotoxicity	Genotoxicity (ISO 10993-3:2014) using: Mouse Lymphoma Assay - 4 and 24 Hr Treatments, 2 extracts	Surface Area: 177.9 cm ²

	Test Report, Biocompatibility, SILQ ClearTract Catheter		Doc No.: TR-SLQ-0012
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 5 of 7

Table 6: Biocompatibility Testing Sample Requirements		
Test / Endpoints	Description (Method)	Sample Size
	Ames Bacterial Reverse Mutation Assay - 2 extracts, 5 Strains	N=1
Bacterial Endotoxins	USP <161> Medical Devices – Bacterial Endotoxin and Pyrogen Tests and USP <85> Bacterial Endotoxins Test: by Kinetic Turbidimetric Assay	N=10 / LOT

7.0 Acceptance Criteria

7.1 Acceptance criteria for each biocompatibility evaluation are defined per ISO 10993 standards and summarized below:

Table 7: Acceptance Criteria	
Evaluation Description	Acceptance Criteria
Genotoxicity	Mouse Lymphoma Assay: The test article will be considered negative (non-mutagenic / non-clastogenic) if the IMF does not exceed the Global Evaluation Factor of 90×10^{-6} at Relative Total Growth (RTG) values $\geq 10\%$. The test article will be considered positive (mutagenic) if the IMF exceeds the Global Evaluation Factor of 90×10^{-6} at Relative Total Growth (RTG) values $\geq 10\%$ for any treatment condition. Ames Bacterial Reverse Mutation Assay: The test article will be considered non-mutagenic in a single dose study if the following are true for all test strains: The calculated fold increase over the negative control is less than 2-fold The number of revertant colonies is within or slightly outside spontaneous negative control reversion ranges and/or the Test Facility historical ranges The background lawn is considered normal or slightly reduced (score of 1 or 2) The test article will be considered mutagenic in a single dose study if the following are true for at least one test strain: The calculated fold increase over the negative control is 2-fold or greater The number of revertant colonies is outside spontaneous negative control reversion ranges The background lawn is considered normal or slightly reduced (score of 1 or 2)
Bacterial Endotoxins	≤ 20 endotoxin units (EU)/Device or ≤ 0.5 endotoxin units (EU)/ml

8.0 Deviations

8.1 There were no deviations.

9.0 Results & Discussion

9.1 The test articles were fully manufactured. Then sent to Steris for EO sterilization per the validated cycle.

9.2 Test articles were shipped to NAMSA (Northwood, Ohio) for the biocompatibility evaluations. The evaluations were conducted under GLP conditions.

	Test Report, Biocompatibility, SILQ ClearTract Catheter	Doc No.: TR-SLQ-0012
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 6 of 7

9.3 The following are the summarized conclusions from each biocompatibility evaluation.

Table 8: – Biocompatibility Testing Summary		
Evaluation Description	Result	Meets Acceptance Criteria
Genotoxicity	Mouse Lymphoma Assay: Under the conditions of this study, the test article was considered negative (non-mutagenic, nonclastogenic) because the test article induced mutant frequency (IMF) did not exceed the Global Evaluation Factor of 90×10^{-6} at Relative Total Growth (RTG) values $\geq 10\%$. Study Outcome: Pass (Non-mutagenic) Ames Bacterial Reverse Mutation Assay: A healthy lawn with no apparent thinning or increase in the size of micro-colonies was observed for all test article plates; therefore these background lawns were considered normal and were given a score of '1'. No precipitate was noted on the test article plates for each strain. Under the conditions of this study and for all strains tested, the mean number of revertant colonies for the test article was less than 2-fold over the mean number of revertant colonies for the concurrent negative control, the number of revertant colonies was within spontaneous negative control reversion ranges or Test Facility historical ranges, and the background lawn was considered normal; therefore, the test article was considered Non-Mutagenic. Study Outcome: Pass (Non-mutagenic)	Yes
Bacterial Endotoxins	N=10 / LOT Lot#1: 12.4 EU/device; 0.31 EU/ml Lot#2: 12.4 EU/device; 0.31 EU/ml Lot#3: 13.2 EU/device; 0.33 EU/ml	Yes

10.0 Conclusion

10.1 The ClearTract Catheter met all biological evaluation requirements indicated in the report, which is in-line with an externally communicating, tissue/bone contacting device with a prolonged contact duration (> 24 hours to 30 days).

11.0 Attachments

- 11.1 NAMS GLP Report No. KMV002-GE12 – Ames Bacterial Reverse Mutation Assay (ISO) – Final Report Mouse
- 11.2 NAMS GLP Report No. KMV001-GE01 –Lymphoma Assay (ISO) – Final Report
- 11.3 LGGS Technical Report Lab No 22-2318-00 – LAL Kinetic Turbidimetric Assay
- 11.4 LGGS Technical Report Lab No 22-2319-00 – LAL Kinetic Turbidimetric Assay
- 11.5 LGGS Technical Report Lab No 22-2320-00 – LAL Kinetic Turbidimetric Assay

	Test Report, Biocompatibility, SILQ ClearTract Catheter		Doc No.: TR-SLQ-0012
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: A
			Page 7 of 7

12.0 Change History

Revision	Change Description
A	Initial release

#3 UPDATED LABELING

Item Number	Item Name
SLQ-4007 rev A	Artwork, IFU SILQ ClearTract SPT
SLQ-LB-4008 rev A	Artwork, Pouch Label
SLQ-LB-4009 rev A	Artwork, Box Label
SLQ-LB-4010 rev A	Artwork, Shipper Box Label



2-Way 100% Silicone ClearTract SPT Catheter

Please read all instructions before using this device

Not made with natural rubber latex

Foley Catheter is MR Safe.

Model	Size
211410SPT	14 FR 10 mL
211605SPT	16 FR 5 mL
211610SPT	16 FR 10 mL
211810SPT	18 FR 10 mL

Contraindications

This device is contraindicated in the presence of conditions which create unacceptable risk during catheterization.

Device Description

The 2-Way 100% Silicone ClearTract Catheter is a standard single-use 2-way Foley catheter that is constructed of medical grade silicone with a surface modification. It incorporates two (2) lumens, one for inflation/deflation of the balloon and the other for drainage of the urinary tract.

The surface of the ClearTract Catheter has a surface modification that reduces surface friction.

Intended Use

The 2-Way 100% Silicone ClearTract Catheter is intended for drainage of the urinary tract. Catheterization is generally accomplished by inserting the catheter through the urethra and into the bladder. However, drainage is sometimes accomplished by suprapubic placement or other placement of the catheter, such as nephrostomy tract. Intended population is adults and pediatrics.

Directions

Installation & Inflation

1. When placement of a urethral catheter is contraindicated or unsuccessful, percutaneous suprapubic urinary bladder catheterization or nephrostomy can be used in accordance with the procedures and policies of the hospital or healthcare facility.
2. Prepare application site in a sterile fashion.
3. Visually inspect the catheter for any signs of damage or deterioration prior to use.
4. If a trial inflation is conducted, follow steps 7 & 8 on the sterile field then evacuate the sterile water with the syringe. Ensure the balloon is empty prior to insertion.
5. Apply medical grade lubricant to distal end of the catheter.
6. Gently advance the catheter through the urethra until distal tip is beyond the bladder neck.
7. Confirm the position of the distal tip. The return of urine typically indicates the distal tip is in the bladder.
8. Fill Luer-tipped syringe with sterile water and connect to the inflation valve.
9. Fill the catheter balloon to the labeled volume. Do NOT exceed labeled volume.
10. Slowly withdraw the catheter until resistance is felt. This is a sign the balloon is resting at the bladder neck.



Deflation & Removal

1. Gently connect Luer-tipped syringe to the inflation valve. Gently aspirate while allowing the pressure in the balloon to assist filling the syringe with water occupying the balloon. Remove and reconnect the syringe in the inflation valve if there is no or slow deflation. Rapid deflation could collapse the lumen and prevent adequate deflation of the balloon.
2. If during the removal process the balloon ruptures, ensure all fragments are carefully removed from patient.
3. Discard after use.

Precautions

1. Sterile technique shall be applied during installation and removal of the catheter.
2. It is recommended to be replaced after 7 days of use, however, shall not exceed 30 days.

Cautions

1. Federal law restricts the device to sale by or on the order of a physician.
2. Never clamp the catheter. Use catheter plug or valve if required.

Warnings

1. Single Use ONLY. Do NOT reuse. Reuse may cause infection, illness, injury or death.
2. Read the instructions carefully before using the device.
3. Do not use this device beyond the labeled expiration date.
4. Do not use the device if the packaging is punctured or damaged.
5. Do not force catheter as injury to the urethra can occur.
6. Do not use needle.
7. Do not exceed recommended capacities.

Potential Adverse Events

With all Foley Catheters, the following potential adverse events could occur:

1. Urinary tract infection
2. Perforation or laceration of the urethra
3. Perforation or laceration of the bladder
4. Abrasions

Symbol Description

	Catalog Number		Do not resterilize
	Lot Number		Do not use if package is damaged
	Use by		Prescription Use Only (US)
	Consult instructions for use		Foley Catheter is MR Safe
	Single Use		Manufacturer Contact Information
	Sterilized by Ethylene Oxide		Device is non-pyrogenic



SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor
Sunny Isles Beach, FL 33160
(424) 309-8523

87654321

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL

1. PRIOR TO PRINTING, VERIFY THE INFORMATION ON THE LABEL IS ACCURATE WHEN COMPARED AGAINST THE WORK ORDER.

2. PRIOR TO PRINTING, VERIFY THE LABEL HAS GONE THROUGH THE INSTALLATION QUALIFICATION (IQ) PROCESS FOR LABELS. REFER TO FM-0044.

3. LABEL CONTENT SHALL BE DEVELOPED PER WI-0009.

2. MATERIAL

1. LABEL: STANDARD WHITE THERMAL TRANSFER

2. RIBBON: WAX / RESIN

3. DIMENSIONS

1. () INDICATES LABEL REFERENCE DIMENSIONS.

4. SAMPLING PLAN

1. IN ADDITION TO THE LABEL GENERATION AND RECONCILIATION FORM, FM-0033, INSPECT N=3 ADDITIONAL SAMPLES AT THE BEGINNING OF THE RUN.

5. GS1 DATA MATRIX VALUES

1. (01): GLOBAL TRADE ITEM NUMBER (GTIN)

2. (10): BATCH OR LOT NUMBER

3. (17): EXPIRATION DATE

6. VISUAL INSPECTION - LABEL TEXT

1. POST-PRINTING, ALL TEXT ON THE LABEL SHALL BE LEGIBLE.

2. ALL VARIABLE LABEL CONTENT (LOT NUMBER, EXPIRATION DATE, ETC.) SHALL MATCH THE WORK ORDER.

3. THE GS1 DATA MATRIX SHALL BE FREE OF MISSING PIXELS AND SHALL NOT HAVE SMUDGES OR STREAKS THROUGH IT. WHEN SCANNED WITH A BARCODE SCANNER, THE TEXT ADJACENT TO THE GS1 DATA MATRIX SHALL MATCH THE TEXT DISPLAYED BY THE SCANNER.

4. (01): VERIFY THE GTIN IS DETAILED IN A HUMAN READABLE FORMAT ADJACENT TO THE GS1 DATA MATRIX. THE GTIN SHALL BE A 14 DIGIT NUMBER.

5. (10): VERIFY THE BATCH OR LOT NUMBER FORMAT IS IN COMPLIANCE WITH QP-019.

6. (17): VERIFY THE EXPIRATION DATE IS IN THE FOLLOWING FORMAT: YYMMDD.

7. VERIFY THE CATHETER STYLE AND CATHETER SIZE ARE REPRESENTATIVE OF A CONFIGURATION FOUND IN THE SKU TABLE.

8. VERIFY THE GTIN IS FOUND IN THE SKU TABLE AND CORRESPONDS WITH THE CORRECT CATHETER SIZE AND STYLE.

9. VERIFY THE REFERENCE NUMBER FOR THE CATHETER CORRESPONDS TO THAT FOUND IN THE SKU TABLE.

10. VERIFY "QUANTITY" IS 1.

7. VISUAL INSPECTION - LABEL SYMBOLS

1. POST-PRINTING, ALL SYMBOLS ON THE LABEL SHALL BE FREE OF SMUDGES / MISSING ARTWORK.

2. THE EXPIRATION DATE NEXT TO ITS CORRESPONDING SYMBOL SHALL BE IN THE FORMAT: YYYY-MM-DD.

REVISIONS

ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

ITEM NO.	PARTNO	DESCRIPTION	QTY.
1	SLQ-AW-4008	LABEL, POUCH, SILQ CLEARTRACT SPT CATHETER	1
2	PWY-5007	RIBBON, THERMAL TRANSFER, WAX/RESIN, 4.3W	1
3	SLQ-5004	LABEL, 2.3125 IN. X 7.50 IN.	1

1

2

3

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH ±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.

FDA
REGISTERED

ISO
CERTIFIED 13485

pathway
New Product Introduction

TITLE:
LABEL, POUCH, SILQ CLEARTRACT SPT CATHETER

SIZE B	DWG. NO. SLQ-LB-4008	REV A
SCALE: 1:1		SHEET 1 OF 2

87654321

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL

1. PRIOR TO PRINTING, VERIFY THE INFORMATION ON THE LABEL IS ACCURATE WHEN COMPARED AGAINST THE WORK ORDER.

2. PRIOR TO PRINTING, VERIFY THE LABEL HAS GONE THROUGH THE INSTALLATION QUALIFICATION (IQ) PROCESS FOR LABELS. REFER TO FM-0044.

3. LABEL CONTENT SHALL BE DEVELOPED PER WI-0009.

2. MATERIAL

1. LABEL: STANDARD WHITE THERMAL TRANSFER

2. RIBBON: WAX / RESIN

3. DIMENSIONS

1. () INDICATES LABEL REFERENCE DIMENSIONS.

4. SAMPLING PLAN

1. IN ADDITION TO THE LABEL GENERATION AND RECONCILIATION FORM, FM-0033, INSPECT N=3 ADDITIONAL SAMPLES AT THE BEGINNING OF THE RUN.

5. GS1 DATA MATRIX VALUES

1. (01): GLOBAL TRADE ITEM NUMBER (GTIN)

2. (10): BATCH OR LOT NUMBER

3. (17): EXPIRATION DATE

6. VISUAL INSPECTION - LABEL TEXT

1. POST-PRINTING, ALL TEXT ON THE LABEL SHALL BE LEGIBLE.

2. ALL VARIABLE LABEL CONTENT (LOT NUMBER, EXPIRATION DATE, ETC.) SHALL MATCH THE WORK ORDER.

3. THE GS1 DATA MATRIX SHALL BE FREE OF MISSING PIXELS AND SHALL NOT HAVE SMUDGES OR STREAKS THROUGH IT. WHEN SCANNED WITH A BARCODE SCANNER, THE TEXT ADJACENT TO THE GS1 DATA MATRIX SHALL MATCH THE TEXT DISPLAYED BY THE SCANNER.

4. (01): VERIFY THE GTIN IS DETAILED IN A HUMAN READABLE FORMAT ADJACENT TO THE GS1 DATA MATRIX. THE GTIN SHALL BE A 14 DIGIT NUMBER.

5. (10): VERIFY THE BATCH OR LOT NUMBER FORMAT IS IN COMPLIANCE WITH QP-019.

6. (17): VERIFY THE EXPIRATION DATE IS IN THE FOLLOWING FORMAT: YYMMDD.

7. VERIFY THE CATHETER STYLE AND CATHETER SIZE ARE REPRESENTATIVE OF A CONFIGURATION FOUND IN THE SKU TABLE.

8. VERIFY THE GTIN IS FOUND IN THE SKU TABLE AND CORRESPONDS WITH THE CORRECT CATHETER SIZE AND STYLE.

9. VERIFY THE REFERENCE NUMBER FOR THE CATHETER CORRESPONDS TO THAT FOUND IN THE SKU TABLE.

10. VERIFY "QUANTITY" IS 10.

7. VISUAL INSPECTION - LABEL SYMBOLS

1. POST-PRINTING, ALL SYMBOLS ON THE LABEL SHALL BE FREE OF SMUDGES / MISSING ARTWORK.

2. THE EXPIRATION DATE NEXT TO ITS CORRESPONDING SYMBOL SHALL BE IN THE FORMAT: YYYY-MM-DD.

REVISIONS

ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

ITEM NO.	PARTNO	DESCRIPTION	QTY.
1	SLQ-AW-4009	LABEL, BOX, SILQ CLEARTRACT SPT CATHETER	1
2	PWY-5007	RIBBON, THERMAL TRANSFER, WAX/RESIN, 4.3W	1
3	SLQ-5004	LABEL, 2.3125 IN. X 7.50 IN.	1

1

2

3

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH ±0.2 BEND ±2.0
NO PLACE DECIMAL ±0.5
ONE PLACE DECIMAL ±0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL

THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.

FDA
REGISTERED


ISO
CERTIFIED 13485

DO NOT SCALE DRAWING

pathway
New Product Introduction

TITLE:
LABEL, BOX, SILQ CLEARTRACT SPT CATHETER

SIZE B	DWG. NO. SLQ-LB-4009	REV A
SCALE: 1:1		SHEET 1 OF 2

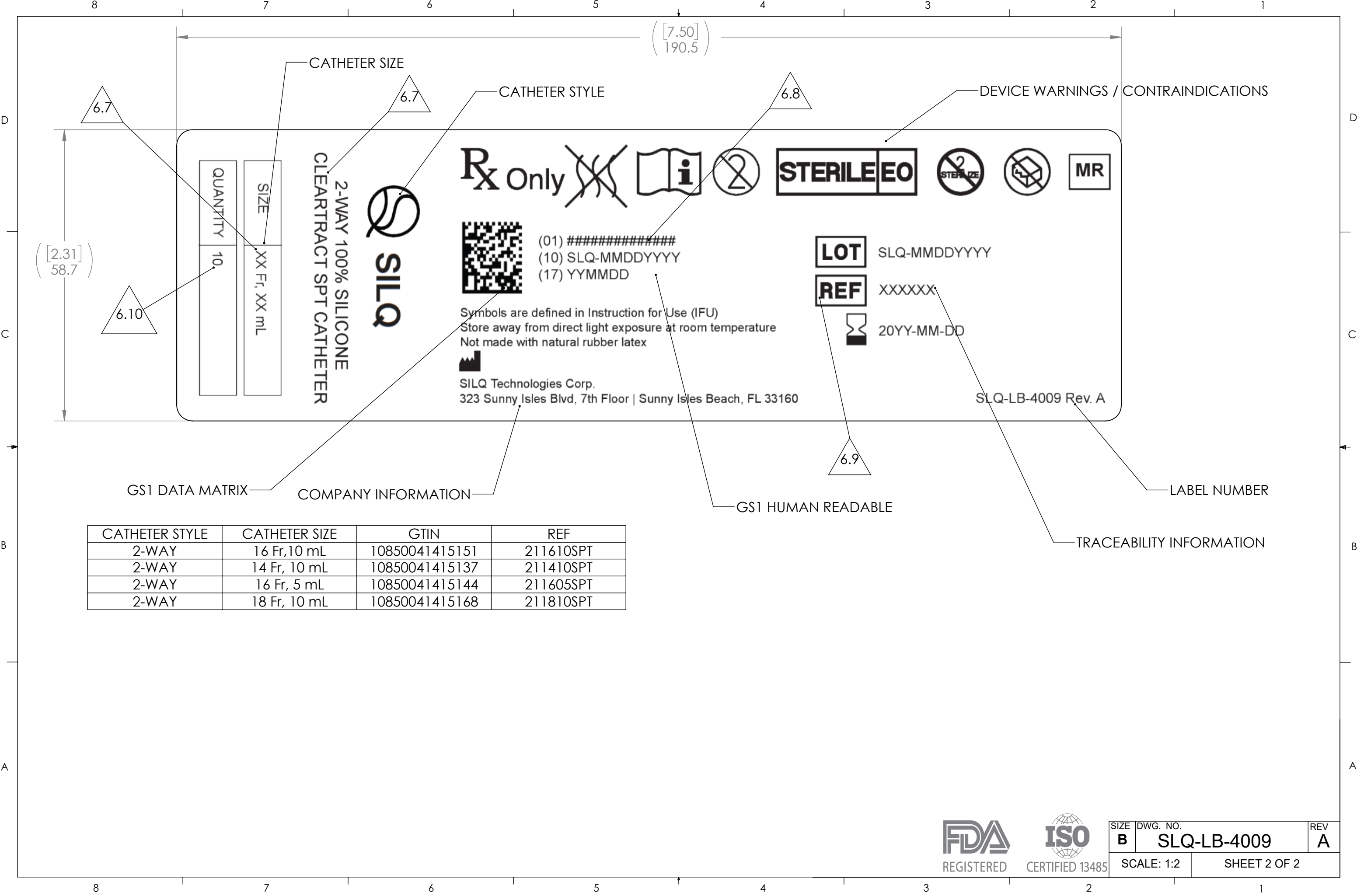

CATHETER STYLE	CATHETER SIZE	GTIN	REF
2-WAY	16 Fr, 10 mL	10850041415151	211610SPT
2-WAY	14 Fr, 10 mL	10850041415137	211410SPT
2-WAY	16 Fr, 5 mL	10850041415144	211605SPT
2-WAY	18 Fr, 10 mL	10850041415168	211810SPT

D

- C

DC

A

pathway
New Product Introduction

SIZE	DWG. NO.	REV
B	SLQ-LB-4010	A


SCALE: 1:2	SHEET 1 OF 2
------------	--------------



;

1

2

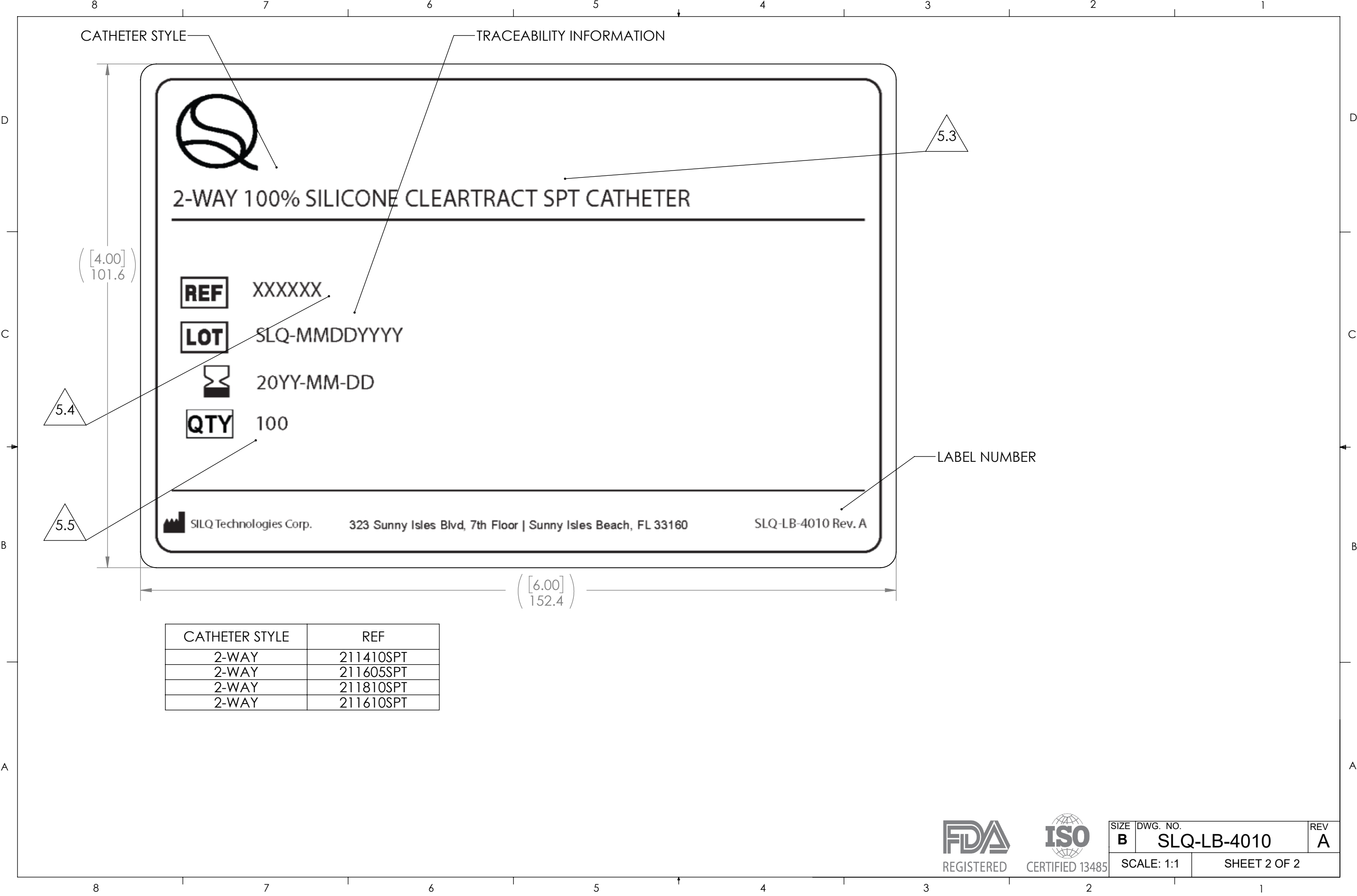

CATHETER STYLE	REF
2-WAY	211410SPT
2-WAY	211605SPT
2-WAY	211810SPT
2-WAY	211610SPT

#4 PRODUCT LABELING CONTENT REVIEW



0005 Rx Only

2-Way 100% Silicone ClearTract SPT Catheter

0084 Please read all instructions before using this device

Not made with natural rubber latex

Foley Catheter is MR Safe.

Model	Size
211410SPT	14 FR 10 mL
211605SPT	16 FR 5 mL
211610SPT	16 FR 10 mL
211810SPT	18 FR 10 mL

0021 0025 Contraindications

This device is contraindicated in the presence of conditions which create unacceptable risk during catheterization.

0006 Device Description

The 2-Way 100% Silicone ClearTract Catheter is a standard single-use 2-way Foley catheter that is constructed of medical grade silicone with a surface modification. It incorporates two (2) lumens, one for inflation/deflation of the balloon and the other for drainage of the urinary tract.

The surface of the ClearTract Catheter has a surface modification that reduces surface friction.

0014 0015 0018 0019 0021 Intended Use

The 2-Way 100% Silicone ClearTract Catheter is intended for drainage of the urinary tract. Catheterization is generally accomplished by inserting the catheter through the urethra and into the bladder. However, drainage is sometimes accomplished by suprapubic placement or other placement of the catheter, such as nephrostomy tract. Intended population is adults and pediatrics.

Directions

0017 0016 0020 Installation & Inflation

1. When placement of a urethral catheter is contraindicated or unsuccessful, percutaneous suprapubic urinary bladder catheterization or nephrostomy can be used in accordance with the procedures and policies of the hospital or healthcare facility.
2. Prepare application site in a sterile fashion.
3. Visually inspect the catheter for any signs of damage or deterioration prior to use.
4. If a trial inflation is conducted, follow steps 7 & 8 on the sterile field then evacuate the sterile water with the syringe. Ensure the balloon is empty prior to insertion.
5. Apply medical grade lubricant to distal end of the catheter.
6. Gently advance the catheter through the urethra until distal tip is beyond the bladder neck.
7. Confirm the position of the distal tip. The return of urine typically indicates the distal tip is in the bladder.
8. Fill Luer-tipped syringe with sterile water and connect to the inflation valve.
9. Fill the catheter balloon to the labeled volume. Do NOT exceed labeled volume.
10. Slowly withdraw the catheter until resistance is felt. This is a sign the balloon is resting at the bladder neck.



Deflation & Removal

1. Gently connect Luer-tipped syringe to the inflation valve. Gently aspirate while allowing the pressure in the balloon to assist filling the syringe with water occupying the balloon. Remove and reconnect the syringe in the inflation valve if there is no or slow deflation. Rapid deflation could collapse the lumen and prevent adequate deflation of the balloon.
2. If during the removal process the balloon ruptures, ensure all fragments are carefully removed from patient.
3. Discard after use.

Precautions

1. Sterile technique shall be applied during installation and removal of the catheter.
2. It is recommended to be replaced after 7 days of use, however, shall not exceed 30 days.

Cautions

1. Federal law restricts the device to sale by or on the order of a physician.
2. Never clamp the catheter. Use catheter plug or valve if required.

Warnings

1. Single Use ONLY. Do NOT reuse. Reuse may cause infection, illness, injury or death.
2. Read the instructions carefully before using the device.
3. Do not use this device beyond the labeled expiration date.
4. Do not use the device if the packaging is punctured or damaged.
5. Do not force catheter as injury to the urethra can occur.
6. Do not use needle.
7. Do not exceed recommended capacities.

Potential Adverse Events

With all Foley Catheters, the following potential adverse events could occur:

1. Urinary tract infection
2. Perforation or laceration of the urethra
3. Perforation or laceration of the bladder
4. Abrasions

Symbol Description

	Catalog Number		Do not resterilize
	Lot Number		Do not use if package is damaged
	Use by		Prescription Use Only (US)
	Consult instructions for use		Foley Catheter is MR Safe
	Single Use		Manufacturer Contact Information
	Sterilized by Ethylene Oxide		Device is non-pyrogenic



SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor
Sunny Isles Beach, FL 33160
(424) 309-8523



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	1

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex



SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

0001
0002

0006



Rx Only



0085



0010 0011 0012



SLQ-MMDDYYYY



XXXXXX



(01) #####
(10) SLQ-MMDDYYYY
(17) YYMMDD

0009



20YY-MM-DD

SLQ-LB-4008 Rev. A

2-WAY 100% SILICONE
CLEARTRACT SPT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	10



Rx Only

0010 0011 0012



(01) #####
(10) SLQ-MMDDYYYY
(17) YYMMDD

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex



0001

SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

0002

0006

0085



STERILE EO



MR

LOT

SLQ-MMDDYYYY

REF


XXXXXX



20YY-MM-DD

0009

SLQ-LB-4009 Rev. A

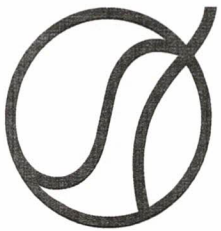

2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF XXXXXX

LOT SLQ-MMDDYYYY

0009  20YY-MM-DD

QTY 100



SILQ Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

0001

0002

	Label Content Design Requirements		Doc No.: WI-0009
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: C
			Page 2 of 7

- 5.1 It is the responsibility of the individual that designs product labels and labeling to comply with this procedure.

6.0 General Requirements

- 6.1 Paper stock and printing materials shall be selected based on its ability to maintain legibility throughout the labeling life (e.g., processing, storage, handling, and distribution).
- 6.2 Adhesive, if applicable, shall be selected based on its ability to maintain adherence throughout the labeling life (e.g., processing, storage, handling, and distribution).
- 6.3 Additional labeling requirements maybe required based on product specific regulations or guidance documents.

7.0 Labeling Requirements

#	Reference	Requirement Description	Device / Label / IFU	Markets
✓ 0001	21CFR801.1 21CFR812.5	Device label shall include the name of manufacturer, packer, or distributor. <i>Abbreviations for "Incorporated", etc. may be used. Qualifying phrase to distinguish relationship "Manufactured for ____", etc. may be used.</i>	Label	US
✓ 0002	21CFR801.1 21CFR812.5	Device label shall include street address, city, state, and zip code of manufacturer, packer, or distributor. <i>Street address may be omitted if listed in current directory.</i>	Label	US
✓ 0003	21CFR801.15	The predominant language shall be English.	Label / IFU	US
✓ 0004	21CFR801.15	Use of symbols shall include either adjacent text in all language(s) or a glossary of symbols.	Label / IFU	US
✓ 0005	21CFR801.15 21CFR801.109	The symbol statement "Rx only" or "[rx] only" or the statement "Caution: Federal law restricts this device to sale by or on the order of a ____" (i.e. "physician", etc.) <i>Except for surgical instruments.</i> "Only" needs to immediately follow "Rx." However, the symbol statement "Rx only" does not necessarily need to be bracketed in quotation marks, and the word "only" may appear in upper or lower case letters, for example, Rx only, Rx Only, or Rx ONLY.	Label / IFU	US
✓ 0006	21CFR801.109	The label of the device includes the method of its application or use. <i>Except for surgical instruments.</i>	Label / IFU	US
✓ 0007	21CFR801.15	If a glossary of symbols is used, it shall contain: - Each symbol used in the labeling - Title and designation number of standard(s) - Title of symbol and reference number (if any) - The symbols meaning	IFU	US
N/A 0008	21CFR801.15	If any representation appears in a foreign language, then all labeling requirements shall also appear in that foreign language.	Label / IFU	US
✓ 0009	21CFR801.18	Dates shall be expressed in the format YYYY-MM-DD (i.e. 4-digit year, 2-digit month, 2-digit day, each separated by a hyphen).	Label / IFU	US
✓ 0010	21CFR801.20 MDR Annex I 23.2	Label shall include a unique device identifier (UDI).	Label	US

	Label Content Design Requirements		Doc No.: WI-0009
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: C
			Page 3 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
✓ 0011	21CFR801.40	UDI must be presented in two forms: - Plain-text - Automatic Identification and Data Capture (AIDC)	Label	US
✓ 0012	21CFR801.40	UDI includes the following information: - Device identifier, which is a fixed portion that identifies the specific version or model of a device and the labeler of the device - Production identifier, which is a variable portion that identifies one or more of the following: lot number, serial number, expiration date, manufacturing date	Label	US
N/A 0013	21CFR801.45	Device meant to be reprocessed must be directly marked with UDI in either or both: - Plain-text - Automatic Identification and Data Capture (AIDC) or alternative technology <i>Except when it interferes with safety and efficacy, or not technically feasible.</i>	Device	US
✓ 0014	21CFR801.5	Instruction for Use includes statements of all conditions, purposes or uses for which such device is intended, prescribed or recommended.	IFU	US
✓ 0015	21CFR801.5	Instruction for Use includes Quantity of dose for: - each intended use - persons of different ages and physical conditions	IFU	US
✓ 0016	21CFR801.5 21CFR801.109	Instruction for Use includes frequency of administration or application.	IFU	US
✓ 0017	21CFR801.5 21CFR801.109	Instruction for Use includes duration of administration or application.	IFU	US
✓ 0018	21CFR801.5 21CFR801.109	Instruction for Use includes time of administration or application in relation to meals, onset of symptoms or other time factors.	IFU	US
✓ 0019	21CFR801.5 21CFR801.109	Instruction for Use includes route or method of administration or application.	IFU	US
✓ 0020	21CFR801.5	Instruction for Use includes preparation for use.	IFU	US
✓ 0021	21CFR801.109	Instruction for Use includes indications for Use, Effects, Hazards, Contraindications, side effects, precautions.	IFU	US
✓ 0022	21CFR801.109	Instruction for Use includes date of issuance or latest revision.	IFU	US
N/A 0023	21CFR812.5	An Investigational Device Exemption (IDE) label includes the quantity of contents.	Label	US
N/A 0024	21CFR812.5	An IDE label includes "CAUTION--Investigational device. Limited by Federal (or United States) law to investigational use."	Label	US
✓ 0025	21CFR812.5	The label or other labeling shall describe all relevant contraindications, hazards, adverse effects, interfering substances or devices, warnings, and precautions.	Label / IFU	US
N/A 0026	21CFR812.5	The labeling of an investigational device shall not bear any statement that is false or misleading in any particular and shall not represent that the device is safe or effective for the purposes for which it is being investigated.	Label / IFU	US

	Label Content Design Requirements		Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 4 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0027	21CFR812.5	An investigational device shipped solely for research on or with laboratory animals shall bear on its label the following statement: "CAUTION--Device for investigational use in laboratory animals or other tests that do not involve human subjects."	Label	US
0028	MDR Annex I 23.1	Information provided (medium, format, content, legibility, location) is appropriate for the device and healthcare practitioner.	Label / IFU	EU
0029	MDR Annex I 23.1	Where practical, information shall appear on the device itself.	Device	EU
0030	MDR Annex I 23.1	Information on device labeling shall take the form of symbols, where appropriate. In areas for which no standards exist, the symbols and colors shall be described in the documentation supplied with the device.	Label / IFU	EU
0031	MDR Annex I 23.1	Labels shall be provided in human-readable format and supplemented with machine-readable information.	Label	EU
0032	MDR Annex I 23.1	Each device shall be accompanied by Instructions for Use (IFU). <i>Except for class I and IIa devices that don't require information for safe use.</i>	IFU	EU
0033	MDR Annex I 23.1, 23.2, 23.3	Device name shall appear on the device label, sterile barrier packaging and IFU.	Label / IFU	EU
0034	MDR Annex I 23.1, 23.2, 23.3	Details needed to identify the device / device description shall appear on the device label, sterile barrier packaging.	Label	EU
0035	MDR Annex I 23.1, 23.2	Packaging contents shall appear on the device label.	Label	EU
0036	MDR Annex I 23.1, 23.2, 23.4	Intended for use statement shall appear on the device label.	Label	EU
0037	MDR Annex I 23.1, 23.2, 23.3, 23.4	Name and address of the manufacturer shall appear on the device label, sterile barrier packaging, IFU.	Label / IFU	EU
0038	MDR Annex I 23.2	Name and Address of authorized representative shall appear on the device label. <i>Except where the manufacturer resides within the EU.</i>	Label	EU
0039	MDR Annex I 23.2, 23.4	Where applicable, device label and IFU shall include an indication device contains: - Medicinal substance, human blood, plasma derivative - Tissues, cells or derivatives of human origin - Tissues, cells or derivatives of animal origin	Label / IFU	EU
0040	MDR Annex I 23.2, 23.4	Where applicable, device label and IFU shall include an indication device contains: - Substances which are carcinogenic, mutagenic or toxic to reproduction - Substances having endocrine-disrupting properties	Label / IFU	EU
0041	MDR Annex I 23.2	The label shall include the lot or serial number of the device along with the appropriate symbol.	Label	EU
0042	MDR Annex I 23.2, 23.3	The sterile barrier packaging shall include the date of manufacture (YYYY-MM).	Label	EU
0043	MDR Annex I 23.2, 23.3	The label and sterile barrier packaging shall include the expiration date (YYYY-MM), if applicable. Otherwise the label shall include the date of manufacture.	Label	EU
0044	MDR Annex I 23.2	The label shall include the UDI number.	Label	EU

	Label Content Design Requirements <i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Doc No.: WI-0009 Rev.: C Page 5 of 7
--	---	--

#	Reference	Requirement Description	Device / Label / IFU	Markets
0045	MDR Annex I 23.2, 23.4	The label and IFU shall indicate any storage and / or handling conditions that apply.	Label / IFU	EU
0046	MDR Annex I 23.3	The sterile barrier packaging shall be indicated it is the sterile barrier.	Label	EU
0047	MDR Annex I 23.2, 23.3, 23.4	The label, sterile barrier packaging, and IFU shall indicate its method and state of sterility that applies, with the appropriate symbol.	Label / IFU	EU
0048	MDR Annex I 23.2	The label shall include any immediate warnings or precautions that apply.	Label	EU
0049	MDR Annex I 23.2, 23.4	The label and IFU shall indicate the device is for single use with the appropriate symbol, if that applies.	Label / IFU	EU
0050	MDR Annex I 23.2	The label shall indicate the device is a reprocessed single-use device, and the number of reprocessing cycles performed any limitations that apply.	Label	EU
0051	MDR Annex I 23.2	The label and sterile barrier packaging shall indicate the device is "custom-made" if applicable.	Label	EU
0052	MDR Annex I 23.2	The label shall indicate the device is a medical device.	Label	EU
0053	MDR Annex I 23.2, 23.3	The label and sterile barrier packaging shall indicate the device is "exclusively for clinical investigation" if applicable.	Label	EU
0054	MDR Annex I 23.2, 23.4	The label and IFU shall indicate the device is composed of / combination of substances intended to be introduced / applied.	Label / IFU	EU
0055	MDR Annex I 23.2	The label shall indicate the serial number for implantable devices, if applicable.	Label	EU
0056	MDR Annex I 23.3	The sterile barrier packaging shall include an instruction to check the IFU if the sterile packaging is damaged or unintentionally opened before use.	Label	EU
0057	MDR Annex I 23.4	Instructions for Use includes device's intended purpose with a clear specification of indications, contra-indications, the patient target group or groups, and of the intended users, as appropriate.	IFU	EU
0058	MDR Annex I 23.4	Instructions for Use includes where applicable, a specification of the clinical benefits to be expected.	IFU	EU
0059	MDR Annex I 23.4	Instructions for Use includes where applicable, links to the summary of safety and clinical performance.	IFU	EU
0060	MDR Annex I 23.1, 23.4	Instructions for Use includes the performance characteristics of the device.	IFU	EU
0061	MDR Annex I 23.4	Instructions for Use includes where applicable, information allowing the healthcare professional to verify if the device is suitable and select the corresponding software and accessories.	IFU	EU
0062	MDR Annex I 23.4	Instructions for Use includes where applicable, information allowing the healthcare professional to verify if the device is suitable and select the corresponding software and accessories.	IFU	EU
0063	MDR Annex I 23.1, 23.4	Instructions for Use includes any residual risks, limitations, contra-indications, precautions, warnings and any undesirable side-effects, including information to be conveyed to the patient in this regard.	IFU	EU

N/A

N/A

	Label Content Design Requirements		Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 6 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0064	MDR Annex I 23.4	Instructions for Use includes specifications the user requires to use the device appropriately (e.g. degree of accuracy).	IFU	EU
0065	MDR Annex I 23.4	Instructions for Use includes details of any preparatory treatment or handling of the device before it is ready for use or during its use, such as sterilisation, final assembly, calibration, etc., including the levels of disinfection required to ensure patient safety and all available methods for achieving those levels of disinfection.	IFU	EU
0066	MDR Annex I 23.4	Instructions for Use includes any requirements for special facilities, or special training, or particular qualifications of the device user and/or other persons.	IFU	EU
0067	MDR Annex I 23.4	Instructions for Use includes the information needed to verify whether the device is properly installed and is ready to perform safely and as intended.	IFU	EU
0068	MDR Annex I 23.4	Instructions for Use includes instructions in the event of the sterile packaging being damaged or unintentionally opened before use, if supplied sterile.	IFU	EU
0069	MDR Annex I 23.4	Instructions for Use includes appropriate instructions for sterilization, if applicable.	IFU	EU
0070	MDR Annex I 23.4	Instructions for Use includes if the device is reusable, information on the appropriate processes for allowing reuse, including cleaning, disinfection, packaging and, where appropriate, the validated method of re-sterilization.	IFU	EU
0071	MDR Annex I 23.4	Instructions for Use includes an indication, if appropriate, that a device can be reused only if it is reconditioned under the responsibility of the manufacturer to comply with the general safety and performance requirements.	IFU	EU
0072	MDR Annex I 23.4	Instructions for Use includes information on known characteristics and technical factors known to the manufacturer that could pose a risk if the device were to be re-used. Applied to single use devices.	IFU	EU
0073	MDR Annex I 23.4	Instructions for Use includes information to identify other devices and equipment in order to obtain a safe combination, and / or any known restrictions. Applies to devices intended for use together with other devices and/or general-purpose equipment.	IFU	EU
0074	MDR Annex I 23.4	Instructions for Use includes detailed information as to the nature, type and where appropriate, the intensity and distribution of the emitted radiation and the means of protecting the patient, user, or other person from unintended radiation during use of the device. Applies to devices that emit radiation for medical purposes.	IFU	EU
0075	MDR Annex I 23.4	Instructions for Use includes information that allows the user and/or patient to be informed of any warnings, precautions, contra- indications, measures to be taken and limitations of use regarding the device.	IFU	EU
0076	MDR Annex I 23.4	Instructions for Use includes the overall qualitative and quantitative information on the materials and substances to which patients can be exposed with implantable devices.	IFU	EU
0077	MDR Annex I 23.4	Instructions for Use includes warnings or precautions to be taken in order to facilitate the safe disposal of the device, its accessories and the consumables used with it, if any.	IFU	EU

	Label Content Design Requirements		Doc No.: WI-0009
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: C
			Page 7 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0078	MDR Annex I 23.4	Instructions for Use includes the circumstances in which the user should consult a healthcare professional, if intended for use by lay persons.	IFU	EU
0079	MDR Annex I 23.4	Instructions for Use includes date of issue of the instructions for use or, if they have been revised, date of issue and identifier of the latest revision of the instructions for use.	IFU	EU
0080	MDR Annex I 23.4	Instructions for Use includes notice to the user and/or patient that any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.	IFU	EU
0081	MDR Annex I 23.4	Instructions for Use includes information to be supplied to the patient with an implanted device.	IFU	EU
0082	MDR Annex I 23.4	Instructions for Use includes, for devices that incorporate electronic programmable systems, including software, or software that are devices in themselves, minimum requirements concerning hardware, IT networks characteristics and IT security measures, including protection against unauthorized access, necessary to run the software as intended.	IFU	EU
0083	MDR Article 20#5, Annex V	CE conformity marking shall consist of the initials "CE" followed by the identification number of the notified body responsible for the device's conformity assessment. The various components of the CE marking must have substantially the same vertical dimension, which may not be less than 5 mm.	Label / IFU	EU
0084	FDA (1768)	If a manufacturer elects to include a statement in medical product labeling indicating that natural rubber latex or synthetic derivatives of natural rubber latex were not used as materials in the manufacture of their medical product and container, FDA recommends the use of the statement "Not made with natural rubber latex." If this statement is made without any qualification, it would apply to the medical product, its container, and any packaging.	Label / IFU	US
0085	ASTM F2503	For devices which are brought into the magnetic resonance (MR) environment, permanently mark the item as MR Safe, MR Conditional, or MR Unsafe using the icons with color or black and white options. Refer to ASTM F2503 TABLE 1 and TABLE 2 for details.	Device / Label / IFU	US

8.0 Records

8.1 Label Content documents will be processed per QP-002. Refer to QP-003 for storage location and retention requirements.

9.0 Flowchart – N/A

10.0 Change History

Revision	Change Description
A	Initial release. This procedure will replace procedure QP-7.09
B	Update header, references, and format of requirements (table).
C	Remove reference to FDA 1150 Guidance

	Label Content Design Requirements	Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: C
		Page 1 of 7

1.0 Purpose

- 1.1 The purpose of this document is to indicate the requirements for labels and labeling associated with medical devices.

2.0 Scope

- 2.1 This applies to labeling that is considered the device labels (primary and secondary packages) and the Instructions for Use (IFU)
- 2.2 The requirements indicated in this document are associated with regulatory requirements associated with FDA and EU.
- 2.3 The procedure does not apply to promotions materials.

3.0 References

- 3.1 21 CFR 801 – Labeling
- 3.2 21 CFR 812 – Investigation Device Exemption
- 3.3 ASTM 2503-20 – Standard Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment
- 3.4 ISO 15223-1 – Medical devices – Symbols to be used with medical device labels, labeling, and information to be supplied – Part 1: General requirements
- 3.5 FDA (1768) Guidance – Recommendations for Labeling Medical Products to Inform Users that the Product or Product Container is not Made with Natural Rubber Latex, December 2, 2014
- 3.6 QP-002 – Document Control
- 3.7 QP-003 – Quality Records

4.0 Definitions

- 4.1 Label – Display of written, printed, or graphic matter upon the immediate container.
- 4.2 Labeling – All labels and other written, printed, or graphic matter upon any article or any of its containers or wrapper, or accompanying such article at any time while a device is held for sale after shipment or delivery for shipment in interstate commerce.
- 4.3 Investigational device – A device that is the object of a clinical investigation or research involving one or more human subjects to determine its safety and effectiveness.
- 4.4 Unique Device Identifier (UDI) – An identifier that adequately identifies a device through its distribution.

5.0 Responsibilities

#5 PRODUCT REQUIREMENTS SPECIFICATION / TRACE MATRIX

Item Number	Item Name
PS-0006 rev E	Product Specification & Trace Matrix, SILQ ClearTract Foley Catheter

	Product Specification / Traceability Matrix – SILQ Foley Products	Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: E
		Page 1 of 9

1.0 Purpose

- 1.1 To document the design requirements for the SILQ ClearTract Foley Catheter products and provide traceability between the risk management reports (i.e, Hazard Analysis, Design Failure Mode Effects Analysis) and Verification & Validation deliverables.

2.0 Scope

- 2.1 SILQ Foley Catheter family of products

3.0 General Information

Table 1: General Information	
Attribute	Description
Device Name	SILQ 2-Way 100% Silicone ClearTract Catheter, 14Fr/10ml SILQ 2-Way 100% Silicone ClearTract SPT Catheter, 14Fr/10ml SILQ 2-Way 100% Silicone ClearTract Catheter, 16Fr/5ml SILQ 2-Way 100% Silicone ClearTract SPT Catheter, 16Fr/5ml SILQ 2-Way 100% Silicone ClearTract Catheter, 16Fr/5-10ml SILQ 2-Way 100% Silicone ClearTract Catheter, 16Fr/10ml SILQ 2-Way 100% Silicone ClearTract SPT Catheter, 16Fr/10ml SILQ 2-Way 100% Silicone ClearTract Catheter, 18Fr/10ml SILQ 2-Way 100% Silicone ClearTract SPTCatheter, 18Fr/10ml
Device's Intended Use / Indications for Use Statement	The 2-Way 100% SiliconeClearTract Catheter is intended for drainage of the urinary tract. Catheterization is generally accomplished by inserting the catheter through the urethra and into the bladder. However, drainage is sometimes accomplished by suprapubic placement or other placement of the catheter, such as nephrostomy tract. Intended population is adults and pediatrics.
User Types	Healthcare practitioners / Clinical user (e.g. nurses, doctors)
User Environment	Hospital
Patient Population	Any patient population
Region or Country	US
Agency	FDA
Device Class	II
Regulation - Description	21 CFR 876.5130 - Urological catheter and accessories
Product Code	EZL - Catheter, Retention Type, Balloon

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: E
								Page 2 of 9

4.0 Requirements and Trace Matrix

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
SECTION 1: Safety Requirements								
Biological	Patient contacting materials are biocompatible.	Patient contacting surfaces and materials shall meet biocompatibility requirement for external communicating device with tissue/bone/dentin contact for a prolong exposure duration (>24 to 30 days) per ISO 10993-1.	Device	1.1	1, 7, 9	HA#8.1, 8.2, 8.3, 9.1	TR-0071 TR-SLQ-0012	N/A
Sterilization	Product is supplied sterile	Device shall be sterile to a sterility assurance level (SAL) of 10 ⁻⁶ when terminally sterilized via Ethylene Oxide (EO).	Device	1.2	2	HA#7.1 DFMEA#5 PFMEA#23,27-32,35-37	TR-0075 TR-HDX-0009	SLQ-21XXXX
Sterilization	Product is non-pyrogenic	Endotoxins limits shall be ≤0.5 EU/mL or ≤20 EU/device for medical devices that directly or indirectly contact the cardiovascular or lymphatic system.	SPT Device	1.3	1, 7, 15, 17	HA#9.2 PFMEA#55	TR-SLQ-0012	SLQ-21XXXXSPT
Biological	EO residual levels meet ISO 10993-7 requirements.	The maximum allowable doses of EO and ECH delivered to patients for a Limited exposure devices: - The average daily dose of EO to patient shall not exceed 4 mg. - The average daily dose of ECH to patient shall not exceed 9 mg	Device	1.4	2, 3, 15, 17	HA#8.2 PFMEA#27-29	TR-0075 TR-HDX-0009	PWI-7374
SECTION 2: Material / Physical Requirements								
Physical	Product is clean	Cleanliness: When inspected by normal or corrected-to-normal vision without magnification under an illuminance of 300 lx to 700 lx, the outer surface of the device shall be smooth and appear free from particles and extraneous matter	Device	2.1	15-17	PFMEA#11	CHR-SLQ-6004	SLQ-6XXX Inspection Form

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: E
								Page 3 of 9

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Material	Product fill valve can be connected to by a syringe	Connector shall be made of rigid or semi-rigid materials.	Connector	2.2	12	HA#10.1, 10.2 DFMEA#4 PFMEA#7-12	N/A	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6000 SLQ-6001 SLQ-6002 SLQ-6004
Material	Product is made of silicone	The device shall be constructed on medical grade silicone material.	Device	2.3	15-17	HA#8.1, 8.2, 8.3, 9.1 DFMEA#1, 2, 3 PFMEA#7-12	TR-0071 TR-SLQ-0012	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6000 SLQ-6001 SLQ-6002 SLQ-6004
Physical	Product meets applicable recognized standards.	Balloon & Shaft Size - The proximal tip, the balloon, and the shaft (for 10 in. distal to the balloon) shall meet the requirements on size and tolerances as is relates to the labeled French size per the tolerances indicated below: Tip = ± 1 Fr Shaft = ± 1 Fr Balloon = + 4 Fr (Uninflated and Deflated after Immersion)	Device	2.4	7, 15-17	HA#10.3 DFMEA#1 PFMEA#7-12	TR-0072 TR-HDX-0001	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6000 SLQ-6001 SLQ-6002 SLQ-6004
Physical	Product is to be offered with a variety of balloon sizes.	The device family includes the following balloon volumes:	Device	2.5	15-17	HA#10.3 DFMEA#1 PFMEA#7-12	TR-0072 TR-HDX-0001	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6000 SLQ-6001 SLQ-6002 SLQ-6004
		Family						
		2-Way						
		Balloon Volume						
		5 ml – 10 ml						

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.							Rev.: E
								Page 4 of 9

Table 2: Requirements & Trace Matrix									
1	2	3		4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description		Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Physical	Product connects to a urine bag or collection device.	The device shall be able to connect into an industry standard or equivalent urine bag to facilitate sterile and easy urine collection.		Device	2.6	N/A	HA#10.1 DFMEA#4 PFMEA#7-12	TR-0072 TR-HDX-0001	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6000 SLQ-6001 SLQ-6002 SLQ-6004
Physical	Product is to be offered with 2-Way and 3-Way lumens.	The device shall be equipped with 2 lumens, which have the following purposes: 1. Drainage 2. Balloon Inflation		Device	2.7	15-17	HA#10.1, 10.2 DFMEA#1, 2 PFMEA#7-12	TR-0072 TR-HDX-0001	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6000 SLQ-6001 SLQ-6002 SLQ-6004
Physical	Product is to be offered with a variety of French sizes.	The device family includes the following range of sizes:		Device	2.8	15-17	HA#10.1 DFMEA#1 PFMEA#7-12	TR-0072 TR-HDX-0001	HDX-6000 HDX-6001 HDX-6002 HDX-6006 SLQ-6000 SLQ-6001 SLQ-6002 SLQ-6004
		Family	French Size						
		2-Way	14Fr – 18Fr						
SECTION 3: Performance Requirements									
Performance	Product meets applicable recognized standards.	Balloon Integrity (Resistance to Rupture) - The balloon shall be inflated easily with distilled or deionized water to labeled volume without showing any evidence of breakage throughout a 24 hour period.		Device	3.1	7	HA#10.1 DFMEA#3 PFMEA#7-12	TR-0072 TR-HDX-0001	N/A
Performance		Balloon Volume Maintenance - The balloon shall maintain its volume throughout a 24 hour period.		Device	3.2	7	HA#10.1 DFMEA#3 PFMEA#7-12	TR-0072 TR-HDX-0001	N/A

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: E
								Page 5 of 9

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Performance		Connector shall be non-interconnectables, which does not appear to provide a secure connection when forcefully assembled and shall easily disengage.	Connector	3.3	12	HA#10.1 DFMEA#4 PFMEA#7-12	TR-0072 TR-HDX-0001	N/A
Performance		Deflation Reliability (Failure to Deflate) - The balloon shall deflate to within 4 French sizes of the label French size within 15 min or be otherwise manipulated to effect drainage within this time period.	Device	3.4	7	HA#10.1, 10.4 DFMEA#4 PFMEA#7-12	TR-0072 TR-HDX-0001	N/A
Performance		Flow Rate though Drainage Lumen - Shall have a minimum average flow rate of 100 cm3/min.	Device	3.5	7	HA#10.2 DFMEA#1 PFMEA#7-12	TR-0072 TR-HDX-0001	N/A
Performance		Inflated Balloon Response to Traction - The entire, inflated balloon of catheters shall not pass into or through the funnel barrel neck (28 Fr x 6 cm) when a 1 lb weight is affixed from the distal catheter end and subjected to static (2 minutes) and dynamic load (2 ft drop).	Device	3.6	7	HA#10.1 DFMEA#3 PFMEA#7-12	TR-0072 TR-HDX-0001	N/A
Performance	Product is intended for urethral or suprapubic placement.	The device shall be easily inserted through the urethra until the deflated balloon passes through the urethral sphincter and in a manner that prevents urethral trauma.	Device	3.7	N/A	HA#10.3 DFMEA#1 PFMEA#7-12	TR-0072 TR-HDX-0001	N/A
SECTION 4: Shelf Life Requirements								
Stability	The product shall have a minimum 1-year shelf life.	The device shall have a final shelf-life period of 3 years, but initially be qualified for shorter shelf-life periods (e.g., 1 year).	Sterile, Packaged Device	4.1	15-16	HA#7.1, 9.1 DFMEA#1, 2, 3, 4, 5	TR-0073	SLQ-21XXXX
SECTION 5: Packaging & Labeling Requirements								
Performance	The product packaging meets applicable recognized standards.	Packaged devices shall maintain its sterility and functionality through transportation	Sterile, Packaged Device	5.1	5, 6, 15-17	HA#14.1 DFMEA#5 PFMEA#23	TR-0074	N/A

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: E
								Page 6 of 9

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
Packaging	The product is to be packaged in a Tyvek/poly pouch.	The peel strength for 1" MFG sealed section of the sterile barrier packaging shall be ≥ 1 lb.	Sterile Barrier	5.2	3, 15-17	HA#7.1, 14.1 DFMEA#5 PFMEA#23	TR-0074	SLQ-8XXXX TM-0003
Packaging		The sterile barrier packaging (include sealed area) shall not exhibit the following visual defects: 1. Unsealed Areas 2. Nonhomogeneous or Undersealed Areas 3. Over-sealed Areas 4. Narrow Seals 5. Channels 6. Wrinkles/Fold-overs/Cracks 7. Tears/Pinholes	Sterile Barrier	5.3	3, 15-17	HA#7.1, 14.1 DFMEA#5 PFMEA#23	TR-0074	SLQ-8XXXX
Packaging		The sterile barrier packaging shall not exhibit a constant stream of bubbles when submerge in water and internal pressurized to a determine test pressure.	Sterile Barrier	5.4	3, 15-17	HA#7.1, 14.1 DFMEA#5	TR-0074	TM-0005
Packaging	Multiple products are packaged in a shelf carton.	The device will be packaged in a tyvek/poly pouch., which 10 pouched devices will be packaged in a box.	Packaged Device	5.5	15-17	HA#7.1, 14.1 DFMEA#5 PFMEA#26	DHR-8XXX	HDX-8000 HDX-8001 HDX-8002 HDX-8006 SLQ-8000 SLQ-8001 SLQ-8002 SLQ-8004
Labeling	Product to be appropriately identified and labeled.	Product labeling shall remain legible and affixed to the primary and secondary packaging throughout the labeling life (e.g. processing, storage, handling, and distribution).	Labeling	5.6	5, 15-17	HA#12.1, 12.2 DFMEA#6,8, 10 PFMEA#18-19, 20-21, 25-26	DHR-8XXX	HDX-8000 HDX-8001 HDX-8002 HDX-8006 SLQ-8000 SLQ-8001 SLQ-8002 SLQ-8004

SECTION 5: User Requirements

	Product Specification / Traceability Matrix – SILQ Foley Products							Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: E
								Page 8 of 9

Table 2: Requirements & Trace Matrix								
1	2	3	4	5	6	7	8	9
Type	User/Customer Requirement	Product Requirement Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference / Design Output
User	N/A	A healthcare provider must be able to easily and quickly implement device.	Device	27	15	HA#13.1	N/A	SLQ-LBL-4000 SLQ-LBL-4001 SLQ-4006 SLQ-4007 SLQ-LBL-4008 SLQ-LBL-4009

5.0 Reference Documents

Table 3: Reference Documents				
Source No	Organization	Document No	Edition	Title
1	ANSI/ AAMI/ ISO	10993-1	2009	Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management process
2	ANSI/ AAMI/ ISO	11135	2014	Sterilization of health care products – Ethylene Oxide – Requirements for development, validation and routine control of a sterilization process for medical devices
3	ISO	11607-1	2019	Packaging for terminally sterilized medical devices – Part 1: Requirements for materials, sterile barrier systems and packaging systems
4	ANSI/ AAMI/ ISO	15223-1	2016	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied
5	ASTM	D4169	2016	Standard Practice For Performance Testing Of Shipping Containers And Systems
6	ASTM	D4332	2014	Standard Practice for Conditioning Containers, Packages, or Packaging Components for Testing.
7	ASTM	F623	2019	Standard Performance Specification for Foley Catheter
8	FDA	21 CFR 801	4/15/2018	Title 21 - Food and Drugs Chapter I – Food and Drug Administration Department of Health Services; Subchapter H – Medical Devices; Part 801 Labeling
9	FDA	FDA-2013-D-0350	6/16/2016	FDA Guidance: Use of International Standard ISO 10993, "Biological Evaluation of Medical Devices Part 1: Evaluation and Testing within a risk management process"
11	ISO	14971	2019	Medical Devices - Application of Risk Management to Medical Devices
12	ISO	80369-1	2018	Small-Bore Connectors For Liquids And Gases In Healthcare Applications - Part 1: General Requirements.

	Product Specification / Traceability Matrix – SILQ Foley Products		Doc No.: PS-0006
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: E
			Page 9 of 9

Table 3: Reference Documents				
Source No	Organization	Document No	Edition	Title
13	FDA	N/A	N/A	Guidance for the Content of Premarket Notifications for Conventional and Antimicrobial Foley Catheters
14	FDA	81 FR 38911	6/15/2016	Use of Symbols in Labeling
15	N/A	RM-0019	Latest	Risk Management Report, SILQ ClearTract Foley Catheter Hazard Analysis (HA)
16	N/A	RM-0020	Latest	Risk Management Report, SILQ ClearTract Foley Catheter Design Failure Mode and Effects Analysis (DFMEA)
17	N/A	RM-0021	Latest	Risk Management Report, SILQ ClearTract Foley Catheter Process Failure Mode and Effects Analysis (PFMEA)

6.0 Change History

Revision	Change Description
A	Initial Release
B	Update format, indications for use
C	Update format; name change Hydrophilix to SILQ; add 16Fr 5-10ml; general updates (spelling); correction to #22 (2-Way); remove #20 (Device is 2-way, not a 3-way), remove ref doc#10
D	Update references to HA and DFMEA, component specs SLQ-6004, labeling SLQ-4000, SLQ-4001, SLQ-4006; add 16Fr/10ml
E	Update TABLE 1 products, Table 2 references to column 6-9. General format/re-organize Table 2 with sections and renumber X.X requirements #s, Replace N/A with User/Customer requirements in column 2. Add requirements # 1.3, 1.4, 2.1, 5.4, 5.6. Clarify requirements # 2.5 and 2.8.

#6 DEVICE MASTER RECORD

BOM & DRAWING

Item Number	Item Name
SLQ-211410 rev A	SILQ 2-way Cleartract SPT 14 Fr x 10 ml Catheter, Sterile, 10-Pack
SLQ-211605 rev A	SILQ 2-way Cleartract SPT 16 Fr x 5 ml Catheter, Sterile, 10-Pack
SLQ-211610 rev A	SILQ 2-way Cleartract SPT 16 Fr x 10 ml Catheter, Sterile, 10-Pack
SLQ-211810 rev A	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack
SLQ-8100 rev A	SILQ 2-way Cleartract SPT 14 Fr x 10 ml Catheter, Non-Sterile, 10-Pack
SLQ-8101 rev A	SILQ 2-way Cleartract SPT 16 Fr x 5 ml Catheter, Non-Sterile, 10-Pack
SLQ-8103 rev A	SILQ 2-way Cleartract SPT 16 Fr x 10 ml Catheter, Non-Sterile, 10-Pack
SLQ-8104 rev A	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack

pathway

Pathway LLC

Printed by Aaron Rogers

Printed on 05/25/2023

Local time zone (GMT-08:00) Pacific Time (America/Los Angeles)

Page Arena > Workspace Items > Bill of Materials > Indented

All information contained in this document is proprietary and confidential.

SLQ-211410SPT

In Design

SILQ 2-way Cleartract SPT 14 Fr x 10 ml Catheter, Sterile, 10-Pack

Revision »A - ECO-001247 Status: Effective, Status Date: 01/12/2023 01:24:56 PM, Shared

Contains 10 first-level Items, 50 line Items, 41 unique Items, 33 of which are shared.

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE		REF	DES
▼	1	SLQ-8100 rev A	SILQ 2-way Cleartract SPT 14 Fr x 10 ml Catheter, Non-Sterile, 10-Pack	Subassembly	In Des	3	1 EA
	1	SLQ-6000 rev A	Catheter, 14 Fr x 10 mL, Silicone, NY Foley, 2-Way	Fabricated Part	In Des	3	10 EA
	2	SLQ-3000 rev B	Solution, SILQ Polymer	Purchased Component	In Prod	4	80 ML
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1HDX-3000 rev A	Solution, Hydrophilix Polymer	80 ML			
	3	SLQ-5000 rev A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	Packaging	In Prod	3	10 EA
	4	SLQ-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	Packaging	In Prod	4	10 EA
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1HDX-5001 rev A	Pouch, Chevron, 3.5 X 21, 10 EA Amcor TLP-058				
▶	5	SLQ-LB-4008 rev A	Label, Pouch, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	10 EA Max 12 month product shelf life
	6	SLQ-5016 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS, Pre-Printed	Packaging	In Prod	3	1 EA
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1SLQ-5015 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS	1 EA			
		2SLQ-5014 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .024 SBS	1 EA			
▶	7	SLQ-LB-4009 rev A	Label, Box, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	2 EA Max 12 month product shelf life
	8	SLQ-4007 rev A	Artwork, IFU, SILQ 2-Way ClearTract SPT Catheter	Labeling	In Prod	2	1 EA
	9	SLQ-5012 rev A	Box, 20 in. L x 18 in. W x 6 in. D, 200 LB Test	Packaging	In Prod	2	0.1 EA
	10	PWY-5000 rev A	Tape, Packaging, 2 Mil, 2 in x 110 yds, Clear	Packaging	In Prod	2	5.2 IN
▶	11	SLQ-LB-4010 rev A	Label, Shipper Box, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	0.2 EA Max 12 month product shelf life

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	    	REF DES
12 MP-SLQ-0002 rev B	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Manufacturing Procedure	In Prod	2	1 EA
13 FM-SLQ-0001 rev A	Form, Run Sheet for SILQ Automation Machine	Form	In Prod	2	1 EA
14 TM-HDX-0001 rev B	Test Method, UV Light Absorbance, Hydrophilix Catheters	Test Method	In Prod	2	1 each
15 TM-HDX-0002 rev A	Test Method, Rinse Test, Hydrophilix Catheter	Test Method	In Prod	2	1 each
16 MP-SLQ-0001 rev B	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	Manufacturing Procedure	In Prod	2	1 EA
2 WI-0006 rev D	Work Order Process	Work Instruction	In Prod	2	1 each
3 FM-0039 rev F	Form, Work Order	Form	In Prod	2	1 each
4 WI-0014 rev A	Line Clearance	Work Instruction	In Prod	2	1 each
5 MP-SLQ-0003 rev A	Manufacturing Process Flow, SILQ Foley Catheter	Manufacturing Procedure	In Prod	2	1 EA
6 SLQ-0001 rev A	Sterilization Load Shipping Pallet Configuration	Subassembly	In Prod	2	1 EA
▶ 7 QP-032 rev F	EO Sterilization Process Control	Quality System Procedure	In Prod	6	1 each
8 WI-0021 rev C	Work Instruction, Pyrogen (LAL) Evaluation Process	Work Instruction	In Prod	2	1 each
▶ 9 QP-010 rev F	Final Product Inspection	Quality System Procedure	In Prod	6	1 each
10 FM-0029 rev A	Form, Certificate of Conformance	Form	In Prod	2	1 each

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
1.

ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY SILQ PRIOR TO PRODUCTION.
2.


FINAL PRODUCT INSPECTION PER QP-010.
3.

PRODUCT EXPIRY DATE IS 12 MONTHS FROM THE DATE OF MANUFACTURE.
2. STERILIZATION
1.

PRODUCT IS STERILIZED VIA ETHYLENE OXIDE (ETO) STERILIZATION PER QP-032 CYCLE PWI-7374.
3. TESTING
1.

PYROGEN EVALUATION PER WI-0021 AT ≤ 0.5 EU/ML OR ≤ 20.0 EU/DEVICE

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

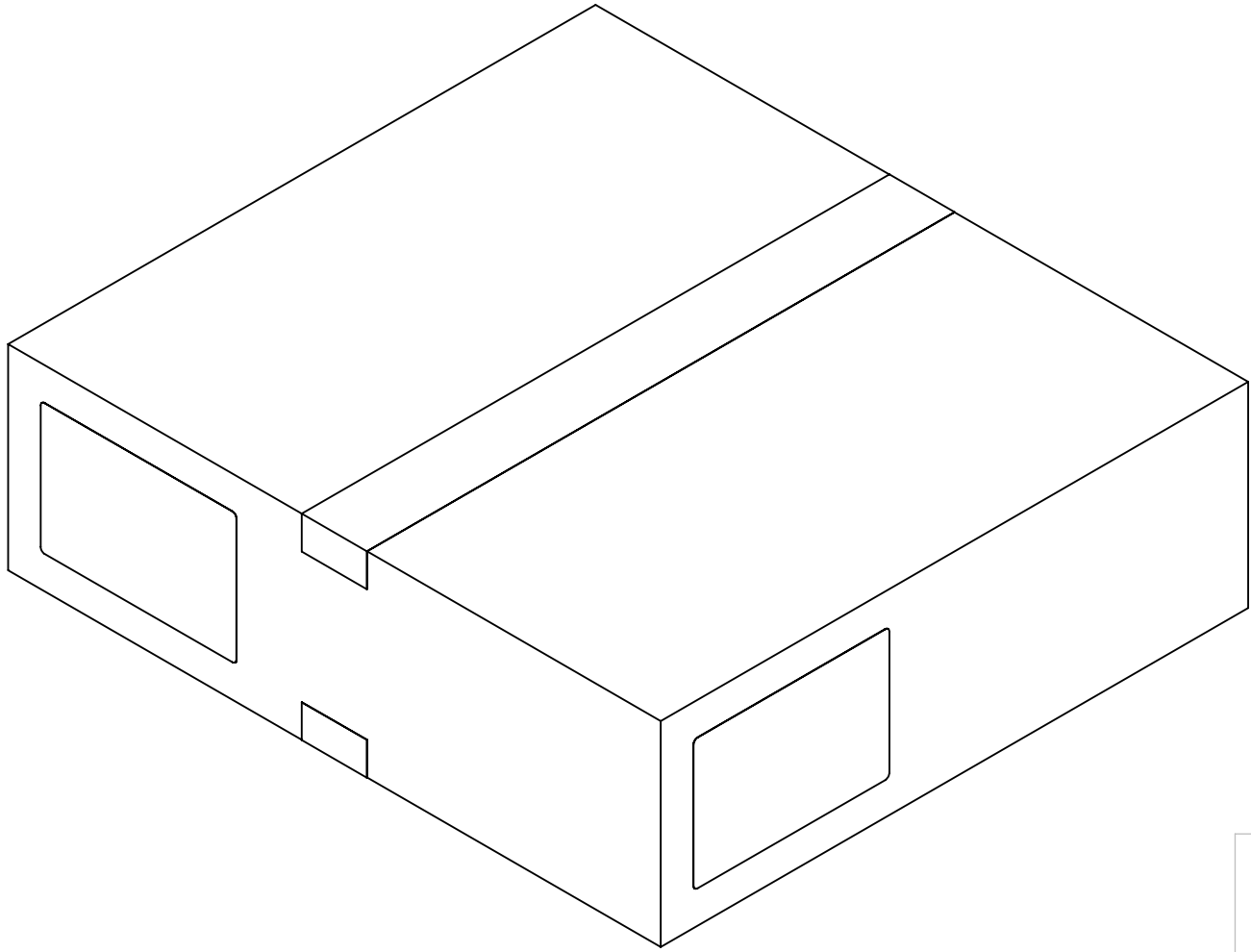

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH ± 0.2 BEND ± 2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ± 0.15
THREE PLACE DECIMAL ± 0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.



TITLE:
2-WAY CLEARTRACT CATHETER, 14 FR X 10 ML, STERILE, 10-PACK

SIZE	DWG. NO.	REV
B	SLQ-211410SPT	A

pathway

Pathway LLC

Printed by Aaron Rogers

Printed on 05/25/2023

Local time zone (GMT-08:00) Pacific Time (America/Los Angeles)

Page Arena > Workspace Items > Bill of Materials > Indented

All information contained in this document is proprietary and confidential.

SLQ-211605SPT

In Production

SILQ 2-way Cleartract SPT 16 Fr x 5 ml Catheter, Sterile, 10-Pack

Revision »A - ECO-001247 Status: Effective, Status Date: 01/12/2023 01:24:56 PM, Shared

Contains 10 first-level Items, 50 line Items, 41 unique Items, 33 of which are shared.

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE		REF	DES
▼	1	SLQ-8101 rev A	SILQ 2-way Cleartract SPT 16 Fr x 5 ml Catheter, Non-Sterile, 10-Pack	Subassembly	In Prod	3	1 EA
	1	SLQ-6001 rev A	Catheter, 16 Fr x 5 mL, Silicone, NY Foley, 2-Way	Fabricated Part	In Prod	3	10 EA
	2	SLQ-3000 rev B	Solution, SILQ Polymer	Purchased Component	In Prod	4	80 ML
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1HDX-3000 rev A	Solution, Hydrophilix Polymer	80 ML			
	3	SLQ-5000 rev A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	Packaging	In Prod	3	10 EA
	4	SLQ-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	Packaging	In Prod	4	10 EA
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1HDX-5001 rev A	Pouch, Chevron, 3.5 X 21, 10 EA Amcor TLP-058				
▶	5	SLQ-LB-4008 rev A	Label, Pouch, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	10 EA Max 36 month product shelf life
	6	SLQ-5016 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS, Pre-Printed	Packaging	In Prod	3	1 EA
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1SLQ-5015 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS	1 EA			
		2SLQ-5014 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .024 SBS	1 EA			
▶	7	SLQ-LB-4009 rev A	Label, Box, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	2 EA Max 36 month product shelf life
	8	SLQ-4007 rev A	Artwork, IFU, SILQ 2-Way ClearTract SPT Catheter	Labeling	In Prod	2	1 EA
	9	SLQ-5012 rev A	Box, 20 in. L x 18 in. W x 6 in. D, 200 LB Test	Packaging	In Prod	2	0.1 EA
	10	PWY-5000 rev A	Tape, Packaging, 2 Mil, 2 in x 110 yds, Clear	Packaging	In Prod	2	5.2 IN
▶	11	SLQ-LB-4010 rev A	Label, Shipper Box, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	0.2 EA Max 36 month product shelf life

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	    	REF DES
12 MP-SLQ-0002 rev B	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Manufacturing Procedure	In Prod	21 EA	
13 FM-SLQ-0001 rev A	Form, Run Sheet for SILQ Automation Machine	Form	In Prod	21 EA	
14 TM-HDX-0001 rev B	Test Method, UV Light Absorbance, Hydrophilix Catheters	Test Method	In Prod	21 each	
15 TM-HDX-0002 rev A	Test Method, Rinse Test, Hydrophilix Catheter	Test Method	In Prod	21 each	
16 MP-SLQ-0001 rev B	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	Manufacturing Procedure	In Prod	21 EA	
2 WI-0006 rev D	Work Order Process	Work Instruction	In Prod	21 each	
3 FM-0039 rev F	Form, Work Order	Form	In Prod	21 each	
4 WI-0014 rev A	Line Clearance	Work Instruction	In Prod	21 each	
5 MP-SLQ-0003 rev A	Manufacturing Process Flow, SILQ Foley Catheter	Manufacturing Procedure	In Prod	21 EA	
6 SLQ-0001 rev A	Sterilization Load Shipping Pallet Configuration	Subassembly	In Prod	21 EA	
7 QP-032 rev F	EO Sterilization Process Control	Quality System Procedure	In Prod	61 each	
8 WI-0021 rev C	Work Instruction, Pyrogen (LAL) Evaluation Process	Work Instruction	In Prod	21 each	
9 QP-010 rev F	Final Product Inspection	Quality System Procedure	In Prod	61 each	
10 FM-0029 rev A	Form, Certificate of Conformance	Form	In Prod	21 each	

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
1.

ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY SILQ PRIOR TO PRODUCTION.
2.


FINAL PRODUCT INSPECTION PER QP-010.
3.

PRODUCT EXPIRY DATE IS 12 MONTHS FROM THE DATE OF MANUFACTURE.
2. STERILIZATION
1.

PRODUCT IS STERILIZED VIA ETHYLENE OXIDE (ETO) STERILIZATION PER QP-032 CYCLE PWI-7374.
3. TESTING
1.

PYROGEN EVALUATION PER WI-0021 AT ≤ 0.5 EU/ML OR ≤ 20.0 EU/DEVICE

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

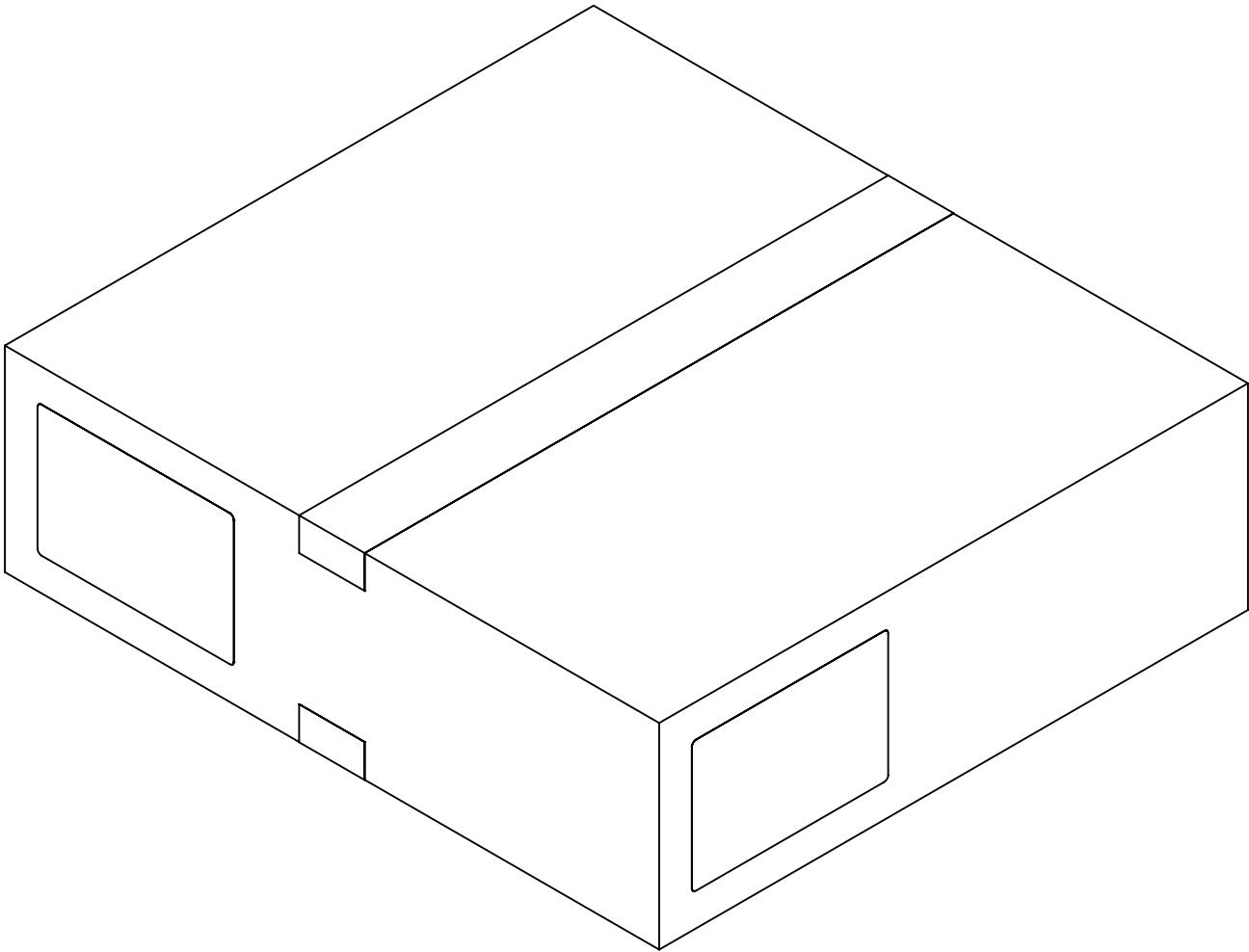

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH ± 0.2 BEND ± 2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ± 0.15
THREE PLACE DECIMAL ± 0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.

pathway

New Product Introduction

TITLE:
2-WAY CLEARTRACT CATHETER, 16 FR X 5 ML, STERILE, 10-PACK

SIZE B	DWG. NO. SLQ-211605SPT	REV A
-----------	---------------------------	----------

DO NOT SCALE DRAWING



pathway

Pathway LLC

Printed by Aaron Rogers

Printed on 05/25/2023

Local time zone (GMT-08:00) Pacific Time (America/Los_Angeles)

Page Arena > Workspace Items > Bill of Materials > Indented

All information contained in this document is proprietary and confidential.

SLQ-211610SPT

In Production

SILQ 2-way Cleartract SPT 16 Fr x 10 ml Catheter, Sterile, 10-Pack

Revision »A - ECO-001247 Status: Effective, Status Date: 01/12/2023 01:24:56 PM, Shared

Contains 10 first-level Items, 50 line Items, 41 unique Items, 33 of which are shared.

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE			
▼	1	SLQ-8104 rev A	SILQ 2-way Cleartract SPT 16 Fr x 10 ml Catheter, Non-Sterile, 10-Pack	Subassembly	In Prod	3	1 EA
	1	SLQ-6004 rev A	Catheter, 16 Fr x 10 mL, Silicone, NY Foley, 2-Way	Fabricated Part	In Prod	3	10 EA
	2	SLQ-3000 rev B	Solution, SILQ Polymer	Purchased Component	In Prod	4	80 ML
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1HDX-3000 rev A	Solution, Hydrophilix Polymer	80 ML			
	3	SLQ-5000 rev A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	Packaging	In Prod	3	10 EA
	4	SLQ-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	Packaging	In Prod	4	10 EA
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1HDX-5001 rev A	Pouch, Chevron, 3.5 X 21, 10 EA Amcor TLP-058				
▶	5	SLQ-LB-4008 rev A	Label, Pouch, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	10 EA Max 36 month product shelf life
	6	SLQ-5016 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS, Pre-Printed	Packaging	In Prod	3	1 EA
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1SLQ-5015 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS	1 EA			
		2SLQ-5014 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .024 SBS	1 EA			
▶	7	SLQ-LB-4009 rev A	Label, Box, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	2 EA Max 36 month product shelf life
	8	SLQ-4007 rev A	Artwork, IFU, SILQ 2-Way ClearTract SPT Catheter	Labeling	In Prod	2	1 EA
	9	SLQ-5012 rev A	Box, 20 in. L x 18 in. W x 6 in. D, 200 LB Test	Packaging	In Prod	2	0.1 EA
	10	PWY-5000 rev A	Tape, Packaging, 2 Mil, 2 in x 110 yds, Clear	Packaging	In Prod	2	5.2 IN
▶	11	SLQ-LB-4010 rev A	Label, Shipper Box, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	0.2 EA Max 36 month product shelf life

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	    	
12 MP-SLQ-0002 rev B	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Manufacturing Procedure	In Prod	2	1 EA
13 FM-SLQ-0001 rev A	Form, Run Sheet for SILQ Automation Machine	Form	In Prod	2	1 EA
14 TM-HDX-0001 rev B	Test Method, UV Light Absorbance, Hydrophilix Catheters	Test Method	In Prod	2	1 each
15 TM-HDX-0002 rev A	Test Method, Rinse Test, Hydrophilix Catheter	Test Method	In Prod	2	1 each
16 MP-SLQ-0001 rev B	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	Manufacturing Procedure	In Prod	2	1 EA
2 WI-0006 rev D	Work Order Process	Work Instruction	In Prod	2	1 each
3 FM-0039 rev F	Form, Work Order	Form	In Prod	2	1 each
4 WI-0014 rev A	Line Clearance	Work Instruction	In Prod	2	1 each
5 MP-SLQ-0003 rev A	Manufacturing Process Flow, SILQ Foley Catheter	Manufacturing Procedure	In Prod	2	1 EA
6 SLQ-0001 rev A	Sterilization Load Shipping Pallet Configuration	Subassembly	In Prod	2	1 EA
▶ 7 QP-032 rev F	EO Sterilization Process Control	Quality System Procedure	In Prod	6	1 each
8 WI-0021 rev C	Work Instruction, Pyrogen (LAL) Evaluation Process	Work Instruction	In Prod	2	1 each
▶ 9 QP-010 rev F	Final Product Inspection	Quality System Procedure	In Prod	6	1 each
10 FM-0029 rev A	Form, Certificate of Conformance	Form	In Prod	2	1 each


NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY SILQ PRIOR TO PRODUCTION.

2. FINAL PRODUCT INSPECTION PER QP-010.

3. PRODUCT EXPIRY DATE IS 12 MONTHS FROM THE DATE OF MANUFACTURE.
2. STERILIZATION
1. PRODUCT IS STERILIZED VIA ETHYLENE OXIDE (ETO) STERILIZATION PER QP-032 CYCLE PWI-7374.
3. TESTING
1. PYROGEN EVALUATION PER WI-0021 AT ≤ 0.5 EU/ML OR ≤ 20.0 EU/DEVICE

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

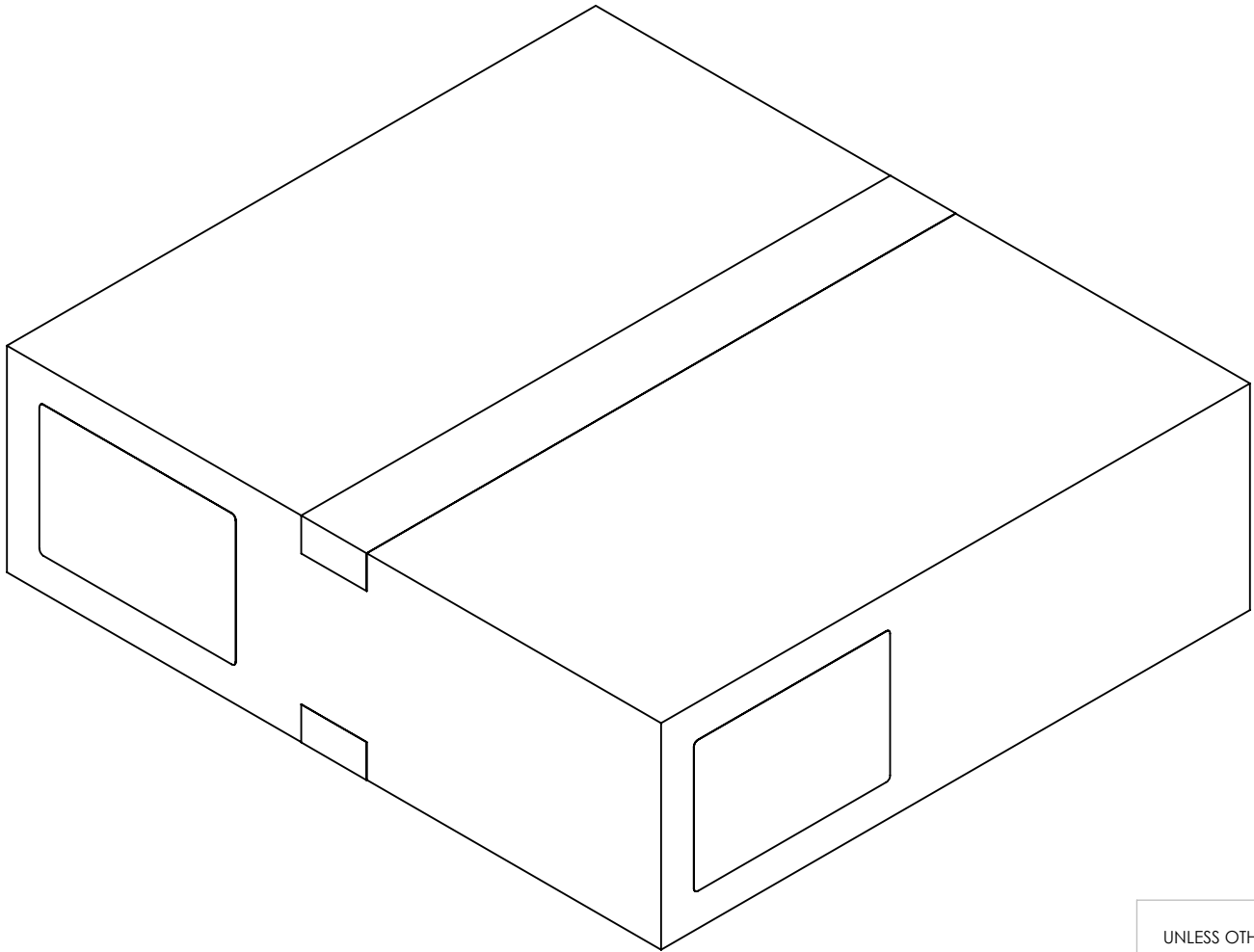

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH ± 0.2 BEND ± 2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ± 0.15
THREE PLACE DECIMAL ± 0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.



New Product Introduction

TITLE:
2-WAY CLEARTRACT CATHETER, 16 FR X 10 ML, STERILE, 10-PACK

SIZE B	DWG. NO. SLQ-211610SPT	REV A
------------------	---------------------------	----------

DO NOT SCALE DRAWING

SCALE: 1:5

SHEET 1 OF 1



pathway

Pathway LLC

Printed by Aaron Rogers

Printed on 05/25/2023

Local time zone (GMT-08:00) Pacific Time (America/Los_Angeles)

Page Arena > Workspace Items > Bill of Materials > Indented

All information contained in this document is proprietary and confidential.

SLQ-211810SPT

In Design

SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack

Revision »A - ECO-001247 Status: Effective, Status Date: 01/12/2023 01:24:56 PM, Shared

Contains 10 first-level Items, 50 line Items, 41 unique Items, 33 of which are shared.

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE		REF	DES
▼	1	SLQ-8102 rev A	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack	Subassembly	In Des	3	1 EA
	1	SLQ-6002 rev A	Catheter, 18 Fr x 10 mL, Silicone, NY Foley, 2-Way	Fabricated Part	In Des	3	10 EA
	2	SLQ-3000 rev B	Solution, SILQ Polymer	Purchased Component	In Prod	4	80 ML
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1HDX-3000 rev A	Solution, Hydrophilix Polymer	80 ML			
	3	SLQ-5000 rev A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	Packaging	In Prod	3	10 EA
	4	SLQ-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	Packaging	In Prod	4	10 EA
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1HDX-5001 rev A	Pouch, Chevron, 3.5 X 21, 10 EA Amcor TLP-058				
▶	5	SLQ-LB-4008 rev A	Label, Pouch, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	10 EA Max 12 month product shelf life
	6	SLQ-5016 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS, Pre-Printed	Packaging	In Prod	3	1 EA
		SUBSTITUTES					
		RANKITEM NUMBER	ITEM NAME		NOTES		
		1SLQ-5015 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS	1 EA			
		2SLQ-5014 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .024 SBS	1 EA			
▶	7	SLQ-LB-4009 rev A	Label, Box, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	2 EA Max 12 month product shelf life
	8	SLQ-4007 rev A	Artwork, IFU, SILQ 2-Way ClearTract SPT Catheter	Labeling	In Prod	2	1 EA
	9	SLQ-5012 rev A	Box, 20 in. L x 18 in. W x 6 in. D, 200 LB Test	Packaging	In Prod	2	0.1 EA
	10	PWY-5000 rev A	Tape, Packaging, 2 Mil, 2 in x 110 yds, Clear	Packaging	In Prod	2	5.2 IN
▶	11	SLQ-LB-4010 rev A	Label, Shipper Box, SILQ ClearTract SPT Catheter	Labeling	In Prod	2	0.2 EA Max 12 month product shelf life

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	    	REF DES
12 MP-SLQ-0002 rev B	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Manufacturing Procedure	In Prod	2	1 EA
13 FM-SLQ-0001 rev A	Form, Run Sheet for SILQ Automation Machine	Form	In Prod	2	1 EA
14 TM-HDX-0001 rev B	Test Method, UV Light Absorbance, Hydrophilix Catheters	Test Method	In Prod	2	1 each
15 TM-HDX-0002 rev A	Test Method, Rinse Test, Hydrophilix Catheter	Test Method	In Prod	2	1 each
16 MP-SLQ-0001 rev B	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	Manufacturing Procedure	In Prod	2	1 EA
2 WI-0006 rev D	Work Order Process	Work Instruction	In Prod	2	1 each
3 FM-0039 rev F	Form, Work Order	Form	In Prod	2	1 each
4 WI-0014 rev A	Line Clearance	Work Instruction	In Prod	2	1 each
5 MP-SLQ-0003 rev A	Manufacturing Process Flow, SILQ Foley Catheter	Manufacturing Procedure	In Prod	2	1 EA
6 SLQ-0001 rev A	Sterilization Load Shipping Pallet Configuration	Subassembly	In Prod	2	1 EA
7 QP-032 rev F	EO Sterilization Process Control	Quality System Procedure	In Prod	6	1 each
8 WI-0021 rev C	Work Instruction, Pyrogen (LAL) Evaluation Process	Work Instruction	In Prod	2	1 each
9 QP-010 rev F	Final Product Inspection	Quality System Procedure	In Prod	6	1 each
10 FM-0029 rev A	Form, Certificate of Conformance	Form	In Prod	2	1 each

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
1.

ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY SILQ PRIOR TO PRODUCTION.
2.


FINAL PRODUCT INSPECTION PER QP-010.
3.

PRODUCT EXPIRY DATE IS 12 MONTHS FROM THE DATE OF MANUFACTURE.
2. STERILIZATION
1.

PRODUCT IS STERILIZED VIA ETHYLENE OXIDE (ETO) STERILIZATION PER QP-032 CYCLE PWI-7374.
3. TESTING
1.

PYROGEN EVALUATION PER WI-0021 AT ≤ 0.5 EU/ML OR ≤ 20.0 EU/DEVICE

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

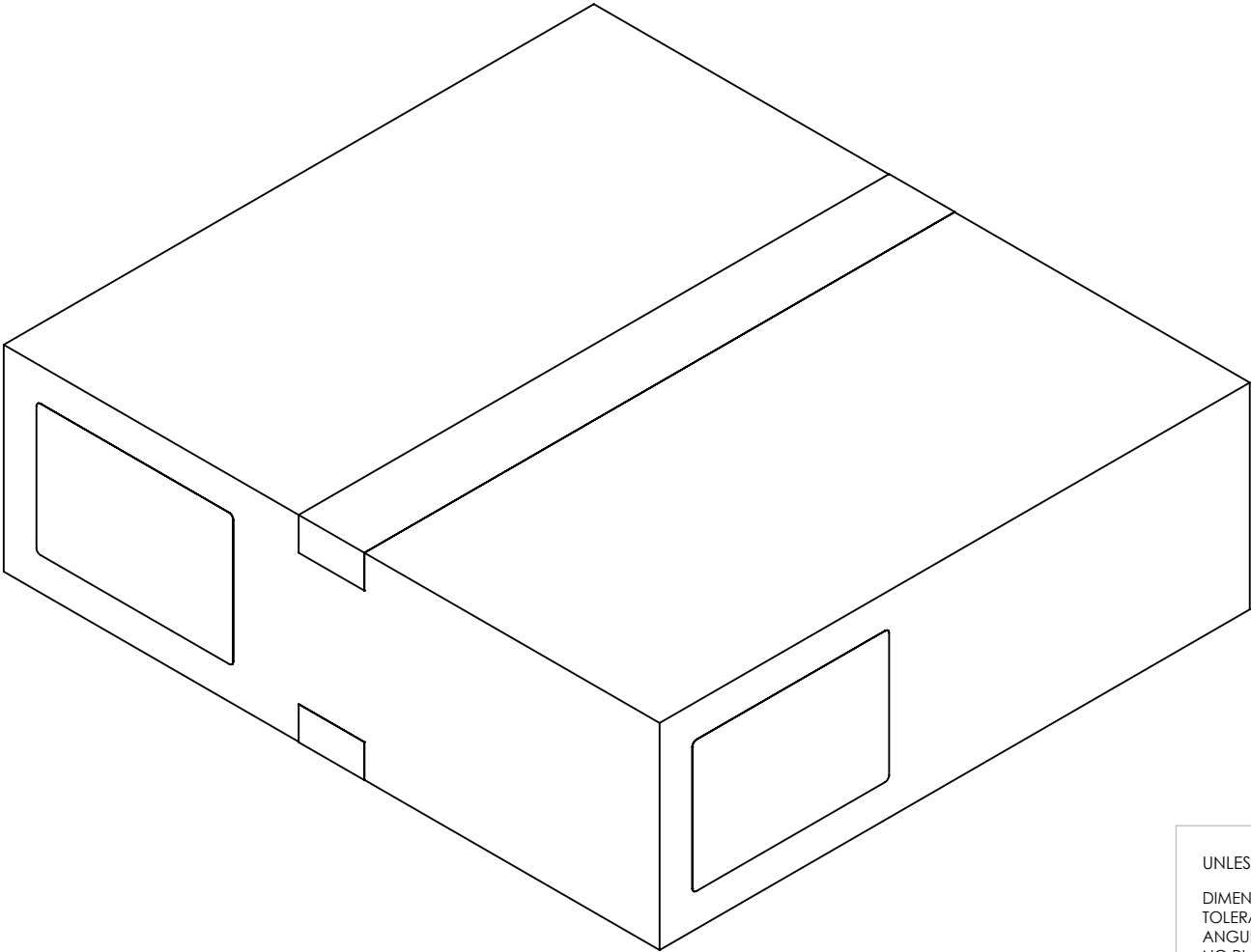


UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH ± 0.2 BEND ± 2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ± 0.15
THREE PLACE DECIMAL ± 0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.



TITLE:
2-WAY CLEARTRACT CATHETER, 18 FR X 10 ML, STERILE, 10-PACK

SIZE	DWG. NO.	REV
B	SLQ-211810SPT	A

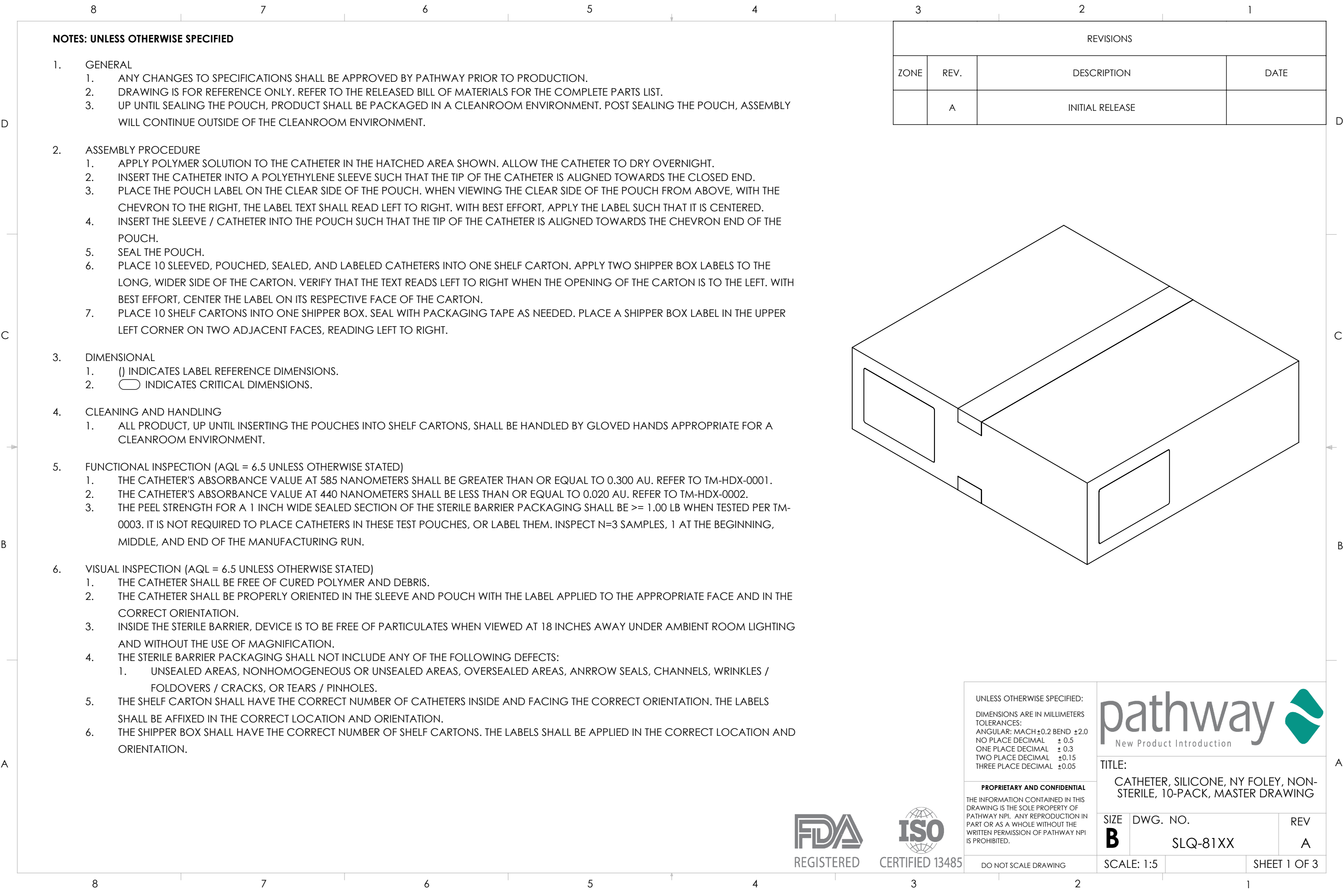

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY PATHWAY PRIOR TO PRODUCTION.

2. DRAWING IS FOR REFERENCE ONLY. REFER TO THE RELEASED BILL OF MATERIALS FOR THE COMPLETE PARTS LIST.

3. UP UNTIL SEALING THE POUCH, PRODUCT SHALL BE PACKAGED IN A CLEANROOM ENVIRONMENT. POST SEALING THE POUCH, ASSEMBLY WILL CONTINUE OUTSIDE OF THE CLEANROOM ENVIRONMENT.
2. ASSEMBLY PROCEDURE
1. APPLY POLYMER SOLUTION TO THE CATHETER IN THE HATCHED AREA SHOWN. ALLOW THE CATHETER TO DRY OVERNIGHT.

2. INSERT THE CATHETER INTO A POLYETHYLENE SLEEVE SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CLOSED END.

3. PLACE THE POUCH LABEL ON THE CLEAR SIDE OF THE POUCH. WHEN VIEWING THE CLEAR SIDE OF THE POUCH FROM ABOVE, WITH THE CHEVRON TO THE RIGHT, THE LABEL TEXT SHALL READ LEFT TO RIGHT. WITH BEST EFFORT, APPLY THE LABEL SUCH THAT IT IS CENTERED.

4. INSERT THE SLEEVE / CATHETER INTO THE POUCH SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CHEVRON END OF THE POUCH.

5. SEAL THE POUCH.

6. PLACE 10 SLEEVED, POUCHED, SEALED, AND LABELED CATHETERS INTO ONE SHELF CARTON. APPLY TWO SHIPPER BOX LABELS TO THE LONG, WIDER SIDE OF THE CARTON. VERIFY THAT THE TEXT READS LEFT TO RIGHT WHEN THE OPENING OF THE CARTON IS TO THE LEFT. WITH BEST EFFORT, CENTER THE LABEL ON ITS RESPECTIVE FACE OF THE CARTON.

7. PLACE 10 SHELF CARTONS INTO ONE SHIPPER BOX. SEAL WITH PACKAGING TAPE AS NEEDED. PLACE A SHIPPER BOX LABEL IN THE UPPER LEFT CORNER ON TWO ADJACENT FACES, READING LEFT TO RIGHT.
3. DIMENSIONAL
1. () INDICATES LABEL REFERENCE DIMENSIONS.

2.  INDICATES CRITICAL DIMENSIONS.
4. CLEANING AND HANDLING
1. ALL PRODUCT, UP UNTIL INSERTING THE POUCHES INTO SHELF CARTONS, SHALL BE HANDLED BY GLOVED HANDS APPROPRIATE FOR A CLEANROOM ENVIRONMENT.
5. FUNCTIONAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
1. THE CATHETER'S ABSORBANCE VALUE AT 585 NANOMETERS SHALL BE GREATER THAN OR EQUAL TO 0.300 AU. REFER TO TM-HDX-0001.

2. THE CATHETER'S ABSORBANCE VALUE AT 440 NANOMETERS SHALL BE LESS THAN OR EQUAL TO 0.020 AU. REFER TO TM-HDX-0002.

3. THE PEEL STRENGTH FOR A 1 INCH WIDE SEALED SECTION OF THE STERILE BARRIER PACKAGING SHALL BE >= 1.00 LB WHEN TESTED PER TM-0003. IT IS NOT REQUIRED TO PLACE CATHETERS IN THESE TEST POUCHES, OR LABEL THEM. INSPECT N=3 SAMPLES, 1 AT THE BEGINNING, MIDDLE, AND END OF THE MANUFACTURING RUN.
6. VISUAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
1. THE CATHETER SHALL BE FREE OF CURED POLYMER AND DEBRIS.

2. THE CATHETER SHALL BE PROPERLY ORIENTED IN THE SLEEVE AND POUCH WITH THE LABEL APPLIED TO THE APPROPRIATE FACE AND IN THE CORRECT ORIENTATION.

3. INSIDE THE STERILE BARRIER, DEVICE IS TO BE FREE OF PARTICULATES WHEN VIEWED AT 18 INCHES AWAY UNDER AMBIENT ROOM LIGHTING AND WITHOUT THE USE OF MAGNIFICATION.


4. THE STERILE BARRIER PACKAGING SHALL NOT INCLUDE ANY OF THE FOLLOWING DEFECTS:

1. UNSEALED AREAS, NONHOMOGENEOUS OR UNSEALED AREAS, OVERSEALED AREAS, ANRRROW SEALS, CHANNELS, WRINKLES / FOLDOVERS / CRACKS, OR TEARS / PINHOLES.

5. THE SHELF CARTON SHALL HAVE THE CORRECT NUMBER OF CATHETERS INSIDE AND FACING THE CORRECT ORIENTATION. THE LABELS SHALL BE AFFIXED IN THE CORRECT LOCATION AND ORIENTATION.

6. THE SHIPPER BOX SHALL HAVE THE CORRECT NUMBER OF SHELF CARTONS. THE LABELS SHALL BE APPLIED IN THE CORRECT LOCATION AND ORIENTATION.

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

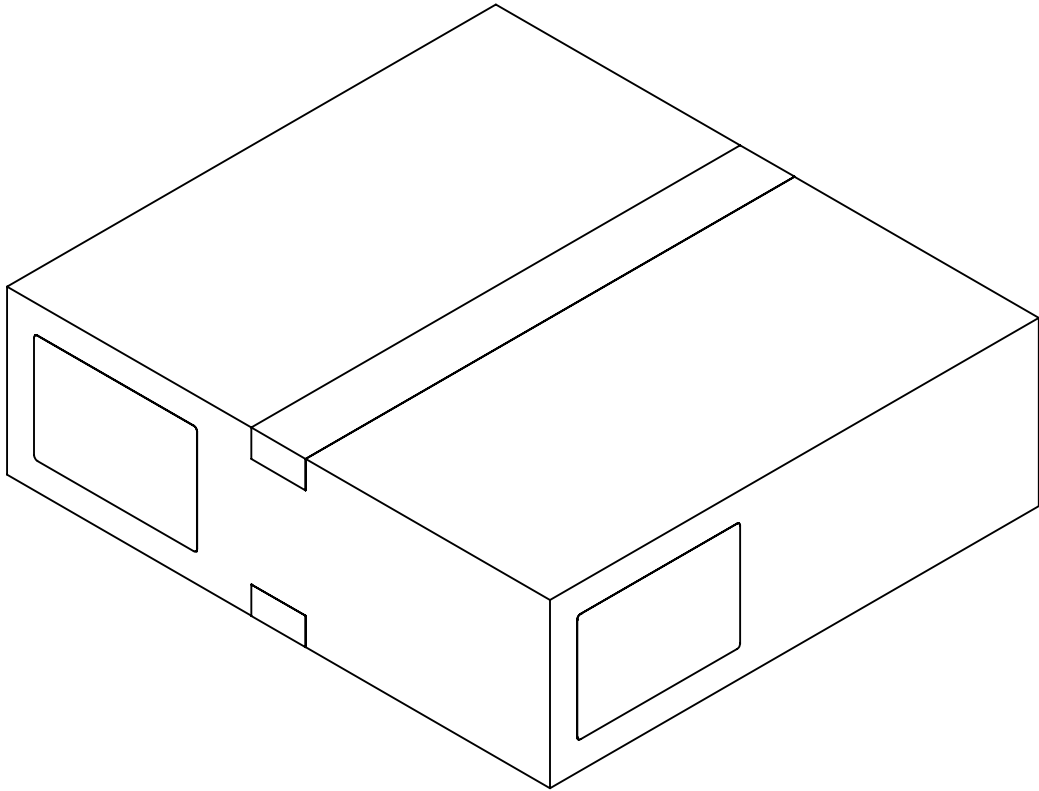

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.

DO NOT SCALE DRAWING


pathway

New Product Introduction

TITLE:
CATHETER, SILICONE, NY FOLEY, NON-STERILE, 10-PACK, MASTER DRAWING


SIZE B	DWG. NO. SLQ-81XX	REV A
SCALE: 1:5		SHEET 1 OF 3



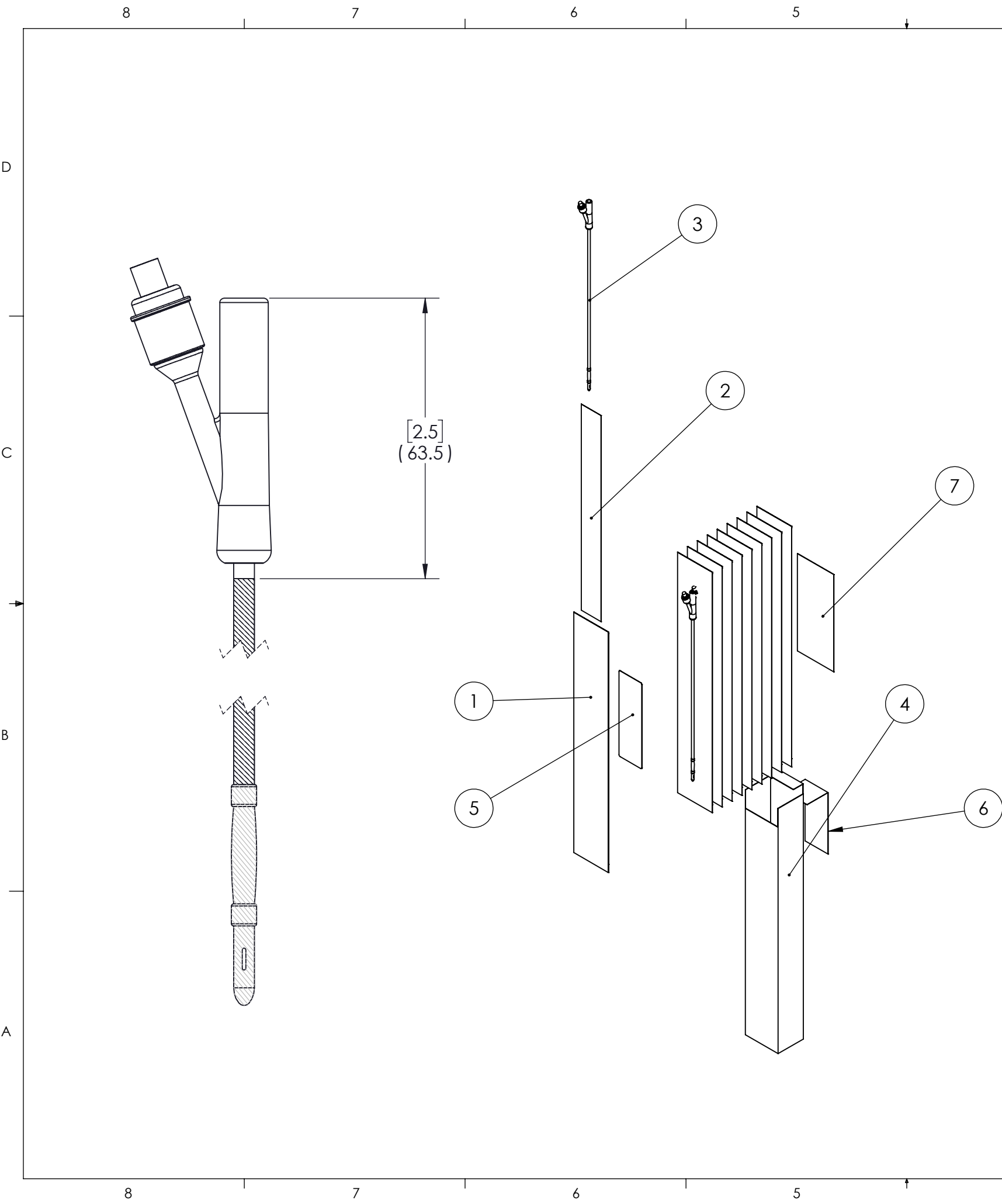

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
1	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	10
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	10
3	SLQ-60XX	SEE TABLE FOR PART SPECIFICATION	10
4	SLQ-5016	SHELF CARTON, 3.25 X 2.5 X 19.75, .018 SBS, PRE-PRINTED	1
5	SLQ-LB-4008	POUCH LABEL, CLEARTRACT SP CATHETER	10
6	SLQ-LB-4009	SHELF CARTON LABEL, CLEARTRACT SP CATHETER	1
7	SLQ-4007	IFU, SILQ 2-WAY CLEARTRACT SP CATHETER	1

PART NUMBER	DESCRIPTION	SLQ-81XX
SLQ-6000	CATHETER, 14 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8100
SLQ-6001	CATHETER, 16 FR X 5 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8101
SLQ-6002	CATHETER, 18 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8102
SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8104

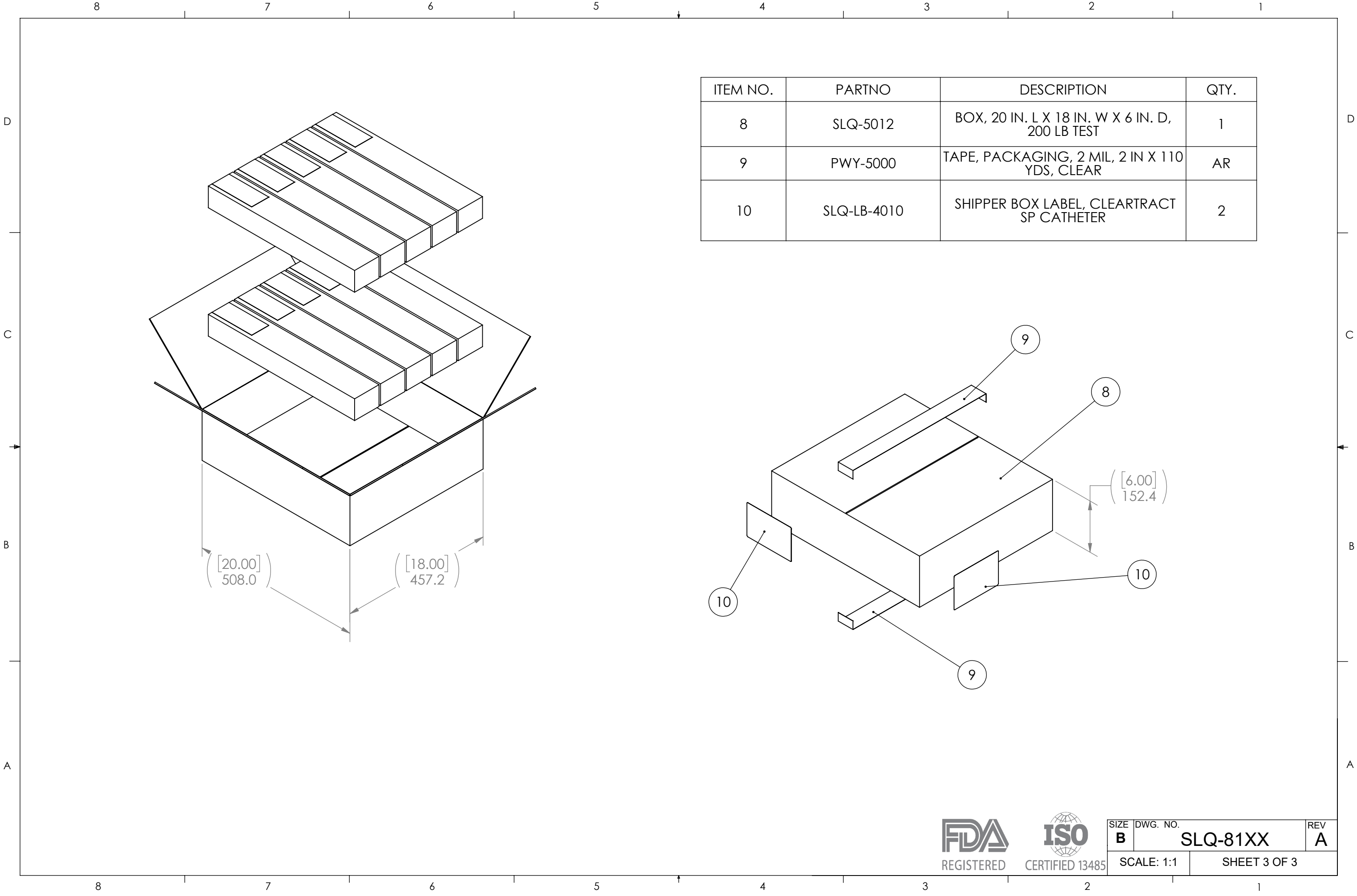

87654321

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL

1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY PATHWAY PRIOR TO PRODUCTION.

2. DRAWING IS FOR REFERENCE ONLY. REFER TO THE RELEASED BILL OF MATERIALS FOR THE COMPLETE PARTS LIST.

3. UP UNTIL SEALING THE POUCH, PRODUCT SHALL BE PACKAGED IN A CLEANROOM ENVIRONMENT. POST SEALING THE POUCH, ASSEMBLY WILL CONTINUE OUTSIDE OF THE CLEANROOM ENVIRONMENT.

2. ASSEMBLY PROCEDURE

1. APPLY POLYMER SOLUTION TO THE CATHETER IN THE HATCHED AREA SHOWN. ALLOW THE CATHETER TO DRY OVERNIGHT.

2. INSERT THE CATHETER INTO A POLYETHYLENE SLEEVE SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CLOSED END.

3. PLACE THE POUCH LABEL ON THE CLEAR SIDE OF THE POUCH. WHEN VIEWING THE CLEAR SIDE OF THE POUCH FROM ABOVE, WITH THE CHEVRON TO THE RIGHT, THE LABEL TEXT SHALL READ LEFT TO RIGHT. WITH BEST EFFORT, APPLY THE LABEL SUCH THAT IT IS CENTERED.

4. INSERT THE SLEEVE / CATHETER INTO THE POUCH SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CHEVRON END OF THE POUCH.

5. SEAL THE POUCH.

6. PLACE 10 SLEEVED, POUCHED, SEALED, AND LABELED CATHETERS INTO ONE SHELF CARTON. APPLY TWO SHIPPER BOX LABELS TO THE LONG, WIDER SIDE OF THE CARTON. VERIFY THAT THE TEXT READS LEFT TO RIGHT WHEN THE OPENING OF THE CARTON IS TO THE LEFT. WITH BEST EFFORT, CENTER THE LABEL ON ITS RESPECTIVE FACE OF THE CARTON.

7. PLACE 10 SHELF CARTONS INTO ONE SHIPPER BOX. SEAL WITH PACKAGING TAPE AS NEEDED. PLACE A SHIPPER BOX LABEL IN THE UPPER LEFT CORNER ON TWO ADJACENT FACES, READING LEFT TO RIGHT.

3. DIMENSIONAL

1. () INDICATES LABEL REFERENCE DIMENSIONS.

2.  INDICATES CRITICAL DIMENSIONS.

4. CLEANING AND HANDLING

1. ALL PRODUCT, UP UNTIL INSERTING THE POUCHES INTO SHELF CARTONS, SHALL BE HANDLED BY GLOVED HANDS APPROPRIATE FOR A CLEANROOM ENVIRONMENT.

5. FUNCTIONAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)

1. THE CATHETER'S ABSORBANCE VALUE AT 585 NANOMETERS SHALL BE GREATER THAN OR EQUAL TO 0.300 AU. REFER TO TM-HDX-0001.

2. THE CATHETER'S ABSORBANCE VALUE AT 440 NANOMETERS SHALL BE LESS THAN OR EQUAL TO 0.020 AU. REFER TO TM-HDX-0002.

3. THE PEEL STRENGTH FOR A 1 INCH WIDE SEALED SECTION OF THE STERILE BARRIER PACKAGING SHALL BE >= 1.00 LB WHEN TESTED PER TM-0003. IT IS NOT REQUIRED TO PLACE CATHETERS IN THESE TEST POUCHES, OR LABEL THEM. INSPECT N=3 SAMPLES, 1 AT THE BEGINNING, MIDDLE, AND END OF THE MANUFACTURING RUN.

6. VISUAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)

1. THE CATHETER SHALL BE FREE OF CURED POLYMER AND DEBRIS.

2. THE CATHETER SHALL BE PROPERLY ORIENTED IN THE SLEEVE AND POUCH WITH THE LABEL APPLIED TO THE APPROPRIATE FACE AND IN THE CORRECT ORIENTATION.

3. INSIDE THE STERILE BARRIER, DEVICE IS TO BE FREE OF PARTICULATES WHEN VIEWED AT 18 INCHES AWAY UNDER AMBIENT ROOM LIGHTING AND WITHOUT THE USE OF MAGNIFICATION.

4. THE STERILE BARRIER PACKAGING SHALL NOT INCLUDE ANY OF THE FOLLOWING DEFECTS:

1. UNSEALED AREAS, NONHOMOGENEOUS OR UNSEALED AREAS, OVERSEALED AREAS, ANRRROW SEALS, CHANNELS, WRINKLES / FOLDOVERS / CRACKS, OR TEARS / PINHOLES.


5. THE SHELF CARTON SHALL HAVE THE CORRECT NUMBER OF CATHETERS INSIDE AND FACING THE CORRECT ORIENTATION. THE LABELS SHALL BE AFFIXED IN THE CORRECT LOCATION AND ORIENTATION.

6. THE SHIPPER BOX SHALL HAVE THE CORRECT NUMBER OF SHELF CARTONS. THE LABELS SHALL BE APPLIED IN THE CORRECT LOCATION AND ORIENTATION.

321

REVISIONS

ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

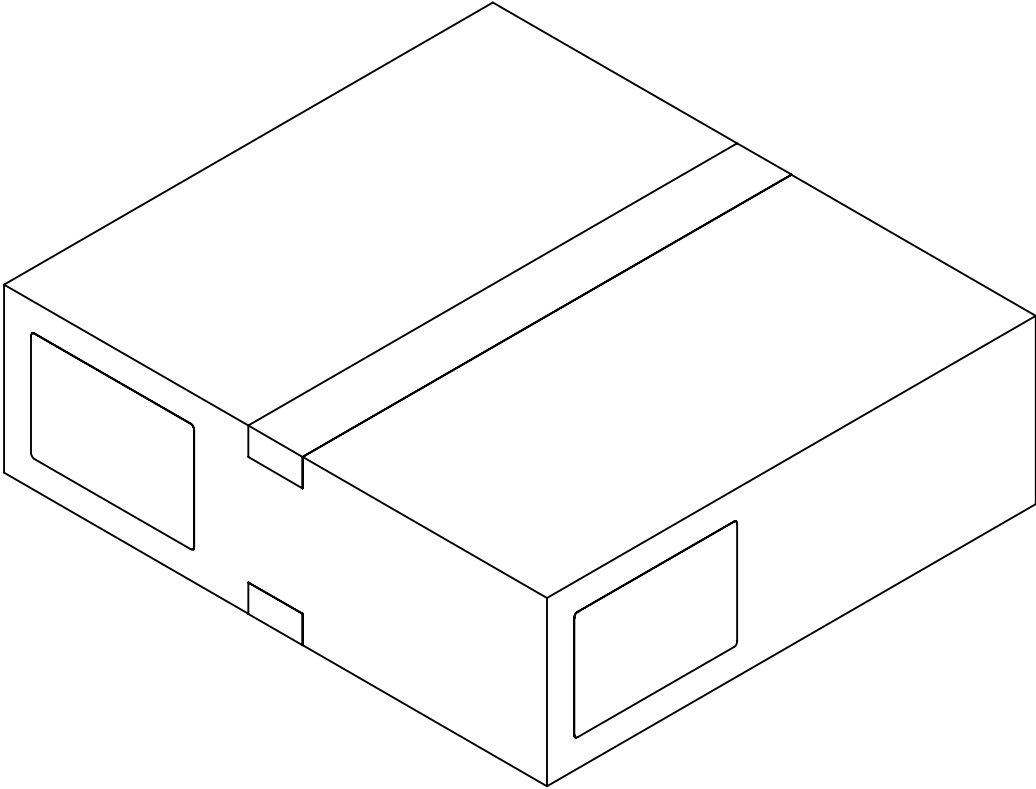

87654321

UNLESS OTHERWISE SPECIFIED:

DIMENSIONS ARE IN MILLIMETERS

TOLERANCES:

ANGULAR: MACH±0.2 BEND ±2.0

NO PLACE DECIMAL ± 0.5

ONE PLACE DECIMAL ± 0.3

TWO PLACE DECIMAL ±0.15

THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL

THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.

DO NOT SCALE DRAWING

pathway

New Product Introduction

TITLE:

CATHETER, SILICONE, NY FOLEY, NON-STERILE, 10-PACK, MASTER DRAWING

SIZE	DWG. NO.	REV
B	SLQ-81XX	A


SCALE: 1:5		SHEET 1 OF 3
------------	--	--------------

FDA

REGISTERED


ISO

CERTIFIED 13485

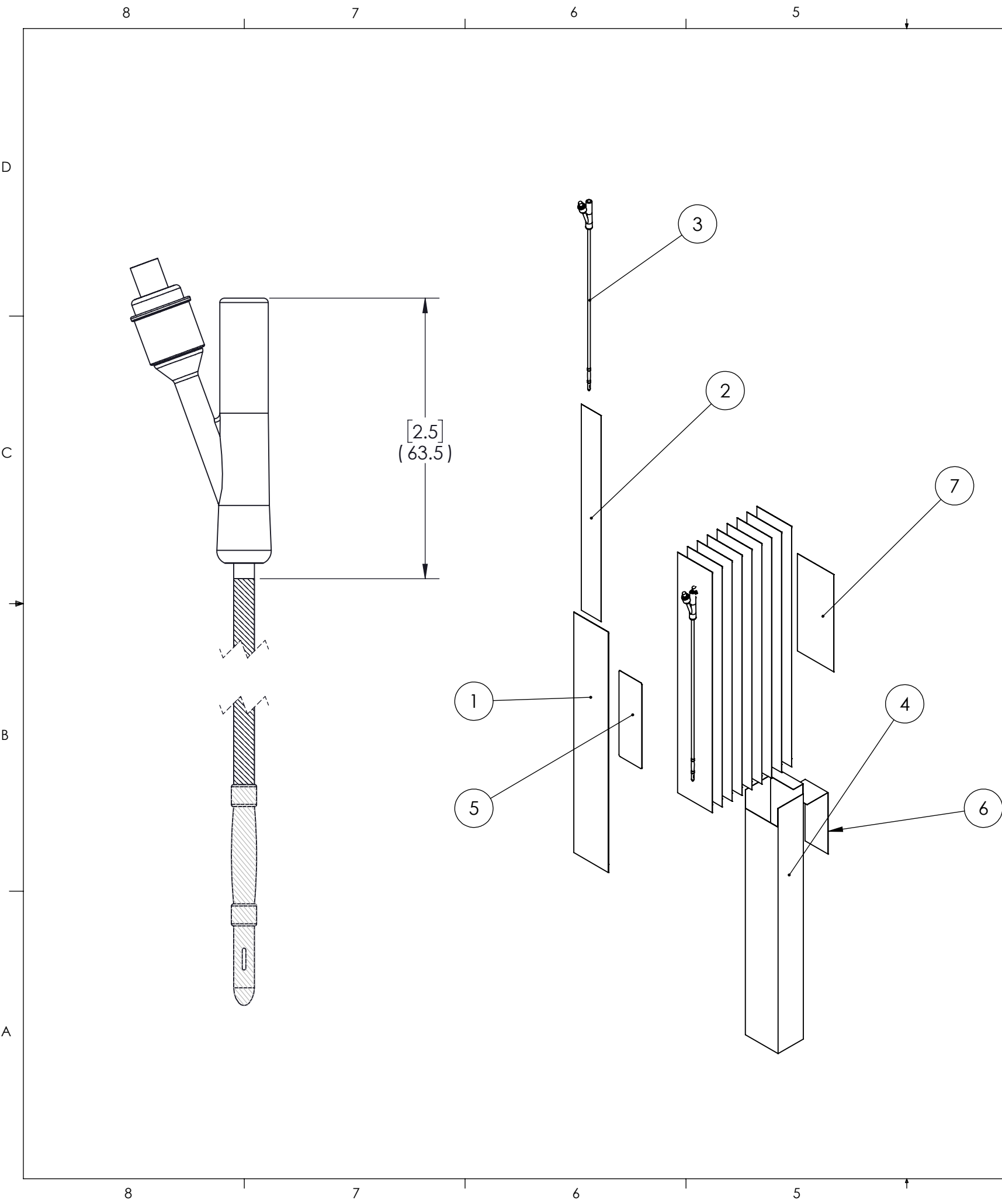

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
1	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	10
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	10
3	SLQ-60XX	SEE TABLE FOR PART SPECIFICATION	10
4	SLQ-5016	SHELF CARTON, 3.25 X 2.5 X 19.75, .018 SBS, PRE-PRINTED	1
5	SLQ-LB-4008	POUCH LABEL, CLEARTRACT SP CATHETER	10
6	SLQ-LB-4009	SHELF CARTON LABEL, CLEARTRACT SP CATHETER	1
7	SLQ-4007	IFU, SILQ 2-WAY CLEARTRACT SP CATHETER	1

PART NUMBER	DESCRIPTION	SLQ-81XX
SLQ-6000	CATHETER, 14 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8100
SLQ-6001	CATHETER, 16 FR X 5 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8101
SLQ-6002	CATHETER, 18 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8102
SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8104

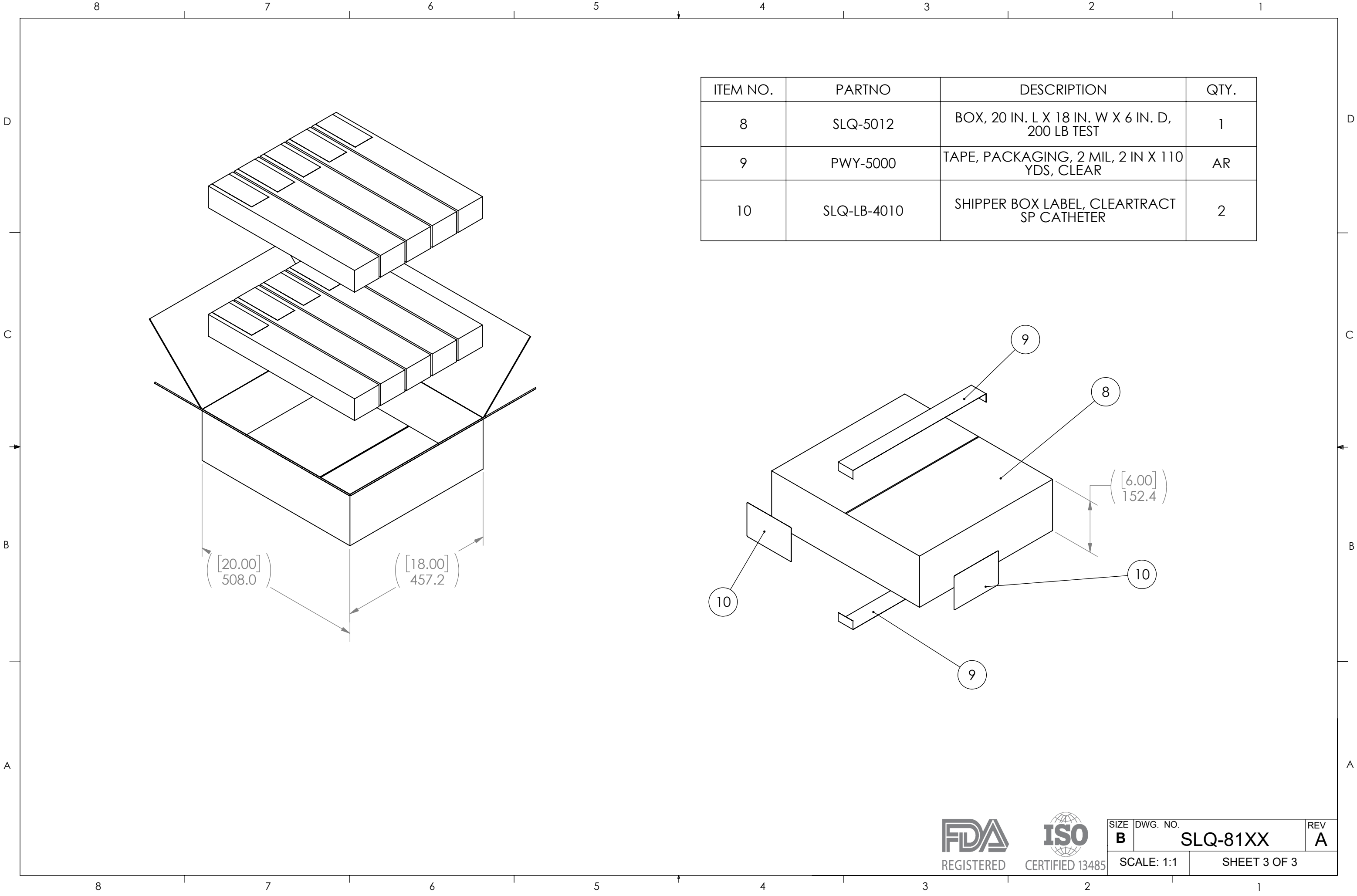

87654321

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL

1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY PATHWAY PRIOR TO PRODUCTION.

2. DRAWING IS FOR REFERENCE ONLY. REFER TO THE RELEASED BILL OF MATERIALS FOR THE COMPLETE PARTS LIST.

3. UP UNTIL SEALING THE POUCH, PRODUCT SHALL BE PACKAGED IN A CLEANROOM ENVIRONMENT. POST SEALING THE POUCH, ASSEMBLY WILL CONTINUE OUTSIDE OF THE CLEANROOM ENVIRONMENT.

2. ASSEMBLY PROCEDURE

1. APPLY POLYMER SOLUTION TO THE CATHETER IN THE HATCHED AREA SHOWN. ALLOW THE CATHETER TO DRY OVERNIGHT.

2. INSERT THE CATHETER INTO A POLYETHYLENE SLEEVE SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CLOSED END.

3. PLACE THE POUCH LABEL ON THE CLEAR SIDE OF THE POUCH. WHEN VIEWING THE CLEAR SIDE OF THE POUCH FROM ABOVE, WITH THE CHEVRON TO THE RIGHT, THE LABEL TEXT SHALL READ LEFT TO RIGHT. WITH BEST EFFORT, APPLY THE LABEL SUCH THAT IT IS CENTERED.

4. INSERT THE SLEEVE / CATHETER INTO THE POUCH SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CHEVRON END OF THE POUCH.

5. SEAL THE POUCH.

6. PLACE 10 SLEEVED, POUCHED, SEALED, AND LABELED CATHETERS INTO ONE SHELF CARTON. APPLY TWO SHIPPER BOX LABELS TO THE LONG, WIDER SIDE OF THE CARTON. VERIFY THAT THE TEXT READS LEFT TO RIGHT WHEN THE OPENING OF THE CARTON IS TO THE LEFT. WITH BEST EFFORT, CENTER THE LABEL ON ITS RESPECTIVE FACE OF THE CARTON.

7. PLACE 10 SHELF CARTONS INTO ONE SHIPPER BOX. SEAL WITH PACKAGING TAPE AS NEEDED. PLACE A SHIPPER BOX LABEL IN THE UPPER LEFT CORNER ON TWO ADJACENT FACES, READING LEFT TO RIGHT.

3. DIMENSIONAL

1. () INDICATES LABEL REFERENCE DIMENSIONS.

2.  INDICATES CRITICAL DIMENSIONS.

4. CLEANING AND HANDLING

1. ALL PRODUCT, UP UNTIL INSERTING THE POUCHES INTO SHELF CARTONS, SHALL BE HANDLED BY GLOVED HANDS APPROPRIATE FOR A CLEANROOM ENVIRONMENT.

5. FUNCTIONAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)

1. THE CATHETER'S ABSORBANCE VALUE AT 585 NANOMETERS SHALL BE GREATER THAN OR EQUAL TO 0.300 AU. REFER TO TM-HDX-0001.

2. THE CATHETER'S ABSORBANCE VALUE AT 440 NANOMETERS SHALL BE LESS THAN OR EQUAL TO 0.020 AU. REFER TO TM-HDX-0002.

3. THE PEEL STRENGTH FOR A 1 INCH WIDE SEALED SECTION OF THE STERILE BARRIER PACKAGING SHALL BE >= 1.00 LB WHEN TESTED PER TM-0003. IT IS NOT REQUIRED TO PLACE CATHETERS IN THESE TEST POUCHES, OR LABEL THEM. INSPECT N=3 SAMPLES, 1 AT THE BEGINNING, MIDDLE, AND END OF THE MANUFACTURING RUN.

6. VISUAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)

1. THE CATHETER SHALL BE FREE OF CURED POLYMER AND DEBRIS.

2. THE CATHETER SHALL BE PROPERLY ORIENTED IN THE SLEEVE AND POUCH WITH THE LABEL APPLIED TO THE APPROPRIATE FACE AND IN THE CORRECT ORIENTATION.

3. INSIDE THE STERILE BARRIER, DEVICE IS TO BE FREE OF PARTICULATES WHEN VIEWED AT 18 INCHES AWAY UNDER AMBIENT ROOM LIGHTING AND WITHOUT THE USE OF MAGNIFICATION.

4. THE STERILE BARRIER PACKAGING SHALL NOT INCLUDE ANY OF THE FOLLOWING DEFECTS:

1. UNSEALED AREAS, NONHOMOGENEOUS OR UNSEALED AREAS, OVERSEALED AREAS, ANRRROW SEALS, CHANNELS, WRINKLES / FOLDOVERS / CRACKS, OR TEARS / PINHOLES.


5. THE SHELF CARTON SHALL HAVE THE CORRECT NUMBER OF CATHETERS INSIDE AND FACING THE CORRECT ORIENTATION. THE LABELS SHALL BE AFFIXED IN THE CORRECT LOCATION AND ORIENTATION.

6. THE SHIPPER BOX SHALL HAVE THE CORRECT NUMBER OF SHELF CARTONS. THE LABELS SHALL BE APPLIED IN THE CORRECT LOCATION AND ORIENTATION.

321

REVISIONS

ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

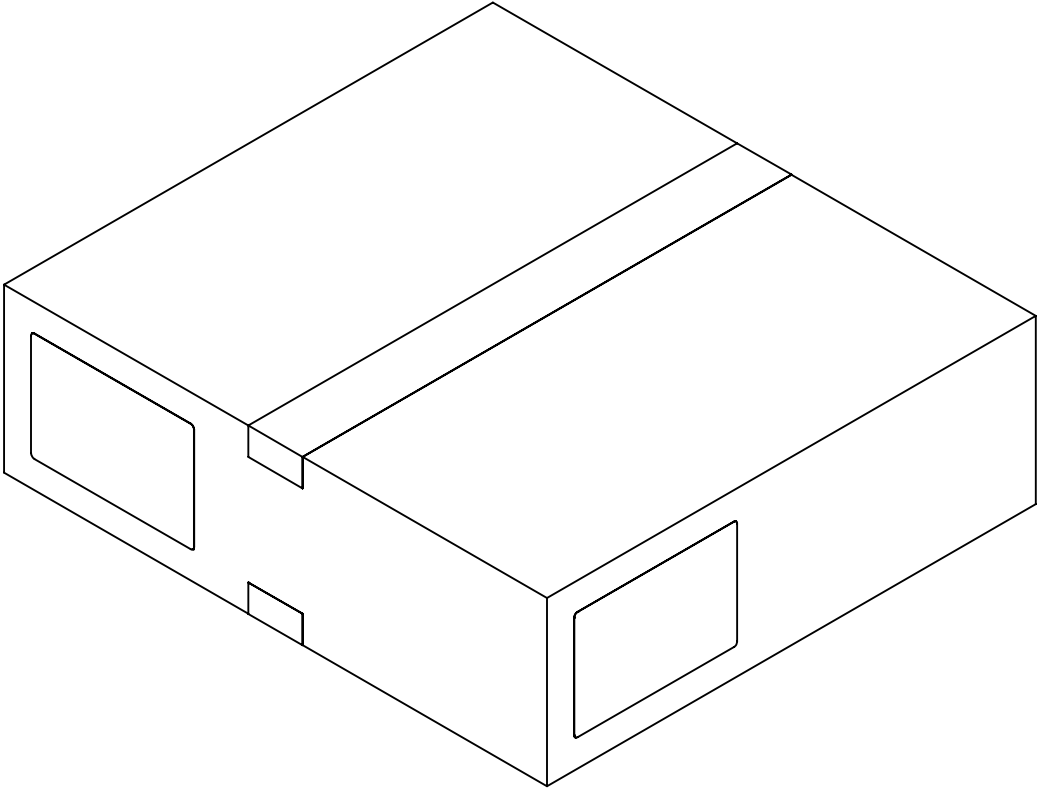

87654321

UNLESS OTHERWISE SPECIFIED:

DIMENSIONS ARE IN MILLIMETERS

TOLERANCES:

ANGULAR: MACH±0.2 BEND ±2.0

NO PLACE DECIMAL ± 0.5

ONE PLACE DECIMAL ± 0.3

TWO PLACE DECIMAL ±0.15

THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL

THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NPI. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NPI IS PROHIBITED.

DO NOT SCALE DRAWING

pathway

New Product Introduction

TITLE:

CATHETER, SILICONE, NY FOLEY, NON-STERILE, 10-PACK, MASTER DRAWING

SIZE	DWG. NO.	REV
B	SLQ-81XX	A


SCALE: 1:5		SHEET 1 OF 3
------------	--	--------------

FDA

REGISTERED


ISO

CERTIFIED 13485

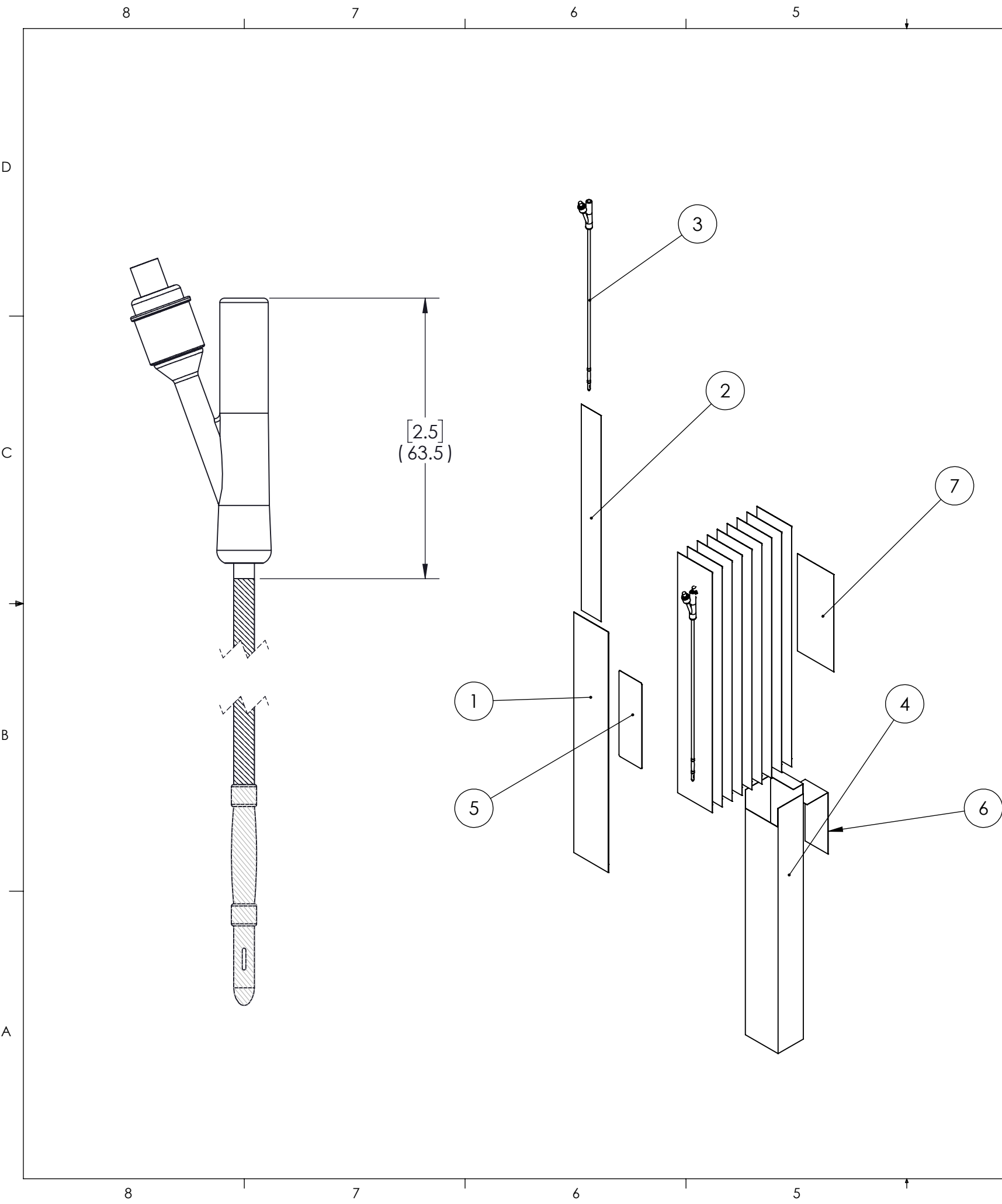

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
1	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	10
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	10
3	SLQ-60XX	SEE TABLE FOR PART SPECIFICATION	10
4	SLQ-5016	SHELF CARTON, 3.25 X 2.5 X 19.75, .018 SBS, PRE-PRINTED	1
5	SLQ-LB-4008	POUCH LABEL, CLEARTRACT SP CATHETER	10
6	SLQ-LB-4009	SHELF CARTON LABEL, CLEARTRACT SP CATHETER	1
7	SLQ-4007	IFU, SILQ 2-WAY CLEARTRACT SP CATHETER	1

PART NUMBER	DESCRIPTION	SLQ-81XX
SLQ-6000	CATHETER, 14 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8100
SLQ-6001	CATHETER, 16 FR X 5 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8101
SLQ-6002	CATHETER, 18 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8102
SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8104

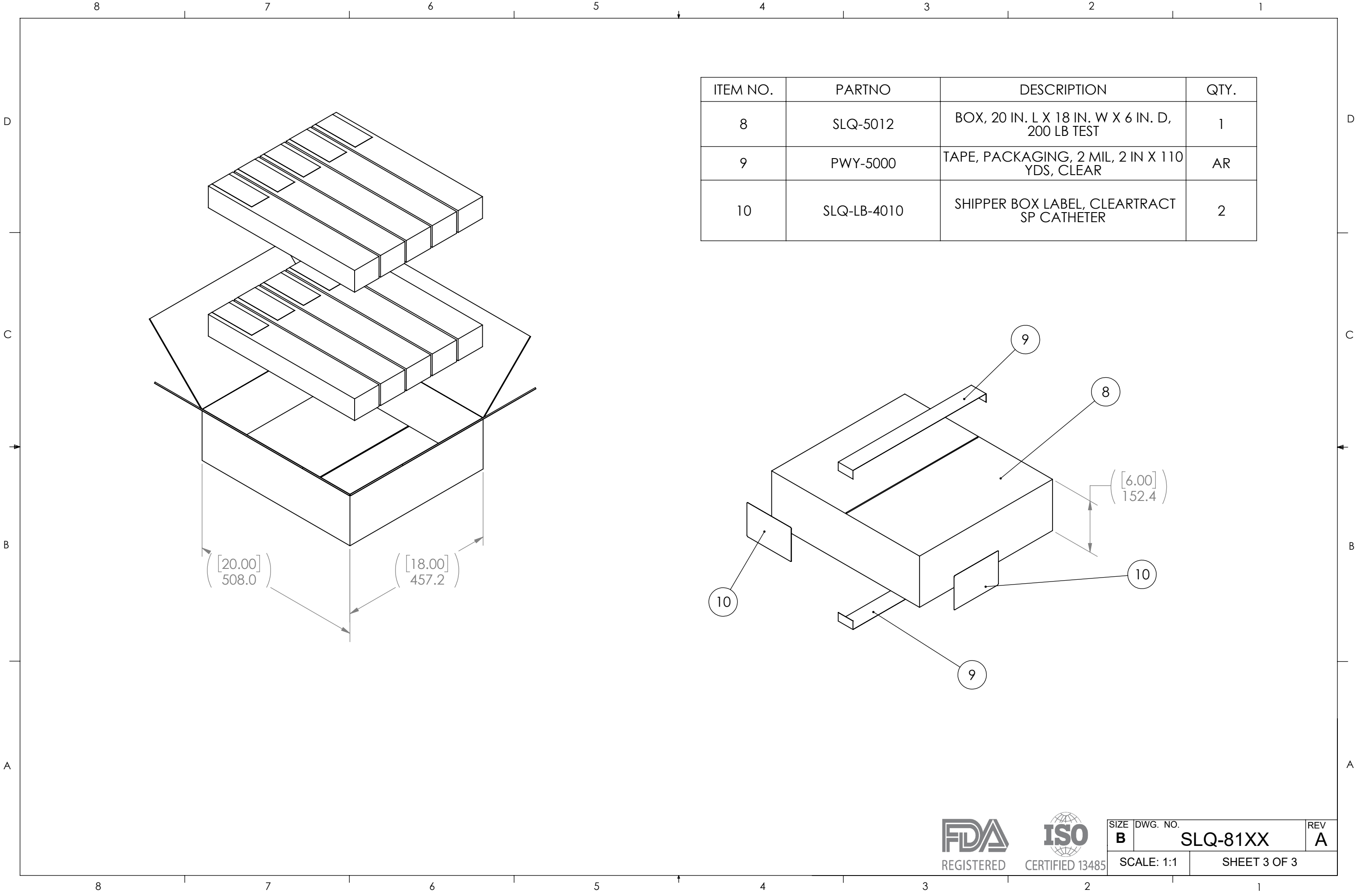

D

- C

C

A


FDA
REGISTERED

pathway
New Product Introduction

A


DO NOT SCALE DRAWING

SCALE: 1:5	SHEET 1 OF 3
------------	--------------

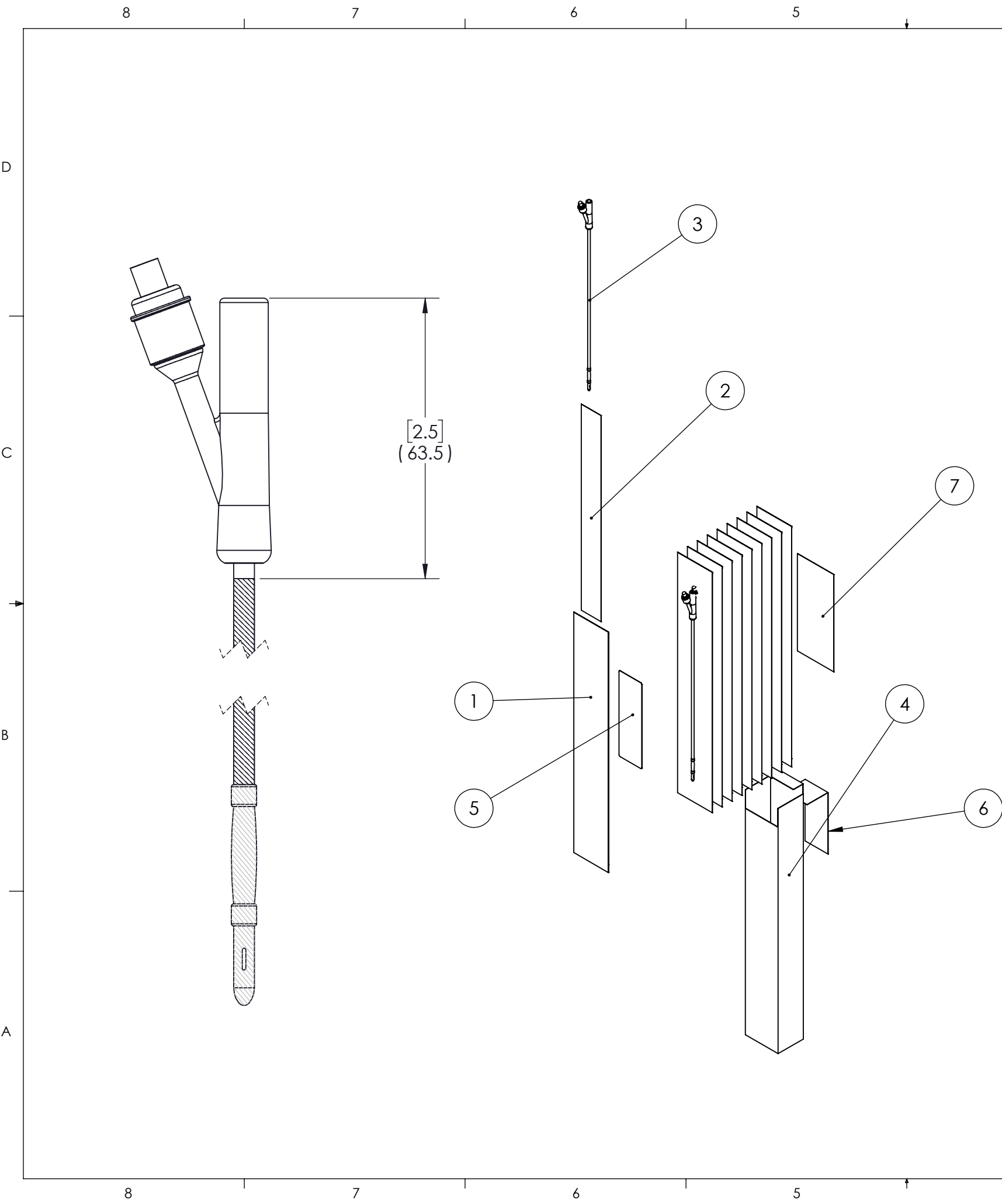

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
1	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	10
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	10
3	SLQ-60XX	SEE TABLE FOR PART SPECIFICATION	10
4	SLQ-5016	SHELF CARTON, 3.25 X 2.5 X 19.75, .018 SBS, PRE-PRINTED	1
5	SLQ-LB-4008	POUCH LABEL, CLEARTRACT SP CATHETER	10
6	SLQ-LB-4009	SHELF CARTON LABEL, CLEARTRACT SP CATHETER	1
7	SLQ-4007	IFU, SILQ 2-WAY CLEARTRACT SP CATHETER	1

PART NUMBER	DESCRIPTION	SLQ-81XX
SLQ-6000	CATHETER, 14 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8100
SLQ-6001	CATHETER, 16 FR X 5 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8101
SLQ-6002	CATHETER, 18 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8102
SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8104

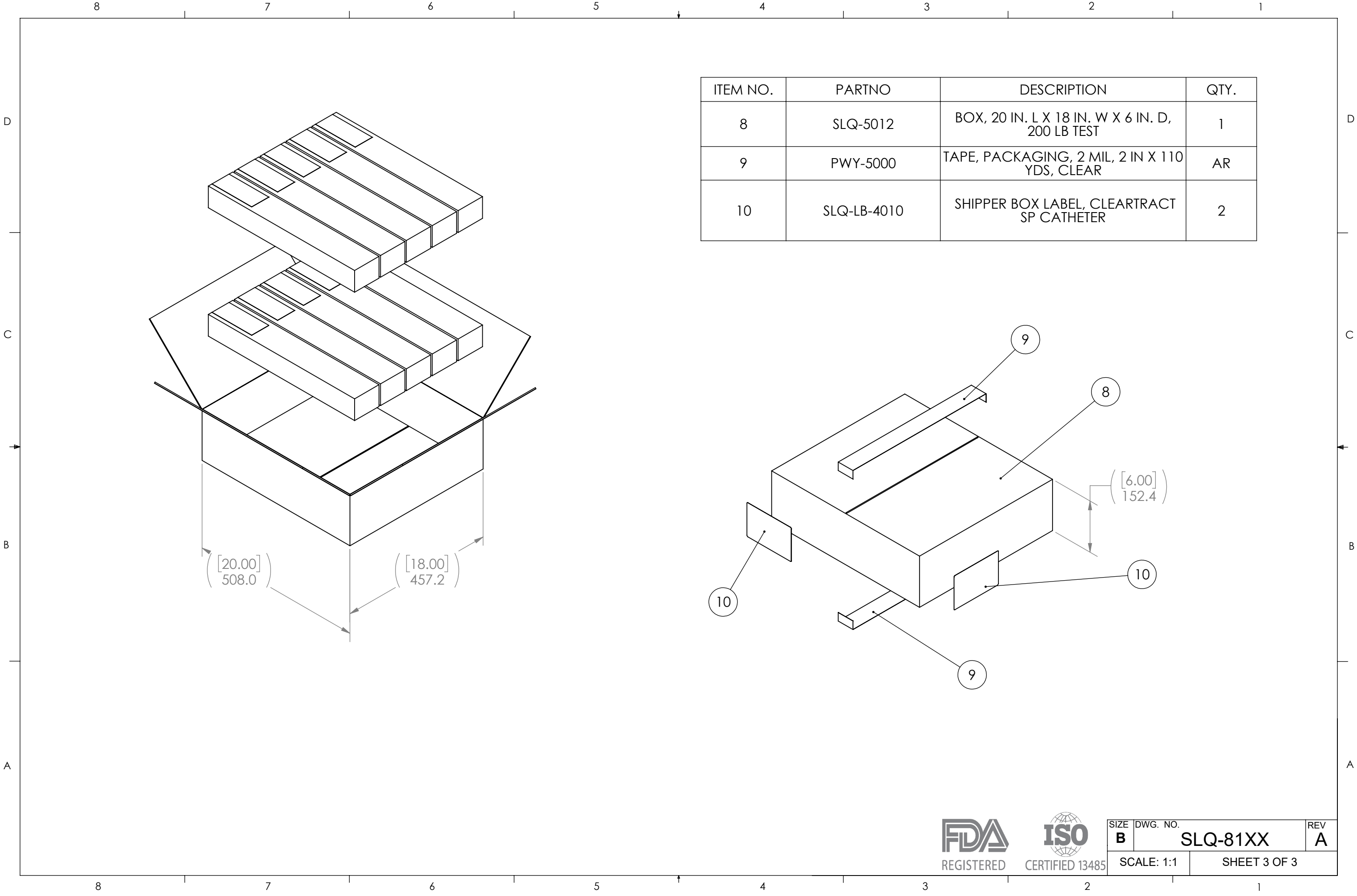

#6 DEVICE MASTER RECORD

INSPECTION FORMS

Item Number	Item Name
SLQ-8100 rev A	SILQ 2-way Cleartract SPT 14 Fr x 10 ml Catheter, Non-Sterile, 10-Pack
SLQ-8101 rev A	SILQ 2-way Cleartract SPT 16 Fr x 5 ml Catheter, Non-Sterile, 10-Pack
SLQ-8103 rev A	SILQ 2-way Cleartract SPT 16 Fr x 10 ml Catheter, Non-Sterile, 10-Pack
SLQ-8104 rev A	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack

Date Inspected		Part No.		Rev.	Description				
				A					
Lot No.		Lot QTY		Approved Supplier(s)					
				Pathway					
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	A	A	A	A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	N/A	0.300	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	N/A	0.020	N/A	N/A	N/A	N/A
	UOM	EA	EA	AU	AU	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Visual	Spectrophotometer	Spectrophotometer	Visual	Visual	Visual	Visual
	Asset No.	N/A	N/A			N/A	N/A	N/A	N/A
	Cal Due Date	N/A	N/A			N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size								
Sample No.		Variable/ Attribute Data							
1									
2									
3									
4									
5									
6									
7									
8									
9									
10									
11									
12									
13									
14									

15									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
Comments									
Disposition:							NCMR No:		
Inspected By:							Date:		
Verified By:							Date:		

Date Inspected	Part No.	Rev.	Description						
		A							
Lot No.	Lot QTY	ES-0059-01 Emplex Model MPS6300-M Continuous Band Sealer - Validated Range: Temperature = 142 °C - 144 °C, Speed = 175 in/min., Pressure = 89 psi ES-0099-01 ACCU-SEAL Model 5300-15-B, Impulse Sealer - Validated Range: Temperature = 270 °F - 280 °F, Dwell = 4 sec., Pressure = 55 psi							
Specification	TM-0003								
Revision									
PRODUCT PREP FOR TEST									
* Seal one sample pouch P/N SLQ-5001 at the start, middle and end of each shift.									
* Cut two 1.00 inch wide x 2.00 inch long strips from the manufacturer sealed end of the pouch.									
* Execute peel strength test per TM-0003 and record results.									
Specification	≥1.0	≥1.0	142 - 144 °C	175 in/min	89 psi	270 - 280 °F	4.0	55.000	
Lower Spec Limit (LSL)	1.0	1.0	142.0	175.0	89	270.0	4.0	55	
Upper Spec Limit (USL)	N/A	N/A	144.0	175.0	89	280.0	4.0	55	
UOM	lb	lb	°C	in/min	psi	°F	seconds	psi	
Measuring Equipment Type	Force Gauge	Force Gauge	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	
Asset No.		0		0	0		0	0	
Cal Due Date		1/0/1900		1/0/1900	1/0/1900		1/0/1900	1/0/1900	
Sample ID	Top Strip	Bottom Strip	1	1	1	1	1	1	
Sample No.	Variable/ Attribute Data								
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									

2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	0.000	0.000	0.000	0.00	0.00	0.00	0.000	0.000	
Max	0.000	0.000	0.000	0.00	0.00	0.00	0.000	0.000	
Mean	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	0.00	#DIV/0!	#DIV/0!	#DIV/0!	
StDev	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	
Comments									

Disposition:		NCMR No:	
Inspected By:		Date:	

Date Inspected		Part No.		Rev.	Description				
				A					
Lot No.		Lot QTY		Approved Supplier(s)					
				Pathway					
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	A	A	A	A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	N/A	0.300	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	N/A	0.020	N/A	N/A	N/A	N/A
	UOM	EA	EA	AU	AU	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Visual	Spectrophotometer	Spectrophotometer	Visual	Visual	Visual	Visual
	Asset No.	N/A	N/A			N/A	N/A	N/A	N/A
	Cal Due Date	N/A	N/A			N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size								
Sample No.		Variable/ Attribute Data							
1									
2									
3									
4									
5									
6									
7									
8									
9									
10									
11									
12									
13									
14									

15									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
Comments									
Disposition:							NCMR No:		
Inspected By:							Date:		
Verified By:							Date:		

Date Inspected	Part No.	Rev.	Description						
		A							
Lot No.	Lot QTY	ES-0059-01 Emplex Model MPS6300-M Continuous Band Sealer - Validated Range: Temperature = 142 °C - 144 °C, Speed = 175 in/min., Pressure = 89 psi ES-0099-01 ACCU-SEAL Model 5300-15-B, Impulse Sealer - Validated Range: Temperature = 270 °F - 280 °F, Dwell = 4 sec., Pressure = 55 psi							
Specification	TM-0003								
Revision									
PRODUCT PREP FOR TEST									
* Seal one sample pouch P/N SLQ-5001 at the start, middle and end of each shift.									
* Cut two 1.00 inch wide x 2.00 inch long strips from the manufacturer sealed end of the pouch.									
* Execute peel strength test per TM-0003 and record results.									
Specification	≥1.0	≥1.0	142 - 144 °C	175 in/min	89 psi	270 - 280 °F	4.0	55.000	
Lower Spec Limit (LSL)	1.0	1.0	142.0	175.0	89	270.0	4.0	55	
Upper Spec Limit (USL)	N/A	N/A	144.0	175.0	89	280.0	4.0	55	
UOM	lb	lb	°C	in/min	psi	°F	seconds	psi	
Measuring Equipment Type	Force Gauge	Force Gauge	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	
Asset No.		0		0	0		0	0	
Cal Due Date		1/0/1900		1/0/1900	1/0/1900		1/0/1900	1/0/1900	
Sample ID	Top Strip	Bottom Strip	1	1	1	1	1	1	
Sample No.	Variable/ Attribute Data								
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									

2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	0.000	0.000	0.000	0.00	0.00	0.00	0.000	0.000	
Max	0.000	0.000	0.000	0.00	0.00	0.00	0.000	0.000	
Mean	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	0.00	#DIV/0!	#DIV/0!	#DIV/0!	
StDev	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	
Comments									

Disposition:		NCMR No:	
Inspected By:		Date:	

Date Inspected		Part No.		Rev.	Description				
				A					
Lot No.		Lot QTY		Approved Supplier(s)					
				Pathway					
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	A	A	A	A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	N/A	0.300	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	N/A	0.020	N/A	N/A	N/A	N/A
	UOM	EA	EA	AU	AU	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Visual	Spectrophotometer	Spectrophotometer	Visual	Visual	Visual	Visual
	Asset No.	N/A	N/A			N/A	N/A	N/A	N/A
	Cal Due Date	N/A	N/A			N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size								
Sample No.		Variable/ Attribute Data							
1									
2									
3									
4									
5									
6									
7									
8									
9									
10									
11									
12									
13									
14									

15									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
Comments									
Disposition:							NCMR No:		
Inspected By:							Date:		
Verified By:							Date:		

Date Inspected	Part No.	Rev.	Description						
		A							
Lot No.	Lot QTY	ES-0059-01 Emplex Model MPS6300-M Continuous Band Sealer - Validated Range: Temperature = 142 °C - 144 °C, Speed = 175 in/min., Pressure = 89 psi ES-0099-01 ACCU-SEAL Model 5300-15-B, Impulse Sealer - Validated Range: Temperature = 270 °F - 280 °F, Dwell = 4 sec., Pressure = 55 psi							
Specification	TM-0003								
Revision									
PRODUCT PREP FOR TEST									
* Seal one sample pouch P/N SLQ-5001 at the start, middle and end of each shift.									
* Cut two 1.00 inch wide x 2.00 inch long strips from the manufacturer sealed end of the pouch.									
* Execute peel strength test per TM-0003 and record results.									
Specification	≥1.0	≥1.0	142 - 144 °C	175 in/min	89 psi	270 - 280 °F	4.0	55.000	
Lower Spec Limit (LSL)	1.0	1.0	142.0	175.0	89	270.0	4.0	55	
Upper Spec Limit (USL)	N/A	N/A	144.0	175.0	89	280.0	4.0	55	
UOM	lb	lb	°C	in/min	psi	°F	seconds	psi	
Measuring Equipment Type	Force Gauge	Force Gauge	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	
Asset No.		0		0	0		0	0	
Cal Due Date		1/0/1900		1/0/1900	1/0/1900		1/0/1900	1/0/1900	
Sample ID	Top Strip	Bottom Strip	1	1	1	1	1	1	
Sample No.	Variable/ Attribute Data								
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									

2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	0.000	0.000	0.000	0.00	0.00	0.00	0.000	0.000	
Max	0.000	0.000	0.000	0.00	0.00	0.00	0.000	0.000	
Mean	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	0.00	#DIV/0!	#DIV/0!	#DIV/0!	
StDev	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	
Comments									

Disposition:		NCMR No:	
Inspected By:		Date:	

Date Inspected		Part No.		Rev.	Description				
				A					
Lot No.		Lot QTY		Approved Supplier(s)					
				Pathway					
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	A	A	A	A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	N/A	0.300	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	N/A	0.020	N/A	N/A	N/A	N/A
	UOM	EA	EA	AU	AU	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Visual	Spectrophotometer	Spectrophotometer	Visual	Visual	Visual	Visual
	Asset No.	N/A	N/A			N/A	N/A	N/A	N/A
	Cal Due Date	N/A	N/A			N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size								
Sample No.		Variable/ Attribute Data							
1									
2									
3									
4									
5									
6									
7									
8									
9									
10									
11									
12									
13									
14									

15									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
Comments									
Disposition:							NCMR No:		
Inspected By:							Date:		
Verified By:							Date:		

Date Inspected	Part No.	Rev.	Description						
		A							
Lot No.	Lot QTY	ES-0059-01 Emplex Model MPS6300-M Continuous Band Sealer - Validated Range: Temperature = 142 °C - 144 °C, Speed = 175 in/min., Pressure = 89 psi ES-0099-01 ACCU-SEAL Model 5300-15-B, Impulse Sealer - Validated Range: Temperature = 270 °F - 280 °F, Dwell = 4 sec., Pressure = 55 psi							
Specification	TM-0003								
Revision									
PRODUCT PREP FOR TEST									
* Seal one sample pouch P/N SLQ-5001 at the start, middle and end of each shift.									
* Cut two 1.00 inch wide x 2.00 inch long strips from the manufacturer sealed end of the pouch.									
* Execute peel strength test per TM-0003 and record results.									
Specification	≥1.0	≥1.0	142 - 144 °C	175 in/min	89 psi	270 - 280 °F	4.0	55.000	
Lower Spec Limit (LSL)	1.0	1.0	142.0	175.0	89	270.0	4.0	55	
Upper Spec Limit (USL)	N/A	N/A	144.0	175.0	89	280.0	4.0	55	
UOM	lb	lb	°C	in/min	psi	°F	seconds	psi	
Measuring Equipment Type	Force Gauge	Force Gauge	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	
Asset No.		0		0	0		0	0	
Cal Due Date		1/0/1900		1/0/1900	1/0/1900		1/0/1900	1/0/1900	
Sample ID	Top Strip	Bottom Strip	1	1	1	1	1	1	
Sample No.	Variable/ Attribute Data								
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									

2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
1: Start (Top Strip)									
2: Start (Bottom Strip)									
3: Middle (Top Strip)									
4: Middle (Bottom Strip)									
5: End (Top Strip)									
6: End (Bottom Strip)									
Result Summary									
Pass QTY									
Fail QTY									
Result (P/F)									
Statistical Summary									
Min	0.000	0.000	0.000	0.00	0.00	0.00	0.000	0.000	
Max	0.000	0.000	0.000	0.00	0.00	0.00	0.000	0.000	
Mean	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	0.00	#DIV/0!	#DIV/0!	#DIV/0!	
StDev	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	#DIV/0!	
Comments									

Disposition:		NCMR No:	
Inspected By:		Date:	

#7 RISK MANAGEMENT REPORT (PFMEA)

Item Number	Item Name
RM-0021 rev D	Risk Management Report, PFMEA, SILQ ClearTract Foley Catheter

1.0 Purpose

- 1.1 The purpose of this document is to identify the associated risk with the manufacturing processes of the SILQ ClearTract Foley Catheter.

2.0 Scope

- 2.1 The product numbers associated with the SILQ ClearTract Foley Catheter is indicated in the following table:

Table 1: SILQ ClearTract Catheter Products	
Part No.	Description
HDX-1410	2-way 100% Silicone ClearTract Catheter, 14 Fr/10ml
HDX-1605	2-way 100% Silicone ClearTract Catheter, 16 Fr/5ml
HDX-1810	2-way 100% Silicone ClearTract Catheter, 18 Fr/10ml
HDX-160510	2-way 100% Silicone ClearTract Catheter, 16 Fr/5-10ml
SLQ-1610	2-Way 100% Silicone ClearTract Catheter, 16Fr/10ml
SLQ-211410	2-way 100% Silicone ClearTract Catheter, 14 Fr/10ml
SLQ-211605	2-way 100% Silicone ClearTract Catheter, 16 Fr/5ml
SLQ-211610	2-way 100% Silicone ClearTract Catheter, 16 Fr/10ml
SLQ-211810	2-way 100% Silicone ClearTract Catheter, 18 Fr/10ml
SLQ-211410	2-way 100% Silicone ClearTract SPT Catheter, 14 Fr/10ml
SLQ-211605	2-way 100% Silicone ClearTract SPT Catheter, 16 Fr/5ml
SLQ-211610	2-way 100% Silicone ClearTract SPT Catheter, 16 Fr/10ml
SLQ-211810	2-way 100% Silicone ClearTract SPT Catheter, 18 Fr/10ml

3.0 References

- 3.1 QP-004 – Risk Management
- 3.2 RM-0018 – Risk Management Plan, SILQ ClearTract Foley Catheter

4.0 Risk Tables Definitions

- 4.1 # – Enumerate each Description with a sequential number
- 4.2 Description – Describes the operation that is being evaluated for risk
- 4.3 Potential Failure Mode (FM) – The potential consequences of the failure
- 4.4 Severity (S) – Measure of the possible consequences of a hazard.

- 4.5 Potential Effect of Failure – All conceivable effects of the potential failure mode
- 4.6 Potential Cause of Failure – All conceivable causes of the potential failure mode
- 4.7 Risk Controls - Process through which decisions are reached and protective measures are implemented for reducing risk to, or maintaining risk within specified levels.
 - 4.7.1 Risk Controls¹ – Initial risk controls that are incorporated.
 - 4.7.2 Risk Controls² – Secondary risk controls, also known as “Mitigation”.
- 4.8 Verification – Evidence source of linkage that initial (i.e., Risk Controls¹) risk controls have implemented, verified, or validated.
- 4.9 Probability – Likelihood of the hazard occurring
 - 4.9.1 Probability¹ (P¹) – The probability assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls¹)
 - 4.9.2 Probability² (P²) – The probability assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls²)
- 4.10 Detection - The ability to screen out a hazard through verification activities (e.g., inspection).
 - 4.10.1 Detection¹ (D¹) – The detection assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls¹)
 - 4.10.2 Detection² (D²) – The detection assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls²)
- 4.11 Frequency – A ‘Frequency Score’ will be calculated using the ratings from probability and detection components. The frequency equation is the following: *Frequency = $\sqrt{Px D}$* **Note:** *The value calculated should be rounded up.*
 - 4.11.1 Frequency¹ (F¹) – The frequency score calculated using Risk Controls¹ and Detection¹
 - 4.11.2 Frequency² (F²) – The frequency score calculated using Risk Controls² and Detection²
- 4.12 Risk Rating – Residual risk determination based on the concept of ‘As Far As Possible’ (AFAP) approach. The AFAP approach provides three (3) possible outcomes, which are 1) the residual risk exceeds the risk acceptability (RED), 2) the residual risk is acceptable to the risk acceptability criteria (GREEN), and 3) the residual risk is between the two (2) states indicated and it can be an option to lower the risk to the lowest practical level without economic considerations (YELLOW).
 - 4.12.1 Risk Rating¹ (R¹) – This calculation is post the implementation of initial risk controls (i.e., Risk Controls¹)
 - 4.12.2 Risk Rating² (R²) – This calculation is post the implementation of the secondary risk controls (i.e., Risk Controls²)

5.0 Risk Acceptability Determination

- 5.1 The criteria of risk acceptability will be determined with a Risk Priority Number (RPN) and on the concept of 'As Far As Possible (AFAP)' approach. Refer to the Risk Rating definition in §4.0.
- 5.2 There are three (3) components used in determining risk acceptability, which are Severity, Probability, and Detection. Refer to the following tables:

Table 2: Severity Ratings	
Index	Description
1	No patient effect
2	Minor injury that resolves without treatment
3	Moderate injury that may require treatment
4	Serious injury requiring medical intervention
5	Life threatening injury, disability, death

Table 3: Probability Ratings	
Index	Description
1	Improbable - <0.001%
2	Remote – 0.001% to .01%
3	Occasional – 0.01% to 0.1%
4	Likely – 0.1% to 1%
5	Frequent - >1%

Table 4: Detection Ratings	
Index	Description
1	Very High – Current controls are almost certain to detect the failure mode.
2	High – High likelihood that current controls will detect the failure mode.
3	Moderate – Moderate likelihood that current controls will detect the failure mode.
4	Low – Low likelihood that current controls will detect the failure mode.
5	Remote – Remote likelihood that current controls will detect the failure mode.

5.3 Risk acceptability will be determined per the following risk index:

Table 5: Risk Rating Table						
Severity	5	5	10	15	20	25
	4	4	8	12	16	20
	3	3	6	9	12	15
	2	2	4	6	8	10
	1	1	2	3	4	5
=Square Root (Probability X Detection)						

Note: The purpose of using the square root of the probability and detection scores is to deemphasize these two (2) components, which places more emphasize on the severity component.

6.0 Risk Assessment Tables

Table 6: SILQ ClearTract Foley Catheter																
#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
1	Receiving Inspection / Quality control	Accept Non-Conforming Material	2	Procedural Complication	Operator Error	Process Instructions	QP-009	1	1	1	2	N/A	N/A	N/A	N/A	N/A
			1		Equipment Error	Calibration Program	QP-020	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Preventative Maintenance Program	QP-015									
		Reject Conforming Material	1	None	Operator Error	Process Instructions	QP-009	1	1	1	1	N/A	N/A	N/A	N/A	N/A
			1		Equipment Error	Calibration Program	QP-020	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Preventative Maintenance Program	QP-015									
		Incorrect / Illegible / No Identification	1	Potential traceability issue	Supplier Error	Purchasing Process	QP-013	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Supplier Qualification	QP-005									
		Damaged product	1	Non-Functional Device	Supplier Error – Incorrect	Purchasing Process	QP-013	1	1	1	1	N/A	N/A	N/A	N/A	N/A

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
					Packaging	Supplier Qualification	QP-005									
					Shipper Error – Mishandling	Process Instructions	QP-009									
		Contamination	2	Cleanroom contamination	Supplier Error	Purchasing Process	QP-013	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Supplier Qualification	QP-005									
					Shipper Error – Mishandling	Process Instructions	QP-009									
		Documentation Error / Incomplete / No CoC/CoA	1	None	Supplier Error	Purchasing Process	QP-013	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Supplier Qualification	QP-005									
2	Line Clearance / Area Preparation	Incorrect Material	1	Non-Functional Device	Operator Error	Process Instructions Process Settings/Line Clearance Form DMR	WI-0014 FM-0039	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Material Expired	1	Non-Functional Device	Operator Error	DMR	FM-0039	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Improper Material Handling	1	Damaged Materials	Operator Error	Process Instructions	WI-0014	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Out of Calibration Equipment	3	Non-Functional Device	Operator Error	Process Instructions	WI-0014	1	1	1	3	N/A	N/A	N/A	N/A	N/A
						Equipment Qualification	IQ Report (TR)									
		Untrained Operators	1	Compromised Device	Inadequate Training	Process Instructions	WI-0014	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Training Program	QP-0017									



Process Failure Mode Effects Analysis (FMEA) – Risk Management Report

SILQ ClearTract Foley Catheter

Doc No.: RM-0021

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: D

Page 6 of 10

Table 6: SILQ ClearTract Foley Catheter

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
3	Catheter tube forming / curing	Incomplete / Malformed Tube	2	Non-Functional Device / Poor flow	Operator Error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	2	2	4	N/A	N/A	N/A	N/A	N/A
4	Balloon forming / dipping	Incomplete / Malformed Balloon	3	Non-Functional Device / Balloon Burst	Operator Error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	2	2	6	N/A	N/A	N/A	N/A	N/A
5	Valve installation	Loose valve / Valve disconnects	1	Non-Functional Device	Operator error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	1	1	1	N/A	N/A	N/A	N/A	N/A
6	Catheter Size Printing	Incorrect print information	2	Urethral damage / over inflation	Operator error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	1	1	2	N/A	N/A	N/A	N/A	N/A
		Illegible	2		Equipment	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	2	1	2	6	N/A	N/A	N/A	N/A	N/A
7	Labeling	Incorrect print information	1	Potential traceability issue	Operator error	Process instruction	MP-0016	2	1	2	2	Label inspection per IP-0002	1	1	1	1
		Illegible	1	Potential traceability issue	Equipment error	Process instruction	MP-0016	2	1	2	2	Label inspection per IP-0002	1	1	1	1



Process Failure Mode Effects Analysis (FMEA) – Risk Management Report

SILQ ClearTract Foley Catheter

Doc No.: RM-0021

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: D

Page 7 of 10

Table 6: SILQ ClearTract Foley Catheter

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
8	Surface Treatment	Incomplete / No Treatment	1	No effect	Equipment error	Process Validation	TR-HDX-0007 (MVP)	1	1	1	1	N/A	N/A	N/A	N/A	N/A
						Start-up Verification Process Instructions	MP-SLQ-0002 TM-HDX-0001 TM-HDX-0002									
						In-Process Inspection	Inspection Plans (Form) HDX-8XXX SLQ-8XXX									
					Operator error	Process instruction	MP-HDX-0004									
						In-Process Inspection	Inspection Plans (Form) HDX-8XXX SLQ-8XXX									
9	Apply Label to Pouch	Misaligned or wrinkled	1	Cosmetic issue	Operator Error	Process Instructions	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	N/A
10	Pouch Loading	Incorrect product configuration	1	Non-Functional Device	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	2	2	2	4	N/A	N/A	N/A	N/A	N/A
11	Tray Sealing (Sterile Barrier)	Weak Seal	4	Breach in sterile barrier	Equipment error	Process validation Inspection process	TR-HDX-0007 (MVP) Inspection Plans (Form) HDX-8XXX SLQ-8XXX	1	1	1	4	N/A	N/A	N/A	N/A	N/A
					Operator error – incorrect sealing parameters	Process instruction Inspection process	TR-0074 Inspection Plans (Form) HDX-8XXX SLQ-8XXX	1	1	1	4	N/A	N/A	N/A	N/A	N/A
					Contamination	Environmental controls	QP-018	1	1	1	4	N/A	N/A	N/A	N/A	N/A

Table 6: SILO ClearTract Foley Catheter

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
12	Assemble Box	Incorrect / Weak box	1	Damage product	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	N/A
	Apply Label to Box	Misaligned or wrinkled	1	Cosmetic issue	Operator Error	Process Instructions	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	N/A
13	Packaging	Incorrect packaging configuration	2	Damage product	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	2	2	2	4	Packaging inspection per IP-0002	1	1	1	2
14	Sterilization	Product not processed / mix-up	4	Potential sterility issue	Operator error	Final Product Inspection	QP-010	1	1	1	4	N/A	N/A	N/A	N/A	N/A
		Cycle Process Excursion	4	Potential sterility issue	Operator error	Supplier;s SOP	WI-00863 GWI-08-004-09	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1
						Supplier Qualification Program	QP-005									
					Equipment error	Process Validation	TR 0075, TR-HDX-0009	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1
						EO Program	QP-032									
		No BI's placed	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1
						Supplier's SOP	WI-00939									
		Incorrect BI's used	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1
						Supplier's SOP	WI-00939									
		Improper Material Handling with Lab	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	1
Supplier's SOP	WI-00939															
Supplier Qualification Program	QP-005															



Process Failure Mode Effects Analysis (FMEA) – Risk Management Report

SILQ ClearTract Foley Catheter

Doc No.: RM-0021

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: D

Page 9 of 10

Table 6: SILQ ClearTract Foley Catheter

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
15	Final Product Inspection	Incomplete DHR	1	Potential traceability issue	Operator error	Process instruction	WI-0006 QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Incorrect DHR	1	Potential traceability issue	Operator error	Process instruction	WI-0006 QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		No CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Incomplete CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A
		Incorrect CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	N/A

	Process Failure Mode Effects Analysis (FMEA) – Risk Management Report SILQ ClearTract Foley Catheter							Doc No.: RM-0021
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>							Rev.: D
								Page 10 of 10

7.0 Risk/ Benefit Assessment

7.1 The following table indicates the risk/benefit analysis for the residual risks that were greater than the risk acceptability limit after all applicable controls/mitigations were applied:

Table 7: Risk/Benefit Assessment for Residual Risk above Acceptability Limit for Process Hazard Assessment									
Table No.	Process Description	Section No.	Potential Failure Mode	Potential Effect of Failure	Potential Cause of Failure	S	F ²	Risk/Benefit Analysis	Does the benefit outweigh the residual risk?
N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A

7.2 There was no residual risk identified that was above the risk acceptability limit per §5.3.

8.0 Conclusion

8.1 Risks have been reduced as low as possible by the following methods:

- Design the product with inherent safety
- Appropriate testing to ensure design criteria were met
- Protective measures/ mechanisms to the product of the device. (critical-to-quality)
- Manufacturing and quality control measures
- Provide information for safety

Based on a review of the risk-management process for the SILQ ClearTract Foley Catheter described herein, the risk management plan based on ISO 14971:2019 has been implemented and the individual, cumulative and overall residual risk have been reduced as low as possible regardless of the magnitude and the benefit of using this device will outweigh risks posed by the device during usage.

9.0 Change History

Revision	Change Description
A	Initial release
B	Add HDX-160510, general updates, add verification references, MP-HDX-0000 to MP-HDX-0004
C	Add SLQ-1610, format updates to tables in §5.0, add sterilization process #14, update verification references
D	Update Scope Table 1, add reference to RM Plan §3.2, add 4 FMs to #1 Receiving Inspection in Table 6

**#8 MARKET CLEARANCE /
K222118 CLEARANCE LETTER**



December 1, 2022

SILQ Technologies, Corp.
% Aaron Rogers
Director of Regulatory and Quality
Pathway LLC
8779 Cottonwood Ave, Suite 105
Santee, California 92071

Re: K222118
Trade/Device Name: 2-Way 100% Silicone Cleartract Catheter
Regulation Number: 21 CFR 876.5130
Regulation Name: Urological Catheter and Accessories
Regulatory Class: II
Product Code: EZL
Dated: October 25, 2022
Received: October 25, 2022

Dear Aaron Rogers:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. Although this letter refers to your product as a device, please be aware that some cleared products may instead be combination products. The 510(k) Premarket Notification Database located at <https://www.accessdata.fda.gov/scripts/cdrh/cfdocs/cfpmn/pmn.cfm> identifies combination product submissions. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration. Please note: CDRH does not evaluate information related to contract liability warranties. We remind you, however, that device labeling must be truthful and not misleading.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's

requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); medical device reporting (reporting of medical device-related adverse events) (21 CFR 803) for devices or postmarketing safety reporting (21 CFR 4, Subpart B) for combination products (see <https://www.fda.gov/combination-products/guidance-regulatory-information/postmarketing-safety-reporting-combination-products>); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820) for devices or current good manufacturing practices (21 CFR 4, Subpart A) for combination products; and, if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR Part 807.97). For questions regarding the reporting of adverse events under the MDR regulation (21 CFR Part 803), please go to <https://www.fda.gov/medical-devices/medical-device-safety/medical-device-reporting-mdr-how-report-medical-device-problems>.

For comprehensive regulatory information about medical devices and radiation-emitting products, including information about labeling regulations, please see Device Advice (<https://www.fda.gov/medical-devices/device-advice-comprehensive-regulatory-assistance>) and CDRH Learn (<https://www.fda.gov/training-and-continuing-education/cdrh-learn>). Additionally, you may contact the Division of Industry and Consumer Education (DICE) to ask a question about a specific regulatory topic. See the DICE website (<https://www.fda.gov/medical-devices/device-advice-comprehensive-regulatory-assistance/contact-us-division-industry-and-consumer-education-dice>) for more information or contact DICE by email (DICE@fda.hhs.gov) or phone (1-800-638-2041 or 301-796-7100).

Sincerely,

Jessica Nguyen, Ph.D.
Assistant Director
DHT3B: Division of Reproductive,
Gynecology and Urology Devices
OHT3: Office of GastroRenal, ObGyn,
General Hospital and Urology Devices
Office of Product Evaluation and Quality
Center for Devices and Radiological Health

Enclosure