

Pathway Innovation LLC.	Template, Design Change Assessment		Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: A
			Page 1 of 2

Project No. / Name	00053 / SILQ ClearTract Foley Catheter Suprapubic Indications for Use				
Product Design Specification / Rev	PS-0006 / D			Date	05-11-2023
Change Type	Product: ☒ Line Extension ☐ Change	Manufacturing: ☐ Change	Rationale	This is a line extension to create new product models with expanded indications for use that includes suprapubic placement of the catheter. New SPT models will include 14Fr10ml, 16Fr5ml, 16Fr10ml, 18Fr10ml.	

Rationale = Describe the justification for the change type and explain the drivers behind the change.

The following table provides an analysis of the potential impact proposed design changes may have with the design requirements described in affected product design specification.

TABLE 1: Change and Requirement Identification						
ID	Type	Component	Change Description	Qualified Part / Rev	Proposed Part / Rev	Requirements Affected
1	A/B/C	IFU	Intended Use changes from: Urological catheter is intended for drainage of the urinary tract. Catheterization is accomplished by inserting the catheter through the urethra and into the bladder. Intended population is adults and pediatrics Intended Use changes to: Urological catheter is intended for drainage of the urinary tract. Catheterization is accomplished by inserting the catheter through the urethra and into the bladder, but also suprapubic placement or by nephrostomy. Intended population is adults and pediatrics	SLQ-4006 / B	SLQ-4007 / A	#1 – Biocompatibility: New endpoint identified (Genotoxicity) #2,3,4 - Labeling
2	A/B/C	Labels / IFU	Update Manufacturer address from: Manufactured for: SILQ Technologies Corp. 12100 Wilshire Blvd, Suite 800 Los Angeles, CA 90025 949-735-2614 Update Manufacturer address to: Manufactured by: SILQ Technologies Corp. 323 Sunny Isles Blvd., 7th Floor Sunny Isles Beach, FL 33160 (424) 309-8523	SLQ-4000 / A SLQ-4001 / A SLQ-4006 / A	SLQ-LBL-4008 / A SLQ-LBL-4009 / A SLQ-4007 / A	#2,3,4 – Labeling

Type: A = Form, B = Fit, C = Function, Performance

The following table provides an impact assessment for the proposed design changes described in Table 1.

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Internal / External Resources	NAMSA – biocompatibility testing
Products Affected	SLQ-211410SPT, SLQ-211605SPT, SLQ-211610SPT, SLQ-211810SPT
Customers Affected	US Market
Distribution Channels	US Market

Pathway Innovation LLC.	Template, Design Change Assessment	Doc No.: T-007
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from.</i>	Rev.: A
		Page 2 of 2

TABLE 2: Change Assessment	
Assessment Topic	Assessment / Rationale
Physical/ Functional (e.g., Bench Testing) Assessment	There are no physical or functional changes to the product. PS-0006 will be updated with reference to design output documents for new model configurations.
Shelf-Life Assessment	All materials that construct the 2-way Catheter have been qualified with a twelve (12) month shelf life. No further action is required.
Dimensional (e.g., 1 st Article) Assessment	There are no dimensional changes to the product. No dimensional verification is required. No further action is required.
Sterilization Assessment	There is no impact to the sterilization validation performed because there were no sterile barrier or pallet configuration changes that affect sterility. No further action is required.
Biocompatibility Assessment	Expanding the indications for use will require genotoxicity testing. The biocompatibility summary report will be amended with additional testing results.
Packaging Assessment	There is no impact to the validated packaging configuration. No further action is required.
Regulatory Assessment	The guidance document “Deciding when to submit a 510(k) for a change to an existing device (1500054)” was used in this assessment, which the attached decision trees were completed. A 510(k) submission is required.
Quality System Assessment	No changes required to the QMS. No further action is required.
Risk Assessment	New risks for modified indications for use and route of administration to be documented and addressed by SILQ in their RMF. SILQ RMF will be updated with reference to new verification activities. RMF is the responsibility of SILQ. No further action is required by Pathway.
Manufacturing Assessment	No changes are required to the manufacturing process. No further action is required.

Contains Nonbinding Recommendations

Note that the first question is whether the change is being made with the intent to significantly improve the safety or effectiveness of the device, for example, to significantly improve clinical outcomes, to mitigate a known risk, in response to adverse events, etc. (Figure 1 – Main Flowchart). If so, the change likely could significantly affect safety or effectiveness and submission of a new 510(k) is likely required. If not, you should continue to follow the logic scheme shown in Figure 1, below.

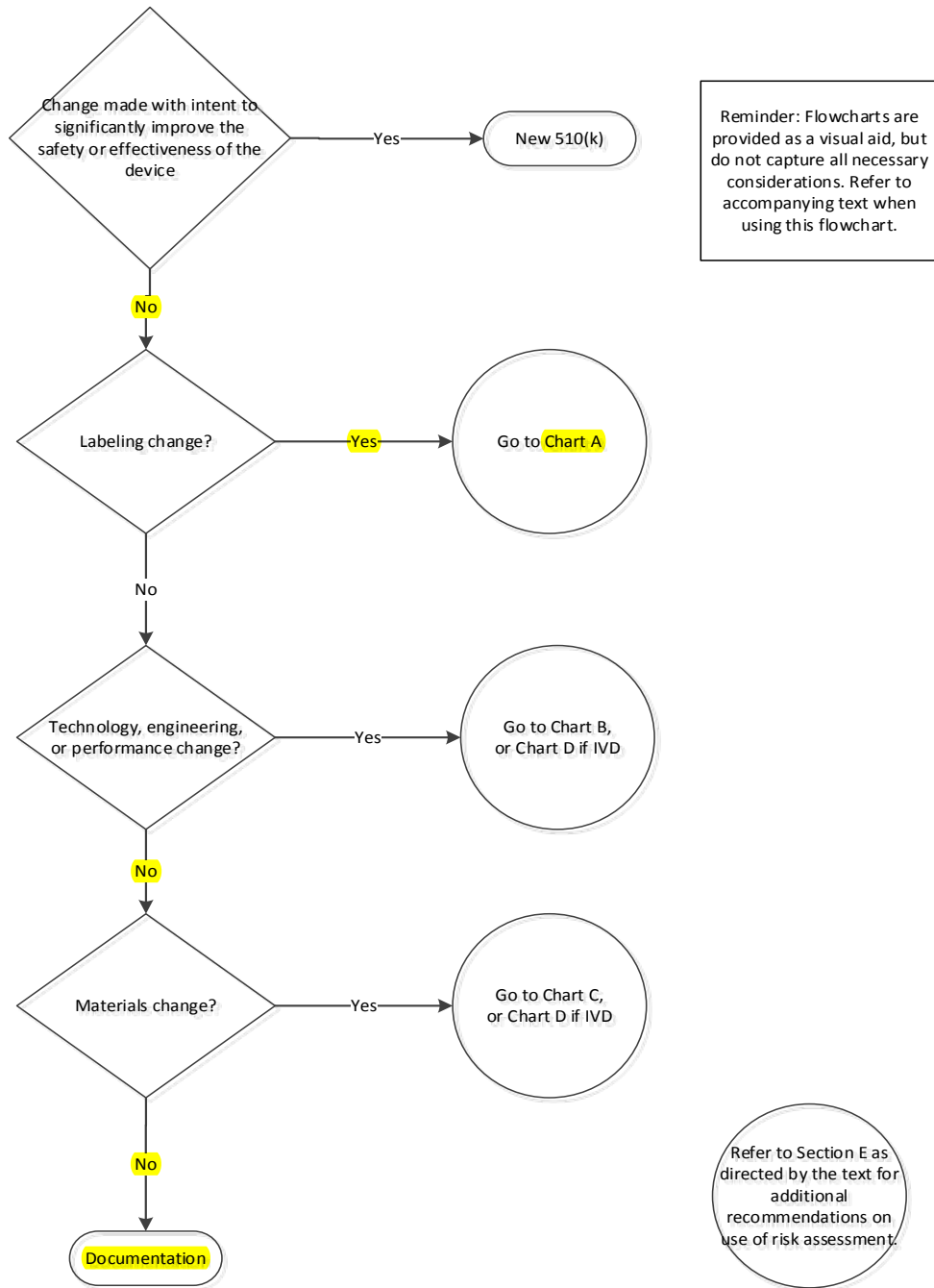


Figure 1 - Main Flowchart

Contains Nonbinding Recommendations

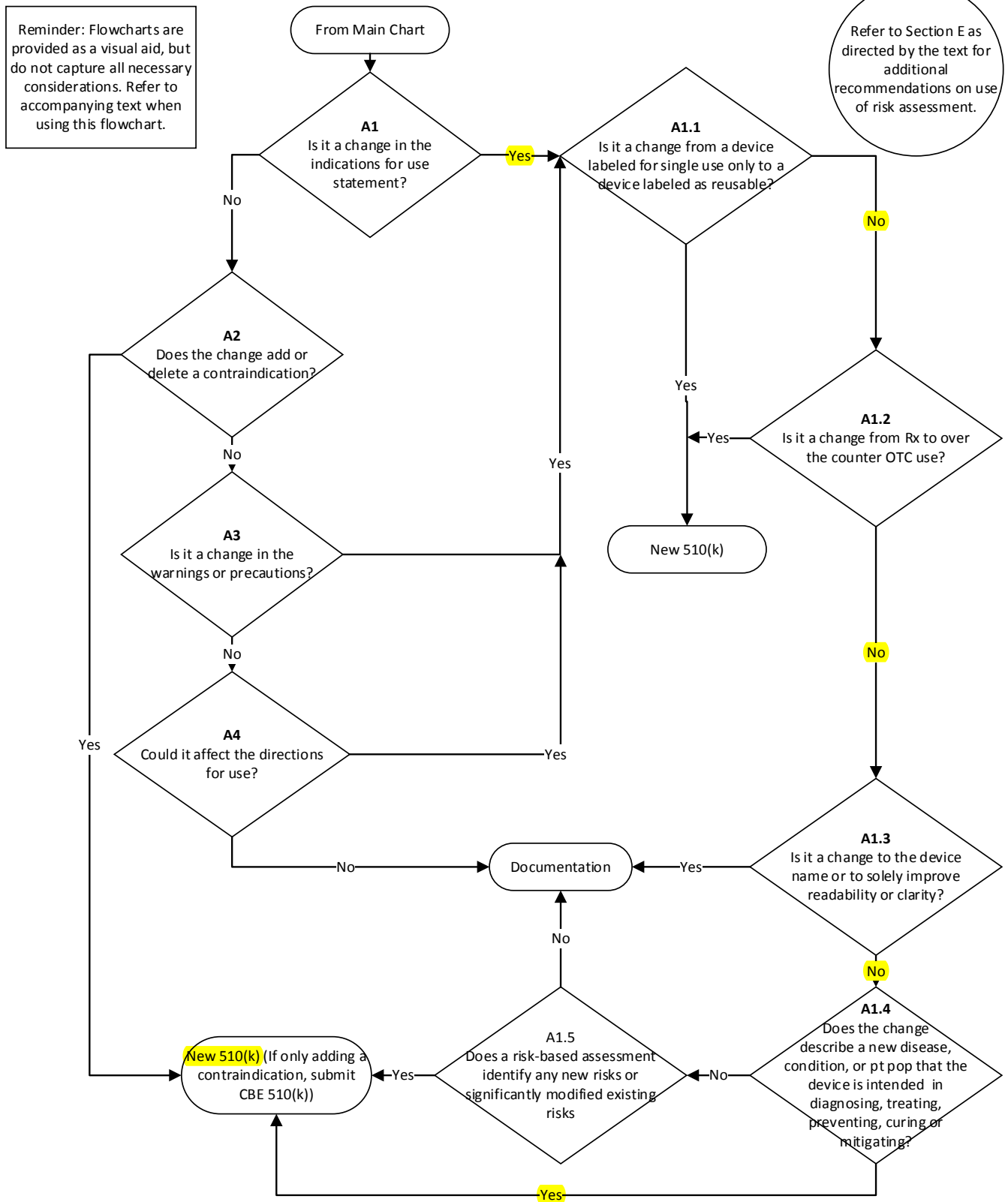


Figure 2 - Flowchart A: Labeling Changes

Contains Nonbinding Recommendations

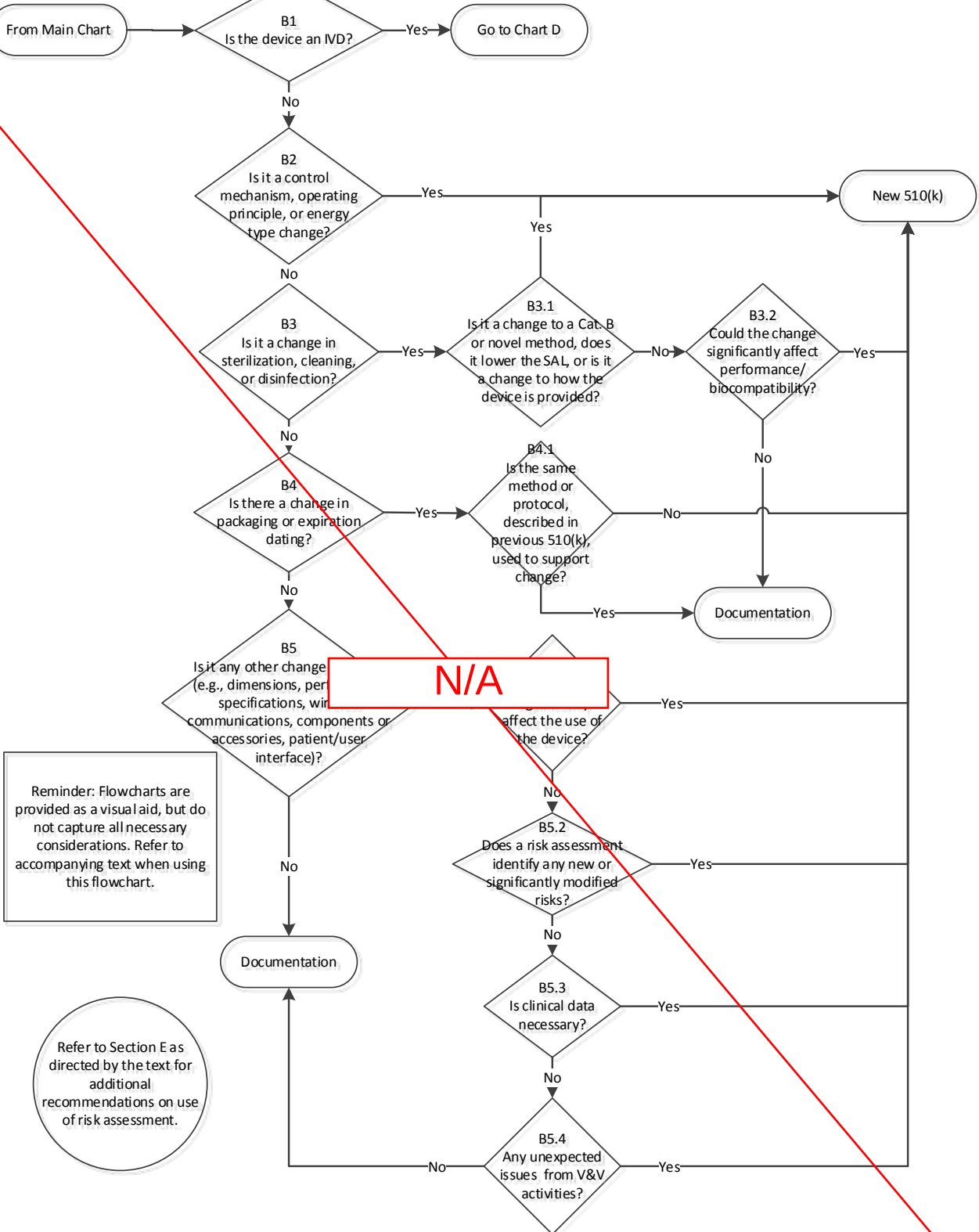


Figure 3 - Flowchart B: Technology, Engineering, and Performance Changes

Contains Nonbinding Recommendations

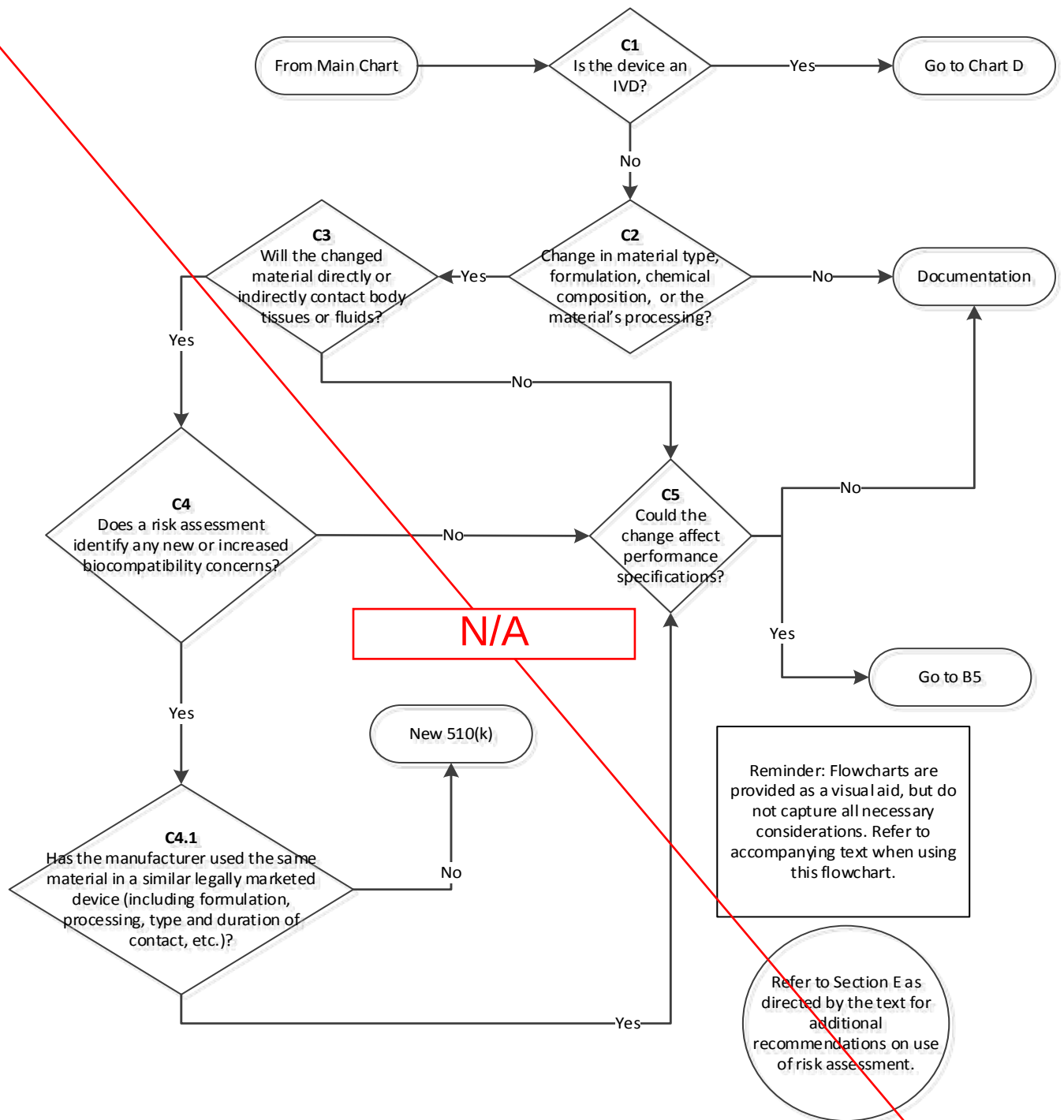


Figure 4 - Flowchart C: Materials Changes