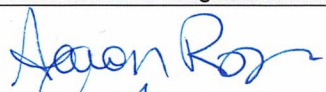
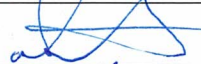
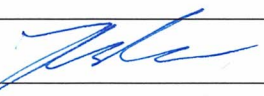


Pathway Innovation LLC.	Design Review		Doc No.: FM-0052
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: B
			Page 1 of 1

Project No.	Product Name	Part No./ Rev	Review Phase	Review Date
00017	Hydrophilix Foley Catheter	N/A	2 (Design Inputs)	7/25/2019

Subject of Review			
#	Topic	Details	Reference Documents
1	Product Specification	Review of document	PS-0006 rev A
2	Hazard Analysis	Review of document. All member participated in developing the risk document.	RM-0019 rev A
3	Project Plan	No changes required to the plan.	PL-0013 rev A
4	Regulatory Plan	No changes required to the plan.	PL-0014 rev A
5	Risk	No new risks discovered, or change required to the RMP or HA.	RM-0018 rev A and #2
6	Review of Phase 2 Deliverables	Each document was reviewed and discussed. Project team in attendance has concluded that the deliverables are complete and satisfy the Design Control requirements in accordance with QP-007 rev D. It is agreed upon to close PHASE 2 and progress into PHASE 3.	Refer to #1 - #5.

Action Items			
#	Description	Assigned To	Due Date
N/A	N/A	N/A	N/A

Reviewers					
Functional Group	Required?	Name	Title	Signature	Date
Quality Assurance / Regulatory Affairs	☒	Aaron Rogers	QA Director		07/25/2019
R&D/ Engineering	☒	Dave Stroup	Managing Partner		07/25/2019
Manufacturing	☒	Jonas Cochran	Managing Partner		7/25/2019
Marketing	☐				
Independent Reviewer	☒	Tom Cassou	Operations / Business Development		7/25/2019
Other	☐				
Other/	☐				

Product Requirement Specification/ Traceability Matrix
Hydrophilix Foley Catheter
General Information Section

Purpose	To document the design requirements for the HDX's Foley Catheter and provide traceability between the risk management reports (i.e, Hazard Analysis, Design Failure Mode Effects Analysis) and Verification & Validation deliverables.
Scope	Foley Catheter (HDX)
Device's Intended Use Statement	Single-use urinary catheter that intended for use in the drainage / irrigation of urinary tract.
FDA Classification	II
FDA Product Code	EZL
FDA Regulation	21 CFR 876.5130 - Urological catheter and accessories
Patient Population	Pediatrics & Adults; Males & Females

Product Requirement Specification/ Traceability Matrix
Hydrophilix Foley Catheter
General Information Section

Source No.	Organization	Document No	Edition	Title
1	ANSI/ AAMI/ ISO	10993-1	2009	Biological Evaluation of Medical Devices Part 1: Evaluation and Testing within a risk management process
2	ANSI/ AAMI/ ISO	11135	2014	Sterilization of health care products – Ethylene Oxide – Requirements for development, validation and routine control of a sterilization process for medical devices
3	ISO	11607-1	2019	Packaging for terminally sterilized medical devices – Part 1: Requirements for materials, sterile barrier systems and packaging systems
4	ANSI/ AAMI/ ISO	15223-1	2016	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied
5	ASTM	D4169	2016	Standard Practice For Performance Testing Of Shipping Containers And Systems
6	ASTM	D4332	2014	Standard Practice for Conditioning Containers, Packages, or Packaging Components for Testing.
7	ASTM	F623	2019	Standard Performance Specification for Foley Catheter
8	FDA	21 CFR 801	4/15/2018	Title 21 - Food and Drugs Chapter I – Food and Drug Administration Department of Health Services; Subchapter H – Medical Devices; Part 801 Labeling
9	FDA	FDA-2013-D-0350	6/16/2016	FDA Guidance: Use of International Standard ISO 10993, “Biological Evaluation of Medical Devices Part 1: Evaluation and Testing within a risk management process”
10	FDA	N/A	1/21/2000	Alternative to Certain Prescription Device Labeling Requirements
11	ISO	14971	2012	Medical Devices - Application of Risk Management to Medical Devices
12	ISO	80369-1	7/10/1905	Small-Bore Connectors For Liquids And Gases In Healthcare Applications - Part 1: General Requirements.
13	FDA	N/A	N/A	Guidance for the Content of Premarket Notifications for Conventional and Antimicrobial Foley Catheters
14	FDA	81 FR 38911	6/15/2016	Use of Symbols in Labeling

Hydrophilix Foley Catheter

General Information Section

Type	Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference	Other Information
Biological	The device shall meet biocompatibility requirement for a surface device with mucosal membrane contact for a prolonged exposure >24 to 30 days, which the following evaluations are required: 1. Cytotoxicity 2. Sensitization 3. Irritation or Intracutaneous Reactivity 4. Material Mediated Pyrogen	Device	1	1, 7, 9		TR-0071	N/A	
Labeling	Information shall be presented in symbol if the symbol is recognized by the regulatory agency(s)	Labeling	2	4, 10, 14		N/A	Refer to product labeling (HDX-4XXX) indicated in BOM for the following finished products: HDX-8000 HDX-8001 HDX-8002	
Labeling	Product labeling (including direction for use) shall comply to FDA regulations on labeling. Specific requirement for devices with antimicrobial agents is indicated in guidance document (14).	Labeling	3	8, 10, 13		N/A	Refer to product labeling (HDX-4XXX) indicated in BOM for the following finished products: HDX-8000 HDX-8001 HDX-8002	
Labeling	Product labeling shall include the following information associated with the device description: 1. Outer Diameter (French Size) 2. Balloon Size (cc) 3. Effective Shaft Size 4. Type of Coating (if applicable) 5. Amount/ concentration of antimicrobial agent (if applicable)	Labeling	4	7, 13		N/A	Refer to product labeling (HDX-4XXX) indicated in BOM for the following finished products: HDX-8000 HDX-8001 HDX-8002	
Material	Connector shall be made of rigid or semi-rigid materials.	Connector	5	12		N/A	HDX-6000 HDX-6001 HDX-6002	
Material	The device shall be constructed on medical grade silicone material.	Device	6	N/A		N/A	HDX-6000 HDX-6001 HDX-6002	
Packaging	The device will be packaged in a tyvek/poly pouch, which 10 pouched devices will be packaged in an box.	Packaged Device	7	N/A		N/A	HDX-8000 HDX-8001 HDX-8002	
Packaging	The peel strength for 1" MFG sealed section of the sterile barrier packaging shall be ≥ 1 lb.	Sterile Barrier	8	3		TR-0074	N/A	
Packaging	The sterile barrier packaging (include sealed area) shall not exhibit the following visual defects: 1. Unsealed Areas 2. Nonhomogeneous or Undersealed Areas 3. Oversealed Areas 4. Narrow Seals 5. Channels 6. Wrinkles/Foldovers/Cracks 7. Tears/Pinholes	Sterile Barrier	9	3		TR-0074	N/A	

Hydrophilix Foley Catheter

General Information Section

Type	Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference	Other Information
Packaging	The sterile barrier packaging shall not exhibit a constant stream of bubbles when submerge in water and internal pressurized to a determine test pressure.	Sterile Barrier	10	3		TR-0074	N/A	
Performance	Balloon Integrity (Resistance to Rupture) - The balloon shall be inflated easily with distilled or deionized water to labeled volume without showing any evidence of breakage throughout a 24 hour period.	Device	11	7		TR-0072	N/A	
Performance	Balloon Volume Maintenance - The balloon shall maintain its volume throughout a 24 hour period.	Device	12	7		TR-0072	N/A	
Performance	Connector shall be non-interconnectables, which does not appear to provide a secure connection when forcefully assembled and shall easily disengage.	Connector	13	12		TR-0072	N/A	
Performance	Deflation Reliability (Failure to Deflate) - The balloon shall deflate to within 4 French sizes of the label French size within 15 min or be otherwise manipulated to effect drainage within this time period.	Device	14	7		TR-0072	N/A	
Performance	Flow Rate through Drainage Lumen - Shall have a minimum average flow rate of 100 cm ³ /min.	Device	15	7		TR-0072	N/A	
Performance	Inflated Balloon Response to Traction - The entire, inflated balloon of catheters shall not pass into or through the funnel barrel neck (28 Fr x 6 cm) when a 1 lb weight is affixed from the distal catheter end and subjected to static (2 minutes) and dynamic load (2 ft drop).	Device	16	7		TR-0072	N/A	
Performance	The device shall be easily inserted through the urethra until the deflated balloon passes through the urethral sphincter and in a manner that prevents urethral trauma.	Device	17	N/A		TR-0072	N/A	
Physical	Balloon & Shaft Size - The proximal tip, the balloon, and the shaft (for 10 in. distal to the balloon) shall meet the requirements on size and tolerances as it relates to the labeled French size per the tolerances indicated below: Tip = ± 1 Fr Shaft = ± 1 Fr Balloon = + 4 Fr (Uninflated and Deflated after Immersion)	Device	18	7		TR-0072	HDX-6000 HDX-6001 HDX-6002	
Physical	The balloon size shall be 10cc.	Device	19	N/A		TR-0072	HDX-6000 HDX-6001 HDX-6002	
Physical	The catheter shall have the ability to flush the outer surface proximal to the inflated balloon.	Device	20	N/A		TR-0072	HDX-6000 HDX-6001 HDX-6002	
Physical	the device shall be able to connect into an industry standard or equivalent urine bag to facilitate sterile and easy urine collection.	Device	21	N/A		TR-0072	HDX-6000 HDX-6001 HDX-6002	
Physical	The device shall be equipment with 3 lumens, which have the following purposes: 1. Drainage 2. Balloon Inflation 3. Irrigation	Device	22	N/A		N/A	HDX-6000 HDX-6001 HDX-6002	

Product Requirement Specification/ Traceability Matrix
Hydrophilix Foley Catheter
General Information Section

Type	Description	Applicable Feature of Device	REQ No.	Regulations/ Standards Reference	Risk Management Reference	Design Verification & Validation Evidence	Specification Reference	Other Information
Physical	The product offering will be French Sizes (Outer Diameter) of 14 Fr (0.184"), 16 Fr (0.210"), and 18 Fr (0.236").	Device	23	N/A		N/A	HDX-6000 HDX-6001 HDX-6002	
Stability	The device shall have a final shelf-life period of 3 years, but initially be qualified for shorter shelf-life periods (e.g., 1 year).	Sterile, Packaged Device	24	N/A		TR-0073	N/A	
Sterilization	Device shall be sterile to a sterility assurance level (SAL) of 10 ⁻⁶ when terminally sterilized via Ethylene Oxide (EO).	Device	25	2		TR-0075	N/A	
Performance	Packaged devices shall maintain it sterility and functionality through transportation	Sterile, Packaged Device	26	5, 6		TR-0074	N/A	
User	A healthcare provider must be able to easily and quickly implement device.	Device	27	N/A		N/A	HDX-4001	

Hydrophilix Foley Catheter

Rev A

General Information Section

Rev	Description of Changes
A	Initial release

Pathway Innovation LLC.	Hazard Analysis (HA) – Risk Management Report, HDX Foley Catheter	Doc No.: RM-0019
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 1 of 9

1.0 Purpose

- 1.1 The purpose of this document is to identify the associated risk with the intended use of the HDX Foley Catheter.

2.0 Scope

- 2.1 The product numbers associated with the HDX Foley Catheter is indicated in the following table:

Table 2.1 – HDX Foley Catheter Products	
Part No.	Description
HDX-8000	14 Fr, 2-Way Foley Catheter
HDX-8001	16 Fr, 2-Way Foley Catheter
HDX-8002	18 Fr, 2-Way Foley Catheter

- 2.2 The intended use of this device is the following:

Single-use urinary catheter that intended for use in the drainage / irrigation of the urinary tract.

3.0 References

- 3.1 QP-004 – Risk Management

4.0 Risk Tables Definitions

- 4.1 Description – Describes the operation that is being evaluated for risk
- 4.2 Verification – Evidence source of linkage that initial (i.e., Risk Controls¹) risk controls have implemented, verified, or validated.
- 4.3 Probability – Likelihood of the hazard occurring
- 4.3.1 Probability¹ (P¹) – The probability assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls¹)
- 4.3.2 Probability² (P²) – The probability assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls²)
- 4.4 Potential Hazard (PH) – Source of harm
- 4.5 Potential Failure Mode (FM) – The potential consequences of the failure
- 4.6 Risk Controls - Process through which decisions are reached and protective measures are implemented for reducing risk to, or maintaining risk within specified levels.
- 4.6.1 Risk Controls¹ – Initial risk controls that are incorporated.
- 4.6.2 Risk Controls² – Secondary risk controls, also known as “Mitigation”.

Pathway Innovation LLC.	Hazard Analysis (HA) – Risk Management Report, HDX Foley Catheter	Doc No.: RM-0019
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 2 of 9

4.7 Risk Rating – Residual risk determination based on the concept of ‘As Far As Possible (AFAP)’ approach. The AFAP approach provides three (3) possible outcomes, which are 1) the residual risk exceeds the risk acceptability (RED), 2) the residual risk is acceptable to the risk acceptability criteria (GREEN), and 3) the residual risk is between the two (2) states indicated and additional controls are suggested to lower the risk to the lowest practical level, without economically considerations, based on management decision (YELLOW).

4.7.1 Risk Rating¹ (R¹) – This calculation is post the implementation of initial risk controls (i.e., Risk Controls¹)

4.7.2 Risk Rating² (R²) – This calculation is post the implementation of the secondary risk controls (i.e., Risk Controls²)

4.8 Severity (S) – Measure of the possible consequences of a hazard.

5.0 Identify Characteristics that could Impact on Safety (Annex C)

C.2.1 What is the intended use and how is the medical device to be used? *Refer to §2.2.*

C.2.2 Is the medical device intended to be implanted? *No*

C.2.3 Is the medical device intended to be in contact with the patient or other persons? *Yes*

C.2.4 What materials or components are utilized in the medical device or are used with, or are in contact with, the medical device? *Silicone*

C.2.5 Is energy delivered to or extracted from the patient? *No*

C.2.6 Are substances delivered to or extracted from the patient? *Yes, substances can be delivery to the urinary tract.*

C.2.7 Are biological materials processed by the medical device for subsequent re-use, transfusion or transplantation? *No*

C.2.8 Is the medical device supplied sterile or intended to be sterilized by the user, or are other microbiological controls applicable? *Yes*

C.2.9 Is the medical device intended to be routinely cleaned and disinfected by the user? *No*

C.2.10 Is the medical device intended to modify the patient environment? *No*

C.2.11 Are measurements taken? *Yes, measurements of urine.*

C.2.12 Is the medical device interpretative? *No*

C.2.13 Is the medical device intended for use in conjunction with other medical devices, medicines or other medical technologies? *Yes, urine collection bag.*

C.2.14 Are there unwanted outputs of energy or substances? *No*

C.2.15 Is the medical device susceptible to environmental influences? *No*

Pathway Innovation LLC.	Hazard Analysis (HA) – Risk Management Report, HDX Foley Catheter	Doc No.: RM-0019
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 3 of 9

- C.2.16 Does the medical device influence the environment? *No*
- C.2.17 Are there essential consumables or accessories associated with the medical device? *No*
- C.2.18 Is maintenance or calibration necessary? *No*
- C.2.19 Does the medical device contain software? *No*
- C.2.20 Does the medical device have a restricted shelf-life? *Yes*
- C.2.21 Are there any delayed or long-term use effects? *No*
- C.2.22 To what mechanical forces will the medical device be subjected? *Tensile force from the urine content in the collection bag; Internal pressure from urine drainage or irrigation of the urinary tract.*
- C.2.23 What determines the lifetime of the medical device? *Product and sterile barrier performance*
- C.2.24 Is the medical device intended for single use? *Yes*
- C.2.25 Is safe decommissioning or disposal of the medical device necessary? *Yes*
- C.2.26 Does installation or use of the medical device require special training or special skills? *No*
- C.2.27 How will information for safe use be provided? *Direction for Use*
- C.2.28 Will new manufacturing processes need to be established or introduced? *No*
- C.2.29 Is successful application of the medical device critically dependent on human factors such as the user interface? *No*
 - C.2.29.1 Can the user interface design features contribute to use error? *No*
 - C.2.29.2 Is the medical device used in an environment where distractions can cause use error? *No*
 - C.2.29.3 Does the medical device have connecting parts or accessories? *Yes, urine collection bag*
 - C.2.29.4 Does the medical device have a control interface? *No*
 - C.2.29.5 Does the medical device display information? *No*
 - C.2.29.6 Is the medical device controlled by a menu? *No*
 - C.2.29.7 Will the medical device be used by persons with special needs? *No*
 - C.2.29.8 Can the user interface be used to initiate user actions? *No*
- C.2.30 Does the medical device use an alarm system? *No*
- C.2.31 In what way(s) might the medical device be deliberately misused? *None*

Pathway Innovation LLC.	Hazard Analysis (HA) – Risk Management Report, HDX Foley Catheter	Doc No.: RM-0019
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 4 of 9

- C.2.32 Does the medical device hold data critical to patient care? *No*
- C.2.33 Is the medical device intended to be mobile or portable? *Yes*
- C.2.34 Does the use of the medical device depend on essential performance? *No*

6.0 Risk Acceptability Determination

- 6.1 The criteria of risk acceptability will be determined with a Risk Priority Number (RPN) and on the concept of 'As Far As Possible' (AFAP) approach. Refer to the Risk Rating definition in §4.0.
- 6.2 There are two (2) components used in determining risk acceptability, which are Severity and Probability. Refer to the following tables:

Table 5.2a - Severity Ratings	
Index	Description
1	No patient effect
2	Minor injury that resolves without treatment
3	Moderate injury that may require treatment
4	Serious injury requiring medical intervention
5	Life threatening injury, disability, death

Table 5.2b – Probability Ratings	
Index	Description
1	Improbable - <0.001%
2	Remote – 0.001% to .01%
3	Occasional – 0.01% to 0.1%
4	Likely – 0.1% to 1%
5	Frequent - >1%

- 6.3 Risk acceptability will be determined per the following risk index:

Table 5.3 – Risk Rating Table						
Severity	5	5	10	15	20	25
	4	4	8	12	16	20
	3	3	6	9	12	15
	2	2	4	6	8	10
	1	1	2	3	4	5
	1	1	2	3	4	5
Probability						

7.0 Risk Assessment Tables

Table 6.0 –Kurin Blood Collection Device

Pathway Innovation LLC.	Hazard Analysis (HA) – Risk Management Report, HDX Foley Catheter										Doc No.: RM-0019
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>										Rev.: A
											Page 5 of 9

#	Hazard Category	PH #	Potential Hazard	S	Potential Failure Mode	Risk Controls ¹	Verification	P ¹	R ¹	Risk Controls ²	P ²	R ²
1	Electromagnetic Energy	The medical device does not have any electromagnetic energy hazards.										
2	Radiation Energy	The medical device does not have any radiation energy hazards.										
3	Thermal Energy	The medical device does not have any thermal energy hazards.										
4	Mechanical Energy	The medical device does not have any mechanical energy hazards.										
5	Biological	5.1	Fever/ infection to patient	3	Device is not provided sterile/ insufficient sterilization process	Device will be terminally sterilized with validated method (Design)	Sterilization validation	1	3	N/A	N/A	N/A
					Device is reused	Warnings statement for single-use (Information)	IFU	2	6	N/A	N/A	N/A
					Inadequate sterile barrier system	Sterile barrier system will be designed for material that are commonly used for single-use medical device (Design)	Design verification - sterile barrier	1	3	N/A	N/A	N/A
						Warning will be incorporated for inspecting packaging prior to use as well not using the device if packaging is damaged (Information)	IFU	2	6	N/A	N/A	N/A
					Use beyond expiration date	Warning statement (Information)	IFU	2	6	N/A	N/A	N/A
		5.2	Fever/infection to clinician	3	Proper handling and disposal is not conducted during clinical application	Instruction and warning statement (Information)	IFU	2	6	N/A	N/A	N/A
6	Chemical	6.1	Allergic reaction	3	Materials (e.g., latex, rubber) that is known to cause allergic response	Material will be selected that do not cause allergic response (Design)	Design verification - biocompatibility	1	3	N/A	N/A	N/A
		6.2	False positives in the blood work evaluation – red blood cell damage	1	Material leach during clinical application	Material selection (Design)	Design verification - biocompatibility	1	1	N/A	N/A	N/A
					Excessive ethylene oxide residuals from sterilization process	Device will be released for distribution with acceptable residual levels (Design)	Sterilization validation	1	1	N/A	N/A	N/A

Pathway Innovation LLC.	Hazard Analysis (HA) – Risk Management Report, HDX Foley Catheter									Doc No.: RM-0019		
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>									Rev.: A		
										Page 6 of 9		

Table 6.0 –Kurin Blood Collection Device												
#	Hazard Category	PH #	Potential Hazard	S	Potential Failure Mode	Risk Controls ¹	Verification	P ¹	R ¹	Risk Controls ²	P ²	R ²
		6.3	Fever	3	Contacting material causes pyrogenic response	Material selection (Design)	Design verification - biocompatibility (i.e., Material mediated pyrogenicity)	1	3	N/A	N/A	N/A
		6.4	Procedure complication – device inoperative	3	Material degrades over usable life	Material selection (Design) Sterilization method selection (Design)	Stability evaluation	1	3	N/A	N/A	N/A
7	Biocompatibility	7.1	Cell Damage	3	Patient contact materials are not biocompatibility	Material selection (Design)	Design verification - biocompatibility	1	3	N/A	N/A	N/A
8	Function (Operational)	8.1	Procedure complication – device inoperative	2	Balloon Burst – Material or dimensional issues (wall thickness)	Material selection (Design) Dimensional Assessment (Design)	Design verification - Balloon Integrity (Resistance to Rupture)	1	2	N/A	N/A	N/A
				1	Balloon deflates or does not inflates – Valve malfunctions or leaks (wall thickness issue)	Valve selection (Design) Dimensional Assessment (Design)	Design verification - Balloon Volume Maintenance	1	1	N/A	N/A	N/A
				1	Balloon does not inflate – Valve malfunctions	Valve selection (Design)	Design verification	1	1	N/A	N/A	N/A
				2	Catheter breaks under normal tensile stresses - Material or dimensional issues	Valve selection (Design) Dimensional Assessment (Design)	Design verification - Inflated Balloon Response to Traction	1	2	N/A	N/A	N/A
				1	Connection issues – Syringe (inflation, irrigation); Collection bag (drainage)	Connection will be standardized to what is incorporated with marketed catheters (Design)	Design verification	1	1	N/A	N/A	N/A
		8.2	Procedure complication – downgraded performance	2	Poor drainage - Insufficient holes and lumen diameter, Lumen collapses (Materials or dimensional issues)	Material selection (Design) Dimensional Assessment (Design)	Design verification - Flow Rate though Drainage Lumen	1	2	N/A	N/A	N/A
				2	Poor irrigation - Insufficient holes and lumen diameter, Lumen collapses (Materials or dimensional issues)	Material selection (Design) Dimensional Assessment (Design)	Design verification	1	2	N/A	N/A	N/A

Pathway Innovation LLC.	Hazard Analysis (HA) – Risk Management Report, HDX Foley Catheter									Doc No.: RM-0019
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>									Rev.: A
										Page 8 of 9

Table 6.0 –Kurin Blood Collection Device												
#	Hazard Category	PH #	Potential Hazard	S	Potential Failure Mode	Risk Controls ¹	Verification	P ¹	R ¹	Risk Controls ²	P ²	R ²
10	Labeling (Product Labeling)	10.1	Loss of traceability - Cannot identify the device or trace back to MFG	1	Device used after expiration date or outside biocompatibility duration limits							
					Occurs during disposal process							
					Incomplete or ambiguous instructions for use	Label and labeling will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review	1	1	N/A	N/A	N/A
11	Operating Instructions (Product Labeling)	11.1	Procedural Complication /Patient discomfort	3	Incomplete or ambiguous instructions for use	IFU will be designed to requirements indicated in regulation for the markets intended to be distributed (Information)	Label content review	1	3	N/A	N/A	N/A
12	Environmental	12.1	Device/ Sterile Barrier Failure /Patient health is compromised	3	Stored outside the prescribe environmental parameters	Warning statement (Information)	IFU	2	6	N/A	N/A	N/A
					Packaging is not adequate to protect the device through transportation	Packaging will be design to protect the device and it sterile barrier through transportation stresses (Design)	Design Verification (Distribution Evaluation)	1	3	N/A	N/A	N/A

8.0 Risk/ Benefit Assessment

8.1 The following table indicates the risk/benefit analysis for the residual risks that were greater than the risk acceptability limit after all applicable controls/mitigations were applied:

Table 7.1 – Risk/Benefit Assessment for Residual Risk above Acceptability Limit for Process Hazard Assessment										
Table No.	Process Description	Section No.	Potential Failure Mode	Potential Effect of Failure	Potential Cause of Failure	S	F ²	Risk/Benefit Analysis		Does the benefit outweigh the residual risk?
N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A		N/A

Pathway Innovation LLC.	Hazard Analysis (HA) – Risk Management Report, HDX Foley Catheter	Doc No.: RM-0019
	<i>This document contains confidential, proprietary information of Pathway Innovation LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: A
		Page 9 of 9

8.2 There was no residual risk identified that was above the risk acceptability limit per §5.3.

9.0 Conclusion

9.1 Risks have been reduced as low as possible by the following methods:

- Design
- Appropriate testing to ensure design criteria were met
- Protective measures (critical-to-quality)
- Manufacturing and quality control measures

9.2 Based on a review of the risk-management process for the HDX Foley Catheter described herein, the risk management plan based on ISO 14971:2012 has been implemented and the individual, cumulative and overall residual risk have been reduced as low as possible regardless of the magnitude and the benefit of using this device will outweigh risks posed by the device during usage.

10.0 Change History

Revision	Change Description
A	Initial release