



Rx Only

2-Way 100% Silicone ClearTract SPT Catheter

Please read all instructions before using this device

Not made with natural rubber latex

Foley Catheter is MR Safe.

Model	Size
211410SPT	14 FR 10 mL
211605SPT	16 FR 5 mL
211610SPT	16 FR 10 mL
211810SPT	18 FR 10 mL

Contraindications

This device is contraindicated in the presence of conditions which create unacceptable risk during catheterization.

Device Description

The 2-Way 100% Silicone ClearTract Catheter is a standard single-use 2-way Foley catheter that is constructed of medical grade silicone with a surface modification. It incorporates two (2) lumens, one for inflation/deflation of the balloon and the other for drainage of the urinary tract.

The surface of the ClearTract Catheter has a surface modification that reduces surface friction.

Intended Use

The 2-Way 100% Silicone ClearTract Catheter is intended for drainage of the urinary tract. Catheterization is generally accomplished by inserting the catheter through the urethra and into the bladder. However, drainage is sometimes accomplished by suprapubic placement or other placement of the catheter, such as nephrostomy tract. Intended population is adults and pediatrics.

Directions

Installation & Inflation

1. When placement of a urethral catheter is contraindicated or unsuccessful, percutaneous suprapubic urinary bladder catheterization or nephrostomy can be used in accordance with the procedures and policies of the hospital or healthcare facility.
2. Prepare application site in a sterile fashion.
3. Visually inspect the catheter for any signs of damage or deterioration prior to use.
4. If a trial inflation is conducted, follow steps 7 & 8 on the sterile field then evacuate the sterile water with the syringe. Ensure the balloon is empty prior to insertion.
5. Apply medical grade lubricant to distal end of the catheter.
6. Gently advance the catheter through the urethra until distal tip is beyond the bladder neck.
7. Confirm the position of the distal tip. The return of urine typically indicates the distal tip is in the bladder.
8. Fill Luer-tipped syringe with sterile water and connect to the inflation valve.
9. Fill the catheter balloon to the labeled volume. Do NOT exceed labeled volume.
10. Slowly withdraw the catheter until resistance is felt. This is a sign the balloon is resting at the bladder neck.



Deflation & Removal

1. Gently connect Luer-tipped syringe to the inflation valve. Gently aspirate while allowing the pressure in the balloon to assist filling the syringe with water occupying the balloon. Remove and reconnect the syringe in the inflation valve if there is no or slow deflation. Rapid deflation could collapse the lumen and prevent adequate deflation of the balloon.
2. If during the removal process the balloon ruptures, ensure all fragments are carefully removed from patient.
3. Discard after use.

Precautions

1. Sterile technique shall be applied during installation and removal of the catheter.
2. It is recommended to be replaced after 7 days of use, however, shall not exceed 30 days.

Cautions

1. Federal law restricts the device to sale by or on the order of a physician.
2. Never clamp the catheter. Use catheter plug or valve if required.

Warnings

1. Single Use ONLY. Do NOT reuse. Reuse may cause infection, illness, injury or death.
2. Read the instructions carefully before using the device.
3. Do not use this device beyond the labeled expiration date.
4. Do not use the device if the packaging is punctured or damaged.
5. Do not force catheter as injury to the urethra can occur.
6. Do not use needle.
7. Do not exceed recommended capacities.

Potential Adverse Events

With all Foley Catheters, the following potential adverse events could occur:

1. Urinary tract infection
2. Perforation or laceration of the urethra
3. Perforation or laceration of the bladder
4. Abrasions

Symbol Description

	Catalog Number		Do not resterilize
	Lot Number		Do not use if package is damaged
	Use by		Prescription Use Only (US)
	Consult instructions for use		Foley Catheter is MR Safe
	Single Use		Manufacturer Contact Information
	Sterilized by Ethylene Oxide		Device is non-pyrogenic



SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor
Sunny Isles Beach, FL 33160
(424) 309-8523



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	1

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex



SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

0001
0002

0006



Rx Only



0085



0010 0011 0012



SLQ-MMDDYYYY



XXXXXX



20YY-MM-DD



(01) #####
(10) SLQ-MMDDYYYY
(17) YYMMDD

SLQ-LB-4008 Rev. A

2-WAY 100% SILICONE
CLEARTRACT SPT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	10



Rx Only



(01) #####
(10) SLQ-MMDDYYYY
(17) YYMMDD

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

0006

0005

0085



STERILE EO




MR

LOT SLQ-MMDDYYYY

REF XXXXXX

20YY-MM-DD

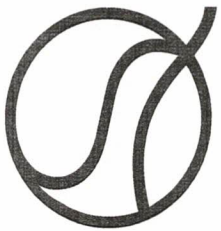

2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF XXXXXX

LOT SLQ-MMDDYYYY

0000  20YY-MM-DD

QTY 100



SILQ Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

0001

0002

	Label Content Design Requirements		Doc No.: WI-0009
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: C
			Page 2 of 7

- 5.1 It is the responsibility of the individual that designs product labels and labeling to comply with this procedure.

6.0 General Requirements

- 6.1 Paper stock and printing materials shall be selected based on its ability to maintain legibility throughout the labeling life (e.g., processing, storage, handling, and distribution).
- 6.2 Adhesive, if applicable, shall be selected based on its ability to maintain adherence throughout the labeling life (e.g., processing, storage, handling, and distribution).
- 6.3 Additional labeling requirements maybe required based on product specific regulations or guidance documents.

7.0 Labeling Requirements

#	Reference	Requirement Description	Device / Label / IFU	Markets
✓ 0001	21CFR801.1 21CFR812.5	Device label shall include the name of manufacturer, packer, or distributor. <i>Abbreviations for "Incorporated", etc. may be used. Qualifying phrase to distinguish relationship "Manufactured for ____", etc. may be used.</i>	Label	US
✓ 0002	21CFR801.1 21CFR812.5	Device label shall include street address, city, state, and zip code of manufacturer, packer, or distributor. <i>Street address may be omitted if listed in current directory.</i>	Label	US
✓ 0003	21CFR801.15	The predominant language shall be English.	Label / IFU	US
✓ 0004	21CFR801.15	Use of symbols shall include either adjacent text in all language(s) or a glossary of symbols.	Label / IFU	US
✓ 0005	21CFR801.15 21CFR801.109	The symbol statement "Rx only" or "[rx] only" or the statement "Caution: Federal law restricts this device to sale by or on the order of a ____" (i.e. "physician", etc.) <i>Except for surgical instruments.</i> "Only" needs to immediately follow "Rx." However, the symbol statement "Rx only" does not necessarily need to be bracketed in quotation marks, and the word "only" may appear in upper or lower case letters, for example, Rx only, Rx Only, or Rx ONLY.	Label / IFU	US
✓ 0006	21CFR801.109	The label of the device includes the method of its application or use. <i>Except for surgical instruments.</i>	Label / IFU	US
✓ 0007	21CFR801.15	If a glossary of symbols is used, it shall contain: - Each symbol used in the labeling - Title and designation number of standard(s) - Title of symbol and reference number (if any) - The symbols meaning	IFU	US
N/A 0008	21CFR801.15	If any representation appears in a foreign language, then all labeling requirements shall also appear in that foreign language.	Label / IFU	US
✓ 0009	21CFR801.18	Dates shall be expressed in the format YYYY-MM-DD (i.e. 4-digit year, 2-digit month, 2-digit day, each separated by a hyphen).	Label / IFU	US
✓ 0010	21CFR801.20 MDR Annex I 23.2	Label shall include a unique device identifier (UDI).	Label	US

	Label Content Design Requirements <i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Doc No.: WI-0009 Rev.: C Page 3 of 7
--	---	--

#	Reference	Requirement Description	Device / Label / IFU	Markets
✓ 0011	21CFR801.40	UDI must be presented in two forms: - Plain-text - Automatic Identification and Data Capture (AIDC)	Label	US
✓ 0012	21CFR801.40	UDI includes the following information: - Device identifier, which is a fixed portion that identifies the specific version or model of a device and the labeler of the device - Production identifier, which is a variable portion that identifies one or more of the following: lot number, serial number, expiration date, manufacturing date	Label	US
N/A 0013	21CFR801.45	Device meant to be reprocessed must be directly marked with UDI in either or both: - Plain-text - Automatic Identification and Data Capture (AIDC) or alternative technology <i>Except when it interferes with safety and efficacy, or not technically feasible.</i>	Device	US
✓ 0014	21CFR801.5	Instruction for Use includes statements of all conditions, purposes or uses for which such device is intended, prescribed or recommended.	IFU	US
✓ 0015	21CFR801.5	Instruction for Use includes Quantity of dose for: - each intended use - persons of different ages and physical conditions	IFU	US
✓ 0016	21CFR801.5 21CFR801.109	Instruction for Use includes frequency of administration or application.	IFU	US
✓ 0017	21CFR801.5 21CFR801.109	Instruction for Use includes duration of administration or application.	IFU	US
✓ 0018	21CFR801.5 21CFR801.109	Instruction for Use includes time of administration or application in relation to meals, onset of symptoms or other time factors.	IFU	US
✓ 0019	21CFR801.5 21CFR801.109	Instruction for Use includes route or method of administration or application.	IFU	US
✓ 0020	21CFR801.5	Instruction for Use includes preparation for use.	IFU	US
✓ 0021	21CFR801.109	Instruction for Use includes indications for Use, Effects, Hazards, Contraindications, side effects, precautions.	IFU	US
✓ 0022	21CFR801.109	Instruction for Use includes date of issuance or latest revision.	IFU	US
N/A 0023	21CFR812.5	An Investigational Device Exemption (IDE) label includes the quantity of contents.	Label	US
N/A 0024	21CFR812.5	An IDE label includes "CAUTION--Investigational device. Limited by Federal (or United States) law to investigational use."	Label	US
✓ 0025	21CFR812.5	The label or other labeling shall describe all relevant contraindications, hazards, adverse effects, interfering substances or devices, warnings, and precautions.	Label / IFU	US
N/A 0026	21CFR812.5	The labeling of an investigational device shall not bear any statement that is false or misleading in any particular and shall not represent that the device is safe or effective for the purposes for which it is being investigated.	Label / IFU	US

	Label Content Design Requirements		Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 4 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0027	21CFR812.5	An investigational device shipped solely for research on or with laboratory animals shall bear on its label the following statement: "CAUTION--Device for investigational use in laboratory animals or other tests that do not involve human subjects."	Label	US
0028	MDR Annex I 23.1	Information provided (medium, format, content, legibility, location) is appropriate for the device and healthcare practitioner.	Label / IFU	EU
0029	MDR Annex I 23.1	Where practical, information shall appear on the device itself.	Device	EU
0030	MDR Annex I 23.1	Information on device labeling shall take the form of symbols, where appropriate. In areas for which no standards exist, the symbols and colors shall be described in the documentation supplied with the device.	Label / IFU	EU
0031	MDR Annex I 23.1	Labels shall be provided in human-readable format and supplemented with machine-readable information.	Label	EU
0032	MDR Annex I 23.1	Each device shall be accompanied by Instructions for Use (IFU). <i>Except for class I and IIa devices that don't require information for safe use.</i>	IFU	EU
0033	MDR Annex I 23.1, 23.2, 23.3	Device name shall appear on the device label, sterile barrier packaging and IFU.	Label / IFU	EU
0034	MDR Annex I 23.1, 23.2, 23.3	Details needed to identify the device / device description shall appear on the device label, sterile barrier packaging.	Label	EU
0035	MDR Annex I 23.1, 23.2	Packaging contents shall appear on the device label.	Label	EU
0036	MDR Annex I 23.1, 23.2, 23.4	Intended for use statement shall appear on the device label.	Label	EU
0037	MDR Annex I 23.1, 23.2, 23.3, 23.4	Name and address of the manufacturer shall appear on the device label, sterile barrier packaging, IFU.	Label / IFU	EU
0038	MDR Annex I 23.2	Name and Address of authorized representative shall appear on the device label. <i>Except where the manufacturer resides within the EU.</i>	Label	EU
0039	MDR Annex I 23.2, 23.4	Where applicable, device label and IFU shall include an indication device contains: - Medicinal substance, human blood, plasma derivative - Tissues, cells or derivatives of human origin - Tissues, cells or derivatives of animal origin	Label / IFU	EU
0040	MDR Annex I 23.2, 23.4	Where applicable, device label and IFU shall include an indication device contains: - Substances which are carcinogenic, mutagenic or toxic to reproduction - Substances having endocrine-disrupting properties	Label / IFU	EU
0041	MDR Annex I 23.2	The label shall include the lot or serial number of the device along with the appropriate symbol.	Label	EU
0042	MDR Annex I 23.2, 23.3	The sterile barrier packaging shall include the date of manufacture (YYYY-MM).	Label	EU
0043	MDR Annex I 23.2, 23.3	The label and sterile barrier packaging shall include the expiration date (YYYY-MM), if applicable. Otherwise the label shall include the date of manufacture.	Label	EU
0044	MDR Annex I 23.2	The label shall include the UDI number.	Label	EU

 pathway	Label Content Design Requirements		Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 5 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0045	MDR Annex I 23.2, 23.4	The label and IFU shall indicate any storage and / or handling conditions that apply.	Label / IFU	EU
0046	MDR Annex I 23.3	The sterile barrier packaging shall be indicated it is the sterile barrier.	Label	EU
0047	MDR Annex I 23.2, 23.3, 23.4	The label, sterile barrier packaging, and IFU shall indicate its method and state of sterility that applies, with the appropriate symbol.	Label / IFU	EU
0048	MDR Annex I 23.2	The label shall include any immediate warnings or precautions that apply.	Label	EU
0049	MDR Annex I 23.2, 23.4	The label and IFU shall indicate the device is for single use with the appropriate symbol, if that applies.	Label / IFU	EU
0050	MDR Annex I 23.2	The label shall indicate the device is a reprocessed single-use device, and the number of reprocessing cycles performed any limitations that apply.	Label	EU
0051	MDR Annex I 23.2	The label and sterile barrier packaging shall indicate the device is "custom-made" if applicable.	Label	EU
0052	MDR Annex I 23.2	The label shall indicate the device is a medical device.	Label	EU
0053	MDR Annex I 23.2, 23.3	The label and sterile barrier packaging shall indicate the device is "exclusively for clinical investigation" if applicable.	Label	EU
0054	MDR Annex I 23.2, 23.4	The label and IFU shall indicate the device is composed of / combination of substances intended to be introduced / applied.	Label / IFU	EU
0055	MDR Annex I 23.2	The label shall indicate the serial number for implantable devices, if applicable.	Label	EU
0056	MDR Annex I 23.3	The sterile barrier packaging shall include an instruction to check the IFU if the sterile packaging is damaged or unintentionally opened before use.	Label	EU
0057	MDR Annex I 23.4	Instructions for Use includes device's intended purpose with a clear specification of indications, contra-indications, the patient target group or groups, and of the intended users, as appropriate.	IFU	EU
0058	MDR Annex I 23.4	Instructions for Use includes where applicable, a specification of the clinical benefits to be expected.	IFU	EU
0059	MDR Annex I 23.4	Instructions for Use includes where applicable, links to the summary of safety and clinical performance.	IFU	EU
0060	MDR Annex I 23.1, 23.4	Instructions for Use includes the performance characteristics of the device.	IFU	EU
0061	MDR Annex I 23.4	Instructions for Use includes where applicable, information allowing the healthcare professional to verify if the device is suitable and select the corresponding software and accessories.	IFU	EU
0062	MDR Annex I 23.4	Instructions for Use includes where applicable, information allowing the healthcare professional to verify if the device is suitable and select the corresponding software and accessories.	IFU	EU
0063	MDR Annex I 23.1, 23.4	Instructions for Use includes any residual risks, limitations, contra-indications, precautions, warnings and any undesirable side-effects, including information to be conveyed to the patient in this regard.	IFU	EU

N/A

N/A

	Label Content Design Requirements		Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>		Rev.: C
			Page 6 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0064	MDR Annex I 23.4	Instructions for Use includes specifications the user requires to use the device appropriately (e.g. degree of accuracy).	IFU	EU
0065	MDR Annex I 23.4	Instructions for Use includes details of any preparatory treatment or handling of the device before it is ready for use or during its use, such as sterilisation, final assembly, calibration, etc., including the levels of disinfection required to ensure patient safety and all available methods for achieving those levels of disinfection.	IFU	EU
0066	MDR Annex I 23.4	Instructions for Use includes any requirements for special facilities, or special training, or particular qualifications of the device user and/or other persons.	IFU	EU
0067	MDR Annex I 23.4	Instructions for Use includes the information needed to verify whether the device is properly installed and is ready to perform safely and as intended.	IFU	EU
0068	MDR Annex I 23.4	Instructions for Use includes instructions in the event of the sterile packaging being damaged or unintentionally opened before use, if supplied sterile.	IFU	EU
0069	MDR Annex I 23.4	Instructions for Use includes appropriate instructions for sterilization, if applicable.	IFU	EU
0070	MDR Annex I 23.4	Instructions for Use includes if the device is reusable, information on the appropriate processes for allowing reuse, including cleaning, disinfection, packaging and, where appropriate, the validated method of re-sterilization.	IFU	EU
0071	MDR Annex I 23.4	Instructions for Use includes an indication, if appropriate, that a device can be reused only if it is reconditioned under the responsibility of the manufacturer to comply with the general safety and performance requirements.	IFU	EU
0072	MDR Annex I 23.4	Instructions for Use includes information on known characteristics and technical factors known to the manufacturer that could pose a risk if the device were to be re-used. Applied to single use devices.	IFU	EU
0073	MDR Annex I 23.4	Instructions for Use includes information to identify other devices and equipment in order to obtain a safe combination, and / or any known restrictions. Applies to devices intended for use together with other devices and/or general-purpose equipment.	IFU	EU
0074	MDR Annex I 23.4	Instructions for Use includes detailed information as to the nature, type and where appropriate, the intensity and distribution of the emitted radiation and the means of protecting the patient, user, or other person from unintended radiation during use of the device. Applies to devices that emit radiation for medical purposes.	IFU	EU
0075	MDR Annex I 23.4	Instructions for Use includes information that allows the user and/or patient to be informed of any warnings, precautions, contra- indications, measures to be taken and limitations of use regarding the device.	IFU	EU
0076	MDR Annex I 23.4	Instructions for Use includes the overall qualitative and quantitative information on the materials and substances to which patients can be exposed with implantable devices.	IFU	EU
0077	MDR Annex I 23.4	Instructions for Use includes warnings or precautions to be taken in order to facilitate the safe disposal of the device, its accessories and the consumables used with it, if any.	IFU	EU

	Label Content Design Requirements		Doc No.: WI-0009
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: C
			Page 7 of 7

#	Reference	Requirement Description	Device / Label / IFU	Markets
0078	MDR Annex I 23.4	Instructions for Use includes the circumstances in which the user should consult a healthcare professional, if intended for use by lay persons.	IFU	EU
0079	MDR Annex I 23.4	Instructions for Use includes date of issue of the instructions for use or, if they have been revised, date of issue and identifier of the latest revision of the instructions for use.	IFU	EU
0080	MDR Annex I 23.4	Instructions for Use includes notice to the user and/or patient that any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.	IFU	EU
0081	MDR Annex I 23.4	Instructions for Use includes information to be supplied to the patient with an implanted device.	IFU	EU
0082	MDR Annex I 23.4	Instructions for Use includes, for devices that incorporate electronic programmable systems, including software, or software that are devices in themselves, minimum requirements concerning hardware, IT networks characteristics and IT security measures, including protection against unauthorized access, necessary to run the software as intended.	IFU	EU
0083	MDR Article 20#5, Annex V	CE conformity marking shall consist of the initials "CE" followed by the identification number of the notified body responsible for the device's conformity assessment. The various components of the CE marking must have substantially the same vertical dimension, which may not be less than 5 mm.	Label / IFU	EU
0084	FDA (1768)	If a manufacturer elects to include a statement in medical product labeling indicating that natural rubber latex or synthetic derivatives of natural rubber latex were not used as materials in the manufacture of their medical product and container, FDA recommends the use of the statement "Not made with natural rubber latex." If this statement is made without any qualification, it would apply to the medical product, its container, and any packaging.	Label / IFU	US
0085	ASTM F2503	For devices which are brought into the magnetic resonance (MR) environment, permanently mark the item as MR Safe, MR Conditional, or MR Unsafe using the icons with color or black and white options. Refer to ASTM F2503 TABLE 1 and TABLE 2 for details.	Device / Label / IFU	US

8.0 Records

8.1 Label Content documents will be processed per QP-002. Refer to QP-003 for storage location and retention requirements.

9.0 Flowchart – N/A

10.0 Change History

Revision	Change Description
A	Initial release. This procedure will replace procedure QP-7.09
B	Update header, references, and format of requirements (table).
C	Remove reference to FDA 1150 Guidance

	Label Content Design Requirements	Doc No.: WI-0009
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.</i>	Rev.: C
		Page 1 of 7

1.0 Purpose

- 1.1 The purpose of this document is to indicate the requirements for labels and labeling associated with medical devices.

2.0 Scope

- 2.1 This applies to labeling that is considered the device labels (primary and secondary packages) and the Instructions for Use (IFU)
- 2.2 The requirements indicated in this document are associated with regulatory requirements associated with FDA and EU.
- 2.3 The procedure does not apply to promotions materials.

3.0 References

- 3.1 21 CFR 801 – Labeling
- 3.2 21 CFR 812 – Investigation Device Exemption
- 3.3 ASTM 2503-20 – Standard Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment
- 3.4 ISO 15223-1 – Medical devices – Symbols to be used with medical device labels, labeling, and information to be supplied – Part 1: General requirements
- 3.5 FDA (1768) Guidance – Recommendations for Labeling Medical Products to Inform Users that the Product or Product Container is not Made with Natural Rubber Latex, December 2, 2014
- 3.6 QP-002 – Document Control
- 3.7 QP-003 – Quality Records

4.0 Definitions

- 4.1 Label – Display of written, printed, or graphic matter upon the immediate container.
- 4.2 Labeling – All labels and other written, printed, or graphic matter upon any article or any of its containers or wrapper, or accompanying such article at any time while a device is held for sale after shipment or delivery for shipment in interstate commerce.
- 4.3 Investigational device – A device that is the object of a clinical investigation or research involving one or more human subjects to determine its safety and effectiveness.
- 4.4 Unique Device Identifier (UDI) – An identifier that adequately identifies a device through its distribution.

5.0 Responsibilities