

Pathway Medtech, LLC.	Type:	Test Report	Doc. No.:	TR-SLQ-0015	Revision:	B
DV Testing, SILQ ClearTract Foley Catheter						

Pathway Medtech, LLC.

Test Report, DV Testing, SILQ
ClearTract Foley Catheter

Prepared For:

SILQ

Pathway Medtech, LLC.	Type:	Test Report	Doc. No.:	TR-SLQ-0015	Revision:	B
DV Testing, SILQ ClearTract Foley Catheter						

SIGNATURE / APPROVAL PAGE:

Approver(s):	Signatures:	Date:
Author / Engineering:	 <u>Will Hail (May 1, 2024 15:07 PDT)</u>	<u>05/01/2024</u>
Quality:	 <u>Aaron Rogers (May 1, 2024 15:36 PDT)</u>	<u>05/01/2024</u>
Customer:	 <u>Ethan Rao (May 1, 2024 15:03 PDT)</u>	<u>05/01/2024</u>

CHANGE HISTORY

Revision	Change Description
A	Initial Release
B	Updated to new documentation template, added data for T=36 testing, small formatting and grammatical changes.

TABLE OF CONTENTS

Signature / Approval Page:.....2

Change History.....2

Table of Contents.....3

1.0 Purpose4

2.0 Scope.....4

3.0 Executive Summary.....4

4.0 Acceptance Criteria.....4

5.0 Material and Equipment6

6.0 Deviations.....6

7.0 Results and Discussion6

8.0 Conclusion15

9.0 References.....16

10.0 Attachments.....16

Pathway Medtech, LLC.	Type:	Test Report	Doc. No.:	TR-SLQ-0015	Revision:	B
DV Testing, SILQ ClearTract Foley Catheter						

1.0 PURPOSE

- 1.1 The purpose of this report is to document the results for the following:
- 1.1.1 Design verification evaluations associated with the SILQ ClearTract Foley Catheter.
 - 1.1.2 Packaging verification evaluations associated with the sterile barrier.
 - 1.1.3 Stability / Shelf Life evaluations.

2.0 SCOPE

- 2.1 This report applies to the SILQ ClearTract Foley Catheters for the following catheter SKUs, pouches, and aging intervals.

Table 1 Report Scope

SKU	Aging Interval
8 Fr, 3-5 mL, 2-Way	T=0, T=12, T=36
24 Fr, 30 mL, 2-Way	T=0, T=12, T=36
SKU	Aging Interval
Steripax (04A 0319F 4203)	T=0, T=12, T=36

- 2.2 Though 3-way and Coude tip styles were identified in the test protocols, they are not featured in this revision of the report.
- 2.3 The 3-year shelf life applies to all sizes and volumes including and between 8 Fr / 3-5 mL to 24 Fr / 30 mL 2-Way ClearTract Catheters.
- 2.4 Test samples were fully manufactured, EO sterilized, and exposed through distribution conditioning and simulated distribution.

3.0 EXECUTIVE SUMMARY

Pathway Medtech generated foley catheters for SILQ to be used for DV testing at the unaged, T=12, and T=36 accelerated aging intervals. Pathway Medtech packed these samples with the new Steripax pouch in order to qualify it for use with the catheters. All units went through the validated sterilization cycle and simulated distribution. Upon returning, Pathway Medtech placed the units designated for aging in the oven, set aside the units designated for real time aging, and started testing T=0 samples. After aging, Pathway Medtech tested the T=12 and T=36 samples. All units, catheters and pouches, satisfied all requirements laid forth in the corresponding protocols. As a result, the catheter SKUs 8 Fr, 3-5 mL, 2-Way and 24 Fr, 30 mL, 2-Way are qualified with a 3-year shelf life. Additionally, the Steripax pouch has also been qualified with a 3-year shelf life.

4.0 ACCEPTANCE CRITERIA

Table 2 Packaging Testing Acceptance Criteria

Term	Meaning
Seal Integrity - Visual Inspection	No wrinkles, creases, voids, channels, or bubbles in seals. Seals must not be over-sealed resulting in melting, brittle material, or material delamination or separation. Seals appear intact and continuous for the entire seal area with a minimum seal width of 9 mm. Note: If no defects are observed from the indicated sample size, then a confidence/ reliability level of 95%/95% will be achieved.

Term	Meaning
Gross Leak (Bubble) Test	No evidence of Gross leaks between the sealed edges or through either flexible barrier when subjected to the established (at or above) test pressure for a minimum of 30 seconds. Note: If no defects are observed from the indicated sample size, then a confidence/ reliability level of 95%/95% will be achieved.
Seal Strength	The peel strength of each 1" section shall be ≥ 1.0 lb. for seals created by Pathway and ≥ 0.6 for seals created by the manufacturer Note: If no defects are observed from the indicated sample size, then a confidence/ reliability level of 95%/95% will be achieved.

Table 3 Design Verification Testing Acceptance Criteria

Term	Meaning																				
Balloon Size and Shaft Size	Tip OD: Nominal ± 1 Fr. Uninflated Balloon OD: Nominal ± 4 Fr. Shaft OD: Nominal ± 1 Fr.																				
Verification of Non-Interconnectable Connector	The catheter connection shall not appear to provide a secure connection when assembled to any surface of the components and shall easily disengage from each of the Small-Bore Connectors.																				
Flow Rate through Drainage Lumen and Irrigation Lumen	<table border="1"> <tr> <td colspan="2">Flow Rate through Drainage Lumen - Shall have a minimum average flow rate, per size, greater than or equal to the values in below.</td></tr> <tr> <td>Size</td><td>Flow Rate [cm^3 / min]</td></tr> <tr> <td>8 Fr</td><td>15</td></tr> <tr> <td>10 Fr</td><td>30</td></tr> <tr> <td>12 Fr</td><td>70</td></tr> <tr> <td>14 Fr – 24 Fr</td><td>100</td></tr> <tr> <td colspan="2">Flow Rate through Irrigation Lumen - Shall have a minimum average flow rate, per size, greater than or equal to the values in below.</td></tr> <tr> <td>Size</td><td>Flow Rate [cm^3 / min]</td></tr> <tr> <td>14 Fr – 20 Fr</td><td>25</td></tr> <tr> <td>22 Fr - 24 Fr</td><td>30</td></tr> </table>	Flow Rate through Drainage Lumen - Shall have a minimum average flow rate, per size, greater than or equal to the values in below.		Size	Flow Rate [cm^3 / min]	8 Fr	15	10 Fr	30	12 Fr	70	14 Fr – 24 Fr	100	Flow Rate through Irrigation Lumen - Shall have a minimum average flow rate, per size, greater than or equal to the values in below.		Size	Flow Rate [cm^3 / min]	14 Fr – 20 Fr	25	22 Fr - 24 Fr	30
Flow Rate through Drainage Lumen - Shall have a minimum average flow rate, per size, greater than or equal to the values in below.																					
Size	Flow Rate [cm^3 / min]																				
8 Fr	15																				
10 Fr	30																				
12 Fr	70																				
14 Fr – 24 Fr	100																				
Flow Rate through Irrigation Lumen - Shall have a minimum average flow rate, per size, greater than or equal to the values in below.																					
Size	Flow Rate [cm^3 / min]																				
14 Fr – 20 Fr	25																				
22 Fr - 24 Fr	30																				
Inflated Balloon Response to Traction	The entire balloon of catheters shall not pass into or through the funnel barrel and shall not rupture.																				
Balloon Integrity (Resistance to Rupture)	The inflation balloon must be inflated easily with distilled or deionized water to its indicated labeled volume (i.e., 5ml, 10ml) without showing any evidence of breakage throughout the 7-day test period.																				
Deflation Reliability (Failure to Deflate)	The balloon shall deflate to within 4 French sizes of the label French size within 15 minutes or be otherwise manipulated to effect drainage after the 7-day immersion period at its labeled fill volume (i.e., 5ml, 10ml).																				
Balloon Volume Maintenance	The catheter shall maintain its labeled volume (i.e., 5ml, 10ml) for 24-hour period.																				

5.0 MATERIAL AND EQUIPMENT

Equipment No.	Description	Calibration Due Date
ES-0012-05	Caliper	06/23/2024
ES-0028-01	Force Gauge, 10 lb.	01/06/2024
*ES-0131-01	Motorized Test Stand	08/09/2023 / 08/09/2024
ES-0003-05	Stopwatch	09/07/2023
ES-0003-06	Stopwatch	06/20/2025
ES-0017-01	Scale, Analytical Balance	4/14/2024
ES-0148-01	Pressure Gauge	11/10/2023
ES-0119-02	Datalogging Thermometer	12/08/2023
ES-0119-03	Datalogging Thermometer	12/08/2023
ES-0149-02	Oven	N/A

*Test stand was calibrated in between T=0 and T=12 tests.

6.0 DEVIATIONS

None.

7.0 RESULTS AND DISCUSSION

7.1 Traceability



7.1.1 The following table details the lot number for each catheter SKU and pouch assessed in this report as well as the accelerated aging interval under which they fall.

Table 4 Traceability Information

Part No.	Lot No.	Qty	Aging Interval	Sterilization Load	Simulated Distribution and Conditioning	Tests Performed
SLQ-210803	SLQ-80140228231E	60	T=0	Cycle ID: PWI-7374 Run ID: 1162-9381	Report Number: WP-14055A	TP-SLQ-0014
SLQ-210803	SLQ-80140228232E	60	T=12			
SLQ-210803	SLQ-80140228233E	60	T=36			
SLQ-212430	SLQ-80080228231E	180	T=0, T=12, T=36			
SLQ-5013	18768	180	T=0, T=12, T=36			TP-SLQ-0014

7.2 Stability Verification and Discussion

7.2.1 Pouched catheters were placed in oven ES-0149-02 in order to undergo accelerated aging. The T=12 set of catheters were in the oven from 6/15/2023 to 07/24/2023 for a total of 39 days. The T=36 set of catheters were in the oven from 6/15/2023 to 10/4/2023 for a total of 110 days. Temperature was recorded every 30 minutes over that duration. Charts of the average daily temperature over that period is shown below.

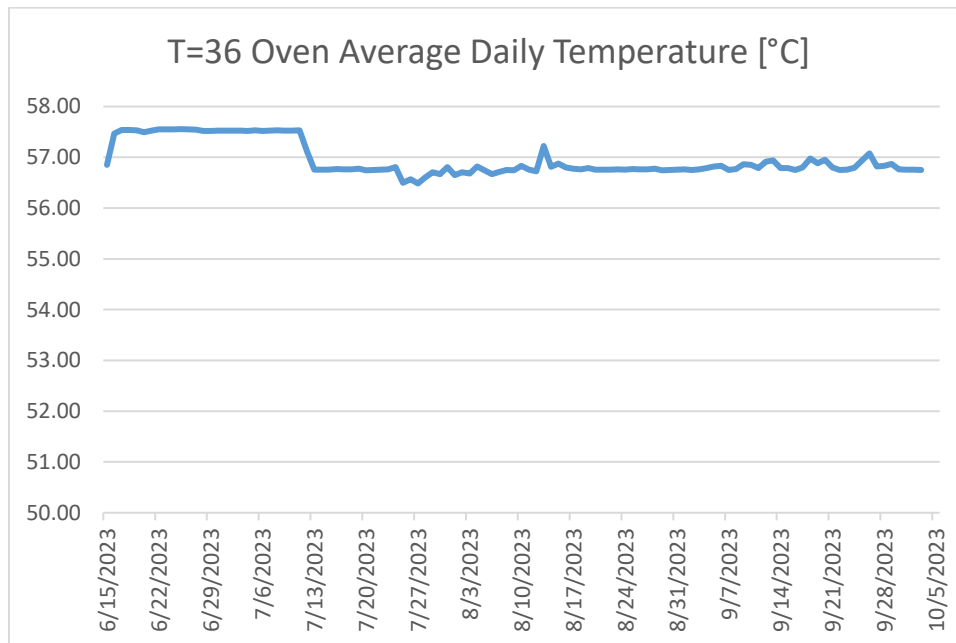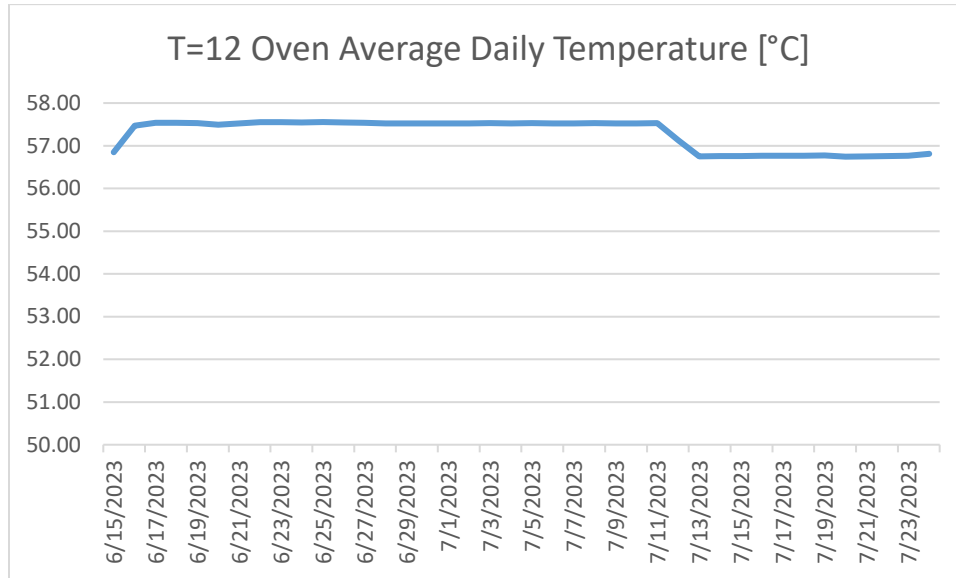

7.2.2 Data collected from ES-0119-03 showed that the product inside of the oven stayed at or above 55°C for 39 days for the T=12 samples and 110 days for the T=36 samples, with the ambient room temperature of 22°C. As a result, there is confidence that product designated as T=12 and T=36 went through sufficient accelerated aging to represent a one year shelf life and three year shelf life respectively.

- 7.2.3 Using the equation below, the accelerated aging time for each data point was found and summed up. As a result, the accelerated aging time for the T=12 samples was calculated to be 461 days, which is equivalent to 15 months. Additionally, the accelerated aging time for the T=36 samples was calculated to be 1,252 days which is equivalent to 41 months.

$$AAF = Q_{10}^{\left\{\frac{(T_{AA}-T_{RT})}{10}\right\}}$$

where:

T_{AA} = accelerated aging temperature (°C)

T_{RT} = ambient temperature (°C)

Q_{10} = reaction rate coefficient (2.0 is a common and conservative estimate)

7.3 Packaging Verification Results and Discussion

7.3.1 Summary

- 7.3.1.1 Prior to testing, the product went through conditioning per ASTM D4332 and distribution simulation D4169 at Westpak, detailed under Westpak report WP-14055A.
- 7.3.1.2 Package testing was conducted at Pathway Medtech in Santee, CA. All packaging had gone through conditioning and simulated distribution prior to testing. There were two subsets of pouches, T=0, T=12, and T=36.
- 7.3.1.3 All test measurement equipment used in this evaluation was calibrated and documented on the appropriate data collection sheet.
- 7.3.1.4 Steripax (04A 0319F 4203) pouches were assessed with 3-Way, 24 Fr, 30 mL (REF: SLQ-312430) inside of them.

7.3.2 Pouch Visual (Ref: ASTM F1886/F1886M-16)

- 7.3.2.1 Visual inspections were performed to ensure the pouch seal created when utilizing ES-0099-01 Accu-Seal Impulse Sealer were complete with no wrinkles, creases, voids, channels, or bubbles in seals. Seals must not be over-sealed resulting in melts, brittle material, or material delamination or separation. Seals appear intact and continuous for the entire seal area with a minimum seal width of 9 mm.
- 7.3.2.2 Recorded results indicate that all samples passed the visual inspection with no evidence of the defects.

Table 5 Pouch Visual Test Results

Pouch	Test Group	Pass Quantity	Fail Quantity	Result
Steripax (04A 0319F 4203)	T=0	58	0	PASS
	T=12	58	0	PASS
	T=36	58	0	PASS

7.3.3 Bubble Test (Ref: ASTM F2096-11 and TM-0005)

Pathway Medtech, LLC.	Type:	Test Report	Doc. No.:	TR-SLQ-0015	Revision:	B
DV Testing, SILQ ClearTract Foley Catheter						

- 7.3.3.1 The testing was conducted per the set up in TM-0005. The resulting test pressure was 0.6 PSI. This value was verified throughout the testing with a calibrated gauge to ensure the pressure remained at or above this value.
- 7.3.3.2 The results of the testing indicate all samples were able to be pressurized to the test pressure with no evidence of gross leaks through any of the seals or through the pouch material itself

Table 6 Bubble Test Results

Pouch	Test Group	Pass Quantity	Fail Quantity	Result
Steripax (04A 0319F 4203)	T=0	58	0	PASS
	T=12	58	0	PASS
	T=36	58	0	PASS

7.3.4 Peel Strength (Ref: ASTM F88/F88M:2015 and TM-0003)

- 7.3.4.1 All Tests were executed utilizing a calibrated force gauge and motorized test stand. The pull speed was set to 12 inches / min for all samples.
- 7.3.4.2 Recorded results indicate that all samples passed the Peel Test.

Table 7 Peel Strength Test Results

Acceptance Criteria	Side "A"	Side "B"	Side "C"	Side "D"
	≥ 1.000 lb.	≥ 0.600 lb.	≥ 0.600 lb.	≥ 0.600 lb.
T=0				
Min	3.682	0.918	1.056	1.000
Max	7.224	2.054	2.544	2.201
Mean	5.645	1.491	1.875	1.545
Standard Deviation	0.933	0.371	0.408	0.361
T=12				
Min	2.264	0.920	1.345	1.038
Max	5.412	2.054	2.411	2.388
Mean	3.119	1.519	1.842	1.399
Standard Deviation	0.774	0.333	0.296	0.315
T=36				
Min	2.114	1.422	1.135	1.510
Max	4.544	5.524	3.510	2.966
Mean	2.846	2.471	1.974	2.252
Standard Deviation	0.523	0.984	0.481	0.365

7.4 Design Verification Results and Discussion

7.4.1 Summary

7.4.1.1 Performance testing of the SILQ ClearTract Foley Catheters was conducted at Pathway Medtech in Santee, CA. All product was treated with SILQ's polymer solution. Test samples were fully manufactured, EO sterilized and exposed through distribution conditioning and simulated distribution. There were two subsets of catheters, T=0, T=12, and T=36.

7.4.1.2 All test measurement equipment used in this evaluation was calibrated and documented on the appropriate data collection sheet.

7.4.1.3 The data collection sheets provide product identification of the sterile catheters that were evaluated as part of the study.

7.4.2 Balloon Size and Shaft Size

7.4.2.1 Thirty (30) units were measured at the catheter tip, balloon, and shaft. The results are detailed in the following table.

Table 8 Balloon and Shaft Size Test Results

Acceptance Criteria	8 Fr, 3-5 mL, 2-Way			24 Fr, 30 mL, 2-Way		
	Tip OD: 0.092 in. – 0.118 in. Uninflated Balloon OD: 0.053 in. – 0.158 in. Shaft OD: 0.092 in. – 0.118 in.			Tip OD: 0.302 in. – 0.328 in. Uninflated Balloon OD: 0.263 in – 0.367 in. Shaft OD: 0.302 in. – 0.328 in.		
	T=0					
T=0	Tip OD [in.]	Balloon OD [in.]	Shaft OD [in.]	Tip OD [in.]	Balloon OD [in.]	Shaft OD [in.]
Min	0.107	0.139	0.107	0.309	0.345	0.311
Max	0.114	0.146	0.115	0.328	0.367	0.327
Mean	0.110	0.143	0.111	0.316	0.361	0.320
Standard Deviation	0.002	0.002	0.002	0.005	0.005	0.004
T=12						
T=12	Tip OD [in.]	Balloon OD [in.]	Shaft OD [in.]	Tip OD [in.]	Balloon OD [in.]	Shaft OD [in.]
Min	0.107	0.141	0.107	0.302	0.352	0.308
Max	0.115	0.147	0.115	0.323	0.367	0.324
Mean	0.111	0.144	0.111	0.314	0.362	0.317
Standard Deviation	0.002	0.002	0.002	0.005	0.004	0.004
T=36						
Min	0.093	0.129	0.099	0.302	0.312	0.302
Max	0.118	0.157	0.118	0.328	0.367	0.327
Mean	0.105	0.141	0.110	0.310	0.354	0.313
Standard Deviation	0.007	0.006	0.006	0.008	0.011	0.007

7.4.2.2 All catheter sizes evaluated were within the specified tolerance as it relates to the labeled size. Acceptance criteria was met. Refer to attached data collection sheets for additional information.

7.4.3 Verification Non-Interconnectable Connector

7.4.3.1 Fifty-eight (58) units from each size was evaluated and the small-bore connector was observed to not provide a secure connection with the catheter and would disengage under its own weight. The results are indicated in the following table:

Table 9 Verification Non-Interconnectable Connector Test Results

SKU	Test Group	Pass Quantity	Fail Quantity	Result
8 Fr, 3-5 mL, 2-Way	T=0	58	0	PASS
	T=12	58	0	PASS
	T=36	58	0	PASS
24 Fr, 30 mL, 2-Way	T=0	58	0	PASS
	T=12	58	0	PASS
	T=36	58	0	PASS

7.4.3.2 Acceptance criteria was met. Refer to attached data collection sheets for additional information.

7.4.4 Flow Rate through Drainage Lumen

7.4.4.1 Thirty (30) units were evaluated to assess flow rate. The volume per minute was measured on each catheter three (3) times, then averaged. The results are indicated in the following table:

Table 10 Flow Rate through Drainage Lumen Test Results

Acceptance Criteria	8 Fr, 3-5 mL, 2-Way	24 Fr, 30 mL, 2-Way
	Flow rate $\geq 15 \text{ cm}^3 / \text{min}$	Flow rate $\geq 100 \text{ cm}^3 / \text{min}$
T=0		
	Flow Rate [cm^3 / min]	Flow Rate [cm^3 / min]
Min	52.0	1003.3
Max	72.3	1043.3
Mean	57.2	1026.6
Standard Deviation	4.3	9.3
T=12		
	Flow Rate [cm^3 / min]	Flow Rate [cm^3 / min]
Min	52.3	1000.0
Max	66.3	1063.3
Mean	58.3	1026.8
Standard Deviation	2.8	18.9
T=36		

Acceptance Criteria	8 Fr, 3-5 mL, 2-Way	24 Fr, 30 mL, 2-Way
	Flow rate $\geq 15 \text{ cm}^3 / \text{min}$	Flow rate $\geq 100 \text{ cm}^3 / \text{min}$
	Flow Rate [cm^3 / min]	Flow Rate [cm^3 / min]
Min	44.3	433.3
Max	60.0	958.3
Mean	52.9	665.8
Standard Deviation	4.3	111.8

7.4.4.2 All catheter evaluated resulted in a flow rate through the drainage lumen greater than $100 \text{ cm}^3/\text{min}$ for the 24 Fr SKU and $15 \text{ cm}^3/\text{min}$ for 8 Fr SKU. Acceptance criteria was met. Refer to attached data collection sheets for additional information.

7.4.5 Inflated Balloon Response to Traction

7.4.5.1 Fifty-eight (58) units were evaluated with a fill volume of 5ml since it is the worst case for this evaluation. All units evaluated did not rupture or pass into or through the funnel barrel in the static or impact load evaluation.

Table 11 Inflated Balloon Response to Traction Test Results

SKU	Test Group	Static Test		Result	Dynamic Test		Result
		Pass Quantity	Fail Quantity		Pass Quantity	Fail Quantity	
8 Fr, 3-5 mL, 2-Way	T=0	58	0	PASS	58	0	PASS
	T=12	58	0	PASS	58	0	PASS
	T=36	58	0	PASS	58	0	PASS
24 Fr, 30 mL, 2-Way	T=0	58	0	PASS	58	0	PASS
	T=12	58	0	PASS	58	0	PASS
	T=36	58	0	PASS	58	0	PASS

7.4.5.2 Acceptance criteria was met. Refer to attached data collection sheets for additional information.

7.4.6 Balloon Integrity (Resistance to Rupture)



7.4.6.1 Fifty-eight (58) units were evaluated. The catheters were filled according to their maximum labeled fill volumes. All units evaluated were not observed to have leaks or breakage for the test duration.


Table 12 Balloon Integrity (Resistance to Rupture) Test Results

SKU	Test Group	Pass Quantity	Fail Quantity	Result
8 Fr, 3-5 mL, 2-Way	T=0	58	0	PASS
	T=12	58	0	PASS
	T=36	58	0	PASS
	T=0	58	0	PASS

SKU	Test Group	Pass Quantity	Fail Quantity	Result
24 Fr, 30 mL, 2-Way	T=12	58	0	PASS
	T=36	58	0	PASS

7.4.6.2 The two charts below detail the average daily temperature of the simulated urine bath of the 7 day period for both the T=0, T=12, and T=36 samples. The data was collected with ES-0119-02. The charts show the average temperature does not fall outside the tolerances given in the test protocol. As a result, there is confidence that the catheter sample sets were exposed to the required temperature bath for the required time duration.

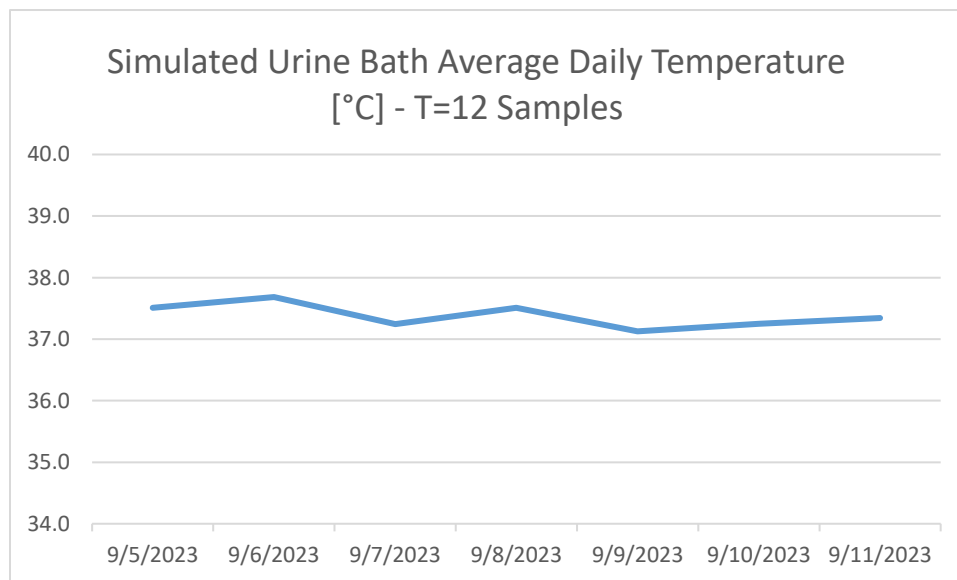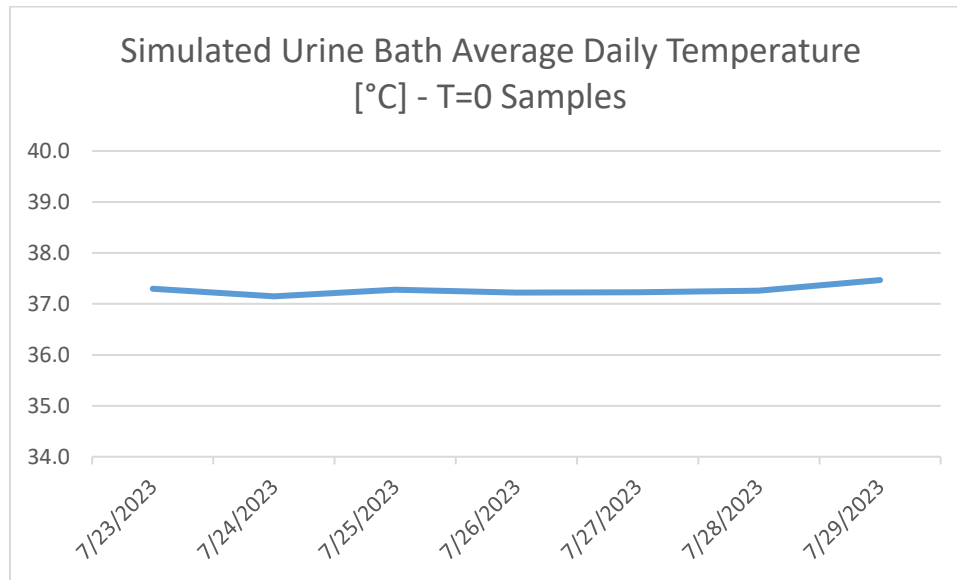

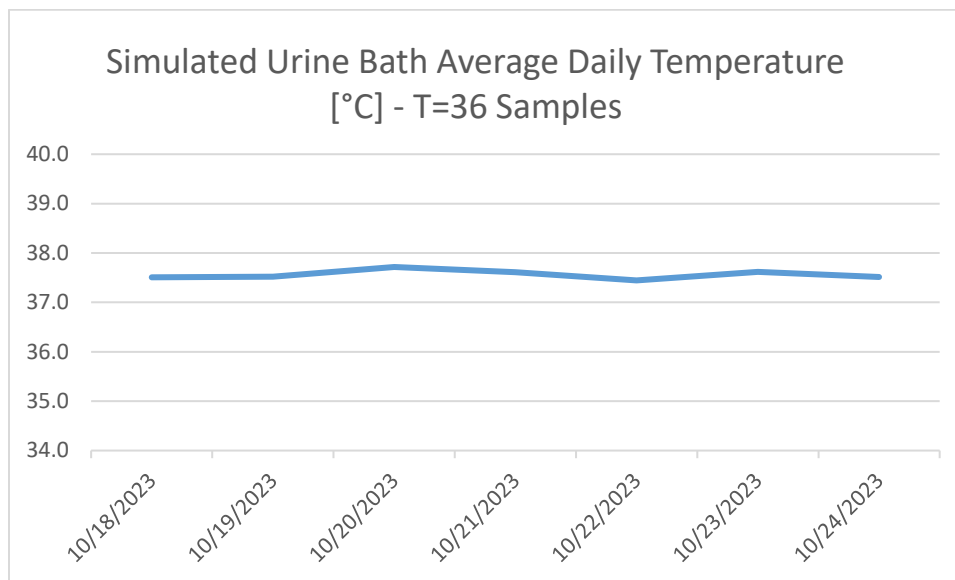

7.4.7 Deflation Reliability (Failure to Deflate)

7.4.7.1 Fifty-eight (58) units were evaluated and filled according to their maximum labeled fill volumes.

Table 13 Deflation Reliability (Failure to Deflate) Test Results

Acceptance Criteria	8 Fr, 3-5 mL, 2-Way	24 Fr, 30 mL, 2-Way
	Deflated Balloon OD: 0.053 in. – 0.158 in.	Deflated Balloon OD: 0.263 in – 0.367 in.
T=0		
	OD [in.]	OD [in.]
Min	0.135	0.359
Max	0.144	0.367
Mean	0.139	0.363
Standard Deviation	0.003	0.002
T=12		
	OD [in.]	OD [in.]
Min	0.141	0.352
Max	0.147	0.367
Mean	0.144	0.362
Standard Deviation	0.001	0.004
T=36		
	OD [in.]	OD [in.]
Min	0.140	0.354
Max	0.147	0.367

Acceptance Criteria	8 Fr, 3-5 mL, 2-Way	24 Fr, 30 mL, 2-Way
	Deflated Balloon OD: 0.053 in. – 0.158 in.	Deflated Balloon OD: 0.263 in – 0.367 in.
Mean	0.143	0.364
Standard Deviation	0.002	0.002

7.4.7.2 None of the catheters evaluated were observed with a rupture, leak, or loss of volume. The deflated catheters all measured within the specified 4Fr from the labeled size. Acceptance criteria was met. Refer to attached data collection sheets for additional information.

7.4.8 Balloon Volume Maintenance

7.4.8.1 Fifty-eight (58) units of each SKU were filled according to their maximum labeled volume. All units evaluated maintained the delivered fill volume for the test duration. Acceptance criteria was met. Refer to attached data collection sheets for additional information.

Table 14 Balloon Volume Maintenance Test Results

SKU	Test Group	Pass Quantity	Fail Quantity	Result
8 Fr, 3-5 mL, 2-Way	T=0	58	0	PASS
	T=12	58	0	PASS
	T=36	58	0	PASS
24 Fr, 30 mL, 2-Way	T=0	58	0	PASS
	T=12	58	0	PASS
	T=36	58	0	PASS

8.0 CONCLUSION

8.1 The 2-Way 8 Fr, 3 – 5 mL and 2-Way 24 Fr, 30 mL SILQ ClearTract Foley Catheter met all the requirements detailed in TP-SLQ-0014 and TP-SLQ-0015 at the unaged (T=0), 1-year (T=12), and 3-year (T=36) accelerated aging time interval.

8.2 The Steripax (04A 0319F 4203) pouches met all requirements detailed in TP-SLQ-0015 and TP-SLQ-0018 at the unaged (T=0), 1-year (T=12), and 3-year (T=36) accelerated aging time interval. These pouches were tested against the most difficult catheter to maintain the sterile barrier, 3-Way 24 Fr. Since the pouches passed with that SKU, there can be confidence it will maintain the sterile barrier with the other two 24 Fr SKUs, 2-Way, and 2-Way Coude Tip, since they are a less challenging packaging configuration. As a result, there can be confidence the Steripax pouches will maintain the sterile barrier with all 24 Fr catheter SKUs.

8.2.1 Based on the above, the SILQ ClearTract Foley Catheter product family is qualified to be placed in the Steripax pouches with a stated shelf life of T=36.

8.3 The ClearTract Catheter has demonstrated a 3-year shelf life by meeting its performance and packaging requirements for 8Fr / 3-5mL to 24Fr / 30mL. Further, these results are applicable to all sizes and volumes including and between 8Fr / 3-5mL to 24Fr / 30mL on the basis they have the same materials of construction and manufacturing processes.

Pathway Medtech, LLC.	Type:	Test Report	Doc. No.:	TR-SLQ-0015	Revision:	B
DV Testing, SILQ ClearTract Foley Catheter						

8.4 With all acceptance criteria met, the results achieve a 95% confidence / 95% reliability level.

9.0 REFERENCES

Document Number	Title
Internal Procedures	
QP-002	Document Control
QP-030	Statistical Techniques
Test Protocols	
TP-SLQ-0014 rev A	Test Protocol, Design Verification, SILQ Foley Catheters
TP-SLQ-0015 rev A	Test Protocol, Stability / Shelf Life, SILQ Foley Catheters
TP-SLQ-0018 Rev B	Test Protocol, Packaging Verification, SILQ Foley Catheters

10.0 ATTACHMENTS

#	Attachment	Edition	Type
1	Data Collection Sheets, 8 Fr	10/26/2023	Testing Data
2	Data Collection Sheets, 24 Fr	10/26/2023	Testing Data
3	Pouch Data Sheets	10/06/2023	Testing Data
4	Engineering Build Record: SLQ-80080228231E	05/31/2023	Device History Record
5	Engineering Build Record: SLQ-80140228231E	05/31/2023	Device History Record
6	Engineering Build Record: SLQ-80140228232E	05/31/2023	Device History Record
7	Engineering Build Record: SLQ-80140228233E	05/31/2023	Device History Record
8	Engineering Build Record: SLQ-80120228231E	05/31/2023	Device History Record
9	Sterilization (Steris) Record	05/04/2023	Sterilization Processing Record
10	Sterilization Packing List	04/21/2023	Packing List
11	Simulated Distribution (Westpak) Record	A	Lab Report
12	Westpak Packing List	05/18/2023	Packing List
13	Training Records	10/02/2023	Training Record

TR-SLQ-0015_B - Test Report, DV Testing, SILQ ClearTract Foley Catheter

Final Audit Report

2024-05-01

Created:	2024-05-01
By:	Will Hail (whail@pathwaymedtech.com)
Status:	Signed
Transaction ID:	CBJCHBCAABAAxHrefNmqKM39n_koGZgUO9IU0zDVC8F0

"TR-SLQ-0015_B - Test Report, DV Testing, SILQ ClearTract Foley Catheter" History

-  Document created by Will Hail (whail@pathwaymedtech.com)
2024-05-01 - 10:02:03 PM GMT
-  Document emailed to Will Hail (whail@pathwaymedtech.com) for signature
2024-05-01 - 10:02:07 PM GMT
-  Document emailed to Aaron Rogers (arogers@pathwaymedtech.com) for signature
2024-05-01 - 10:02:07 PM GMT
-  Document emailed to Ethan Rao (ethanr@silq.tech) for signature
2024-05-01 - 10:02:07 PM GMT
-  Email viewed by Ethan Rao (ethanr@silq.tech)
2024-05-01 - 10:02:53 PM GMT
-  Document e-signed by Ethan Rao (ethanr@silq.tech)
Signature Date: 2024-05-01 - 10:03:51 PM GMT - Time Source: server
-  Email viewed by Will Hail (whail@pathwaymedtech.com)
2024-05-01 - 10:07:18 PM GMT
-  Document e-signed by Will Hail (whail@pathwaymedtech.com)
Signature Date: 2024-05-01 - 10:07:26 PM GMT - Time Source: server
-  Email viewed by Aaron Rogers (arogers@pathwaymedtech.com)
2024-05-01 - 10:35:13 PM GMT
-  Document e-signed by Aaron Rogers (arogers@pathwaymedtech.com)
Signature Date: 2024-05-01 - 10:36:01 PM GMT - Time Source: server

✔ Agreement completed.

2024-05-01 - 10:36:01 PM GMT