


Certificate of Conformance

Pathway Medtech LLC. certifies that the products or services indicated below have been manufactured in accordance with the requirements specified on referenced documents.

Part Number:	SLQ-211810SPT
Description:	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack
Revision:	B
Lot Number(s):	SLQ-01242025
Expiration Date:	2028-01-31
Quantity Released:	369/36.9 packs of 10

The product indicated can be transferred to a "release" inventory location and is cleared for distribution. The individual signing below authorizes the release of the product.

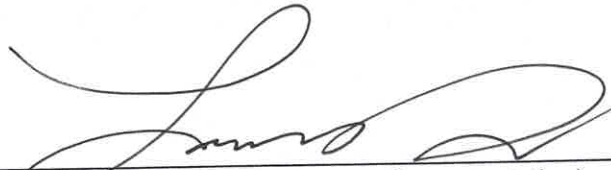


Signature (Quality Assurance Representative)

02/19/2025

Date

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
 1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY PATHWAY PRIOR TO PRODUCTION.
 2. DRAWING IS FOR REFERENCE ONLY. REFER TO THE RELEASED BILL OF MATERIALS FOR THE COMPLETE PARTS LIST.
 3. UP UNTIL SEALING THE POUCH, PRODUCT SHALL BE PACKAGED IN A CLEANROOM ENVIRONMENT. POST SEALING THE POUCH, ASSEMBLY WILL CONTINUE OUTSIDE OF THE CLEANROOM ENVIRONMENT.
2. ASSEMBLY PROCEDURE
 1. APPLY POLYMER SOLUTION TO THE CATHETER IN THE HATCHED AREA SHOWN. ALLOW THE CATHETER TO DRY OVERNIGHT.
 2. INSERT THE CATHETER INTO A POLYETHYLENE SLEEVE SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CLOSED END.
 3. PLACE THE POUCH LABEL ON THE CLEAR SIDE OF THE POUCH. WHEN VIEWING THE CLEAR SIDE OF THE POUCH FROM ABOVE, WITH THE CHEVRON TO THE RIGHT, THE LABEL TEXT SHALL READ LEFT TO RIGHT. WITH BEST EFFORT, APPLY THE LABEL SUCH THAT IT IS CENTERED.
 4. INSERT THE SLEEVE / CATHETER INTO THE POUCH SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CHEVRON END OF THE POUCH.
 5. SEAL THE POUCH.
 6. PLACE 10 SLEEVED, POUCHED, SEALED, AND LABELED CATHETERS INTO ONE SHELF CARTON. APPLY TWO SHIPPER BOX LABELS TO THE LONG, WIDER SIDE OF THE CARTON. VERIFY THAT THE TEXT READS LEFT TO RIGHT WHEN THE OPENING OF THE CARTON IS TO THE LEFT. WITH BEST EFFORT, CENTER THE LABEL ON ITS RESPECTIVE FACE OF THE CARTON.
 7. PLACE 10 SHELF CARTONS INTO ONE SHIPPER BOX. SEAL WITH PACKAGING TAPE AS NEEDED. PLACE A SHIPPER BOX LABEL IN THE UPPER LEFT CORNER ON TWO ADJACENT FACES, READING LEFT TO RIGHT.
3. DIMENSIONAL
 1. $\square$ INDICATES LABEL REFERENCE DIMENSIONS.
 2. $\square$ INDICATES CRITICAL DIMENSIONS.
4. CLEANING AND HANDLING
 1. ALL PRODUCT, UP UNTIL INSERTING THE POUCHES INTO SHELF CARTONS, SHALL BE HANDLED BY GLOVED HANDS APPROPRIATE FOR A CLEANROOM ENVIRONMENT.
5. FUNCTIONAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
 1. THE CATHETER'S ABSORBANCE VALUE AT 585 NANOMETERS SHALL BE GREATER THAN OR EQUAL TO 0.300 AU. REFER TO TM-HDX-0001.
 2. THE CATHETER'S ABSORBANCE VALUE AT 440 NANOMETERS SHALL BE LESS THAN OR EQUAL TO 0.020 AU. REFER TO TM-HDX-0002.
 3. THE PEEL STRENGTH FOR A 1 INCH WIDE SEALED SECTION OF THE STERILE BARRIER PACKAGING SHALL BE ≥ 1.00 LB WHEN TESTED PER TM-0003. IT IS NOT REQUIRED TO PLACE CATHETERS IN THESE TEST POUCHES, OR LABEL THEM. INSPECT N=3 SAMPLES, 1 AT THE BEGINNING, MIDDLE, AND END OF THE MANUFACTURING RUN.
6. VISUAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
 1. THE CATHETER SHALL BE FREE OF CURED POLYMER AND DEBRIS.
 2. THE CATHETER SHALL BE PROPERLY ORIENTED IN THE SLEEVE AND POUCH WITH THE LABEL APPLIED TO THE APPROPRIATE FACE AND IN THE CORRECT ORIENTATION.
 3. INSIDE THE STERILE BARRIER, DEVICE IS TO BE FREE OF PARTICULATES WHEN VIEWED AT 18 INCHES AWAY UNDER AMBIENT ROOM LIGHTING AND WITHOUT THE USE OF MAGNIFICATION.
 4. THE STERILE BARRIER PACKAGING SHALL NOT INCLUDE ANY OF THE FOLLOWING DEFECTS:
 1. UNSEALED AREAS, NONHOMOGENEOUS OR UNSEALED AREAS, OVERSEALED AREAS, ANROW SEALS, CHANNELS, WRINKLES / FOLDOVERS / CRACKS, OR TEARS / PINHOLES.
 5. THE SHELF CARTON SHALL HAVE THE CORRECT NUMBER OF CATHETERS INSIDE AND FACING THE CORRECT ORIENTATION. THE LABELS SHALL BE AFFIXED IN THE CORRECT LOCATION AND ORIENTATION.
 6. THE SHIPPER BOX SHALL HAVE THE CORRECT NUMBER OF SHELF CARTONS. THE LABELS SHALL BE APPLIED IN THE CORRECT LOCATION AND ORIENTATION.

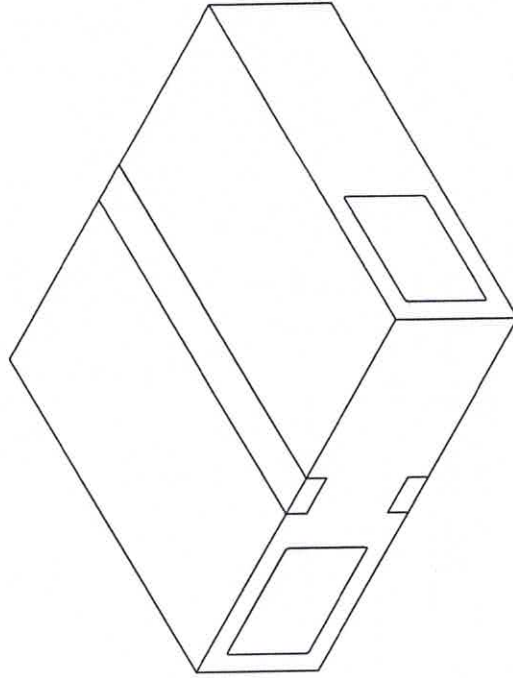

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	



TITLE:

CATHETER, SILICONE NY FOLEY, NON-STERILE, 10-PACK, MASTER DRAWING

SIZE DWG. NO. B SLQ-81XX

REV A

SCALE: 1:5

SHEET 1 OF 3


UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ±0.5
HATCHED AREAS: ±0.5
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY INC. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY INC. IS PROHIBITED.



CERTIFIED 13485

REGISTERED

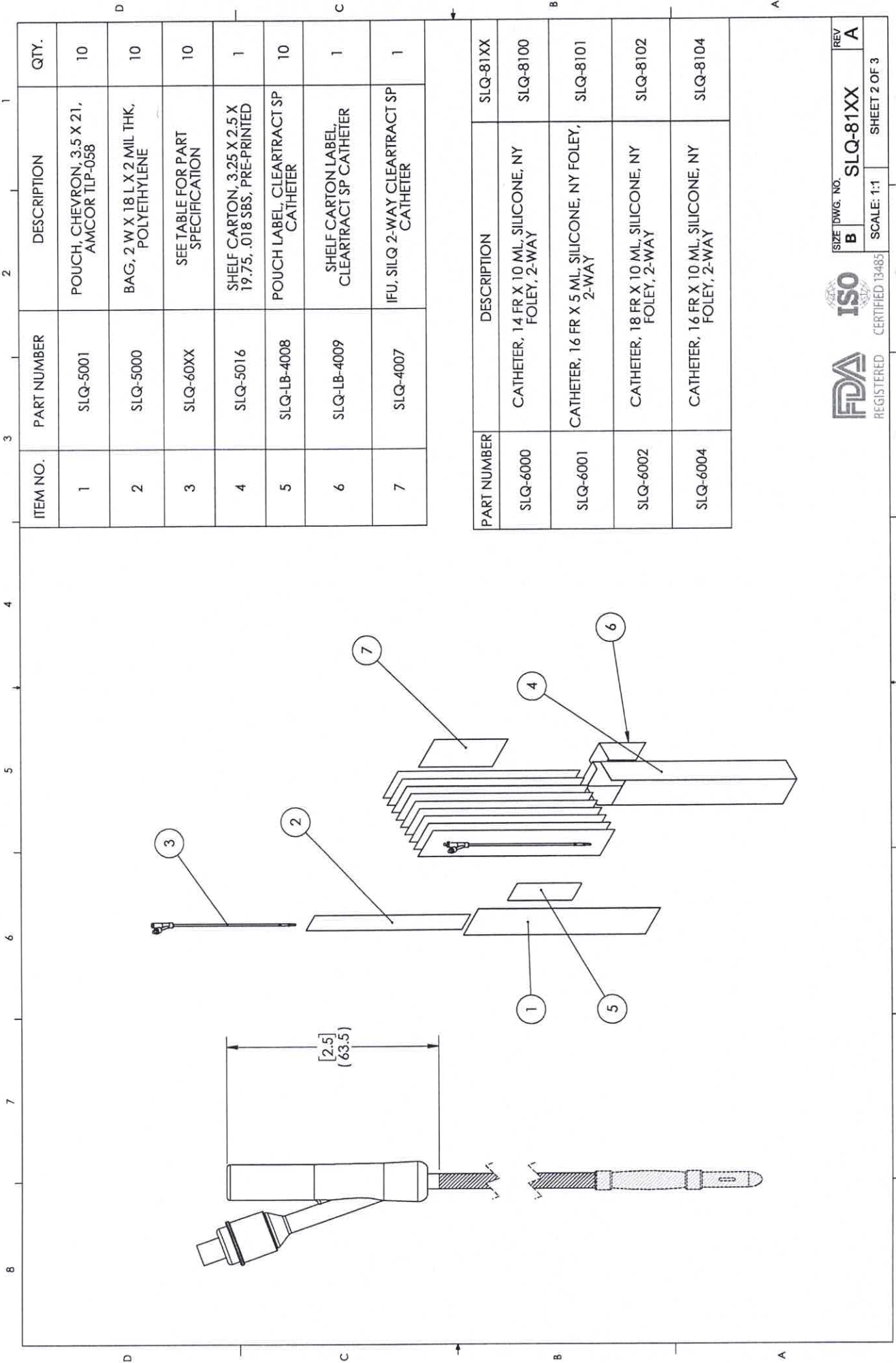

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
1	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	10
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	10
3	SLQ-60XX	SEE TABLE FOR PART SPECIFICATION	10
4	SLQ-5016	SHELF CARTON, 3.25 X 2.5 X 19.75, .018 SBS, PRE-PRINTED	1
5	SLQ-LB-4008	POUCH LABEL, CLEARTRACT SP CATHETER	10
6	SLQ-LB-4009	SHELF CARTON LABEL, CLEARTRACT SP CATHETER	1
7	SLQ-4007	IFU, SILQ 2-WAY CLEARTRACT SP CATHETER	1

PART NUMBER	DESCRIPTION	SLQ-81XX
SLQ-6000	CATHETER, 14 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8100
SLQ-6001	CATHETER, 16 FR X 5 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8101
SLQ-6002	CATHETER, 18 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8102
SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8104



SIZE DWG. NO. **B**



REV **A**

SLQ-81XX

SCALE: 1:1

SHEET 2 OF 3

ITEM NO.	PARTNO	DESCRIPTION	QTY.
8	SLQ-5012	BOX, 20 IN. L X 18 IN. W X 6 IN. D, 200 LB TEST	1
9	PWY-5000	TAPE, PACKAGING, 2 MIL, 2 IN X 110 YDS, CLEAR	AR
10	SLQ-LB-4010	SHIPPER BOX LABEL, CLEARTRACT SP CATHETER	2

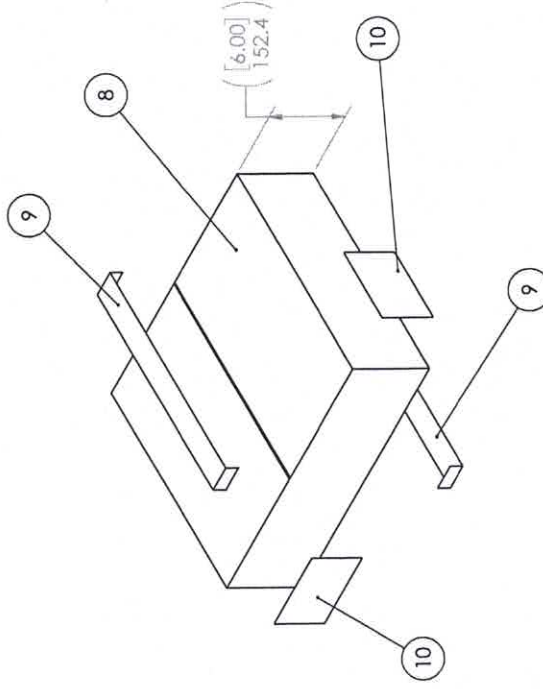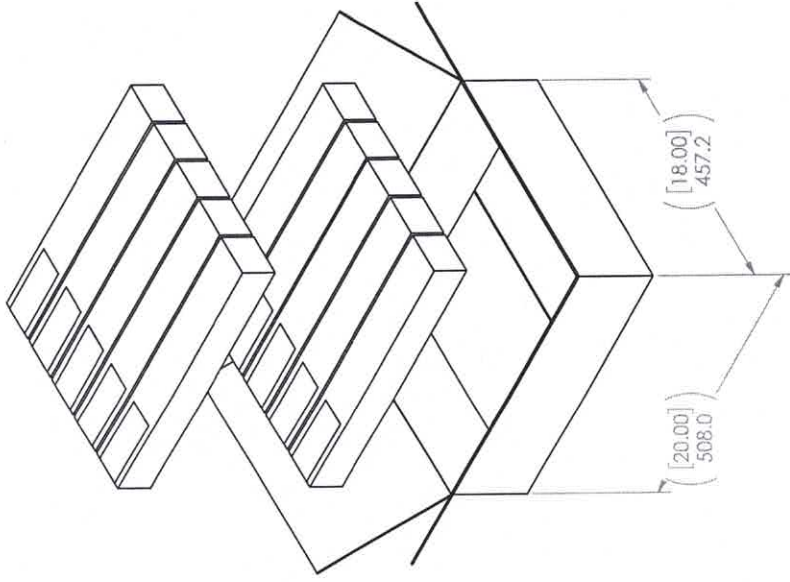



Printed by Sean Tamayo
Printed on 01/23/2025
Local time zone (GMT-08:00) Pacific Time (America/Los Angeles)
Page Arena > Workspace Items > Bill of Materials > Indented
All information contained in this document is proprietary and confidential.

SLQ-211810SPT

In Production

SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack

Revision »B - ECO-001669 Status: Effective, Status Date: 05/14/2024 03:24:42 PM, Shared

Contains 10 first-level Items, 50 line Items, 41 unique Items, 33 of which are shared.

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	△	🔍	🔄	📄	📁	📧
▶ 1 SLQ-8102 rev A	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack	📦 Subassembly	In Prod	📦	3	1	EA		
2 WI-0006 rev E	Work Order Process	📄 Work Instruction	In Prod	📄	2	1	each		
3 FM-0039 rev G	Form, Work Order	📄 Form	In Prod	📄	2	1	each		
4 WI-0014 rev B	Line Clearance Process	📄 Work Instruction	In Prod	📄	2	1	each		
5 MP-SLQ-0003 rev A	Manufacturing Process Flow, SILQ Foley Catheter	📄 Manufacturing Procedure	In Prod	📄	2	1	EA		
6 SLQ-0001 rev A	Sterilization Load Shipping Pallet Configuration	📦 Subassembly	In Prod	📦	2	1	EA		
▶ 7 QP-032 rev F	EO Sterilization Process Control	^{SOP} 📄 Quality System Procedure	In Prod	📄	7	1	each		
8 WI-0021 rev C	Work Instruction, Pyrogen (LAL) Evaluation Process	📄 Work Instruction	In Prod	📄	2	1	each		
▶ 9 QP-010 rev F	Final Product Inspection	^{SOP} 📄 Quality System Procedure	In Prod	📄	6	1	each		
10 FM-0029 rev B	Form, Certificate of Conformance	📄 Form	In Prod	📄	2	1	each		

NOTES: UNLESS OTHERWISE SPECIFIED

1.

GENERAL

1.

ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY SILQ PRIOR TO PRODUCTION.

2.

FINAL PRODUCT INSPECTION PER QP-010.

3.

PRODUCT EXPIRY DATE IS 36 MONTHS FROM THE DATE OF MANUFACTURE.
2.

STERILIZATION


1.

PRODUCT IS STERILIZED VIA ETHYLENE OXIDE (ETO) STERILIZATION PER QP-032 CYCLE PWI-7374.
3.

TESTING

1.

PYROGEN EVALUATION PER WI-0021 AT ≤ 0.5 EU/ML OR ≤ 20.0 EU/DEVICE

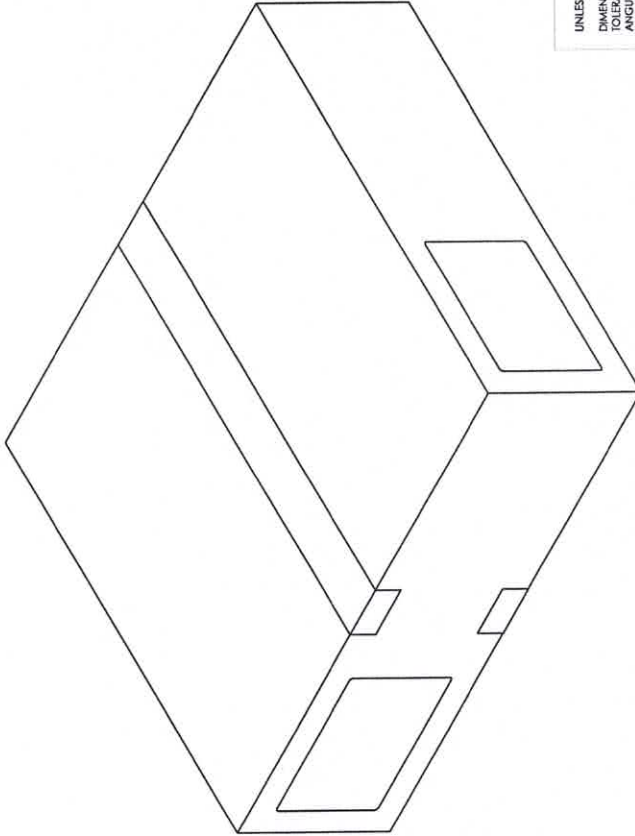

REV. A		1
DESCRIPTION		
INITIAL RELEASE		
CHANGE PRODUCT EXPIRY DATE FROM 12 TO 36 MONTHS, REVISED FORMAT TO MATCH T-018		



TITLE:

2-WAY CLEARTRACT CATHETER,
18 FR X 10 ML, STERILE, 10-PACK

SIZE DWG. NO. REV

B SLQ-211810SPT B

1-018 Rev. A SHEET 1 OF 1

DRAWN BY:

B. SANDOVAL

DRAWING NOT TO SCALE

1-018 Rev. A

SHEET 1 OF 1



Pathway Medtech, LLC.

Printed by Sean Tamayo
Printed on 01/23/2025
Local time zone (GMT-08:00) Pacific Time (America/Los Angeles)
Page Arena > Workspace Items > Bill of Materials > Indented
All information contained in this document is proprietary and confidential.

SLQ-8102

In Production

SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack

Revision »A - ECO-001669 Status: Effective, Status Date: 05/14/2024 03:24:42 PM, Shared

Contains 16 first-level Items, 31 line Items, 24 unique Items, 23 of which are shared.

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE
1 SLQ-6002 rev A	Catheter, 18 Fr x 10 mL, Silicone, NY Foley, 2-Way	Fabricated Part	In Prod
2 SLQ-3000 rev B	Solution, SILQ Polymer	Purchased Component	In Prod
SUBSTITUTES			
RANK#ITEM NUMBER		ITEM NAME	NOTES
1 HDX-3000 rev A		Solution, Hydrophilix Polymer	80 ML
3 SLQ-5000 rev A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	Packaging	In Prod
4 SLQ-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	Packaging	In Prod
SUBSTITUTES			
RANK#ITEM NUMBER		ITEM NAME	NOTES
1 HDX-5001 rev A		Pouch, Chevron, 3.5 X 21, Amcor TLP-058	10 EA
5 SLQ-LB-4008 rev A	Label, Pouch, SILQ ClearTract SPT Catheter	Labeling	In Prod

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE			
6	SLQ-5016 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS, Pre-Printed	Packaging	In Prod		3	1 EA
SUBSTITUTES							
		RANKQ ITEM NUMBER	ITEM NAME	NOTES			
		1	SLQ-5015 rev A	Shelf Carton, 3.25 1 EA x 2.5 x 19.75, .018 SBS			
		2	SLQ-5014 rev A	Shelf Carton, 3.25 1 EA x 2.5 x 19.75, .024 SBS			
7	SLQ-LB-4009 rev A	Label, Box, SILQ ClearTract SPT Catheter	Labeling	In Prod		2	2 EA
8	SLQ-4007 rev A	Artwork, IFU, SILQ 2-Way ClearTract SPT Catheter	Labeling	In Prod		2	1 EA
9	SLQ-5012 rev A	Box, 20 in. L x 18 in. W x 6 in. D, 200 LB Test	Packaging	In Prod		2	0.1 EA
10	PWY-5000 rev A	Tape, Packaging, 2 Mil, 2 in x 110 yds, Clear	Packaging	In Prod		2	5.2 IN
11	SLQ-LB-4010 rev A	Label, Shipper Box, SILQ ClearTract SPT Catheter	Labeling	In Prod		2	0.2 EA
12	MP-SLQ-0002 rev B	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	Manufacturing Procedure	In Prod		2	1 EA
13	FM-SLQ-0001 rev A	Form, Run Sheet for SILQ Automation Machine	Form	In Prod		2	1 EA
14	TM-HDX-0001 rev B	Test Method, UV Light Absorbance, Hydrophilix Catheters	Test Method	In Prod		2	1 each
15	TM-HDX-0002 rev A	Test Method, Rinse Test, Hydrophilix Catheter	Test Method	In Prod		2	1 each
16	MP-SLQ-0001 rev B	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	Manufacturing Procedure	In Prod		2	1 EA

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 1 of 7


Production Planning – Document the products identification information with the planned build quantity and any deviation(s) that may apply. Attach a copy of the device's specification/drawing/BOM and any applicable deviation(s) to the work order form prior to transfer to manufacturing.

Part No.	Rev	Lot No./ Serial No.	Expiration Date	Quantity	Deviation No.
SLQ-211810SPT	B	SLQ-01242025	2028-01-30	50/ packs of 10	N/A
Comments					
SLQ 2-way Cleartract SPT 18 Fr x 10ml Catheter, Sterile, 10-Pack					
Only u/s catheters were coated due to material shortages. ST 02/06/2025					
Name			Signature		Date
Sean Tamayo					01/24/2025

Lot No/Serial Number	SLQ-01242025
----------------------	--------------

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 2 of 7

Line Clearance - The individual performing the line clearance certifies that line clearance activities were conducted per WI-0014 prior to starting manufacturing activities. Quality management or designee verifies the line clearance was performed. *Note: The date line clearance occurred is considered the start date of manufacturing.*

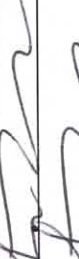Actions	Completed.	Performed by	Date
1. Clear Area of previous WO Raw Materials.	☒ Completed.		01/24/2025
2. Verify kitted Raw Materials (RM) for: a. RM meet BOM revision level requirements. b. RM have not Expired. c. RM have no physical damage	☒ Completed. ☒ Completed. ☒ Completed.		01/24/2025
3. Equipment is calibrated and properly functioning. <i>Note: Equipment identification and calibration due date will be documented in the work order.</i>	☒ Completed. ☐ N/A		01/24/2025
4. Proper documentation is present and available.	☒ Completed.		01/24/2025
5. Production Associates are trained and qualified to perform the process.	☒ Completed. ☐ N/A		01/24/2025
6. All work surfaces cleaned with IPA.	☒ Completed. ☐ N/A		01/24/2025
7. Labels printed per Label Generation and Reconciliation Form.	☐ Completed. <i>1</i> ☒ N/A		01/24/2025
8. Verify validated processes have been set to critical points within specification.	☒ Completed. ☐ N/A		01/24/2025
9. Verify Start-up activities are within specification	☒ Completed. ☐ N/A		01/24/2025
Conclusion: Line clearance performed	☒ Completed		01/24/2025
Quality Management			
Name		Signature	
Larry Dacera			
		Date	
		01/24/2025	

① labels to be printed at a later date. ST 01/24/2025

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 3 of 7

Material/ Procedure Traceability (Check appropriate Box):

- ☐ – Identification information for the procedures and materials used correspond to the attached product's BOM and are captured on the attached Production Pick List (ref. P/N FM-0059).
- ☒ - Identification information for the procedures and materials used correspond to the attached product's BOM and are captured in the table below (N/A or delete table if not applicable).

BOM Item No.	Type	Rev	Material Designation	Lot No./Serial No.	Expiration Date	Rev	Material Designation	Lot No./Serial No.	Expiration Date
SLQ-8102	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-01242025	2028-01-31		☐ Primary ☐ Substitute		
1	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	20241204	N/A		☐ Primary ☐ Substitute		
2	☐ Document ☒ Material	B	☒ Primary ☐ Substitute	2302-01-1	02/22/2024, REF NCR 24-027		☐ Primary ☐ Substitute		
3	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
4	☐ Document ☒ Material	A	☐ Primary ☒ Substitute	A199054	N/A		☐ Primary ☐ Substitute	N/A To 02/19/2025	
5	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-01242025	2028-01-31		☐ Primary ☐ Substitute		
6	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	000072 011223	N/A		☐ Primary ☐ Substitute		
7	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-01242025	2028-01-31		☐ Primary ☐ Substitute		
8	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
9	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
10	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
11	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
12	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
13	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		



Work Order Form

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Doc No.: FM-0039
 Rev.: G
 Page 4 of 7

BOM Item No.	Type	Rev	Material Designation	Lot No./Serial No.	Expiration Date	Rev	Material Designation	Lot No./Serial No.	Expiration Date
14	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
15	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
16	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
SLQ-211810SPT	☐ Document ☒ Material	B	☒ Primary ☐ Substitute	SLQ-01242025	2028-01-31		☐ Primary ☐ Substitute		
1	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
2	☒ Document ☐ Material	E	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
3	☒ Document ☐ Material	G	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
4	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
5	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
6	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
7	☒ Document ☐ Material	F	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
8	☒ Document ☐ Material	C	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
9	☒ Document ☐ Material	F	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
10	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 5 of 7

Training Traceability – All individuals indicated in the following table have been trained to the referenced document(s) which support the manufacturing, inspection and test activities of the product identified in the work order. Leadership representative has certified that the individual(s) have received proper training prior to conducting manufacturing activities.

Document No.	Rev	Print Name	Signature
MP-SLQ-0001	B	Saverdy Pen	Saverdy Pen
		James Romero	James Romero
		Aiden Emerick	Aiden Emerick
		N/A	
MP-SLQ-0002	B	Saverdy Pen	Saverdy Pen
		James Romero	James Romero
		Aiden Emerick	Aiden Emerick
		N/A	
Leadership Representative ①			
Name		Signature	Date
Sean Lammey		Sean Lammey	01/25/2025

① ST 02/06/2025

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 6 of 7

Equipment Traceability – Identify any equipment that is used in the manufacturing/processing of the product in the following table.

Equipment Description	Asset No.	Cal Due Date (If applicable)	Equipment Description	Asset No.	Cal Due Date (if applicable)
Accu Seal Impulse Sealer	ES-0099-01	08/09/2025	Catheter Automation Machine	ES-0091-01	N/A

Production Reporting – Document scrapped information, date manufacturing completed, and final quantity manufactured.

BOM Item No.	Defect No. 1		Defect No. 2		Defect No. 3		Defect No. 4		Defect No. 5		Comments
	Code	QTY	Code	QTY	Code	QTY	Code	QTY	Code	QTY	
SLQ-6002	3	58		10		28	N/A				Code 10: QTY 18 Consumed for client testing performed by Ethan Rao QTY 10 Sent out for testing at LGGs
Device Assembly							N/A				

Defect Codes: 1=Damage, 2=Dimensional, 3=Visual, 4=Voids/Shorts, 5=Discoloration, 6=Contamination, 7=Operator Error, 8=Occlusion, 9=Leak, 10 =Miscellaneous (Indicate in Comment Section)

Manufacturing Completion Date	02/06/2025		Quantity MFG	369/36.9 10 packs	
Name	Sean Tamayo		Signature		
			Date	02/06/2025	

① ST 02/06/2025


Lot No/Serial Number	SLQ-01242025
----------------------	--------------

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 7 of 7


Documentation Checklist – If applicable, the following documentation shall be included upon completion of the manufacturing build. Document packet shall be organized in the order indicated below. Ensure that good documentation practices are applied per QP-003 "Quality Records". Check box(es) ☐ that apply to each requirement. The individual who signs has confirmed that all requirements have been met.

1. Work Order (WO) Record	☒ Present and Completed	
2. Product's Specification/Drawing & BOM (Arena PLM)	☒ Present	
3. Deviations	☐ Present	☒ N/A
4. Product Labeling & Generation/ Reconciliation Record	☒ Present, Completed and Results Meet Specification	☐ N/A
5. Non-Conforming Reports (NCR)	☐ Present and Completed	☒ N/A
6. Inspection Records	☒ Present, Completed and Results Meet Specification	
7. Critical Process Parameter Record	☐ Present, Completed and Results Meet Specification	☐ N/A
Name		Date
Larry Dacoron		02/06/2025

Release Authorization for Further Processing (e.g., terminal sterilization) – The individual signing below is certifying that they independently reviewed the Device History Record (DHR) and confirmed that all required documentation is included, records are complete, good documentation practices were followed, and specifications have been met (if not NCR is included)

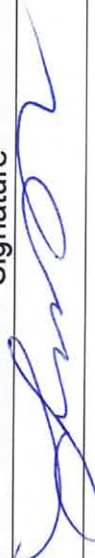Name	Signature	Date
Larry Dacoron		02/06/2025

Lot No/Serial Number

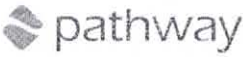	Label Generation & Reconciliation					Doc No.: FM-0033		
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>							Rev.: D
								Page 1 of 2

Label Information							
Part No.	SLQ-LB-4009	Rev.	A	Lot No.	SLQ-01242025 ☐ N/A	Label QTY	10 ☐ N/A
Description	Label, Box, SILQ ClearTract SPT Catheter	UDI No.	10850041415168 ☐ N/A		Expiration Date:	2028-01-31 ☐ N/A	

Label Quantity Calculation	
Required Label Quantity (A)	50
2% Overage ($A \times 0.02 = B$)	2
QC Retain Labels	2
Quantity of Labels Required ($A + B + 2 = C$)	54
Quantity of Labels Printed ($C = D$)	54
Production's Initial	ST
Date	01/31/2025

Label Verification	
What is the revision of the label per the BOM?	A
What is the revision on the printed label?	A
Is the printed information correct per specifications?	☒ Yes ☐ No
First and last label printed is affixed to back of form	☒ Yes ☐ No
Quantity of labels released to production (E)	52
QC's Initials	GS
Date	02/06/2025

Label Reconciliation	
Quantity of labels used by production (F)	37
Quantity scrapped by production (G)	0
Quantity of unused labels returned to QC and scrapped (H)	15
($E - F - G - H = 0$)	0
Are all labels accounted for?	☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

- Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
- Review the appropriate workstations and surrounding areas for the missing labels.

Investigational Results	
☒ N/A	

	Label Generation & Reconciliation	Doc No.: FM-0033
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>	Rev.: D
		Page 2 of 2
QC's Initials		Date
		02/04/2025

SIZE	18 Fr, 10 mL
QUANTITY	10

2-WAY 100% SILICONE
CLEARTRACT SPT CATHETER



Rx Only    **STERILE EO**   **MR**



(01) 10850041415168
(10) SLQ-01242025
(17) 280131

LOT SLQ-01242025
REF 211810SPT
2028-01-31

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

 SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

First

SIZE	18 Fr, 10 mL
QUANTITY	10

2-WAY 100% SILICONE
CLEARTRACT SPT CATHETER



Rx Only    **STERILE EO**   **MR**



(01) 10850041415168
(10) SLQ-01242025
(17) 280131

LOT SLQ-01242025
REF 211810SPT
2028-01-31

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

 SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

last

NOTES: UNLESS OTHERWISE SPECIFIED

- GENERAL
 - PRIOR TO PRINTING, VERIFY THE INFORMATION ON THE LABEL IS ACCURATE WHEN COMPARED AGAINST THE WORK ORDER.
 - PRIOR TO PRINTING, VERIFY THE LABEL HAS GONE THROUGH THE INSTALLATION QUALIFICATION (IQ) PROCESS FOR LABELS. REFER TO FM-0044.
 - LABEL CONTENT SHALL BE DEVELOPED PER WI-0009.

- MATERIAL
 - LABEL: STANDARD WHITE THERMAL TRANSFER
 - RIBBON: WAX / RESIN


- DIMENSIONS
 - () INDICATES LABEL REFERENCE DIMENSIONS.

- SAMPLING PLAN
 - IN ADDITION TO THE LABEL GENERATION AND RECONCILIATION FORM, FM-0033, INSPECT N=3 ADDITIONAL SAMPLES AT THE BEGINNING OF THE RUN.

- GS1 DATA MATRIX VALUES
 - (01): GLOBAL TRADE ITEM NUMBER (GTIN)
 - (10): BATCH OR LOT NUMBER
 - (17): EXPIRATION DATE

- VISUAL INSPECTION - LABEL TEXT
 - POST-PRINTING, ALL TEXT ON THE LABEL SHALL BE LEGIBLE.
 - ALL VARIABLE LABEL CONTENT (LOT NUMBER, EXPIRATION DATE, ETC.) SHALL MATCH THE WORK ORDER.
 - THE GS1 DATA MATRIX SHALL BE FREE OF MISSING PIXELS AND SHALL NOT HAVE SMUDGES OR STREAKS THROUGH IT. WHEN SCANNED WITH A BARCODE SCANNER, THE TEXT ADJACENT TO THE GS1 DATA MATRIX SHALL MATCH THE TEXT DISPLAYED BY THE SCANNER.
 - THE GS1 DATA MATRIX. THE GTIN SHALL BE A 14 DIGIT NUMBER.
 - (10): VERIFY THE BATCH OR LOT NUMBER FORMAT IS IN COMPLIANCE WITH QP-019.
 - (17): VERIFY THE EXPIRATION DATE IS IN THE FOLLOWING FORMAT: YYMMDD.
 - VERIFY THE CATHETER STYLE AND CATHETER SIZE ARE REPRESENTATIVE OF A CONFIGURATION FOUND IN THE SKU TABLE.
 - VERIFY THE GTIN IS FOUND IN THE SKU TABLE AND CORRESPONDS WITH THE CORRECT CATHETER SIZE AND STYLE.
 - VERIFY THE REFERENCE NUMBER FOR THE CATHETER CORRESPONDS TO THAT FOUND IN THE SKU TABLE.
 - VERIFY "QUANTITY" IS 10.

- VISUAL INSPECTION - LABEL SYMBOLS
 - POST-PRINTING, ALL SYMBOLS ON THE LABEL SHALL BE FREE OF SMUDGES / MISSING ARTWORK.
 - THE EXPIRATION DATE NEXT TO ITS CORRESPONDING SYMBOL SHALL BE IN THE FORMAT: YYYY-MM-DD.

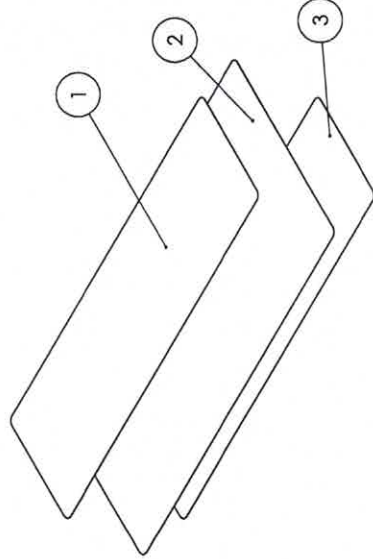

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

ITEM NO.	PARTNO	DESCRIPTION	QTY.
1	SLQ-AW-4009	LABEL, BOX, SILQ CLEARTRACT SPT CATHETER	1
2	PWY-5007	RIBBON, THERMAL TRANSFER, WAX/RESIN, 4.3W	1
3	SLQ-5004	LABEL, 2.3125 IN. X 7.50 IN.	1

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH+0.2 BEND ±2.0
NO PLACE DECIMAL ±0.5
ONE PLACE DECIMAL ±0.3
TWO PLACE DECIMAL ±0.1
THREE PLACE DECIMAL ±0.05

THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NT. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NT IS PROHIBITED.


pathway
New Product Introduction

TITLE:
LABEL, BOX, SILQ CLEARTRACT SPT CATHETER

SIZE DWG. NO. REV
B SLQ-LB-4009 A

SCALE: 1:1 SHEET 1 OF 2



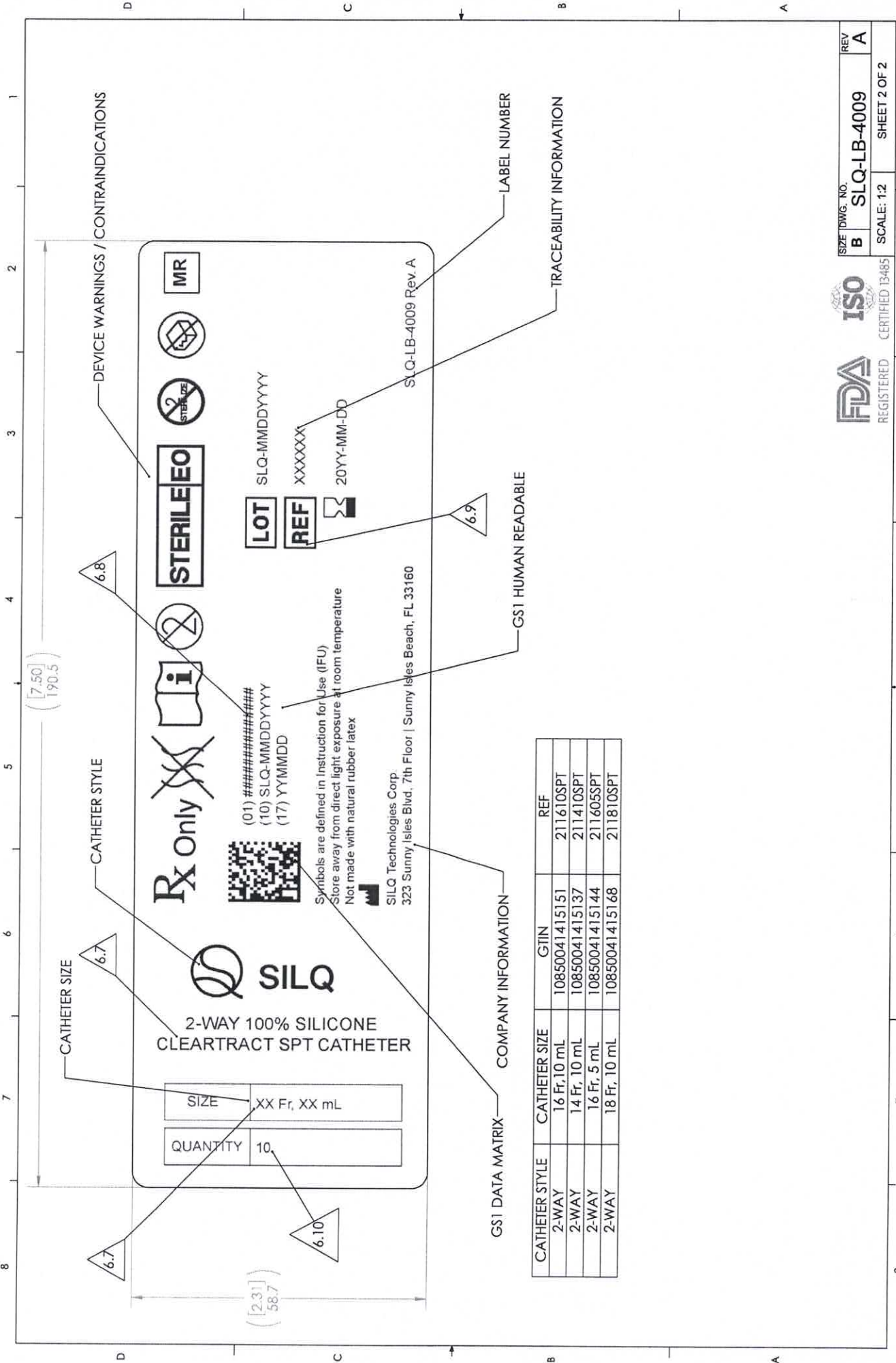

2-WAY 100% SILICONE
CLEARTRACT SPT CATHETER

SIZE	XX Fr, XX mL
QUANTITY	10

Rx Only

(01) #####
(10) SLQ-MMDDYY
(17) YYMMDD

STERILE

LOT SLQ-MMDDYY
REF XXXXX
20YY-MM-DD

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

CATHETER STYLE	CATHETER SIZE	GTIN	REF
2-WAY	16 Fr, 10 mL	10850041415151	211610SPT
2-WAY	14 Fr, 10 mL	10850041415137	211410SPT
2-WAY	16 Fr, 5 mL	10850041415144	211605SPT
2-WAY	18 Fr, 10 mL	10850041415168	211810SPT

	Label Generation & Reconciliation					Doc No.: FM-0033		
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.							Rev.: D
								Page 1 of 2

Label Information							
Part No.	SLQ-LB-4008	Rev.	A	Lot No.	SLQ-01242025 ☐ N/A	Label QTY	1 ☐ N/A
Description	Label, Pouch, SILQ ClearTract SPT Catheter	UDI No.	00850041415161 ☐ N/A		Expiration Date:	2028-01-31 ☐ N/A	

Label Quantity Calculation	
Required Label Quantity (A)	500
2% Overage (A×0.02=B)	10
QC Retain Labels	2
Quantity of Labels Required (A+B+2=C)	512
Quantity of Labels Printed (C=D)	512
Production's Initial	ST
Date	01/24/2025

Label Verification	
What is the revision of the label per the BOM?	A
What is the revision on the printed label?	A
Is the printed information correct per specifications?	☒ Yes ☐ No
First and last label printed is affixed to back of form	☒ Yes ☐ No
Quantity of labels released to production (E)	510
QC's Initials	AS
Date	02/06/2026

Label Reconciliation	
Quantity of labels used by production (F)	369
Quantity scrapped by production (G)	0
Quantity of unused labels returned to QC and scrapped (H)	141
(E-F-G-H=0)	0
Are all labels accounted for?	☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

1. Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
2. Review the appropriate workstations and surrounding areas for the missing labels.

Investigational Results
☒ N/A

	Label Generation & Reconciliation		Doc No.: FM-0033
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: D
			Page 2 of 2
QC's Initials			Date
			02/06/2025



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	18 Fr, 10 mL
QUANTITY	1

R_x Only



STERILE EO



MR

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

 SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

LOT SLQ-01242025

REF 211810SPT

 2028-01-31



(01) 00850041415161
(10) SLQ-01242025
(17) 280131

SLQ-LB-4008 Rev. A

First



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	18 Fr, 10 mL
QUANTITY	1

R_x Only



STERILE EO



MR

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

 SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

LOT SLQ-01242025

REF 211810SPT

 2028-01-31



(01) 00850041415161
(10) SLQ-01242025
(17) 280131

SLQ-LB-4008 Rev. A

Last

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
 1. PRIOR TO PRINTING, VERIFY THE INFORMATION ON THE LABEL IS ACCURATE WHEN COMPARED AGAINST THE WORK ORDER.
 2. PRIOR TO PRINTING, VERIFY THE LABEL HAS GONE THROUGH THE INSTALLATION QUALIFICATION (IQ) PROCESS FOR LABELS. REFER TO FM-0044.
 3. LABEL CONTENT SHALL BE DEVELOPED PER WI-0009.

2. MATERIAL
 1. LABEL: STANDARD WHITE THERMAL TRANSFER
 2. RIBBON: WAX / RESIN

3. DIMENSIONS
 1. () INDICATES LABEL REFERENCE DIMENSIONS.

4. SAMPLING PLAN
 1. IN ADDITION TO THE LABEL GENERATION AND RECONCILIATION FORM, FM-0033, INSPECT N=3 ADDITIONAL SAMPLES AT THE BEGINNING OF THE RUN.


5. GS1 DATA MATRIX VALUES
 1. (01): GLOBAL TRADE ITEM NUMBER (GTIN)
 2. (10): BATCH OR LOT NUMBER
 3. (17): EXPIRATION DATE

VISUAL INSPECTION - LABEL TEXT

1. POST-PRINTING, ALL TEXT ON THE LABEL SHALL BE LEGIBLE.
2. ALL VARIABLE LABEL CONTENT (LOT NUMBER, EXPIRATION DATE, ETC.) SHALL MATCH THE WORK ORDER.
3. THE GS1 DATA MATRIX SHALL BE FREE OF MISSING PIXELS AND SHALL NOT HAVE SMUDGES OR STREAKS THROUGH IT. WHEN SCANNED WITH A BARCODE SCANNER, THE TEXT ADJACENT TO THE GS1 DATA MATRIX SHALL MATCH THE TEXT DISPLAYED BY THE SCANNER.
4. (01): VERIFY THE GTIN IS DETAILED IN A HUMAN READABLE FORMAT ADJACENT TO THE GS1 DATA MATRIX. THE GTIN SHALL BE A 14 DIGIT NUMBER.
5. (10): VERIFY THE BATCH OR LOT NUMBER FORMAT IS IN COMPLIANCE WITH QP-019.
6. (17): VERIFY THE EXPIRATION DATE IS IN THE FOLLOWING FORMAT: YYMMDD. VERIFY THE CATHETER STYLE AND CATHETER SIZE ARE REPRESENTATIVE OF A CONFIGURATION FOUND IN THE SKU TABLE.
7. VERIFY THE GTIN IS FOUND IN THE SKU TABLE AND CORRESPONDS WITH THE CORRECT CATHETER SIZE AND STYLE.
8. VERIFY THE REFERENCE NUMBER FOR THE CATHETER CORRESPONDS TO THAT FOUND IN THE SKU TABLE.
9. VERIFY "QUANTITY" IS 1.

VISUAL INSPECTION - LABEL SYMBOLS

1. POST-PRINTING, ALL SYMBOLS ON THE LABEL SHALL BE FREE OF SMUDGES / MISSING ARTWORK.
2. THE EXPIRATION DATE NEXT TO ITS CORRESPONDING SYMBOL SHALL BE IN THE FORMAT: YYYY-MM-DD.

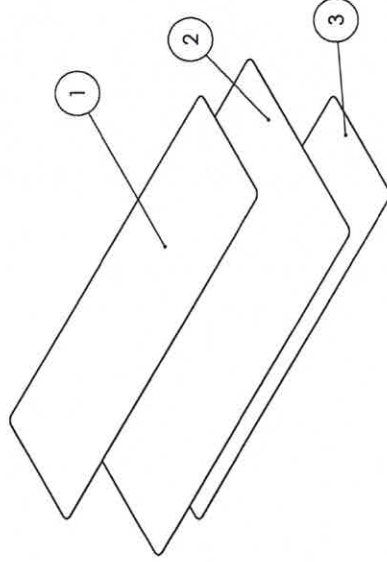

ITEM NO.	PARTNO	DESCRIPTION	QTY.
1	SLQ-AW-4008	LABEL, POUCH, SILQ CLEARTRACT SPT CATHETER	1
2	PWY-5007	RIBBON, THERMAL TRANSFER, WAX/RESIN, 4.3W	1
3	SLQ-5004	LABEL, 2.3125 IN. X 7.50 IN.	1

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

pathway
New Product Introduction

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH±0.2 BEND ±2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ± 0.1
THREE PLACE DECIMAL ± 0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF PATHWAY NP. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF PATHWAY NP IS PROHIBITED.

TITLE:
LABEL, POUCH, SILQ CLEARTRACT SPT CATHETER


SIZE DWG. NO. **B SLQ-LB-4008**

SCALE: 1:1

DO NOT SCALE DRAWING



REGISTERED CERTIFIED 13485

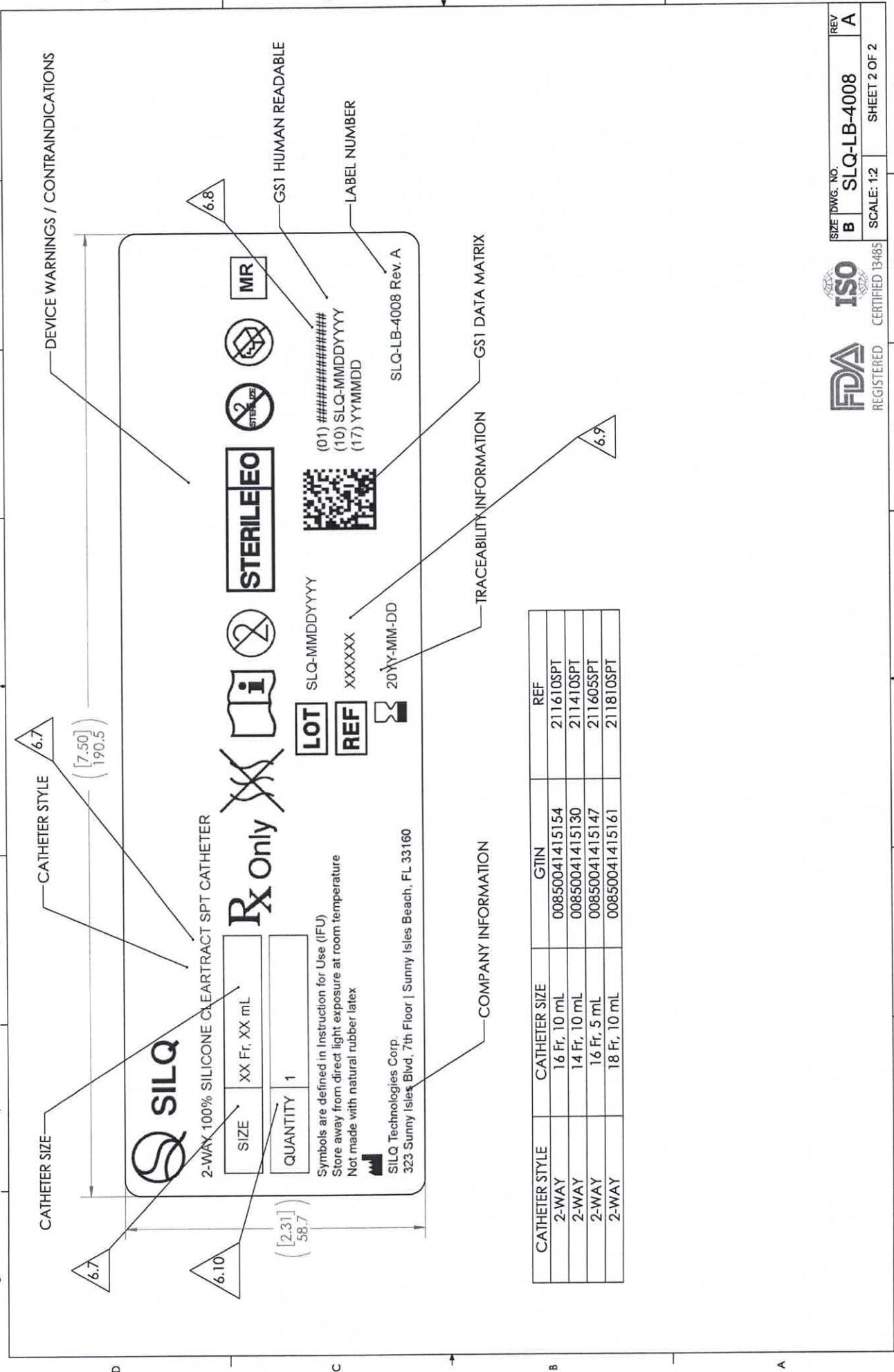

CATHETER STYLE	CATHETER SIZE	GTIN	REF
2-WAY	16 Fr, 10 mL	00850041415154	211610SPT
2-WAY	14 Fr, 10 mL	00850041415130	211410SPT
2-WAY	16 Fr, 5 mL	00850041415147	211605SPT
2-WAY	18 Fr, 10 mL	00850041415161	211810SPT

	Label Generation & Reconciliation				Doc No.: FM-0033	
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.					
	Rev.: D					
						Page 1 of 2

Label Information						
Part No.	SLQ-LB-4010	Rev.	A	Lot No.	SLQ-01242025 ☐ N/A	Label QTY
						100 ☐ N/A
Description	Label, Shipper Box, SILQ ClearTract SPT Catheter	UDI No.	10850041415168 ☐ N/A		Expiration Date:	2028-01-31 ☐ N/A

Label Quantity Calculation	
Required Label Quantity (A)	10
2% Overage (A×0.02=B)	1
QC Retain Labels	2
Quantity of Labels Required (A+B+2=C)	13
Quantity of Labels Printed (C=D)	13
Production's Initial	ST
Date	01/31/2025

Label Verification	
What is the revision of the label per the BOM?	A
What is the revision on the printed label?	A
Is the printed information correct per specifications?	☒ Yes ☐ No
First and last label printed is affixed to back of form	☒ Yes ☐ No
Quantity of labels released to production (E)	11
QC's Initials	BS
Date	02/04/2025

Label Reconciliation	
Quantity of labels used by production (F)	4
Quantity scrapped by production (G)	0
Quantity of unused labels returned to QC and scrapped (H)	7
(E-F-G-H=0)	0
Are all labels accounted for?	☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

1. Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
2. Review the appropriate workstations and surrounding areas for the missing labels.

Investigational Results
☒ N/A

	Label Generation & Reconciliation	Doc No.: FM-0033
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>	Rev.: D
		Page 2 of 2
QC's Initials		Date 



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF

211810SPT

LOT

SLQ-01242025



2028-01-31

QTY

100

 SILQ Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

First



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF

211810SPT

LOT

SLQ-01242025



2028-01-31

QTY

100

 SILQ Technologies Corp.


323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

Last

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
 1. PRIOR TO PRINTING, VERIFY THE INFORMATION ON THE LABEL IS ACCURATE WHEN COMPARED AGAINST THE WORK ORDER.
 2. PRIOR TO PRINTING, VERIFY THE LABEL HAS GONE THROUGH THE INSTALLATION QUALIFICATION (IQ) PROCESS FOR LABELS. REFER TO FM-0044.
 3. LABEL CONTENT SHALL BE DEVELOPED PER WI-0009.
2. MATERIAL
 1. LABEL: STANDARD WHITE THERMAL TRANSFER
 2. RIBBON: WAX / RESIN
3. DIMENSIONS
 1. () INDICATES LABEL REFERENCE DIMENSIONS.
4. SAMPLING PLAN
 1. IN ADDITION TO THE LABEL GENERATION AND RECONCILIATION FORM, FM-0033, INSPECT N=3 ADDITIONAL SAMPLES AT THE BEGINNING OF THE RUN.
5. VISUAL INSPECTION - LABEL TEXT
 1. POST-PRINTING, ALL TEXT ON THE LABEL SHALL BE LEGIBLE.
 2. ALL VARIABLE LABEL CONTENT (LOT NUMBER, EXPIRATION DATE, ETC.) SHALL MATCH THE WORK ORDER.
 3. VERIFY THE CATHETER STYLE IS REPRESENTATIVE OF A CONFIGURATION FOUND IN THE SKU TABLE.
 4. VERIFY THE REFERENCE NUMBER FOR THE CATHETER CORRESPONDS TO THAT FOUND IN THE SKU TABLE.
 5. VERIFY "QUANTITY" IS 100.
6. VISUAL INSPECTION - LABEL SYMBOLS
 1. POST-PRINTING, ALL SYMBOLS ON THE LABEL SHALL BE FREE OF SMUDGES / MISSING ARTWORK.
 2. THE EXPIRATION DATE NEXT TO ITS CORRESPONDING SYMBOL SHALL BE IN THE FORMAT: YYYY-MM-DD.

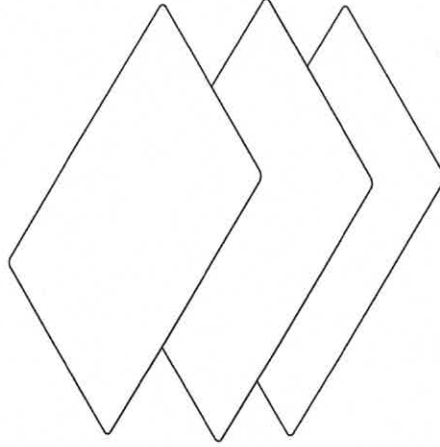

ITEM NO.	PARTNO	DESCRIPTION	QTY.
1	SLQ-AW-4010	ARTWORK, LABEL, SHIPPER BOX, SILQ CLEARTRACT SPT CATHETER	1
2	PWY-5004	LABEL, 4 X 6	1
3	PWY-5007	RIBBON, THERMAL TRANSFER, WAX/RESIN, 4.3W	1

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
	A	INITIAL RELEASE	

UNLESS OTHERWISE SPECIFIED:
 DIMENSIONS ARE IN MILLIMETERS
 TOLERANCES:
 ANGULAR: MACH±0.2 BEND ±2.0
 NO PLACE DECIMAL ±0.3
 ONE PLACE DECIMAL ±0.3
 TWO PLACE DECIMAL ±0.15
 THREE PLACE DECIMAL ±0.05

TITLE:

LABEL, SHIPPER BOX, SILQ
 CLEARTRACT SPT CATHETER

SIZE DWG. NO. REV

B SLQ-LB-4010 A

SCALE 1:2

DO NOT SCALE DRAWING

CERTIFIED 13485



pathway
 New Product Introduction

PROPERTY AND CONFIDENTIAL
 THE INFORMATION CONTAINED IN THIS
 DRAWING IS THE PROPERTY OF PATHWAY
 MEDICAL, INC. ANY REPRODUCTION IN
 PART OR AS A WHOLE WITHOUT THE
 WRITTEN PERMISSION OF PATHWAY MEDICAL
 IS PROHIBITED.

REGISTERED

4

5

6

7

8

9

10

11

12

13

14

15

16

17

18

19

20

21

TRACEABILITY INFORMATION

CATHETER STYLE



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF

XXXXXX

LOT

SLQ-MMDDYYYY



20YY-MM-DD

QTY

100



SLQ Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

LABEL NUMBER

([6.00])
152.4

CATHETER STYLE	REF
2-WAY	211410SPT
2-WAY	211605SPT
2-WAY	211810SPT
2-WAY	211610SPT



SIZE DWG. NO.

B SLQ-LB-4010

REV

A


SCALE: 1:1

SHEET 2 OF 2

Date Inspected	Part No.	Rev.	Description						
2/5/2025	SLQ-8102	A	SILQ 2-Way Cleartract SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack						
Lot No.		Lot QTY	ES-0059-01 Emplex Model MPS6300-M Continuous Band Sealer - Validated Range: Temperature = 142 °C - 144 °C, Speed = 175 in/min., Pressure = 89 psi ES-0099-01 ACCU-SEAL Model 5300-15-B, Impulse Sealer - Validated Range: Temperature = 270 °F - 280 °F, Dwell = 4 sec., Pressure = 55 psi						
Specification		TM-0003							
Revision		C							
PRODUCT PREP FOR TEST									
* Seal one sample pouch P/N SLQ-5001 at the start, middle and end of each shift.									
* Cut two 1.00 inch wide x 2.00 inch long strips from the manufacturer sealed end of the pouch.									
* Execute peel strength test per TM-0003 and record results.									
Specification	≥1.0	≥1.0	142-144 °C	175 in/min	89 psi	270 - 280 °F	4.0	55.000	
Lower Spec Limit (LSL)	1.0	1.0	142.0	175.0	89	270.0	4.0	55	
Upper Spec Limit (USL)	N/A	N/A	144.0	175.0	89	280.0	4.0	55	
UOM	lb	lb	°C	in/min	psi	°F	seconds	psi	
Measuring Equipment Type	Force Gauge	Force Gauge	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	
Asset No.	ES-0028-01	ES-0028-01	N/A	N/A	N/A	ES-0099-01	ES-0099-01	ES-0099-01	
Cal Due Date			N/A	N/A	N/A	8/9/2024	8/9/2024	8/9/2024	
Sample ID	Top Strip	Bottom Strip	N/A	N/A	N/A	1	1	1	
Sample No.	Variable/ Attribute Data								
1: Start (Top Strip)	6.270	N/A	N/A	N/A	N/A	270.0	4.0	55.0	
2: Start (Bottom Strip)	N/A	6.598	N/A	N/A	N/A	N/A	N/A	N/A	
3: Middle (Top Strip)	7.796	N/A	N/A	N/A	N/A	N/A	N/A	N/A	
4: Middle (Bottom Strip)	N/A	6.494	N/A	N/A	N/A	N/A	N/A	N/A	
5: End (Top Strip)	6.930	N/A	N/A	N/A	N/A	N/A	N/A	N/A	
6: End (Bottom Strip)	N/A	6.434	N/A	N/A	N/A	N/A	N/A	N/A	
Result Summary									
Pass QTY	3	3	N/A	N/A	N/A	1	1	1	
Fail QTY	0	0	N/A	N/A	N/A	0	0	0	
Result (P/F)	P	P	N/A	N/A	N/A	P	P	P	
Statistical Summary									
Min	6.270	6.434	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	
Max	7.796	6.598	#REF!	#REF!	0.00	#REF!	#REF!	#REF!	
Mean	6.999	6.509	#REF!	#REF!	0.00	#REF!	#REF!	#REF!	
StDev	0.765	0.083	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	
Comments									
N/A						NCMR No: N/A			
Disposition: ACCEPT						Date: 02/05/2025			
Inspected By: [Signature]									



Date Inspected		Part No.		Rev.	Description				
02/06/2025		SLQ-8102		A	SILQ 2-Way Cleartract SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack				
Lot No.		Lot QTY		Approved Supplier(s)					
SLQ-01242025		500/50 packs of 10		Pathway					
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	A	A	A	A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	N/A	0.300	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
	UOM	EA	EA	AU	EA	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Visual	Spectrophotometer	Visual	Visual	Visual	Visual	Visual
	Asset No.	N/A	N/A	See Comments	N/A	N/A	N/A	N/A	N/A
	Cal Due Date	N/A	N/A	See Comments	N/A	N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size	9	9	9	9	9	9	9	5
Sample No.		Variable/ Attribute Data							
1	N/A	N/A	0.423	N/A	N/A	N/A	N/A	N/A	N/A
2	N/A	N/A	0.420	N/A	N/A	N/A	N/A	N/A	N/A
3	N/A	N/A	0.464	N/A	N/A	N/A	N/A	N/A	N/A
4	N/A	N/A	0.499	N/A	N/A	N/A	N/A	N/A	N/A
5	N/A	N/A	0.408	N/A	N/A	N/A	N/A	N/A	N/A
6	N/A	N/A	0.451	N/A	N/A	N/A	N/A	N/A	N/A
7	N/A	N/A	0.453	N/A	N/A	N/A	N/A	N/A	N/A

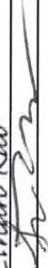8		N/A	N/A	0.418	0.000	N/A	N/A	N/A	N/A	N/A	N/A
9		N/A	N/A	0.475	0.004	N/A	N/A	N/A	N/A	N/A	N/A
Result Summary											
Pass QTY		N/A	N/A	9	9	N/A	N/A	N/A	N/A	N/A	N/A
Fail QTY		N/A	N/A	0	0	N/A	N/A	N/A	N/A	N/A	N/A
Result (P/F)		N/A	N/A	P	P	N/A	N/A	N/A	N/A	N/A	N/A
Statistical Summary											
Min		N/A	N/A	0.408	0.000	N/A	N/A	N/A	N/A	N/A	N/A
Max		N/A	N/A	0.499	0.007	N/A	N/A	N/A	N/A	N/A	N/A
Mean		N/A	N/A	0.450	0.003	N/A	N/A	N/A	N/A	N/A	N/A
StDev		N/A	N/A	0.030	0.002	N/A	N/A	N/A	N/A	N/A	N/A
Comments											
Rinse test and UV Absorbance Test were performed at SILQ. The rest of the visual inspections were performed at Pathway.											
Disposition:	ACCEPT					NCMR No:	N/A				
Inspected By:	Ethan Rao					Date:	2/6/2025				
Verified By:						Date:	02/06/2025				

Date Inspected		Part No.		Rev.	Description				
01/24/2025-02/06/2025		SLQ-8102		A	SILQ 2-Way Cleartract SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack				
Lot No.		Lot QTY		Approved Supplier(s)					
SLQ-01242025		500/50 packs of 10		Pathway					
In-Process Inspection	Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box. IP-0002
	Revision	A	A	A	A	A	A	A	A
	Lower Spec Limit (LSL)	N/A	N/A	0.300	N/A	N/A	N/A	N/A	N/A
	Upper Spec Limit (USL)	N/A	N/A	N/A	0.020	N/A	N/A	N/A	N/A
	UOM	EA	EA	AU	AU	EA	EA	EA	EA
	Measuring Equipment Type	Visual	Visual	Spectrophotometer	Spectrophotometer	Visual	Visual	Visual	Visual
	Asset No.	N/A	N/A	See Comments	See Comments	N/A	N/A	N/A	N/A
	Cal Due Date	N/A	N/A	See Comments	See Comments	N/A	N/A	N/A	N/A
	AQL	6.5	6.5	6.5	6.5	6.5	6.5	6.5	6.5
	Sample Size	9	9	9	9	9	9	9	4
Sample No.		Variable/ Attribute Data							
1	P	P	N/A	N/A	P	P	P	P	P
2	P	P	N/A	N/A	P	P	P	P	P
3	P	P	N/A	N/A	P	P	P	P	P
4	P	P	N/A	N/A	P	P	P	P	P
5	P	P	N/A	N/A	P	P	P	P	N/A
6	P	P	N/A	N/A	P	P	P	P	N/A
7	P	P	N/A	N/A	P	P	P	P	N/A
8	P	P	N/A	N/A	P	P	P	P	N/A
9	P	P	N/A	N/A	P	P	P	P	N/A

Result Summary

Pass QTY	9	9	N/A	N/A	9	9	9	9	4
Fail QTY	0	0	N/A	N/A	0	0	0	0	0
Result (P/F)	P	P	N/A	N/A	P	P	P	P	P
Statistical Summary									
Min	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	0.000	0.000	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	#DIV/0!	#DIV/0!	N/A	N/A	N/A	N/A	N/A

Rinse test and UV Absorbance Test were performed at SILQ. The rest of the visual inspections were performed at Pathway.

Disposition:	ACCEPT	NCMR No:	N/A
Inspected By:	<i>Richie Willis / James Korman / S. Lopez</i>	Date:	<i>11/6/23</i>
Verified By:	<i>[Signature]</i>	Date:	<i>11/6/23</i>

	Form, Run Sheet for SILQ Automation Machine		Doc No.: FM-SQL-0001
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: A
			Page 1 of 1

Part Number	Revision	Description	Machine Asset No.
SLQ-8102	A	SILQ 2-Way ClearTract SPT 18Fr x 10ml Catheter, Non-Sterile, 10-Pack	ES-0091-01

Equipment Setup / Startup		Adjustments during Production		
Setting [Units]	Value	Value	Value	Value
Setup / Adjustment Date	01/24/2025	01/24/2025		
Setup / Adjustment Time	10:35	11:15		
Fill Time [seconds]	9 seconds	8.5 seconds		
Confirm Cylinder Actuation All 4 pairs [P/F]	P			
UV Cure Time [seconds]	120 seconds	N/A ST 02/06/2025		
Drain Time [seconds]	30 seconds			
Rinse Time [seconds]	6 seconds	5.5		
Final Drain [seconds]	30 seconds			
Number of Drain / Rinse Cycles [Count]	3	N/A ST 02/06/2025		
Confirm Cylinder Actuation All 4 pairs [P/F]	P			
Size Setting	18 Fr			
Performed By Name	James Romero	James Romero		
Performed By Signature	James Romero	James Romero		
Verified By Name	Giovanna Siqueiros	Giovanna Siqueiros		
Verified By Signature	G.S.	G.S.		

Engineering / Lead Supervisor

QA

Date Inspected	Part No.		Rev.	Description							
2/19/25	SLQ-211810SPT		B	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack							
Lot No.		Lot QTY									
SLQ-01242025		369/36.9 packs of 10									
Identification of the product is traceable between the DHR, sterilization records, and laboratory report	DHR is completed and meets good documentation practices per QP-003	Sterilization processing record meets parameters of validated process	Sterility laboratory report meets acceptance criteria indicated in QP-010	Endotoxin laboratory report (LAL) meets acceptance criteria per WI-0021	No Damage to Outer Packing on Shipment back from Sterilizer	Verify Quantity Received back from Sterilizer matches Packing List					
Specification	QP-010 Revision	F	F	F	F	F	F				
Lower Spec Limit (LSL)	N/A	N/A	N/A	N/A	20.0	N/A	N/A				
Upper Spec Limit (USL)	N/A	N/A	N/A	N/A	N/A	N/A	N/A				
UOM	N/A	N/A	N/A	N/A	EU	N/A	N/A				
Measuring Equipment Type	Visual	Visual	Visual	Visual	Visual	Visual	Visual				
Asset No.	N/A	N/A	N/A	N/A	N/A	N/A	N/A				
Cal Due Date	N/A	N/A	N/A	N/A	N/A	N/A	N/A				
AQL	See Below	See Below	See Below	See Below	See Below	See Below	See Below				
DHR											
Sterilization Processing Record	DHR	Sterilization Processing Record	Sterility Laboratory Report	Endotoxin (LAL) Laboratory Report	Shipment	Shipment					
Laboratory Report											
Sample Size											
Sample No.	Variable/ Attribute Data										
1	P	P	P	P	P	P	P				
Result Summary											
Pass QTY	1	1	1	1	1	1	1				
Fail QTY	0	0	0	0	0	0	0				
Result (P/F)	P	P	P	P	P	P	P				
Statistical Summary											
Min	N/A	N/A	N/A	N/A	N/A	N/A	N/A				
Max	N/A	N/A	N/A	N/A	N/A	N/A	N/A				
Mean	N/A	N/A	N/A	N/A	N/A	N/A	N/A				
StDev	N/A	N/A	N/A	N/A	N/A	N/A	N/A				
Comments											
N/A ZO 2/19/25											
Disposition:	Accept						NCR No:	N/A			
Inspected By:	Lary Dacron						Date:	02/19/25			
Verified By:	Gordon Jacques						Date:	02/19/25			

Certificate of Analysis



Customer: Pathway Medtech LLC
Address: 8779 Cottonwood Avenue, Ste 105
Santee
California
USA
92071
Purchase Order No: P2029
Product: PWI 1162-16739 / Product Lot #: 250808 (BI Lot: 0394)
Description: NA004
Standard Method: USP
Internal Work Instruction: Q08-WI-001720 (LAB-003) Sterility Test Procedure
Customer Protocol: NA

Job No: BI-USA-BRO-202500512
Report Number: 146068
Date Issued: 19-Feb-2025
Page 1 of 1

Date Sample(s) Received: 12-Feb-2025
Date On Test: 12-Feb-2025

Analytical Results

Label ID / Designation: NA

Lot / Sterilization Batch: NA

Microbiology

Result

Biological Indicators

BI Type	Strip
BI Lot Number	0394
BI Expiration Date	08-08-2025
PCD Type	NA
PCD Lot Number	NA
Number of samples	7
Check 7	Negative
Positive Control Check 7	Positive
Number of positives	0
Negative Controls	Negative
Comment	NA

The results presented relate only to the items tested above as supplied by the Customer and received by the laboratory.
This test report should not be reproduced except in full, without prior written permission.

Released by:

Carrie Stinson

Quality Technician

Authorisation Date: 19-February-2025

Testing Facility
USA - Brooklyn Park, MN
STERIS
9303 West Broadway Avenue

Brooklyn Park
Minnesota
55445

STERIS

Ethylene Oxide Customer Specification

Customer: Pathway Innovation LLC Cycle/Specification #: PWI-7374 Rev: 1

Customer Contact: Aaron Rogers Phone Number: (619)415-0103

After-Hours Customer Contact: Aaron Rogers Phone Number: (619)415-0103

Approved Chambers: 6 Applicable Product/Item: PSM & HDX Foley Catheter

Product Classification (Check all that apply):

- ☒ Medical Device ☐ Pharmaceutical ☐ Active Pharmaceutical Ingredient ☐ Labware ☐ Botanicals ☐ Food / Spice
☐ Food Packaging ☐ Cosmetics ☐ Animal Food / Pet Treats ☐ Other (define) _____
☐ Check here if any of the above listed codes is an implantable medical device.
☐ Check here if any of the above listed codes include a battery or similar energy source.

1.0 General (any special comments): **Full/Routine Cycle**

2.0 Handling of Damaged Product and Count Discrepancies

Damage:

- ☒ Process all damaged cases (customer will be notified when damage occurs)
☐ Do not process damaged cases (damage to be returned to customer non-processed, customer will be notified when damage occurs)
☐ Disposition each damaged case (damages will be placed on hold until a customer disposition is obtained)
Note: Processing and/or shipping delays may result)
☐ Other (Specify): _____

Count Discrepancies:

- ☐ Continue processing at STERIS count (customer will be notified of discrepancy)
☒ Disposition discrepancy prior to processing (when discrepancy is noted upon receipt, product will be placed on hold until customer disposition is obtained)
☐ Disposition discrepancy prior to shipping (when discrepancy is noted during processing, product will be placed on hold until customer disposition is obtained)
☐ Other (Specify): _____

3.0 Sample Placement

- Samples Placed by: ☐ Customer ☒ STERIS
- Number and type of Samples: 7 Mail BIs, NA004 Spore Strip in a letter envelope.
- Sample Placement Instructions: See Placement Diagram

Front	Left	Right	Back
XX	XX	XX	X

Ethylene Oxide Customer Specification

Customer: Pathway Innovation LLC Cycle/Specification #: PWI-7374 Rev: 1

4.0 Product Staging Requirements

- First Pallet into Chamber: 1
- Load Temperature requirements: N/A
- Pallets per Load: Minimum 1 Maximum: 2
- Other: N/A
- Other: N/A

5.0 EO Processing

Preconditioning Specification (if applicable)			
Process variables	Units	Set point	Tolerance
Temperature	°F	115	±10
Relative Humidity	%RH	50	+30/-10
Dwell Time	Hours	18	16 - 32
Transfer Time to Chamber	Minutes	ASAP	≤60

Preconditioning Special Instructions

Extend the preconditioning dwell time to compensate for any time that the room is below the minimum specified limit for temperature and/ or relative humidity.

Handling of PCR Nonconformities

Temperature/RH low;

- ☒ Add time to meet minimum specified time in specification
☐ Other: _____

Temperature/RH high;

- ☐ Continue process and evaluate after processing
☐ Stop process and transfer product to warehouse
☒ Other See Special Instructions

Ethylene Oxide Customer Specification

Customer: Pathway Innovation LLC Cycle/Specification #: PWI-7374 Rev: 1

Cycle Parameters				
Process Variables	Unit	Set point	Tolerance	Rate
EO Concentration (Calculated)	mg/L	623.5	590.9 – 656.7	
Humidity (Calculated)	%	25.3	11.0 – 43.5	N/A
Chamber Temperature				
Cycle Start	°F	125	±10	N/A
Exposure	°F	125	±5	N/A
Initial Evacuation A				
Evacuation Pressure	InHgA	2.0	±0.5	≤2.5/min
Stabilization				
Hold Pressure	InHgA	2.0	<0.3	N/A
Dwell	Minutes	5	>1	N/A
Leak Test A				
Hold Pressure	InHgA	2.0	<0.3	N/A
Dwell	Minutes	5	±1	N/A
Nitrogen Dilution				
Inject Pressure	InHgA	28.0	±0.5	≤2.5/min
Evac Pressure	InHgA	2.0	±0.5	≤2.5/min
Repetitions	#	1	N/A	N/A
Humidity Injection				
Inject Steam	InHgA	3.0	±0.5	N/A
Humidity Dwell				
Dwell Pressure	InHgA	3.0	±0.5	N/A
Dwell Time	Minutes	60	55 - 75	N/A
Sterilant Injection				
Inject Pressure	InHgA	14.3	±0.5	≤2.5/min
Exposure				
Pressure	InHgA	14.3	±0.5	N/A
Temperature	°F	125	±5	N/A
Make up gas	EO	N/A	N/A	N/A
Dwell Time	Minutes	240	235 – 255	N/A
Sterilant Removal				
Evacuation Pressure	InHgA	2.0	±0.5	≤2.5/min
Nitrogen Wash				
Inject Pressure	InHgA	28.0	±0.5	≤2.5/min
Evacuation Pressure	InHgA	2.0	±0.5	≤2.5/min
Repetitions	#	3	N/A	N/A
Final Air Brake				
End Pressure	InHgA	27.9	Atmosphere	≤2.5/min

Add rows or attach additional sheet as necessary

Ethylene Oxide Customer Specification

Customer: Pathway Innovation LLC Cycle/Specification #: PWI-7374 Rev: 1**Cycle Special Instructions**

N/A

Heated Aeration (if applicable)

Process Variables	Units	Set point	Tolerance
Temperature	°F	113	±10
Dwell Time	Hours	48	48 - 72

Approved Aeration Cells

N/A

Aeration Special Instructions

Extend the dwell time to compensate for the time the temperature is out of tolerance.

Handling of Aeration Nonconformities**Time Short or Temperature Low**☒ *Add time to meet minimum specified in heated aeration*☐ *Other:* _____**Temperature High**☐ *Continue process and evaluate after processing*☒ *Other* See Special Instructions

Ethylene Oxide Customer Specification

Customer: Pathway Innovation LLC Cycle/Specification #: PWI-7374 Rev: 1

6.0 Sample Retrieval:

- Samples Retrieved by: ☐ Customer ☒ STERIS
- Number and type of Samples: 7 Mail BI's
- Sample Retrieval Instructions: Retrieve within 8 hours aeration

7.0 Sample Shipping Instructions:

Ship Via: FEDEX

Charge to Customer Account: ☒ Yes ☐ No

Customer Account Number: 681054936

Ship to: STERIS Labs

Additional Information: Enter Run/Load number into the Reference field

8.0 Document Transmittal:

- ☒ Electronic (check if yes) E-mail addresses: arogers@pathwaynpi.com, btamayo@pathwaynpi.com
- ☐ Hard Copy (check if yes)

Ship Via: N/A

Charge to Customer Account: ☐ Yes ☒ No

Customer Account Number: N/A

Ship to: N/A

Additional Information: N/A

Documentation Required (check if yes)

Record	Electronically	Hard Copy
Cycle Printout	☒	☐
Customer Specification and Process Record	☒	☐
Product Receiving & Accountability	☒	☐
PCR Charts	☒	☐
Cycle Printout	☒	☐
Aeration Charts	☒	☐
Certificate of Processing	☒	☐

Ethylene Oxide Customer Specification

Customer: Pathway Innovation LLC Cycle/Specification #: PWI-7374 Rev: 1

9.0 Shipping Instructions

- ☒ Customer to arrange shipping
☐ STERIS to arrange shipping

Trucking Company #1

Preferred Carrier	Contact Name	Phone Number	Freight Terms	If Third Party, please provide billing Information:
N/A	N/A	N/A	☐ Collect ☐ Third Party	Address: N/A Contact Name: N/A Phone Number: N/A

Trucking Company #2 (if applicable)

Preferred Carrier	Contact Name	Phone Number	Freight Terms	If Third Party, please provide billing Information:
N/A	N/A	N/A	☐ Collect ☐ Third Party	Address: N/A Contact Name: N/A Phone Number: N/A

Special Instructions (Include all special loading, processing or shipping instructions. Attach separate sheet if necessary):

N/A

Ethylene Oxide Customer Specification

Customer: Pathway Innovation LLC Cycle/Specification #: PWI-7374 Rev: 1

10.0 Approvals

CUSTOMER Approval

Customer Signature: Aaron Rogers

Date: 08/06/2020

Print Name: AARON ROGERS

Title: DIR RA/BA

Note: If form is completed electronically, please print out and sign in ink.

STERIS Approval

Production Signature: [Signature]

Date: 9/2/2020

QS/RC Signature: [Signature]

Date: 09/03/2020

11.0 Change Description

Revision	Description of Changes	Effective Date
0	New Full Cycle in accordance with TP-0017 Rev B Pathway Innovation LLC	05/01/2018
1	Update Contact Name, Sample Placement, Retrieval and Shipping Instructions, Add Routine to Process Name	08/06/2020 <u>09/03/2020</u> (i)

[Signature]
 09/03/2020
 Entry error.
 Effective date once finalized.

STERIS: Ethylene Oxide Certificate Of Processing

Prepared PATHWAY INNOVATION LLC
For
Ethylene Oxide Run ID 1162-16739

<u>Product Code</u>	<u>Lot Number</u>	<u>Quantity (EA)</u>
PWI-7374 Full Cycle	SLQ-01242025	3
<u>Customer ItemID:</u>	SLQ-211810SPT	

Preconditioning

<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>
Preconditioning Room	09-Feb-2025 1:00 PM	10-Feb-2025 6:30 PM	1770

		<u>Min</u>	<u>Max</u>
Duration (min)	Spec:	960	1920
Relative Humidity (%RH)	Spec:	40.0	80.0
	Actual:	45.0	63.0
Temperature (°F)	Spec:	105.0	125.0
	Actual:	112.0	121.0

Vessel Processing

<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>
Chamber 6	10-Feb-2025 6:35 PM	11-Feb-2025 3:31 AM	536
Cycle ID: PWI-7374	Check Value: N/A		

		<u>Min</u>	<u>Max</u>
Duration (min)	Spec:	1	1440

EO Drum Information

<u>Drum ID</u>	<u>Drum Serial Number</u>	<u>Drum Lot Number</u>	<u>EO Usage (LB)</u>
2222	E001547	UTLX930340M24	17.0
Total EO Usage:			17.0

EO Drum Information (Spec)

	<u>EO Usage (LB)</u>
Min	1.0
Max	100.0

<u>Product Code</u>	<u>Lot Number</u>	<u>Quantity (EA)</u>	
Aeration			
<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>
Aeration Room	11-Feb-2025 3:41 AM	13-Feb-2025 4:36 AM	2935
	Total Aeration Time:		2935
		<u>Min</u>	<u>Max</u>
	Duration (min) Spec:	2880	4320
	Aeration in Room Temperature (°F) Spec:	103.0	123.0
	Actual:	111.0	115.0

EO Spec ID/Rev/Description ODMS Spec #5 Rev 1 (PWI-7374 Rev. 1 Medical Product)

Product meets Customer specifications; zero nonconformities occurred during this run.

Customer Specification ID: 769

Signature Manifest

Reviewed and E-Signed By: **Antonio Garner (QS & RC Technician I)**

DateTime Esigned 13-Feb-2025 11:06 AM

Document Content Revision 1

Operating facilities are in compliance with applicable state and federal regulations (FDA, NRC, EPA, and OSHA) and provide services under a quality system which meets the requirements of FDA QSR, EN/ISO 13485, and in alignment with EN ANSI/AAMI/ISO 11135.

Processing Location

43425 Business Park Dr
Temecula CA 92590
Phone: 951-694-9340

*** End of Ethylene Oxide Certificate of Processing Document ***

STERIS: Receiving Report
PATHWAY INNOVATION LLC (8130)
Ethylene Oxide

Product Received From:

PATHWAY INNOVATION LLC
8779 Cottonwood Ave
Suite 105
Santee CA 92071
325

Carrier: NorthStar
Trailer: N/A
Bill of Lading: 8590930
Seal: None
PO(s): P2028

Product counts provided by Customer

Receipt Date/Time: 09-Feb-2025 9:11 AM
Pallet Count: 1
Unloaded By: Alex Salazar
Counted By: Alex Salazar

Container Unload: ☐
Promise Date: 11-Mar-2025 10:11 AM

Item [Pkg]

PWI-7374 Full Cycle [1]

Customer Item ID: SLQ-211810SPT

Lot Number
SLQ-01242025

Quantity

Customer STERIS CD UOM
3 3 0 EA

Receipt Totals: 3 3 0

Inventory by Received Pallet

Pallet (SN) Item ID [Pkg]

1 (95880) PWI-7374 Full Cycle [1]
Customer Item ID: SLQ-211810SPT

Lot Number
SLQ-01242025

Quantity

Customer STERIS UOM
3 3 EA

Warehouse Location Pallet Total
N/A N/A 1

Receipt ID
1161-12293



STERIS: Receiving Report
PATHWAY INNOVATION LLC (8130)
Ethylene Oxide

Product Received From:

PATHWAY INNOVATION LLC
8779 Cottonwood Ave
Suite 105
Santee CA 92071
325

Carrier: **NorthStar**

Trailer: **N/A**

Bill of Lading: **8590930**

Seal: **None**

PO(s): **P2028**

Product counts provided by Customer

Receipt Date/Time: **09-Feb-2025 9:11 AM**

Pallet Count: **1**

Unloaded By: **Alex Salazar**

Counted By: **Alex Salazar**

Container Unload: ☐

Promise Date: **11-Mar-2025 10:11 AM**

Receipt ID
1161-12293



Signature Manifest

Reviewed By: **Alex Salazar (Material Handler Warehouse II)**

Date/Time Signed **09-Feb-2025 9:38 AM**

Document Content Revision 1

Processing Location

43425 Business Park Dr
Temecula CA 92590
Phone: 951-694-9340

GW1-08-004-01

Last Revised in Rel 2.0.0.0

Rel Date: **31-Mar-2017**

Document ID: 57775

Page 2 of 2

PATHWAY INNOVATION LLC

Packing List No: PL-SLQ-0057

Page No: 1 of 1


Company: STERIS AST
Attn: Diana Burns
Address: 43425 Business Park Drive
Temecula, CA 92590


Ship Date: 2/7/2025

Tracking No: TBD

Delivery Company: NorthStar Courier
Delivery Method: Courier
FOB: Santee, CA

Note: Routine Processing under PWI-7374

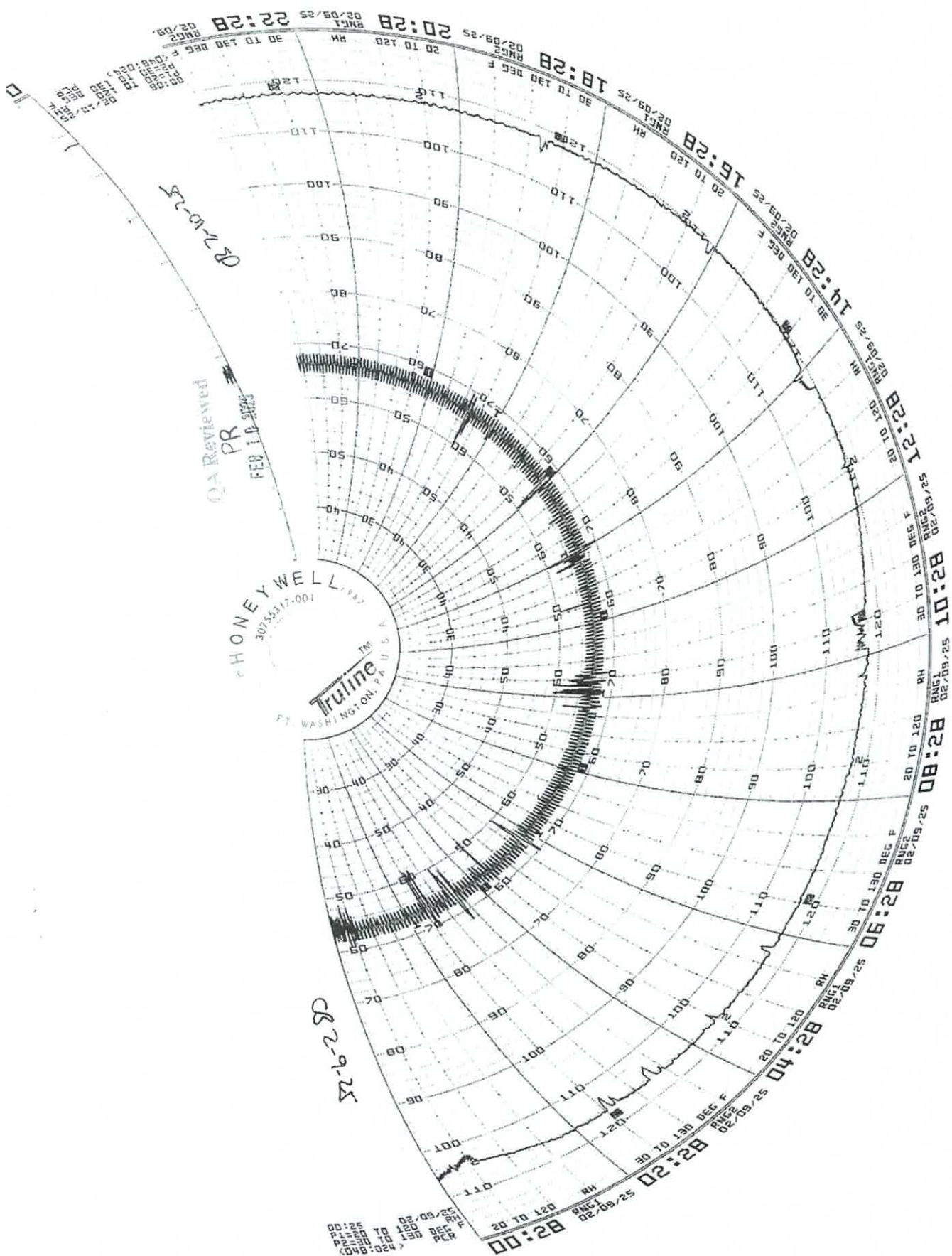

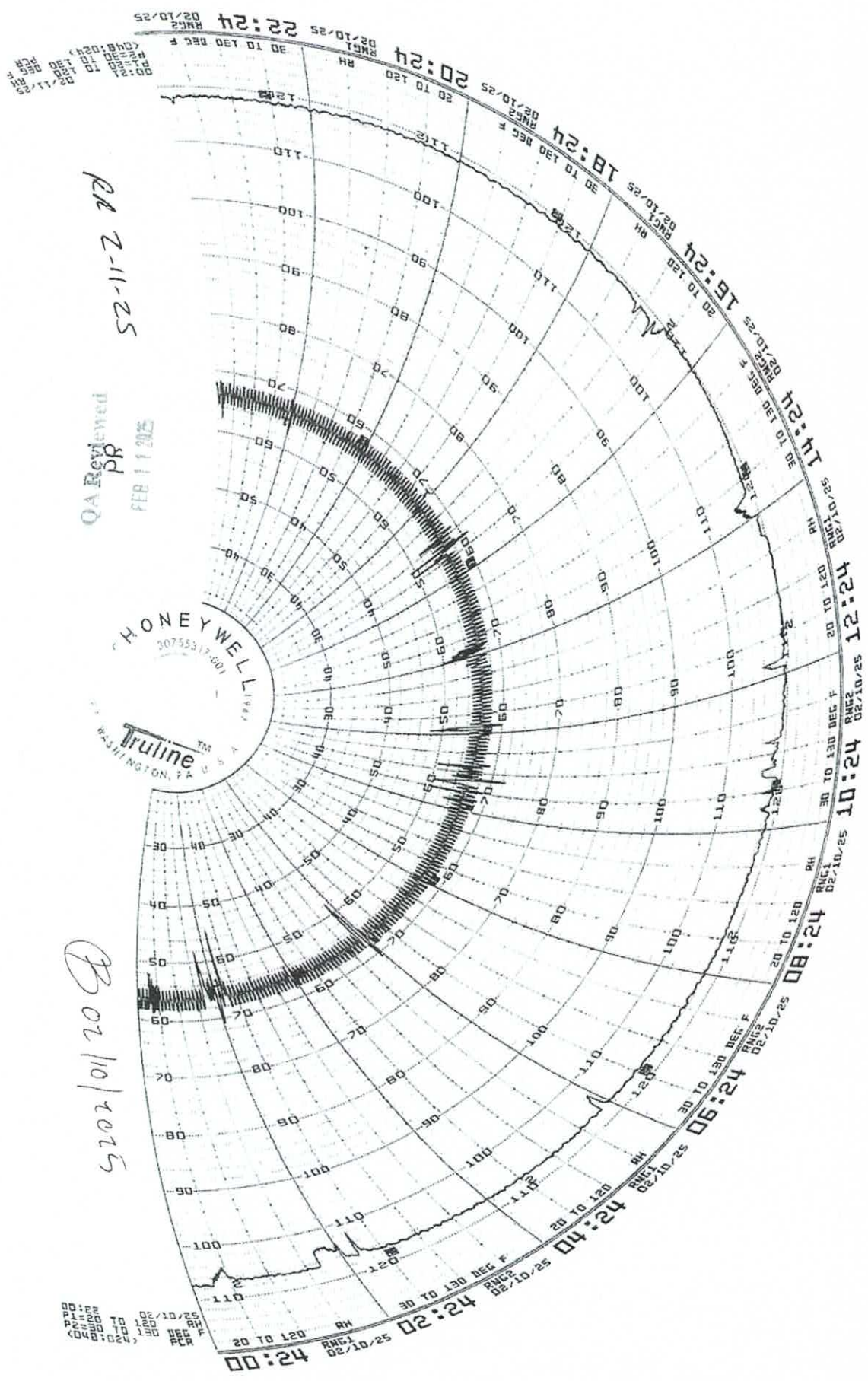

RA 2-11-25

QA Reviewed
PR
FEB 11 2025

02/10/2025

**STERIS**

Applied Sterilization Technologies

EO Process Report

STERIS

Applied Sterilization Technologies

43425 Business Park Drive

Temecula, CA 92590

Tel: (951) 694-9340

Fax: (951) 699-5981

Customer Name Pathway FC**Load #** PWI-1162-16739**Date Processed** 02/10/2025 18:35:50**Pressure Units** "HgA**Cycle Number** 7374**Cycle Last Modified:** 4/30/2018 1:56:49PM**Chamber #** 6**Cycle Counter #** 4,853**Temperature Units** Fahrenheit**Phase Name****Initial Evacuation**

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/10/2025 18:35:50	00:00:00	28.9	127
02/10/2025 18:40:51	00:05:00	21.5	126
02/10/2025 18:45:51	00:10:00	14.8	125
02/10/2025 18:50:51	00:15:00	8.0	124
02/10/2025 18:55:51	00:20:00	3.0	124
02/10/2025 18:57:54	00:22:03	2.0	125

Summary of Initial Evacuation

Pressure Start	28.9	Min Pressure	2.0	Min Temp	124	Pressure Change	-26.9
Pressure End	2.0	Max Pressure	28.9	Max Temp	127	Phase Time	00:22:03

Stabilization

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/10/2025 18:57:54	00:22:04	2.0	125
02/10/2025 18:59:54	00:24:04	2.1	126
02/10/2025 18:59:54	00:24:04	2.1	126

Summary of Stabilization

Pressure Start	2.0	Min Pressure	2.0	Min Temp	125	Pressure Change	0.1
Pressure End	2.1	Max Pressure	2.1	Max Temp	126	Phase Time	00:02:00

Leak Check

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/10/2025 18:59:54	00:24:04	2.1	126
02/10/2025 19:00:54	00:25:04	2.1	126
02/10/2025 19:01:54	00:26:04	2.1	126
02/10/2025 19:02:55	00:27:04	2.2	126
02/10/2025 19:03:54	00:28:04	2.2	126
02/10/2025 19:04:54	00:29:04	2.2	126

Summary of Leak Check

Pressure Start	2.1	Min Pressure	2.1	Min Temp	126	Pressure Change	0.1
Pressure End	2.2	Max Pressure	2.2	Max Temp	126	Phase Time	00:05:00

Initial Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/10/2025 19:04:54	00:29:04	2.2	126
02/10/2025 19:09:55	00:34:04	9.3	128
02/10/2025 19:14:54	00:39:04	16.8	128
02/10/2025 19:19:54	00:44:04	24.3	128
02/10/2025 19:22:24	00:46:34	28.0	128

Customer Name Pathway FC
Load # PWI-1162-16739

STERIS
Applied Sterilization Technologies
Temecula, CA

Summary of Initial Nitrogen Wash

Pressure Start	2.2	Min Pressure	2.2	Min Temp	126	Pressure Change	25.8
Pressure End	28.0	Max Pressure	28.0	Max Temp	128	Phase Time	00:17:30

Initial Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/10/2025 19:22:24	00:46:34	28.0	128
02/10/2025 19:27:24	00:51:34	21.3	124
02/10/2025 19:32:25	00:56:34	14.3	124
02/10/2025 19:37:25	01:01:34	7.3	123
02/10/2025 19:42:25	01:06:34	2.7	123
02/10/2025 19:43:54	01:08:03	2.0	123

Summary of Initial Evacuation

Pressure Start	28.0	Min Pressure	2.0	Min Temp	123	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	128	Phase Time	00:21:29

Humidification

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	Avg. % RH
02/10/2025 19:43:54	01:08:03	2.0	123	8
02/10/2025 19:44:54	01:09:03	2.6	124	25
02/10/2025 19:45:54	01:10:03	3.0	125	35
02/10/2025 19:46:00	01:10:09	3.0	125	36

Summary of Humidification

Pressure Start	2.0	Min Pressure	2.0	Min Temp	123	Pressure Change	1.0	Min RH	8
Pressure End	3.0	Max Pressure	3.0	Max Temp	125	Phase Time	00:02:06	Max RH	36

Humidity Dwell

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	Avg. % RH
02/10/2025 19:46:00	01:10:09	3.0	125	36
02/10/2025 19:50:59	01:15:09	3.0	126	37
02/10/2025 19:55:59	01:20:09	3.0	125	38
02/10/2025 20:01:00	01:25:09	3.0	125	38
02/10/2025 20:06:00	01:30:09	3.0	125	38
02/10/2025 20:11:00	01:35:09	3.1	125	38
02/10/2025 20:16:00	01:40:10	3.1	125	38
02/10/2025 20:21:00	01:45:10	3.1	125	38
02/10/2025 20:26:00	01:50:10	3.1	125	38
02/10/2025 20:31:00	01:55:10	3.1	126	38
02/10/2025 20:36:00	02:00:10	3.1	125	38
02/10/2025 20:41:01	02:05:10	3.1	125	39
02/10/2025 20:46:00	02:10:09	3.1	125	39

Summary of Humidity Dwell

Pressure Start	3.0	Min Pressure	3.0	Min Temp	125	Pressure Change	0.1	Min RH	35
Pressure End	3.1	Max Pressure	3.1	Max Temp	126	Phase Time	01:00:00	Max RH	40

Sterilant Injection

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	EO Drum 1 (lb)	EO Drum 2 (lb)	EO CONC
02/10/2025 20:46:00	02:10:09	3.1	125	337	0	0
02/10/2025 20:51:00	02:15:09	7.7	126	333	0	225
02/10/2025 20:56:00	02:20:09	12.8	126	326	0	527
02/10/2025 20:57:31	02:21:40	14.3	126	323	0	611

Customer Name Pathway FC
Load # PWI-1162-16739

STERIS
Applied Sterilization Technologies
Temecula, CA

Summary of Sterilant Injection

Pressure Start 3.1	Min Pressure 3.1	Min Temp 125	Drum 1 Start 337	Drum 2 Start 0
Pressure End 14.3	Max Pressure 14.3	Max Temp 126	Drum 1 End 323	Drum 2 End 0
Pressure Change 11.2	MIN EO CONC 0.00		EtO Used 1 14	EtO Used 2 0
Phase Time 00:11:31	MAX EO CONC 611.00			
			Total EO Weight 14	

Exposure

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	EO Drum 1 (lb)	EO Drum 2 (lb)	EO CONC
02/10/2025 20:57:31	02:21:40	14.3	126	323	0	611
02/10/2025 21:02:31	02:26:40	14.2	125	324	0	632
02/10/2025 21:07:31	02:31:41	14.1	125	324	0	628
02/10/2025 21:12:32	02:36:41	14.1	125	324	0	627
02/10/2025 21:17:31	02:41:41	14.3	125	323	0	635
02/10/2025 21:22:31	02:46:41	14.3	125	323	0	633
02/10/2025 21:27:31	02:51:41	14.3	125	323	0	632
02/10/2025 21:32:32	02:56:41	14.3	125	323	0	628
02/10/2025 21:37:32	03:01:41	14.3	125	323	0	628
02/10/2025 21:42:32	03:06:41	14.3	125	323	0	627
02/10/2025 21:47:32	03:11:41	14.3	125	323	0	624
02/10/2025 21:52:31	03:16:41	14.3	125	323	0	624
02/10/2025 21:57:32	03:21:41	14.3	125	322	0	622
02/10/2025 22:02:32	03:26:41	14.3	125	322	0	621
02/10/2025 22:07:32	03:31:42	14.3	125	322	0	621
02/10/2025 22:12:32	03:36:42	14.3	125	322	0	618
02/10/2025 22:17:33	03:41:42	14.3	125	322	0	618
02/10/2025 22:22:32	03:46:42	14.3	125	322	0	618
02/10/2025 22:27:32	03:51:42	14.3	125	322	0	615
02/10/2025 22:32:32	03:56:42	14.3	125	322	0	615
02/10/2025 22:37:32	04:01:42	14.3	125	322	0	616
02/10/2025 22:42:33	04:06:42	14.3	125	322	0	614
02/10/2025 22:47:33	04:11:42	14.3	125	322	0	615
02/10/2025 22:52:33	04:16:42	14.3	125	322	0	612
02/10/2025 22:57:32	04:21:42	14.3	125	322	0	612
02/10/2025 23:02:33	04:26:42	14.3	125	321	0	612
02/10/2025 23:07:33	04:31:43	14.3	125	321	0	609
02/10/2025 23:12:33	04:36:43	14.3	125	321	0	610
02/10/2025 23:17:33	04:41:43	14.3	125	321	0	610
02/10/2025 23:22:34	04:46:43	14.3	125	321	0	608
02/10/2025 23:27:34	04:51:43	14.3	125	321	0	609
02/10/2025 23:32:33	04:56:43	14.3	125	321	0	608
02/10/2025 23:37:33	05:01:43	14.3	125	321	0	607
02/10/2025 23:42:33	05:06:43	14.3	125	321	0	608
02/10/2025 23:47:34	05:11:43	14.3	125	321	0	606
02/10/2025 23:52:34	05:16:43	14.3	125	321	0	606
02/10/2025 23:57:34	05:21:43	14.3	125	321	0	606
02/11/2025 00:02:34	05:26:43	14.3	125	321	0	606
02/11/2025 00:07:34	05:31:44	14.3	125	321	0	606
02/11/2025 00:12:34	05:36:44	14.3	125	321	0	605
02/11/2025 00:17:34	05:41:44	14.3	125	320	0	604
02/11/2025 00:22:34	05:46:44	14.3	125	320	0	605
02/11/2025 00:27:34	05:51:44	14.3	125	321	0	603
02/11/2025 00:32:35	05:56:44	14.3	125	320	0	604
02/11/2025 00:37:35	06:01:44	14.3	125	320	0	604
02/11/2025 00:42:34	06:06:44	14.3	125	320	0	603
02/11/2025 00:47:34	06:11:44	14.3	125	320	0	603
02/11/2025 00:52:35	06:16:44	14.3	125	320	0	603
02/11/2025 00:57:31	06:21:40	14.3	125	320	0	602

Customer Name Pathway FC
Load # PWI-1162-16739

STERIS
Applied Sterilization Technologies
Temecula, CA

Summary of Exposure

Pressure Start 14.3	Min Pressure 14.1	Min Temp 125	Drum 1 Start 323	Drum 2 Start 0
Pressure End 14.3	Max Pressure 14.4	Max Temp 126	Drum 1 End 320	Drum 2 End 0
Pressure Change 0.0	MIN EO CONC 602.00		EtO Used 1 3	EtO Used 2 0
Phase Time 04:00:00	MAX EO CONC 637.00		Total EO Weight 3	

Sterilant Removal

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/11/2025 00:57:31	06:21:41	14.3	125
02/11/2025 01:02:31	06:26:41	9.9	123
02/11/2025 01:07:31	06:31:41	5.3	123
02/11/2025 01:12:25	06:36:34	2.0	123

Summary of Sterilant Removal

Pressure Start 14.3	Min Pressure 2.0	Min Temp 123	Pressure Change -12.3
Pressure End 2.0	Max Pressure 14.3	Max Temp 125	Phase Time 00:14:53

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/11/2025 01:12:25	06:36:34	2.0	123
02/11/2025 01:17:24	06:41:34	9.8	126
02/11/2025 01:22:24	06:46:34	17.0	127
02/11/2025 01:27:24	06:51:34	24.3	127
02/11/2025 01:29:59	06:54:08	28.0	127

Summary of Post Nitrogen Wash

Pressure Start 2.0	Min Pressure 2.0	Min Temp 123	Pressure Change 26.0
Pressure End 28.0	Max Pressure 28.0	Max Temp 127	Phase Time 00:17:35

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/11/2025 01:29:59	06:54:08	28.0	127
02/11/2025 01:34:59	06:59:09	21.3	124
02/11/2025 01:39:59	07:04:09	14.3	123
02/11/2025 01:45:00	07:09:09	7.4	122
02/11/2025 01:49:59	07:14:09	2.8	122
02/11/2025 01:51:41	07:15:50	2.0	122

Summary of Post Evacuation

Pressure Start 28.0	Min Pressure 2.0	Min Temp 122	Pressure Change -26.0
Pressure End 2.0	Max Pressure 28.0	Max Temp 127	Phase Time 00:21:42

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/11/2025 01:51:41	07:15:50	2.0	122
02/11/2025 01:56:40	07:20:51	9.1	126
02/11/2025 02:01:41	07:25:51	16.6	127
02/11/2025 02:06:41	07:30:51	24.1	128
02/11/2025 02:09:20	07:33:29	28.0	128

Summary of Post Nitrogen Wash

Pressure Start 2.0	Min Pressure 2.0	Min Temp 122	Pressure Change 26.0
Pressure End 28.0	Max Pressure 28.0	Max Temp 128	Phase Time 00:17:39

Customer Name Pathway FC
Load # PWI-1162-16739

STERIS
Applied Sterilization Technologies
Temecula, CA

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/11/2025 02:09:20	07:33:29	28.0	128
02/11/2025 02:14:20	07:38:29	21.3	124
02/11/2025 02:19:20	07:43:29	14.3	123
02/11/2025 02:24:19	07:48:29	7.3	122
02/11/2025 02:29:20	07:53:29	2.7	122
02/11/2025 02:30:47	07:54:57	2.0	122

Summary of Post Evacuation

Pressure Start	28.0	Min Pressure	2.0	Min Temp	121	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	128	Phase Time	00:21:27

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/11/2025 02:30:47	07:54:57	2.0	122
02/11/2025 02:35:47	07:59:57	9.1	126
02/11/2025 02:40:48	08:04:57	16.6	127
02/11/2025 02:45:47	08:09:57	24.1	128
02/11/2025 02:48:25	08:12:35	28.0	128

Summary of Post Nitrogen Wash

Pressure Start	2.0	Min Pressure	2.0	Min Temp	122	Pressure Change	26.0
Pressure End	28.0	Max Pressure	28.0	Max Temp	128	Phase Time	00:17:39

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/11/2025 02:48:25	08:12:35	28.0	128
02/11/2025 02:53:26	08:17:35	21.3	124
02/11/2025 02:58:26	08:22:36	14.3	123
02/11/2025 03:03:26	08:27:36	7.3	122
02/11/2025 03:08:26	08:32:36	2.6	121
02/11/2025 03:09:47	08:33:57	2.0	122

Summary of Post Evacuation

Pressure Start	28.0	Min Pressure	2.0	Min Temp	121	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	128	Phase Time	00:21:21

Air Backfill

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
02/11/2025 03:09:47	08:33:57	2.0	122
02/11/2025 03:14:47	08:38:57	7.8	126
02/11/2025 03:19:47	08:43:57	14.1	126
02/11/2025 03:24:47	08:48:57	20.3	125
02/11/2025 03:29:47	08:53:57	26.5	125
02/11/2025 03:31:02	08:55:12	27.9	125

Summary of Air Backfill


Pressure Start	2.0	Min Pressure	2.0	Min Temp	122	Pressure Change	25.9
Pressure End	27.9	Max Pressure	27.9	Max Temp	126	Phase Time	00:21:15


Customer Name Pathway FC
Load # PWI-1162-16739


STERIS
Applied Sterilization Technologies
Temecula, CA

End of Cycle

Production Review  Date 07/11/2025
Signature

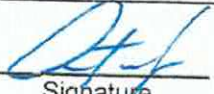Quality Review  Date 2.13.25
Signature

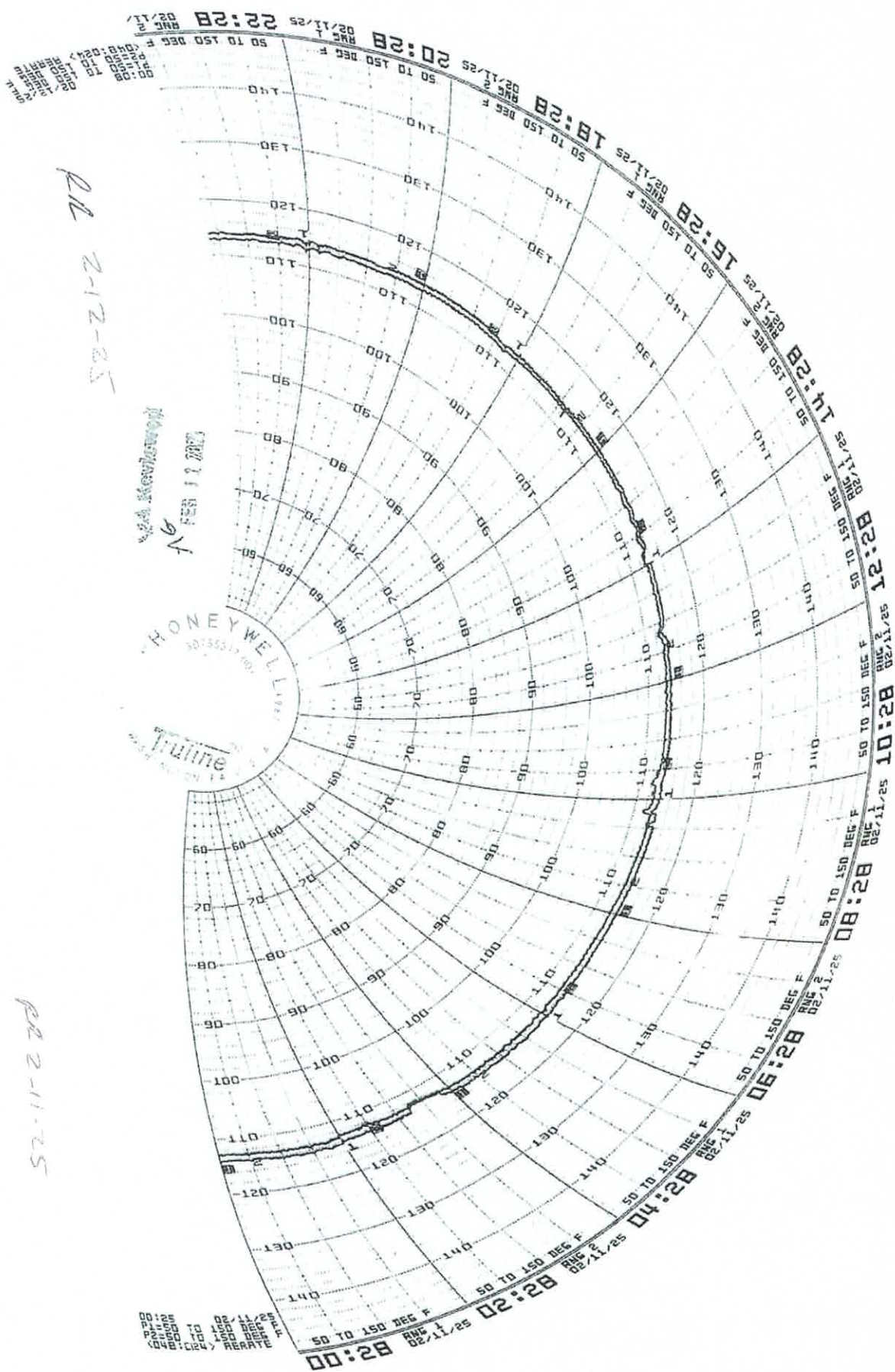

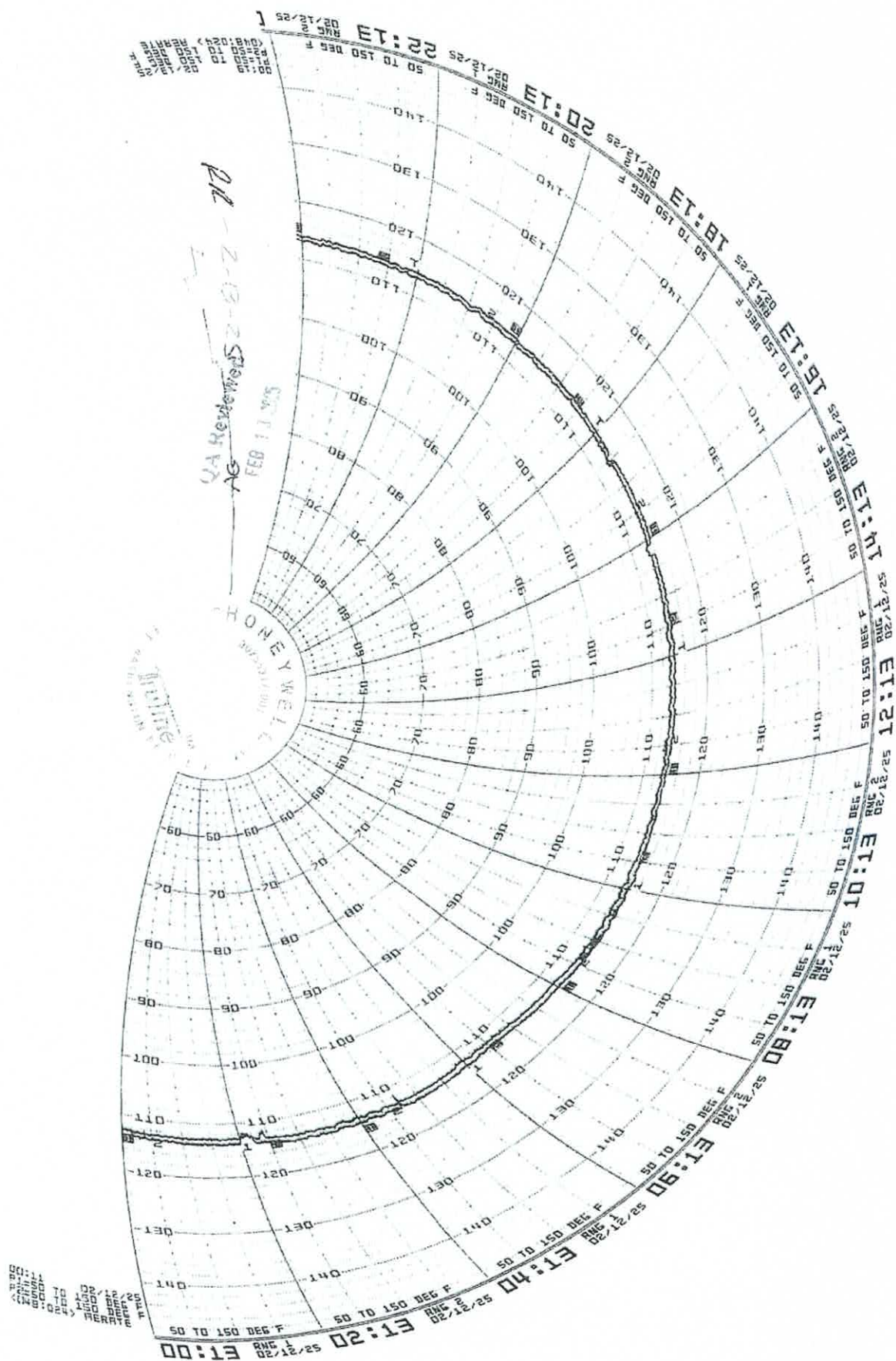

RR 2-12-25

**LGGs, Inc. Technical Report**

SPONSOR: Pathway Medtech LLC
8779 Cottonwood Ave, Suite 105
Santee, CA 92071

Lab No. 25-0269-00
P.O.No. P2013
Lot No. SLQ-01242025
Date Received: 01-29-25

Analysis: Limulus Amebocyte Lysate (LAL) Kinetic Turbidimetric Assay

Test Article: SLQ-8102

Condition of Test Article at Receipt: Good

Procedure: Ten (10) products were extracted following procedures outlined in test specification number, BEKT-PTW-001-00. One-tenth ml of the test and/or control solution and one-tenth of lysate solution, reconstituted per suppliers directions, were placed in sterile microplate wells and incubated at 37°C for a minimum of one hour in the Kinetic QCL Reader. The reaction time with the LAL/substrate was determined turbidimetrically and compared to a standard curve to determine endotoxin concentration. Positive Product Control (PPC) was simultaneously prepared to demonstrate that the test article did not significantly interfere with the LAL/substrate reaction.

Results:	Test Article Extract Dilution:	<u>2.0</u>
	PPC Testing:	<u>Acceptable</u>
	Test Article Extract:	<u><0.02 EU/ml; <0.8 EU/device</u>
	Date Tested:	<u>01-29-25</u>

Key to Results: Type of Product	Current FDA Requirements
1) Medical Device	1) Less than or equal to 0.5 EU/ml 2) Less than or equal to 20.0 EU/device

Comments: Test article as submitted conforms to FDA and USP requirements for end-product release of medical devices.

Deviations: None

Completed: 1/30/2025 Tech RR / GS Approved by: 
Technical Manager (or deputy)

The conclusions and data herein apply only to the items tested and are not indicative of the quality or condition of apparently identical product. The report is submitted for the exclusive use of the sponsor. No use or mention of LGGs, Inc. in connection with any form of advertising or public announcement or reproduction of the report except in full may be made without written approval of LGGs, Inc.