

Certificate of Conformance

Pathway LLC certifies that the product indicated below meets all requirements documented in its approved specification and the company's Quality Management System.

Part Number: **SLQ-211810SPT**
Description: **SILQ 2-way Clearteact SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack**
Revision: **B**
Lot Number: /
Serial Number(s): **SLQ-81020515241**
Expiration Date: **2027-05-15**
Quantity Released
(Case Quantity): **38 Packs of 10**

The product indicated can be transferred to a "release" inventory location and is cleared for distribution. The individual signing below authorizes the release of the product.



Signature (Quality Assurance Representative)


02/29/2024

Date

Date Inspected	Part No.	Rev.	Description
07/29/2024	SLQ-211810SPT	B	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack
Lot No.	Lot QTY		
SLQ-81020515241	38 Packs of 10		
Identification of the product is traceable between the DHR, sterilization records, and laboratory report	DHR is completed and meets good documentation practices per QP-003.	Sterilization processing record meets parameters of validated process	Sterility laboratory report meets acceptance criteria indicated in QP-010
Endotoxin laboratory report (LAL) meets acceptance criteria per WI-0021	No Damage to Outer Packing on Shipment back from Sterilizer	Verify Quantity Received back from Sterilizer matches Packing List	
Specification	QP-010 Revision	F	F
Lower Spec Limit (LSL)	N/A	N/A	N/A
Upper Spec Limit (USL)	N/A	N/A	N/A
UOM	N/A	N/A	N/A
Measuring Equipment Type	Visual	Visual	Visual
Asset No.	N/A	N/A	N/A
Cal Due Date	N/A	N/A	N/A
AQL	See Below	See Below	See Below
Sample Size	DHR	Sterilization Processing Record	Laboratory Report
	DHR	Sterilization Processing Record	Sterility Laboratory Report
	Endotoxin (LAL) Laboratory Report	Shipment	Shipment
Sample No.	Variable/ Attribute Data		
1	P	P	P
Pass QTY	1	1	1
Fail QTY	0	0	0
Result (P/F)	P	P	P
Min	N/A	N/A	N/A
Max	N/A	N/A	N/A
Mean	N/A	N/A	N/A
StDev	N/A	N/A	N/A
Comments			
ONIA mly 07/29/2024			
Disposition:	ACCEPT	NCR No:	24-027
Inspected By:	Melissa Garcia	Date:	07/29/2024
Verified By:	Larry Racusin	Date:	07/29/2024

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
 1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY PATHWAY PRIOR TO PRODUCTION.
 2. DRAWING IS FOR REFERENCE ONLY. REFER TO THE RELEASED BILL OF MATERIALS FOR THE COMPLETE PARTS LIST.
 3. UP UNTIL SEALING THE POUCH, PRODUCT SHALL BE PACKAGED IN A CLEANROOM ENVIRONMENT. POST SEALING THE POUCH, ASSEMBLY WILL CONTINUE OUTSIDE OF THE CLEANROOM ENVIRONMENT.
2. ASSEMBLY PROCEDURE
 1. APPLY POLYMER SOLUTION TO THE CATHETER IN THE HATCHED AREA SHOWN. ALLOW THE CATHETER TO DRY OVERNIGHT.
 2. INSERT THE CATHETER INTO A POLYETHYLENE SLEEVE SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CLOSED END.
 3. PLACE THE POUCH LABEL ON THE CLEAR SIDE OF THE POUCH. WHEN VIEWING THE CLEAR SIDE OF THE POUCH FROM ABOVE, WITH THE CHEVRON TO THE RIGHT, THE LABEL TEXT SHALL READ LEFT TO RIGHT, WITH BEST EFFORT, APPLY THE LABEL SUCH THAT IT IS CENTERED.
 4. INSERT THE SLEEVE / CATHETER INTO THE POUCH SUCH THAT THE TIP OF THE CATHETER IS ALIGNED TOWARDS THE CHEVRON END OF THE POUCH.
 5. SEAL THE POUCH.
 6. PLACE TO SLEEVED, POUCHED, SEALED, AND LABELED CATHETERS INTO ONE SHELF CARTON. APPLY TWO SHIPPER BOX LABELS TO THE LONG, WIDER SIDE OF THE CARTON. VERIFY THAT THE TEXT READS LEFT TO RIGHT WHEN THE OPENING OF THE CARTON IS TO THE LEFT. WITH BEST EFFORT, CENTER THE LABEL ON ITS RESPECTIVE FACE OF THE CARTON.
 7. PLACE TO SHELF CARTONS INTO ONE SHIPPER BOX. SEAL WITH PACKAGING TAPE AS NEEDED. PLACE A SHIPPER BOX LABEL IN THE UPPER LEFT CORNER ON TWO ADJACENT FACES, READING LEFT TO RIGHT.
3. DIMENSIONAL
 1. (I) INDICATES LABEL REFERENCE DIMENSIONS.
 2. (O) INDICATES CRITICAL DIMENSIONS.
4. CLEANING AND HANDLING
 1. ALL PRODUCT, UP UNTIL, INSERTING THE POUCHES INTO SHELF CARTONS, SHALL BE HANDLED BY GLOVED HANDS APPROPRIATE FOR A CLEANROOM ENVIRONMENT.
5. FUNCTIONAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
 1. THE CATHETER'S ABSORBANCE VALUE AT 565 NANOMETERS SHALL BE GREATER THAN OR EQUAL TO 0.300 AU. REFER TO TM-HDX-0001.
 2. THE CATHETER'S ABSORBANCE VALUE AT 440 NANOMETERS SHALL BE LESS THAN OR EQUAL TO 0.020 AU. REFER TO TM-HDX-0002.
 3. THE PEEL STRENGTH FOR A 1 INCH WIDE SEALED SECTION OF THE STERILE BARRIER PACKAGING SHALL BE ≥ 1.00 LB WHEN TESTED PER TM-0003. IT IS NOT REQUIRED TO PLACE CATHETERS IN THESE TEST POUCHES, OR LABEL THEM. INSPECT N=3 SAMPLES, 1 AT THE BEGINNING, MIDDLE, AND END OF THE MANUFACTURING RUN.
6. VISUAL INSPECTION (AQL = 6.5 UNLESS OTHERWISE STATED)
 1. THE CATHETER SHALL BE FREE OF CURED POLYMER AND DEBRIS.
 2. THE CATHETER SHALL BE PROPERLY ORIENTED IN THE SLEEVE AND POUCH WITH THE LABEL APPLIED TO THE APPROPRIATE FACE AND IN THE CORRECT ORIENTATION.
 3. INSIDE THE STERILE BARRIER, DEVICE IS TO BE FREE OF PARTICULATES WHEN VIEWED AT 18 INCHES AWAY UNDER AMBIENT ROOM LIGHTING AND WITHOUT THE USE OF MAGNIFICATION.
 4. THE STERILE BARRIER PACKAGING SHALL NOT INCLUDE ANY OF THE FOLLOWING DEFECTS:
 1. UNSEALED AREAS, NONHOMOGENEOUS OR UNSEALED AREAS, OVERSEALED AREAS, ANROW SEALS, CHANNELS, WRINKLES / FOLDOVERS / CRACKS, OR TEARS / PINHOLES.
 5. THE SHELF CARTON SHALL HAVE THE CORRECT NUMBER OF CATHETERS INSIDE AND FACING THE CORRECT ORIENTATION. THE LABELS SHALL BE AFFIXED IN THE CORRECT LOCATION AND ORIENTATION.
 6. THE SHIPPER BOX SHALL HAVE THE CORRECT NUMBER OF SHELF CARTONS. THE LABELS SHALL BE APPLIED IN THE CORRECT LOCATION AND ORIENTATION.

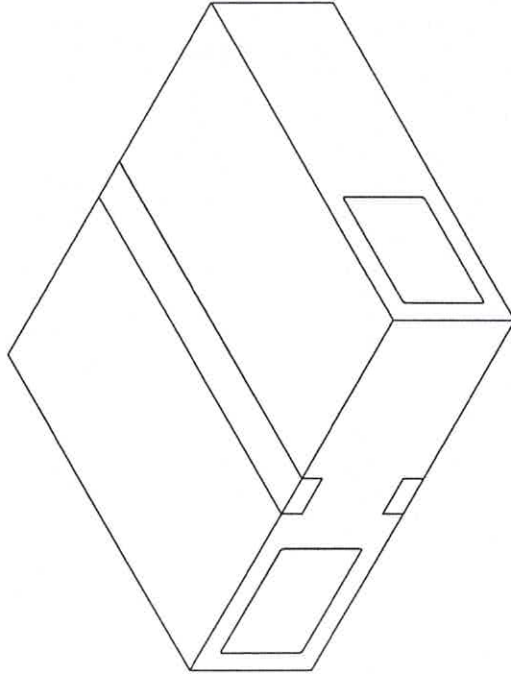

REVISIONS			
ZONE	REV.	DESCRIPTION	DATE
A		INITIAL RELEASE	

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
ANGULAR: MACH+0.2 BEND ± 2.0
NO PLACE DECIMAL ± 0.5
ONE PLACE DECIMAL ± 0.3
TWO PLACE DECIMAL ± 0.15
THREE PLACE DECIMAL ± 0.05

PROPERTY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS
DRAWING IS THE SOLE PROPERTY OF
PATHWAY NP. ANY REPRODUCTION IN
PART OR AS A WHOLE WITHOUT THE
WRITTEN PERMISSION OF PATHWAY NP
IS PROHIBITED.

pathway
New Product Introduction

TITLE:
CATHETER, SILICONE, NY FOLEY, NON-
STERILE, 10-PACK, MASTER DRAWING

SIZE DWG. NO. REV
B SLQ-81XX A

DO NOT SCALE DRAWING


SCALE: 1:5

SHEET 1 OF 3



REGISTERED



CERTIFIED 13485

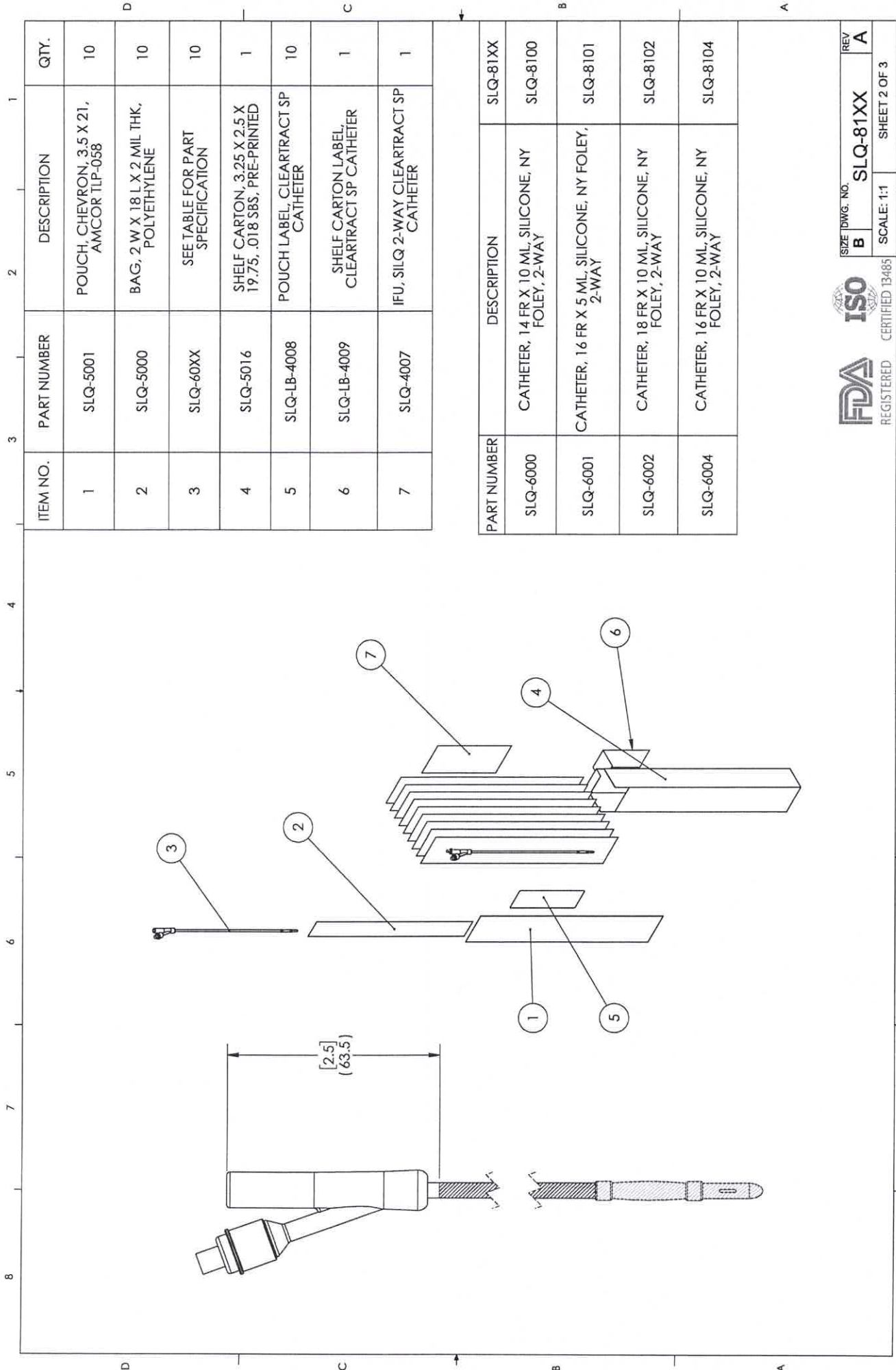

ITEM NO.	PART NUMBER	DESCRIPTION	QTY.
1	SLQ-5001	POUCH, CHEVRON, 3.5 X 21, AMCOR TLP-058	10
2	SLQ-5000	BAG, 2 W X 18 L X 2 MIL THK, POLYETHYLENE	10
3	SLQ-60XX	SEE TABLE FOR PART SPECIFICATION	10
4	SLQ-5016	SHELF CARTON, 3.25 X 2.5 X 19.75, .018 SBS, PRE-PRINTED	1
5	SLQ-LB-4008	POUCH LABEL, CLEARTRACT SP CATHETER	10
6	SLQ-LB-4009	SHELF CARTON LABEL, CLEARTRACT SP CATHETER	1
7	SLQ-4007	IFU, SILQ 2-WAY CLEARTRACT SP CATHETER	1

PART NUMBER	DESCRIPTION	
SLQ-6000	CATHETER, 14 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-81XX
SLQ-6001	CATHETER, 16 FR X 5 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8100
SLQ-6002	CATHETER, 18 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8101
SLQ-6004	CATHETER, 16 FR X 10 ML, SILICONE, NY FOLEY, 2-WAY	SLQ-8102
		SLQ-8104

ITEM NO.	PARTNO	DESCRIPTION	QTY.
8	SLQ-5012	BOX, 20 IN. L X 18 IN. W X 6 IN. D, 200 LB TEST	1
9	PWY-5000	TAPE, PACKAGING, 2 MIL, 2 IN X 110 YDS, CLEAR	AR
10	SLQ-LB-4010	SHIPPER BOX LABEL, CLEARTRACT SP CATHETER	2

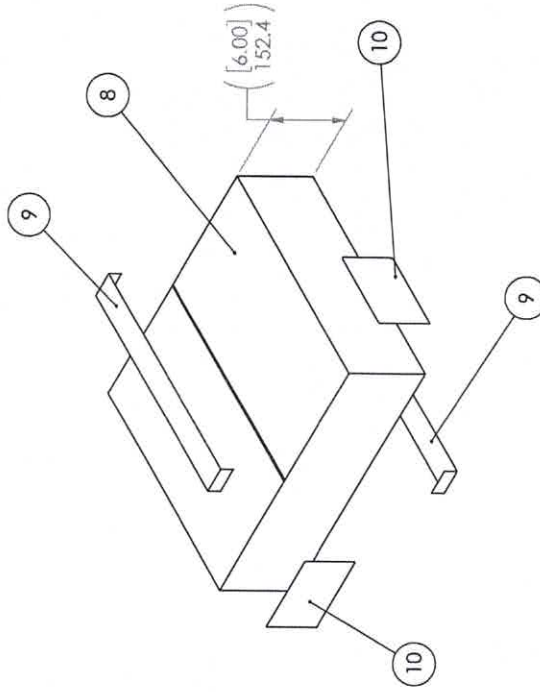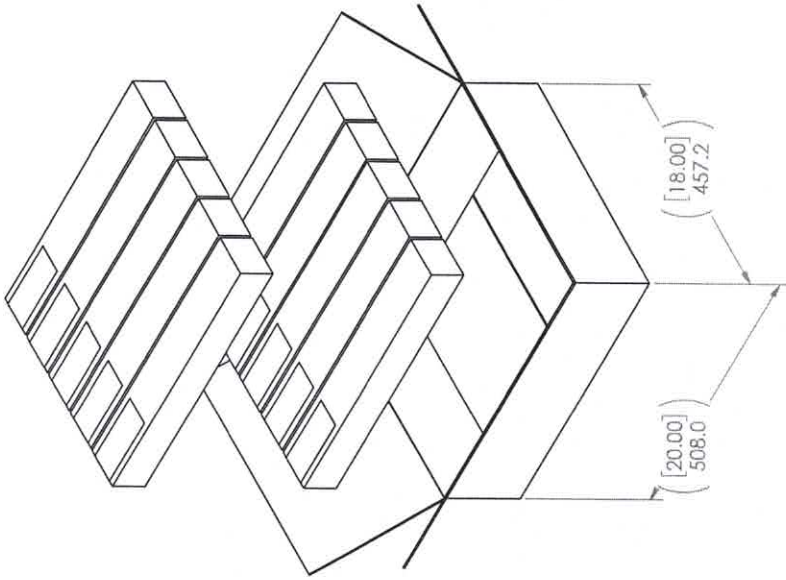

SIZE DWG. NO.	REV
B	A
SLQ-81XX	
SCALE: 1:1	
SHEET 3 OF 3	



Pathway Medtech, LLC.

Printed by Sean Tamayo

Printed on 07/19/2024

Local time zone (GMT-08:00) Pacific Time (America/Los Angeles)

Page Arena > Workspace Items > Bill of Materials > Indented

All information contained in this document is proprietary and confidential.

SLQ-8102

In Production

SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack



Revision »A - ECO-001669 Status: Effective, Status Date: 05/14/2024 03:24:42 PM, Shared

Contains 16 first-level Items, 31 line Items, 24 unique Items, 23 of which are shared.

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE
1 SLQ-6002 rev A	Catheter, 18 Fr x 10 mL, Silicone, NY Foley, 2-Way	Fabricated Part	In Prod
2 SLQ-3000 rev B	Solution, SILQ Polymer	Purchased Component	In Prod
SUBSTITUTES			
	RANK#ITEM NUMBER	ITEM NAME	NOTES
	1 HDX-3000 rev A	Solution, 80 ML Hydrophilix Polymer	
3 SLQ-5000 rev A	Bag, 2 W x 18 L x 2 Mil Thk, Polyethylene	Packaging	In Prod
4 SLQ-5001 rev A	Pouch, Chevron, 3.5 X 21, Amcor TLP-058	Packaging	In Prod
SUBSTITUTES			
	RANK#ITEM NUMBER	ITEM NAME	NOTES
	1 HDX-5001 rev A	Pouch, Chevron, 10 EA 3.5 X 21, Amcor TLP-058	
5 SLQ-LB-4008 rev A	Label, Pouch, SILQ ClearTract SPT Catheter	Labeling	In Prod

#	ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	     		
6	SLQ-5016 rev A	Shelf Carton, 3.25 x 2.5 x 19.75, .018 SBS, Pre-Printed	 Packaging	In Prod		3	1 EA
SUBSTITUTES							
		RANK# ITEM NUMBER	ITEM NAME		NOTES		
	1	SLQ-5015 rev A	Shelf Carton, 3.25 1 EA x 2.5 x 19.75, .018 SBS				
	2	SLQ-5014 rev A	Shelf Carton, 3.25 1 EA x 2.5 x 19.75, .024 SBS				
7	SLQ-LB-4009 rev A	Label, Box, SILQ ClearTract SPT Catheter	 Labeling	In Prod		2	2 EA
8	SLQ-4007 rev A	Artwork, IFU, SILQ 2-Way ClearTract SPT Catheter	 Labeling	In Prod		2	1 EA
9	SLQ-5012 rev A	Box, 20 in. L x 18 in. W x 6 in. D, 200 LB Test	 Packaging	In Prod		2	0.1 EA
10	PWY-5000 rev A	Tape, Packaging, 2 Mil, 2 in x 110 yds, Clear	 Packaging	In Prod		2	5.2 IN
11	SLQ-LB-4010 rev A	Label, Shipper Box, SILQ ClearTract SPT Catheter	 Labeling	In Prod		2	0.2 EA
12	MP-SLQ-0002 rev B	Manufacturing Procedure, SILQ ClearTract Catheter Surface Treatment Using the Automation Machine	 Manufacturing Procedure	In Prod		2	1 EA
13	FM-SLQ-0001 rev A	Form, Run Sheet for SILQ Automation Machine	 Form	In Prod		2	1 EA
14	TM-HDX-0001 rev B	Test Method, UV Light Absorbance, Hydrophilix Catheters	 Test Method	In Prod		2	1 each
15	TM-HDX-0002 rev A	Test Method, Rinse Test, Hydrophilix Catheter	 Test Method	In Prod		2	1 each
16	MP-SLQ-0001 rev B	Manufacturing Procedure, Labeling and Packaging, SILQ ClearTract Catheters	 Manufacturing Procedure	In Prod		2	1 EA

ITEM NO.	PARTNO	DESCRIPTION	QTY.
8	SLQ-5012	BOX, 20 IN. L X 18 IN. W X 6 IN. D, 200 LB TEST	1
9	PWY-5000	TAPE, PACKAGING, 2 MIL, 2 IN X 110 YDS, CLEAR	AR
10	SLQ-LB-4010	SHIPPER BOX LABEL, CLEARTRACT SP CATHETER	2

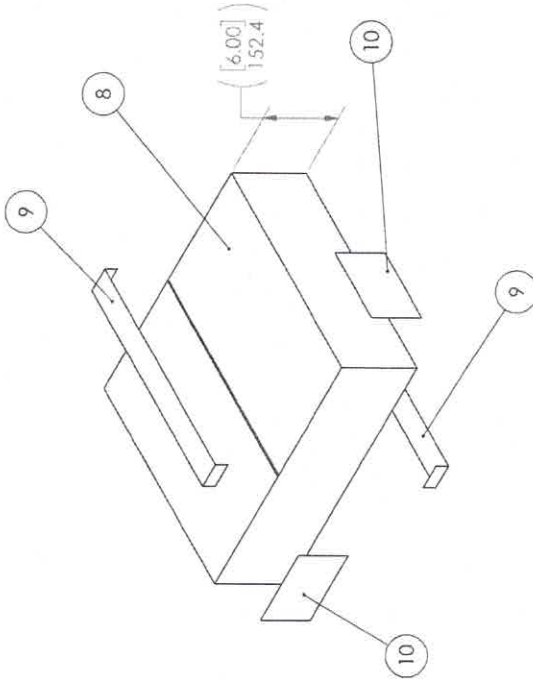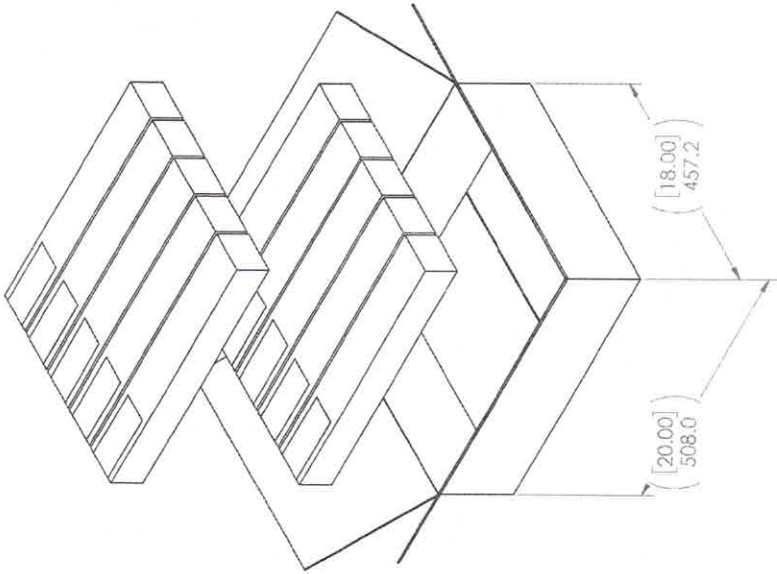


SIZE	DWG. NO.	REV
B	SLQ-81XX	A

SHEET 3 OF 3

NOTES: UNLESS OTHERWISE SPECIFIED

1. GENERAL
 1. ANY CHANGES TO SPECIFICATIONS SHALL BE APPROVED BY SILQ PRIOR TO PRODUCTION.
 2. FINAL PRODUCT INSPECTION PER QP-010.
 3. PRODUCT EXPIRY DATE IS 36 MONTHS FROM THE DATE OF MANUFACTURE.
2. STERILIZATION
 1. PRODUCT IS STERILIZED VIA ETHYLENE OXIDE (ETO) STERILIZATION PER QP-032 CYCLE PWI-7374.
3. TESTING
 1. PYROGEN EVALUATION PER WI-0021 AT ≤ 0.5 EU/ML OR ≤ 20.0 EU/DEVICE

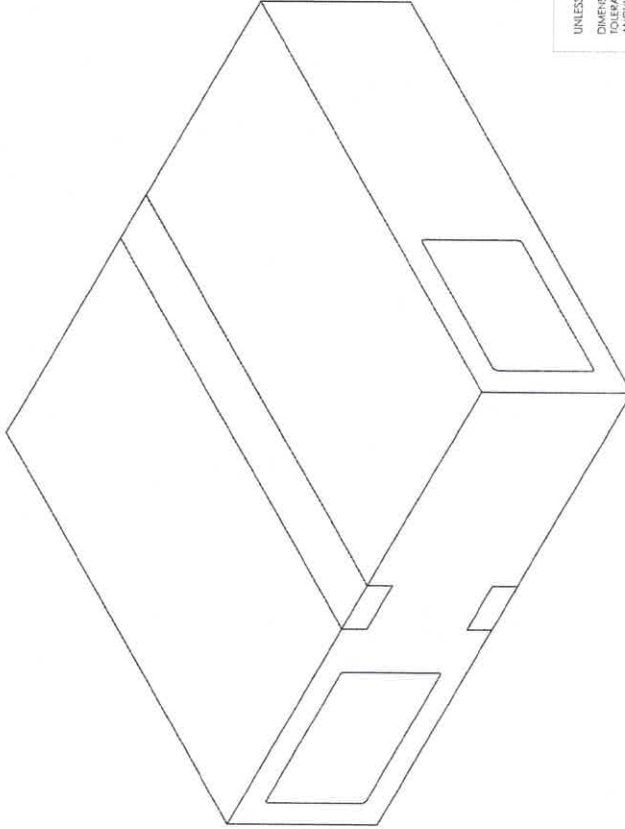

UNLESS OTHERWISE SPECIFIED:
DIMENSIONS ARE IN MILLIMETERS
TOLERANCES:
FRACTIONS: ±0.25
DECIMALS: ±0.25
NO PLACE DECIMAL ±0.3
ONE PLACE DECIMAL ±0.3
TWO PLACE DECIMAL ±0.15
THREE PLACE DECIMAL ±0.05

PROPRIETARY AND CONFIDENTIAL
THE INFORMATION CONTAINED IN THIS
DRAWING IS THE SOLE PROPERTY OF
PATHWAY MEDTECH. NO REPRODUCTION
OR TRANSMISSION OF THIS INFORMATION
IN ANY FORM OR BY ANY MEANS
WITHOUT THE WRITTEN PERMISSION OF
PATHWAY MEDTECH IS PROHIBITED.



TITLE:

2-WAY CLEARTRACT CATHETER,
18 FR X 10 ML, STERILE, 10-PACK

SIZE DWG. NO.

B SLQ-211810SPT

REV

B

DRAWN BY:

B. SANDOVAL

DRAWING NOT TO SCALE

T-018 REV. A

SHEET 1 OF 1

SLQ-211810SPT

In Production

SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Sterile, 10-Pack

Revision »B - ECO-001669 Status: Effective, Status Date: 05/14/2024 03:24:42 PM, Shared

Contains **10** first-level items, **50** line items, **41** unique items, **33** of which are shared.

# ITEM NUMBER	ITEM NAME	CATEGORY	PHASE	△ ↕ ↺ ↻ □ ☑ ☒
▶ 1 SLQ-8102 rev A	SILQ 2-way Cleartract SPT 18 Fr x 10 ml Catheter, Non-Sterile, 10-Pack	 Subassembly	In Prod	3 1 EA
2 WI-0006 rev D	Work Order Process	 Work Instruction	In Prod	2 1 each
3 FM-0039 rev G	Form, Work Order	 Form	In Prod	2 1 each
4 WI-0014 rev A	Line Clearance	 Work Instruction	In Prod	2 1 each
5 MP-SLQ-0003 rev A	Manufacturing Process Flow, SILQ Foley Catheter	 Manufacturing Procedure	In Prod	2 1 EA
6 SLQ-0001 rev A	Sterilization Load Shipping Pallet Configuration	 Subassembly	In Prod	2 1 EA
▶ 7 QP-032 rev F	EO Sterilization Process Control	 SOP Quality System Procedure	In Prod	7 1 each
8 WI-0021 rev C	Work Instruction, Pyrogen (LAL) Evaluation Process	 Work Instruction	In Prod	2 1 each
▶ 9 QP-010 rev F	Final Product Inspection	 SOP Quality System Procedure	In Prod	6 1 each
10 FM-0029 rev A	Form, Certificate of Conformance	 Form	In Prod	2 1 each

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 1 of 7

Production Planning – Document the products identification information with the planned build quantity and any deviation(s) that may apply. Attach a copy of the device's specification/drawing/BOM and any applicable deviation(s) to the work order form prior to transfer to manufacturing.

Part No.	Rev	Lot No./ Serial No.	Expiration Date	Quantity	Deviation No.
SLQ-211810SPT	B	SLQ-81020515241	2027-05-15	405	NA
Comments					
SLQ 2-way Cleartract SPT 18 Fr x 10ml Catheter, Sterile, 10-Pack					
Name			Signature	Date	
Sean Tamarzo				05/15/2024	

Lot No/Serial Number	SLQ-81020515241
----------------------	-----------------



Work Order Form


Doc No.: FM-0039

Rev.: G

Page 2 of 7

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Line Clearance - The individual performing the line clearance certifies that line clearance activities were conducted per WI-0014 prior to starting manufacturing activities. Quality management or designee verifies the line clearance was performed. *Note: The date line clearance occurred is considered the start date of manufacturing.*

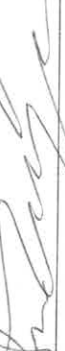Actions		Performed by	Date
1. Clear Area of previous WO Raw Materials.	☒ Completed.		05/15/2024
2. Verify kitted Raw Materials (RM) for: a. RM meet BOM revision level requirements. b. RM have not Expired. c. RM have no physical damage	a. ☒ Completed. b. ☒ Completed. c. ☒ Completed.		05/15/2024
3. Equipment is calibrated and properly functioning.	☒ Completed. ☐ N/A		05/15/2024
<i>Note: Equipment identification and calibration due date will be documented in the work order.</i>			
4. Proper documentation is present and available.	☒ Completed.		05/15/2024
5. Production Associates are trained and qualified to perform the process.	☒ Completed. ☐ N/A		05/15/2024
6. All work surfaces cleaned with IPA.	☒ Completed. ☐ N/A		05/15/2024
7. Labels printed per Label Generation and Reconciliation Form.	☒ Completed. ☐ N/A		05/15/2024
8. Verify validated processes have been set to critical points within specification.	☒ Completed. ☐ N/A		05/15/2024
9. Verify Start-up activities are within specification	☐ Completed. ☒ N/A		05/15/2024
Conclusion: Line clearance performed	☒ Completed		05/15/2024
Quality Management			
Name B. TAMAYO 05/15/2024		Signature B. Tamayo B.T.	Date 05/15/2024

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 3 of 7

Material/ Procedure Traceability (Check appropriate Box):

- ☐ – Identification information for the procedures and materials used correspond to the attached product's BOM and are captured on the attached Production Pick List (ref. P/N FM-0059).
- ☒ - Identification information for the procedures and materials used correspond to the attached product's BOM and are captured in the table below (N/A or delete table if not applicable).

BOM Item No.	Type	Rev	Material Designation	Lot No./Serial No.	Expiration Date	Rev	Material Designation	Lot No./Serial No.	Expiration Date
SLQ-8102	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-81020515241	2027-05-15		☐ Primary ☐ Substitute		
1	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	20211122	N/A		☐ Primary ☒ Substitute		
2	☐ Document ☒ Material	B	☒ Primary ☐ Substitute	2302-01-2	02/23/2024, REF NCR-24-027		☐ Primary ☐ Substitute		
3	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
4	☐ Document ☒ Material	A	☐ Primary ☒ Substitute	A199054	N/A		☐ Primary ☐ Substitute	N/A ST 05/20/2024	
5	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-81020515241	2027-05-15		☐ Primary ☐ Substitute		
6	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	000072 011223	N/A		☐ Primary ☐ Substitute		
7	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	SLQ-81020515241	2027-05-15		☐ Primary ☐ Substitute		
8	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
9	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
10	☐ Document ☒ Material	A	☐ Primary ☒ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
11	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
12	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		

Lot No./Serial Number	SLQ-81020515241
-----------------------	-----------------



Work Order Form

Doc No.: FM-0039

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.

Rev.: G

Page 4 of 7

BOM Item No.	Type	Rev	Material Designation	Lot No./Serial No.	Expiration Date	Rev	Material Designation	Lot No./Serial No.	Expiration Date
13	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
14	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
15	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
16	☒ Document ☐ Material	B	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
SLQ-211810SPT	☐ Document ☒ Material	B	☒ Primary ☐ Substitute	SLQ-81020515241	2027-05-15		☐ Primary ☐ Substitute		
1	☐ Document ☒ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
2	☒ Document ☐ Material	D	☐ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
3	☒ Document ☐ Material	G	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
4	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
5	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
6	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
7	☒ Document ☐ Material	F	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
8	☒ Document ☐ Material	C	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
9	☒ Document ☐ Material	F	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		
10	☒ Document ☐ Material	A	☒ Primary ☐ Substitute	N/A	N/A		☐ Primary ☐ Substitute		

N/A STCS/20120224

Lot No./Serial Number

SLQ-81020515241

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 5 of 7

Training Traceability – All individuals indicated in the following table have been trained to the referenced document(s) which support the manufacturing, inspection and test activities of the product identified in the work order. Leadership representative has certified that the individual(s) have received proper training prior to conducting manufacturing activities.

Document No.	Rev	Print Name	Signature
MP-SLQ-0001	B	EVA TAMAYO	<i>Eva Tamayo</i>
		Aiden Emerick	<i>Aiden Emerick</i>
		N/A	ST 05/20/2024
MP-SLQ-0002	B	EVA TAMAYO	<i>Eva Tamayo</i>
		Aiden Emerick	<i>Aiden Emerick</i>
		N/A	ST 05/20/2024
Leadership Representative			
Name		Signature	Date
Sean Tamayo		<i>Sean Tamayo</i>	05/15/2024

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 6 of 7

Equipment Traceability – Identify any equipment that is used in the manufacturing/processing of the product in the following table.

Equipment Description	Asset No.	Cal Due Date (If applicable)	Equipment Description	Asset No.	Cal Due Date (if applicable)
Accu Seal Impulse Sealer	ES-0099-01	08/09/2024	Catheter Automation Machine	ES-0091-01	N/A

Production Reporting – Document scrapped information, date manufacturing completed, and final quantity manufactured.

BOM Item No.	Defect No. 1		Defect No. 2		Defect No. 3		Defect No. 4		Defect No. 5		Comments
	Code	QTY	Code	QTY	Code	QTY	Code	QTY	Code	QTY	
SLQ-8102	3	7									Inspected/Dispositioned by customer Sent to SILQ for testing.
SLQ-8102	10	18									
Device Assembly											

Defect Codes: 1=Damage, 2=Dimensional, 3=Visual, 4=Voids/Shorts, 5=Discoloration, 6=Contamination, 7=Operator Error, 8=Occlusion, 9=Leak, 10 =Miscellaneous (Indicate in Comment Section)

Manufacturing Completion Date	05/20/2024	Quantity MFG	38 packs of 10
Name	EVA TAMAYO	Signature	
		Date	05/20/2024

	Work Order Form		Doc No.: FM-0039
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: G
			Page 7 of 7

Documentation Checklist – If applicable, the following documentation shall be included upon completion of the manufacturing build. Document packet shall be organized in the order indicated below. Ensure that good documentation practices are applied per QP-003 "Quality Records". Check box(es) ☐ that apply to each requirement. The individual who signs has confirmed that all requirements have been met.

1. Work Order (WO) Record	☒ Present and Completed	
2. Product's Specification/Drawing & BOM (Arena PLM)	☒ Present	
3. Deviations	☐ Present	☒ N/A
4. Product Labeling & Generation/ Reconciliation Record	☒ Present, Completed and Results Meet Specification	☐ N/A
5. Non-Conforming Reports (NCR)	☒ Present and Completed	☐ N/A
6. Inspection Records	☒ Present, Completed and Results Meet Specification	
7. Critical Process Parameter Record	☐ Present, Completed and Results Meet Specification	☒ N/A
Name	Signature	Date
Melissa Garcia	Melissa Garcia	07/29/2024

Release Authorization for Further Processing (e.g., terminal sterilization) – The individual signing below is certifying that they independently reviewed the Device History Record (DHR) and confirmed that all required documentation is included, records are complete, good documentation practices were followed, and specifications have been met (if not NCR is included)

Name	Signature	Date
Melissa Garcia	Melissa Garcia	07/29/2024

	Form, Run Sheet for SILQ Automation Machine		Doc No.: FM-SQL-0001
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as an uncontrolled "For Reference Only", unless accompanied with a controlled stamp. Users of this document must verify that they are using the most current revision of the controlled document before performing work.		Rev.: A
			Page 1 of 1

Part Number	Revision	Description	Machine Asset No.
SLQ-8102	A	SILQ 2-Way ClearTract SPT 18Fr x 10ml Catheter, Non-Sterile, 10-Pack	ES-0091-01


Equipment Setup / Startup		Adjustments during Production		
Setting [Units]	Value	Value	Value	Value
Setup / Adjustment Date	05/15/2024			
Setup / Adjustment Time	10:10 AM			
Fill Time [seconds]	28.0			
Confirm Cylinder Actuation All 4 pairs [P/F]	P			
UV Cure Time [seconds]	120.0			
Drain Time [seconds]	30.0			
Rinse Time [seconds]	6.3			
Final Drain [seconds]	30.0			
Number of Drain / Rinse Cycles [Count]	3.0			
Confirm Cylinder Actuation All 4 pairs [P/F]	P			
Size Setting	18Fr			
Performed By Name	Sean Tamayo			
Performed By Signature				
Verified By Name	B. Tamayo			
Verified By Signature				

Engineering / Lead Supervisor

QA

VIA ST 05/20/2024

Pathway Medtech, LLC.	Nonconforming Record	NCR No.:	24-027
-----------------------	-----------------------------	----------	---------------

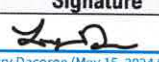Source:	☐ Receiving ☒ Start-Up Verification (QC) ☐ In-Process (QC) ☐ In-Process (Production)				Date:	05/14/2024
	☐ Final Product ☐ RMA (No.) ☐ Miscellaneous:					
Part No.:	SLQ-3000	Revision:	B	Lot No. / SN:	2302-01-2, 2302-01-1	
Description	Solution, SILQ Polymer		Supplier	SILQ Industries	Inspector	Larry Dacoron
NONCONFORMITY(S) DETAIL(S)						
#	Description	QTY Rejected	QTY Inspected	AQL	Specification	
					Procedure	Revision
1	Expired solution (02/22/2024 and 02/23/2024)	2	2	N/A	SLQ-3000	B
CONTAINMENT ACTIVITY						
Contained QTY		Comments				
2 barrels		N/A				
INITIAL DISPOSITION						
#	Disposition					
1	☐ Increase Sampling ☐ Sort (See instructions below) ☐ Rework (See instructions below) ☒ Hold/Quarantine ☐ N/A					
Note: Applies to Rework Dispositions Only		Risk Analysis Required?:		Risk Analysis:		
		☐ Yes ☒ No		N/A		
Additional Instructions:		N/A				
APPROVALS						
Signature		Date	Signature		Date	
☒ QA	 Larry Dacoron (May 15, 2024 08:30 PDT)	05/15/24	☐ MFG			
☒ CUSTOMER	 Ethan Rao (May 15, 2024 08:48 PDT)	05/15/24	☐			
SORT / REWORK COMPLETED BY						☐ N/A
Signature		Date	Signature		Date	
N/A		N/A	N/A		N/A	
RE-INSPECTION						
Description	QTY Rejected	QTY Inspected	AQL	Specification		
				Procedure	Revision	
N/A	N/A	N/A	N/A	N/A	N/A	
Inspected By:						Date
Signature	N/A					N/A
FINAL DISPOSITION						
#	Disposition					QTY
1	Meets Specification after Initial Disposition					0
	Use As Is (UAI)					2
	Return to Vendor (RTV)					0
	Scrap					0
					Total:	2
Note: Applies to UAI Dispositions Only		Risk Analysis Required?:		Risk Analysis:		
		☒ Yes ☐ No		Expired solutions may undergo chemical changes over time, leading to potential degradation of the product. Using expired solutions may yield inaccurate or unreliable data, undermining the validity of the study results.		
Additional Instructions:		The suspension will be used in the manufacturing process of the product. Any product that will use the expired solution will be placed in quarantine after sterilization. After manufacturing of the product, the expired suspension will be subject to client's full chemical analysis QC process analogous to what was done upon initial validation of the suspension. An additional 50 16Fr catheters were manufactured using the same expired suspension, all of which will be tested after sterilization to further validate the suspension's acceptability. These 18Fr catheters will be subject to the standard sampling plan & QC process. If the results of the chemical analysis of the expired suspension, the testing of the 50 16Fr catheters, and the standard testing of this product are all acceptable, the product will be released from quarantine.				
CAPA No. ☒ N/A		N/A				
APPROVALS						
Signature		Date	Signature		Date	
☒ QA	 Larry Dacoron (May 15, 2024 08:30 PDT)	05/15/24	☒ MFG	 Sean Tamayo (May 15, 2024 08:35 PDT)	05/15/24	

Pathway Medtech, LLC.		Nonconforming Record		NCR No.:	24-027
☒ CUSTOMER	<u>Ethan Rao</u> Ethan Rao (May 15, 2024 08:48 PDT)	05/15/24	☐ N/A	N/A	N/A
Customer Comments / Special Disposition Instructions					
MATERIAL DESTRUCTION (APPLIES TO SCRAP ONLY)					☒ N/A
Performed by / Signature		Date	QTY Destroyed		
N/A		N/A	N/A		

NCR 24-027 - Expired solution

Final Audit Report

2024-05-15

Created:	2024-05-15
By:	Larry Dacoron (ldacoron@pathwaynpi.com)
Status:	Signed
Transaction ID:	CBJCHBCAABAAJBfWGX-TJUPOJiKfnYn3DZGHllz0g4N4

"NCR 24-027 - Expired solution" History

-  Document created by Larry Dacoron (ldacoron@pathwaynpi.com)
2024-05-15 - 3:29:41 PM GMT
-  Document emailed to Ethan (ethanr@silq.tech) for signature
2024-05-15 - 3:29:45 PM GMT
-  Document emailed to Sean (stamayo@pathwaynpi.com) for signature
2024-05-15 - 3:29:45 PM GMT
-  Document emailed to Larry Dacoron (ldacoron@pathwaynpi.com) for signature
2024-05-15 - 3:29:45 PM GMT
-  Email viewed by Larry Dacoron (ldacoron@pathwaynpi.com)
2024-05-15 - 3:30:12 PM GMT
-  Document e-signed by Larry Dacoron (ldacoron@pathwaynpi.com)
Signature Date: 2024-05-15 - 3:30:28 PM GMT - Time Source: server
-  Email viewed by Sean (stamayo@pathwaynpi.com)
2024-05-15 - 3:35:08 PM GMT
-  Signer Sean (stamayo@pathwaynpi.com) entered name at signing as Sean Tamayo
2024-05-15 - 3:35:40 PM GMT
-  Document e-signed by Sean Tamayo (stamayo@pathwaynpi.com)
Signature Date: 2024-05-15 - 3:35:42 PM GMT - Time Source: server
-  Email viewed by Ethan (ethanr@silq.tech)
2024-05-15 - 3:47:49 PM GMT
-  Signer Ethan (ethanr@silq.tech) entered name at signing as Ethan Rao
2024-05-15 - 3:48:19 PM GMT



Document e-signed by Ethan Rao (ethanr@silq.tech)

Signature Date: 2024-05-15 - 3:48:21 PM GMT - Time Source: server



Agreement completed.

2024-05-15 - 3:48:21 PM GMT



Adobe Acrobat Sign

	DOCUMENT TITLE: Quality Control Procedure, Crystal Violet Absorbance, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ003	Effective Date: 06Jun2024



Crystal Violet Absorbance, Silq Surface Treatment

Quality Control Procedure

	DOCUMENT TITLE: Quality Control Procedure, Crystal Violet Absorbance, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ003	Effective Date: 06Jun2024

1 Purpose

- 1.1 The purpose of this document is to detail the quality control procedure used to verify the presence of Silq's Surface Treatment coating on catheters using a crystal violet dye.

2 Scope

- 2.1 This procedure applies to catheters treated with the Silq Surface Treatment process.

3 Documents

- 3.1 Reference Documents
 - 3.1.1 QM.SLQ049 Workstation Practices SOP
 - 3.1.2 QC-E.SLQ015 Quality Control Qualification Procedure, Jenway 7205 Spectrophotometer

4 Definitions

- 4.1 Abbreviations / Acronyms
 - 4.1.1 **cm:** Centimeter
 - 4.1.2 **DI:** Deionized
 - 4.1.3 **in:** Inches
 - 4.1.4 **mg:** Milligram
 - 4.1.5 **mL:** Milliliter
 - 4.1.6 **nm:** Nanometer

5 Responsibilities

- 5.1 Pathway MedTech (Pathway) is responsible for:
 - 5.1.1 Manufacturing, including final packaging and sterilization, of catheters in accordance with established manufacturing procedures
- 5.2 Silq is responsible for:
 - 5.2.1 Performing crystal violet absorbance test in accordance with this procedure
 - 5.2.2 Documenting the results of the crystal violet absorption test

6 Safety Information

	DOCUMENT TITLE: Quality Control Procedure, Crystal Violet Absorbance, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ003	Effective Date: 06Jun2024

- 6.1 Gloves (non-latex) and safety eye wear are to be worn at all times during the execution of this procedure.

7 Equipment / Materials / Supplies

- 7.1 Spectrophotometer (E.SLQ015)
- 7.2 UV Cuvette (S.SLQ011)
- 7.3 Crystal Violet Dye Solution (S.SLQ008)
- 7.4 DI Water
- 7.5 Acetic Acid (S.SLQ009)
- 7.6 Centrifuge tubes, 50 mL, Polypropylene, Flat Cap (S.SLQ002)
- 7.7 Glass test tubes (S.SLQ010)
- 7.8 X-Acto Knife or Equivalent
- 7.9 Inch Ruler
- 7.10 Timer
- 7.11 Non-Latex Gloves
- 7.12 Metal Tray

8 Line Clearance and Setup

- 8.1 Perform line clearance of the work area in accordance with Workstation Practices QM.SLQ049.
- 8.2 Power up the spectrophotometer and allow to warm up for at least 15 minutes. Set the wavelength to 585nm.
- 8.3 The spectrophotometer shall be qualified daily before use. Once the spectrophotometer is warmed up, qualify the spectrophotometer in accordance with QC-E.SLQ015.

9 Procedure

- 9.1 Secure 3 glass tubes (S.SLQ010).
- 9.2 Fill one tube with a well stirred crystal violet solution (S.SLQ008) to the demarcated line.
- 9.3 Fill the other two tubes with DI Water to the demarcated line.
- 9.4 Take a dried, surface treated catheter sample and immerse it in the tube containing the crystal violet solution. Allow the sample to soak for 1 minute.

	DOCUMENT TITLE: Quality Control Procedure, Crystal Violet Absorbance, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ003	Effective Date: 06Jun2024

- 9.5 Slowly remove the catheter sample from the tube allowing for all crystal violet solution to drain. Immediately dip the catheter sample into a tube containing DI water and again slowly remove the catheter sample allowing for all residual liquid to drain. Repeat with the 2nd tube of DI water.
- 9.6 Place the catheter sample on the clean, dry metal tray.
- 9.7 Using an X-Acto knife (or equivalent), cut off the tip of the catheter sample just above the islet at the bottom.
- 9.8 Using a ruler, measure 10" from the end of the catheter sample and make another cut.
- 9.9 Cut the 10" catheter sample into five 2" segments. Use tweezers to place all segments into a 50mL centrifuge tube (S.SLQ002) containing 25 mL of the 2:1 DI water: Acetic acid solution.
- 9.10 Shake the 50mL centrifuge tube for 30 seconds until all dye has been removed from the catheter segments.
- 9.11 Pour 2 mL of the solution into a clean UV cuvette (S.SLQ011).
- 9.12 Separately, prepare a cuvette containing 2 mL of the 2:1 DI water:acetic acid solution. Use this cuvette to blank the spectrophotometer at 585nm.
- 9.13 Measure the absorbance value of the sample cuvette and document results in Attachment 1.

10 Acceptance Criteria

- 10.1 The absorbance value shall be a minimum of 0.400.

11 Attachments

- 11.1 Attachment 1: Crystal Violet Absorbance Datasheet

	DOCUMENT TITLE: Quality Control Procedure, Crystal Violet Absorbance, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ003	Effective Date: 06Jun2024

Attachment 1: Crystal Violet Absorbance Datasheet

Tester (Print name & Sign): Ethan Rao *Ethan Rao*

Date of Test: 18 July 2024

General Product Information

Part Number:	SLQ-8104	Part Description:	SILQ 2-Way Cleartract SPT 16 Fr x 10mL Catheter
Part Revision:	B	Lot #:	SLQ-80040411241E

Supplies Traceability

Part Number	Description	Lot Number	Expiration Date
S.SLQ008	Crystal violet dye	SHBR0031	01 July 2027
S.SLQ009	Acetic acid	231511	01 May 2028

Equipment Traceability

Part Number	Description	Equipment ID	Last qualified on*
E.SLQ015	Jenway 7205 Spectrophotometer	E.SLQ015	18 July 2024

* Attach QC Qualification Datasheet

Summary of Results

Catheter Sample #	Result	Acceptance criteria	Meets Criteria?
1	0.458	$\geq .400$	☒ Yes ☐ No
2	0.467	$\geq .400$	☒ Yes ☐ No
3	0.415	$\geq .400$	☒ Yes ☐ No
4	0.493	$\geq .400$	☒ Yes ☐ No
5	0.416	$\geq .400$	☒ Yes ☐ No
6	0.429	$\geq .400$	☒ Yes ☐ No
7	0.444	$\geq .400$	☒ Yes ☐ No
8	0.433	$\geq .400$	☒ Yes ☐ No
9	0.499	$\geq .400$	☒ Yes ☐ No
10	0.481	$\geq .400$	☒ Yes ☐ No
11	0.496	$\geq .400$	☒ Yes ☐ No
12	0.452	$\geq .400$	☒ Yes ☐ No
13	0.469	$\geq .400$	☒ Yes ☐ No
14	0.412	$\geq .400$	☒ Yes ☐ No
15	0.449	$\geq .400$	☒ Yes ☐ No
16	0.487	$\geq .400$	☒ Yes ☐ No
17	0.467	$\geq .400$	☒ Yes ☐ No
18	0.493	$\geq .400$	☒ Yes ☐ No
19	0.514	$\geq .400$	☒ Yes ☐ No

	DOCUMENT TITLE: Quality Control Procedure, Crystal Violet Absorbance, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ003	Effective Date: 06Jun2024

Catheter Sample #	Result	Acceptance criteria	Meets Criteria?
20	0.465	$\geq .400$	☒ Yes ☐ No
21	0.438	$\geq .400$	☒ Yes ☐ No
22	0.482	$\geq .400$	☒ Yes ☐ No

OVERALL TEST RESULTS (circle one): PASS / FAIL

Reviewed By (Sign and Date): Brian McVerry 19 July 2024



DOCUMENT TITLE:
Test Protocol, Suspension Shelf Life Evaluation,
Silq Surface Treatment

Revision:
A

DOCUMENT #:
VV.SLQ028

Effective Date:
06Jun2024



Suspension Shelf Life Evaluation, Silq Surface Treatment

Test Protocol

	DOCUMENT TITLE: Test Protocol, Suspension Shelf Life Evaluation, Silq Surface Treatment	Revision: A
	DOCUMENT #: VV.SLQ028	Effective Date: 06Jun2024

1 Purpose

- 1.1 The purpose of this protocol is to verify the Silq Suspension, part number C.SLQ001, meets specified performance requirements for a defined shelf life (e.g., 2 years).
- 1.2 Performance testing consists of:
 - 1.2.1 Verifying the quality attributes of the Suspension (i.e., appearance, residue on evaporation, pH, ¹H NMR, and ¹⁹F NMR)
 - 1.2.2 Confirming the presence and adherence of the Suspension to the catheter

2 Scope

- 2.1 This test protocol applies to the Silq Suspension, part number C.SLQ001
- 2.2 The Silq Surface Treatment process may be applied to a variety of polymeric substrates. A representative catheter (16 Fr Foley catheter) will be chosen to evaluate the performance of aged Suspension solution.

3 Documents

- 3.1 Reference Documents
 - 3.1.1 QC-C.SLQ001 Quality Control Procedure, Suspension Analysis
 - 3.1.2 SP-C.SLQ001 Specification, Suspension
 - 3.1.3 OM.SLQ003 Quality Control Procedure, Crystal Violet Absorbance
 - 3.1.4 OM.SLQ004 Quality Control Procedure, Coating Adherence Verification

4 Background

- 4.1 The C.SLQ001 Suspension solution is a water-based coating solution for polymeric and elastomeric surfaces. The solution contains a three-component methacrylate copolymer, polymerized via free radical polymerization from polysulfobetaine methacrylate, perfluorophenylazide methacrylate, and methacrylic acid. The material is dispersed at 1% w/v% in high purity water.

5 Training

- 5.1 Applicable personnel will be trained to this protocol (Training Record Form in Appendix 1).

6 Test Configuration

- 6.1 Aged Suspension will be tested to confirm its attributes meets specification.

**DOCUMENT TITLE:**Test Protocol, Suspension Shelf Life Evaluation,
Silq Surface Treatment**Revision:**

A

DOCUMENT #:

VV.SLQ028

Effective Date:

06Jun2024

- 6.2 Finished 16 Fr Foley catheters built to the approved manufacturing procedure by Pathway will be used for the evaluation of the aged Suspension on catheters. 16 Fr Foley catheters are representative samples of the ClearTract Foley catheter product family.
- 6.3 Manufacturing lot numbers will be recorded for traceability.
- 6.4 All test samples of the finished 16 Fr Foley catheters will be packaged and subjected to the validated ethylene oxide sterilization cycle.
- 6.5 Sample size
- 6.5.1 Shelf life evaluation will consist of the following tests with the listed sample sizes:

Test Group	Test Description	Sample Size
Quality attributes (Suspension)	Verify the attributes of the Suspension (i.e., appearance, residue on evaporation, pH, ¹ H NMR, and ¹⁹ F NMR)	100%
Crystal Violet Absorbance (Finished catheter)	Verify the presence of Suspension using a crystal violet dye solution	A sufficient number of test samples to achieve 90%/90% confidence/reliability interval using variable data analysis
Coating Adherence Verification (Finished catheter)	Verify the Suspension remains adhered to the catheter	A sufficient number of test samples to achieve 90%/90% confidence/reliability interval using variable data analysis

- 6.5.2 For quality attribute testing of the Suspension, each lot of Suspension will be tested for a sample size of 100%.
- 6.5.3 For crystal absorbance and coating adherence verification testing of finished catheters, a sufficient number of test samples will be selected to perform variable data analysis at a 90/90% confidence/reliability interval. These statistical intervals were determined using a risk-based approach in compliance with the statistical techniques work instruction QM.SLQ011.

7 Test Equipment and Materials / Supplies

- 7.1 See the equipment, materials, and supplies required for the execution of this test protocol defined in documents QC-C.SLQ001 Quality Control Procedure, Suspension Analysis, OM.SLQ003 Quality Control Procedure, Crystal Violet Absorbance, and OM.SLQ004 Quality Control Procedure, Coating Adherence Verification.

	DOCUMENT TITLE: Test Protocol, Suspension Shelf Life Evaluation, Silq Surface Treatment	Revision: A
	DOCUMENT #: VV.SLQ028	Effective Date: 06Jun2024

8 Test Procedure: Suspension

- 8.1.1 Prepare and test the aged C.SLQ001 Suspension in accordance with QC-C.SLQ001 Quality Control Procedure, Suspension Analysis.
- 8.1.2 Document the test results on the Suspension Quality Control Report template (TMP1-QC-C.SLQ001).
- 8.1.3 Record the summary of the quality control test in Appendix 2. Attach a copy of the Suspension Quality Control Report to the final test report.

9 Test Procedure: Finished Catheter

- 9.1.1 Using the aged Suspension, Pathway shall manufacture 16 Fr ClearTract Foley catheter test samples in accordance with established manufacturing procedures. Following sterilization, the test samples will be delivered to Silq for testing.
- 9.1.2 The presence of Suspension on the test samples shall be verified in accordance with OM.SLQ003 Quality Control Procedure, Crystal Violet Absorbance. Record the results of the crystal violet absorbance test in Appendix 2.
- 9.1.3 Adherence of the Suspension on the test samples shall be verified in accordance with OM.SLQ004 Quality Control Procedure, Coating Adherence Verification. Record the results of the coating adherence verification test in Appendix 2.

10 Analysis of Data

- 10.1 Testing will collect variable data. The sample size corresponds to a 90% confidence / 90% reliability interval.
- 10.2 Variable data analysis techniques using a tolerance interval ("k" factor) will be used to determine if the quantitative results fall within the required confidence and reliability intervals.

11 Acceptance Criteria

- 11.1 Suspension
 - 11.1.1 Appearance – Color: Opaque white
 - 11.1.2 Appearance – Form: Inviscid liquid
 - 11.1.3 Residue on Evaporation: 0.5% – 1.5%
 - 11.1.4 pH: 5 – 8
 - 11.1.5 ¹H NMR: A successful result is one that includes broad peaks within 0.2ppm of 0.99, 2.2, 2.9, 3.1, 3.5, 3.7, 4.4, and 4.75 ppm indicating

	DOCUMENT TITLE: Test Protocol, Suspension Shelf Life Evaluation, Silq Surface Treatment	Revision: A
	DOCUMENT #: VV.SLQ028	Effective Date: 06Jun2024

successful incorporation of the zwitterion moiety into the polymer. The peak at 4.75 ppm is the 'solvent peak' resulting from the deuterium oxide solvent. All peaks should be easily identifiable and have an amplitude of at least 5% of the solvent peak.

11.1.6 ¹⁹F NMR: A successful result is one that includes broad peaks within 5ppm of -140 and -150ppm indicating successful incorporation of the PFPA moiety into the polymer. These peaks should be broad with a small amplitude but still easily identifiable.

11.2 Crystal Violet Absorbance: The absorbance value shall be a minimum of 0.400.

11.3 Coating Adherence Verification: The absorbance value shall be a less than or equal to 0.02.

12 Attachments

12.1 Attachment 1: Protocol Training Record Form

12.2 Attachment 2: Shelf Life Evaluation Data Sheet

	DOCUMENT TITLE: Test Protocol, Suspension Shelf Life Evaluation, Silq Surface Treatment	Revision: A
	DOCUMENT #: VV.SLQ028	Effective Date: 06Jun2024

Attachment 1 – Protocol Training Record Form

By signing this form, I understand the requirements and responsibilities to executing this protocol.

Trainee Name	Initials	Signature	Date	Department / Company
Ethan Rao	ER	Ethan Rao	07 June 2024	Silq Technologies

	DOCUMENT TITLE: Test Protocol, Suspension Shelf Life Evaluation, Silq Surface Treatment	Revision: A
	DOCUMENT #: VV.SLQ028	Effective Date: 06Jun2024

Attachment 2 – Suspension Shelf Life Evaluation Data Sheet

Aged Suspension Information

Part Number:	C.SLQ001	Manufacture & Expiration Date:	Mfg date: 25 Aug 2021 Exp date: 25 Aug 2022
Part Revision:	A	Batch / Lot #:	HYX-825021

Suspension Quality Control Test Results

Attribute	Result	Meets Criteria?
Appearance (color)	See attached test report*	☒ PASS ☐ FAIL
Appearance (form)	See attached test report*	☒ PASS ☐ FAIL
Residue on Evaporation	See attached test report*	☒ PASS ☐ FAIL
pH	See attached test report*	☒ PASS ☐ FAIL
¹ H NMR	See attached test report*	☒ PASS ☐ FAIL
¹⁹ F NMR	See attached test report*	☒ PASS ☐ FAIL

* Attach the Suspension Quality Control test report to this data sheet

Finished Catheter Testing

Attribute	Result	Meets Criteria?
Crystal violet absorbance	See attached test report**	☒ PASS ☐ FAIL
Coating adherence	See attached test report***	☒ PASS ☐ FAIL

** Attach the Crystal Violet Absorbance Quality Control test results to this data sheet

*** Attach the Coating Adherence Quality Control test results to this data sheet

OVERALL TEST RESULTS (circle one): **PASS** / FAIL

Completed By (Print name, Sign and Date): Ethan Rao Ethan Rao 18 July 2024

Reviewed By (Print name, Sign and Date): Brian McVerry Brian McVerry 19 July 2024

	DOCUMENT TITLE: Quality Control Procedure, Coating Adherence Verification, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ004	Effective Date: 06Jun2024



Coating Adherence Verification, Silq Surface Treatment

Quality Control Procedure

	DOCUMENT TITLE: Quality Control Procedure, Coating Adherence Verification, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ004	Effective Date: 06Jun2024

1 Purpose

- 1.1 The purpose of this document is to detail the quality control procedure used to verify the coating using the Silq Surface Treatment process remains adhered to the catheter. The coated catheter is rinsed with deionized water and the rinse is examined for the presence of coating using spectroscopy.

2 Scope

- 2.1 This procedure applies to catheters treated with the Silq Surface Treatment process.

3 Documents

- 3.1 Reference Documents
 - 3.1.1 QM.SLQ049 Workstation Practices SOP
 - 3.1.2 QC-E.SLQ015 Quality Control Qualification Procedure, Jenway 7205 Spectrophotometer

4 Definitions

- 4.1 Abbreviations / Acronyms
 - 4.1.1 **cm:** Centimeter
 - 4.1.2 **DI:** Deionized
 - 4.1.3 **in:** Inches
 - 4.1.4 **mL:** Milliliter
 - 4.1.5 **nm:** Nanometer

5 Responsibilities

- 5.1 Pathway MedTech (Pathway) is responsible for:
 - 5.1.1 Manufacturing, including final packaging and sterilization, of catheters in accordance with established manufacturing procedures
- 5.2 Silq is responsible for:
 - 5.2.1 Performing crystal violet absorbance test in accordance with this procedure
 - 5.2.2 Documenting the results of the crystal violet absorption test

6 Safety Information

- 6.1 Gloves (non-latex) and safety eye wear are to be worn at all times during the execution of this procedure.

	DOCUMENT TITLE: Quality Control Procedure, Coating Adherence Verification, Silq Surface Treatment	Revision: A
	DOCUMENT #: OM.SLQ004	Effective Date: 06Jun2024

7 Equipment / Materials / Supplies

- 7.1 Spectrophotometer (E.SLQ015)
- 7.2 UV Cuvette (S.SLQ011)
- 7.3 DI Water
- 7.4 Glass test tubes (S.SLQ010)
- 7.5 Non-Latex Gloves

8 Line Clearance and Setup

- 8.1 Perform line clearance of the work area in accordance with Workstation Practices QM.SLQ049.
- 8.2 Power up the spectrophotometer and allow to warm up for at least 15 minutes. Set the wavelength to 585nm.
- 8.3 The spectrophotometer shall be qualified daily before use. Once the spectrophotometer is warmed up, qualify the spectrophotometer in accordance with QC-E.SLQ015.

9 Procedure

- 9.1 Secure a glass tube (S.SLQ010) and fill with DI water to the demarcated line.
- 9.2 Take a dried, surface treated catheter and immerse it in the tube. Gently lift the catheter out of the tube and solution and then resubmerge it for a total of five times.
- 9.3 Remove the catheter from the tube.
- 9.4 Using a gloved finger to plug the top of the tube, invert the tube at least 3 times vigorously to ensure that any materials present in solution will be homogenously distributed.
- 9.5 Pour out 2mL of the solution into a clean cuvette.
- 9.6 Separately, pour 2mL of pure DI water into a UV cuvette (S.SLQ011) and use this cuvette to blank the spectrophotometer at 440 nm.
- 9.7 Measure the absorbance value of the sample cuvette and document the results in Attachment 1.

10 Acceptance Criteria

- 10.1 The absorbance value shall be a less than or equal to 0.02.

11 Attachments

- 11.1 Attachment 1: Coating Adherence Verification Datasheet



DOCUMENT TITLE:
Quality Control Procedure, Coating Adherence
Verification, Silq Surface Treatment

DOCUMENT #:
OM.SLQ004

Revision:
A

Effective Date:
06Jun2024

Attachment 1: Coating Adherence Verification Datasheet

Tester (Print name & Sign): Ethan Rao

Ethan Rao

Date of Test: 18 July 2024

General Product Information

Part Number:	SLQ-8104	Part Description:	SILQ 2-Way Cleartract SPT 16 Fr x 10mL Catheter
Part Revision:	B	Lot #:	SLQ-80040411241E

Equipment Traceability

Part Number	Description	Equipment ID	Last qualified on*
E.SLQ015	Jenway 7205 Spectrophotometer	E.SLQ015	18 July 2024

* Attach QC Qualification Datasheet

Summary of Results

Catheter Sample #	Result	Acceptance criteria	Meets Criteria?
1	0.002	$\leq .02$	☒ Yes ☐ No
2	0.006	$\leq .02$	☒ Yes ☐ No
3	0.000	$\leq .02$	☒ Yes ☐ No
4	0.003	$\leq .02$	☒ Yes ☐ No
5	0.009	$\leq .02$	☒ Yes ☐ No
6	0.002	$\leq .02$	☒ Yes ☐ No
7	0.006	$\leq .02$	☒ Yes ☐ No
8	0.001	$\leq .02$	☒ Yes ☐ No
9	0.101	$\leq .02$	☒ Yes ☐ No
10	0.010	$\leq .02$	☒ Yes ☐ No
11	0.006	$\leq .02$	☒ Yes ☐ No
12	0.005	$\leq .02$	☒ Yes ☐ No
13	0.003	$\leq .02$	☒ Yes ☐ No
14	0.000	$\leq .02$	☒ Yes ☐ No
15	0.004	$\leq .02$	☒ Yes ☐ No
16	0.009	$\leq .02$	☒ Yes ☐ No
17	0.001	$\leq .02$	☒ Yes ☐ No
18	0.008	$\leq .02$	☒ Yes ☐ No
19	0.009	$\leq .02$	☒ Yes ☐ No
20	0.010	$\leq .02$	☒ Yes ☐ No
21	0.012	$\leq .02$	☒ Yes ☐ No
22	0.008	$\leq .02$	☒ Yes ☐ No

OVERALL TEST RESULTS (circle one): **PASS** / FAIL

Reviewed By (Sign and Date): Brian McVerry 19 July 2024

**LGGS, Inc. Technical Report**

SPONSOR: Pathway LLC
8779 Cottonwood Ave Suite 105
Santee, CA 92071

Lab No. 24-3127-00
P.O.No. P1839
Lot No. SLQ-81020515241
Date Received: 07-17-24

Analysis: Limulus Amebocyte Lysate (LAL) Kinetic Turbidimetric Assay

Test Article: PN: SLQ-211810SPT

Condition of Test Article at Receipt: Good

Procedure: Three (3) products were extracted following procedures outlined in test specification number, BEKT-PTW-001-00. One-tenth ml of the test and/or control solution and one-tenth of lysate solution, reconstituted per suppliers directions, were placed in sterile microplate wells and incubated at 37°C for a minimum of one hour in the Kinetic QCL Reader. The reaction time with the LAL/substrate was determined turbidimetrically and compared to a standard curve to determine endotoxin concentration. Positive Product Control (PPC) was simultaneously prepared to demonstrate that the test article did not significantly interfere with the LAL/substrate reaction.

Results:	Test Article Extract Dilution:	<u>2.0</u>
	PPC Testing:	<u>Acceptable</u>
	Test Article Extract:	<u>0.20 EU/ml; 8.0 EU/device</u>
	Date Tested:	<u>07-17-24</u>

Key to Results: Type of Product


Current FDA Requirements

1) Medical Device

1) Less than or equal to 0.5 EU/ml
Less than or equal to 20.0 EU/device

Comments: Test article as submitted conforms to FDA and USP requirements for end-product release of medical devices.

Deviations: None

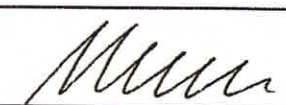Completed: 7/18/2024 Tech RR / GS Approved by: 
Technical Manager (or deputy)

The conclusions and data herein apply only to the items tested and are not indicative of the quality or condition of apparently identical product. The report is submitted for the exclusive use of the sponsor. No use or mention of LGGS, Inc. in connection with any form of advertising or public announcement or reproduction of the report except in full may be made without written approval of LGGS, Inc.

**LGGS, Inc. Technical Report**

SPONSOR: Pathway LLC
8779 Cottonwood Ave Suite 105
Santee, CA 92071

Lab No. 24-3128-00
P.O.No. P1839
Lot No. Not Supplied
Date Received: 07-17-24

Analysis: Limulus Amebocyte Lysate (LAL) Kinetic Turbidimetric Assay

Test Article: Ningbo Luke 16Fr 5-15ml Catheter Samples

Condition of Test Article at Receipt: Good

Procedure: Ten (10) products were extracted following procedures outlined in test specification number, BEKT-PTW-017-00. One-tenth ml of the test and/or control solution and one-tenth of lysate solution, reconstituted per suppliers directions, were placed in sterile microplate wells and incubated at 37°C for a minimum of one hour in the Kinetic QCL Reader. The reaction time with the LAL/substrate was determined turbidimetrically and compared to a standard curve to determine endotoxin concentration. Positive Product Control (PPC) was simultaneously prepared to demonstrate that the test article did not significantly interfere with the LAL/substrate reaction.

Results:	Test Article Extract Dilution:	<u>2.0</u>
	PPC Testing:	<u>Acceptable</u>
	Test Article Extract:	<u><0.02 EU/ml; <0.8 EU/device</u>
	Date Tested:	<u>07-17-24</u>

Key to Results: Type of Product


Current FDA Requirements

1) Medical Device

1) Less than or equal to 0.5 EU/ml
Less than or equal to 20.0 EU/device

Comments: Test article as submitted conforms to FDA and USP requirements for end-product release of medical devices.

Deviations: None

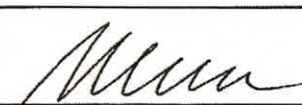Completed: 7/18/2024 Tech RR / GS Approved by: 
Technical Manager (or deputy)

The conclusions and data herein apply only to the items tested and are not indicative of the quality or condition of apparently identical product. The report is submitted for the exclusive use of the sponsor. No use or mention of LGGS, Inc. in connection with any form of advertising or public announcement or reproduction of the report except in full may be made without written approval of LGGS, Inc.

Date Inspected	Part No.	Rev.	Description						
05/17/2024	SLQ-8102	B	SILQ 2-Way Cleartract SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack						
Lot No.	Lot QTY	ES-0059-01 Emplex Model MPS6300-M Continuous Band Sealer - Validated Range: Temperature = 142 °C - 144 °C, Speed = 175 in/min., Pressure = 89 psi ES-0099-01 ACCU-SEAL Model 5300-15-B, Impulse Sealer - Validated Range: Temperature = 270 °F - 280 °F, Dwell = 4 sec., Pressure = 55 psi							
SLQ-81020515241	405								
Specification	TM-0003								
Revision	C								
PRODUCT PREP FOR TEST									
* Seal one sample pouch P/N SLQ-5001 at the start, middle and end of each shift.									
* Cut two 1.00 inch wide x 2.00 inch long strips from the manufacturer sealed end of the pouch.									
* Execute peel strength test per TM-0003 and record results.									
Specification	≥1.0	≥1.0	142-144 °C	175 in/min	89 psi	270 - 280 °F	4.0	55.000	
Lower Spec Limit (LSL)	1.0	1.0	142.0	175.0	89	270.0	4.0	55	
Upper Spec Limit (USL)	N/A	N/A	144.0	175.0	89	280.0	4.0	55	
UOM	lb	lb	°C	in/min	psi	°F	seconds	psi	
Measuring Equipment Type	Force Gauge	Force Gauge	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual BandSealer Gauge ONLY ES-0059-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	Visual ImpulseSealer Gauge ONLY ES-0099-01	
Asset No.	ES-0028-01	ES-0028-01	N/A	N/A	N/A	ES-0099-01	ES-0099-01	ES-0099-01	
Cal Due Date	12/18/2024	12/18/2024	N/A	N/A	N/A	8/9/2024	8/9/2024	8/9/2024	
Sample ID	Top Strip	Bottom Strip	N/A	N/A	N/A	1	1	1	
Sample No.	Variable/ Attribute Data								
1: Start (Top Strip)	6.522	N/A	N/A	N/A	N/A	270.0	4.0	55.0	
2: Start (Bottom Strip)	N/A	5.925	N/A	N/A	N/A	270.0	4.0	55.0	
3: Middle (Top Strip)	6.260	N/A	N/A	N/A	N/A	270.0	4.0	55.0	
4: Middle (Bottom Strip)	N/A	5.162	N/A	N/A	N/A	270.0	4.0	55.0	
5: End (Top Strip)	6.012	N/A	N/A	N/A	N/A	270.0	4.0	55.0	
6: End (Bottom Strip)	N/A	6.245	N/A	N/A	N/A	270.0	4.0	55.0	
Result Summary									
Pass QTY	3	3	N/A	N/A	N/A	6	6	6	
Fail QTY	0	0	N/A	N/A	N/A	0	0	0	
Result (P/F)	P	P	N/A	N/A	N/A	P	P	P	
Statistical Summary									
Min	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	
Max	#REF!	#REF!	#REF!	#REF!	0.00	#REF!	#REF!	#REF!	
Mean	#REF!	#REF!	#REF!	#REF!	0.00	#REF!	#REF!	#REF!	
StDev	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	#REF!	
Comments									
N/A									
Disposition: ACCEPT							NCMR No: N/A		
Inspected By:							Date: 05/17/2024		

Date Inspected		Part No.		Rev.	Description						
05/15/2024		SLQ-8102		B	SILQ 2-Way Clearfract SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack						
Lot No.		Lot QTY		Approved Supplier(s)							
SLQ-81020515241		405		Pathway							
In-Process Inspection		Catheter has been coated with SILQ's polymer solution (SLQ-3000)	Catheter has been rinsed of cured polymer solution and is free of debris	Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001	Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002	Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".	Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.	Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton.	Shipper Box Inspection: 1) Correct quantity in Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the inner carton, or shipper box.	
		Revision	B	MH 07/29/24							
		Lower Spec Limit (LSL)	N/A	A-B A-B							
		Upper Spec Limit (USL)	N/A	0.300							
		UOM	EA	N/A							
		Measuring Equipment Type	Visual	Spectrophotometer							
		Asset No.	N/A	See Comments							
		Cal Due Date	N/A	See Comments							
		AQL	6.5	6.5							
		Sample Size	9	9							
Sample No.		Variable/ Attribute Data									
1	p		p		p		p		p		
2	p		p		p		p		p		
3	p		p		p		p		p		
4	p		p		p		p		p		
5	p		p		p		p		p		
6	p		p		p		p		p		
7	p		p		p		p		p		
8	p		p		p		p		p		
9	p		p		p		p		p		
		Result Summary									

[illegible]

Date Inspected		Part No.		Rev.		Description				
06/01/2024		SLQ-8102		B		SILQ 2-Way Cleartract SPT 18 Fr x 10 mL Catheter, Non-Sterile, 10-Pack				
Lot No.		Lot QTY		Rinsed catheter		Approved Supplier(s)				
SLQ-81020515241		405		Coated catheter absorbance @585nm ≥0.300 AU Ref TM-HDX-0001		Pathway				
Catheter has been coated with SILQ's polymer solution (SLQ-3000)		Catheter has been rinsed of cured polymer solution and is free of debris		Pouch Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 1, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".		Carton Label Inspection: 1) Label content matches DMR 2) P/N, Lot No., QTY 10, and Exp date on labels match Work Order requirements. 3) Labels are legible, not missing text, 4) Labels with no stains or discoloration. 5) Non HRI Text exp date format is "YYYY-MM-DD".				
In-Process Inspection		Rinsed catheter absorbance @440nm ≤0.020 AU Ref TM-HDX-0002		Pouch Inspection: 1) Catheter is properly oriented in the sleeve and pouch. 2) Pouch seal is continuous with no wrinkles, gaps, folds, bubbles or burn/overheating clear spots. Seal width at least 1/8 inch. 3) Pouch is with no tears, punctures or pin holes. 4) Pouch observed with no particulate/loose debris inside pouch. 5) Labels affixed & positioned in correct location, is legible.		Inner Carton Inspection: 1) Correct quantity in 10-Pack Box. 2) Product configuration matches specification. 3) Labels affixed & positioned in correct location. 4) No stains or discoloration that affects appearance. 5) No damage observed to the device, its sterile barrier, or inner carton. IP-0002				
Revision		B		B		B				
Lower Spec Limit (LSL)		N/A		0.300		N/A				
Upper Spec Limit (USL)		N/A		N/A		N/A				
UOM		EA		AU		EA				
Measuring Equipment Type		Visual		Spectrophotometer		Visual				
Asset No.		N/A		See Comments		N/A				
Cal Due Date		N/A		See Comments		N/A				
AQL		6.5		6.5		6.5				
Sample Size		9		9		9				
Sample No.						Variable/ Attribute Data				
1		N/A		0.413		N/A				
2		N/A		0.482		N/A				
3		N/A		0.446		N/A				
4		N/A		0.501		N/A				
5		N/A		0.423		N/A				
6		N/A		0.476		N/A				
7		N/A		0.491		N/A				
8		N/A		0.422		N/A				
9		N/A		0.472		N/A				



Result Summary										
Pass QTY	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Fail QTY	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Result (P/F)	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Statistical Summary										
Min	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Max	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Mean	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
StDev	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
Comments										
Rinse test and UV Absorbance Test were performed at SLQ. The rest of the visual inspections were performed at Pathway.										
Disposition:	ACCEPT							NCMR No:	N/A	
Inspected By:	Ethan Rao							Date:	02 June 2024	
Verified By:	[Signature]							Date:	08/02/2024	

	Label Generation & Reconciliation				Doc No.: FM-0033	
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>					
	Rev.: D Page 1 of 2					

Label Information						
Part No.	SLQ-LB-4008	Rev.	A	Lot No.	SLQ-81020515241 ☐ N/A	Label QTY ☐ N/A
Description	Label, Pouch, SILQ ClearTract SPT Catheter	UDI No.	00850041415161 ☐ N/A		Expiration Date:	2027-05-15 ☐ N/A

Label Quantity Calculation	
Required Label Quantity (A)	405
2% Overage ($A \times 0.02 = B$)	8
QC Retain Labels	2
Quantity of Labels Required ($A + B + 2 = C$)	415
Quantity of Labels Printed ($C = D$)	415
Production's Initial	 Date

Label Verification	
What is the revision of the label per the BOM?	A
What is the revision on the printed label?	A
Is the printed information correct per specifications?	☒ Yes ☐ No
First and last label printed is affixed to back of form	☒ Yes ☐ No
Quantity of labels released to production (E)	413
QC's Initials	 Date

Label Reconciliation	
Quantity of labels used by production (F)	413
Quantity scrapped by production (G)	0
Quantity of unused labels returned to QC and scrapped (H)	0
($E - F - G - H = 0$)	0
Are all labels accounted for?	☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

1. Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
2. Review the appropriate workstations and surrounding areas for the missing labels.

Investigational Results	
☒ N/A	

	Label Generation & Reconciliation		Doc No.: FM-0033
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>		Rev.: D
			Page 2 of 2
QC's Initials	MH	Date	07/29/2024



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	18 Fr, 10 mL
QUANTITY	1

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

R_x Only



STERILE EO



MR

LOT SLQ-81020515241

REF 211810SPT

 2027-05-15



(01) 00850041415161
(10) SLQ-81020515241
(17) 270515

SLQ-LB-4008 Rev. A

First



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

SIZE	18 Fr, 10 mL
QUANTITY	1

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

R_x Only



STERILE EO



MR

LOT SLQ-81020515241

REF 211810SPT

 2027-05-15



(01) 00850041415161
(10) SLQ-81020515241
(17) 270515

SLQ-LB-4008 Rev. A

Last



Label Generation & Reconciliation

Doc No.: FM-0033

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.

Rev.: D

Page 1 of 2

Label Information

Part No.	SLQ-LB-4009	Rev.	A	Lot No.	SLQ-81020515241 ☐ N/A	Label QTY	10 ☐ N/A
Description	Label, Box, SILQ ClearTract SPT Catheter	UDI No.	10850041415168 ☐ N/A	Expiration Date:	2027-05-15 ☐ N/A		

Label Quantity Calculation

Required Label Quantity (A)		41
2% Overage (A×0.02=B)		1
QC Retain Labels		2
Quantity of Labels Required (A+B+2=C)		44
Quantity of Labels Printed (C=D)		44
Production's Initial		Date 05/15/2024

Label Verification

What is the revision of the label per the BOM?		A
What is the revision on the printed label?		A
Is the printed information correct per specifications?		☒ Yes ☐ No
First and last label printed is affixed to back of form		☒ Yes ☐ No
Quantity of labels released to production (E)		42
QC's Initials	MH	Date 07/29/2024

Label Reconciliation

Quantity of labels used by production (F)		42
Quantity scrapped by production (G)		0
Quantity of unused labels returned to QC and scrapped (H)		0
(E-F-G-H=0)		0
Are all labels accounted for?		☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

1. Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
2. Review the appropriate workstations and surrounding areas for the missing labels.

Investigational Results

☒ N/A

	Label Generation & Reconciliation	Doc No.: FM-0033	
	This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.		Rev.: D
			Page 2 of 2
QC's Initials	ME	Date	07/29/2024

SIZE	18 Fr, 10 mL
QUANTITY	10

2-WAY 100% SILICONE
CLEARTRACT SPT CATHETER



R_x Only    **STERILE EO**   **MR**



(01) 10850041415168
(10) SLQ-81020515241
(17) 270515

LOT SLQ-81020515241
REF 211810SPT
2027-05-15

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex


 SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

First

SIZE	18 Fr, 10 mL
QUANTITY	10

2-WAY 100% SILICONE
CLEARTRACT SPT CATHETER

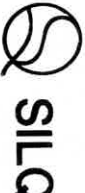

R_x Only    **STERILE EO**   **MR**



(01) 10850041415168
(10) SLQ-81020515241
(17) 270515

LOT SLQ-81020515241
REF 211810SPT
2027-05-15

Symbols are defined in Instruction for Use (IFU)
Store away from direct light exposure at room temperature
Not made with natural rubber latex

 SILQ Technologies Corp.
323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4009 Rev. A

Last



Label Generation & Reconciliation

Doc No.: FM-0033

This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.

Rev.: D

Page 1 of 2

Label Information

Part No.	SLQ-LB-4010	Rev.	A	Lot No.	SLQ-81020515241 ☐ N/A	Label QTY	100 ☐ N/A
Description	Label, Shipper Box, SILQ ClearTract SPT Catheter	UDI No.	10850041415168 ☐ N/A		Expiration Date:	2027-05-15 ☐ N/A	

Label Quantity Calculation

Required Label Quantity (A)	8		
2% Overage ($A \times 0.02 = B$)	1		
QC Retain Labels	2		
Quantity of Labels Required ($A + B + 2 = C$)	11		
Quantity of Labels Printed ($C = D$)	11		
Production's Initial		Date	05/15/2024

Label Verification

What is the revision of the label per the BOM?	A		
What is the revision on the printed label?	A		
Is the printed information correct per specifications?	☒ Yes ☐ No		
First and last label printed is affixed to back of form	☒ Yes ☐ No		
Quantity of labels released to production (E)	9		
QC's Initials		Date	05/15/2024

Label Reconciliation

Quantity of labels used by production (F)	9
Quantity scrapped by production (G)	0
Quantity of unused labels returned to QC and scrapped (H)	0
($E - F - G - H = 0$)	0
Are all labels accounted for?	☒ Yes ☐ No

If the counts cannot be reconciled, perform the following activities:

1. Verify that all products are labeled with the information recorded in the "Label Generation" section. Document the product quantity verified in designated field below.
2. Review the appropriate workstations and surrounding areas for the missing labels.

Investigational Results

☒ N/A

	Label Generation & Reconciliation	Doc No.: FM-0033	
	<i>This document contains confidential, proprietary information of Pathway LLC and shall not be copied or reproduced without prior written permission there from.</i>	Rev.: D	
		Page 2 of 2	
QC's Initials	MH	Date	07/29/2024



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF

211810SPT

LOT

SLQ-81020515241



2027-05-15

QTY

100



SILQ Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

First



2-WAY 100% SILICONE CLEARTRACT SPT CATHETER

REF

211810SPT

LOT

SLQ-81020515241



2027-05-15

QTY

100



SILQ Technologies Corp.

323 Sunny Isles Blvd, 7th Floor | Sunny Isles Beach, FL 33160

SLQ-LB-4010 Rev. A

Last

Bioburden Report



Pathway Medtech LLC
8779 Cottonwood Ave, STE 105
Santee, CA 92071

SAMPLE #: 24061001
DATE RECEIVED: 26-Jun-2024
TEST INITIATED: 03-Jul-2024
TEST COMPLETED: 08-Jul-2024

SAMPLE NAME/DESCRIPTION: SLQ-8102

UNIQUE IDENTIFIER: SLQ-81020515241

ALERT LIMIT: NG CFU/sample

ACTION LIMIT: NG CFU/sample

TEST RESULTS CFU/sample

Sample ID	Sample Designation	Aerobes	Aerobic Sporeformers	Fungi
1	NG	17	6	< 6
2	NG	17	20	< 6
3	NG	56	8	11
	Average	29.9	11.2	< 7.5
	Range	17 to 56	6 to 20	0 to 11
	Standard Deviation	18.5	6	~ 4.6

Sample with "<" numbers indicate 0 CFU (Colony Forming Units) were recovered and "~" indicates approximately. The stated value reflects any dilutions or Sample Item Portions used during the testing process.

CORRECTION FACTOR APPLIED: 1.40

PD #: 002075

REFERENCE: Processed according to LAB-002: Bioburden Product Test

COMMENTS: NA

AUTHORIZED

BY: 

DATE: 16 Jul 2024

Specific test results may not be indicative of the characteristics of any other samples from the same lot or similar lots. This test report shall not be reproduced, except in full, without prior written approval. Liability is limited to the costs of the tests. The uncertainty of measurement associated with the measurement result reported in this certificate is available from the organization upon request. Results apply to samples as received.

STERIS: Ethylene Oxide Certificate Of Processing

Prepared PATHWAY INNOVATION LLC

For

Ethylene Oxide Run ID 1162-14225

<u>Product Code</u>	<u>Lot Number</u>	<u>Quantity (EA)</u>
PWI-7374 Full Cycle	Unique Run ID: 81020515241	123

Preconditioning

<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>
Preconditioning Room	30-Jun-2024 4:52 AM	30-Jun-2024 10:45 PM	1073

		<u>Min</u>	<u>Max</u>
Duration (min)	Spec:	960	1920
Relative Humidity (%RH)	Spec:	40.0	80.0
	Actual:	53.0	69.0
Temperature (°F)	Spec:	105.0	125.0
	Actual:	111.0	116.0

Vessel Processing

<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>
Chamber 6	30-Jun-2024 10:53 PM	01-Jul-2024 8:20 AM	567
Cycle ID: PWI-7374	Check Value: N/A		

		<u>Min</u>	<u>Max</u>
Duration (min)	Spec:	1	1440

EO Drum Information

<u>Drum ID</u>	<u>Drum Serial Number</u>	<u>Drum Lot Number</u>	<u>EO Usage (LB)</u>
1993	E000269	UTLX902993E24	18.0

Total EO Usage: 18.0

EO Drum Information (Spec)

	<u>EO Usage (LB)</u>
Min	1.0
Max	100.0

Product Code

Lot Number

Quantity (EA)

Aeration

<u>Equipment Name</u>	<u>Start Date</u>	<u>End Date</u>	<u>Duration (min)</u>
Aeration Room	01-Jul-2024 8:29 AM	03-Jul-2024 8:55 AM	2906
		Total Aeration Time:	2906
		<u>Min</u>	<u>Max</u>
	Duration (min) Spec:	2880	4320
	Aeration in Room Temperature (°F) Spec:	103.0	123.0
	Actual:	109.0	114.0

EO Spec ID/Rev/Description ODMS Spec #5 Rev 1 (PWI-7374 Rev. 1 Medical Product)

Product meets Customer specifications; zero nonconformities occurred during this run.

Customer Specification ID: 769

Signature Manifest

Reviewed and E-Signed By: **Shirley Carter (Senior Quality Technician)**

DateTime Esigned 05-Jul-2024 2:39 PM

Document Content Revision 1

Operating facilities are in compliance with applicable state and federal regulations (FDA, NRC, EPA, and OSHA) and provide services under a quality system which meets the requirements of FDA QSR, EN/ISO 13485, and in alignment with EN ANSI/AAMI/ISO 11135.

Processing Location

43425 Business Park Dr
Temecula CA 92590
Phone: 951-694-9340

*** End of Ethylene Oxide Certificate of Processing Document ***

STERIS: Receiving Report
PATHWAY INNOVATION LLC (8130)
Ethylene Oxide

Receipt ID
1161-10583

Product Received From:

PATHWAY INNOVATION LLC
8779 Cottonwood Ave
Suite 105
Santee CA 92071
325

Carrier: NorthStar

Trailer: N/A

Bill of Lading: N/A

Seal: None

PO(s): P1800

Product counts provided by Customer

Receipt Date/Time: 26-Jun-2024 1:22 PM

Pallet Count: 1

Unloaded By: Dave Smith

Counted By: Diana Burns

Container Unload: ☐

Promise Date: 26-Jul-2024 1:22 PM

Item [Pkg]

PWL-7374 Full Cycle [1]

Lot Number

Unique Run ID: 81020515241

Quantity

Customer STERIS CD UOM

123 123 0 EA

Receipt Totals:

123 123 0

Inventory by Received Pallet

Pallet (SN) Item ID [Pkg]

1 (81357) PWL-7374 Full Cycle [1]

Lot Number

Unique Run ID: 81020515241

Quantity

Customer STERIS UOM

123 123 EA

Warehouse

Location

N/A

Pallet

Total

1

Signature Manifest

Reviewed By: Diana Burns (Manager, Production)

Date/Time Signed 26-Jun-2024 1:45 PM

Document Content Revision 1

Processing Location

43425 Business Park Dr
Temecula CA 92590
Phone: 951-694-9340

GW-08-004-01

Last Revised in Rel 2.0.0.0

Rel Date: 31-Mar-2017

Document ID: 50788

Page 1 of 1

PATHWAY INNOVATION LLC

**8779 Cottonwood Ave
Suite 105
Santee, CA 92071
(619)415-0103**

Packing List No: PL-SLQ-0048

Page No: 1 of 1

Ship To:

Company: STERIS AST
Attn: Diana Burns
Address: 43425 Business Park Drive
Temecula, CA 92590

Phone: 951-694-9340

Ship Date: 6/17/2024

PO# P1800


Tracking No: TBD

Delivery Company: NorthStar Courier

Delivery Method: Courier

FOB: Santee, CA

[illegible]

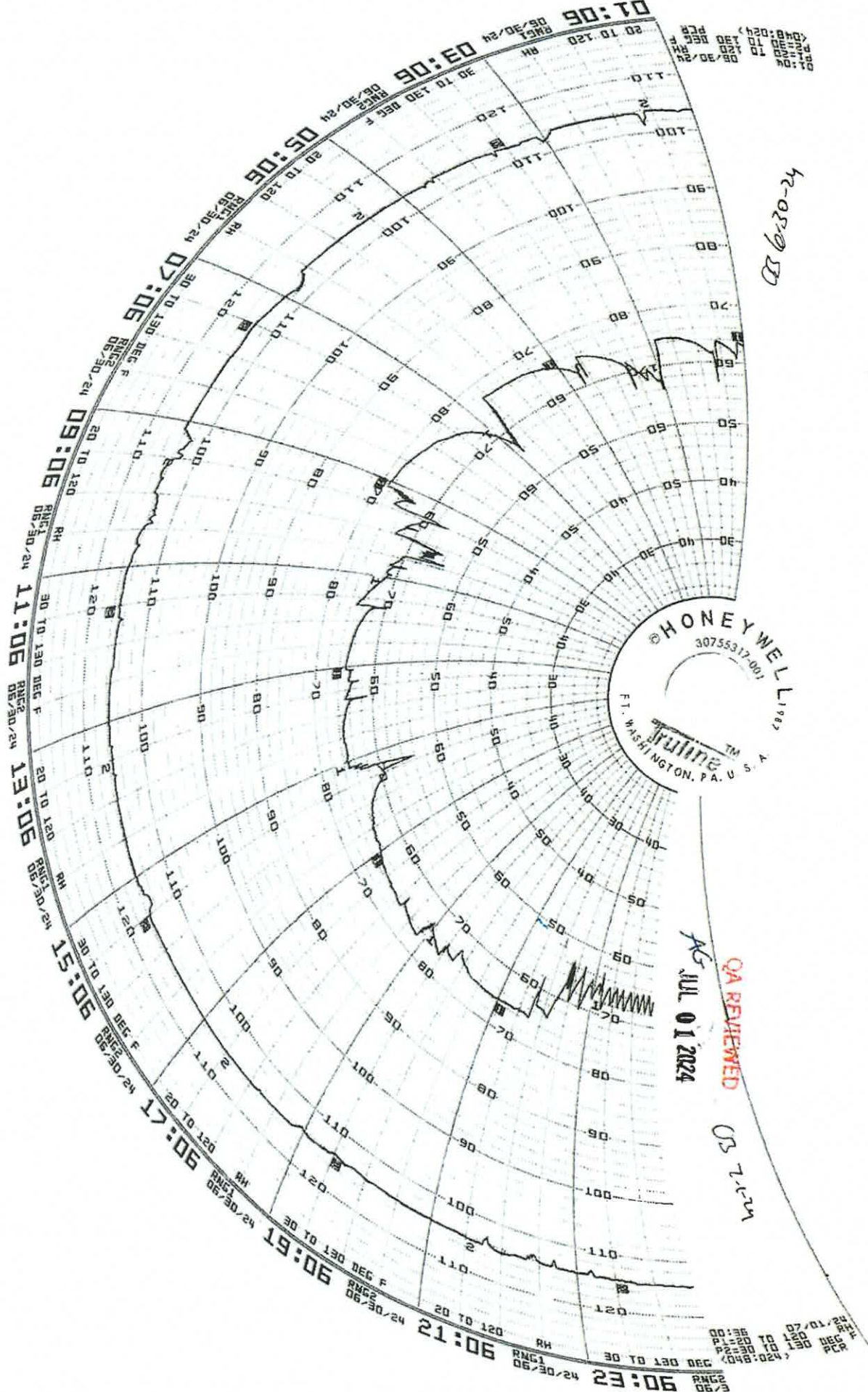

Customer Name Pathway FC

Cycle Number 7374

Load # PWI 1162-14225

Cycle Last Modified: 4/30/2018 1:56:49PM

Date Processed 06/30/2024 22:53:15

Chamber # 6

Cycle Counter # 4,420

Pressure Units "HgA

Temperature Units Fahrenheit

Phase Name

Initial Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
06/30/2024 22:53:15	00:00:00	28.9	127
06/30/2024 22:58:14	00:05:00	21.5	126
06/30/2024 23:03:14	00:10:00	14.7	126
06/30/2024 23:08:15	00:15:00	8.2	125
06/30/2024 23:13:15	00:20:00	4.2	125
06/30/2024 23:18:15	00:25:00	2.3	126
06/30/2024 23:19:44	00:26:29	2.0	126

Summary of Initial Evacuation

Pressure Start	28.9	Min Pressure	2.0	Min Temp	125	Pressure Change	-26.9
Pressure End	2.0	Max Pressure	28.9	Max Temp	127	Phase Time	00:26:29

Stabilization

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
06/30/2024 23:19:44	00:26:29	2.0	126
06/30/2024 23:21:43	00:28:29	2.1	127
06/30/2024 23:21:43	00:28:29	2.1	127

Summary of Stabilization

Pressure Start	2.0	Min Pressure	2.0	Min Temp	126	Pressure Change	0.1
Pressure End	2.1	Max Pressure	2.1	Max Temp	127	Phase Time	00:02:00

Leak Check

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
06/30/2024 23:21:43	00:28:29	2.1	127
06/30/2024 23:22:44	00:29:29	2.1	127
06/30/2024 23:23:43	00:30:29	2.1	127
06/30/2024 23:24:44	00:31:29	2.1	127
06/30/2024 23:25:44	00:32:30	2.2	127
06/30/2024 23:26:44	00:33:29	2.2	127

Summary of Leak Check

Pressure Start	2.1	Min Pressure	2.1	Min Temp	127	Pressure Change	0.1
Pressure End	2.2	Max Pressure	2.2	Max Temp	127	Phase Time	00:05:00

Initial Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
06/30/2024 23:26:44	00:33:29	2.2	127
06/30/2024 23:31:44	00:38:29	9.3	129
06/30/2024 23:36:44	00:43:29	16.7	128
06/30/2024 23:41:44	00:48:30	24.2	128

06/30/2024 23:44:16 00:51:01 28.0 128

Summary of Initial Nitrogen Wash

Pressure Start	2.2	Min Pressure	2.2	Min Temp	127	Pressure Change	25.8
Pressure End	28.0	Max Pressure	28.0	Max Temp	129	Phase Time	00:17:32

Initial Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
06/30/2024 23:44:16	00:51:01	28.0	128
06/30/2024 23:49:16	00:56:01	21.3	124
06/30/2024 23:54:16	01:01:01	14.3	124
06/30/2024 23:59:16	01:06:01	7.8	123
07/01/2024 00:04:15	01:11:02	4.2	124
07/01/2024 00:09:16	01:16:02	2.4	125
07/01/2024 00:11:03	01:17:49	2.0	125

Summary of Initial Evacuation

Pressure Start	28.0	Min Pressure	2.0	Min Temp	123	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	128	Phase Time	00:26:47

Humidification

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	Avg. % RH
07/01/2024 00:11:03	01:17:49	2.0	125	3
07/01/2024 00:12:03	01:18:49	2.4	126	10
07/01/2024 00:13:03	01:19:49	2.7	126	21
07/01/2024 00:14:03	01:20:49	2.9	126	29
07/01/2024 00:14:12	01:20:58	3.0	126	30

Summary of Humidification

Pressure Start	2.0	Min Pressure	2.0	Min Temp	125	Pressure Change	1.0	Min RH	3
Pressure End	3.0	Max Pressure	3.0	Max Temp	126	Phase Time	00:03:09	Max RH	30

Humidity Dwell

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	Avg. % RH
07/01/2024 00:14:12	01:20:58	3.0	126	30
07/01/2024 00:19:12	01:25:58	3.0	126	32
07/01/2024 00:24:12	01:30:58	3.0	125	32
07/01/2024 00:29:12	01:35:58	3.0	125	32
07/01/2024 00:34:13	01:40:58	3.0	125	33
07/01/2024 00:39:13	01:45:58	3.1	125	33
07/01/2024 00:44:13	01:50:58	3.1	125	32
07/01/2024 00:49:12	01:55:58	3.1	125	33
07/01/2024 00:54:12	02:00:58	3.1	125	33
07/01/2024 00:59:13	02:05:58	3.1	125	33
07/01/2024 01:04:13	02:10:59	3.1	125	34
07/01/2024 01:09:13	02:15:59	3.1	125	33
07/01/2024 01:14:12	02:20:58	3.2	125	34

Summary of Humidity Dwell

Pressure Start	3.0	Min Pressure	3.0	Min Temp	125	Pressure Change	0.2	Min RH	29
Pressure End	3.2	Max Pressure	3.2	Max Temp	126	Phase Time	01:00:00	Max RH	34

Sterilant Injection

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	EO Drum 1 (lb)	EO Drum 2 (lb)	EO CONC
07/01/2024 01:14:12	02:20:58	3.2	125	355	0	1
07/01/2024 01:19:13	02:25:58	7.7	126	350	0	232
07/01/2024 01:24:13	02:30:58	12.8	126	342	0	537
07/01/2024 01:25:40	02:32:26	14.3	126	340	0	619

Summary of Sterilant Injection

Pressure Start	3.2	Min Pressure	3.2	Min Temp	125	Drum 1 Start	355	Drum 2 Start	0
Pressure End	14.3	Max Pressure	14.3	Max Temp	126	Drum 1 End	340	Drum 2 End	0
Pressure Change	11.2	MIN EO CONC	1.00	EtO Used 1	15	EtO Used 2	0		
Phase Time	00:11:28	MAX EO CONC	619.00	Total EO Weight	15				

Exposure

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp	EO Drum 1 (lb)	EO Drum 2 (lb)	EO CONC
07/01/2024 01:25:40	02:32:26	14.3	126	340	0	619
07/01/2024 01:30:41	02:37:26	14.3	125	341	0	642
07/01/2024 01:35:41	02:42:26	14.2	125	340	0	636
07/01/2024 01:40:41	02:47:26	14.2	125	340	0	629
07/01/2024 01:45:40	02:52:26	14.1	125	340	0	623
07/01/2024 01:50:40	02:57:26	14.3	125	340	0	630
07/01/2024 01:55:41	03:02:26	14.3	125	339	0	627
07/01/2024 02:00:41	03:07:27	14.3	125	339	0	622
07/01/2024 02:05:41	03:12:27	14.3	125	339	0	619
07/01/2024 02:10:41	03:17:27	14.3	125	339	0	616
07/01/2024 02:15:42	03:22:27	14.2	125	339	0	612
07/01/2024 02:20:41	03:27:27	14.2	125	339	0	610
07/01/2024 02:25:41	03:32:27	14.2	125	339	0	607
07/01/2024 02:30:41	03:37:27	14.2	125	339	0	604
07/01/2024 02:35:41	03:42:27	14.2	125	339	0	603
07/01/2024 02:40:42	03:47:27	14.2	125	339	0	600
07/01/2024 02:45:42	03:52:27	14.2	125	339	0	599
07/01/2024 02:50:42	03:57:27	14.2	125	338	0	597
07/01/2024 02:55:41	04:02:27	14.2	125	338	0	594
07/01/2024 03:00:42	04:07:28	14.2	125	338	0	594
07/01/2024 03:05:42	04:12:28	14.2	125	338	0	592
07/01/2024 03:10:42	04:17:28	14.2	125	338	0	591
07/01/2024 03:15:42	04:22:28	14.2	125	338	0	591
07/01/2024 03:20:43	04:27:28	14.2	125	338	0	588
07/01/2024 03:25:42	04:32:28	14.2	125	338	0	587
07/01/2024 03:30:42	04:37:28	14.2	125	338	0	588
07/01/2024 03:35:42	04:42:28	14.2	125	338	0	585
07/01/2024 03:40:42	04:47:28	14.2	125	338	0	583
07/01/2024 03:45:43	04:52:28	14.2	125	338	0	583
07/01/2024 03:50:43	04:57:28	14.2	125	338	0	582
07/01/2024 03:55:43	05:02:28	14.2	125	338	0	581
07/01/2024 04:00:42	05:07:28	14.2	125	338	0	581
07/01/2024 04:05:43	05:12:29	14.2	125	338	0	580
07/01/2024 04:10:43	05:17:29	14.2	125	338	0	579
07/01/2024 04:15:43	05:22:29	14.3	125	337	0	579
07/01/2024 04:20:43	05:27:29	14.3	125	337	0	579
07/01/2024 04:25:43	05:32:29	14.3	125	337	0	577
07/01/2024 04:30:44	05:37:29	14.3	125	337	0	577
07/01/2024 04:35:43	05:42:29	14.3	125	337	0	576
07/01/2024 04:40:43	05:47:29	14.3	125	337	0	575
07/01/2024 04:45:43	05:52:29	14.3	125	337	0	576
07/01/2024 04:50:44	05:57:29	14.3	125	337	0	575
07/01/2024 04:55:44	06:02:29	14.3	125	337	0	573
07/01/2024 05:00:44	06:07:29	14.3	125	337	0	573
07/01/2024 05:05:44	06:12:30	14.3	125	337	0	573
07/01/2024 05:10:43	06:17:30	14.3	125	337	0	572
07/01/2024 05:15:44	06:22:30	14.3	125	337	0	570
07/01/2024 05:20:44	06:27:30	14.3	125	337	0	572
07/01/2024 05:25:40	06:32:26	14.3	125	337	0	571

Summary of Exposure

Pressure Start	14.3	Min Pressure	14.1	Min Temp	125	Drum 1 Start	340	Drum 2 Start	0
Pressure End	14.3	Max Pressure	14.4	Max Temp	126	Drum 1 End	337	Drum 2 End	0
Pressure Change	0.0	MIN EO CONC	570.00	EtO Used 1	3	EtO Used 2	0		
Phase Time	04:00:00	MAX EO CONC	642.00						
						Total EO Weight	3		

Sterilant Removal

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
07/01/2024 05:25:40	06:32:26	14.3	125
07/01/2024 05:30:40	06:37:26	9.9	123
07/01/2024 05:35:40	06:42:26	5.5	123
07/01/2024 05:40:41	06:47:26	2.9	123
07/01/2024 05:44:26	06:51:12	2.0	123

Summary of Sterilant Removal

Pressure Start	14.3	Min Pressure	2.0	Min Temp	123	Pressure Change	-12.3
Pressure End	2.0	Max Pressure	14.3	Max Temp	125	Phase Time	00:18:46

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
07/01/2024 05:44:26	06:51:12	2.0	123
07/01/2024 05:49:27	06:56:12	9.8	126
07/01/2024 05:54:27	07:01:12	16.9	126
07/01/2024 05:59:27	07:06:12	24.2	126
07/01/2024 06:02:01	07:08:47	28.0	127

Summary of Post Nitrogen Wash

Pressure Start	2.0	Min Pressure	2.0	Min Temp	123	Pressure Change	26.0
Pressure End	28.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:17:35

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
07/01/2024 06:02:01	07:08:47	28.0	127
07/01/2024 06:07:02	07:13:47	21.3	123
07/01/2024 06:12:02	07:18:47	14.3	122
07/01/2024 06:17:02	07:23:47	7.8	121
07/01/2024 06:22:01	07:28:48	4.3	122
07/01/2024 06:27:02	07:33:48	2.5	123
07/01/2024 06:29:12	07:35:58	2.0	123

Summary of Post Evacuation

Pressure Start	28.0	Min Pressure	2.0	Min Temp	121	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:27:11

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
07/01/2024 06:29:12	07:35:58	2.0	123
07/01/2024 06:34:12	07:40:58	9.0	126
07/01/2024 06:39:13	07:45:58	16.6	126
07/01/2024 06:44:12	07:50:58	24.1	127
07/01/2024 06:46:52	07:53:37	28.0	127

Summary of Post Nitrogen Wash

Pressure Start	2.0	Min Pressure	2.0	Min Temp	123	Pressure Change	26.0
Pressure End	28.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:17:39

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
07/01/2024 06:46:52	07:53:37	28.0	127
07/01/2024 06:51:52	07:58:37	21.3	123
07/01/2024 06:56:51	08:03:37	14.3	122
07/01/2024 07:01:51	08:08:37	7.8	121
07/01/2024 07:06:52	08:13:38	4.3	122
07/01/2024 07:11:52	08:18:38	2.5	123
07/01/2024 07:14:05	08:20:50	2.0	123

Summary of Post Evacuation

Pressure Start	28.0	Min Pressure	2.0	Min Temp	121	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:27:13

Post Nitrogen Wash

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
07/01/2024 07:14:05	08:20:50	2.0	123
07/01/2024 07:19:04	08:25:50	9.1	126
07/01/2024 07:24:04	08:30:50	16.6	127
07/01/2024 07:29:04	08:35:50	24.1	127
07/01/2024 07:31:44	08:38:29	28.0	127

Summary of Post Nitrogen Wash

Pressure Start	2.0	Min Pressure	2.0	Min Temp	123	Pressure Change	26.0
Pressure End	28.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:17:39

Post Evacuation

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
07/01/2024 07:31:44	08:38:29	28.0	127
07/01/2024 07:36:43	08:43:29	21.3	123
07/01/2024 07:41:44	08:48:30	14.3	122
07/01/2024 07:46:44	08:53:30	7.9	121
07/01/2024 07:51:44	08:58:30	4.3	122
07/01/2024 07:56:44	09:03:30	2.5	123
07/01/2024 07:58:56	09:05:42	2.0	123

Summary of Post Evacuation

Pressure Start	28.0	Min Pressure	2.0	Min Temp	121	Pressure Change	-26.0
Pressure End	2.0	Max Pressure	28.0	Max Temp	127	Phase Time	00:27:13

Air Backfill

Date / Time	Cycle Time	Avg. CH Press	Avg. CH Temp
07/01/2024 07:58:56	09:05:42	2.0	123
07/01/2024 08:03:56	09:10:42	7.8	125
07/01/2024 08:08:57	09:15:42	14.1	125
07/01/2024 08:13:57	09:20:42	20.3	125
07/01/2024 08:18:57	09:25:43	26.5	124
07/01/2024 08:20:12	09:26:57	27.9	124

Summary of Air Backfill

Pressure Start	2.0	Min Pressure	2.0	Min Temp	123	Pressure Change	25.9
Pressure End	27.9	Max Pressure	27.9	Max Temp	125	Phase Time	00:21:15


Customer Name Pathway FC
Load # PWI 1162-14225

STENS
Applied Sterilization Technologies
Temecula, CA


End of Cycle

Production Review SEUS Date 7/1/24
Signature

Quality Review B.S. Carter Date 7/5/24
Signature

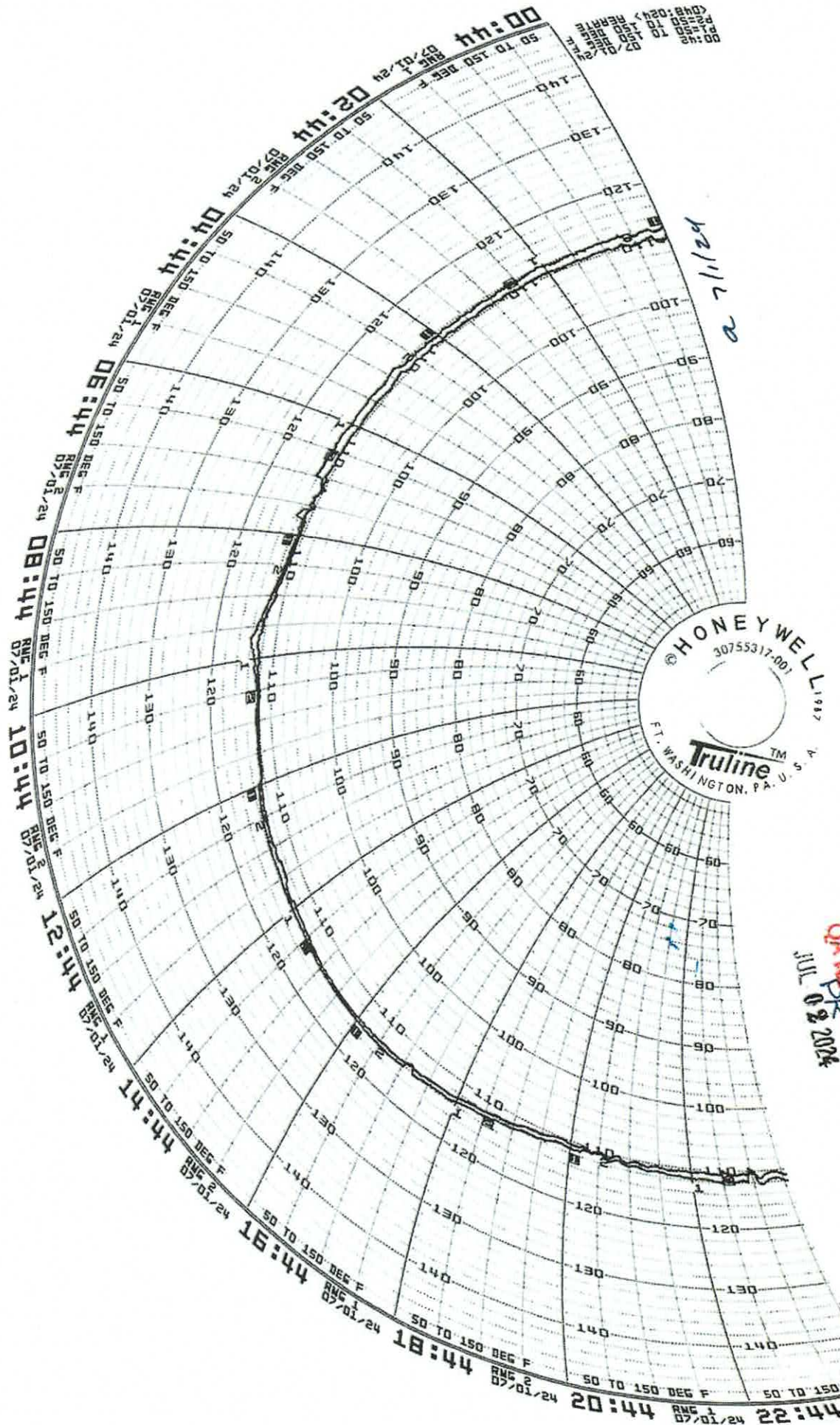

QA REVIEWED
JUL 02 2024

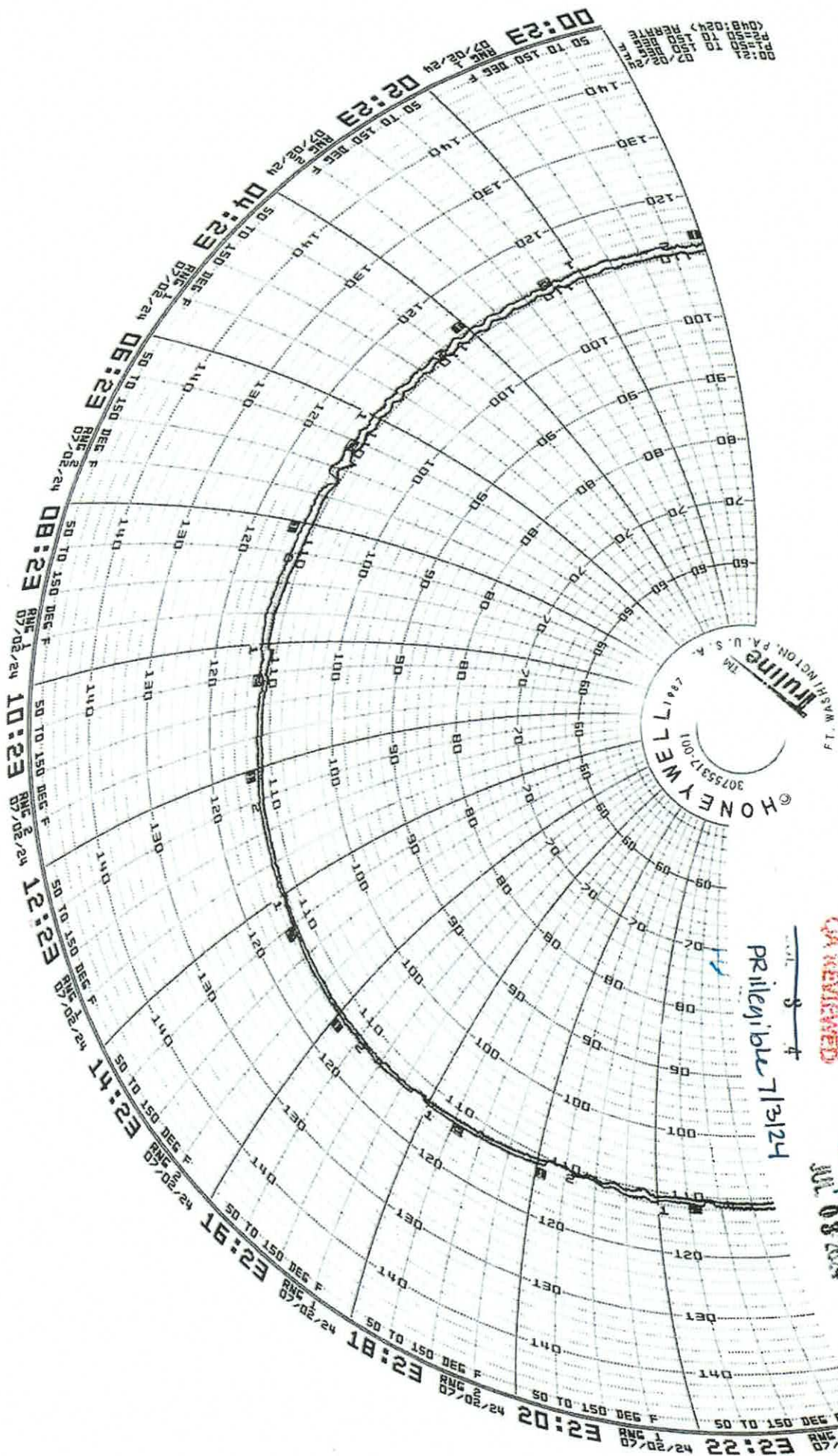


Y2-2-L 20

Prilegible 7/3/24

QA REVIEWED

Y2-3-L 20
QA REVIEWED
JUL 08 2024

00:11
01:11
02:11
03:11
04:11
05:11
06:11
07:11
08:11
09:11
10:11
11:11
12:11
13:11
14:11
15:11
16:11
17:11
18:11
19:11
20:11
21:11
22:11
23:11
00:11
01:11
02:11
03:11
04:11
05:11
06:11
07:11
08:11
09:11
10:11
11:11
12:11
13:11
14:11
15:11
16:11
17:11
18:11
19:11
20:11
21:11
22:11
23:11

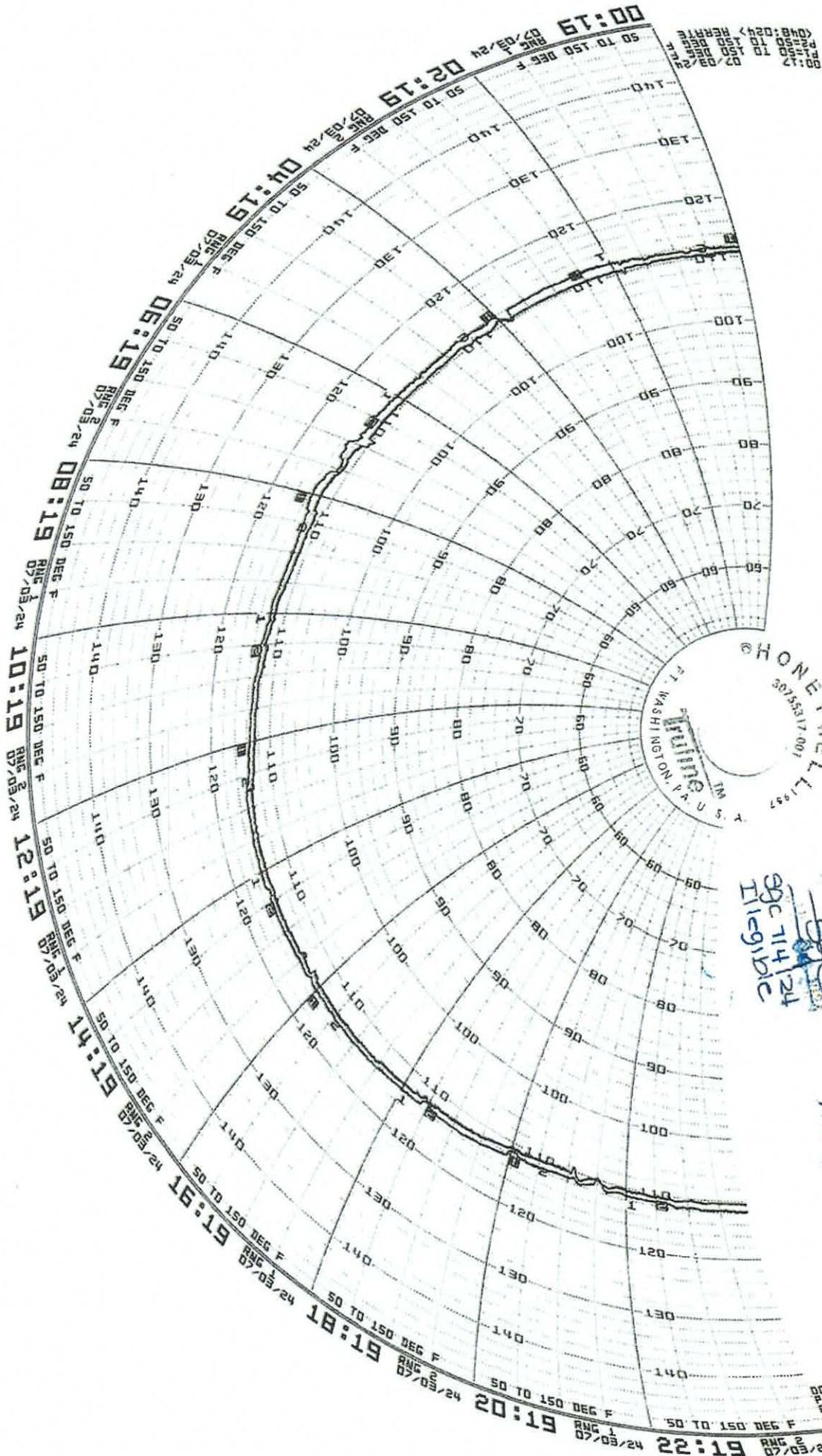

RR 7-3-24

HONEYWELL
30755317001
F. WASHINGTON, D.C. U.S.A.

30714124
IICIGIC

h2-4-L 02

00:19
01:19
02:19
03:19
04:19
05:19
06:19
07:19
08:19
09:19
10:19
11:19
12:19
13:19
14:19
15:19
16:19
17:19
18:19
19:19
20:19
21:19
22:19
23:19
24:19
25:19
26:19
27:19
28:19
29:19
30:19
31:19
32:19
33:19
34:19
35:19
36:19
37:19
38:19
39:19
40:19
41:19
42:19
43:19
44:19
45:19
46:19
47:19
48:19
49:19
50:19
51:19
52:19
53:19
54:19
55:19
56:19
57:19
58:19
59:19
60:19
61:19
62:19
63:19
64:19
65:19
66:19
67:19
68:19
69:19
70:19
71:19
72:19
73:19
74:19
75:19
76:19
77:19
78:19
79:19
80:19
81:19
82:19
83:19
84:19
85:19
86:19
87:19
88:19
89:19
90:19
91:19
92:19
93:19
94:19
95:19
96:19
97:19
98:19
99:19
100:19