

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

Pathway Medtech, LLC.

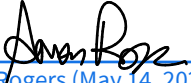
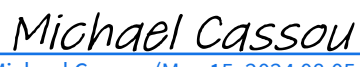
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report , ClearTrack Catheter

Prepared For:

SILQ Technologies, Corp.

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

SIGNATURE / APPROVAL PAGE / RISK MANAGEMENT TEAM:

Approver(s):	Signatures:	Date:
Author / Quality:	 Aaron Rogers (May 14, 2024 16:49 PDT)	05/14/2024
Customer:	 Ethan Rao (May 14, 2024 17:03 PDT)	05/14/2024
Engineering:	 Will Hail (May 16, 2024 10:24 PDT)	05/16/2024
Manufacturing:	 Michael Cassou (May 15, 2024 08:05 PDT)	05/15/2024

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

CHANGE HISTORY

Revision	Change Description
A	Initial Release
B	Add HDX-160510, general updates, add verification references, MP-HDX-0000 to MP-HDX-0004
C	Add SLQ-1610, format updates to tables in §5.0, add sterilization process #14, update verification references
D	Update Scope Table 1, add reference to RM Plan §3.2, add 4 FMs to #1 Receiving Inspection in Table 6
E	Update format, Update Table 1 for SPT, Update #s in Table 6, Add new risks #12, 21, 22, 24, 27-29, 31, 32, 36, 43, 44
F	Add reference to TR-SLQ-0007

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

TABLE OF CONTENTS

Signature / Approval Page / Risk Management Team:..... 2

Change History..... 3

Table of Contents..... 4

1.0 Purpose 5

2.0 Scope..... 5

3.0 Definitions 6

4.0 Risk Acceptability Determination 7

5.0 Risk Assessment Tables 8

6.0 Risk / Benefit Assessment 15

7.0 Conclusion 16

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

1.0 PURPOSE

The purpose of this risk report is to identify the associated risk with the manufacturing process of the ClearTract Catheter.

2.0 SCOPE

2.1 This risk report applies to the products identified with the associated part number(s) in the following table.

Table 1 Product Part Numbers

REF / Part Number	Description
HDX-1410	2-way 100% Silicone ClearTract Catheter, 14 Fr/10ml
HDX-1605	2-way 100% Silicone ClearTract Catheter, 16 Fr/5ml
HDX-1810	2-way 100% Silicone ClearTract Catheter, 18 Fr/10ml
HDX-160510	2-way 100% Silicone ClearTract Catheter, 16 Fr/5-10ml
SLQ-1610	2-Way 100% Silicone ClearTract Catheter, 16Fr/10ml
SLQ-211410	2-way 100% Silicone ClearTract Catheter, 14 Fr/10ml
SLQ-211605	2-way 100% Silicone ClearTract Catheter, 16 Fr/5ml
SLQ-211610	2-way 100% Silicone ClearTract Catheter, 16 Fr/10ml
SLQ-211810	2-way 100% Silicone ClearTract Catheter, 18 Fr/10ml
SLQ-212430	2-way 100% Silicone ClearTract Catheter, 24 Fr/30ml
SLQ-211410SPT	2-way 100% Silicone ClearTract SPT Catheter, 14 Fr/10ml
SLQ-211605SPT	2-way 100% Silicone ClearTract SPT Catheter, 16 Fr/5ml
SLQ-211610SPT	2-way 100% Silicone ClearTract SPT Catheter, 16 Fr/10ml
SLQ-211810SPT	2-way 100% Silicone ClearTract SPT Catheter, 18 Fr/10ml
SLQ-212430SPT	2-way 100% Silicone ClearTract SPT Catheter, 24 Fr/30ml

2.2 Members of the Risk Management Team that participated in developing this report are indicated on Page 2.

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

3.0 DEFINITIONS

Term	Meaning
#	Enumerate each Process with a sequential number
Component / Process Description	Describes the component / operation that is being evaluated for risk
PH#	Enumerate each Potential Failure Mode with a sequential number
Potential Failure Mode (FM)	The potential consequences of the failure
Severity (S)	Measure of the possible consequences of a hazard
Potential Effect of Failure	All conceivable effects of the potential failure mode
Potential Cause of Failure	All conceivable causes of the potential failure mode
Risk Controls	Process through which decisions are reached and protective measures are implemented for reduce risk to or maintaining risk within specified levels.
Risk Controls ¹	Initial risk controls that are incorporated.
Risk Controls ²	Secondary risk controls, also known as "Mitigation".
Verification	Evidence source of linkage that initial (i.e., Risk Controls ¹) risk controls have implemented, verified, or validated.
Probability	Likelihood of the hazard occurring
Probability ¹ (P ¹)	The probability assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls ¹)
Probability ² (P ²)	The probability assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls ²)
Detection	The ability to screen out a hazard through verification activities (e.g., inspection).
Detection ¹ (D ¹)	The detection assessment is post the implementation/ verification of initial risk controls (i.e., Risk Controls ¹)
Detection ² (D ²)	The detection assessment is post the implementation/ verification of the secondary risk controls (i.e., Risk Controls ²)
Essential Performance (EP)	Identifying design feature that is essential for proper functioning of the device.
Frequency	A 'Frequency Score' will be calculated using the ratings from probability and detection components. The frequency equation is the following: $\text{Frequency} = \sqrt{P \times D}$ Note: The value calculated should be rounded up.
Frequency ¹ (F ¹)	The frequency score calculated using Risk Controls ¹ and Detection ¹
Frequency ² (F ²)	The frequency score calculated using Risk Controls ² and Detection ²
Risk Rating	Residual risk determination based on the concept of 'As Far As Possible (AFAP)' approach. The AFAP approach provides three (3) possible outcomes, which are 1) the residual risk exceeds the risk acceptability (RED), 2) the residual risk is acceptable to the risk acceptability criteria (GREEN), and 3) the residual risk is between the two (2) states indicated and it can be an option to lower the risk to the lowest practical level without economic considerations (YELLOW).
Risk Rating ¹ (R ¹)	This calculation is post the implementation of initial risk controls (i.e., Risk Controls ¹)
Risk Rating ² (R ²)	This calculation is post the implementation of the secondary risk controls (i.e., Risk Controls ²)

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

4.0 RISK ACCEPTABILITY DETERMINATION

- 4.1 The criteria of risk acceptability will be determined with a Risk Priority Number (RPN) and on the concept of ‘As Far As Possible’ approach. Refer to the Risk Rating definition in §3.0.
- 4.2 There are three (3) components used in determining risk acceptability, which are Severity, Probability, and Detection. Refer to the following tables:

Table 2 Severity Ratings

Index	Description
1	No Patient effect
2	Minor injury that resolves without treatment
3	Moderate injury that may require treatment
4	Serious injury requiring medical intervention
5	Life threatening injury, disability, death

Table 3 Probability Ratings

Index	Description
1	Improbable - <0.001%
2	Remote – 0.001% to .01%
3	Occasional – 0.01% to 0.1%
4	Likely – 0.1% to 1%
5	Frequent - >1%

Table 4 Detection Ratings

Index	Description
1	Very High – Current controls are almost certain to detect the failure mode.
2	High – High likelihood that current controls will detect the failure mode.
3	Moderate – Moderate likelihood that current controls will detect the failure mode.
4	Low – Low likelihood that current controls will detect the failure mode.
5	Remote – Remote likelihood that current controls will detect the failure mode.

- 4.3 Risk acceptability will be determined per the following risk index:

Table 5 Risk Rating

5	5	10	15	20	25	Legend
---	---	----	----	----	----	--------

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

	4	4	8	12	16	20	Risk is not acceptable. Risk Mitigation is required.
	3	3	6	9	12	15	
	2	2	4	6	8	10	Risk is acceptable (AFAP). Risk Mitigation is recommended.
	1	1	2	3	4	5	
		1	2	3	4	5	Risk is acceptable. Risk Mitigation not required.
		Frequency					

Note: The purpose of using the square root of the probability and detection scores is to deemphasize these two (2) components, which places more emphasize on the severity component.

5.0 RISK ASSESSMENT TABLES

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
1	Line Clearance / Area Preparation	Incorrect Material	1	Non-Functional Device	Operator Error	Process Instructions Process Settings/Line Clearance Form DMR	WI-0014 FM-0039	1	1	1	1	N/A	N/A	N/A	N/A	1
2		Material Expired	1	Non-Functional Device	Operator Error	DMR	FM-0039	1	1	1	1	N/A	N/A	N/A	N/A	1
3		Improper Material Handling	1	Damaged Materials	Operator Error	Process Instructions	WI-0014	1	1	1	1	N/A	N/A	N/A	N/A	1
4		Out of Calibration Equipment	3	Non-Functional Device	Operator Error	Process Instructions	WI-0014	1	1	1	3	N/A	N/A	N/A	N/A	3
						Equipment Qualification	IQ Report (TR)	1	1	1	3	N/A	N/A	N/A	N/A	3
5		Untrained Operators	1	Compromised Device	Inadequate Training	Process Instructions	WI-0014	1	1	1	1	N/A	N/A	N/A	N/A	1
6	Receiving Inspection / Quality control	Accept Non-Conforming Material	2	Procedural Complication	Operator Error	Process Instructions	QP-009	1	1	1	2	N/A	N/A	N/A	N/A	2
			1		Equipment Error	Calibration Program	QP-020	1	1	1	1	N/A	N/A	N/A	N/A	1

Pathway Medtech, LLC.		Type:	Risk Management (RM)			Doc. No.:		RM-0021				Revision:		F		
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter																
#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
						Preventative Maintenance Program	QP-015	1	1	1	1	N/A	N/A	N/A	N/A	1
7		Reject Conforming Material	1	None	Operator Error	Process Instructions	QP-009	1	1	1	1	N/A	N/A	N/A	N/A	1
					Equipment Error	Calibration Program	QP-020	1	1	1	1	N/A	N/A	N/A	1	
						Preventative Maintenance Program	QP-015	1	1	1	1	N/A	N/A	N/A	1	
8		Incorrect / Illegible / No Identification	1	Potential traceability issue	Supplier Error	Purchasing Process	QP-013	1	1	1	1	N/A	N/A	N/A	N/A	1
						Supplier Qualification	QP-005	1	1	1	1	N/A	N/A	N/A	N/A	1
9		Damaged product	1	Non-Functional Device	Supplier Error – Incorrect Packaging	Purchasing Process	QP-013	1	1	1	1	N/A	N/A	N/A	N/A	1
						Supplier Qualification	QP-005	1	1	1	1	N/A	N/A	N/A	N/A	1
					Shipper Error – Mishandling	Process Instructions	QP-009	1	1	1	1	N/A	N/A	N/A	N/A	1
10		Contamination	2	Cleanroom contamination	Supplier Error	Purchasing Process	QP-013	1	1	1	1	N/A	N/A	N/A	N/A	1
						Supplier Qualification	QP-005	1	1	1	1	N/A	N/A	N/A	N/A	1
					Shipper Error – Mishandling	Process Instructions	QP-009	1	1	1	1	N/A	N/A	N/A	N/A	1
11		Documentation Error / Incomplete / No CoC/CoA	1	None	Supplier Error	Purchasing Process	QP-013	1	1	1	1	N/A	N/A	N/A	N/A	1
						Supplier Qualification	QP-005	1	1	1	1	N/A	N/A	N/A	N/A	1
						Training Program	QP-0017	1	1	1	1	N/A	N/A	N/A	N/A	1
12	Environmental Controls /		3	Microbial Contamination	Operator Error	Process Instructions	WI-0012	2	3	3	9	Product Sterilization	1	1	1	3

2021 Rev 1 - ECO-001003 / ECO-001003 / Int. Production / Effective: 03/10/2024 10:33:24 AM

Pathway Medtech, LLC.			Type:	Risk Management (RM)		Doc. No.:		RM-0021				Revision:		F		
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter																
#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
	Controlled Environment Area (CEA)	Increased particulate levels			Equipment Error	Preventative Maintenance Program	QP-015	1	3	2	6	Product Sterilization	1	1	1	3
13	Catheter tube forming / curing	Incomplete / Malformed Tube	2	Non-Functional Device / Poor flow	Operator Error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	2	2	4	N/A	N/A	N/A	N/A	4
14	Balloon forming / dipping	Incomplete / Malformed Balloon	3	Non-Functional Device / Balloon Burst	Operator Error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	2	2	6	N/A	N/A	N/A	N/A	6
15	Valve installation	Loose valve / Valve disconnects	1	Non-Functional Device	Operator error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	1	1	1	N/A	N/A	N/A	N/A	1
16	Catheter Size Printing	Incorrect print information	2	Urethral damage / over inflation	Operator error	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form) HDX-6XXX SLQ-6XXX	1	1	1	2	N/A	N/A	N/A	N/A	2
17		Illegible	2		Equipment	Supplier Qualification, Receiving Inspection In-Process Inspection	Ningbo Yingmed Catheter Inspection Plans (Form)	2	1	2	4	N/A	N/A	N/A	N/A	4

This document contains confidential, proprietary information of Pathway Medtech and shall not be copied, reproduced, or otherwise disclosed to any third party without prior written permission there from. If any portion of this document is in hardcopy form, it must be considered and treated as uncontrolled "For Reference Only".

T-005 rev E

Pathway Medtech, LLC.		Type:	Risk Management (RM)		Doc. No.:		RM-0021		Revision:		F					
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter																
#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
							HDX-6XXX SLQ-6XXX									
18	Label Printing	Incorrect print information	2	Potential traceability issue	Operator error	Process instruction	MP-0016	2	2	2	4	Label inspection per QP-031/FM-0033	2	1	2	4
19		Illegible	2	Potential traceability issue	Equipment error	Process instruction	MP-0016	2	2	2	4	Label inspection per QP-031/FM-0033	2	1	2	4
20	Surface Treatment	Incomplete / No Treatment	1	No effect	Equipment error	Process Validation	TR-HDX-0007 (MVR) TR-SLQ-0007 (MVR)	1	1	1	1	N/A	N/A	N/A	N/A	1
						Start-up Verification Process Instructions	MP-SLQ-0002 TM-HDX-0001 TM-HDX-0002	1	1	1	1	N/A	N/A	N/A	N/A	1
						In-Process Inspection	Inspection Plans (Form) HDX-8XXX SLQ-8XXX	1	1	1	1	N/A	N/A	N/A	N/A	1
					Operator error	Process instruction	MP-HDX-0004	1	1	1	1	N/A	N/A	N/A	N/A	1
						In-Process Inspection	Inspection Plans (Form) HDX-8XXX SLQ-8XXX	1	1	1	1	N/A	N/A	N/A	N/A	1
21	Visual check pouch label	Accept Non-Conforming Material	1	Potential traceability issue	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	1
22		Reject Conforming Material	1	None	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	1
23	Apply Label to Pouch	Misaligned or wrinkled	1	Cosmetic issue	Operator Error	Process Instructions	MP-HDX-0001	1	1	1	1	N/A	N/A	N/A	N/A	1

Pathway Medtech, LLC.			Type:	Risk Management (RM)		Doc. No.:		RM-0021				Revision:		F		
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter																
#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
							MP-SLQ-0001									
24		Incorrect packaging configuration	1	Potential misbranding / labeling issue	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	2	2	2	2	N/A	N/A	N/A	N/A	2
25	Pouch Loading	Incorrect product configuration	1	Non-Functional Device	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	2	2	2	2	N/A	N/A	N/A	N/A	2
26	Pouch Sealing (Sterile Barrier)	Weak Seal	4	Breach in sterile barrier	Equipment error	Process validation Inspection process	TR-HDX-0007 (MVR) TR-SLQ-0007 (MVR) Inspection Plans (Form) HDX-8XXX SLQ-8XXX	1	1	1	4	N/A	N/A	N/A	N/A	4
					Operator error – incorrect sealing parameters	Process instruction Inspection process	TR-0074 Inspection Plans (Form) HDX-8XXX SLQ-8XXX	1	1	1	4	N/A	N/A	N/A	N/A	4
					Contamination	Environmental controls	QP-018	1	1	1	4	N/A	N/A	N/A	N/A	4
27		Airborne / Surface microbial contamination	4	Microbial contamination / Non-sterile device	Operator error	Gowning instructions	WI-0012	2	3	3	12	Product Sterilization	1	1	1	4
					Equipment error	Preventative Maintenance Program	QP-015	1	3	2	8	Product Sterilization	1	1	1	4
28		Material microbial contamination	4	Microbial contamination / Non-sterile device	Operator error	Material Handling and Storage	QP-021	1	1	1	4	Product Sterilization	1	1	1	4
29			4	Material degradation	Operator error	Material Handling and Storage	QP-021	1	1	1	4	N/A	N/A	N/A	N/A	4

0021 Rev 1 - ECO-001003 / ECO-001003 / In Production / Effective: 03/10/2024 10:33:24 AM

Pathway Medtech, LLC.			Type:	Risk Management (RM)		Doc. No.:		RM-0021				Revision:		F		
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter																
#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
		Adverse environmental conditions				Environmental Controls	QP-018									
					Equipment error	Calibration Program	QP-020	1	1	1	4	N/A	N/A	N/A	N/A	4
				Preventive Maintenance Program		QP-015										
30	Assemble Box	Incorrect / Weak box	1	Damage product	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	1
31	Visual check carton label	Accept Non-Conforming Material	1	Potential traceability issue	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	1
32		Reject Conforming Material	1	None	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	1
33	Apply Label to Box	Misaligned or wrinkled	1	Cosmetic issue	Operator Error	Process Instructions	MP-HDX-0001 MP-SLQ-0001	1	1	1	1	N/A	N/A	N/A	N/A	1
34	Packaging	Incorrect packaging configuration	2	Damage product	Operator error	Process instruction	MP-HDX-0001 MP-SLQ-0001	2	2	2	4	N/A	N/A	N/A	N/A	4
35	Sterilization	Product not processed / mix-up	4	Potential sterility issue	Operator error	Final Product Inspection	QP-010	1	1	1	4	N/A	N/A	N/A	N/A	4
36		Product duplicate process mix-up	1	Damaged / Non-functional / Degraded performance	Operator error	Process Instructions	QP-021	1	1	1	1	N/A	N/A	N/A	N/A	1
						EO Program	QP-032	1	1	1	1	N/A	N/A	N/A	N/A	1
37		Cycle Process Excursion	4	Potential sterility issue	Operator error	Supplier's SOP	WI-00863 GWI-08-004-09	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4

Pathway Medtech, LLC.			Type:	Risk Management (RM)			Doc. No.:		RM-0021			Revision:		F		
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter																
#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
						Cycle Specification	PWI-7374	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
						Supplier Qualification Program	QP-005	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
					Equipment error	Process Validation	TR 0075, TR-HDX-0009	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
						EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
38		No BI's placed	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
						Supplier's SOP	WI-00939	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
39		Incorrect BI's used	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
						Supplier's SOP	WI-00939	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
40		Improper Material Handling with Lab	4	Potential sterility issue	Operator error	EO Program	QP-032	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
						Supplier's SOP	WI-00939	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
						Supplier Qualification Program	QP-005	1	1	1	4	Final Product Inspection per QP-010	1	1	1	4
41		Final Product Inspection	Incomplete DHR	1	Potential traceability issue	Operator error	Process instruction	WI-0006 QP-010	1	1	1	1	N/A	N/A	N/A	N/A

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

#	Process Description / Function	Potential Failure Mode	S	Potential Effects of Failure	Potential Causes of Failure	Risk Controls ¹	Verification	P ¹	D ¹	F ¹	R ¹	Risk Controls ²	P ²	D ²	F ²	R ²
42		Incorrect DHR	1	Potential traceability issue	Operator error	Process instruction	WI-0006 QP-010	1	1	1	1	N/A	N/A	N/A	N/A	1
43		Accept Nonconforming Sterilization Process Record	4	Potential sterility issue	Operator error	Process instruction	QP-010	1	1	1	4	N/A	N/A	N/A	N/A	4
44		Accept Nonconforming Sterility Report	4	Potential sterility issue	Operator error	Process instruction	QP-010	1	1	1	4	N/A	N/A	N/A	N/A	4
45		Accept Nonconforming LAL/BET Report	4	Potential sterility issue	Operator error	Process instruction	QP-010	1	1	1	4	N/A	N/A	N/A	N/A	4
46		No CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	1
47		Incomplete CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	1
48		Incorrect CofC	1	N/A	Operator error	Process instruction	QP-010	1	1	1	1	N/A	N/A	N/A	N/A	1

6.0 RISK / BENEFIT ASSESSMENT

6.1 The following table indicates the risk / benefit analysis for the residual risks that were greater than the risk acceptability limit after all applicable controls / mitigations were applied.

Table 6 Risk / Benefit Assessment for Residual Risk above Acceptability Limit for Process Risk Assessment

#	Process Description	Potential Failure Mode	Potential Effect of Failure	Potential Cause of Failure	S	F ²	Risk/Benefit Analysis	Does the benefit outweigh the residual risk?
N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A

6.2 There was no residual risk identified that was above the risk acceptability limit.

Pathway Medtech, LLC.	Type:	Risk Management (RM)	Doc. No.:	RM-0021	Revision:	F
Process Failure Mode and Effects Analysis (FMEA) Risk Management Report, SILQ ClearTract Catheter						

7.0 CONCLUSION

- 7.1 All Identified risks regardless of acceptability were reduced as far as possible (AFAP). All risks, including residual risks and regardless of their dimension, were reduced “as far as possible” (AFAP). Risks have been reduced as low as possible by the following methods:
- Design the product with inherent safety
 - Appropriate testing to ensure design criteria were met
 - Protective measures (critical-to-quality)
 - Manufacturing and quality control measures
 - Provide information for safety
- 7.2 Based on a review of the risk-management process for the device described herein, the risk management plan based on ISO 14971:2019 has been implemented and the individual and overall residual risk have been reduced as low as possible regardless of the magnitude and the benefit of using this device will outweigh risks posed by the device during usage.

RM-0021_F_SILQ PFMEA

Final Audit Report

2024-05-16

Created:	2024-05-14
By:	Aaron Rogers (arogers@pathwaymedtech.com)
Status:	Signed
Transaction ID:	CBJCHBCAABAALVgmaFbLXBasx4z3U4QnMhaw14RdoOYa

"RM-0021_F_SILQ PFMEA" History

-  Document created by Aaron Rogers (arogers@pathwaymedtech.com)
2024-05-14 - 11:49:33 PM GMT
-  Document emailed to Aaron Rogers (arogers@pathwaymedtech.com) for signature
2024-05-14 - 11:49:38 PM GMT
-  Document emailed to Ethan Rao (ethanr@silq.tech) for signature
2024-05-14 - 11:49:39 PM GMT
-  Document emailed to Will Hail (whail@pathwaymedtech.com) for signature
2024-05-14 - 11:49:39 PM GMT
-  Document emailed to Michael Cassou (mcassou@pathwaymedtech.com) for signature
2024-05-14 - 11:49:39 PM GMT
-  Document e-signed by Aaron Rogers (arogers@pathwaymedtech.com)
Signature Date: 2024-05-14 - 11:49:45 PM GMT - Time Source: server
-  Email viewed by Ethan Rao (ethanr@silq.tech)
2024-05-15 - 0:02:46 AM GMT
-  Document e-signed by Ethan Rao (ethanr@silq.tech)
Signature Date: 2024-05-15 - 0:03:05 AM GMT - Time Source: server
-  Email viewed by Michael Cassou (mcassou@pathwaymedtech.com)
2024-05-15 - 3:04:21 PM GMT
-  Document e-signed by Michael Cassou (mcassou@pathwaymedtech.com)
Signature Date: 2024-05-15 - 3:05:57 PM GMT - Time Source: server
-  Email viewed by Will Hail (whail@pathwaymedtech.com)
2024-05-16 - 5:24:32 PM GMT



Adobe Acrobat Sign

 Document e-signed by Will Hail (whail@pathwaymedtech.com)

Signature Date: 2024-05-16 - 5:24:43 PM GMT - Time Source: server

 Agreement completed.

2024-05-16 - 5:24:43 PM GMT